



VistA Blood Establishment Computer Software (VBECS) – Vision Interface Configuration and Setup Guide, Version 3.0

July 2017

Department of Veterans Affairs
Enterprise Project Management Office

Revision History

Date	Revision	Description	Author
01-25-17	1.0	Initial version. (Defect 424731)	BBM team
		(Defect 497919) Updated Fig. 6, Fig. 7, Fig. 10, Fig. 22 and Fig. 30 to reflect changed configuration file names and additional rules. Updated Appendices A and B with modified test codes for phenotyping. Updated Appendices C and D with modified rules. Updated section Instrument Manager version with information that the interface requires Specimen Management Module activated in IM.	
05-05-17	2.0	Added Chapters 4 and 5 to the Set Up Automated Instrument section.	BBM team
06-30-17	3.0	(Defect 546927) Fixed typo in phenotyping codes in Appendix A and Appendix B.	BBM team

Table of Contents

REVISION HISTORY.....	2
TABLE OF CONTENTS.....	3
INTRODUCTION	5
RELATED MANUALS AND REFERENCE MATERIALS	6
SET UP AUTOMATED INSTRUMENT	7
1 Network Connectivity Setup	7
2 Testing Profiles Setup.....	7
3 User ID Setup	7
4 Manual versus Automated review	7
5 Limitations of the crossmatch testing.....	8
SET UP INSTRUMENT MANAGER	9
1 Instrument Manger Version.....	9
2 Instrument Manager to Automated Instrument Connectivity	9
3 Installing Instrument Driver	10
4 Set up Instrument Side Configuration	10
4.1 Download instrument configuration files	10
4.2 Import Instrument side configuration	12
4.3 Verify test code mapping for instrument side configuration	15
4.4 Verify rules for instrument side configuration.....	17
4.5 Configure rules for instrument side configuration.....	18
4.6 Import VBECS (HL7) side configuration.....	21
4.7 Verify test code mapping for VBECS side configuration.....	24
4.8 Verify rules for VBECS (HL7) side configuration.....	25
5 Set Up HL7 Connection to VBECS TEST	29
6 Set Up Instrument Connection	31
7 Test New Connections.....	33
8 Validate Instrument connectivity to VBECS TEST	34
9 Set up HL7 Connection to VBECS PROD.....	34
GLOSSARY	39
APPENDICES	40
APPENDIX A: INSTRUMENT SIDE MAPPING	40
APPENDIX B: HL7 (VBECS) SIDE MAPPING.....	44
APPENDIX C: INSTRUMENT SIDE RULES.....	48
APPENDIX D: VBECS (HL7) SIDE RULES	53
APPENDIX E: VISION TESTING PROFILES.....	54
APPENDIX F: CA SDM SUPPORT TICKET TEMPLATE.....	55
APPENDIX G: VALIDATION PLANNING AND EXAMPLE TEST SCENARIOS	56
TEST GROUP ONE: AI INTERFACE DISABLED	57
TEST GROUP TWO: VERIFY AI INDIVIDUAL TEST(S).....	58
TEST GROUP 2 SCENARIO 1: VERIFY AI TAS TEST	59
TEST GROUP 2 SCENARIO 2: VERIFY AI SEROLOGIC CROSSMATCH TEST	60
TEST GROUP 2 SCENARIO 3: VERIFY AI PATIENT DIAGNOSTIC TESTS.....	61


This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.


TEST GROUP 2 SCENARIO 4: VERIFY AI BLOOD UNIT TESTS	62
---	----


Introduction

VBECS is a Blood Bank application that facilitates ongoing compliance with Food and Drug Administration (FDA) standards for medical devices and enhances the VA's ability to produce high-quality blood products and services to veterans. The system follows blood bank standards, standards of national accrediting agencies, FDA regulations and VA policies.

VBECS 2.2.0 introduced a new interface for blood bank testing performed by blood bank instrumentation to VBECS. The implementation of the interface and its associated validation are described in this guide.

 *Unauthorized access or misuse of this system and/or its data is a federal crime. Use of all data must be in accordance with VA security and privacy policies.*


 *The U.S. FDA classifies this software as a medical device. Unauthorized modifications will render this device an adulterated medical device under Section 501 of the Medical Device Amendments to the Federal Food, Drug, and Cosmetic Act. Acquiring and implementing this software through the Freedom of Information Act requires the implementer to assume total responsibility for the software and become a registered manufacturer of a medical device, subject to FDA regulations. Adding to or updating VBECS software without permission is prohibited.*

 *Instructions in this Setup Guide must be followed for the interface to deliver information to VBECS. Local validation is required to confirm proper operation before use. Validation and verification is required to ensure connectivity to VBECS.*

This guide is provided to assist you with the multi-faceted required setup of your local blood bank testing instrument(s), Data Innovations Instrument Manager (DI IM) and VBECS to electronically transmit instrument test results to VBECS for use in the transfusion service.

There are specific setup requirements to test and transmit those testing results to VBECS for review using DI IM (**Figure 1**).

Figure 1: Hardware and Interface Configuration

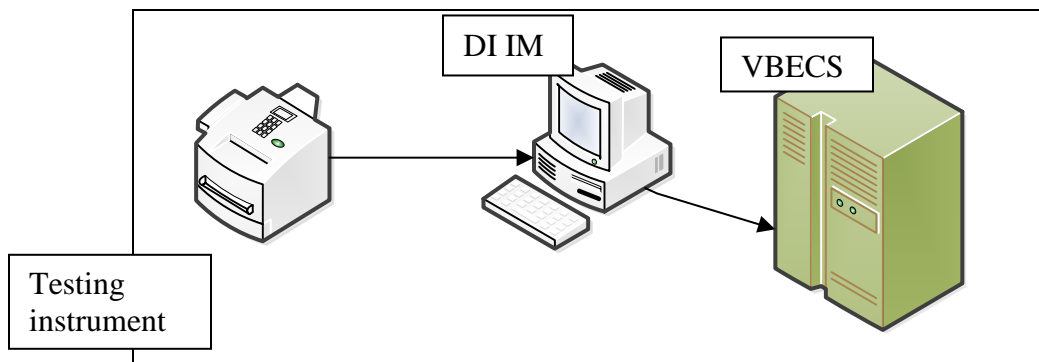


Your local testing instrument(s) communicates with DI IM via an instrument specific driver provided by DI that must be downloaded from DI and installed, locally.

DI IM communicates directly with VBECS via a generic HL7 interface driver that must be downloaded from DI and installed, locally. This driver is then customized for VBECS by downloading and installing the driver configuration file.

VBECS has an interface that must be configured in VBECS Administrator to receive messages from DI IM (**Figure 2**).

Figure 2: Setup Path

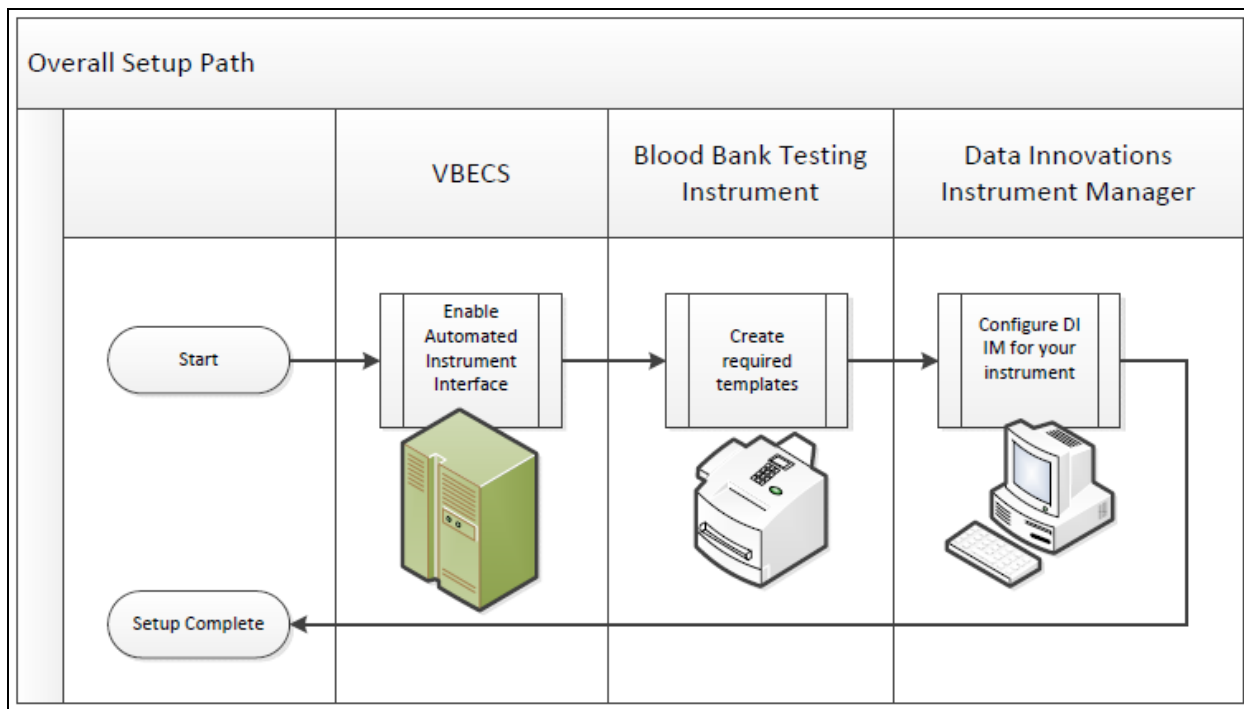

Related *Manuals and Reference Materials*

- *VistA Blood Establishment Computer Software (VBECS) 2.2.0 Technical Manual-Security Guide*
- *Data Innovations Instrument Manager Manual*
- *Blood Bank Analyzer User's Guide (Instrument Manual)*

Set Up Automated Instrument

1 Network Connectivity Setup

In order to ensure a proper functioning of an interface between an Automated Instrument and Instrument Manager the Instrument needs to be connected to the VA network. The static IP Address and Port number has to be assigned to the Instrument (further referred in this document as **<Instrument IP>** and **<Instrument Port>**). Please refer to the Instrument Manual or contact your vendor for the instructions about how to perform this setup. Check with Local IT staff to establish the connection to Instrument Manager.

2 Testing Profiles Setup



Test results sent using profiles that are not supported will not be accepted in VBECS.

Please refer to **Appendix E: Vision Testing Profiles** for the list of supported instrument profiles. It is up to each blood bank to choose which test profiles will be used as long as they are selected from the list in **Appendix E**.



AABB Standards require a serological XM to detect ABO incompatibility (5.16.1) and a local policy must be in place if the site is not performing an IS AHG as part of their serologic crossmatch test, manual or using an instrument.

3 User ID Setup



Failure to set up proper user IDs on an instrument will prevent instrument results from being accepted in VBECS.

In order for VBECS to properly recognize the person who performed testing on an Instrument, all users accessing the instrument must have their user IDs set up to match their network user IDs (e.g., VHAISHBURNISK).

Please refer to the Instrument Manual or contact your vendor for the instructions about how to perform setup of user IDs.

4 Manual versus Automated review

Vision can be configured to automatically review and send test results to LIS. However, in this configuration the test results are not associated with the name of the person who performed testing on the instrument. VBECS requires that all results sent from an instrument are accompanied by the ID of the user who initiated testing, therefore for the interface to function properly Vision MAY NOT be configured for the full automated review of the results.

The manufacturer recommends that the following configuration is used to ensure that test results are always accompanied by user ID:

Enforce Manual Review: No

Enforce Manual Review for Cards with Error Grades: Yes

All Maintenance Tasks must be completed successfully: Yes

All column reactions must be without errors: Yes

Valid Test Result Interpretation (incl. non-discrepant columns): Yes

Above/below positive reaction threshold: Yes

Entire order must be completed: No



Please consult user's manual for the Vision or contact Ortho Diagnostics support if you need help in changing the manual review settings on the instrument.

5 Limitations of the crossmatch testing

Due to the limitation of the Data Innovation's driver for Vision it is not possible to successfully send results of the crossmatch test for multiple donor samples from one crossmatch order. In order to perform crossmatch for multiple donor samples users have to create multiple crossmatch orders on the Vision. Each crossmatch order may be associated with only one donor sample. Such results may then be transmitted to VBECS, one crossmatch order/donor sample at a time.

Set Up Instrument Manager

1 Instrument Manger Version

Please verify that you are using Instrument Manager Version 8.13 or greater (Figure 3). In Instrument Manager, navigate to **Help -> About Instrument Manager**.

Figure 3: Instrument Manager Version Screen



*If your version of Instrument Manger is older than 8.13 please STOP executing this guide and update the software first. **Do not proceed until the issue is resolved.***



*If your Instrument Manager is greater than 8.13 you may proceed. The user must execute the instructions and validate functionality on newer version. Discrepancies in the instructions must be reported as a CA-SDM ticket. See **Appendix F**.*



*Instrument Manager must have the **Specimen Management Module** licensed and activated. **Do not proceed until this issue is resolved.***

2 Instrument Manager to Automated Instrument Connectivity

Please contact your local network administrative staff and ensure that your local network allows two-way TCP/IP connectivity between **<Instrument Manager IP>** address and **<Instrument IP>** address on **<Instrument Port>**.

3 Installing Instrument Driver

Two drivers are required for the correct operation of the Automated Instrument interface to VBECS (Figure 4):

Figure 4: Required Drivers

Driver	Name	Version	Description
Ortho Vision	ocdvisii	8.00.0003	Instrument to DI IM driver
Data Innovations LLC, Configurable HL7	diihl7ml	8.00.0049	DI IM to VBECS (HL7) driver

Please refer to the user's manual for Instrument Manager or contact Data Innovations for the instructions on how to install drivers for Instrument Manager.

After successful installation of drivers please go to **Report -> Available Drivers** menu option in Instrument Manager and verify that the drivers listed in Figure 4 are present.

4 Set up Instrument Side Configuration



Execute instructions in this section for each Vision instrument that will be connecting to VBECS.



Modifying rules or test code mappings as imported using this Instrument Manager configuration Setup Guide may lead to malfunction of the Automated Instrument to VBECS interface.

Prerequisites for the Instrument Manager Configuration files download:

- You must be an administrator on *<Instrument Manager Server>*.
- Once the above prerequisite is met you may proceed.


4.1 Download instrument configuration files



Use local procedures for copying the instrument configuration files to the Instrument Manager server.


1. Navigate to **<http://vaww.oed.portal.va.gov/projects/vbecs/Pages/Instrument-Manager-Configuration-Files-Release.aspx>** (Figure 5).

Figure 5: Connecting to the VBECS SharePoint

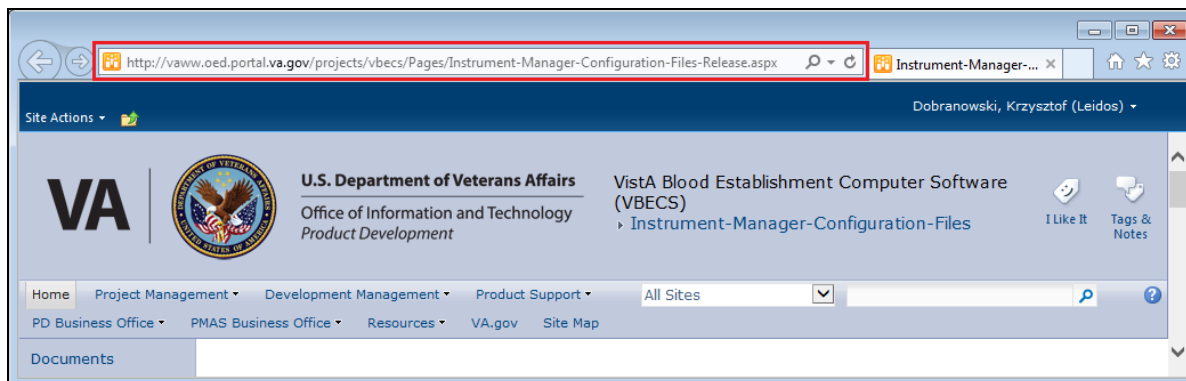


2. To download a file from the SharePoint, right-click on it and select **Save target as** (Figure 6).

Figure 6: Example of Save target as...

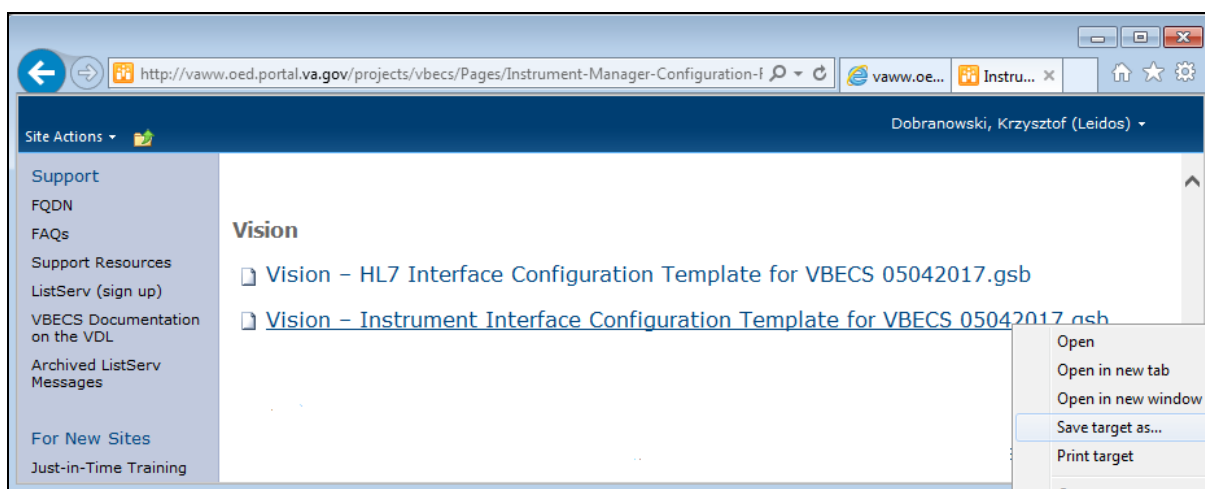

3. In the next screen, specify the directory to save (Figure 7).

Figure 7: Example of Save As

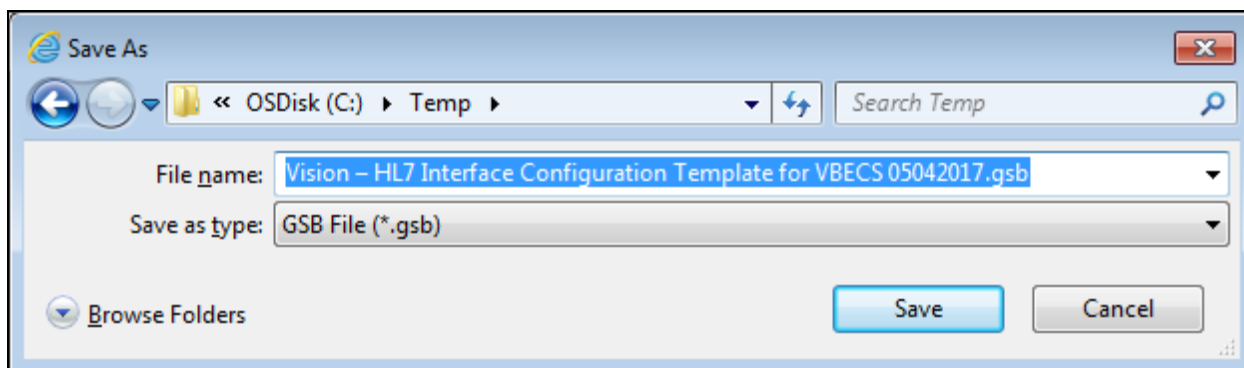

4. Save both the **Vision – HL7 Interface Configuration Template for VBECS 05042017.gsb** and the **Vision – Instrument Interface Configuration Template for VBECS 05042017.gsb** files to the **C:\temp** directory.

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

5. Per local procedures, copy both files to **C:\temp** on the Instrument Manager server
6. On the Instrument Manager server; Click **Start**, and in the “**Search programs and files**” box type **Run** and hit enter. Type **powershell** and click **OK** to launch PowerShell
7. Copy and paste or type the following commands to generate checksums for configuration files:

```
certutil -hashfile "C:\Temp\Vision – HL7 Interface Configuration Template for VBECS  
05042017.gsb" MD5 <press Enter>
```

```
certutil -hashfile "C:\Temp\Vision – Instrument Interface Configuration Template for VBECS  
05042017.gsb" MD5 <press Enter>
```



To copy, highlight the lines in grey and enter CTRL-C.

To insert the copied line into a PowerShell window, right click in the PowerShell window and select “Paste”.

8. Verify that checksums for both files match those shown in Figure 8.

Figure 8: Instrument Manager Configuration File checksums

```
Administrator: Windows PowerShell
PS C:\> certutil -hashfile "C:\Temp\Vision – HL7 Interface Configuration Template for VBECS 05042017.gsb" MD5
MD5 hash of file C:\Temp\Vision – HL7 Interface Configuration Template for VBECS 05042017.gsb:
6b c3 7c 1b fd c6 cb 23 31 01 0f dc 3c de 11 ee
CertUtil: -hashfile command completed successfully.
PS C:\> certutil -hashfile "C:\Temp\Vision – Instrument Interface Configuration Template for VBECS 05042017.gsb" MD5
MD5 hash of file C:\Temp\Vision – Instrument Interface Configuration Template for VBECS 05042017.gsb:
88 a1 7a fb 87 2f 56 f1 bb 49 0a 45 bb eb 29 d2
CertUtil: -hashfile command completed successfully.
PS C:\>
```




If the checksums do not match, stop and file a national CA SDM ticket to coordinate assistance with installation using the template in Appendix F. Do not proceed until the issue is resolved.

9. Close the **PowerShell** window.

4.2 Import Instrument side configuration


1. After downloading and verifying configuration files, open Instrument Manager and navigate to **Configuration -> Configuration Editor**.
2. Click on the **Import** button (Figure 9).

Figure 9: Example of Configuration Editor

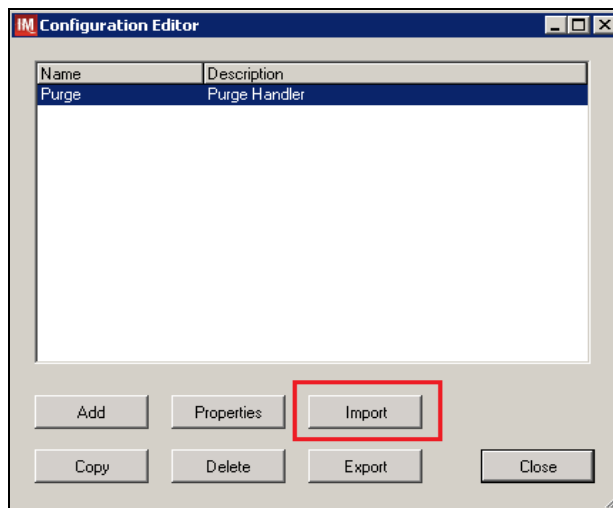

3. Once the import screen opens, click the **Browse** button and select **C:\Temp** folder (Figure 10).
4. Click OK.

Figure 10: Example of Configuration Import Screen

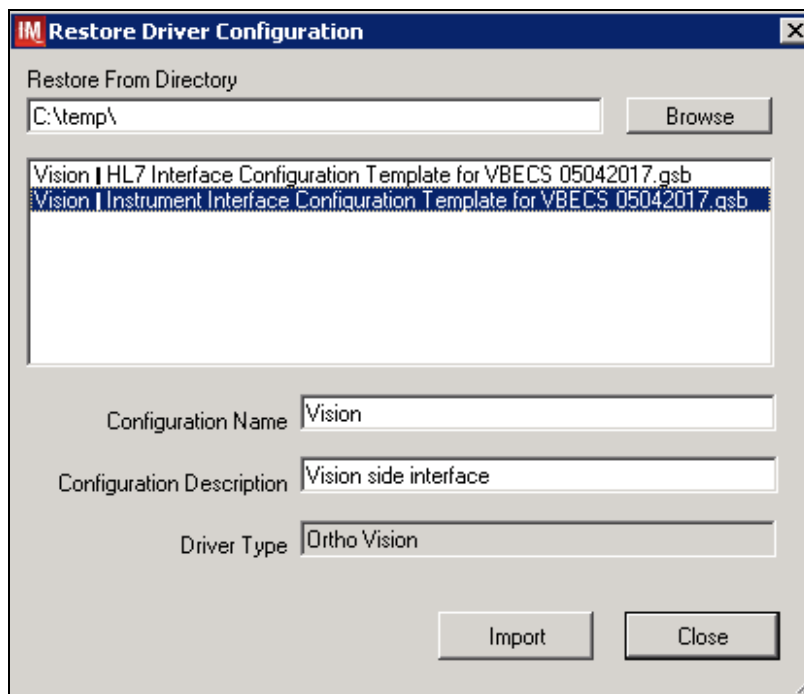

5. Select **Vision – Instrument Interface Configuration Template for VBECS 05042017.gsb** file from the list.
6. Enter **Configuration Name** that contains 3 letter location code of the instrument (e.g. HIN for Hines VAMC), word Vision and sequence number (1 for the first instrument, 2 for the second etc.).
Example **Configuration Name** for instrument configuration located at Hines would be:


This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

HIN_Vision_1

This configuration name will be further referred in this document as *<Instrument Side Configuration>*.

7. Enter **Configuration Description** and click **Import** button. Verify that the following confirmation window displays (Figure 11).


Figure 11: Example of Successful Configuration Import

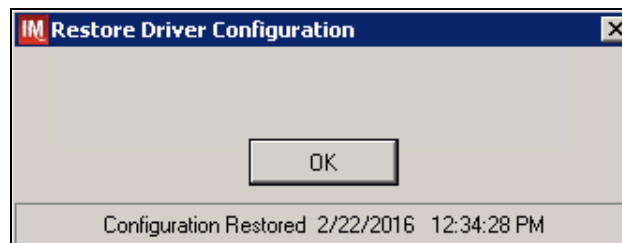

If you are using newer version of the driver than the one mentioned in section 3, the Instrument Manager will warn you about the discrepancy in driver versions. Please acknowledge this warning and continue.

8. Click **OK** and then close the **Restore Driver Configuration** window.
9. Verify that **Configuration Editor** shows the new configuration on the list (Figure 12).

Figure 12: Example of Configuration Editor Window Showing Newly Imported Configuration

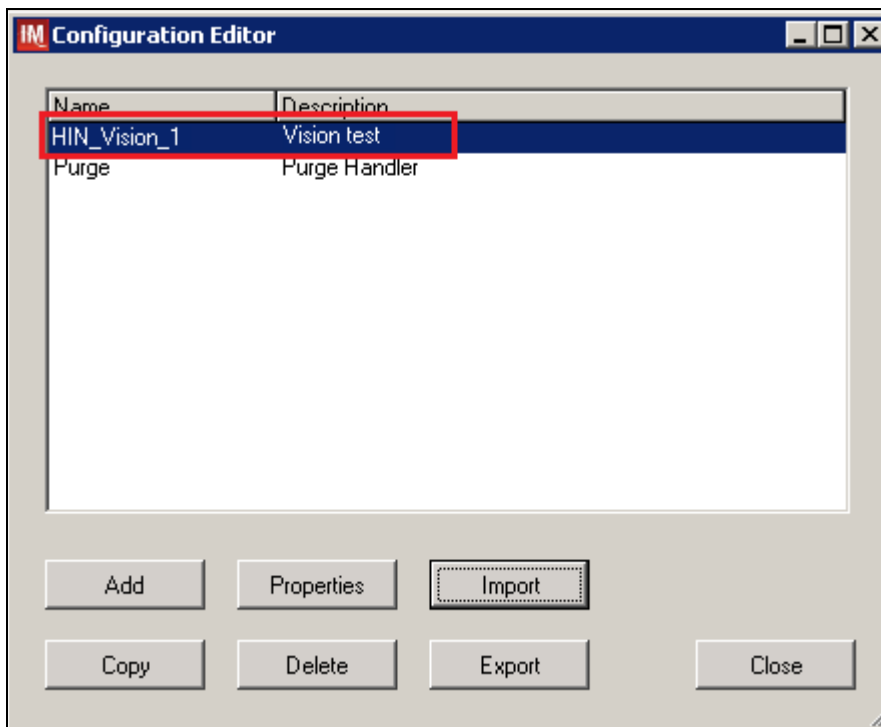

10. Close the Configuration Editor window.

4.3 Verify test code mapping for instrument side configuration

1. Navigate to **Reports -> Configuration Options and Mappings**.
2. Select the **<Instrument Side Configuration Name>** from the pull down menu (Figure 13).
3. Select the **Results Test Code Mapping Tab** (Figure 13).
4. Verify the Results Test Code Mapping Report Tab matches the list in- **Appendix A: Instrument Side mapping**.

Figure 13: Example of Results Test Code Mapping Tab


Instrument Test Code	IM Test Code	Fluid	Fluid Description	Test Resulting Option	Display Name	Default Test Code	LOINC
VBECS ABDRev-ABO	VBECS ABDRev-ABO	CENTBLOOD	CENTBLOOD	OD	ABOInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO	VBECS ABDRev-ABO	PACKEDCELLS	PACKEDCELLS	OD	ABOInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-A1-Cells	VBECS ABDRev-ABO-A1-Cells	CENTBLOOD	CENTBLOOD	OD	A1Cells	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-A1-Cells	VBECS ABDRev-ABO-A1-Cells	PACKEDCELLS	PACKEDCELLS	OD	A1Cells	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-A	VBECS ABDRev-ABO-Anti-A	CENTBLOOD	CENTBLOOD	OD	AntiA	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-A	VBECS ABDRev-ABO-Anti-A	PACKEDCELLS	PACKEDCELLS	OD	AntiA	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-B	VBECS ABDRev-ABO-Anti-B	CENTBLOOD	CENTBLOOD	OD	AntiB	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-B	VBECS ABDRev-ABO-Anti-B	PACKEDCELLS	PACKEDCELLS	OD	AntiB	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-B-Cells	VBECS ABDRev-ABO-B-Cells	CENTBLOOD	CENTBLOOD	OD	BCells	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-B-Cells	VBECS ABDRev-ABO-B-Cells	PACKEDCELLS	PACKEDCELLS	OD	BCells	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh	VBECS ABDRev-Rh	CENTBLOOD	CENTBLOOD	OD	RhInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh	VBECS ABDRev-Rh	PACKEDCELLS	PACKEDCELLS	OD	RhInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Anti-D	VBECS ABDRev-Rh-Anti-D	CENTBLOOD	CENTBLOOD	OD	AntiD	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Anti-D	VBECS ABDRev-Rh-Anti-D	PACKEDCELLS	PACKEDCELLS	OD	AntiD	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Ctrl	VBECS ABDRev-Rh-Ctrl	CENTBLOOD	CENTBLOOD	OD	DControl	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Ctrl	VBECS ABDRev-Rh-Ctrl	PACKEDCELLS	PACKEDCELLS	OD	DControl	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr	VBECS ABS 2 Cell-ABScr	CENTBLOOD	CENTBLOOD	OD	ABScInterp	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr	VBECS ABS 2 Cell-ABScr	PACKEDCELLS	PACKEDCELLS	OD	ABScInterp	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel I	VBECS ABS 2 Cell-ABScr-0.8-Sel I	CENTBLOOD	CENTBLOOD	OD	SC1	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel I	VBECS ABS 2 Cell-ABScr-0.8-Sel I	PACKEDCELLS	PACKEDCELLS	OD	SC1	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel II	VBECS ABS 2 Cell-ABScr-0.8-Sel II	CENTBLOOD	CENTBLOOD	OD	SC2	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel II	VBECS ABS 2 Cell-ABScr-0.8-Sel II	PACKEDCELLS	PACKEDCELLS	OD	SC2	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr	VBECS ABS 3 Cell-ABScr	CENTBLOOD	CENTBLOOD	OD	ABScInterp	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr	VBECS ABS 3 Cell-ABScr	PACKEDCELLS	PACKEDCELLS	OD	ABScInterp	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr-0.8-Surg1	VBECS ABS 3 Cell-ABScr-0.8-Surg1	CENTBLOOD	CENTBLOOD	OD	SC1	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr-0.8-Surg1	VBECS ABS 3 Cell-ABScr-0.8-Surg1	PACKEDCELLS	PACKEDCELLS	OD	SC1	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr-0.8-Surg2	VBECS ABS 3 Cell-ABScr-0.8-Surg2	CENTBLOOD	CENTBLOOD	OD	SC2	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr-0.8-Surg2	VBECS ABS 3 Cell-ABScr-0.8-Surg2	PACKEDCELLS	PACKEDCELLS	OD	SC2	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr-0.8-Surg3	VBECS ABS 3 Cell-ABScr-0.8-Surg3	CENTBLOOD	CENTBLOOD	OD	SC3	<input checked="" type="checkbox"/>	

STOP *If mismatches in Test Codes names, missing or extra Test Codes are encountered, file a national CA SDM ticket to coordinate assistance with installation using the template in Appendix F. Do not proceed until the issue is resolved.*

4.4 Verify rules for instrument side configuration

1. Remaining in the Configuration Options and Mapping Window, select the **Rules Tab** (Figure 14).
2. Verify the Rules Tab matches the list in *Appendix C: Instrument Side Rules*.

Figure 14: Example of Rules Tab

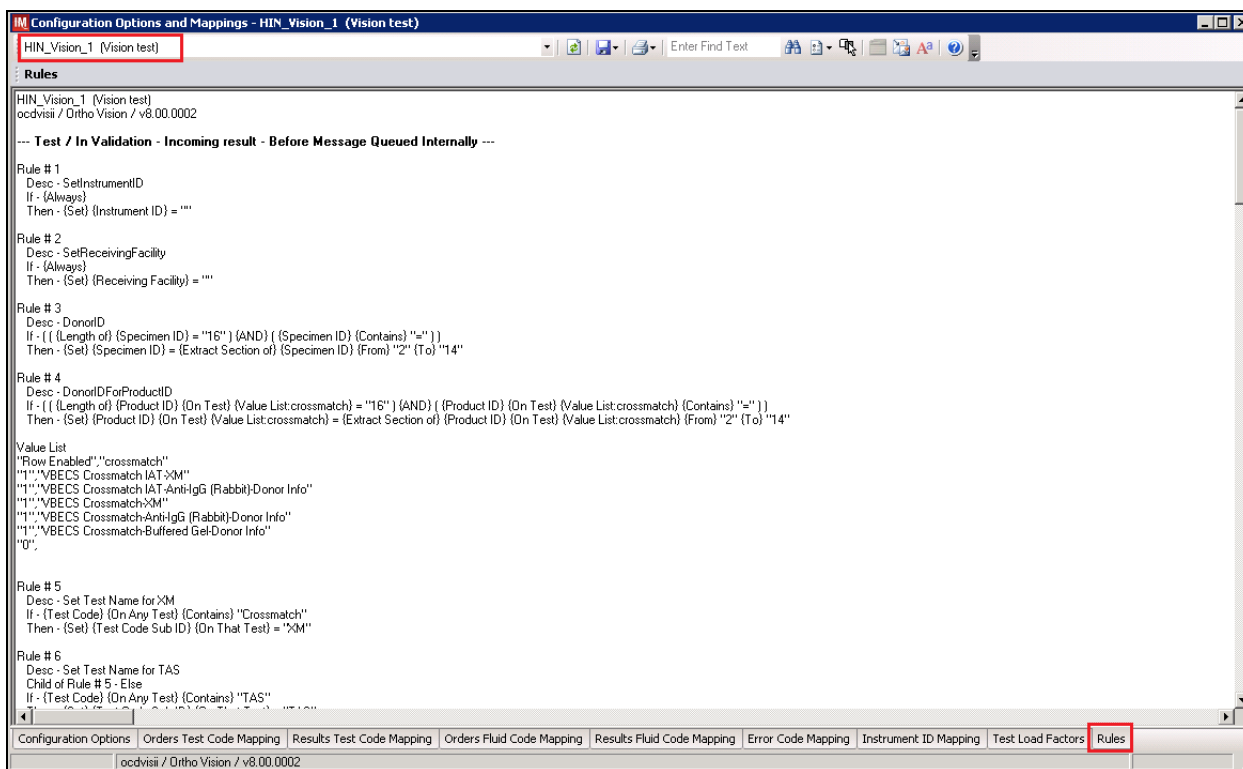

3. Close the Configuration Options and Mappings window.


STOP *If problems are encountered, file a national CA SDM ticket to coordinate assistance with installation using the template in Appendix F. Do not proceed until the issue is resolved.*

4.5 Configure rules for instrument side configuration

In this section you will establish the name of the instrument associated with test results for VBECS reports. *If you are setting up multiple instruments, make sure that each has a unique name.*

1. Navigate to Configuration -> Specimen Management Configuration -> Rules Processing.
2. Select **<Instrument Side Configuration Name>** from the pull down menu (Figure 15).
3. Locate rule **SetInstrumentID**.
4. Click on the **Then** line in the rule.
5. Modify the rule by typing **<Instrument Name>** between the quotation marks in the lower box in (Figure 15).

Figure 15: Example of Instrument Name Setup

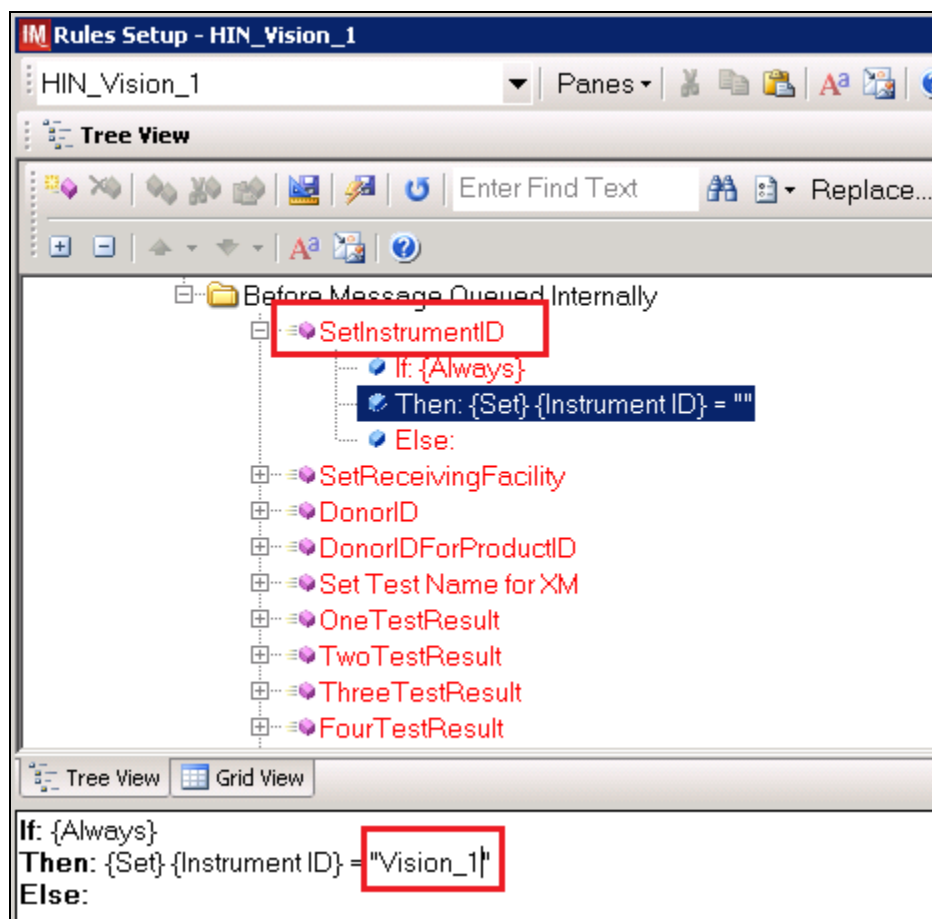

6. Locate rule **SetReceivingFacility**. Click on **Then** line in the rule. Modify the rule by typing **<Division Code>** between the quotation marks as shown in (Figure 16).


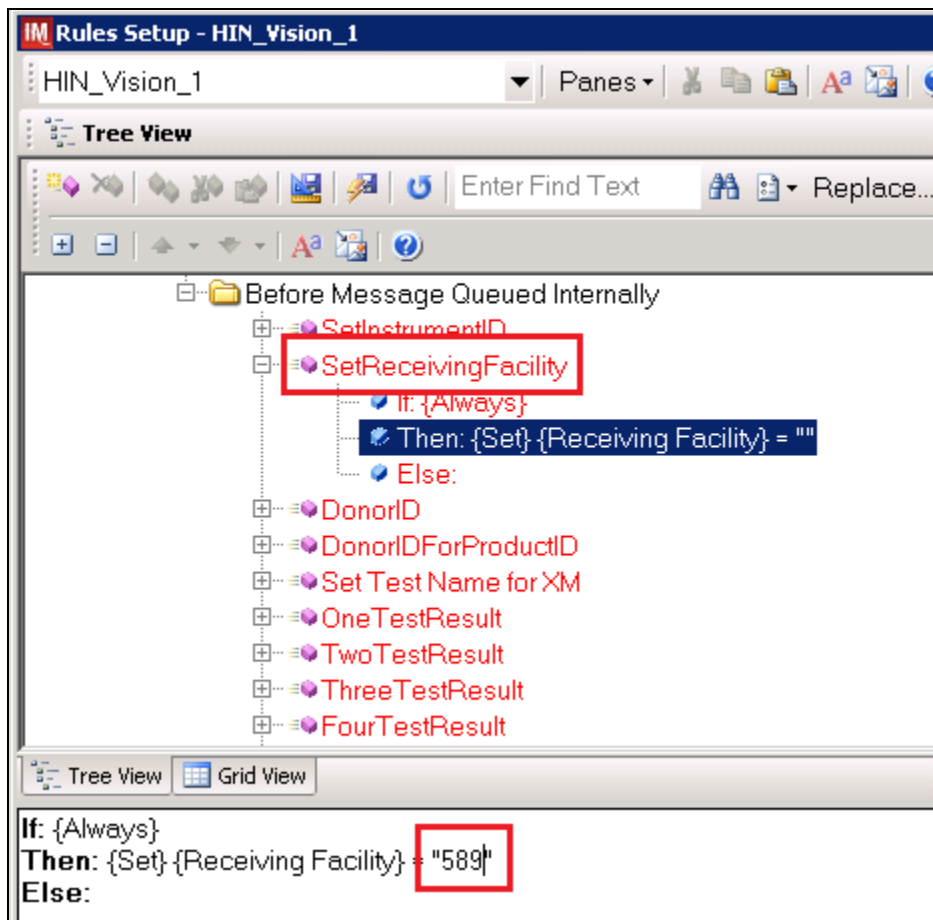
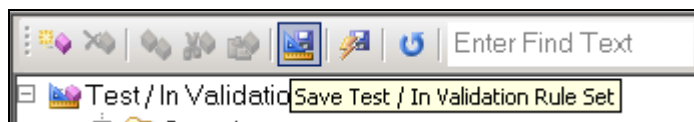
 **<Division Code>** also known as Station Number in VistA is a unique alphanumeric code that is associated with each hospital (e.g., **589** for VA Heartland West VAMC). This code must match the division code configured in **VBECS Administrator** application for a given site.

Figure 16: Example of Division Code Setup




7. Click on the **Save Test / In Validation Rule Set** button located in the toolbar above the rules (Figure 17).


Figure 17: Example of Save Rules Button



8. Verify that the message in Figure 18 is received:

Figure 18: Example of Rule Save

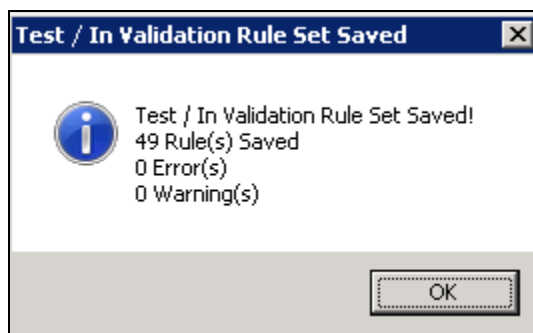

 ***If problems are encountered, file a national CA SDM ticket to coordinate assistance with installation using the template in Appendix F. Do not proceed until the issue is resolved.***


9. Click on the **Save Live Rule Set** button located in the upper toolbar and click **Yes** to confirm (Figure 19).

Figure 19: Example of Save Rules in Live Set Button



10. Navigate to Reports->Configuration Options and Mappings and select the Rules Tab. Scroll down to ***Live – Incoming result – Before Message Queued Internally*** (Figure 20).
11. Verify that the Live Rule Set matches ***Appendix C: Instrument Side Rules*** and includes changes made to **SetReceivingFacility** (Figure 20) and **SetInstrumentID**.

Figure 20: Example of Live Rules Set View

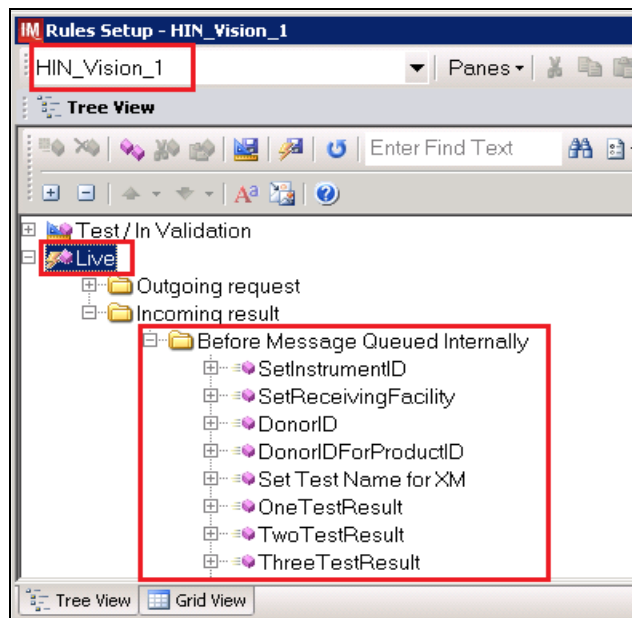

STOP If problems are encountered, file a national CA SDM ticket to coordinate assistance with installation using the template in Appendix F. Do not proceed until the issue is resolved.

12. Close the **Configuration Options and Mappings** window.


4.6 Import VBECS (HL7) side configuration

i Only one HL7 configuration is needed even if you use multiple Vision instruments. All Visions will share the same HL7 configuration.

STOP Modifying rules or test code mappings in the Instrument Manager configuration outside of this Setup Guide is not allowed and may lead to malfunction of the Automated Instrument to VBECS interface.

1. Navigate to **Configuration -> Configuration Editor** (Figure 21).
2. Click on the **Import** button.

Figure 21: Example of Configuration Editor Window

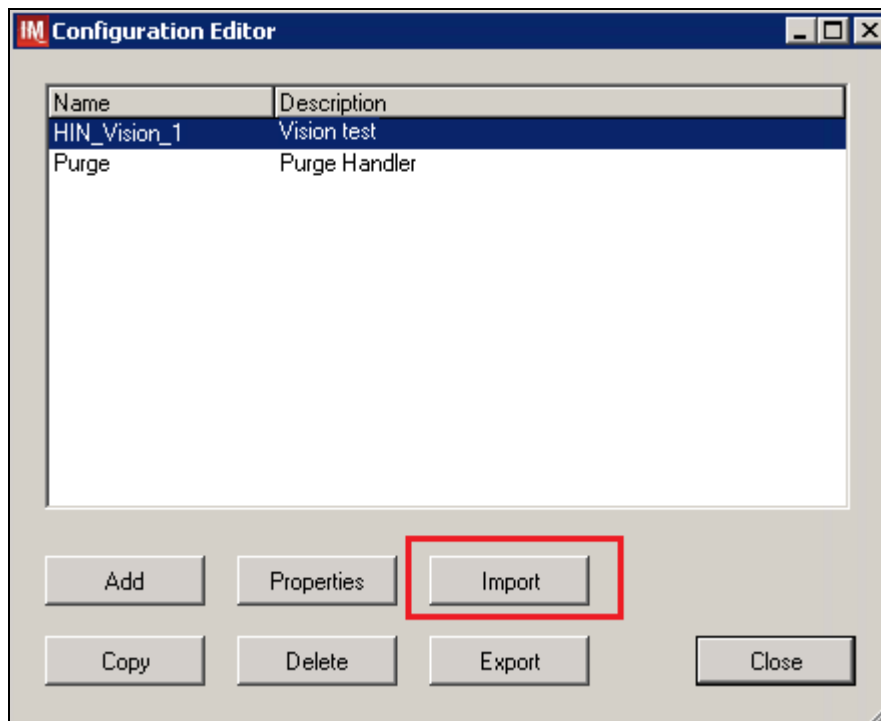


3. Once the **Restore Driver Configuration** window opens, click **Browse** the button and select **C:\Temp** folder (Figure 22).
4. Select **Vision – HL7 Interface Configuration Template for VBECS 05042017.gsb** file from the list. Enter **Configuration Name** that contains 3 letter location code of the instrument (e.g. **HIN** for Hines VAMC), word **VBECS** and sequence number (**1** for the first configuration, **2** for the second etc.). Example Configuration Name for VBECS side configuration located at Hines would be:

HIN_VBECS_1


This configuration name will be further referred in this document as *<HL7 Side Configuration>*.


Figure 22: Example of Configuration Import Window

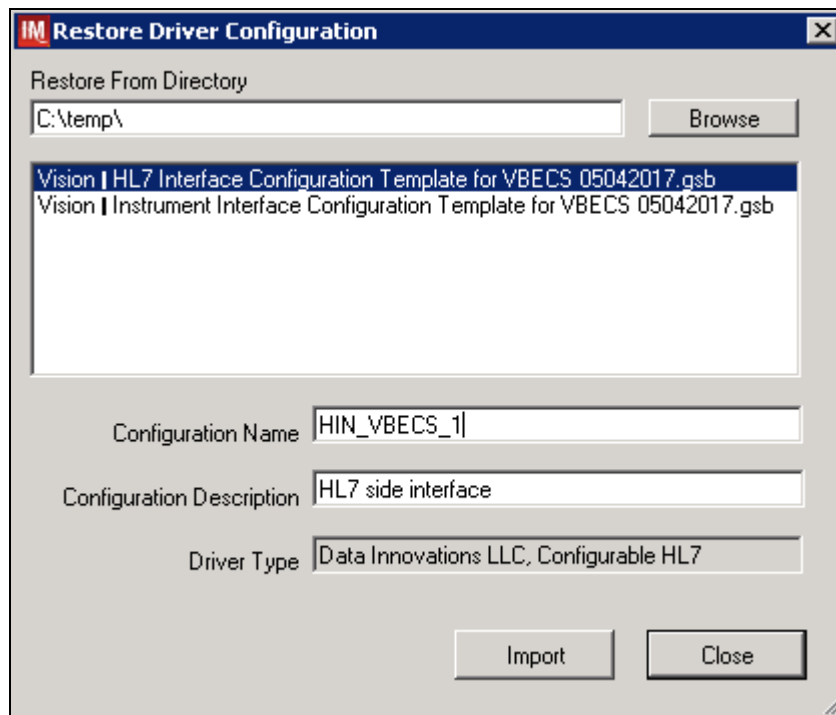

5. Enter **Configuration Description** and click **Import** button. Verify that the confirmation window in Figure 23 displays.


Figure 23: Example of Successful Configuration Import

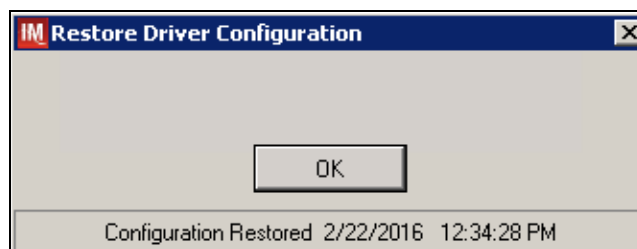

 *If you are using newer version of the driver than the one mentioned in section 3, the Instrument Manager will warn you about the discrepancy in driver versions. Please acknowledge this warning and continue.*

6. Click **OK** and close the **Restore Driver Configuration** window.
7. Verify that **Configuration Editor** shows the new configuration on the list (Figure 24).

Figure 24: Example of Newly Imported HL7 Configuration

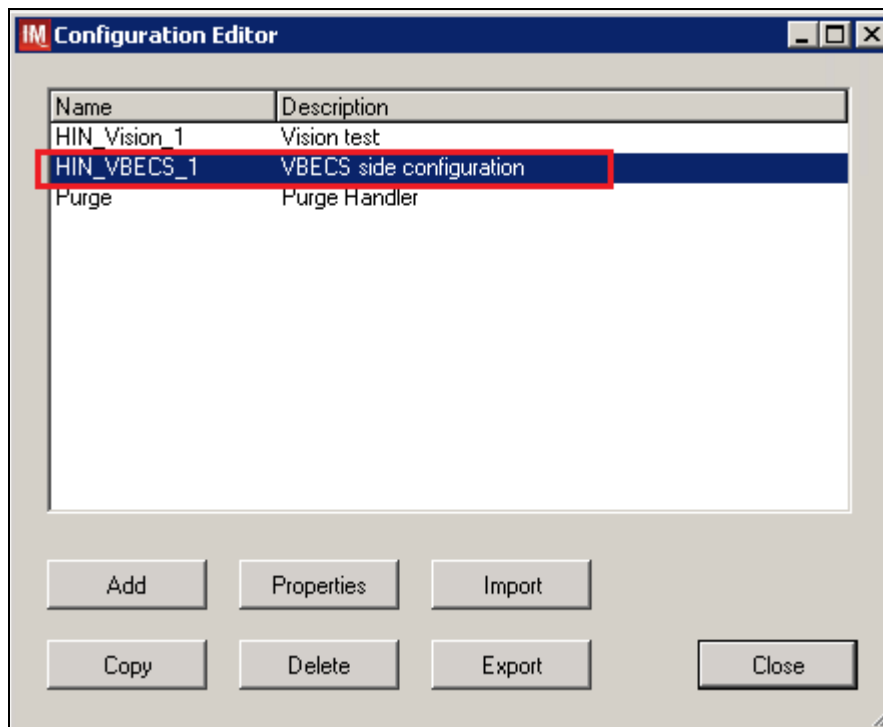

4.7 Verify test code mapping for VBECS side configuration

5. Navigate to **Reports -> Configuration Options and Mappings**.
6. Select the **<HL7 Side Configuration Name>** from the pull down menu (Figure 25).
7. Select the **Results Test Code Mapping Tab** (Figure 25).
8. Verify the Results Test Code Mapping Report Tab matches the list in **Appendix B: HL7 (VBECS) Side Mapping**.

Figure 25: Example of HL7 Configuration Report Window

The screenshot shows a window titled "IM Configuration Options and Mappings - HIN_VBECS_1 (VBECS side configuration)". Below the title bar is a search bar with the text "HIN_VBECS_1 (VBECS side configuration)". The main area is titled "Results Test Code Mapping" and contains a table with the following columns: Instrument Test Code, IM Test Code, Fluid, Fluid Description, Test Resulting Option, Display Name, Default Test Code, and LOINC. The table lists various test codes and their corresponding mappings, with checkboxes in the "Default Test Code" column.

Instrument Test Code	IM Test Code	Fluid	Fluid Description	Test Resulting Option	Display Name	Default Test Code	LOINC
▶ AHGInterp	AHGInterp	CENTBLOOD	CENTBLOOD	D	AHGInterp	<input checked="" type="checkbox"/>	
AHGInterp	AHGInterp	PACKEDCELLS	PACKEDCELLS	D	AHGInterp	<input checked="" type="checkbox"/>	
ISInterp	ISInterp	CENTBLOOD	CENTBLOOD	D	ISInterp	<input checked="" type="checkbox"/>	
ISInterp	ISInterp	PACKEDCELLS	PACKEDCELLS	D	ISInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO	ABOInterp	CENTBLOOD	CENTBLOOD	D	ABOInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO	ABOInterp	PACKEDCELLS	PACKEDCELLS	D	ABOInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-A1-Cells	A1Cells	CENTBLOOD	CENTBLOOD	D	A1Cells	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-A1-Cells	A1Cells	PACKEDCELLS	PACKEDCELLS	D	A1Cells	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-A	AntiA	CENTBLOOD	CENTBLOOD	D	AntiA	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-A	AntiA	PACKEDCELLS	PACKEDCELLS	D	AntiA	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-B	AntiB	CENTBLOOD	CENTBLOOD	D	AntiB	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-Anti-B	AntiB	PACKEDCELLS	PACKEDCELLS	D	AntiB	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-B-Cells	BCells	CENTBLOOD	CENTBLOOD	D	BCells	<input checked="" type="checkbox"/>	
VBECS ABDRev-ABO-B-Cells	BCells	PACKEDCELLS	PACKEDCELLS	D	BCells	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh	RhInterp	CENTBLOOD	CENTBLOOD	D	RhInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh	RhInterp	PACKEDCELLS	PACKEDCELLS	D	RhInterp	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Anti-D	AntiD	CENTBLOOD	CENTBLOOD	D	AntiD	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Anti-D	AntiD	PACKEDCELLS	PACKEDCELLS	D	AntiD	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Ctrl	DControl	CENTBLOOD	CENTBLOOD	D	DControl	<input checked="" type="checkbox"/>	
VBECS ABDRev-Rh-Ctrl	DControl	PACKEDCELLS	PACKEDCELLS	D	DControl	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr	ABSInterp	CENTBLOOD	CENTBLOOD	D	ABSInterp	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr	ABSInterp	PACKEDCELLS	PACKEDCELLS	D	ABSInterp	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel I	SC1	CENTBLOOD	CENTBLOOD	D	SC1	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel I	SC1	PACKEDCELLS	PACKEDCELLS	D	SC1	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel II	SC2	CENTBLOOD	CENTBLOOD	D	SC2	<input checked="" type="checkbox"/>	
VBECS ABS 2 Cell-ABScr-0.8-Sel II	SC2	PACKEDCELLS	PACKEDCELLS	D	SC2	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr	ABSInterp	CENTBLOOD	CENTBLOOD	D	ABSInterp	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr	ABSInterp	PACKEDCELLS	PACKEDCELLS	D	ABSInterp	<input checked="" type="checkbox"/>	
VBECS ABS 3 Cell-ABScr-0.8-Sura1	SC1	CENTBLOOD	CENTBLOOD	D	SC1	<input checked="" type="checkbox"/>	

At the bottom of the window, there are tabs for "Configuration Options", "Orders Test Code Mapping", "Results Test Code Mapping" (which is selected), "Orders Fluid Code Mapping", "Results Fluid Code Mapping", and "Rules". The status bar at the bottom indicates "dihi7ml / Data Innovations LLC, Configurable HL7 / v8.00.0049".

9. Close **Configuration Options and Mappings** window.



*If mismatches in Test Codes names, missing or extra Test Codes are encountered, **file a national CA SDM ticket to coordinate assistance with installation using template in Appendix F.** Do not proceed until the issue is resolved.*

4.8 Verify rules for VBECS (HL7) side configuration




Modifying rules or test code mappings as imported using this Instrument Manager configuration Setup Guide may lead to malfunction of the Automated Instrument to VBECS interface.

1. Navigate to **Configuration -> Specimen Management Configuration -> Rules Processing**.

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.


2. On the **Rule Processing** window, select **<HL7 Side Configuration>** from the drop-down in the upper left corner.
3. Next, click on the **Quick Filter** icon located on the right side of the **Configuration Name** (Figure 26)

Figure 26: Example of Rules Processing Toolbar

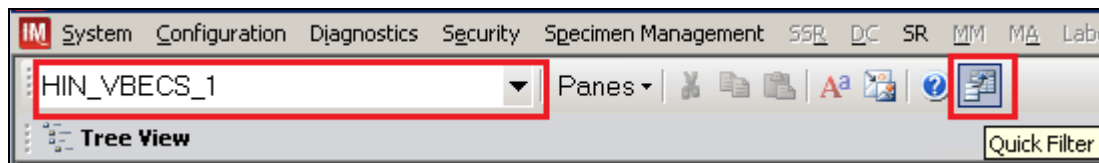


4. On the **Quick Filter** window, check the box next to the **Outgoing Result** and click **Close** button (Figure 27).

Figure 27: Example of Quick Filter Window

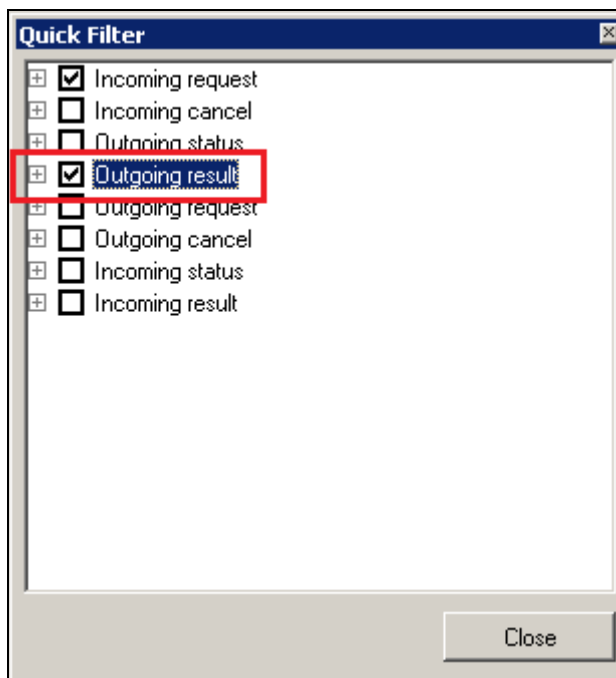

5. Expand **Outgoing Result** section and **Before Message Sent to this Connection** sections (Figure 28).

Figure 28: Example of Rules Window

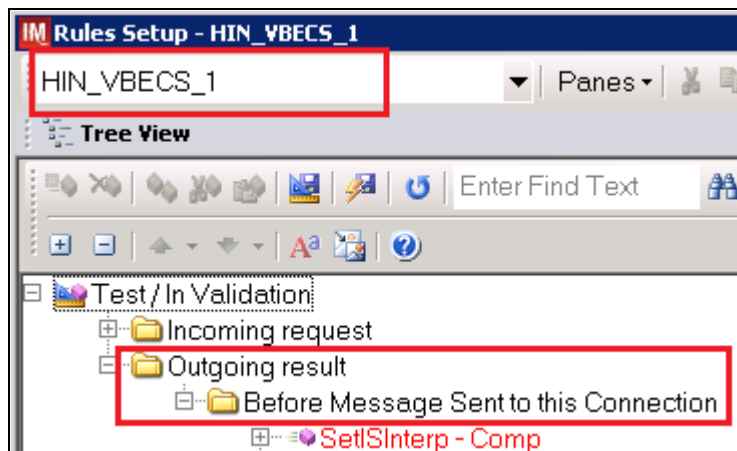

6. Expand all rules and verify that they match rules listed in *Appendix D: VBECS (HL7) Side Rules*.

STOP If mismatches in rules are encountered, file a national CA SDM ticket to coordinate assistance with installation using template in Appendix F. Do not proceed until the issue is resolved.


7. Click on the **Save Live Rule Set** button located in the upper toolbar and click **Yes** to confirm (Figure 29).

Figure 29: Example of Save Rules in Live Set Button



8. Expand the **Live Rule Set** and verify that rule text matches *Appendix D: VBECS (HL7) Side Rules*. (Figure 30).

Figure 30: Example of Live Rules Set View

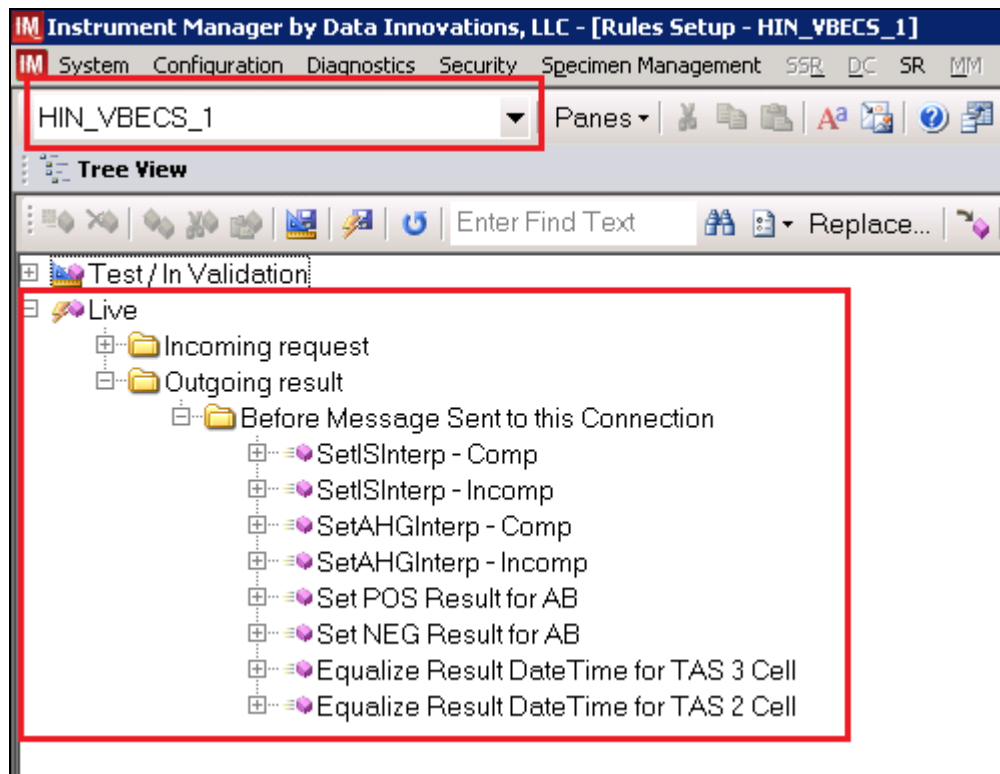

*If problems are encountered, **file a national CA SDM ticket to coordinate assistance with installation using the template in Appendix F.** Do not proceed until the issue is resolved.*

9. Close the **Rules Setup** window.

5 Set Up HL7 Connection to VBECS TEST



If you are using multiple Vision instruments they can all share one connection to VBECS.



Modifying connection settings in the Instrument Manager outside of this Setup Guide is not allowed and may lead to malfunction of the Automated Instrument to VBECS interface.

1. Navigate to **Configuration -> Connection Assignment**.
2. On the **Connection Assignment** window, click the **Add** button.
3. On the **Connection Properties** window, enter the **Connection Name**. Enter a name that contains **<HL7 Side Configuration>** and word **Connection**. For example:

HIN_VBECS_1_Connection

4. Select **<HL7 Side Configuration>** as **Configuration Name**.
5. Check **Include in Specimen Management** check box.
6. Select **TCP/IP** connection in Device. (Figure 31).

Figure 31: Example of Connection Properties Window

IM Connection Properties - HIN_VBECS_1_Connection

Connection Name:

Configuration Name:

Site:

Location:

☒ Start on System Start

Destination Line(s):

Number of Days to Keep:

Incoming Messages:

Outgoing Messages:

Communications Trace:

Error Messages:

Driver Data:

Advanced Options:

☐ Explode Coded Entries for this Connection

☒ Include in Specimen Management

Default user:

☐ Update Specimen Management on Status Messages

☐ Include in Specimen Storage and Retrieval

Device:

☐ NULL ☐ COM ☒ TCP/IP

Device Parameters ...


Close

7. Click on **Device Parameters** button.

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

8. Enter **TCP/IP Address** and **TCP/IP Port Number** that matches **VBECS TEST Application IP Address and IP Port Number** configured in **VBECS TEST Administrator** application for **Auto Instrument Interface** (Figure 32). Please refer to *VistA Blood Establishment Computer Software (VBECS) 2.2.0 Technical Manual Security Guide* for instruction on how to configure interfaces for VBECS.


Figure 32: Example of TCP/IP Configuration Window

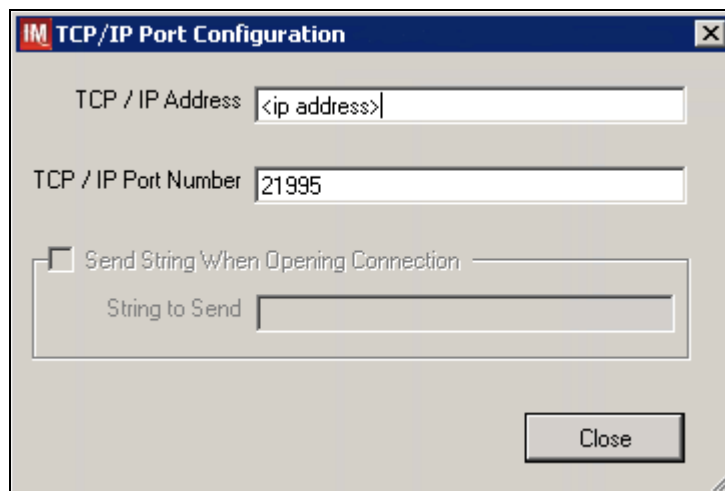

The screenshot shows a window titled "IM TCP/IP Port Configuration". It contains two text input fields: "TCP / IP Address" with the placeholder text "<ip address>" and "TCP / IP Port Number" with the value "21995". Below these fields is a checkbox labeled "Send String When Opening Connection". Under the checkbox is a text input field labeled "String to Send". A "Close" button is located at the bottom right of the window.

9. Click close on the **TCP/IP Port Configuration** window and click **Yes** to save.
10. Click close on the **Connection Properties** window and click **Yes** to save.
11. Verify the newly created connection shows on **Connection Assignment** window (Figure 33).


Figure 33: Example of Newly Created Connection Figure

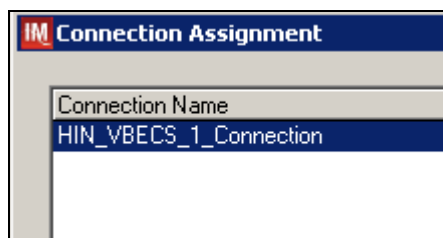


The screenshot shows a window titled "IM Connection Assignment". It contains a table with the following content:

Connection Name
HIN_VBECS_1_Connection

6 Set Up Instrument Connection



Execute instructions in this section for each instrument that will be connecting to VBECS.




Modifying connection settings in the Instrument Manager outside of this Setup Guide is not allowed and may lead to malfunction of the Automated Instrument to VBECS interface.

1. On the **Connection Assignment** window, click the **Add** button.
2. On the **Connection Properties** window, enter the **Connection Name**. Enter a name that contains *<Instrument Side Configuration>* and word **Connection**. For example:

HIN_Vision_1_Connection

3. Select *<Instrument Side Configuration>* as **Configuration Name**.
4. Check **Include in Specimen Management** check box.
5. Select **TCP/IP** connection.
6. On the **Destination Lines** list check the box next to connection configured in previous section (Figure 34).


Figure 34: Example of Connection Properties Window

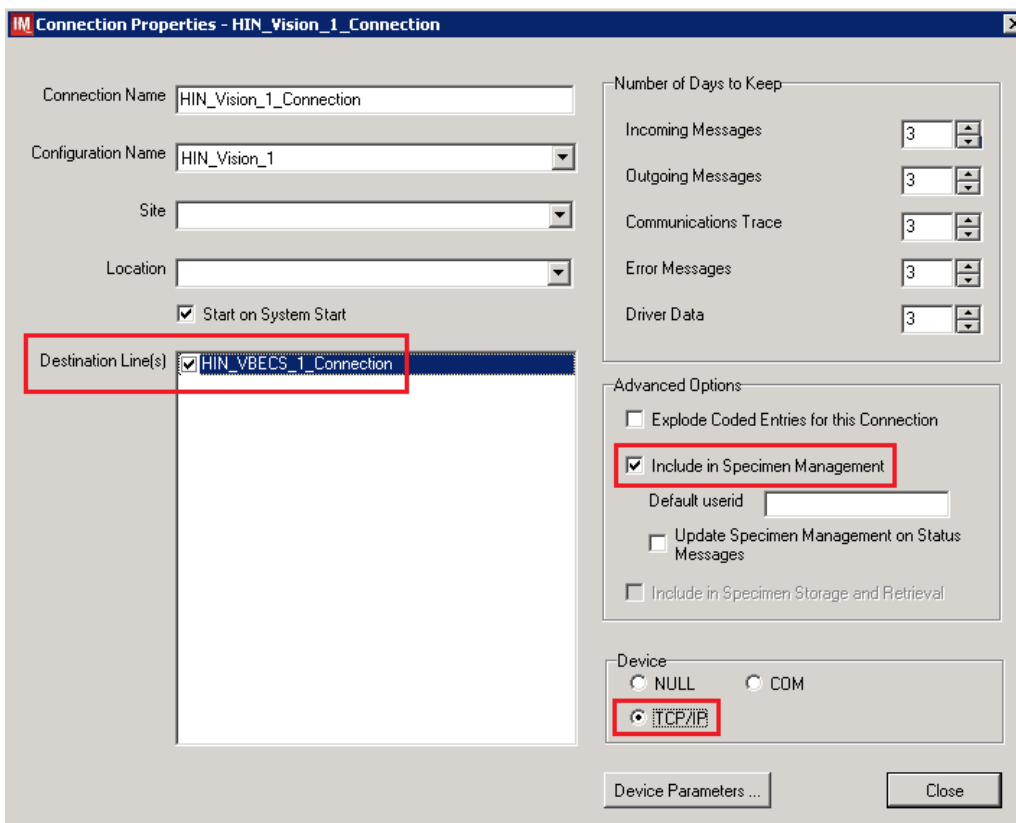

The screenshot shows the 'Connection Properties - HIN_Vision_1_Connection' window. The 'Connection Name' field is 'HIN_Vision_1_Connection'. The 'Configuration Name' dropdown is 'HIN_Vision_1'. The 'Site' and 'Location' dropdowns are empty. The 'Start on System Start' checkbox is checked. The 'Destination Line(s)' list has 'HIN_VBECS_1_Connection' selected. The 'Number of Days to Keep' section has 'Incoming Messages', 'Outgoing Messages', 'Communications Trace', 'Error Messages', and 'Driver Data' all set to 3. The 'Advanced Options' section has 'Include in Specimen Management' checked. The 'Device' section has 'TCP/IP' selected. The 'Close' button is at the bottom right.


7. Click on **Device Parameters** button.
8. **TCP/IP Address** remains blank.
9. Set **TCP/IP Port Number** to *<Instrument Port>* discussed in Section **Network Connectivity Setup**.
10. Uncheck **Send String When Opening Connection** check box (Figure 35).

Figure 35: Example of TCP/IP Configuration Window

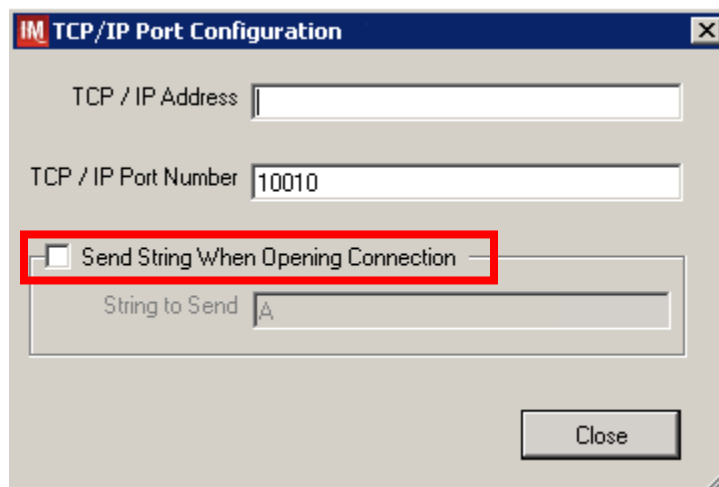

11. Click close on the **TCP/IP Port Configuration** window and click **Yes** to save.
12. Click close on the **Connection Properties** window and click **Yes** to save.
13. Verify the newly created connection shows on **Connection Assignment** window (Figure 36).
14. Close the **Connection Assignment** window.


Figure 36: Example of Newly Created Connection

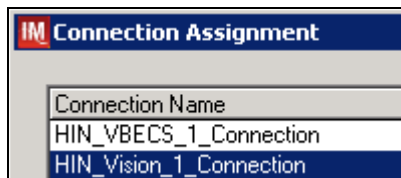

7 Test New Connections


1. Navigate to **System -> Status**.
2. Right-click each newly created connection and choose option to **Start Selected Connections** (Figure 37).

Figure 37: Example of Connection Status Window

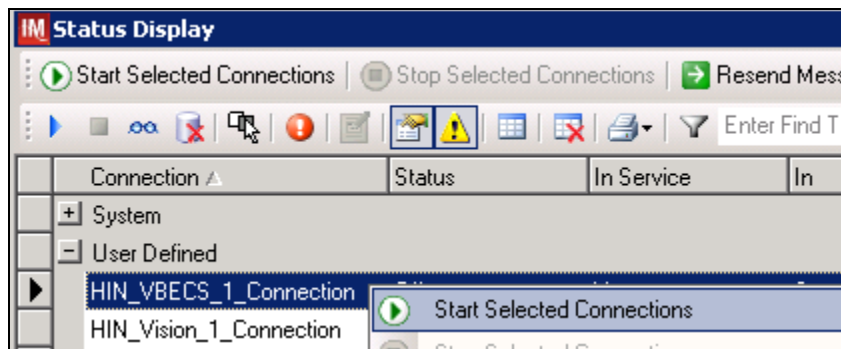

3. Verify that all newly created connections are showing **Status** of **On** after about a minute or so. (Figure 38).

Figure 38: Example of Successful Connection Test

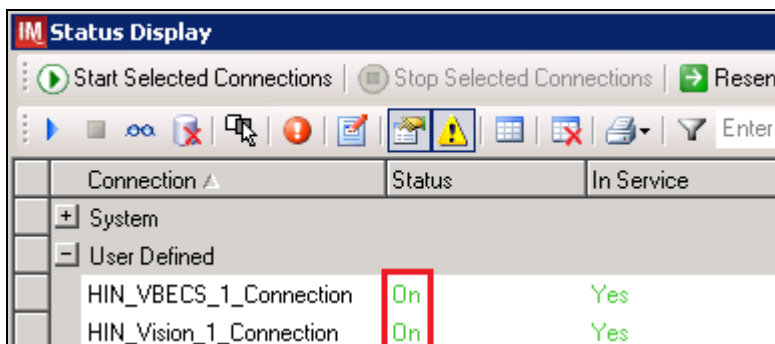


STOP If one or more connections fail to start, file a national CA SDM ticket to coordinate assistance with installation using template in Appendix F. Do not proceed until the issue is resolved.

4. Close the **Status Display**.

8 Validate Instrument connectivity to VBECS TEST

Execute validation instructions from *Appendix G* to verify that Instrument is properly interfacing with VBECS.




If one or more validation scenarios fail, file a national CA SDM ticket to coordinate assistance with installation using template in Appendix F. Do not proceed until the issue is resolved.

9 Set up HL7 Connection to VBECS PROD


1. Navigate to **System -> Status**.
2. Right-click <VBECS_Connection> and choose option to **Stop Selected Connections** (Figure 39).

Figure 39: Example of Stop Selected Connections

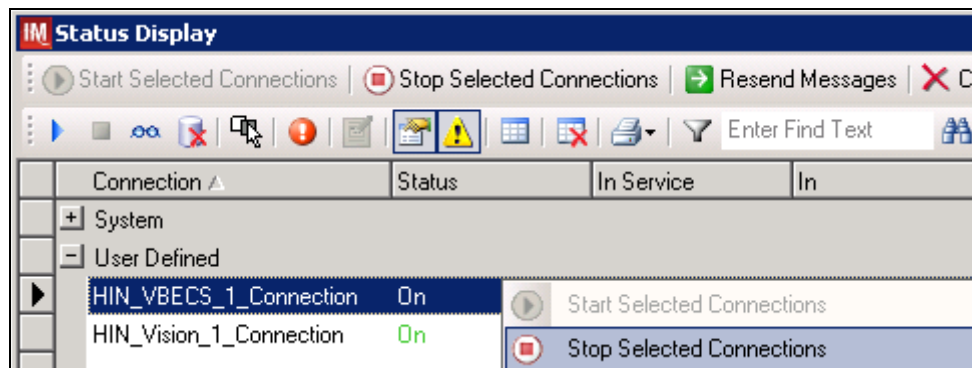

3. Wait until connection status changes to **Off**. Navigate to **Configuration -> Connection Assignment**.
4. Select <VBECS_Connection> and click **Properties** (Figure 40).

Figure 40: Example of Connection Assignment

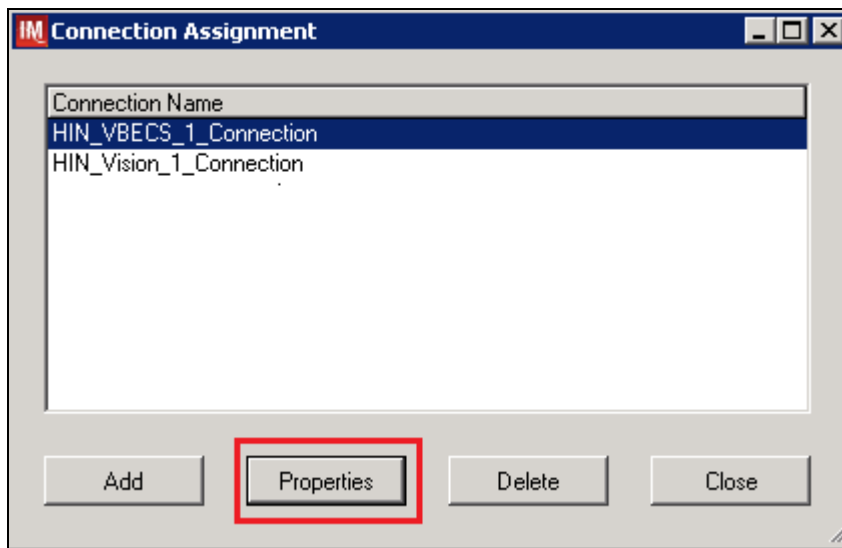

5. On **Connection Properties** window click **Device Parameters** (Figure 41).

Figure 41: Example of Connection Properties

Connection Name: HIN_VBECS_1_Connection

Configuration Name: HIN_VBECS_1

Site:

Location:

☒ Start on System Start

Destination Line(s):
☐ HIN_VBECS_1_Connection
☐ HIN_Vision_1_Connection

Number of Days to Keep:

Incoming Messages: 3

Outgoing Messages: 3

Communications Trace: 3

Error Messages: 3

Driver Data: 3

Advanced Options:

☐ Explode Coded Entries for this Connection

☒ Include in Specimen Management

Default userid:

☐ Update Specimen Management on Status Messages

☐ Include in Specimen Storage and Retrieval

Device:


☐ NULL ☐ COM ☒ TCP/IP

Device Parameters ...

Close


6. Enter **TCP/IP Address** and **TCP/IP Port Number** that matches **VBECS PROD Application IP Address and IP Port Number** configured in **VBECS PROD Administrator** application for **Auto Instrument Interface** (Figure 42). Please refer to *VistA Blood Establishment Computer Software (VBECS) Technical Manual Security Guide* for instruction on how to configure interfaces for VBECS.

Figure 42: Example of TCP/IP Configuration Window

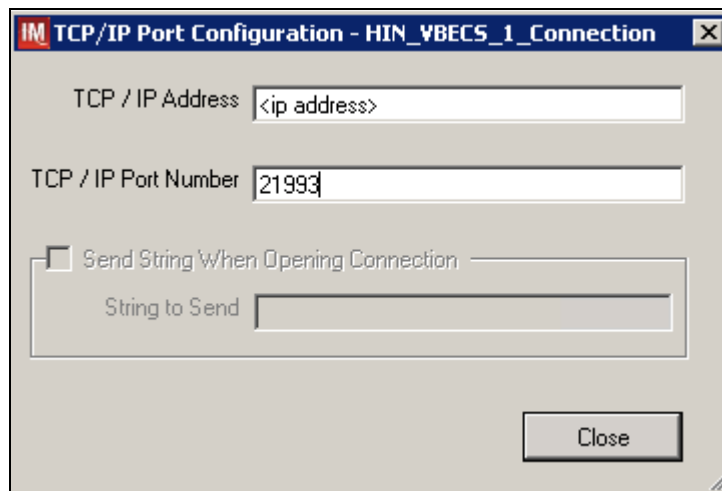


7. Close the **TCP/IP Port Configuration** and click **Yes** to confirm changes.
8. Close the **Connection Properties window** and click **Yes** to confirm changes.
9. Close the **Connection Assignment** window.
10. Navigate to **System -> Status**.
11. Right-click **<VBECS_Connection>** and choose option to **Start Selected Connections** (Figure 43).

Figure 43: Example of Connection Status Window

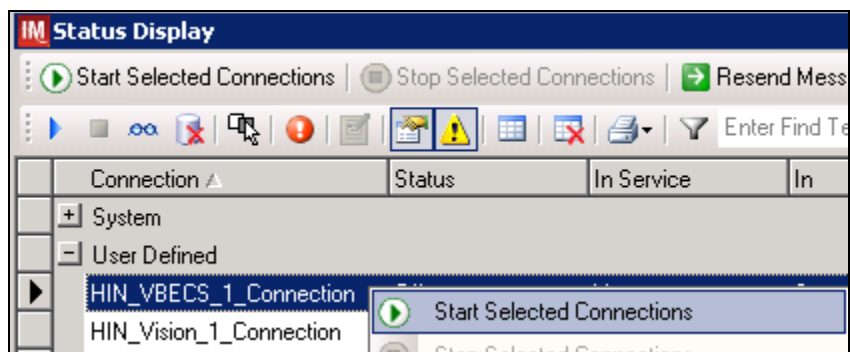

12. Verify that **<VBECS_Connection>** is showing Status of **On** after about a minute or so. (Figure 44).

Figure 44: Example of Successful Connection Test Figure 44

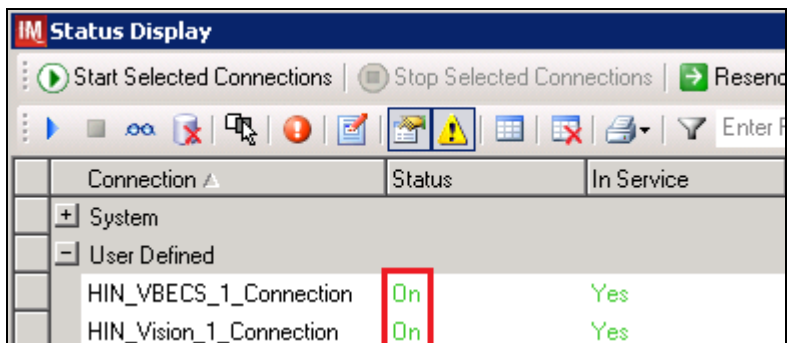

*If connection fails to start, **file a national CA SDM ticket to coordinate assistance with installation using template in Appendix F. Do not proceed until the issue is resolved.***

13. Close the **Status Display** and Log Off the system.

Glossary

Acronym, Term	Definition
Automated Instrument	Blood Bank Analyzer that performs blood testing.
Division Code	Also known as Station Number in Vista is the unique alphanumeric code that is associated with each hospital (e.g. 589 for VA Heartland West).
Instrument Manager (IM)	Software created by Data Innovations that serves as a middleware between Automated Instrument and VBECS. It translates messages containing test results sent from an instrument into HL7 messages that are then parsed into VBECS.
VA	Department of Veterans Affairs.
VBECS	VistA Blood Establishment Computer Software.
VistA	Veterans Health Information Systems and Technology Architecture. VistA software, developed by the VA, is used to support clinical and administrative functions at VA Medical Centers nationwide. VistA is composed of packages that undergo a verification process to ensure conformity with name spacing and other VistA standards and conventions.

Appendices

Appendix A: Instrument Side Mapping

Table 1- Appendix A: Instrument Side mapping

Instrument Test Code	IM Test Code	Fluid	Display Name
VBECS ABDRev-ABO	VBECS ABDRev-ABO	CENTBLOOD	ABOInterp
VBECS ABDRev-ABO	VBECS ABDRev-ABO	PACKEDCELLS	ABOInterp
VBECS ABDRev-ABO-A1-Cells	VBECS ABDRev-ABO-A1-Cells	CENTBLOOD	A1Cells
VBECS ABDRev-ABO-A1-Cells	VBECS ABDRev-ABO-A1-Cells	PACKEDCELLS	A1Cells
VBECS ABDRev-ABO-Anti-A	VBECS ABDRev-ABO-Anti-A	CENTBLOOD	AntiA
VBECS ABDRev-ABO-Anti-A	VBECS ABDRev-ABO-Anti-A	PACKEDCELLS	AntiA
VBECS ABDRev-ABO-Anti-B	VBECS ABDRev-ABO-Anti-B	CENTBLOOD	AntiB
VBECS ABDRev-ABO-Anti-B	VBECS ABDRev-ABO-Anti-B	PACKEDCELLS	AntiB
VBECS ABDRev-ABO-B-Cells	VBECS ABDRev-ABO-B-Cells	CENTBLOOD	BCells
VBECS ABDRev-ABO-B-Cells	VBECS ABDRev-ABO-B-Cells	PACKEDCELLS	BCells
VBECS ABDRev-Rh	VBECS ABDRev-Rh	CENTBLOOD	RhInterp
VBECS ABDRev-Rh	VBECS ABDRev-Rh	PACKEDCELLS	RhInterp
VBECS ABDRev-Rh-Anti-D	VBECS ABDRev-Rh-Anti-D	CENTBLOOD	AntiD
VBECS ABDRev-Rh-Anti-D	VBECS ABDRev-Rh-Anti-D	PACKEDCELLS	AntiD
VBECS ABDRev-Rh-Ctrl	VBECS ABDRev-Rh-Ctrl	CENTBLOOD	DControl
VBECS ABDRev-Rh-Ctrl	VBECS ABDRev-Rh-Ctrl	PACKEDCELLS	DControl
VBECS ABS 2 Cell-ABScr	VBECS ABS 2 Cell-ABScr	CENTBLOOD	ABSInterp
VBECS ABS 2 Cell-ABScr	VBECS ABS 2 Cell-ABScr	PACKEDCELLS	ABSInterp
VBECS ABS 2 Cell-ABScr-0.8-Sel I	VBECS ABS 2 Cell-ABScr-0.8-Sel I	CENTBLOOD	SC1
VBECS ABS 2 Cell-ABScr-0.8-Sel I	VBECS ABS 2 Cell-ABScr-0.8-Sel I	PACKEDCELLS	SC1
VBECS ABS 2 Cell-ABScr-0.8-Sel II	VBECS ABS 2 Cell-ABScr-0.8-Sel II	CENTBLOOD	SC2
VBECS ABS 2 Cell-ABScr-0.8-Sel II	VBECS ABS 2 Cell-ABScr-0.8-Sel II	PACKEDCELLS	SC2
VBECS ABS 3 Cell-ABScr	VBECS ABS 3 Cell-ABScr	CENTBLOOD	ABSInterp
VBECS ABS 3 Cell-ABScr	VBECS ABS 3 Cell-ABScr	PACKEDCELLS	ABSInterp
VBECS ABS 3 Cell-ABScr-0.8-Surg1	VBECS ABS 3 Cell-ABScr-0.8-Surg1	CENTBLOOD	SC1
VBECS ABS 3 Cell-ABScr-0.8-Surg1	VBECS ABS 3 Cell-ABScr-0.8-Surg1	PACKEDCELLS	SC1
VBECS ABS 3 Cell-ABScr-0.8-Surg2	VBECS ABS 3 Cell-ABScr-0.8-Surg2	CENTBLOOD	SC2
VBECS ABS 3 Cell-ABScr-0.8-Surg2	VBECS ABS 3 Cell-ABScr-0.8-Surg2	PACKEDCELLS	SC2
VBECS ABS 3 Cell-ABScr-0.8-Surg3	VBECS ABS 3 Cell-ABScr-0.8-Surg3	CENTBLOOD	SC3
VBECS ABS 3 Cell-ABScr-0.8-Surg3	VBECS ABS 3 Cell-ABScr-0.8-Surg3	PACKEDCELLS	SC3
VBECS Confirm AB-ABO	VBECS Confirm AB-ABO	CENTBLOOD	ABOInterp
VBECS Confirm AB-ABO	VBECS Confirm AB-ABO	PACKEDCELLS	ABOInterp
VBECS Confirm AB-ABO-Anti-A	VBECS Confirm AB-ABO-Anti-A	CENTBLOOD	AntiA
VBECS Confirm AB-ABO-Anti-A	VBECS Confirm AB-ABO-Anti-A	PACKEDCELLS	AntiA

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Instrument Test Code	IM Test Code	Fluid	Display Name
VBECS Confirm AB-ABO-Anti-B	VBECS Confirm AB-ABO-Anti-B	CENTBLOOD	AntiB
VBECS Confirm AB-ABO-Anti-B	VBECS Confirm AB-ABO-Anti-B	PACKEDCELLS	AntiB
VBECS Confirm ABD-ABO	VBECS Confirm ABD-ABO	CENTBLOOD	ABOInterp
VBECS Confirm ABD-ABO	VBECS Confirm ABD-ABO	PACKEDCELLS	ABOInterp
VBECS Confirm ABD-ABO-Anti-A	VBECS Confirm ABD-ABO-Anti-A	CENTBLOOD	AntiA
VBECS Confirm ABD-ABO-Anti-A	VBECS Confirm ABD-ABO-Anti-A	PACKEDCELLS	AntiA
VBECS Confirm ABD-ABO-Anti-B	VBECS Confirm ABD-ABO-Anti-B	CENTBLOOD	AntiB
VBECS Confirm ABD-ABO-Anti-B	VBECS Confirm ABD-ABO-Anti-B	PACKEDCELLS	AntiB
VBECS Confirm ABD-Rh	VBECS Confirm ABD-Rh	CENTBLOOD	RhInterp
VBECS Confirm ABD-Rh	VBECS Confirm ABD-Rh	PACKEDCELLS	RhInterp
VBECS Confirm ABD-Rh-Anti-D	VBECS Confirm ABD-Rh-Anti-D	CENTBLOOD	AntiD
VBECS Confirm ABD-Rh-Anti-D	VBECS Confirm ABD-Rh-Anti-D	PACKEDCELLS	AntiD
VBECS Confirm O Neg-Anti-AB	VBECS Confirm O Neg-Anti-AB	CENTBLOOD	ABOInterp
VBECS Confirm O Neg-Anti-AB	VBECS Confirm O Neg-Anti-AB	PACKEDCELLS	ABOInterp
VBECS Confirm O Neg-Anti-AB-Anti-A+B	VBECS Confirm O Neg-Anti-AB-Anti-A+B	CENTBLOOD	AntiAB
VBECS Confirm O Neg-Anti-AB-Anti-A+B	VBECS Confirm O Neg-Anti-AB-Anti-A+B	PACKEDCELLS	AntiAB
VBECS Confirm O Neg-Rh	VBECS Confirm O Neg-Rh	CENTBLOOD	RhInterp
VBECS Confirm O Neg-Rh	VBECS Confirm O Neg-Rh	PACKEDCELLS	RhInterp
VBECS Confirm O Neg-Rh-Anti-D	VBECS Confirm O Neg-Rh-Anti-D	CENTBLOOD	AntiD
VBECS Confirm O Neg-Rh-Anti-D	VBECS Confirm O Neg-Rh-Anti-D	PACKEDCELLS	AntiD
VBECS Confirm O Pos-Anti-AB	VBECS Confirm O Pos-Anti-AB	CENTBLOOD	ABOInterp
VBECS Confirm O Pos-Anti-AB	VBECS Confirm O Pos-Anti-AB	PACKEDCELLS	ABOInterp
VBECS Confirm O Pos-Anti-AB-Anti-A+B	VBECS Confirm O Pos-Anti-AB-Anti-A+B	CENTBLOOD	AntiAB
VBECS Confirm O Pos-Anti-AB-Anti-A+B	VBECS Confirm O Pos-Anti-AB-Anti-A+B	PACKEDCELLS	AntiAB
VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info	VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info	CENTBLOOD	AHG
VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info	VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info	PACKEDCELLS	AHG
VBECS Crossmatch IAT-XM	VBECS Crossmatch IAT-XM	CENTBLOOD	XMInterp
VBECS Crossmatch IAT-XM	VBECS Crossmatch IAT-XM	PACKEDCELLS	XMInterp
VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info	VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info	CENTBLOOD	AHG
VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info	VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info	PACKEDCELLS	AHG
VBECS Crossmatch-Buffered Gel-Donor Info	VBECS Crossmatch-Buffered Gel-Donor Info	CENTBLOOD	IS
VBECS Crossmatch-Buffered Gel-Donor Info	VBECS Crossmatch-Buffered Gel-Donor Info	PACKEDCELLS	IS
VBECS Crossmatch-XM	VBECS Crossmatch-XM	CENTBLOOD	XMInterp
VBECS Crossmatch-XM	VBECS Crossmatch-XM	PACKEDCELLS	XMInterp
VBECS DAT IgG-IgG	VBECS DAT IgG-IgG	CENTBLOOD	IgGInterp
VBECS DAT IgG-IgG	VBECS DAT IgG-IgG	PACKEDCELLS	IgGInterp
VBECS DAT IgG-IgG-IgG	VBECS DAT IgG-IgG-IgG	CENTBLOOD	IgG

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Instrument Test Code	IM Test Code	Fluid	Display Name
VBECS DAT IgG-IgG-IgG	VBECS DAT IgG-IgG-IgG	PACKEDCELLS	IgG
VBECS DAT Poly+IgG-IgG	VBECS DAT Poly+IgG-IgG	CENTBLOOD	IgGInterp
VBECS DAT Poly+IgG-IgG	VBECS DAT Poly+IgG-IgG	PACKEDCELLS	IgGInterp
VBECS DAT Poly+IgG-IgG-IgG	VBECS DAT Poly+IgG-IgG-IgG	CENTBLOOD	IgG
VBECS DAT Poly+IgG-IgG-IgG	VBECS DAT Poly+IgG-IgG-IgG	PACKEDCELLS	IgG
VBECS DAT Poly+IgG-Poly	VBECS DAT Poly+IgG-Poly	CENTBLOOD	DATInterp
VBECS DAT Poly+IgG-Poly	VBECS DAT Poly+IgG-Poly	PACKEDCELLS	DATInterp
VBECS DAT Poly+IgG-Poly-Poly	VBECS DAT Poly+IgG-Poly-Poly	CENTBLOOD	AHG
VBECS DAT Poly+IgG-Poly-Poly	VBECS DAT Poly+IgG-Poly-Poly	PACKEDCELLS	AHG
VBECS DAT-Poly	VBECS DAT-Poly	CENTBLOOD	DATInterp
VBECS DAT-Poly	VBECS DAT-Poly	PACKEDCELLS	DATInterp
VBECS DAT-Poly-Poly	VBECS DAT-Poly-Poly	CENTBLOOD	AHG
VBECS DAT-Poly-Poly	VBECS DAT-Poly-Poly	PACKEDCELLS	AHG
VBECS Phenotyping C-Anti-C (RH2)	VBECS Phenotyping C-Anti-C (RH2)	CENTBLOOD	CInterp
VBECS Phenotyping C-Anti-C (RH2)	VBECS Phenotyping C-Anti-C (RH2)	PACKEDCELLS	CInterp
VBECS Phenotyping C-Anti-C (RH2)-Anti-C (RH2)	VBECS Phenotyping C-Anti-C (RH2)-Anti-C (RH2)	CENTBLOOD	Anti-C
VBECS Phenotyping C-Anti-C (RH2)-Anti-C (RH2)	VBECS Phenotyping C-Anti-C (RH2)-Anti-C (RH2)	PACKEDCELLS	Anti-C
VBECS Phenotyping E-Anti-E (RH3)	VBECS Phenotyping E-Anti-E (RH3)	CENTBLOOD	EInterp
VBECS Phenotyping E-Anti-E (RH3)	VBECS Phenotyping E-Anti-E (RH3)	PACKEDCELLS	EInterp
VBECS Phenotyping E-Anti-E (RH3)-Anti-E (RH3)	VBECS Phenotyping E-Anti-E (RH3)-Anti-E (RH3)	CENTBLOOD	Anti-E
VBECS Phenotyping E-Anti-E (RH3)-Anti-E (RH3)	VBECS Phenotyping E-Anti-E (RH3)-Anti-E (RH3)	PACKEDCELLS	Anti-E
VBECS Phenotyping smc-Anti-c (RH4)	VBECS Phenotyping smc-Anti-c (RH4)	CENTBLOOD	cInterp
VBECS Phenotyping smc-Anti-c (RH4)	VBECS Phenotyping smc-Anti-c (RH4)	PACKEDCELLS	cInterp
VBECS Phenotyping smc-Anti-c (RH4)-Anti-c (RH4)	VBECS Phenotyping smc-Anti-c (RH4)-Anti-c (RH4)	CENTBLOOD	Anti-c
VBECS Phenotyping smc-Anti-c (RH4)-Anti-c (RH4)	VBECS Phenotyping smc-Anti-c (RH4)-Anti-c (RH4)	PACKEDCELLS	Anti-c
VBECS Phenotyping sme-Anti-e (RH5)	VBECS Phenotyping sme-Anti-e (RH5)	CENTBLOOD	eInterp
VBECS Phenotyping sme-Anti-e (RH5)	VBECS Phenotyping sme-Anti-e (RH5)	PACKEDCELLS	eInterp
VBECS Phenotyping sme-Anti-e (RH5)-Anti-e (RH5)	VBECS Phenotyping sme-Anti-e (RH5)-Anti-e (RH5)	CENTBLOOD	Anti-e
VBECS Phenotyping sme-Anti-e (RH5)-Anti-e (RH5)	VBECS Phenotyping sme-Anti-e (RH5)-Anti-e (RH5)	PACKEDCELLS	Anti-e
VBECS TAS 2 Cell-ABO	VBECS TAS 2 Cell-ABO	CENTBLOOD	ABOInterp
VBECS TAS 2 Cell-ABO	VBECS TAS 2 Cell-ABO	PACKEDCELLS	ABOInterp
VBECS TAS 2 Cell-ABO-A1-Cells	VBECS TAS 2 Cell-ABO-A1-Cells	CENTBLOOD	A1Cells
VBECS TAS 2 Cell-ABO-A1-Cells	VBECS TAS 2 Cell-ABO-A1-Cells	PACKEDCELLS	A1Cells
VBECS TAS 2 Cell-ABO-Anti-A	VBECS TAS 2 Cell-ABO-Anti-A	CENTBLOOD	AntiA
VBECS TAS 2 Cell-ABO-Anti-A	VBECS TAS 2 Cell-ABO-Anti-A	PACKEDCELLS	AntiA
VBECS TAS 2 Cell-ABO-Anti-B	VBECS TAS 2 Cell-ABO-Anti-B	CENTBLOOD	AntiB
VBECS TAS 2 Cell-ABO-Anti-B	VBECS TAS 2 Cell-ABO-Anti-B	PACKEDCELLS	AntiB

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Instrument Test Code	IM Test Code	Fluid	Display Name
VBECs TAS 2 Cell-ABO-B-Cells	VBECs TAS 2 Cell-ABO-B-Cells	CENTBLOOD	BCells
VBECs TAS 2 Cell-ABO-B-Cells	VBECs TAS 2 Cell-ABO-B-Cells	PACKEDCELLS	BCells
VBECs TAS 2 Cell-ABScr	VBECs TAS 2 Cell-ABScr	CENTBLOOD	ABSInterp
VBECs TAS 2 Cell-ABScr	VBECs TAS 2 Cell-ABScr	PACKEDCELLS	ABSInterp
VBECs TAS 2 Cell-ABScr-0.8-Sel I	VBECs TAS 2 Cell-ABScr-0.8-Sel I	CENTBLOOD	SC1
VBECs TAS 2 Cell-ABScr-0.8-Sel I	VBECs TAS 2 Cell-ABScr-0.8-Sel I	PACKEDCELLS	SC1
VBECs TAS 2 Cell-ABScr-0.8-Sel II	VBECs TAS 2 Cell-ABScr-0.8-Sel II	CENTBLOOD	SC2
VBECs TAS 2 Cell-ABScr-0.8-Sel II	VBECs TAS 2 Cell-ABScr-0.8-Sel II	PACKEDCELLS	SC2
VBECs TAS 2 Cell-Rh	VBECs TAS 2 Cell-Rh	CENTBLOOD	RhInterp
VBECs TAS 2 Cell-Rh	VBECs TAS 2 Cell-Rh	PACKEDCELLS	RhInterp
VBECs TAS 2 Cell-Rh-Anti-D	VBECs TAS 2 Cell-Rh-Anti-D	CENTBLOOD	AntiD
VBECs TAS 2 Cell-Rh-Anti-D	VBECs TAS 2 Cell-Rh-Anti-D	PACKEDCELLS	AntiD
VBECs TAS 2 Cell-Rh-Ctrl	VBECs TAS 2 Cell-Rh-Ctrl	CENTBLOOD	DControl
VBECs TAS 2 Cell-Rh-Ctrl	VBECs TAS 2 Cell-Rh-Ctrl	PACKEDCELLS	DControl
VBECs TAS 3 Cell-ABO	VBECs TAS 3 Cell-ABO	CENTBLOOD	ABOInterp
VBECs TAS 3 Cell-ABO	VBECs TAS 3 Cell-ABO	PACKEDCELLS	ABOInterp
VBECs TAS 3 Cell-ABO-A1-Cells	VBECs TAS 3 Cell-ABO-A1-Cells	CENTBLOOD	A1Cells
VBECs TAS 3 Cell-ABO-A1-Cells	VBECs TAS 3 Cell-ABO-A1-Cells	PACKEDCELLS	A1Cells
VBECs TAS 3 Cell-ABO-Anti-A	VBECs TAS 3 Cell-ABO-Anti-A	CENTBLOOD	AntiA
VBECs TAS 3 Cell-ABO-Anti-A	VBECs TAS 3 Cell-ABO-Anti-A	PACKEDCELLS	AntiA
VBECs TAS 3 Cell-ABO-Anti-B	VBECs TAS 3 Cell-ABO-Anti-B	CENTBLOOD	AntiB
VBECs TAS 3 Cell-ABO-Anti-B	VBECs TAS 3 Cell-ABO-Anti-B	PACKEDCELLS	AntiB
VBECs TAS 3 Cell-ABO-B-Cells	VBECs TAS 3 Cell-ABO-B-Cells	CENTBLOOD	BCells
VBECs TAS 3 Cell-ABO-B-Cells	VBECs TAS 3 Cell-ABO-B-Cells	PACKEDCELLS	BCells
VBECs TAS 3 Cell-ABScr	VBECs TAS 3 Cell-ABScr	CENTBLOOD	ABSInterp
VBECs TAS 3 Cell-ABScr	VBECs TAS 3 Cell-ABScr	PACKEDCELLS	ABSInterp
VBECs TAS 3 Cell-ABScr-0.8-Surg1	VBECs TAS 3 Cell-ABScr-0.8-Surg1	CENTBLOOD	SC1
VBECs TAS 3 Cell-ABScr-0.8-Surg1	VBECs TAS 3 Cell-ABScr-0.8-Surg1	PACKEDCELLS	SC1
VBECs TAS 3 Cell-ABScr-0.8-Surg2	VBECs TAS 3 Cell-ABScr-0.8-Surg2	CENTBLOOD	SC2
VBECs TAS 3 Cell-ABScr-0.8-Surg2	VBECs TAS 3 Cell-ABScr-0.8-Surg2	PACKEDCELLS	SC2
VBECs TAS 3 Cell-ABScr-0.8-Surg3	VBECs TAS 3 Cell-ABScr-0.8-Surg3	CENTBLOOD	SC3
VBECs TAS 3 Cell-ABScr-0.8-Surg3	VBECs TAS 3 Cell-ABScr-0.8-Surg3	PACKEDCELLS	SC3
VBECs TAS 3 Cell-Rh	VBECs TAS 3 Cell-Rh	CENTBLOOD	RhInterp
VBECs TAS 3 Cell-Rh	VBECs TAS 3 Cell-Rh	PACKEDCELLS	RhInterp
VBECs TAS 3 Cell-Rh-Anti-D	VBECs TAS 3 Cell-Rh-Anti-D	CENTBLOOD	AntiD
VBECs TAS 3 Cell-Rh-Anti-D	VBECs TAS 3 Cell-Rh-Anti-D	PACKEDCELLS	AntiD
VBECs TAS 3 Cell-Rh-Ctrl	VBECs TAS 3 Cell-Rh-Ctrl	CENTBLOOD	DControl
VBECs TAS 3 Cell-Rh-Ctrl	VBECs TAS 3 Cell-Rh-Ctrl	PACKEDCELLS	DControl

Appendix B: HL7 (VBECS) Side Mapping

Table 2- Appendix B: HL7 (VBECS) Side Mapping

Instrument Test Code	IM Test Code	Fluid
AHGInterp	AHGInterp	CENTBLOOD
AHGInterp	AHGInterp	PACKEDCELLS
ISInterp	ISInterp	CENTBLOOD
ISInterp	ISInterp	PACKEDCELLS
VBECS ABDRev-ABO	ABOInterp	CENTBLOOD
VBECS ABDRev-ABO	ABOInterp	PACKEDCELLS
VBECS ABDRev-ABO-A1-Cells	A1Cells	CENTBLOOD
VBECS ABDRev-ABO-A1-Cells	A1Cells	PACKEDCELLS
VBECS ABDRev-ABO-Anti-A	AntiA	CENTBLOOD
VBECS ABDRev-ABO-Anti-A	AntiA	PACKEDCELLS
VBECS ABDRev-ABO-Anti-B	AntiB	CENTBLOOD
VBECS ABDRev-ABO-Anti-B	AntiB	PACKEDCELLS
VBECS ABDRev-ABO-B-Cells	BCells	CENTBLOOD
VBECS ABDRev-ABO-B-Cells	BCells	PACKEDCELLS
VBECS ABDRev-Rh	RhInterp	CENTBLOOD
VBECS ABDRev-Rh	RhInterp	PACKEDCELLS
VBECS ABDRev-Rh-Anti-D	AntiD	CENTBLOOD
VBECS ABDRev-Rh-Anti-D	AntiD	PACKEDCELLS
VBECS ABDRev-Rh-Ctrl	DControl	CENTBLOOD
VBECS ABDRev-Rh-Ctrl	DControl	PACKEDCELLS
VBECS ABS 2 Cell-ABScr	ABSInterp	CENTBLOOD
VBECS ABS 2 Cell-ABScr	ABSInterp	PACKEDCELLS
VBECS ABS 2 Cell-ABScr-0.8-Sel I	SC1	CENTBLOOD
VBECS ABS 2 Cell-ABScr-0.8-Sel I	SC1	PACKEDCELLS
VBECS ABS 2 Cell-ABScr-0.8-Sel II	SC2	CENTBLOOD
VBECS ABS 2 Cell-ABScr-0.8-Sel II	SC2	PACKEDCELLS
VBECS ABS 3 Cell-ABScr	ABSInterp	CENTBLOOD
VBECS ABS 3 Cell-ABScr	ABSInterp	PACKEDCELLS
VBECS ABS 3 Cell-ABScr-0.8-Surg1	SC1	CENTBLOOD
VBECS ABS 3 Cell-ABScr-0.8-Surg1	SC1	PACKEDCELLS
VBECS ABS 3 Cell-ABScr-0.8-Surg2	SC2	CENTBLOOD
VBECS ABS 3 Cell-ABScr-0.8-Surg2	SC2	PACKEDCELLS
VBECS ABS 3 Cell-ABScr-0.8-Surg3	SC3	CENTBLOOD
VBECS ABS 3 Cell-ABScr-0.8-Surg3	SC3	PACKEDCELLS
VBECS Confirm AB-ABO	ABOInterp	CENTBLOOD
VBECS Confirm AB-ABO	ABOInterp	PACKEDCELLS

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Instrument Test Code	IM Test Code	Fluid
VBECS Confirm AB-ABO-Anti-A	AntiA	CENTBLOOD
VBECS Confirm AB-ABO-Anti-A	AntiA	PACKEDCELLS
VBECS Confirm AB-ABO-Anti-B	AntiB	CENTBLOOD
VBECS Confirm AB-ABO-Anti-B	AntiB	PACKEDCELLS
VBECS Confirm ABD-ABO	ABOInterp	CENTBLOOD
VBECS Confirm ABD-ABO	ABOInterp	PACKEDCELLS
VBECS Confirm ABD-ABO-Anti-A	AntiA	CENTBLOOD
VBECS Confirm ABD-ABO-Anti-A	AntiA	PACKEDCELLS
VBECS Confirm ABD-ABO-Anti-B	AntiB	CENTBLOOD
VBECS Confirm ABD-ABO-Anti-B	AntiB	PACKEDCELLS
VBECS Confirm ABD-Rh	RhInterp	CENTBLOOD
VBECS Confirm ABD-Rh	RhInterp	PACKEDCELLS
VBECS Confirm ABD-Rh-Anti-D	AntiD	CENTBLOOD
VBECS Confirm ABD-Rh-Anti-D	AntiD	PACKEDCELLS
VBECS Confirm O Neg-Anti-AB	ABOInterp	CENTBLOOD
VBECS Confirm O Neg-Anti-AB	ABOInterp	PACKEDCELLS
VBECS Confirm O Neg-Anti-AB-Anti-A+B	AntiAB	CENTBLOOD
VBECS Confirm O Neg-Anti-AB-Anti-A+B	AntiAB	PACKEDCELLS
VBECS Confirm O Neg-Rh	RhInterp	CENTBLOOD
VBECS Confirm O Neg-Rh	RhInterp	PACKEDCELLS
VBECS Confirm O Neg-Rh-Anti-D	AntiD	CENTBLOOD
VBECS Confirm O Neg-Rh-Anti-D	AntiD	PACKEDCELLS
VBECS Confirm O Pos-Anti-AB	ABOInterp	CENTBLOOD
VBECS Confirm O Pos-Anti-AB	ABOInterp	PACKEDCELLS
VBECS Confirm O Pos-Anti-AB-Anti-A+B	AntiAB	CENTBLOOD
VBECS Confirm O Pos-Anti-AB-Anti-A+B	AntiAB	PACKEDCELLS
VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info	AHG	CENTBLOOD
VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info	AHG	PACKEDCELLS
VBECS Crossmatch IAT-XM	XMInterp	CENTBLOOD
VBECS Crossmatch IAT-XM	XMInterp	PACKEDCELLS
VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info	AHG	CENTBLOOD
VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info	AHG	PACKEDCELLS
VBECS Crossmatch-Buffered Gel-Donor Info	IS	CENTBLOOD
VBECS Crossmatch-Buffered Gel-Donor Info	IS	PACKEDCELLS
VBECS DAT IgG-IgG	IgGInterp	CENTBLOOD
VBECS DAT IgG-IgG	IgGInterp	PACKEDCELLS
VBECS DAT IgG-IgG-IgG	IgG	CENTBLOOD
VBECS DAT IgG-IgG-IgG	IgG	PACKEDCELLS
VBECS DAT Poly+IgG-IgG	IgGInterp	CENTBLOOD
VBECS DAT Poly+IgG-IgG	IgGInterp	PACKEDCELLS

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Instrument Test Code	IM Test Code	Fluid
VBECS DAT Poly+IgG-IgG-IgG	IgG	CENTBLOOD
VBECS DAT Poly+IgG-IgG-IgG	IgG	PACKEDCELLS
VBECS DAT Poly+IgG-Poly	DATInterp	CENTBLOOD
VBECS DAT Poly+IgG-Poly	DATInterp	PACKEDCELLS
VBECS DAT Poly+IgG-Poly-Poly	AHG	CENTBLOOD
VBECS DAT Poly+IgG-Poly-Poly	AHG	PACKEDCELLS
VBECS DAT-Poly	DATInterp	CENTBLOOD
VBECS DAT-Poly	DATInterp	PACKEDCELLS
VBECS DAT-Poly-Poly	AHG	CENTBLOOD
VBECS DAT-Poly-Poly	AHG	PACKEDCELLS
VBECS Phenotyping C-Anti-C (RH2)-Anti-C (RH2)	AntiC	CENTBLOOD
VBECS Phenotyping C-Anti-C (RH2)-Anti-C (RH2)	AntiC	PACKEDCELLS
VBECS Phenotyping E-Anti-E (RH3)-Anti-E (RH3)	AntiE	CENTBLOOD
VBECS Phenotyping E-Anti-E (RH3)-Anti-E (RH3)	AntiE	PACKEDCELLS
VBECS Phenotyping smc-Anti-c (RH4)-Anti-c (RH4)	Antic	CENTBLOOD
VBECS Phenotyping smc-Anti-c (RH4)-Anti-c (RH4)	Antic	PACKEDCELLS
VBECS Phenotyping sme-Anti-e (RH5)-Anti-e (RH5)	Antie	CENTBLOOD
VBECS Phenotyping sme-Anti-e (RH5)-Anti-e (RH5)	Antie	PACKEDCELLS
VBECS TAS 2 Cell-ABO	ABOInterp	CENTBLOOD
VBECS TAS 2 Cell-ABO	ABOInterp	PACKEDCELLS
VBECS TAS 2 Cell-ABO-A1-Cells	A1Cells	CENTBLOOD
VBECS TAS 2 Cell-ABO-A1-Cells	A1Cells	PACKEDCELLS
VBECS TAS 2 Cell-ABO-Anti-A	AntiA	CENTBLOOD
VBECS TAS 2 Cell-ABO-Anti-A	AntiA	PACKEDCELLS
VBECS TAS 2 Cell-ABO-Anti-B	AntiB	CENTBLOOD
VBECS TAS 2 Cell-ABO-Anti-B	AntiB	PACKEDCELLS
VBECS TAS 2 Cell-ABO-B-Cells	BCells	CENTBLOOD
VBECS TAS 2 Cell-ABO-B-Cells	BCells	PACKEDCELLS
VBECS TAS 2 Cell-ABScr	ABScrInterp	CENTBLOOD
VBECS TAS 2 Cell-ABScr	ABScrInterp	PACKEDCELLS
VBECS TAS 2 Cell-ABScr-0.8-Sel I	SC1	CENTBLOOD
VBECS TAS 2 Cell-ABScr-0.8-Sel I	SC1	PACKEDCELLS
VBECS TAS 2 Cell-ABScr-0.8-Sel II	SC2	CENTBLOOD
VBECS TAS 2 Cell-ABScr-0.8-Sel II	SC2	PACKEDCELLS
VBECS TAS 2 Cell-Rh	RhInterp	CENTBLOOD
VBECS TAS 2 Cell-Rh	RhInterp	PACKEDCELLS
VBECS TAS 2 Cell-Rh-Anti-D	AntiD	CENTBLOOD
VBECS TAS 2 Cell-Rh-Anti-D	AntiD	PACKEDCELLS
VBECS TAS 2 Cell-Rh-Ctrl	DControl	CENTBLOOD
VBECS TAS 2 Cell-Rh-Ctrl	DControl	PACKEDCELLS

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Instrument Test Code	IM Test Code	Fluid
VBECS TAS 3 Cell-ABO	ABOInterp	CENTBLOOD
VBECS TAS 3 Cell-ABO	ABOInterp	PACKEDCELLS
VBECS TAS 3 Cell-ABO-A1-Cells	A1Cells	CENTBLOOD
VBECS TAS 3 Cell-ABO-A1-Cells	A1Cells	PACKEDCELLS
VBECS TAS 3 Cell-ABO-Anti-A	AntiA	CENTBLOOD
VBECS TAS 3 Cell-ABO-Anti-A	AntiA	PACKEDCELLS
VBECS TAS 3 Cell-ABO-Anti-B	AntiB	CENTBLOOD
VBECS TAS 3 Cell-ABO-Anti-B	AntiB	PACKEDCELLS
VBECS TAS 3 Cell-ABO-B-Cells	BCells	CENTBLOOD
VBECS TAS 3 Cell-ABO-B-Cells	BCells	PACKEDCELLS
VBECS TAS 3 Cell-ABScr	ABSInterp	CENTBLOOD
VBECS TAS 3 Cell-ABScr	ABSInterp	PACKEDCELLS
VBECS TAS 3 Cell-ABScr-0.8-Surg1	SC1	CENTBLOOD
VBECS TAS 3 Cell-ABScr-0.8-Surg1	SC1	PACKEDCELLS
VBECS TAS 3 Cell-ABScr-0.8-Surg2	SC2	CENTBLOOD
VBECS TAS 3 Cell-ABScr-0.8-Surg2	SC2	PACKEDCELLS
VBECS TAS 3 Cell-ABScr-0.8-Surg3	SC3	CENTBLOOD
VBECS TAS 3 Cell-ABScr-0.8-Surg3	SC3	PACKEDCELLS
VBECS TAS 3 Cell-Rh	RhInterp	CENTBLOOD
VBECS TAS 3 Cell-Rh	RhInterp	PACKEDCELLS
VBECS TAS 3 Cell-Rh-Anti-D	AntiD	CENTBLOOD
VBECS TAS 3 Cell-Rh-Anti-D	AntiD	PACKEDCELLS
VBECS TAS 3 Cell-Rh-Ctrl	DControl	CENTBLOOD
VBECS TAS 3 Cell-Rh-Ctrl	DControl	PACKEDCELLS

Appendix C: Instrument Side Rules

Vision (Vision side interface)
ocdvisii / Ortho Vision / v8.00.0002

--- Test / In Validation - Incoming result - Before Message Queued Internally ---

Rule # 1

Desc - SetInstrumentID
If - {Always}
Then - {Set} {Instrument ID} = ""

Rule # 2

Desc - SetReceivingFacility
If - {Always}
Then - {Set} {Receiving Facility} = ""

Rule # 3

Desc - DonorID
If - (({Length of} {Specimen ID} = "16") {AND} ({Specimen ID} {Contains} "="))
Then - {Set} {Specimen ID} = {Extract Section of} {Specimen ID} {From} "2" {To} "14"

Rule # 4

Desc - DonorIDForProductID
If - (({Length of} {Product ID} {On Test} {Value List:crossmatch} = "16") {AND} ({Product ID} {On Test} {Value List:crossmatch} {Contains} "="))
Then - {Set} {Product ID} {On Test} {Value List:crossmatch} = {Extract Section of} {Product ID} {On Test} {Value List:crossmatch} {From} "2" {To} "14"

Value List

"Row Enabled","crossmatch"
"1","VBECS Crossmatch IAT-XM"
"1","VBECS Crossmatch IAT-Anti-IgG (Rabbit)-Donor Info"
"1","VBECS Crossmatch-XM"
"1","VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info"
"1","VBECS Crossmatch-Buffered Gel-Donor Info"
"0",

Rule # 5

Desc - Set Test Name for XM
If - {Test Code} {On Any Test} {Contains} "Crossmatch"
Then - {Set} {Test Code Sub ID} {On That Test} = "XM"

Rule # 6

Desc - Set Test Name for TAS
Child of Rule # 5 - Else
If - {Test Code} {On Any Test} {Contains} "TAS"
Then - {Set} {Test Code Sub ID} {On That Test} = "TAS"

Rule # 7

Desc - Set Test Name for DAT
Child of Rule # 6 - Else
If - {Test Code} {On Any Test} {Contains} "DAT"
Then - {Set} {Test Code Sub ID} {On That Test} = "DAT"

Rule # 8

Desc - Set Test Name for Patient ABORh
Child of Rule # 7 - Else
If - {Test Code} {On Any Test} {Contains} "ABDRev"
Then - {Set} {Test Code Sub ID} {On That Test} = "Patient ABORh"

This document is version controlled via approval signature. Do not assume that a printed, copied, or downloaded document is current. Consult the Electronic Document Control System to verify use of the current version.

Rule # 9

Desc - Set Test Name for ABS
Child of Rule # 8 - Else
If - {Test Code} {On Any Test} {Contains} "ABS"
Then - {Set} {Test Code Sub ID} {On That Test} = "ABS"

Rule # 10

Desc - Set Test Name for Unit ABORh
Child of Rule # 9 - Else
If - {Test Code} {On Any Test} {Contains} "Confirm"
Then - {Set} {Test Code Sub ID} {On That Test} = "Unit ABORh"

Rule # 11

Desc - Set Test Name for Antigen Typing
Child of Rule # 10 - Else
If - {Test Code} {On Any Test} {Contains} "Phenotyp"
Then - {Set} {Test Code Sub ID} {On That Test} = "Antigen Typing"

Rule # 12

Desc - OneTestResult
If - {Result} {On Any Test} = "10"
Then - {Set} {Result} {On That Test} = "1+"

Rule # 13

Desc - TwoTestResult
If - {Result} {On Any Test} = "20"
Then - {Set} {Result} {On That Test} = "2+"

Rule # 14

Desc - ThreeTestResult
If - {Result} {On Any Test} = "30"
Then - {Set} {Result} {On That Test} = "3+"

Rule # 15

Desc - FourTestResult
If - {Result} {On Any Test} = "40"
Then - {Set} {Result} {On That Test} = "4+"

Rule # 16

Desc - CMPXMResult
If - {Result} {On Any Test} = "CMP"
Then - {Set} {Result} {On That Test} = "Compatible"

Rule # 17

Desc - INCMXPXMResult
If - {Result} {On Any Test} = "INCM"
Then - {Set} {Result} {On That Test} = "Incompatible"

Rule # 18

Desc - POSResult
If - {Result} {On Any Test} = "POS"
Then - {Set} {Result} {On That Test} = "Pos"

Rule # 19

Desc - NEGRResult
If - {Result} {On Any Test} = "NEG"
Then - {Set} {Result} {On That Test} = "Neg"

Rule # 20

Desc - Set Result Status R
If - {Result Status} {On Any Test} = "R"
Then - {Set} {Result Status} {On That Test} = "F"

Rule # 21

Desc - Set Result Status Empty
If - {Result Status} {On Any Test} = ""
Then - {Set} {Result Status} {On That Test} = "F"

Rule # 22

Desc - Set Result DateTime for TAS-ABO 3 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS TAS 3 Cell-ABO"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 3 Cell-ABO"

Rule # 23

Desc - Set Result DateTime for TAS-Rh 3 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS TAS 3 Cell-Rh"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 3 Cell-Rh"

Rule # 24

Desc - Set Result DateTime for TAS-ABS 3 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS TAS 3 Cell-ABScr"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 3 Cell-ABScr"

Rule # 25

Desc - Set Result DateTime for TAS-ABO 2 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS TAS 2 Cell-ABO"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 2 Cell-ABO"

Rule # 26

Desc - Set Result DateTime for TAS-Rh 2 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS TAS 2 Cell-Rh"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 2 Cell-Rh"

Rule # 27

Desc - Set Result DateTime for TAS-ABS 2 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS TAS 2 Cell-ABScr"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 2 Cell-ABScr"

Rule # 28

Desc - Set Result DateTime for Patient ABORh - ABO
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS ABDRev-ABO"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS ABDRev-ABO"

Rule # 29

Desc - Set Result DateTime for Patient ABORh - Rh
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS ABDRev-Rh"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS ABDRev-Rh"

Rule # 30

Desc - Set Result DateTime for ABS 3 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS ABS 3 Cell-ABScr"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS ABS 3 Cell-ABScr"

Rule # 31

Desc - Set Result DateTime for ABS 2 Cell
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS ABS 2 Cell-ABScr"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS ABS 2 Cell-ABScr"

Rule # 32

Desc - Set Result DateTime for AGC
If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Phenotyping C-Anti-C (RH2)"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Phenotyping C-Anti-C (RH2)"

Rule # 33

Desc - Set Result DateTime for AGE

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Phenotyping E-Anti-E (RH3)"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Phenotyping E-Anti-E (RH3)"

Rule # 34

Desc - Set Result DateTime for AGc

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Phenotyping smc-Anti-c (RH4)"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Phenotyping smc-Anti-c (RH4)"

Rule # 35

Desc - Set Result DateTime for AGE

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Phenotyping sme-Anti-e (RH5)"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Phenotyping sme-Anti-e (RH5)"

Rule # 36

Desc - Set Result DateTime for Crossmatch IAT

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Crossmatch IAT"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Crossmatch IAT-XM"

Rule # 37

Desc - Set Result DateTime for Crossmatch

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Crossmatch" {AND} {Test Code} {On That Test} {NOT} {Contains} "IAT"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Crossmatch-XM"

Rule # 38

Desc - Set Result DateTime for DAT IgG

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS DAT IgG-IgG"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS DAT IgG-IgG"

Rule # 39

Desc - Set Result DateTime for DAT

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS DAT-Poly"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS DAT-Poly"

Rule # 40

Desc - Set Result DateTime for Unit ABO

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Confirm AB-ABO"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Confirm AB-ABO"

Rule # 41

Desc - Set Result DateTime for Unit ABORh - ABO

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Confirm ABD-ABO"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Confirm ABD-ABO"

Rule # 42

Desc - Set Result DateTime for Unit ABORh - Rh

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Confirm ABD-Rh"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Confirm ABD-Rh"

Rule # 43

Desc - Set Result DateTime for Unit ABORh - O Pos

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Confirm O Pos-Anti-AB"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Confirm O Pos-Anti-AB"

Rule # 44

Desc - Set Result DateTime for Unit ABORh - O Neg - ABO

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Confirm O Neg-Anti-AB"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Confirm O Neg-Anti-AB"

Rule # 45

Desc - Set Result DateTime for Unit ABORh - O Neg - Rh

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS Confirm O Neg-Rh"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS Confirm O Neg-Rh"

Rule # 46

Desc - Set Result DateTime for DAT+IgG IgG

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS DAT Poly+IgG-IgG"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS DAT Poly+IgG-IgG"

Rule # 47

Desc - Set Result DateTime for DAT+IgG Poly

If - {Result Date/Time} {On Any Test} = "" {AND} {Test Code} {On That Test} {Contains} "VBECS DAT Poly+IgG-Poly"
Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS DAT Poly+IgG-Poly"

Rule # 48

Desc - SetISInterp

If - {Test Code} {On Any Test} = "VBECS Crossmatch-Buffered Gel-Donor Info"

Then - {Add Test} "ISInterp" {AND} {Set} {Product ID} {On Test} "ISInterp" = {Product ID} {On That Test} {AND} {Set} {Test Code Sub ID} {On Test} "ISInterp" = {Test Code Sub ID} {On That Test}

Rule # 49

Desc - SetAHGInterp

If - {Test Code} {On Any Test} = "VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info"

Then - {Add Test} "AHGInterp" {AND} {Set} {Product ID} {On Test} "AHGInterp" = {Product ID} {On That Test} {AND} {Set} {Test Code Sub ID} {On Test} "AHGInterp" = {Test Code Sub ID} {On That Test}

Appendix D: VBECS (HL7) Side Rules

HIN_VBECS_1 (HL7 side interface)

diuhl7ml / Data Innovations LLC, Configurable HL7 / v8.00.0049

--- Test / In Validation - Outgoing result - Before Message Sent to this Connection ---

Rule # 1

Desc - SetISInterp - Comp

If - {Test Code} {On Any Test} = "VBECS Crossmatch-Buffered Gel-Donor Info" {AND} {Result} {On That Test} = "0"

Then - {Set} {Result} {On Test} "ISInterp" = "Compatible" {AND} {Set} {Result Status} {On Test} "ISInterp" = {Result Status} {On That Test} {AND} {Set} {Result Date/Time} {On Test} "ISInterp" = {Result Date/Time} {On That Test}

Rule # 2

Desc - SetISInterp - Incomp

If - {Test Code} {On Any Test} = "VBECS Crossmatch-Buffered Gel-Donor Info" {AND} {Result} {On That Test} {Contains} "+"

Then - {Set} {Result} {On Test} "ISInterp" = "Incompatible" {AND} {Set} {Result Status} {On Test} "ISInterp" = {Result Status} {On That Test} {AND} {Set} {Result Date/Time} {On Test} "ISInterp" = {Result Date/Time} {On That Test}

Rule # 3

Desc - SetAHGInterp - Comp

If - {Test Code} {On Any Test} = "VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info" {AND} {Result} {On That Test} = "0"

Then - {Set} {Result} {On Test} "AHGInterp" = "Compatible" {AND} {Set} {Result Status} {On Test} "AHGInterp" = {Result Status} {On That Test} {AND} {Set} {Result Date/Time} {On Test} "AHGInterp" = {Result Date/Time} {On That Test}

Rule # 4

Desc - SetAHGInterp - Incomp

If - {Test Code} {On Any Test} = "VBECS Crossmatch-Anti-IgG (Rabbit)-Donor Info" {AND} {Result} {On That Test} {Contains} "+"

Then - {Set} {Result} {On Test} "AHGInterp" = "Incompatible" {AND} {Set} {Result Status} {On Test} "AHGInterp" = {Result Status} {On That Test} {AND} {Set} {Result Date/Time} {On Test} "AHGInterp" = {Result Date/Time} {On That Test}

Rule # 5

Desc - Set POS Result for AB

If - {Result} {On Any Test} = "Pos" {AND} {Test Code} {On That Test} {Contains} "Anti-AB"

Then - {Set} {Result} {On That Test} = "AB"

Rule # 6

Desc - Set NEG Result for AB

If - {Result} {On Any Test} = "Neg" {AND} {Test Code} {On That Test} {Contains} "Anti-AB"

Then - {Set} {Result} {On That Test} = "O"

Rule # 7

Desc - Equalize Result DateTime for TAS 3 Cell

If - {Test Code} {On Any Test} {Contains} "VBECS TAS 3 Cell"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 3 Cell-ABScr"

Rule # 8

Desc - Equalize Result DateTime for TAS 2 Cell

If - {Test Code} {On Any Test} {Contains} "VBECS TAS 2 Cell"

Then - {Set} {Result Date/Time} {On That Test} = {Result Date/Time} {On Test} "VBECS TAS 2 Cell-ABScr"

Appendix E: Vision Testing Profiles

Table 3: Vision Testing Profiles

Profile Name	Card Type(s)	Test Name	VBECS Test
VBECS TAS 3 Cell	A/B/D Mono and Reverse Anti-IgG	A/B/D and Reverse Group 3 Cell Screen	TAS
VBECS TAS 2 Cell	A/B/D Mono and Reverse Anti-IgG	A/B/D and Reverse Group 2 Cell Screen	TAS
VBECS ABDRev	A/B/D Mono and Reverse	A/B/D and Reverse Group	Patient ABO/Rh
VBECS ABS 3 Cell	Anti-IgG	3 Cell Screen	Antibody Screen
VBECS ABS 2 Cell	Anti-IgG	2 Cell Screen	Antibody Screen
VBECS DAT	Anti-IgG, -C3d; Polyspecific (Rabbit)	DAT (IgG C3d)	DAT (Polyspecific)
VBECS DAT IgG	Anti-IgG	DAT (IgG)	DAT IgG
VBECS DAT Poly+IgG	Anti-IgG, -C3d; Polyspecific (Rabbit)	DAT Specificity	DAT (Polyspecific) and DAT IgG
VBECS Phenotyping C	Anti-C	Anti-C	Patient and Unit AGC
VBECS Phenotyping smc	Anti-c	Anti-c	Patient and Unit AGc
VBECS Phenotyping E	Anti-E	Anti-E	Patient and Unit AGE
VBECS Phenotyping sme	Anti-e	Anti-e	Patient and Unit AGe
VBECS Crossmatch	Anti-IgG Buffered Gel	Crossmatch-IAT Crossmatch-IS	Full Crossmatch (IS + IAT)
VBECS Crossmatch IAT	Anti-IgG	Crossmatch-IAT	IAT only Crossmatch
VBECS Confirm ABD	A/B/D Mono Grouping	A/B/D Grouping	Unit ABO/Rh Confirmation
VBECS Confirm AB	A/B Mono Grouping	A/B Group	Unit ABO Confirmation
VBECS Confirm O Pos	Anti-A,B	Anti-A,B	Unit ABO Confirmation for O Pos
VBECS Confirm O Neg	Anti-A,B	Anti-A,B and Anti-D	Unit ABO Confirmation for O Neg


Appendix F: CA SDM Support Ticket Template

Please use the following to complete your ticket:

- **Incident area:** NTL.APP.HealtheVet VistA.VBECS 2_0
- **Summary:** VBECS 2.2.0 Automated Instrument configuration assistance request
- **Description:**
 - Enter contact information for the main blood bank and the name of the responsible individual(s) for all communications.
 - Enter preferred date and time.
- **Property Questions:** Respond NO to the required questions.
- **Group:** automatically fills in correctly do not change it from. *NTL Alert Blood Bank & VBECS*
- **Add additional information as required.**

Appendix G: Validation Planning and Example Test Scenarios

The following is a flowchart to help assess any one change and plan accordingly.

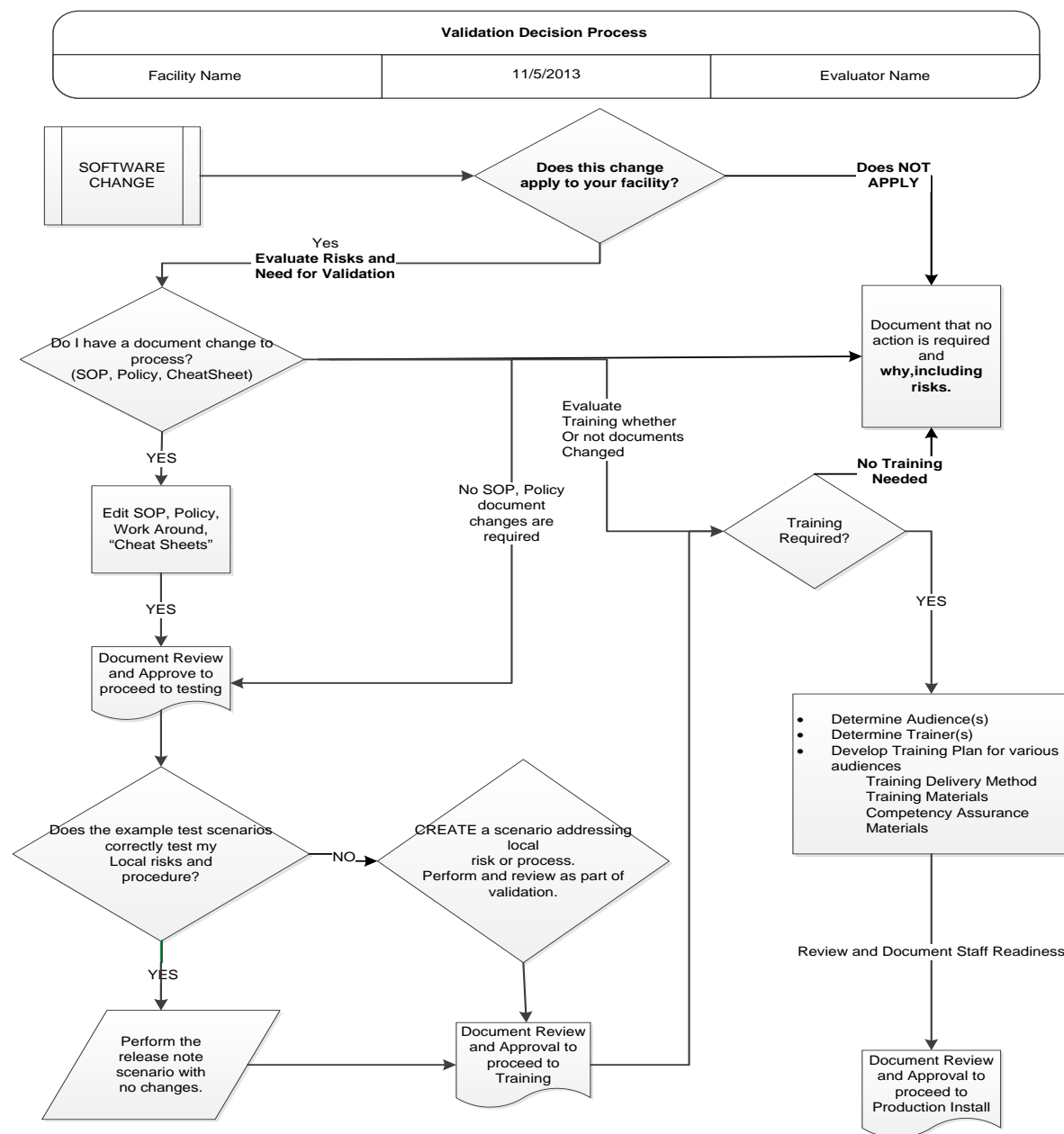



These are examples of test scenarios. Each site is responsible for evaluating changes for their intended use, assess risk, and for establishing additional validation test scenarios.

All test scenarios have a suggested user role, this may require adjustment to align with the patient or unit data selected to execute the scenario. Process any overrides as well, based on patient or unit selection.

**The Expected Outcome numbering corresponds to the Step number where the change verification appears.*

Test Group One: AI interface disabled

Test Objective: Demonstrate that the system will reject test results sent from automated instrument if the Automated Instrument interface is disabled.

Option: VBECS Administrator, Automated Instrument Interface

At least one instrument available in any division and configured for connection with VBECS via Instrument Manager.

Note: This may be executed at only one division if VBECS is used in a multidivisional configuration as the interface is enabled/disabled for all or none.

Data	<p>Before beginning, verify current configuration, activity, and status of the existing interfaces, VistALink, and CPRS in your test account. Make sure that the VBECS-OERR HL7 link in VistA is NOT shutdown.</p> <ol style="list-style-type: none"> 1. Create RBC and TAS orders in CPRS and accept them in VBECS.
User	<p>VBECS Administrator access is required to configure the AUTOMATED TESTING interface. No specific user role is required to process TAS and RBC orders in VBECS.</p>
Steps	<ol style="list-style-type: none"> 1. Log into VBECS Administrator: 2. Disable the AUTOMATED TESTING interface using Configure Interfaces option. 3. A message appears (You are about to disable Auto Instrument Interface. It will cause VBECS to stop sending and receiving messages via that interface. Continue?), Click Yes. 4. Close the window as all fields are disabled. 5. Perform TAS test on the instrument and send results of testing to VBECS
Expected Outcome	<ol style="list-style-type: none"> 5. Verify that after sending test results from an instrument fails it is possible to complete TAS test using manual testing grids and subsequently issue and transfuse blood unit to a patient.
Reports:	<ol style="list-style-type: none"> 6. Review the Audit Trail report for changes to the interface.

Test Group Two: Verify AI individual test(s)

Test Objective: Demonstrate that the system will allow the selected user role to perform normal workflow activities using your local blood bank testing instrument.

Note: Validate all tests that are performed at your site using an automated instrument.



Prior to initiating any testing of the automated instrument interface, see the VistA Blood Establishment Computer Software (VBECS) 2.2.0 Technical Manual-Security Guide for instructions for configuring your local instrument and Data Innovations Instrument Manager.



Configure your Test UID to clearly differentiate from a Production UID during your validation, for example, production is 1234567890, with 123 being your consistent prefix, change 123 or 12 to TS or TST or 999. Please change the Numeric Identifier for the Blood Bank Accession Area in your site Test Account.

Enable Automated Instrument Interface in VBECS Admin if it was disabled while performing testing in Test Group One.

Verify your VBECS processes with the Automated Instrument interface.

Your local test plan will demonstrate that the system will perform normal daily work per your local policies, procedures and local validation plan that may include:

Blood Units: Automated Instrument or via the short cut key 

- ABO-Rh confirmation
- Unit Antigen Typing or Repeat

Patients: Automated Instrument or via the short cut key 

- Perform a Type and Screen test (ABS and ABS with Auto-control only)
- Perform crossmatch tests: Serologic (*Selected in Blood Units: Select Units before testing.**)
- Patient Tests:
 - ABO/Rh or Repeat
 - Antibody Screen or Repeat
 - Direct Antiglobulin Test or Repeat
 - Antigen Typing or Repeat
- Overrides, (ABO/Rh Discrepancy, Crossmatch incompatible: Give with Medical Director Approval)
- Report data from these actions is available for retrieval.
 - Testing Worklist Report
 - Patient History Report
 - Unit History Report
 - Exception Report (ABO/Rh Discrepancy, Give with Medical Director Approval)

*Remember to click NO to proceeding to the serologic crossmatch when selecting units for automated instrument testing.

Test Group 2 Scenario 1: Verify AI TAS test

Note: the Type and Screen (TAS) is a combination of two tests and may be saved as individual tests without completing the TAS as a whole.

Data	VBECS: Select a patient. CPRS: Place a Type and Screen (TAS) order for the patient LAB: Accession the order VBECS: Accept the order. (Orders, Accept Orders)
User	No specific user role is required to process TAS order in VBECS.
Steps	<ol style="list-style-type: none"> 1. User checks the Patient testing list to make sure the order is accepted in VBECS and appears on the Pending Task List (PTL). (Patient, Patient Testing, Diagnostic Tests). 2. Close the PTL. 3. Process the specimen on the instrument using the recommended TAS profile for that instrument. Complete all work needed to transmit the information to VBECS. 4. Select Patients, Automated Instrument to review TAS results. 5. Select the specimen UID, scanning the UID is preferred. 6. Review TAS test results. 7. Accept only the ABS or ABO/Rh test. Close the window. 8. Open the PTL. 9. Try to select TAS on PTL and complete testing on it manually. 10. Open the Automated instrument window and accept the second part of TAS. 11. Close the automated instrument window. 12. Check Reports
Expected Outcome	<ol style="list-style-type: none"> 5. Verify that the specimen UID is selectable by scanning, entry or patient selection 6. Verify that the correct test results appear on the Automated Instrument review list. 8. Verify that the TAS appears on the PTL with a status of "Instrument Results Pending Review". 9. Verify that system prevents user from completing TAS since they are still pending results from an instrument for it. 12. Verify that the results and comments appear as expected on the reports: <ul style="list-style-type: none"> • Testing Worklist Report • Patient History Report • Exception Report

Test Group 2 Scenario 2: Verify AI serologic crossmatch test	
Data	<p>VBECS: Select a patient. CPRS: Place a Type and Screen (TAS) order and Red Blood Cell order for the patient LAB: Accession the orders</p> <p>VBECS: Accept the order. (Orders, Accept Orders) Process the TAS to completion.</p> <p>Select a blood unit for a selected patient*:</p> <ul style="list-style-type: none"> • Previously entered into the division's inventory (Blood Units, Incoming Shipment) • ABO compatible • May be available or selected for another patient (available, selected, crossmatched) status. • May or may not trigger selection overrides <p>*Remember to click NO to proceeding to the serologic crossmatch when selecting units for automated instrument testing.</p>
User	No specific user role is required to process crossmatch test in VBECS.
Steps	<ol style="list-style-type: none"> 1. Process the component unit's specimen on the instrument using the recommended profile for that instrument. Complete all work needed to transmit the information to VBECS. 2. Select Patients, Automated Instrument to review crossmatch results. 3. Select the specimen UID, scanning the UID is preferred. 4. Review results, select compatibility and accept crossmatch result. 5. Print or do not print tag as desired. 6. Accept the test and close the automated instrument window. 7. Check Reports.
Expected Outcome	<ol style="list-style-type: none"> 4. Verify that crossmatch results sent from an instrument show correctly on the screen. 7. Verify that the results and comments appear as expected on the reports: <ul style="list-style-type: none"> • Testing Worklist Report • Patient History Report (interpretations only). • Exception Report

Test Group 2 Scenario 3: Verify AI patient diagnostic tests (ABO/Rh, Antibody Screen, Direct Antiglobulin Test, Patient Antigen Typing, and the reflex test)	
Data	VBECS: Select a patient. CPRS: Place a diagnostic test order for the patient LAB: Accession the order VBECS: Accept the order. (Orders, Accept Orders)
User	No specific user role is required to process diagnostic tests in VBECS.
Steps	<ol style="list-style-type: none"> 1. User checks the Patient testing list to make sure the order is accepted in VBECS and appears on the Pending Task List (PTL). (Patients, Patient Testing, Diagnostic Tests). 2. Close the PTL. 3. Process the specimen on the instrument using the recommended testing profile for that instrument. Complete all work needed to transmit the information to VBECS. 4. Select Patients, Automated Instrument to review test results. 5. Select the specimen UID, scanning the UID is preferred. 6. Review and accept test results 7. Close the automated instrument window. 8. Check reports.
Expected Outcome	<ol style="list-style-type: none"> 5. Verify that the specimen UID is selectable by scanning, entry or patient selection. 6. Verify that patient test results sent from an instrument show correctly on the screen. 8. Verify that the results and comments appear as expected on the reports: <ul style="list-style-type: none"> • Testing Worklist Report • Patient History Report • Exception Report

Test Group 2 Scenario 4: Verify AI blood unit tests	
Data	<p>VBECS: Select a blood unit previously entered into the division's inventory (Blood Units, Incoming Shipment).</p> <p>For ABO/Rh Confirmation testing on the instrument, unit should be in a Limited status.</p> <p>For Unit Antigen Typing, the unit may or may not have been confirmed.</p>
User	No specific user role is required to process blood unit tests in VBECS.
Steps	<ol style="list-style-type: none"> 1. Process the component unit's specimen on the instrument using the recommended profile for that instrument. Complete all work needed to transmit the information to VBECS. 2. Select Blood Units, Automated Instrument to review test results. 3. Select the blood component unit's Donor Identification Number (DIN), scanning the DIN is preferred. 4. Select the product code (if there are multiple blood units with the same product code) 5. Review the transmitted blood unit test. 6. Accept the test and close the automated instrument window. 7. Check Reports
Expected Outcome	<ol style="list-style-type: none"> 5. Verify that blood unit test results show correctly on the screen. 7. Verify that the results and comments appear as expected on the reports: <ul style="list-style-type: none"> • Testing Worklist Report • Unit History Report (interpretations only). • Exception Report