

**Laboratory**  
**Version 5.2**  
**Search/Extract Technical and User Guide**



**September 2015**

**Department of Veterans Affairs (VA)**  
**Office of Information and Technology (OI&T)**  
**Product Development (PD)**

## Revision History

Date	Version	Description	Author
9/2015	LR*5.2*442	<p>Updates for LR*5.2*442, ICD-10 PTF Modifications:</p> <p>Updated title page and footers, updated Preface (p.iv), added new DBIA #6130 (p.15), reformatted and updated Routine List (p.16), updated patches required (p.17), updated Installation Instructions (p.33), updated Post Installation Instructions (p.38).</p>	<p>T.B., PM</p> <p>K.R., Technical Writer</p>
07/2014	LR*5.2*421	<p>Updated for ICD-10 Remediation:</p> <p>Added Revision History.</p> <p>Updated TOC.</p> <p>p. iv Updated information in Preface.</p> <p>pp. <a href="#">14</a>, <a href="#">20</a>, <a href="#">53</a>, <a href="#">58</a>, <a href="#">87</a>, <a href="#">109</a>, <a href="#">125</a>, General Changes from ICD-9 to ICD</p> <p>p. <a href="#">16</a> Updated Routine List and updated Software Requirements</p> <p>p. <a href="#">27</a> Updated File (#69.5) file fields for Help.</p> <p>p. <a href="#">32</a> Added fields to Data Dictionary.</p> <p>pp. <a href="#">17</a>, <a href="#">33</a>, <a href="#">38</a> Added notes that refer users to the LR*5.2*421 Release Notes for Required Patches, Installation and Post Installation instructions.</p> <p>p. <a href="#">45</a> Updated DG1 Segment of HL7 message to utilize ~ rather than ^.</p> <p>pp. <a href="#">46</a>, <a href="#">47</a>, <a href="#">116</a>, <a href="#">135</a>, <a href="#">136</a>, <a href="#">139</a> Removed Personally Identifiable Information (PII).</p> <p>pp. <a href="#">62</a>, <a href="#">66</a>, <a href="#">69-70</a>, <a href="#">72-73</a>, <a href="#">75-76</a>, <a href="#">79-80</a>, <a href="#">82-83</a>, <a href="#">85-85</a>, <a href="#">88-89</a>, <a href="#">91-92</a>, <a href="#">94-95</a>, <a href="#">97-98</a>, <a href="#">100-101</a>, <a href="#">104</a>, <a href="#">129-130</a>, <a href="#">132-133</a>, <a href="#">145</a> Updated Primary Menu and examples for addition of Code System prompt.</p>	<p>K.K., VA PM;</p> <p>D.G., ICD-10 PM;</p> <p>E. P., Technical Writer</p>

Date	Version	Description	Author
	Multiple pages	Changed ICDM-9 to ICD-10.	
	Multiple pages	Austin Automation Center (AAC) changed to Austin Information Technology Center (AITC)	
05/2008	LR*5.2*175	Initial release.	VA OI&T

## Preface

The Veterans Health Information Systems and Architecture (*VISTA*) Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide provides assistance for installing, implementing, and maintaining the *VISTA* Laboratory Search/Extract Patch LR\*5.2\*175 software application.

The International Classification of Diseases, Tenth Revision (ICD-10) Remediation patch LR\* 5.2\*421 makes the following changes to the Emerging Pathogens Initiative (EPI) application:

- The following fields and screens have been updated to refer to “ICD” rather than “ICD9”:
  - Laboratory Search/Extract Parameters Input screens
  - Enter/Edit Local Pathogens screens
  - Detailed Verification Report
  - Help text
- Within the Enter/Edit Local Pathogens and Laboratory Search/Extract Parameters Input screens, users are prompted to specify a code set on which to search prior to entering an ICD code. Based on this input, the system will only allow ICD-9 entry or ICD-10 entry.
- The Pathogen Inquiry option has been modified to list both ICD-9 and ICD-10 codes.
- The Generate Local Report/Spreadsheet option has been modified to include both the Diagnosis Code Set Designation and the Diagnosis Code.
- Health Level 7 (HL7) Reports that are sent to Austin Information Technology Center (AITC) have been modified to include ICD-10 Codes and Descriptions, which are included in the DG1 HL7 Segments.

**The ICD-10 PTF Modifications patch LR\*5.2\*442 made changes to accommodate the expanded number of ICD-10 codes that can now be entered in a patient record.**

### Recommended Users:

Department of Veterans Affairs Medical Center (DVAMC) Information Resource Management (IRM) staff

Laboratory Information Manager (LIM)

Representative from the Microbiology section for the Emerging Pathogens Initiative (i.e., director, supervisor, or technologist)

Total Quality Improvement/Quality Improvement/Quality Assurance (TQI/QI/QA)  
staff (or person at the facility with similar function)

## **VISTA Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide Sections:**

**Pre-Installation Information:** This section contains pre-installation information.

**Installation Instructions:** This section of the guide provides an example of the VISTA Laboratory Search/Extract Patch LR\*5.2\*175 installation process.

**Post Installation Instructions:** This section provides all the necessary information required for the IRM and LIM staff to implement the software application.

**Appendix A - EPI User Guide:** The EPI User Guide provides the necessary information for implementing and maintaining the EPI search/extract criteria.

**Appendix B - NCH User Guide:** The NCH User Guide provides the necessary information for implementing and maintaining the NCH search/extract criteria.

**Appendix C - Editing Files, Input Screens, Linking Data, and Request Form:** This section provides instructions for editing files, input screens, linking data, and a Workload and Suffixes Codes Request Form.

# Table of Contents

VISTA LABORATORY SEARCH/EXTRACT PATCH LR*5.2*175 TECHNICAL AND USER GUIDE SECTIONS: -----	VI
<b>OVERVIEW -----</b>	<b>1</b>
MANDATE -----	1
FUNCTIONALITY -----	2
<i>Mail Messages</i> -----	2
EPI Mail Messages:-----	3
EPI HL7 formatted Mail Messages-----	3
EPI Verification Report Mail Messages -----	3
EPI Confirmation Mail Messages-----	3
EPI Processing Report Mail Message-----	3
NCH Mail Messages:-----	4
NCH HL7 Formatted Mail Messages -----	4
NCH Verification Report mail messages -----	4
HL7 Message Status List-----	4
NCH HL7 formatted Acknowledgment mail message-----	4
NCH Acknowledgment mail message -----	5
NCH VA Alert mail message -----	5
ENHANCEMENTS -----	6
<i>Cytology Search</i> -----	6
<i>New Option</i> -----	6
<i>Input Screens</i> -----	6
<i>Mail Groups</i> -----	6
<i>LRNCH Protocol</i> -----	7
<i>New mail messages</i> -----	7
<i>Files</i> -----	7
LAB SEARCH/EXTRACT PROTOCOL file (#69.4) New Fields -----	7
LAB SEARCH/EXTRACT file (#69.5) New Fields -----	7
MODIFICATIONS -----	9
<i>Files Renamed</i> -----	9
<i>Menu Renamed</i> -----	9
<i>Options Renamed</i> -----	9
<i>Option Replaced</i> -----	9
ORIENTATION -----	11
<i>Screen Displays</i> -----	11
<i>Computer Dialogue</i> -----	11
<i>User Response</i> -----	11
<i>Return Symbol</i> -----	11
<i>Tab Symbol</i> -----	11
<i>Technical and User Guide Distributions</i> -----	12
<i>Hard Copy</i> -----	12
<i>Electronic Distributions</i> -----	12
<i>ANONYMOUS.SOFTWARE Accounts</i> -----	12
<i>References</i> -----	12
<b>PRE-INSTALLATION INFORMATION-----</b>	<b>13</b>
HARDWARE AND OPERATING SYSTEM REQUIREMENTS-----	13
<i>Digital Equipment Corporation (DEC) Alpha Series</i> -----	13
<i>Personal Computer (PC) System</i> -----	13
SYSTEM PERFORMANCE CAPACITY-----	13
INSTALLATION TIME -----	13
USERS ON THE SYSTEM -----	14
STAFFING REQUIREMENTS-----	14
<i>IRM Staff</i> -----	14

## Table of Contents

<i>Laboratory Staff</i> -----	14
PERIODIC REVIEWS -----	14
BACKUP ROUTINES -----	15
TEST SITES -----	15
KERNEL INSTALLATION AND DISTRIBUTION SYSTEM (KIDS)-----	15
DATABASE INTEGRATION AGREEMENTS (DBIA)-----	15
NAMESPACE -----	16
ROUTINE LIST -----	16
VISTA SOFTWARE REQUIREMENTS -----	17
PATCHES REQUIRED -----	17
HEALTH LEVEL SEVEN (HL7)-----	17
DOMAIN-----	17
PROTOCOLS-----	18
MAIL GROUPS-----	18
LABORATORY SEARCH/EXTRACT MENU AND OPTIONS -----	19
LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN -----	20
PROMPTS DEFINITIONS-----	20
VISTA LAB SEARCH/EXTRACT FILES DATA DICTIONARIES-----	22
<i>LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4)</i> -----	22
<i>LAB SEARCH/EXTRACT file (#69.5)</i> -----	25
<b>INSTALLATION INSTRUCTIONS -----</b>	<b>33</b>
<b>POST INSTALLATION INSTRUCTIONS-----</b>	<b>38</b>
<i>DSM/Alpha and Open M Sites</i> -----	38
TROUBLE SHOOTING -----	42
<b>HEALTH LEVEL SEVEN (HL7) PROTOCOL -----</b>	<b>44</b>
APPENDIX - A -----	52
<b>EPI USER GUIDE -----</b>	<b>54</b>
PURPOSE -----	54
EPI OBJECTIVE:-----	54
EPI DATA TRANSMISSION -----	55
<i>EPI HL7 formatted Mail Messages</i> -----	55
<i>EPI Verification Report Mail Messages</i> -----	55
<i>EPI Confirmation Mail Messages</i> -----	55
<i>EPI Processing Report Mail Message</i> -----	55
AUSTIN AUTOMATION CENTER DATABASE PROCESSING -----	56
<i>Numerator file:</i> -----	56
<i>Denominator file:</i> -----	57
LAB SEARCH/EXTRACT PROTOCOL EDIT [LREPI PROTOCOL EDIT] OPTION -----	58
EPI DESCRIPTIONS AND INPUT EXAMPLES -----	59
LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN -----	60
<i>Candida (Reference #8)</i> -----	61
<i>Clostridium difficile (Reference #4)</i> -----	65
<i>Creutzfeldt-Jakob Disease (CJD) (Reference #13)</i> -----	69
<i>Cryptosporidium (Reference #9)</i> -----	72
<i>Dengue (Reference #12)</i> -----	75
<i>E. coli O157:H7 (Reference #10)</i> -----	78
<i>Hepatitis C Antibody Positive (Reference #2)</i> -----	81
<i>Legionella (Reference #7)</i> -----	84
<i>Leishmaniasis (Reference #14)</i> -----	88
<i>Malaria (Reference #11)</i> -----	91
<i>Penicillin- Resistant Pneumococcus (Reference #3)</i> -----	94
<i>Streptococcus-Group A (Reference #6)</i> -----	97



## Table of Contents

<i>Tuberculosis (Reference #5)</i> -----	100
<i>Vancomycin-Resistant Enterococcus (VRE) (Reference #1)</i> -----	103
CONCLUSION -----	107
EPI HELPFUL HINTS: -----	108
<i>Clostridium difficile</i> -----	108
<i>Validating EPI Data Capture</i> -----	110
<i>Emerging Pathogens Verification Report Message</i> -----	111
<i>Protocols</i> -----	113
<i>Domain</i> -----	113
<i>EPI Mail Groups</i> -----	113
<i>Office of the Director (00)</i> -----	113
<i>EPI Mail Groups</i> -----	114
<i>Adding Mail Groups</i> -----	115
<i>Starting the Lower Level Protocol of the HL7 V. 1.6 Background Job</i> -----	115
<i>EPI HL7 Format Mail Message</i> -----	116
<i>EPI Confirmation Message</i> -----	116
<i>EPI Processing Report</i> -----	117
<i>EPI Table of Reject and Warning Codes</i> -----	118
<b>NCH USER GUIDE</b> -----	<b>124</b>
OVERVIEW-----	124
<i>Mandate</i> -----	125
<i>NCH Database Access</i> -----	125
<i>Impact</i> -----	125
<i>National Roll-Up</i> -----	125
<i>NCH Search and Extract Criteria</i> -----	126
<i>Recommended Users</i> -----	126
<i>Periodic Reviews</i> -----	126
<i>NCH Data Transmission</i> -----	127
NCH HL7 formatted mail messages-----	127
NCH Verification Report mail message -----	127
NCH Acknowledgment mail message -----	127
NCH VA View Alert mail message -----	127
<i>Cholesterol Screening</i> -----	128
Cholesterol Screening for Hyperlipidemia -----	128
<i>Papanicolaou (Pap) Screening</i> -----	132
Pap Smear for Cervical Cancer -----	132
NCH MAIL MESSAGES-----	135
NCH HL7 Formatted Mail Message -----	135
<i>NCH Verification Report mail message</i> -----	136
<i>NCH HL7 Mail Message Status List</i> -----	137
<i>NCH HL7 Formatted (Acknowledgment) mail message</i> -----	138
<i>NCH Acknowledgment mail message</i> -----	139
NCH VA ALERT-----	140
<b>EDITING FILES/SCREENS, LINKING DATA, REQUEST FORM</b> -----	<b>144</b>
EDITING TOPOGRAPHY FILE (#61) -----	144
HOW TO LINK ANTIMICROBIAL ENTRIES TO WORKLOAD CODES ENTRIES -----	145
<i>Using the Antimicrobial Link Update [LREPILK] options</i> -----	145
How to link entries using the AUTO option-----	145
How to add and delete entries to a file using the MANUAL option -----	145
How to add entries using the SEMI-AUTO option. -----	146
<i>How to Delete an Entry from the Laboratory Search/Extract Parameters Input Screen</i> -----	147
<i>How to add an entry using the Laboratory Search/Extract Parameters Input Screen</i> -----	148
<i>Additional Workload and Suffixes Codes Request Form</i> -----	149

## Overview

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 is an enhancement to the **VISTA** Laboratory Emerging Pathogens Initiative Patch LR\*5.2\*132 software application. The Laboratory Search/Extract enhancement software supports the following two national initiatives:

1. **Emerging Pathogens Initiative (EPI):** The Department of Veterans Affairs (DVA) Headquarters Infectious Disease Program Office assists with identifying new antibiotic-resistant, otherwise problematic pathogens within the DVA. Using this objective information, plans may be formulated on a national level for intervention strategies and resource needs. Results from the aggregate data may also be shared with the appropriate public health authorities for planning on the national level (i.e., private health care sectors).
2. **National Center for Health Promotion Cholesterol and Pap Screening (NCH):** Pursuant to the Congressional mandates stipulated in its enabling legislation (PL. 102-585, U.S.C. 17, 1704), the Department of Veterans Affairs (DVA) National Center for Health Promotion (NCHP) is tasked with monitoring and improving the prevalence of health promotion screening activities provided to veterans VA-wide. The purpose of the NCH Cholesterol and Pap Screening database is to monitor cholesterol and Pap screening activities at a national level, with the ultimate goal for improving detection and treatment of hyperlipidemia and cervical cancer. The NCH Cholesterol and Pap Screening database will also provide a valuable resource for clinical and health services researchers on screening activities and health outcomes, particularly for high-risk and special emphasis group such as veterans with hyperlipidemia, older veterans, and female veterans.

## Mandate

The installation of the **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software should be completed in accordance with PL. 102-585, U.S.C. 17, 1704. It is recommended that all VAMCs installed and implemented the software as soon as possible to help facilitate these important initiatives.

## Functionality

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 enhancement software application tracks the EPI and NCH ongoing initiatives within the VAMCs without the necessity of any additional data entry. The software will automatically search and extract data from several **VISTA** Software databases (e.g., Laboratory, PIMS, and Social Work). The LAB DATA file (#63) (i.e., contains verified lab data results), PTF file (#45), and the PATIENT file (#2) data are used for the defined search/extract criteria.

The **VISTA** Laboratory Search/Extract software application automatically generates the data as defined by the EPI and NCH criteria's. The EPI search/extract data is processed on the 15<sup>th</sup> of each month. The NCH search/extract data is processed daily.

## Mail Messages

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 enhancement software application automatically processes the EPI and NCH data producing HL7 formatted mail messages. A Verification Report mail message is then produced for each HL7 formatted mail message. The HL7 formatted mail messages are used for sending the EPI and NCH data to the Austin Information Technology Center (AITC) database.

The Verification Report mail message is a copy of the HL7 formatted mail messages in a human readable format. The Verification Report allows the users to review the EPI and NCH data and make corrections (e.g., complete social security numbers, valid Date of Births, and Period of Services, etc.) as deemed necessary.

The new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option enables the user to generate a Verification Report mail message (i.e., human readable format) as often as needed for reviewing the EPI and NCH data. This option will automatically send the NCH data to the AITC database each time the option is manually run.

## EPI Mail Messages:

**NOTE:** See Appendix A for examples of the EPI mail messages.

### EPI HL7 formatted Mail Messages

The **VISTA** Laboratory Search/Extract software automatically processes the EPI data on the 15<sup>th</sup> of each month producing HL7 formatted mail messages. A Verification Report mail message is then produced (i.e., in a human readable format) for each HL7 mail message. The EPI HL7 formatted mail messages **must** be **manually** forwarded on the 15<sup>th</sup> of each month to the AITC database via the Q-EPI.MED.VA.GOV domain. Use the VA MailMan software to manually forward the EPI data to XXX@ Q-EPI.MED.VA.GOV domain.

### EPI Verification Report Mail Messages

The EPI Verification Report mail messages are sent to the EPI-REPORT mail group on the 15<sup>th</sup> of each month. The members of the EPI-REPORT mail group should review the EPI Verification Report and make data corrections (e.g., social security number, date of birth, period of service, etc.) as deemed necessary. The software automatically processes the EPI data corrections on the 15<sup>th</sup> of each month. **Note:** The EPI Verification Report may be generated as often as needed using the new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option.

### EPI Confirmation Mail Messages

EPI Confirmation mail messages are sent to the sending facility EPI-REPORT mail group **after** the EPI HL7 mail messages has been received by the AITC database.

### EPI Processing Report Mail Message

EPI Processing Report mail messages are sent to the sending facility EPI-REPORT mail group at the end of the AITC processing cycle (i.e., the 25<sup>th</sup> of each month). The EPI Processing Report mail messages confirms that the EPI data has been processed and lists any errors and/or warning codes requiring corrections.

## **NCH Mail Messages:**

**NOTE:** See Appendix B for examples of the NCH mail messages.

### NCH HL7 Formatted Mail Messages

The **VISTA** Laboratory Search/Extract software automatically processes the NCH data daily producing HL7 formatted mail messages. A Verification Report mail message is then produced (i.e., in a human readable form) for each HL7 formatted mail message. The software automatically transmits the NCH HL7 formatted mail messages daily to the AITC database.

### NCH Verification Report mail messages

The NCH Verification Report mail messages are generated daily and sent to the LR NCH-Report mail group. The members of the LR NCH-Report mail group may use the Verification Report to review the NCH data and make corrections (e.g., social security number, date of birth, period of service, etc.) as deemed necessary. The new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option enables the users to generate the NCH Verification Report and to transmit the NCH data corrections to the AITC database as often as needed.

### HL7 Message Status List

The HL7 Message Status List contains the status of each transmission sent by the site to the AITC database.

### NCH HL7 formatted Acknowledgment mail message

The NCH HL7 formatted Acknowledgment mail message is sent by the AITC to the sending site. This message is translated in a human readable format and sent to the LR NCH Report mail group.

### NCH Acknowledgment mail message

The NCH Acknowledgment mail message displays the status of the NCH HL7 formatted mail message transmission to the AITC. The NCH Acknowledgment mail message is sent to the LR NCH-Report mail group.

### NCH VA Alert mail message

A NCH VA Alert mail message is sent to the LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4), Report Mail Group field (#1), after the NCH data has been processed by the AITC.

## Enhancements

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software consists of the following enhancements:

### Cytology Search

The Cytology Search/Extract searches the LAB DATA file (#63) based on the criteria defined by the NCH Cholesterol and NCH Pap search criteria in the LAB SEARCH/EXTRACT file (#69.5).

### New Option

Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option: This option is used to run the EPI or NCH extracts manually.

Lab Search/Extract Protocol Edit [LREPI PROTOCOL EDIT] option: Use this option to edit the LAB SEARCH/EXTRACT PROTOCOL file (#69.4).

### Input Screens

The Emerging Pathogens Parameters Input Screen is **renamed** Laboratory Search/Extract Parameters Input Screen. The Laboratory Search/Extract Parameters Input Screens were enhanced by adding a fifth screen to accommodate the following four new prompts (i.e., Lag Days, Sex, Before Date of Birth, and After Date of Birth).

### Mail Groups

Patch LR\*5.2\*175 automatically installs the LR NCH and LR NCH-Report mail groups during the installation process. The LR NCH and LR NCH-Report mail groups Coordinator(s) may be entered during the installation process. Other mail groups members **must** be added to the mail groups after the installation process.

## LRNCH Protocol

The LRNCH protocol is added to the LAB SEARCH/EXTRACT file (69.5), Protocol field (#12) during patch installation process.

### New mail messages

- NCH HL7 formatted mail messages
- NCH Verification Report mail messages
- NCH HL7 formatted Status List
- NCH HL7 formatted (Acknowledgment) mail messages
- NCH Acknowledgment messages
- CH VA Alert messages

N

## Files

### **LAB SEARCH/EXTRACT PROTOCOL file (#69.4) New Fields**

**Send Alert field (#3):** Enter 1 for YES to receive an alert when a search is run automatically.

**Receive Alert field (#5):** Select the user (s) or group responsible for receiving alerts when the search has been executed automatically.

**Title field (#4):** This field will over-ride the default title for the Verification Report. If this field is left, blank the default title Verification Report will be used.

### **LAB SEARCH/EXTRACT file (#69.5) New Fields**

**Lag Days field (#5):** This is the number of days that the search should wait after the test was ordered to ensure that the results have been entered in the system. Enter the number of lag days desired. The amounts allowed will vary depending on the type.

**Sex field (#16): FOR FUTURE USE ONLY**

**Before Date of Birth field (#17):** Patients born after the date entered will not be included in the report. Enter a date to screen out patients born after the date entered.



**After Date of Birth field (#18):** Persons born before the entered will not be included in the report.

**SNOMED Codes field (#8):** This is a list of SNOMED codes to be included in the search.

**PTF field (#14):** This field is used to follow the inpatient PTF information.

## Modifications

Listed below are the modifications made to the **VISTA** Laboratory Emerging Pathogens Initiative (EPI) Patch LR\*5.2\*132 software application files, menu, and options. These modifications were made to accomplish the enhancements for the **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software application.

### Files Renamed

The EMERGING PATHOGENS PROTOCOL file (#69.4) is **renamed** the LAB SEARCH/EXTRACT PROTOCOL file (#69.4).

The EMERGING PATHOGENS file (#69.5) is **renamed** the LAB SEARCH/EXTRACT file (#69.5).

### Menu Renamed

The Emerging Pathogens Primary Menu [LREPI EMERGING PATHOGENS MENU] is **renamed** the Lab Search/Extract Primary Menu [LREPI SEARCH EXTRACT MENU].

### Options Renamed


The Emerging Pathogens Parameter Setup [LREPI PARAMETER UPDATE] option is **renamed** the Lab Search/Extract Parameter Setup [LREPI PARAMETER SETUP] option.

The Emerging Pathogens Nightly Task [LREPI NIGHTLY TASK] is **renamed** the Lab Search/Extract Nightly Task [LREPI NIGHTLY TASK].

### Option Replaced

The Emerging Pathogens Manual Run [LREPI (EPI) MANUAL RUN] option is **replaced** with the new Lab Search/Extract Manual Run (Enhanced) [LREPI ENCHANCE MANUAL RUN] option.

## Laboratory Search/Extract Process Flowchart

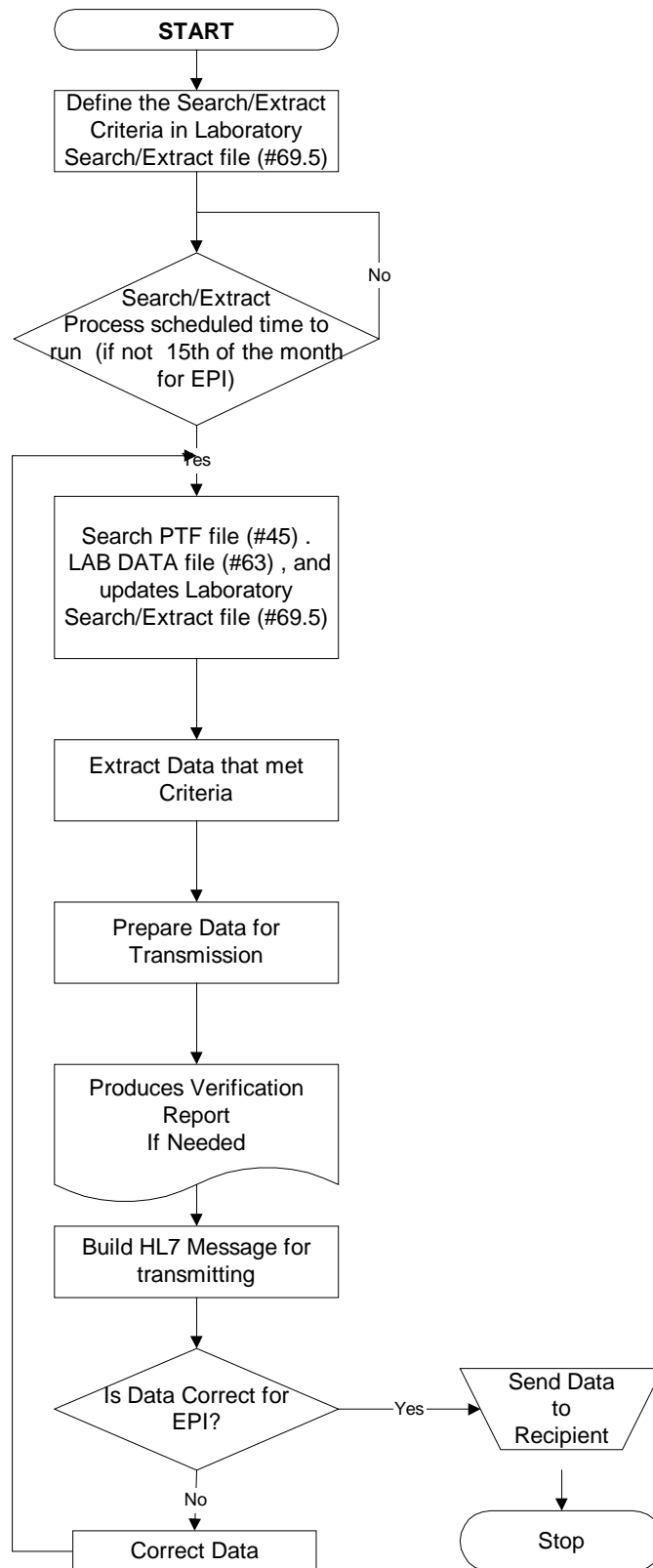

## Orientation

This section addresses the Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide computer screen dialogue, symbols, reference manuals, and the electronic and hard copy distributions.

### Screen Displays

The Lab Search/Extract Primary [LREPI SEARCH EXTRACT MENU] Menu options are using VA FileMan/ScreenMan forms for editing and displaying data. For detailed instructions using ScreenMan forms please refer to the VA FileMan V. 21.0 User Manual, Section 6 - ScreenMan.

### Computer Dialogue

The computer dialogue appears in Courier font, no larger than 10 points.

**Example:** Courier font 10 points

### User Response

User entry response appears in boldface type Courier font, no larger than 10 points.

**Example:** Boldface type

### Return Symbol

User response to computer dialogue is followed by the <RET> symbol which appears in Courier font, no larger than 10 points, and bolded.

**Example:** <RET>

### Tab Symbol

User response to computer dialogue is followed by the <Tab> symbol that appears in Courier font, no larger than 10 points, and bolded.

**Example:** <Tab>

## Technical and User Guide Distributions

The following information states how to obtain the Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide electronic and hard copy distributions.

### Hard Copy

The Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide hard copies are distributed by National Center for Documentation (NCD) to all Veteran Affairs Medical Centers (VAMCs).

### Electronic Distributions

The Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide is available on the Intranet in Hyper Text Markup Language (HTML) and Portable Document Format (PDF) at the following Intranet address:  
[http://152.127.1.95/softserv/clin\\_nar.row/lab/](http://152.127.1.95/softserv/clin_nar.row/lab/)

## ANONYMOUS.SOFTWARE Accounts

The Laboratory Search/Extract Patch LR\*5.2\*175 Technical and User Guide is also available on the ANONYMOUS SOFTWARE accounts at the Albany, Hines, and Salt Lake City Chief Information Officer Field Offices (CIOFOs) in PDF (i.e., P175TUG.PDF) at the following FTP addresses:

<b>CIOFOs</b>	<b>FTP Address</b>
Albany	152.127.1.5 - anonymous.software
Hines	152.129.1.110 - anonymous.software
Salt Lake City	152.131.2.1 - anonymous.software

## References

Kernel V. 8.0 Systems Manual  
HL7 V. 1.6 Manuals  
PIMS V. 5.3 Manuals  
VA FileMan V. 21.0 User Manual, Section 6 - ScreenMan  
MailMan V. 7.1 Manuals

## Pre-Installation Information

### Hardware and Operating System Requirements

**VISTA** software operates on two hardware platforms. The hardware platforms are listed in the mini-computer category, which provides multi-tasking and multi-user capabilities. The hardware platforms systems used are:

#### Digital Equipment Corporation (DEC) Alpha Series

Digital Equipment Corporation (DEC) Alpha series is using the DEC Open Virtual Memory System (VMS), Version 6.1 or greater, operating system. This platform uses the DEC System Mumps (DSM), Version 6.3 or greater, of American National Standards Institutes (ANSI) of Massachusetts General Hospital Utility Multi-Programming System (MUMPS) also known as 'M' language. MUMPS is a Federal Information Processing Standard (FIPS) language.

#### Personal Computer (PC) System

Personal Computer (PC) System with 486 or Pentium computer processor chip is using the Microsoft Disk Operating System (MS-DOS). The platform uses Open-M, of the American National Standards Institutes (ANSI) of Massachusetts General Hospital Utility Multi-Programming System (MUMPS) also known as 'M' language. MUMPS is a Federal Information Processing Standard (FIPS) language.

### System Performance Capacity

There are no changes in the performance of the system once Patch LR\*5.2\*175 is installed.

### Installation Time

Installation time is less than 2 minutes during off peak hours and less than 5 minutes during peak hours.

## Users on the System

Users may remain on system and no options need to be placed out of service.

## Staffing Requirements

The following staff is required to successfully install and implement Patch LR\*5.2\*175 software:

### IRM Staff

An IRM staff is required for installing Patch LR\*5.2\*175, setting up the domains, mail groups, and menu assignments.

### Laboratory Staff

It is **highly recommended** that the following person (s) jointly participate in reviewing the parameter descriptions:

Laboratory Information Manager (LIM)

Representative from the Microbiology section for the Emerging Pathogens Initiative (i.e., director, supervisor, or technologist)

Total Quality Improvement/Quality Improvement/Quality Assurance (TQI/QI/QA) staff (or person at the facility with similar function)

## Periodic Reviews

The Lab Search/Extract parameter descriptions require an ongoing review process (as directed by (VAHQ) Infectious Disease Program and NHCP Offices). The person(s) participating in the ongoing review process is responsible for ensuring the following requirements are kept current:

Periodic reviews of the ICD codes

Periodic reviews of the Lab Search/Extract Parameter Setup

Annual review of the NCHP parameter descriptions

## Backup Routines

It is highly recommended that a backup of the transport global be performed before installing Patch LR\*5.2\*175.

## Test Sites

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software application was tested by the following sites:

Test Sites	Test Phases	Platform
Asheville, NC VAMC	Alpha	DEC/Alpha
Brockton/West Roxbury, MA VAMC	Alpha	DEC/Alpha
Durham, NC VAMC	Beta	DEC/Alpha
Lexington, KY VAMC	Beta	DEC/Alpha
Long Beach, CA VAMC	Beta	DEC/Alpha
Muskogee, OK VAMC	Beta	Open M

## Kernel Installation and Distribution System (KIDS)

The Laboratory Search/Extract Patch LR\*5.2\*175 distribution is using KIDS. For further instructions on using KIDS, please refer to the Kernel Version 8.0 Systems Manual.

## Database Integration Agreements (DBIA)

The following new DBIA was approved for VistA Laboratory ICD-10 PTF Modifications patch LR\*5.2\*442:

- **DBIA#6130**

The following DBIA were approved for **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175:

DBIA #418      DBIA #1372      DBIA #1888      DBIA #2488



## Namespace

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 namespace is Laboratory's LR.

## Routine List

**NOTE: The below routine list includes three routines (LRAPQAT1, LREPI3 and LREPI5) that have been modified for patch LR\*5.2\*442.**

**NOTE: The below routine list includes three routines (LR421P, LREPICD, and LRESPIXDG) that have been added for Patch LR\*5.2\*421.**

- LR132
- LR132P
- LR175
- LR175P
- LR421P
- LRAPQAT1
- LREPI
- LREPI1A
- LREPI2
- LREPI3
- LREPI4
- LREPI5
- LREPIAK
- LREPICD
- LREPICY
- LREPILK
- LREPIPH
- LREPIPI
- LREPIRM
- LREPIRN
- LREPIRP
- LREPIRP3
- LREPIRP5
- LREPIRP7
- LREPIRS1
- LREPIRS3
- LRESPIXDG

## VISTA Software Requirements

Packages	Versions (or Greater)
VA FileMan	21 (with patches installed)
Kernel	8.0 (with patches installed)
Laboratory	5.2 (with patches installed)
PIMS	5.3 (with patches installed)
Health Level 7	1.6 (with patches installed)
Social Work	3.0 (with patches installed)
MailMan	7.1 (with patches installed)

## Patches Required

**NOTE:** For a list of patches required to install patch LR\*5.2\*442, please refer to the ICD-10 PTF Modifications Installation Guide:  
<http://www.va.gov/vdl/application.asp?appid=118>

**NOTE:** For a list of patches required to install patch LR\*5.2\*421, please refer to the ICD-10 Release Notes.

Prior to the installation of Patch LR\*5.2\*175, the following patches **MUST** be installed:

Packages	Patches
MailMan V.7.1	XM*DBA*116 (NCH-Lab Domain)
Health Level Seven V. 1.6	HL*1.6*34
Laboratory V. 5.2	LR*5.2*132
	LR*5.2*157

## Health Level Seven (HL7)

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software uses the **VISTA** HL7 software to transmit data to the AITC database. The data resides in the AITC database. The data is use for assisting public health surveillance activities with new health care initiatives.

## Domain

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software application Domain name is XXX@Q-NCH. Sites **must** install Patch XM\*DBA\*116. See VA MailMan V. 7.1 Manual for instruction on how to set up the domain after the patch has been installed.

## Protocols

**LRNCH** is the Event Driver protocol used for sending the NCH HL7 messages to the Austin Automation Center.

**LRNCH SEND CLIENT** is the Client server used for sending NCH HL7 mail messages to the Austin Automation Center.

**LRNCHAITC** is the Subscriber server used for processing the NCH Acknowledgment mail messages received from the Austin Automation Center.

## Mail Groups

The Office of the Director (00) at each facility should determine responsible individual(s)/function(s) for the mail groups. The Laboratory Search/Extract Patch LR\*5.2\*175 creates the following two mail groups during the installation process.

**LR NCH:** This mail group is used to transmit the NCH HL7 messages derived from the parameters defined in the LAB SEARCH/EXTRACT file (#69.5) to the AITC database.

**LR NCH-REPORT:** This mail group delivers the Verification Report mail message and Acknowledgment mail messages (i.e., human readable format) to the recipients assigned to this mail group.

## Laboratory Search/Extract Menu and Options

The Lab Search/Extract Primary Menu [LREPI SEARCH EXTRACT MENU] is the only menu included with the VISTA Laboratory Search/Extract software. The menu contains four options. The Lab Search/Extract Nightly Task [LREPI NIGHTLY TASK] is a stand-alone option that is included with the software. There are no locks or security keys associated with the menu or options.

**NOTE:** The Lab Search/Extract Primary Menu [LREPI SEARCH EXTRACT MENU] options are using VA FileMan screen displays, referred to as ScreenMan. For detailed instructions on how to use the screens please review the VA FileMan V. 21.0 User Manual, Section 6 ScreenMan.

**Lab Search/Extract Primary Menu [LREPI SEARCH EXTRACT MENU]:** This is the primary menu that contains the Lab Search/Extract options.

**Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCE MANUAL RUN] option:** This option is used to run the EPI and NCH extracts manually.

**Antimicrobial Link Update [LREPILK] option:** This option allows the user to link the ANTIMICROBIAL SUSCEPIBILTY file (#62.06) with the WKLD CODE file (#64).

**NOTE:** Please see the Appendix C section of this guide for instructions on “How to Link Antimicrobial Entries to Workload Code Entries” using the Antimicrobial Link Update [LREPILK] option.

**Lab Search/Extract Parameter Setup [LREPI PARAMETER SETUP]option:** This option allows the users to setup the EPI and NCH parameter search/extract criteria.

**Lab Search/Extract Protocol Edit [LREPI PROTOCOL EDIT] option:** Use this option to edit the LAB SEARCH/EXTRACT PROTOCOL file (#69.4).

**Lab Search/Extract Nightly Task [LREPI NIGHTLY TASK] option:** The Lab Search/Extract Nightly Task [LREPI NIGHTLY TASK] option **must** be scheduled to run each night by TaskMan. This option will build a HL7 Message and send it to the defined locations specified by the EPI and NCH protocols.

## Laboratory Search/Extract Parameters Input Screen

### Prompts Definitions

Laboratory Search/Extract Parameters Input Screen Prompts	Laboratory Search/Extract Parameters Input Screen Prompt Definitions
Laboratory Test (s)	Consider these synonymous with, chemistry, serology, hematology, and “blood/serum” tests. Results anticipated to be found here would have had a test done under the chemistry/hematology accession areas, even if physically performed in microbiology and other areas. Select tests from the LABORATORY TEST file (#60).
Indicator	Select the code that will determine how to match lab results. ‘1’ FOR Use Reference Ranges ‘2’ FOR Contains ‘3’ FOR Greater Than ‘4’ FOR Less Than ‘5’ FOR Equal To
Value	Positive, etc. Answer must be 1-15 characters in length. This is a Free Text field.
ICD-9	ICD-9 standardized code used nationwide in federal and non-federal/private health care facilities. Select from the ICD-9 DIAGNOSIS file (#80).
ICD-10	ICD-10 standardized code used nationwide in federal and non-federal/private health care facilities. Select from the ICD DIAGNOSIS file (#80).
ICD Description	Title of ICDM diagnosis
Selected Etiology	Consider synonymous with organism, final microbial diagnosis/isolate. Select from the ETIOLOGY FIELD file (#61.2)
Selected SNOMED codes	Answer with SNOMED CODES You may enter a new SNOMED CODE, if you wish. Answer must be 1-15 characters in length.
Antimicrobial Susceptibility	Enter the Antimicrobial that will be used in screening out sensitive Etiologies (e.g., “Vancomycin” for Vancomycin Resistant Enterococcus). Select from the ANTIMICROBIAL SUSCEPTIBILITY file (#62.6).
NLT Code	Displays the associated NLT code if linked. If no NLT Code is displayed use the Antimicrobial Link Update option.
NLT Description	Displays the Description of the linked NLT code.
Topography Selection	Enter a date to screen out patients born before the date entered. Examples of Valid Dates: JAN 20 1957 or 20 JAN 57 or 1/20/57 or 012057 T (for TODAY), T+1 (for TOMORROW), T+2, T+7, etc. T-1 (for YESTERDAY), T-3W (for 3 WEEKS AGO), etc.
Include	Selection of Topography screens all others out except the ones selected. For “ALL” leave blank. Not to be used in conjunction with the exclude Topography selection. Select from the TOPOGRAPHY file (#61).

Laboratory Search/Extract Parameters Input Screen Prompts	Laboratory Search/Extract Parameters Input Screen Prompt Definitions
Exclude	Select the Topography to screen out. Not to be used in conjunction with the Include Topography selection. Select from the TOPOGRAPHY file (#61).
First Encounter:	Limits the output to the first encounter for the patient. Otherwise list all encounters. Choose: '1' FOR YES '0' FOR NO
Follow PTF:	Indicates if the PTF record will be followed until a discharge has been entered. Choose: '1' FOR YES '0' FOR NO
Before Date Of Birth:	Enter a date to screen out patients born before the date entered. Examples of Valid Dates: JAN 20 1957 or 20 JAN 57 or 1/20/57 or 012057 T (for TODAY), T+1 (for TOMORROW), T+2, T+7, etc. T-1 (for YESTERDAY), T-3W (for 3 WEEKS AGO), etc.
After Date Of Birth:	A birthrate to screen patients (i.e., patients DOB after 1/1/1950).
Select SEX:	<b>FOR FUTURE USE ONLY.</b>
Run Date:	Date that the last Auto Search/Extract processed.
Protocol:	Defines the protocol used to define the output messages. Select from the LAB SEARCH/EXTRACT PROTOCOL file (#69.4).
Run Cycle:	Enter the date that the last Auto Search/Extract processed.
Lag Days:	Enter the numbers of lag days desired. The number allowed will vary depending on the cycle type.
General Description:	To review or edit the General Description prompt use the <b>&lt;Tab&gt;</b> key.

## VISTA Lab Search/Extract Files Data Dictionaries

The following are examples of the LABORATORY SEARCH/EXTRACT files standard Data Dictionaries:

### LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4)

STANDARD DATA DICTIONARY #69.4 -- LAB SEARCH/EXTRACT PROTOCOL FILE 06/11/98  
 PAGE 1  
 STORED IN ^LAB(69.4, (2 ENTRIES) SITE: BROCKTON/W.ROXBURY UCI: VAH,ROU  
 (VERSION 5.2)

DATA ELEMENT	NAME TITLE	GLOBAL LOCATION	DATA TYPE
-----------------	---------------	--------------------	--------------

-----  
 This file contains additional search and extract parameters to be used by the Laboratory Search/Extract software. These parameters are not specific to entries in LAB SEARCH/EXTRACT file (#69.5), but are specific to the Protocol (refer to Protocol file #101) being used.

POINTED TO BY: PROTOCOL field (#12) of the LAB SEARCH/EXTRACT File (#69.5)

CROSS  
 REFERENCED BY: PROTOCOL(B)

69.4,.01	PROTOCOL	0;1 POINTER TO PROTOCOL FILE (#101) (Required)
	INPUT TRANSFORM:	S DINUM=X
	LAST EDITED:	NOV 08, 1996
	HELP-PROMPT:	Select the Protocol that will be used to build the message.
	DESCRIPTION:	This is the entry in the Protocol file (#101) that will be used to build the HL7 Message. This allows additional parameters to be associated with the protocol.
	NOTES:	XXXX--CAN'T BE ALTERED EXCEPT BY PROGRAMMER
	CROSS-REFERENCE:	69.4^B 1)= S ^LAB(69.4,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.4,"B",\$E(X,1,30),DA)

69.4,1	Report Mail Group	0;2 POINTER TO MAIL GROUP FILE (#3.8)
	LAST EDITED:	NOV 08, 1996
	HELP-PROMPT:	Select which mail group will receive the verification report.
	DESCRIPTION:	This is the mail group that will receive the verification report.

69.4,2	Message Size	0;3 NUMBER
	INPUT TRANSFORM:	K:+X'=X!(X>999999)!(X<100)!(X?.E1"."1N.N) X

## Installation Instructions

```
LAST EDITED:      DEC 04, 1996
HELP-PROMPT:      Enter the maximum number of lines a message
                   may contain (100 to 999999) before it breaks
                   into another message.
DESCRIPTION:       This determines how big the HL7 message will
                   be before it breaks into another message.
```

```
69.4,3      Send Alert      0;4 SET
```

```

LAST EDITED:      '1' FOR YES;
HELP-PROMPT:      JUL 31, 1997
                   Enter '1' or 'YES' to send an alert when a
                   search is run automatically. Leave this field
                   blank if no alert is desired.
DESCRIPTION:      When this field is set to '1' or 'YES' an
                   alert is transmitted to the recipients listed
                   in the RECEIVE ALERTS field when a search is
                   run automatically for this protocol.

```

69.4,4	Title	0;5 FREE TEXT
--------	-------	---------------

```

INPUT TRANSFORM:  K:$L(X)>40!($L(X)<2) X
LAST EDITED:      AUG 01, 1997
HELP-PROMPT:      Enter a new title (2 to 40 characters) for the
                   Verification report. If this field is left
                   blank the default title "Verification Report"
                   will be used.
DESCRIPTION:      This field will over-ride the default title
                   for the Verification Report.

```

69.4,5	Receive Alerts	1;0	VARIABLE POINTER	Multiple #69.45
	DESCRIPTION:	This is the User(s) or Group that will receive alerts stating the extract has been run. automatically.		

69.45,.01	Receive Alerts	0;1	VARIABLE POINTER
		(Multiply asked)	

FILE	ORDER	PREFIX	LAYGO	MESSAGE
200	1	U	n	Responsible User
3.8	2	G	n	Responsible Group

```

LAST EDITED:      AUG 05, 1997

HELP-PROMPT:      Select the User(s) or Group responsible for
                   receiving alerts that the search has been
                   run automatically.
DESCRIPTION:       This is the user(s) or Group that will
                   receive alerts that the search has been run
                   automatically.
CROSS-REFERENCE:   69.45^B
                   1)= S
                   ^LAB(69.4,DA(1),1,"B",$E(X,1,30),DA)="
                   2)= K ^LAB(69.4,DA(1),1,"B",$E(X,1,30),DA)
                   3)= Required Index for Variable Pointer

```

FILES POINTED TO

FIELDS

MAIL GROUP (#3.8)

Report Mail Group (#1)



## Installation Instructions

```

PROTOCOL (#101)                                PROTOCOL (#.01)

INPUT TEMPLATE(S):

PRINT TEMPLATE(S):
CAPTIONED                                         USER #0

SORT TEMPLATE(S):

FORM(S)/BLOCK(S):
LREPIPROT                                         SEP 05, 1997@13:34  USER #0
  LRPROTHEAD                                     DD #69.4
  LRPROT1                                       DD #69.4
  LRPROT2                                       DD #69.45
```

LAB SEARCH/EXTRACT file (#69.5)

STANDARD DATA DICTIONARY #69.5 -- LAB SEARCH/EXTRACT FILE 06/11/98 PAGE 1  
 STORED IN ^LAB(69.5, (17 ENTRIES) SITE: BROCKTON/W.ROXBURY UCI: VAH,ROU  
 (VERSION 5.2)

DATA ELEMENT	NAME TITLE	GLOBAL LOCATION	DATA TYPE
-----------------	---------------	--------------------	--------------

-----  
 This file contains search criteria used by the Laboratory Search/Extract software. This file should only be edited using Lab Search/Extract Parameter Setup option [LREPI PARAMETER SETUP] provided with this software.

## CROSS

REFERENCED BY: NAME(B), REFERENCE NUMBER(C)

69.5,.01	NAME	0;1 FREE TEXT (Required)
	INPUT TRANSFORM:	K:\$L(X)>50!(\$L(X)<3)!'(X'?1P.E)!(X'?.ANP) X
	LAST EDITED:	DEC 17, 1996
	HELP-PROMPT:	Enter a Name (3 to 50 characters) for the Search/Extract parameter you are defining.
	DESCRIPTION:	This is the name of the Search/Extract parameter you are defining.
	NOTES:	XXXX--CAN'T BE ALTERED EXCEPT BY PROGRAMMER
	CROSS-REFERENCE:	69.5^B 1)= S ^LAB(69.5,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,"B",\$E(X,1,30),DA)
69.5,.05	REFERENCE NUMBER	0;9 NUMBER
	INPUT TRANSFORM:	K:+X'=X!(X>999)!(X<1)!(X?.E1"."1N.N)!(X'>99)!\$D(^LAB(69.5,"C",X)) X
	LAST EDITED:	NOV 29, 1996
	HELP-PROMPT:	Type a unique number between 100 and 999 to identify this entry. Numbers from 1 to 99 are reserved for future use.
	DESCRIPTION:	This is a unique number used to identify this entry.
	NOTES:	UNEDITABLE XXXX--CAN'T BE ALTERED EXCEPT BY PROGRAMMER
	CROSS-REFERENCE:	69.5^C 1)= S ^LAB(69.5,"C",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,"C",\$E(X,1,30),DA)
69.5,1	ACTIVE	0;2 SET
		'1' FOR YES; '0' FOR NO;
	LAST EDITED:	JUN 09, 1998
	HELP-PROMPT:	'1' or 'YES' indicates that this is an active entry.
	DESCRIPTION:	This defines if this entry is active or not.

## Installation Instructions

69.5,2	LAB TEST	1;0 POINTER Multiple #69.52
	DESCRIPTION:	This is the test that is to be searched for and retrieved.
69.52,.01	LAB TEST	0;1 POINTER TO LABORATORY TEST FILE (#60) (Multiply asked)
	INPUT TRANSFORM:	S DIC("S")="I \$P(\$G(^ (0)),U,4)="CH"!"(\$P(\$G( ^ (0)),U,4)="CY"" D ^DIC K DIC S DIC=DIE,X= +Y K:Y<0 X
	LAST EDITED:	FEB 04, 1998
	HELP-PROMPT:	Enter the lab test that is to be searched for and retrieved. Consider this synonymous with chemistry, serology, hematology "blood/serum" or Cytology tests.
	DESCRIPTION:	This is the lab test that is to be searched for and retrieved.
	SCREEN:	S DIC("S")="I \$P(\$G(^ (0)),U,4)="CH"!"(\$P(\$G( ^ (0)),U,4)="CY""
	EXPLANATION:	Only CH subscripts are selectable.
	NOTES:	XXXX--CAN'T BE ALTERED EXCEPT BY PROGRAMMER
	CROSS-REFERENCE:	69.52^B 1)= S ^LAB(69.5,DA(1),1,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,DA(1),1,"B",\$E(X,1,30),DA)
69.52,1	INDICATOR	0;2 SET
		'1' FOR Use Reference Ranges; '2' FOR Contains; '3' FOR Greater Than; '4' FOR Less Than; '5' FOR Equal To;
	LAST EDITED:	FEB 05, 1997
	HELP-PROMPT:	Select the code that will determine how to match lab results.
	DESCRIPTION:	This indicates if the search for the lab test is conditional.
69.52,2	INDICATED VALUE	0;3 FREE TEXT
	INPUT TRANSFORM:	K:\$L(X)>30!(\$L(X)<1) X
	LAST EDITED:	FEB 05, 1998
	HELP-PROMPT:	Enter the data to be compared using the INDICATOR field.
	DESCRIPTION:	If the search is conditional this defines the criteria.
69.5,3	ETIOLOGY	2;0 POINTER Multiple #69.53
	DESCRIPTION:	This defines the Etiology to search for.
69.53,.01	ETIOLOGY	0;1 POINTER TO ETIOLOGY FIELD FILE (#61.2) (Multiply asked)
LAST EDITED:	AUG 29, 1996	

## Installation Instructions

	HELP-PROMPT:	Select the Etiology to search for.
	DESCRIPTION:	This defines the Etiology to search for. Select the appropriate Etiology.
	CROSS-REFERENCE:	69.53^B 1)= S ^LAB(69.5,DA(1),2,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,DA(1),2,"B",\$E(X,1,30),DA)
69.5,4	ICD DIAGNOSIS	3;0 POINTER Multiple #69.54
	LAST EDITED:	JUN 12, 2012
	DESCRIPTION:	This defines the ICD to search for.
69.54,.01	ICD DIAGNOSIS	0;1 POINTER TO ICD DIAGNOSIS FILE (#80) (Multiply asked)
	LAST EDITED:	SEP 17, 2012
	HELP-PROMPT:	Select the ICD standardized code, used nationwide in federal and non-federal private health care facilities, to be included in the search.
	DESCRIPTION:	This defines an ICD for use in emerging pathogens data search/extract.
	CROSS-REFERENCE:	69.54^B 1)= S LAB(69.5,DA(1),3,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,DA(1),3,"B",\$E(X,1,30),DA)
69.54,1	CODING SYSTEM	0;2 POINTER TO ICD CODING SYSTEMS FILE (#80.4)
	LAST EDITED:	AUG 10, 2012
	HELP-PROMPT:	Enter the applicable ICD coding system for this ICD.
69.5,5	ANTIMICROBIAL SUSCEPTIBILITY	4;0 POINTER Multiple #69.55
	LAST EDITED:	JAN 22, 1997
	DESCRIPTION:	This determines that if any of the Etiologies selected are to be resistant to any Antimicrobials.
69.55,.01	ANTIMICROBIAL SUSCEPTIBILITY	0;1 POINTER TO ANTIMICROBIAL SUSCEPTIBILITY FILE (#62.06) (Multiply asked)
	LAST EDITED:	JAN 22, 1997
	HELP-PROMPT:	Enter the Antimicrobial that will be used in screening out sensitive Etiologies.
	DESCRIPTION:	This determines that if any of the selected are to be resistant to Antimicrobials. Select the appropriate Antimicrobials to screen out the Etiologies.
	CROSS-REFERENCE:	69.55^B 1)= S ^LAB(69.5,DA(1),4,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,DA(1),4,"B",\$E(X,1,30),DA)
69.5,6	INCLUDED SITES	5;0 POINTER Multiple #69.56

## Installation Instructions

```

                LAST EDITED:      OCT 04, 1996
                DESCRIPTION:      This determines what Topography to screen for.

69.56,.01      TOPOGRAPHY        0;1 POINTER TO TOPOGRAPHY FIELD FILE (#61)
                                (Multiply asked)

                LAST EDITED:      OCT 04, 1996
                HELP-PROMPT:      Selection of a Topography screens all others
                                out except the ones selected. For "ALL"
                                leave blank. Not to be used in conjunction
                                with the exclude Topography selection.
                DESCRIPTION:      This determines what Topography to screen
                                for. Select the appropriate Topography to
                                include in the extract.

                CROSS-REFERENCE:  69.56^B
                                1)= S
                                ^LAB(69.5,DA(1),5,"B",$E(X,1,30),DA)=" "
                                2)= K ^LAB(69.5,DA(1),5,"B",$E(X,1,30),DA)

69.5,7         EXCLUDED SITES    6;0 POINTER Multiple #69.57

                DESCRIPTION:      This determines what Topography to screen
                                out.

69.57,.01      TOPOGRAPHY        0;1 POINTER TO TOPOGRAPHY FIELD FILE (#61)
                                (Multiply asked)

                LAST EDITED:      OCT 04, 1996
                HELP-PROMPT:      Select the Topography to screen out. Not to
                                be used in conjunction with the Include
                                Topography selection.
                DESCRIPTION:      This determines what Topography to screen
                                out. Select the appropriate Topography to be
                                excluded from the extract.

                CROSS-REFERENCE:  69.57^B
                                1)= S
                                ^LAB(69.5,DA(1),6,"B",$E(X,1,30),DA)=" "
                                2)= K ^LAB(69.5,DA(1),6,"B",$E(X,1,30),DA)

69.5,8         SNOMED CODES      9;0 Multiple #69.58

                DESCRIPTION:      This is a list of SNOMED codes to be included
                                in the search.

69.58,.01      SNOMED CODES      0;1 FREE TEXT (Multiply asked)

                INPUT TRANSFORM:  K:$L(X)>15!($L(X)<1) X
                LAST EDITED:      AUG 20, 1997
                HELP-PROMPT:      Enter any SNOMED codes to be included in the
                                search. Enter one code per line.
                DESCRIPTION:      This is a list of SNOMED code to be included
                                in the search.

                CROSS-REFERENCE:  69.58^B
                                1)= S
                                ^LAB(69.5,DA(1),9,"B",$E(X,1,30),DA)=" "
                                2)= K ^LAB(69.5,DA(1),9,"B",$E(X,1,30),DA)

69.5,9         RUN DATE          0;4 DATE

                INPUT TRANSFORM:  S %DT="ESTX" D ^%DT S X=Y K:Y<1 X
                LAST EDITED:      FEB 06, 1997
                HELP-PROMPT:      Enter the date that the last Auto

```

## Installation Instructions

	DESCRIPTION:	Search/Extract processed. The date that the last Auto Search/Extract processed
69.5,10	CYCLE	0;5 SET  'M' FOR MONTHLY; 'D' FOR DAILY; LAST EDITED: DEC 03, 1997 HELP-PROMPT: Select whether this parameter is to be used 'D'aily or 'M'onthly DESCRIPTION: This field defines how often the this entry is acted upon by the nightly task.  CROSS-REFERENCE: ^^TRIGGER^69.5^10.5 1)= K DIV S DIV=X,D0=DA,DIV(0)=D0 S Y(1)=\$S(\$D( ^LAB(69.5,D0,0)):^(0),1:"") S X=\$P(Y(1),U,3),X=X S DIU=X K Y S X="" X ^DD(69.5,10,1,1,1.4) 1.4)= S DIH=\$S(\$D(^LAB(69.5,DIV(0),0)):^(0),1:" "),DIV=X S \$P(^^(0),U,3)=DIV,DIH=69.5,DIG=10.5 D^DICR:\$O(^DD(DIH,DIG,1,0))>0  2)= K DIV S DIV=X,D0=DA,DIV(0)=D0 S Y(1)=\$S(\$D(^LAB(69.5,D0,0)):^(0),1:"") S X=\$P(Y(1),U,3),X=X S DIU=X K Y S X="" X ^DD(69.5,10,1,1,2.4) 2.4)= S DIH=\$S(\$D(^LAB(69.5,DIV(0),0)):^(0),1: " "),DIV=X S \$P(^^(0),U,3)=DIV,DIH=69.5,DIG=10.5 D ^DICR:\$O(^DD(DIH,DIG,1,0))>0 CREATE VALUE)= "" DELETE VALUE)= "" FIELD)= LAG DAY
69.5,10.5	LAG DAYS	0;3 NUMBER  INPUT TRANSFORM: S CYC=\$P(^LAB(69.5,DA,0),U,5),MAX=\$S(CYC="D":99 ,CYC="M":26,1:0) K:+X'=X!(X>MAX)!(X<1)!(X?.E1".  "1N.N) X K MAX,CYC LAST EDITED: SEP 04, 1997 HELP-PROMPT: Enter the numbers of lag days desired. The amount allowed will vary depending on the cycle type. DESCRIPTION: This is the number of days that the search should wait after the test was ordered to ensure that the results have been entered in the system.  NOTES: XXXX--CAN'T BE ALTERED EXCEPT BY PROGRAMMER TRIGGERED by the CYCLE field of the LAB SEARCH/EXTRACT File
69.5,11	FIRST ENCOUNTER	0;6 SET  '1' FOR YES; '0' FOR NO; LAST EDITED: DEC 30, 1996 HELP-PROMPT: Enter '1' to limit the output to the first encounter for the patient. '0' will list all encounters.

## Installation Instructions

	DESCRIPTION:	This determines if after the first encounter is found and extracted should sequential encounters be extracted.
69.5,12	PROTOCOL	0;7 POINTER TO LAB SEARCH/EXTRACT PROTOCOL FILE (#69.4)
	LAST EDITED:	NOV 08, 1996
	HELP-PROMPT:	Select the Protocol to be used to define the output messages.
	DESCRIPTION:	This defines what protocol is associated with the parameters.
69.5,13	FOLLOW PTF	0;8 SET
		'1' FOR YES; '0' FOR NO;
	LAST EDITED:	OCT 17, 1996
	HELP-PROMPT:	'YES' or '1' indicates that the PTF record will be followed until a discharge has been entered.
	DESCRIPTION:	This determines whether the discharge information should be updated upon discharge if a inpatient encounter does not have an existing discharge.
69.5,14	PTF	7;0 POINTER Multiple #69.514 (Add New Entry without Asking)
	DESCRIPTION:	This is the PTF file entry to be followed.
69.514,.01	PTF	0;1 POINTER TO PTF FILE (#45)
	LAST EDITED:	OCT 17, 1996
	HELP-PROMPT:	Enter the PTF file entry to be followed.
	DESCRIPTION:	This is the PTF file entry to be followed.
	CROSS-REFERENCE:	69.514^B 1)= S ^LAB(69.5,DA(1),7,"B",\$E(X,1,30),DA)=" " 2)= K ^LAB(69.5,DA(1),7,"B",\$E(X,1,30),DA)
69.514,1	DATE	0;2 DATE
	INPUT TRANSFORM:	S %DT="E" D ^%DT S X=Y K:Y<1 X
	LAST EDITED:	DEC 31, 1996
	HELP-PROMPT:	Enter the date that the Inpatient discharge information was included in the report as an update.
	DESCRIPTION:	This is the date that the Inpatient discharge information was included in the report as a update.
69.5,15	Description	8;0 WORD-PROCESSING #69.515
	DESCRIPTION:	This is the general description for the entry.
69.5,16	SEX	0;10 SET
		'M' FOR MALE; 'F' FOR FEMALE; 'O' FOR OTHER;
	LAST EDITED:	NOV 17, 1997
	HELP-PROMPT:	Enter the sex code to be included in the

## Installation Instructions

```
search. Leave blank if all types are to be
included.
DESCRIPTION: This is the sex code to be included in the
search. All types will be included if this
field is left blank.

69.5,17 BEFORE DATE OF BIRTH 0;11 DATE

INPUT TRANSFORM: S %DT="EX" D ^%DT S X=Y K:Y<1 X
LAST EDITED: NOV 17, 1997
HELP-PROMPT: Enter a date to screen out patients born after
the date entered.
DESCRIPTION: Patients born after the date entered will not
be included in the report.

69.5,18 AFTER DATE OF BIRTH 0;12 DATE

INPUT TRANSFORM: S %DT="EX" D ^%DT S X=Y K:Y<1 X
LAST EDITED: NOV 17, 1997
HELP-PROMPT: A brithrate to screen patients. i.e. patients
DOB after 1/1/1950.
before the date entered.
DESCRIPTION: Persons born before the entered will not be
included in the report.
```

FILES POINTED TO	FIELDS
ANTIMICROBIAL SUSCEPTIBILITY (#62.06)	ANTIMICROBIAL
SUSCEPTIBILITY:ANTIMICROBIAL SUSCEPTIBILITY (#.01)	
ETIOLOGY FIELD (#61.2)	ETIOLOGY:ETIOLOGY (#.01)
ICD CODING SYSTEMS (#80.4)	ICD DIAGNOSIS:CODING SYSTEM (#1)
ICD DIAGNOSIS (#80)	ICD DIAGNOSIS:ICD DIAGNOSIS (#.01)
LAB SEARCH/EXTRACT PROTOCOL(#69.4)PROTOCOL (#12)	
LABORATORY TEST (#60)	LAB TEST:LAB TEST (#.01)
PTF (#45)	PTF:PTF (#.01)
TOPOGRAPHY FIELD (#61)	INCLUDED SITES:TOPOGRAPHY (#.01)
	EXCLUDED SITES:TOPOGRAPHY (#.01)

INPUT TEMPLATE(S):

PRINT TEMPLATE(S):  
CAPTIONED

USER #0

SORT TEMPLATE(S):

FORM(S)/BLOCK(S):

LREPI	OCT 07, 1996@10:13	USER #0
LREPIHEAD	DD #69.5	
LREPI2	DD #69.52	
LREPI3	DD #69.54	
LREPI1	DD #69.5	
LREPI11	DD #69.5	



## Installation Instructions

LREPI4	DD #69.53
LREPI5	DD #69.55
LREPI6	DD #69.5
LREPI12	DD #69.58
LREPI7	DD #69.5
LREPI8	DD #69.56
LREPI9	DD #69.57
LREPI13	DD #69.5
LREPI10	DD #69.5

## Installation Instructions

**NOTE:** For patch LR\*5.2\*442 installation instructions, please refer to the ICD-10 PTF Modifications Installation Guide:

<http://www.va.gov/vdl/application.asp?appid=118>

**NOTE:** For Patch LR\*5.2\*421 installation instructions, please refer to the ICD-10 Release Notes.

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 is using the Kernel Installation and Distribution System (KIDS).

**NOTE:** For further instructions on using KIDS, please refer to the Kernel V. 8.0 Systems Manual, Chapter 26, pages 393-409.

1. Use the 'INSTALL/CHECK MESSAGE' option on the PackMan menu. This option will load the KIDS package in this message onto your system.
2. Review your mapped set. If the routines are mapped, they should be removed from the mapped set at this time.
3. The patch has now been loaded into a Transport global on your system. You now need to use KIDS to install the Transport global. On the KIDS menu, under the Installation menu, use the following options:
  - Print Transport Global
  - Compare Transport Global to Current System
  - Verify Checksums in Transport Global
  - Backup a Transport Global
4. Users may remain on system and no options need to be placed out of service.
5. Installation time is less than 2 minutes during off peak hours and less than 5 minutes during peak hours.
6. Installation of this patch requires no additional memory space.
7. From the Installation Menu of the KIDS menu, run the option Install Package(s) Select the package LR\*5.2\*175 and proceed with the install.
8. If any routines were unmapped as part of step 2, they should be returned to the

mapped set once the installation has run to completion.

9. Delete routine LR175 and LR175P.
10. Remove the Emerging Pathogens Manual Run [LREPI (EPI) MANUAL RUN] option.

## Installation Instructions

**Example:** Terminal screen dialogue seen during the KIDS install (dates shown in the example will **not** be the same as those in the released version).

Select Programmer Options Option: **KID** Kernel Installation & Distribution System<RET>

Select Kernel Installation & Distribution System Option: **IN**stallation<RET>

Select Installation Option: **MAIL**Man Menu<RET>

VA MailMan 7.1 service for DOE.DAVID\_R@JDV.ISC-DALLAS.VA.GOV

You last used MailMan: 01 Jun 98 14:12

You have 1 new message.

Select MailMan Menu Option: **NEW** messages and responses<RET>

Subj: LR\*5.2\*175 [15320] 9 Jun 1998 16:50 EST 6046 Lines

From: <"NPM [#27200849]"@FORUM.VA.GOV> in 'IN' basket. Page 1 \*\*NEW\*\*

-----  
\$TXT Created by DOE DAVID R at BROCKTON.VA.GOV (KIDS) on TUESDAY, 06/09/98 at 16:15

=====

Run Date: JUN 09, 1998	Designation: LR*5.2*175
------------------------	-------------------------

Package : LR - LAB SERVICE	Priority: Mandatory
----------------------------	---------------------

Version : 5.2	Status: Under Development
---------------	---------------------------

=====

Subject: HEALTH PROMOTIONS ENHANCE TO EPI

Category:

- Routine

- Enhancement (Mandatory)

- Other

Select MESSAGE Action: IGNORE (in IN basket)//X<RET>

Select PackMan function: **6** INSTALL/CHECK MESSAGE<RET>

Line 3 Message #13928486 Unloading KIDS Distribution LR\*5.2\*175

LR\*5.2\*175

## Installation Instructions

Want to Continue with Load? YES//<RET>

Loading Distribution.

LR\*5.2\*175

Will first run the Environment Check Routine, LR175

Select PackMan function:<RET>

Select MESSAGE Action: IGNORE (in IN basket)// ^ <RET>

Ignored

Select MailMan Menu Option:<RET>

Select Installation Option: **IN**stall Package(s)<RET>

Select INSTALL NAME:       **LR\*5.2\*175**               Loaded from Distribution  
7/10/98@14:31:58

=> LR\*5.2\*175

This Distribution was loaded on Jul 10, 1998@14:31:58 with header of

LR\*5.2\*175

It consisted of the following Install(s):

LR\*5.2\*175

LR\*5.2\*175

Will first run the Environment Check Routine, LR175

Environment Check is Ok ---

Install Questions for LR\*5.2\*175

Incoming Files:

69.4           LAB SEARCH/EXTRACT PROTOCOL  
\*BUT YOU ALREADY HAVE 'EMERGING PATH' PROTOCOL' AS File #69.4!  
Shall I write over your EMERGING PATH' PROTOCOL File? YES//<RET>

69.5           LAB SEARCH/EXTRACT  
\*BUT YOU ALREADY HAVE 'EMERGING PATHOGENSS'AS File #69.5!  
Shall I write over your EMERGING PATHOGENSS File? YES//<RET>

Incoming Mail Groups:

## Installation Instructions

Enter the Coordinator for Mail Group 'LR NCH': **DOE BLOW** BD

Enter the Coordinator for Mail Group 'LR NCH-REPORT': **DOE BLOW** BD

Want to DISABLE Scheduled Options, Menu Options, and Protocols? YES//**NO<RET>**

Enter the Device you want to print the Install messages.

You can queue the install by enter a 'Q' at the device prompt.  
Enter a '^' to abort the install.

DEVICE: HOME// [**Enter your site printer's name**]**<RET>**

Adding Protocal 'LRNCH' to the LAB SEARCH/EXTRACT File (69.5)

Updating Routine file

Updating KIDS files

LR\*5.2\*175 Installed.

Jul 10, 1998@14:33:47.

Install Message sent #13928998

-----

100%	25	50	75
Complete			

Install Completed

## Post Installation Instructions

**NOTE:** There are no post installation instructions for LR\*5.2\*442.

**NOTE:** For Patch LR\*5.2\*421 post installation instructions, please refer to the ICD-10 Release Notes for LR\*5.2\*421.

It is highly recommended that the Laboratory Search/Extract Patch LR\*5.2\*175 post installation instructions are done in the following order:

### DSM/Alpha and Open M Sites

If you have disabled journaling, you may now re-enable journaling. If using a mapped system, rebuild the map set now.

**IRM - Step #1:** Ensure that **XXX@Q-NCH** domain is defined as instructed by VA MailMan V.7.1, Patch XM\*DBA\*116.

**IRM – Step #2:** Add designated members to the LR NCH and LR NCH-Report mail groups.

**NOTE:** The Office of the Director should designate the recipients of these mail groups.

### **Example:** LR NCH mail group setup

```
NAME: LRNCH          TYPE: public
ALLOW SELF ENROLLMENT?: NO
LAST REFERENCED: MAY 01, 1998 RESTRICTIONS: UNRESTRICTED
COORDINATOR: DOE BLOW
MEMBER: [Enter Member's Name]
MEMBER: [Enter Member's Name]
DESCRIPTION: LRNCH Mail group
REMOTE MEMBERS:XXX@Q-NCH
```

### **Example:** LR NCH-Report mail group setup

```
NAME: LRNCH-REPORT   TYPE: public
ALLOW SELF ENROLLMENT?: NO REFERENCE COUNT: 56
LAST REFERENCED: MAY 01, 1998 RESTRICTIONS: UNRESTRICTED
COORDINATOR: DOE BLOW
MEMBER: [Enter Member's Name]
MEMBER: [Enter Member's Name]
DESCRIPTION: LR NCH-REPORT Mail group
```

**IRM - Step #3:** Assign the Lab Search/Extract Primary [LREPI SEARCH EXTRACT MENU] menu to the designated user(s).

**NOTE:** It is highly recommended that a LIM, TQI/QA/QI, and a representative from the Microbiology section for EPI (e.g., director, supervisor, or technologist) be assigned the Lab Search/Extract Primary [LREPI SEARCH EXTRACT MENU] menu. These individuals should be responsible for initially setting up the NCH parameters and performing periodic reviews to ensure that the parameters remain current.

**LIM - Step #4:** Use the Lab Search/Extract Parameter Setup [LREPI PARAMETER SETUP] option to set up the NCH Cholesterol and Pap smear parameters (e.g., in Appendix B). Lab test names for Cholesterol and Pap smears may differ at each site.

**Note:** The user assigned the Lab Search/Extract Primary Menu [LREPI SEARCH EXTRACT MENU] should be designated to set up the parameters.

To verify that the fields in the LAB SEARCH/EXTRACT file (#69.5) are correctly defined, do a VA FileMan inquiry of the file. The fields should be defined in File (#69.5) as shown in the example below.

**Example: NCH Cholesterol setup**

NAME: NCH CHOLESTEROL	ACTIVE: <b>NO&lt;RET&gt;</b>
LAG DAYS: <b>10&lt;RET&gt;</b>	RUN DATE: JUN 22, 1998
CYCLE: <b>DAILY&lt;RET&gt;</b>	PROTOCOL: <b>LRNCH&lt;RET&gt;</b>
FOLLOW PTF: <b>YES&lt;RET&gt;</b>	REFERENCE NUMBER: 50
LAB TEST: <b>CHOLESTEROL&lt;RET&gt;</b>	
LAB TEST: <b>HDL CHOLESTEROL&lt;RET&gt;</b>	
LAB TEST: <b>TRIGLYCERIDE&lt;RET&gt;</b>	
LAB TEST: <b>LDL CHOLESTEROL&lt;RET&gt;</b>	
LAB TEST: <b>HEALTH FAIR-CHOLESTEK PANEL&lt;RET&gt;</b>	

**Example: NCH Pap Smear setup**

NAME: NCH PAP SMEAR	ACTIVE: NO
LAG DAYS: 10	RUN DATE: JUN 22, 1998
CYCLE: DAILY	PROTOCOL: LRNCH
FOLLOW PTF: YES	REFERENCE NUMBER: 51
SEX: FEMALE	
LAB TEST: PAP SMEAR	
TOPOGRAPHY: VAGINAL CYTOLOGIC MATERIAL	
TOPOGRAPHY: CERVICAL CYTOLOGIC MATERIAL	



**IRM - Step #5:** Set up the EPI and LRNCH protocols using the Lab Search/Extract Protocol Edit [LREPI PROTOCOL EDIT] option.

**Example: EPI Protocol Setup**

Protocol Parameters Setup Definition

---

PROTOCOL: LREPI

Title: Emerging Pathogens Initiative (EPI)

Message Size: **32000**

Report Mail Group: **EPI-REPORT**

Send Alert: **YES**

Send Alert To

**DOE,Jane**

**Example: LRNCH Protocol Setup**

Protocol Parameters Setup Definition

---

PROTOCOL: LRNCH

Title: National Center for Health Promotion

Message Size: **32000**

Report Mail Group: **LR NCH-REPORT**

Send Alert: **YES**

Send Alert To

**DOE,Jane**

**IRM - Step #6:** Start the Lower Level Protocol of the HL7 V. 1.6 background job.

**Example:**

Select Systems Manager Menu Option: **HL7 Main<RET>Menu**

- 1 V1.5 OPTIONS
- 2 V1.6 OPTIONS ...
- 3 Activate/Inactivate Application
- 4 Print/Display Menu ...
- 5 Purge Message Text File Entries

Select HL7 Main Menu Option: **2<RET>V1.6 OPTIONS**

- 1 Communications Server
- 2 Interface Workbench
- 3 Message Requeuer

Select V1.6 OPTIONS Option: **1<RET> Communications Server**

## Post-Installation Instructions

- 1 Edit Communication Server parameters
- 2 Manage incoming & outgoing filers ...
- 3 Monitor incoming & outgoing filers
- 4 Start LLP
- 5 Stop LLP
- 6 Systems Link Monitor
- 7 Logical Link Queue Management ...
- 8 Report

Select Communications Server Option: **4<RET>** Start LLP

This option is used to launch the lower level protocol for the appropriate device. Please select the node with which you want to communicate.

Select HL LOGICAL LINK NODE: **L RNCH<RET>**

The LLP was last shutdown on JAN 30, 1997 12:06:19.

Select one of the following:

- |   |            |
|---|------------|
| F | FOREGROUND |
| B | BACKGROUND |
| Q | QUIT       |

Method for running the receiver: **B//<RET>** ACKGROUND

Job was queued as 131225.

### **IRM - Step #7:** Schedule the Lab Search/Extract Nightly Task [LREPI NIGHTLY TASK] option to run nightly.

**Example:** Select OPTION to schedule or reschedule: **LREPI NIGHTLY TASK<RET>**

Lab Search/Extract Nightly Task

...OK? Yes//**<RET>** (Yes)

(R)

Edit Option Schedule

Option Name: LREPI NIGHTLY TASK

Menu Text: Lab Search/Extract Nightly Task TASK ID: 2298854

---

QUEUED TO RUN AT WHAT TIME: **APR 29,1998@01:45<RET>**

DEVICE FOR QUEUED JOB OUTPUT:**<RET>**

QUEUED TO RUN ON VOLUME SET: **ROU<RET>**

RESCHEDULING FREQUENCY: **1D<RET>**

TASK PARAMETERS:**<RET>**

SPECIAL QUEUEING:**<RET>**

## Trouble Shooting

1. If a NCH HL7 mail message is sent to the AITC but a NCH Verification Report mail message was **not** sent to the users please check the following:
  - Check the LR NCH and LR NCH-Report mail group memberships.
  - Check the member's name in the MAIL GROUP file (#3.8), Remote field (#12) in the LR NCH mail group for incorrect member names.
  - Ensure that the LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4), Report Mail Group field (#1) entry contains LR NCH-Report.
2. If no VA Alert mail messages are sent, please check for the following:
  - Ensure that the LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4), Send Alert field (#3) is set to "YES".
3. If the Lab Search/Extract Nightly Task [LREPI NIGHTLY TASK] option job does not generate a NCH HL7 mail message and Verification Report mail message but the new Lab Search/Extract Manual Run (Enhanced) [LREPI (EPI) ENHANCED MANUAL RUN] option does, please check the following:
  - Ensure that the LABORATORY SEARCH/EXTRACT file (#69.5), Active field (#1) entry for the NCH Cholesterol and NCH Pap smear Search/Extract is set to "NO".



## Health Level Seven (HL7) Protocol

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software application uses the **VISTA** Laboratory, PIMS, and Social Work databases for the EPI and NCH search/extract criteria. The **VISTA** HL7 software is used to transmit the data to the AITC database.

### 3. General Specifications

#### 3.1 Communication Protocol

The electronic **VISTA** MailMan software application is used as the communications protocol for sending the EPI and NCH HL7 mail messages between **VISTA** database and AITC database.

#### 3.2 Application Processing Rules

The HL7 protocol itself describes the basic rules for application processing by the sending and receiving systems. The HL7 Version 2.2 protocol is used. The ORU message is sent using the HL7 batch protocol.

#### 3.3 Message

The following HL7 mail message is used to support the exchange of data.

**ORU----- OBSERVATIONAL RESULTS UNSOLICITED**

#### 3.4 Segments

The following HL7 segments are used to support the exchange of data.

DG1	Diagnosis	OBX	Observation Results
<b>MSH----- MESSAGE HEADER</b>			
<b>----- PID</b>			
<b>-----PATIENT IDENTIFICATION</b>			
NTE	Notes and Comments	PV1	Patient Visit
OBR	Observation Request		

#### 3.5 Fields

The HL7 fields on the following page are used to support the exchange of data for each of the segments listed in the 3.4 Segments.

SEGMENT	FIELD SEQUENCE NUMBER	FIELD ELEMENT NAME	USER/HL7 DEFINED
<b>DG1</b>	1	Set ID-Diagnosis (Sequence #)	HL7
	3	Diagnosis Code (Code(id) ~Text (St.) ~ Name of coding system (st)	HL7
<b>MSH</b>	1	Field Separator	HL7
	2	Encoding Characters	HL7
	3	Sending Application	HL7
	4	Sending Facility	HL7
	5	Receiving Application	HL7
	6	Receiving Facility	HL7
	7	Date/Time of Message	HL7
	8	Security	HL7
	9	Message Type	HL7
	10	Message Control ID	HL7
	11	Processing ID	HL7
	12	Version ID	HL7
<b>OBR</b>	1	Set ID-Observation Request (Seq #)	HL7
	4	Universal Service ID (identifier ^ text ^ name of coding system ^ alt id ^ alt text ^ alt coding system)	HL7
	7	Observation Date/Time	HL7
	15	Specimen Source (Specimen source code (CE) ^^ text (TX) )	HL7 (Table 0070)
	26	Parent Results (OBX observation id of parent ^OBX sub ID	HL7
<b>NTE</b>	1	Set ID Notes and Comments (Seq #)	HL7
	3	Comment	HL7
<b>OBX</b>	1	Set Id-Observational Simple (seq. #)	HL7
	2	Value Type	HL7
	3	Observation Identifier (identifier ^ text ^ name of coding system ^ alt id ^ alt text ^ alt coding system)	HL7
	4	Sub Id	HL7
	5	Observation Value (Result)	HL7
	6	Units (Units)	HL7
	8	Abnormal Flags	HL7 (Table 0078)
	15	Date/Time of the Observation (Verified Date/Time)	HL7
<b>PID</b>	2	Patient ID (External ID)	HL7
	3	Patient ID (Internal ID)	HL7
	5	Patient Name	HL7
	7	Date of Birth	HL7
	8	Sex	HL7 (Table 0001)
	10	Race	HL7 (Table VA07)
	11	Address (Homeless)	HL7
	19	SSN	HL7
	27	Veteran's Military Status Period Of Service	HL7 (Table Va011)
<b>PV1</b>	1	Set ID - Patient Visit	HL7
	2	Patient Class	HL7
	36	Discharge Disposition	HL7
	44	Admit Date/Time (Event Date/Time)	HL7
	45	Discharge Date/Time	HL7

## 4.0 Transaction Specifications

### 4.1 General

The **VISTA** software sends an ORU observation result type HL7 message whenever one or more of the defined initiatives are identified.

### 4.2 Specific Transaction

#### A. Identified Encounter

When EPI and NCHP-CPSdb data are identified an ORU message is sent to the AITC. These ORU messages consist of the following segments:

#### Example: EPI ORU Message

```

ORU  OBSERVATIONAL RESULT UNSOLICITED
MSH  Message Header
NTE  Notes and Comments
PID  Patient Identification
PV1  Patient Visit
NTE  Notes and Comments
DG1  Diagnosis
OBR  Observation Report
OBX  Results

```

```

MSH |~|&|EPI-XXX|170|EPI-XXX|170|19961018113521|||ORU~R01|107|P|2.2|||USA
NTE |REPORTING DATE FROM 19850101 TO 19961018
PID |1|000-16-7946~0~M10|5~5~M10|EPI~PATONE||19220912|M||7|||||||000167946
PV1 |1|O|||||||||||||||||||||||||||||19950315151907
NTE |1|Vanc-Res Enterococcus
DG1 |1|451.19^DEEP PHLEBITIS-LEG NEC^I9
DG1 |2|511.9~PLEURAL EFFUSION NOS~I9
DG1 |3|670.02~MAJOR PUERP INF-DEL P/P~I9
DG1 |4|331.0~ALZHEIMER'S DISEASE~I9
DG1 |5|500.~COAL WORKERS' PNEUMOCON~I9
OBR |1||^CHEMISTRY TEST^VANLT||19950315151907|||||SER^^SERUM
OBX |1|ST|84330.0000^Glucose Quant^VANLT^175^GLUCOSE1^VA60||25|mg/dL|70-125|L*
NTE |2|2^Hepatitis C antibody
OBR |2||^CHEMISTRY TEST^VANLT||19950315151907|||||SER^^SERUM
OBX |1|ST|84330.0000^Glucose Quant^VANLT^175^GLUCOSE1^VA60||25|mg/dL|70-125|L*
PID |2|000-45-6666~8~M10|7~7~M10|EPI~PATTWO||19591229|F||7||||||000456666
PV1 |1|O|||||||||||||||||||||||||||||19950315152721
NTE |1|1^Vanc-Res Enterococcus
OBR |1||87999.0000^MICRO CULTURE^VANLT||198612100835|||||^BLOOD
OBX |1|CE|87993.0000^BACTERIOLOGY CULTURE^VANLT|1|^ESCHERICHIA COLI
OBR |2|^ANTIBIOTIC MIC^VANLT||198612100835|||||^BLOOD|||||87993.0000^1
OBX |1|ST|81812.0000^Neomycin^VANLT^18^NEOMYCN^VA62.06||||R
OBX |2|ST|^35^BACTRCN^VA62.06||||R
OBX |3|ST|81852.0000^Penicillin^VANLT^23^PENICLN^VA62.06||||R
OBX |4|ST|81676.0000^Clindamycin^VANLT^33^CLINDAM^VA62.06||||S
OBX |5|ST|81307.0000^Gentamicin^VANLT^33^GENTMCN^VA62.06||||R
OBX |6|ST|81656.0000^Chloramphenicol^VANLT^10^CHLORAM^VA62.06||||R
OBX |7|ST|81946.0000^Tetracycline NOS^VANLT^27^TETRCLN^VA62.06||||R

```

**Example: NCH ORU Mail Message**

```

MSH^~|\&^NCH-LAB^525^NCH-
LAB^525^19980416021403^^ORU~R01^1741793^P^2.2^^NE^AL^USA|
NTE^^~REPORTING DATE FROM 19980415 TO 19980415~2
PID^50^666-16-1010~4~M10^69836~2~M10^^TESTONE,
PATIENT^^19230609^M^^6^~02132^^^^ ^^^666161010^^^^^^2
PV1^1^O^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^19980414075625
NTE^1^50~NCH CHOLESTEROL
OBR^1^^^81121.0000~CHEMISTRY
TEST~VANLT^^^19980414075625^^^^^^SER~~SERUM^^^R/
CH 0414 333
OBX^1^ST^82466.0000~Cholesterol Total~VANLT~183~CHOLESTEROL~VA60^^227^mg/dl^"
S
EE TEST DESC"-^^^^^^19980414100054
OBX^2^ST^83013.0000~Cholesterol HDL~VANLT~244~HDL
CHOLESTEROL~VA60^^pending^mg/
dl^30-65^^^^^^19980414100054
OBX^3^ST^82350.0000~Calculation~VANLT~901~LDL
CHOLESTEROL~VA60^^pending^mg/dl^"
SEE TEST DESC"-^^^^^^19980414100054
OBX^4^ST^84480.0000~Triglycerides w o
extract~VANLT~205~TRIGLYCERIDE~VA60^^230^
mg/dl^35-200^H^^^^^^19980414100054
PID^51^000-54-0927~4~M10^73453~2~M10^^EPI~PATTHREE~E.^19620728^F^^7^~02332
^^^^^^000540927^^^^^^X
PV1^1^O^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^19980414101934
NTE^1^50~NCH CHOLESTEROL
OBR^1^^^81121.0000~CHEMISTRY
TEST~VANLT^^^19980414101934^^^^^^SER~~SERUM^^^R/
CH 0414 439
OBX^1^ST^82466.0000~Cholesterol Total~VANLT~183~CHOLESTEROL~VA60^^169^mg/dl^"
SEE TEST DESC -^^^^^^19980414115804
OBX^2^ST^83013.0000~Cholesterol HDL~VANLT~244~HDL
CHOLESTEROL~VA60^^pending^mg/
dl^30-65^^^^^^19980414115804
OBX^3^ST^82350.0000~Calculation~VANLT~901~LDL
CHOLESTEROL~VA60^^pending^mg/dl^"
SEE TEST DESC"-^^^^^^19980414115804
OBX^4^ST^84480.0000~Triglycerides w o
extract~VANLT~205~TRIGLYCERIDE~VA60^^72^m
g/dl^35-200^^^^^^19980414115804
PID^52^000-30-6640~4~M10^74913~4~M10^^EPI~PATSEVEN.^19400125^M^^7^~02021^^^^
^^000306640^^^^^^5
PV1^1^O^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^19980414112550
NTE^1^50~NCH CHOLESTEROL
OBR^1^^^81121.0000~CHEMISTRY
TEST~VANLT^^^19980414112550^^^^^^SER~~SERUM^^^R/
CH 0414 530
OBX^1^ST^82466.0000~Cholesterol Total~VANLT~183~CHOLESTEROL~VA60^^257^mg/dl^"
S
EE TEST DESC"-^^^^^^19980414141943
PID^53^000-42-9046~8~M10^75677~2~M10^^EPI JR.~PATFOUR.^19490720^M^^6^~145
85^^^^^^000429046^^^^^^7
PV1^1^I^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^^19980413164159
NTE^1^50~NCH CHOLESTEROL
OBR^1^^^81121.0000~CHEMISTRY
TEST~VANLT^^^19980414060010^^^^^^SER~~SERUM^^^R/
CH 0414 306
OBX^1^ST^82466.0000~Cholesterol Total~VANLT~183~CHOLESTEROL~VA60^^169^mg/dl^"
SEE TEST DESC-^^^^^^19980414094610
OBX^2^ST^83013.0000~Cholesterol HDL~VANLT~244~HDL
CHOLESTEROL~VA60^^pending^mg/
dl^30-65^^^^^^19980414094610

```



```
OBX^3^ST^82350.0000~Calculation~VANLT~901~LDL
CHOLESTEROL~VA60^^pending^mg/dl^"
SEE TEST DESC"-^^^^^^19980414094610
OBX^4^ST^84480.0000~Triglycerides w o
extract~VANLT~205~TRIGLYCERIDE~VA60^^184^
mg/dl^35-200^^^^^^19980414094610
```

**Table VA011 - Period of Service**

Value	Description
0	KOREAN
1	WORLD WAR I
2	WORLD WAR II
3	SPANISH AMERICAN
4	PRE-KOREAN
5	POST-KOREAN
6	OPERATION DESERT SHIELD
7	VIETNAM ERA
8	POST-VIETNAM
9	OTHER OR NONE
A	ARMY--ACTIVE DUTY
B	NAVY, MARINE--ACTIVE DUTY
C	AIR FORCE--ACTIVE DUTY
D	COAST GUARD--ACTIVE DUTY
E	RETIRED, UNIFORMED FORCES
F	MEDICAL REMEDIAL ENLIST
G	MERCHANT SEAMEN--USPHS
H	OTHER USPHS BENEFICIARIES
I	OBSERVATION/EXAMINATION
J	OFFICE OF WORKERS COMP.
K	JOB CORPS/PEACE CORPS
L	RAILROAD RETIREMENT
M	BENEFICIARIES-FOREIGN GOV
N	HUMANITARIAN (NON-VET)
O	CHAMPUS RESTORE
P	OTHER REIMBURS. (NON-VET)
Q	OTHER FEDERAL - DEPENDENT
R	DONORS (NON-VET)
S	SPECIAL STUDIES (NON-VET)
T	OTHER NON-VETERANS
U	CHAMPVA--SPOUSE, CHILD
V	CHAMPUS
W	CZECHOSLOVAKIA/POLAND SVC
X	PERSIAN GULF WAR
Y	CAV/NPS
Z	MERCHANT MARINE

**Table 0070 - Specimen Source Codes**

Abbreviations	Descriptions	Abbreviations	Descriptions	Abbreviations	Descriptions
ABS	Abscess	FLU	Body fluid, unsp	SER	Serum
AMN	Amniotic fluid	GAS	Gas	SKN	Skin
ASP	Aspirate	GAST	Gastric fluid/contents	SKM	Skeletal muscle
BPH	Basophils	GEN	Genital	SPRM	Spermatozoa
BIFL	Bile fluid	GENC	Genital cervix	SPT	Sputum
BBL	Blood bag	GENV	Genital vaginal	SPTT	Sputum tracheal aspirate
BLDC	Blood capillary	HAR	Hair	STON	Stone (use CALC)
BPV	Blood product unit	IHG	Inhaled Gas	STL	Stool = Fecal
BLDV	Blood venous	IT	Intubation tube	SWT	Sweat
BON	Bone	ISLT	Isolate	SNV	Synovial fluid (Joint fluid)
BRTH	Breath (use EXHLD)	LAM	Lamella	TEAR	Tears
BRO	Bronchial	WBC	Leukocytes	THRT	Throat
BRN	Burn	LN	Line	THRB	Thrombocyte (platelet)
CALC	Calculus (=Stone)	LNA	Line arterial	TISS	Tissue
CDM	Cardiac muscle	LVN	Line venous	TISG	Tissue gall bladder
CNL	Cannula	LIQ	Liquid NOS	TLGI	Tissue large intestine
CTP	Catheter tip	LYM	Lymphocytes	TLNG	Tissue lung
CSF	Cerebral spinal fluid	MAC	Macrophages	TISPL	Tissue placenta
CVM	Cervical mucus	MAR	Marrow	TSMI	Tissue small intestine
CVX	Cervix	MEC	Meconium	TISU	Tissue ulcer
COL	Colostrum	MBLD	Menstrual blood	TUB	Tube NOS
CBLD	Cord blood	MLK	Milk	ULC	Ulcer
CNJT	Conjunctiva	MILK	Breast milk	UMB	Umbilical blood
CUR	Curettage	NAIL	Nail	UMED	Unknown medicine
CYST	Cyst	NOS	Nose (nasal passage)	URTH	Urethra
DIAF	Dialysis fluid	ORH	Other	UR	Urine
DOSE	Dose med or substance	PAFL	Pancreatic fluid	URC	Urine clean catch
DRN	Drain	PAT	Patient	URT	Urine catheter
DUFL	Duodenal fluid	PRT	Peritoneal fluid ascites	URNS	Urine sediment
EAR	Ear	PLC	Placenta	USUB	Unknown substance
EARW	Ear wax (cerumen)	PLAS	Plasma	VOM	Vomit
ELT	Electrode	PLB	Plasma bag	BLD	Whole blood
ENDC	Endocardium	PLR	Pleural fluid (thoracentesis fld)	BDY	Whole body
ENDM	Endometrium	PMN	Polymorphonuclear neutrophils	WAT	Water
EOS	Eosinophils	PPP	Platelet poor plasma	WICK	Wick
RBC	Erythrocytes	PRP	Platelet rich plasma	WND	Wound
EYE	Eye	PUS	Pus	WNDA	Wound abscess
EXHLD	Exhaled gas (breath)	RT	Route of medicine	WNDE	Wound exudate
FIB	Fibroblasts	SAL	Saliva	WNDD	Wound drainage
FLT	Filter	SEM	Seminal fluid	XXX	To be specified in another part of the message
FIST	Fistula				

**Table VA07 - Race**

<b>Value</b>	<b>Description</b>
1	HISPANIC, WHITE
2	HISPANIC, BLACK
3	AMERICAN INDIAN OR ALASKA NATIVE
4	BLACK, NOT OF HISPANIC ORIGIN
5	ASIAN OR PACIFIC ISLANDER
6	WHITE NOT OF HISPANIC ORIGIN
7	UNKNOWN

**Table 0001 - Sex**

<b>Value</b>	<b>Description</b>
F	FEMALE
M	MALE
O	OTHER

**Table 0078 - Abnormal flags**

<b>Value</b>	<b>Description</b>
L	Below low normal
H	Above high normal
LL	Below lower panic limits
HH	Above upper panic limits
<b>For microbiology sensitivities only</b>	
S	Sensitive
R	Resistant
I	Intermediate
MS	Moderately sensitive
VS	Very sensitive



# APPENDIX - A

## EPI USER GUIDE



# EPI User Guide

The EPI User Guide provides all the necessary information, instructions, illustrations, and examples required for the EPI coordinators, Laboratory personnel, and other users to implement and maintain the Laboratory Search/Extract software. This information **must** be adhered to as recommended to assure a successful implementation and utilization of the software.

## Purpose

The VAHQ Infectious Disease Program Office Emerging Pathogens Initiative is to identify new antibiotic-resistant, otherwise problematic pathogens. Using this objective information, plans may be formulated on a national level for intervention strategies and resource needs. Results of aggregate data may also be shared with appropriate public health authorities for planning on the national level for the non-VA and private health care sectors. The **VISTA** Laboratory Search/Extract software tracks Emerging Pathogens for DVA and the national level without the necessity for additional site data entry.

## EPI Objective:

Emerging Pathogens (as defined by VAHQ) act as triggers for data acquisition for the **VISTA** Laboratory Search/Extract software. The software then retrieves relevant, predetermined, and patient-specific data for transmission to the AITC database repository. Once at that location, the data are analyzed using a SAS based statistical software. VAHQ Reports may then be generated for appropriate use and distribution at the national level.

- Identify Emerging Pathogens
- Extract specific data associated with the Emerging Pathogen
- Transmit data to AITC
- Create national Statistical Analysis System (SAS) data sets for Infectious Diseases Program Office access
- Intermittent periodic review **must** be done as determined by the sites. The Verification Report is **not** used to determine Laboratory or ICD collected data totals or validation (i.e., X numbers of cases of positive tests for Hepatitis C or X positive culture results for Streptococcus, Group A).

## EPI Data Transmission

### EPI HL7 formatted Mail Messages

The VISTA Laboratory Search/Extract software automatically processes the EPI data on the 15<sup>th</sup> of each month producing HL7 formatted mail messages. A Verification Report mail message is then produced (i.e., in a human readable format) for each HL7 mail message. The EPI HL7 formatted mail messages **must** be **manually** forwarded on the 15<sup>th</sup> of each month to the AITC database via the Q-EPI.MED.VA.GOV domain. Use the VA MailMan software to manually forward the EPI data to XXX@ Q-EPI.MED.VA.GOV domain.

### EPI Verification Report Mail Messages

The EPI Verification Report mail messages are sent to the EPI-REPORT mail group on the 15<sup>th</sup> of each month. The members of the EPI-REPORT mail group should review the EPI Verification Report and make data corrections (e.g., social security number, date of birth, period of service, etc.) as deemed necessary. The software automatically processes the EPI data corrections on the 15<sup>th</sup> of each month. **Note:** The EPI Verification Report may be generated as often as needed using the new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option.

### EPI Confirmation Mail Messages

EPI Confirmation mail messages are sent to the sending facility EPI-REPORT mail group **after** the EPI HL7 mail messages has been received by the AITC database.

### EPI Processing Report Mail Message

EPI Processing Report mail messages are sent to the sending facility EPI-REPORT mail group at the end of the AITC processing cycle (i.e., the 25<sup>th</sup> of each month). The EPI Processing Report mail messages confirms that the EPI data has been processed and lists any errors and/or warning codes requiring corrections.



## Austin Automation Center Database Processing

The Austin Information Technology Center (AITC) created two file structures, both in Statistical Analysis System (SAS) file format. These two file structures are used primarily as a source of data for the Infectious Diseases Program Office. The data is available to the Infectious Diseases Program Office to be manipulated and used for analysis and reporting.

The two file structures are referred to as the “Numerators” and “Denominators” because of their planned utilization.

### Numerator file:

The Numerator file contains an accumulation of data sent by all DVA medical centers. It contains twelve individual months worth of data and is updated monthly. Each month the oldest month data is removed from the Numerator file and the latest month’s data is added.

Upon receipt of the DVA medical centers monthly input, AITC transmits an acknowledgment message back to the sending DVA medical centers. The acknowledgment message identifies any “problem transmissions”. The “problem transmissions” records any discrepancies in the field format or the actual field value. AITC identifies the data as invalid records or just warning the DVA medical center that the record has some discrepancies. However, the data is **not** being rejected. Both the “problem transmissions” and the “accepted records” are documented on a Processing Report that is transmitted from the AITC to the DVA medical centers. The “Processing Report” itemizes all transmissions received by AITC, document the records status as either being accepted or rejected (with the reason code identified), and a warning stating that something is unusual about the value of one or more fields (with the warning reason code identified). An example of the “Tables of Reject and Warning codes” are located in the Appendix - B section of this guide.

The Numerator file information is specific to unique patients with a VAHQ designated Emerging Pathogen. The Emerging Pathogens are flagged through the **VISTA** software process. The Numerator file data are collected and transmitted to AITC monthly by the DVA facilities medical centers.

## Denominator file:

The Denominator file provides the Infectious Diseases Program Office, data elements for each DVA medical center. The source of these data elements is the corporate medical database residing in AITC. The individual files that these data elements are extracted from are the National Patient Care (NPC), Inpatient Treatment File (PTF), Automated Management Information System (AMIS), and Cost Distribution Report (CDR) systems.

The data elements are:

- Unique SSN served (inpatient and outpatient together)
- Total # of discharges
- Total unique SSN discharges
- Inpatient hospital days
- Inpatient ICU days
- Unique SSN encounters for both inpatient and outpatient

A “running 12 month” accumulation is required (i.e., there will always be one year’s worth of monthly counts) with the oldest month dropped off each cycle and a new one added.

**NOTE:** The need to track individual station data and to consolidate by parent station has not been specified. At this time we are only gathering by individual station number.

## Lab Search/Extract Protocol Edit [LREPI PROTOCOL EDIT] option

The new Lab Search/Extract Protocol Edit [LREPI PROTOCOL EDIT] option is located on the Lab Search/Extract Primary Menu [LREPI SEARCH EXTRACT MENU]. This new option is used for editing the EPI protocol.

### Example: EPI Protocol Parameters Setup Definition

Protocol Parameters Setup Definition

---

PROTOCOL: LREPI

Title: Emerging Pathogens Initiative (EPI)

Message Size: **32000**

Report Mail Group: **EPI-REPORT**

Send Alert: **YES**

Send Alert To

**DOE,Jane**

## EPI Descriptions and Input Examples

This section contains the 14 Emerging Pathogens descriptions and input examples that the **VISTA** Laboratory Search/Extract software is tracking. The following EPI descriptions **must** be reviewed for compliance (as specified by the VAHQ Infectious Disease Program Office) before defining the Lab Search/Extract Parameters Setup.

<i>Candida</i>	<i>Legionella</i>
<i>Clostridium difficile</i>	Leishmaniasis
Creutzfeldt-Jakob Disease	Malaria
<i>Cryptosporidium</i>	Pen- Res Pneumococcus
Dengue	<i>Streptococcus</i> -Group A
<i>E. coli</i> O157:H7	Tuberculosis
<i>Hepatitis C Antibody Pos</i>	<i>Vanc-Res Enterococcus</i>

**NOTE:** It is **highly recommended** that the following person(s) jointly participate in the review process for the 14 EPI descriptions:

Laboratory Information Manager (LIM)

Total Quality Improvement/Quality Improvement/Quality Assurance (TQI/QI/QA) staff (or person at the site with similar function)

Representative from the Microbiology section for the Emerging Pathogens Initiative (i.e., director, supervisor, or technologist)

The 14 Emerging Pathogens parameter descriptions will require an ongoing review process (as specified by the VAHQ Infectious Disease Program Office). The person(s) participating in the ongoing review process is responsible for ensuring the following requirements are kept current:

Periodic reviews of the ICD codes

Periodic reviews of the Lab Search/Extract Parameter Setup [LREPI PARAMETER SETUP] for the EPI

Annual review of the 14 Emerging Pathogens descriptions (as specified by the VAHQ Infectious Disease Program Office)

## Laboratory Search/Extract Parameters Input Screen

The following examples displays how to use the Laboratory Search/Extract Parameters Input Screen for defining the emerging pathogens parameters. Some of the Input Screen examples show partially pre-populated fields. The ETIOLOGY FIELD file (#61.2) site specific data is used to partially pre-populate the fields in the LAB SEARCH/EXTRACT file (#69.5). However, further data entries are required for site specific data. Additional data entries can be added or deleted to meet your site-specific needs.

**NOTE:** The Lab Search/Extract Primary [LREPI SEARCH EXTRACT MENU] Menu options are using VA FileMan screens displays, referred to as ScreenMan. For detailed instructions on how to use the screens please review the VA FileMan V. 21.0 User Manual, Section 6 ScreenMan.

**NOTE:** There are a number of different ways that sites have chosen to enter results into the **VISTA** database. As long as the results are in a retrievable format (straight from the **VISTA** database without additional manual input needed), how it is entered is **not** of significance to the Emerging Pathogen Initiative. However, two preferred methods make it easy to capture the data. Please reference the Helpful Hints section of this guide for the two methods.

**NOTE:** Please be consistent with site specific data spelling or alternate spelling to assure accurate data capture.

## Candida (Reference #8)

Fungal infections are rising in significance especially in severely ill patients. The same is true for bloodstream infections acquired in the hospital, especially those associated with intravenous lines. Fungal bloodstream infections are increasing in prevalence.

As a marker of bloodstream infections, the fungus *Candida* (and *Torulopsis*) has been chosen as an initial indicator organism. This organism may **not** be a prevalent or significant entity at your site, however, its presence is more likely to be indicative of serious or true infection than other organisms. The fungus *Candida* (and *Torulopsis*) may commonly be isolated from the blood in association with IV lines. Additionally, this yeast is more likely to be associated with nosocomial acquisition than other organisms (i.e., *Staphylococcus aureus* and coagulase negative *Staphylococcus*), which can cause a number of community acquired syndromes **not** at all related to IV lines.

All episodes of *Candida* (*Torulopsis*, yeast) isolation from blood or a blood source (central line, IV catheter tip, etc.) are being tracked. The **VISTA** Laboratory Search/Extract software has provided a partial pre-populated list of (etiologies/organisms) that fit the description for *Candida* (*Torulopsis*, yeast) to choose. These (etiologies/organisms) should be used, in addition to any site specific (etiologies/organisms) which may also fit the description.

**Example:**

```

                                Lab Search/Extract Primary Menu

    ENH    Lab Search/Extract Manual Run (Enhanced)
    VR     Print Detailed Verification Report
    LO     Local Pathogen Menu ...
    PI     Pathogen Inquiry
    UP     Lab EPI Parameter Setup
           Lab EPI Protocol Edit
    LK     Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
    CANDIDA
    CLOSTRIDIUM DIFFICILE
    CREUTZFELDT-JAKOB DISEASE
    CRYPTOSPORIDIUM
    DENGUE
    E. COLI 0157:H7
    HEPATITIS C ANTIBODY POS
    LEGIONELLA
    LEISHMANIASIS
    MALARIA
    NCH CHOLESTEROL
    NCH PAP SMEAR
    PEN-RES PNEUMOCOCCUS
    STREPTOCOCCUS GROUP A
    TUBERCULOSIS
    VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: Candida <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN				Page 1 of 5
NAME: Candida		ACTIVE: YES		
<hr/>				
Laboratory Test(s)	Indicator		Value	
<b>&lt;RET&gt;</b>				
ICD Coding System [ICD-9 or ICD-10]? (9/10): <b>&lt;RET&gt;</b>				
ICD Code	Cd Set	ICD Description		
<b>&lt;RET&gt;</b>				
<hr/>				
Exit	Save	Next Page	Refresh	
COMMAND: <b>N&lt;RET&gt;</b>		Press <PF1>H for help		Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: CANDIDA		ACTIVE: YES
<hr/>		
Selected Etiology	Selected Snomed Codes	
<b>Examples:</b> CANDIDA		
CANDIDA GUILLIERMONDII		

CANDIDA KRUSEI		
CANDIDA PARAPSILOSIS		
CANDIDA PSEUDOTROPICALIS		
CANDIDA SKIN TEST ANTIGEN		
CANDIDA STELLATOIDEA		
CANDIDA TROPICALIS		
CANDIDA, NOS		

**<RET>**

**Note:** During the post Init, the ETIOLOGY FIELD file (#61.2) was searched to pre-populate the Etiology field (#3) in the EMERGING PATHOGENS file (#69.5). Listed above are examples of etiology entries which may have been populated from your site's file. Additional etiologies may be added or deleted at the Selected Etiology prompt to meet your site specific needs.

**Note:** If spelling differences occur within your ETIOLOGY FIELD file (#61.2), be consistent with your local file and spell the results here, as it is spelled in your file (even if it is spelled differently in the example). We are concerned more importantly with data recovery.

Antimicrobial Susceptibility	NLT Code	NLT Description
------------------------------	----------	-----------------

**<RET>**

---

Exit	Save	Next Page	Refresh
------	------	-----------	---------

COMMAND: **N<RET>** Press <PF1>H for help    Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
---	--	-------------

NAME: Candida	ACTIVE:	YES
---------------	---------	-----

---

Topography Selection

Include	Exclude
Blood<RET>	<RET>
Bloodstream<RET>	
Catheter Tip<RET>	

**Note:** These are only suggestions. Please add accordingly to your site definition.

---

Exit	Save	Next Page	Refresh
------	------	-----------	---------

COMMAND: **N<RET>** Press <PF1>H for help    Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
---	--	-------------

NAME: Candida	ACTIVE:	YES
---------------	---------	-----

---

FIRST ENCOUNTER:<RET>	Follow PTF:YES<RET>
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>
Select SEX:<RET>	

---

Exit	Save	Refresh
------	------	---------



COMMAND: **E<RET>** Press <PF1>H for help Insert  
Save changes before leaving form (Y/N)?**Y<RET>**

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN Page 5 of 5  
NAME:Candida ACTIVE: YES

Run Date:<RET> Protocol:LREPI<RET>

Run Cycle:MONTHLY<RET> Lag Days:15<RET>

General Description:<TAB>

Exit Save Refresh

COMMAND: **E<RET>** Press <PF1>H for help Insert  
Save changes before leaving form (Y/N)?**Y<RET>**

## Clostridium difficile (Reference #4)

Disease associated with the presence of *Clostridium difficile* enterotoxin A can cause significant morbidity, as well as mortality. It is of importance, as its predominant acquisition seems to occur nosocomially. Presence of Clostridial toxin (either enterotoxin A or cytotoxin L) by assay (whether it be EIA, latex agglutination, cytotoxicity of cell culture  $\pm$  neutralization, or culture of organism with subsequent colony testing) is the best indicator that an inflammatory diarrheal disease is due to presence of *Clostridium difficile*.

Laboratory Services are quite varied as to how they identify the presence of *Clostridium difficile*. Some labs are set up to identify *C. difficile* as the final microbiological (bacterial) etiology of a culture, even if a culture method was not used. Other labs use a final etiology of “see comment” and then enter the results in a free text format. Still others enter the text under a hematology or chemistry format where a reference range and “positive” and “negative” result values can be entered. Wherever the facility lab places the results which are used to demonstrate the presence of toxin-producing *C. difficile*, we need to be able to track them (that means it **must** occur as a retrievable “positive” or “negative” result, or as a “bacterial etiology”). Results in a “Comments” or “Free-text” section are **not** acceptable.

There are a number of different ways that sites have chosen to enter *Clostridium difficile* toxin assay results into the **VISTA** database. As long as the toxin assay results are in a retrievable format (straight from the **VISTA** database without additional manual input needed), how it is entered is **not** of significance to the Emerging Pathogen Initiative. However, there are two preferred methods that make it easy to capture the data. Please reference the Appendix - B section of this guide for the two methods.

**Example:**

```

                                Lab Search/Extract Primary Menu

ENH    Lab Search/Extract Manual Run (Enhanced)
VR      Print Detailed Verification Report
LO      Local Pathogen Menu ...
PI      Pathogen Inquiry
UP      Lab EPI Parameter Setup
        Lab EPI Protocol Edit
LK      Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup
Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS
Select LAB SEARCH/EXTRACT NAME: CLOSTRIDIUM DIFFICILE <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: CLOSTRIDIUM DIFFICILE	ACTIVE:	YES
Laboratory Test(s)	Indicator	Value
<b>Clostridium&lt;RET&gt;</b> difficile toxin	<b>Contains&lt;RET&gt;</b>	<b>Pos&lt;RET&gt;</b>
<p><b>Note:</b> This example is only a suggestion. Please add accordingly to your site definition.</p> <p>ICD Coding System [ICD-9 or ICD-10]? (9/10): 9</p> <p>ICD Code Cd Set ICD Description</p> <p><b>&lt;RET&gt;</b></p>		
Exit	Save	Next Page
Refresh		
COMMAND: <b>N&lt;RET&gt;</b>		Press <PF1>H for help
		Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: CLOSTRIDIUM DIFFICILE	ACTIVE:	YES
Selected Etiology <b>Clostridium difficile toxin positive&lt;RET&gt;</b>		Selected Snomed Codes
<b>Note:</b> This is only a suggestion. Please add accordingly to your site definition.		
Antimicrobial Susceptibility <RET>	NLT Code	NLT Description
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		Press <PF1>H for help    Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: CLOSTRIDIUM DIFFICILE	ACTIVE:	YES
Topography Selection		
Include <RET>	Exclude <RET>	
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		Press <PF1>H for help    Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: CLOSTRIDIUM DIFFICILE	ACTIVE: YES	
First Encounter:<RET>	Follow PTF: <b>YES</b> <RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Selected SEX:<RET>		
Exit      Save      Refresh		
COMMAND: <b>E</b> <RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: CLOSTRIDIUM DIFFICILE	ACTIVE:	YES
Run Date:<RET>	Protocol: <b>LREPI</b> <RET>	
Run Cycle: <b>MONTHLY</b> <RET>	Lag Days: <b>15</b> <RET>	
General Description:<TAB>		
Exit      Save      Refresh		
COMMAND: <b>E</b> <RET>		Press <PF1>H for help    Insert
Save changes before leaving form (Y/N)? <b>Y</b> <RET>		

## Creutzfeldt-Jakob Disease (CJD) (Reference #13)

*Creutzfeldt-Jakob Disease (CJD)* disease is a rare illness associated with prions. The DVA has chosen to follow this entity because of historic problems with certain blood products used in the private and public health care sectors. The data will be one of a number of ways used to identify changes in trends of incidence of this illness. This task is remarkably complex because of the long incubation period of CJD. There are no specific tests for diagnosis other than central nervous system histology combined with clinical presentation. As such, this entity is followed through ICD coding.

### Example:

```

                                Lab Search/Extract Primary Menu

ENH    Lab Search/Extract Manual Run (Enhanced)
VR      Print Detailed Verification Report
LO      Local Pathogen Menu ...
PI      Pathogen Inquiry
UP      Lab EPI Parameter Setup
        Lab EPI Protocol Edit
LK      Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: CREUTZFELDT-JAKOB DISEASE <RET>

                                LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN                                Page 1 of 5
NAME:CREUTZFELDT-JAKOB DISEASE                                ACTIVE:    YES

Laboratory Test(s)                                Indicator                                Value
<RET>

```

ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>

ICD Code	Cd Set	ICD Description
046.1	ICD-9	JAKOB-CREUTZFELDT DIS

<RET>

Exit      Save      Next Page      Refresh

COMMAND: **N<RET>**      Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN      Page 2 of 5

NAME: CREUTZFELDT-JAKOB DISEASE      ACTIVE:      YES

Selected Etiology      Selected Snomed Codes<RET>

Antimicrobial Susceptibility	NLT Code	NLT Description
<RET>		

Exit      Save      Next Page      Refresh

COMMAND: **N<RET>**      Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN      Page 3 of 5

NAME: CREUTZFELDT-JAKOB DISEASE      ACTIVE:      YES

Topography Selection

Include	Exclude
<RET>	<RET>

Exit      Save      Next Page      Refresh

COMMAND: **N<RET>**      Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: CREUTZFELDT-JAKOB DISEASE	ACTIVE:	YES
<hr/>		
First Encounter: <RET>	Follow PTF: YES<RET>	
BEFORE DATE OF BIRTH: <RET>	AFTER DATE OF BIRTH: <RET>	
Select SEX: <RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>	Press <PF1>H for help	Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: CREUTZFELDT-JAKOB DISEASE	ACTIVE:	YES
<hr/>		
Run Date:<RET>	Protocol:LREPI<RET>	
Run Cycle:MONTHLY<RET>	Lag Days:15<RET>	
General Description:<TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>	Press <PF1>H for help	Insert
Save changes before leaving form (Y/N)?Y<RET>		



## Cryptosporidium (Reference #9)

The parasite *Cryptosporidium parvum* is a cause of water-borne diarrheal disease. It has gained recent prominence after evaluation of the outbreak in the greater Milwaukee area in 1993 which is estimated to have affected <400,000 persons. In addition to affecting HIV-infected persons and young children, information exists which demonstrates that the chronically ill, elderly are also a higher risk group than the general population. Microbiology laboratory data (parasitology for most laboratories) as well as ICD coding is used to track this disease, both are narrowly defined parameters.

**NOTE:** Microsporidiosis is a similar disease, however, the EPI do **not** currently wish to follow this disease process. Microsporidian etiologies should **not** be entered.

### **Example:**

```

                                Lab Search/Extract Primary Menu

ENH      Lab Search/Extract Manual Run (Enhanced)
VR        Print Detailed Verification Report
LO        Local Pathogen Menu ...
PI        Pathogen Inquiry
UP        Lab EPI Parameter Setup
          Lab EPI Protocol Edit
LK        Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup
Select LAB SEARCH/EXTRACT NAME: ?<RET>
  Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
  Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y  (Yes)<RET>
Choose from:
  CANDIDA
  CLOSTRIDIUM DIFFICILE
  CREUTZFELDT-JAKOB DISEASE
  CRYPTOSPORIDIUM
  DENGUE
  E. COLI 0157:H7
  HEPATITIS C ANTIBODY POS
  LEGIONELLA
  LEISHMANIASIS
  MALARIA
  NCH CHOLESTEROL
  NCH PAP SMEAR
  PEN-RES PNEUMOCOCCUS
  STREPTOCOCCUS GROUP A
  TUBERCULOSIS
  VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: CRYPTOSPORIDIUM <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 1 of 5						
NAME: CRYPTOSPORIDIUM		ACTIVE: YES							
<table style="width: 100%; border-collapse: collapse;"> <tr> <th style="text-align: left; width: 50%;">Laboratory Test(s)</th> <th style="text-align: left; width: 30%;">Indicator</th> <th style="text-align: left; width: 20%;">Value</th> </tr> <tr> <td colspan="3" style="padding-top: 5px;"><b>&lt;RET&gt;</b></td> </tr> </table>				Laboratory Test(s)	Indicator	Value	<b>&lt;RET&gt;</b>		
Laboratory Test(s)	Indicator	Value							
<b>&lt;RET&gt;</b>									
ICD Coding System [ICD-9 or ICD-10]? (9/10): <b>&lt;RET&gt;</b>									
ICD Code 007.8 <b>&lt;RET&gt;</b>	Cd Set ICD-9	ICD Description PROTOZOAL INTEST DIS N							
<table style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 20%;">Exit</td> <td style="width: 20%;">Save</td> <td style="width: 20%;">Next Page</td> <td style="width: 40%;">Refresh</td> </tr> </table>				Exit	Save	Next Page	Refresh		
Exit	Save	Next Page	Refresh						
COMMAND: <b>N&lt;RET&gt;</b>		Press <PF1>H for help    Insert							

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5						
NAME: CRYPTOSPORIDIUM		ACTIVE: YES							
<table style="width: 100%; border-collapse: collapse;"> <tr> <th style="text-align: left; width: 50%;">Selected Etiology</th> <th style="text-align: left; width: 50%;">Selected Snomed Codes</th> </tr> <tr> <td colspan="2" style="padding-top: 5px;"><b>Cryptosporidium&lt;RET&gt;</b></td> </tr> </table>				Selected Etiology	Selected Snomed Codes	<b>Cryptosporidium&lt;RET&gt;</b>			
Selected Etiology	Selected Snomed Codes								
<b>Cryptosporidium&lt;RET&gt;</b>									
<b>Note: If Cryptosporidium is reported under parasitology, add Cryptosporidium species at the Etiology prompt.</b>									
<table style="width: 100%; border-collapse: collapse;"> <tr> <th style="text-align: left; width: 40%;">Antimicrobial Susceptibility</th> <th style="text-align: left; width: 20%;">NLT Code</th> <th style="text-align: left; width: 40%;">NLT Description</th> </tr> <tr> <td colspan="3" style="padding-top: 5px;"><b>&lt;RET&gt;</b></td> </tr> </table>				Antimicrobial Susceptibility	NLT Code	NLT Description	<b>&lt;RET&gt;</b>		
Antimicrobial Susceptibility	NLT Code	NLT Description							
<b>&lt;RET&gt;</b>									
<table style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 20%;">Exit</td> <td style="width: 20%;">Save</td> <td style="width: 20%;">Next Page</td> <td style="width: 40%;">Refresh</td> </tr> </table>				Exit	Save	Next Page	Refresh		
Exit	Save	Next Page	Refresh						
COMMAND: <b>N&lt;RET&gt;</b>		Press <PF1>H for help    Insert							

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: CRYPTOSPORIDIUM	ACTIVE:	YES
Topography Selection		
Include <RET>	Exclude <RET>	
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: CRYPTOSPORIDIUM	ACTIVE:	YES
First Encounter:<RET>	Follow PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: CRYPTOSPORIDIUM	ACTIVE:	YES
Run Date:	Protocol: LREPI<RET>	
Run Cycle:MONTHLY<RET>	Lag Days:15<RET>	
General Description:<TAB>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help      Insert
Save changes before leaving form (Y/N)?Y<RET>		

## Dengue (Reference #12)

The mosquito-borne disease of Dengue Hemorrhagic Fever is a rare but re-emerging infection, especially in the Caribbean. The VA has seen cases of Dengue Hemorrhagic Fever over the last several years. Most of these cases have been in Dengue endemic areas served by the VA. However, as our society becomes more mobile, and the area of Dengue endemity expands, more cases are likely to occur. Because microbiologic culture is not routinely done and serology can be difficult to track, initially ICD coded diagnoses is used to track this entity.

### Example:

```

                                Lab Search/Extract Primary Menu

ENH   Lab Search/Extract Manual Run (Enhanced)
VR     Print Detailed Verification Report
LO     Local Pathogen Menu ...
PI     Pathogen Inquiry
UP     Lab EPI Parameter Setup
       Lab EPI Protocol Edit
LK     Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: DENGUE <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: DENGUE	ACTIVE: YES	

Laboratory Test(s)	Indicator	Value
<RET>		
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>		
ICD Code	Cd Set	ICD Description
061.	ICD-9	DENGUE
065.4	ICD-9	MOSQUITO-BORNE HEM FEVER
<RET>		
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		
Press <PF1>H for help		
Insert		

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: DENGUE		ACTIVE: YES	
Etiology		Selected Snomed Codes	
<RET>			
Antimicrobial Susceptibility	NLT Code	NLT Description	
<RET>			
Exit	Save	Next Page	Refresh
COMMAND: N<RET>			
Press <PF1>H for help			
Insert			

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: DENGUE		ACTIVE: YES
<hr/>		
Topography Selection		
Include <RET>	Exclude <RET>	
<hr/>		
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: DENGUE		ACTIVE: YES
<hr/>		
First Encounter: <RET>	Follow PTF: YES <RET>	
BEFORE DATE OF BIRTH: <RET>	AFTER DATE OF BIRTH: <RET>	
Select SEX: <RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: DENGUE		ACTIVE      YES
<hr/>		
Run Date: <RET>	Protocol: LREPI <RET>	
Run Cycle: MONTHLY <RET>	Lag Days: 15 <RET>	
General Description: <TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help      Insert
Save changes before leaving form (Y/N)? Y<RET>		

## E. coli O157:H7 (Reference #10)

*Escherichia coli* serotype O157 (*E. coli* O157) has gained prominence as a food-borne illness with potentially life threatening complications coming from the associated Hemolytic Uremic Syndrome. Not all sites routinely culture for the presence of *E. coli* O157 in stool specimens submitted for culture. In addition, *E. coli* O157 is not a microbiologic (bacterial) etiology pre-existing in the most recent - national microbiology lab package. In order to nationally track cultures positive for this organism, each site will need to make an etiology specific for *E. coli* O157 (e.g. *Escherichia coli* O157, *E. coli* O157, *E. coli* serotype O157, etc.). Some sites have already done this and will **not** need to generate a new entry.

**NOTE:** Entering *Escherichia coli* or *E. coli* from the bacterial etiology and then entering “serotype O157” or “O157”, under the “Comments” or “Free Text” section is **not** acceptable, as it will **not** allow the data to be retrieved nationally.

All subsequent positive cultures for this organism **must** then be entered under the new etiology.

Other serotypes of *E. coli* will also cause disease, but we will not currently track these as O157 causes by far, the majority of cases of interest for the national database.

The EPI criteria is dependent on your site. If your site already has an etiology that will select positive cultures for *E. coli* O157, then enter that etiology. However, if your site had to enter a new etiology to accommodate the EPI criteria, be sure to enter this new etiology here.

**Example:**

```

                                Lab Search/Extract Primary Menu

    ENH    Lab Search/Extract Manual Run (Enhanced)
    VR      Print Detailed Verification Report
    LO      Local Pathogen Menu ...
    PI      Pathogen Inquiry
    UP      Lab EPI Parameter Setup
           Lab EPI Protocol Edit
    LK      Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
  Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
  Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
  CANDIDA
  CLOSTRIDIUM DIFFICILE
  CREUTZFELDT-JAKOB DISEASE
  CRYPTOSPORIDIUM
  DENGUE
  E. COLI 0157:H7
  HEPATITIS C ANTIBODY POS
  LEGIONELLA
  LEISHMANIASIS
  MALARIA
  NCH CHOLESTEROL
  NCH PAP SMEAR
  PEN-RES PNEUMOCOCCUS
  STREPTOCOCCUS GROUP A
  TUBERCULOSIS
  VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: E. COLI 0157:H7 <RET>

```



LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 1 of 5
NAME: E. COLI 0157:H7		ACTIVE: YES	
Laboratory Test(s)	Indicator	Value	
<RET>			
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>			
ICD Code	Cd Set	ICD Description	
<RET>			
Exit	Save	Next Page	Refresh
COMMAND: N<RET>		Press <PF1>H for help	

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: E. COLI 0157:H7		ACTIVE: YES	
Selected Etiology			
<b>Example: Escherichia coli 0157</b> <RET>			
<b>Note:</b> Entering <i>Escherichia coli</i> or <i>E. coli</i> from the bacterial etiology and then entering "serotype 0157" or "0157", under the Comments section or in free text is <b>not</b> acceptable as it will <b>not</b> allow the data to be retrieved nationally).			
Antimicrobial Susceptibility	NLT Code	NLT Description	
<RET>			
Exit	Save	Next Page	Refresh
COMMAND: N<RET>		Press <PF1>H for help	

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: E. COLI 0157:H7		ACTIVE: YES
<hr/> <p style="text-align: center;">Topography Selection</p>		
Include <RET>	Exclude <RET>	
<hr/> <p>Exit      Save      Next Page      Refresh</p>		
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: E. COLI 0157:H7		ACTIVE: YES
<hr/>		
First Encounter:<RET>	Follow PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>		
<hr/> <p>Exit      Save      Refresh</p>		
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: E. COLI 0157:H7		ACTIVE      YES
<hr/>		
Run Date:<RET>	Protocol: LREPI<RET>	
Run Cycle: MONTHLY<RET>	Lag Days: 15<RET>	
General Description:<TAB>		
<hr/> <p>Exit      Save      Refresh</p>		
COMMAND: E<RET>		Press <PF1>H for help      Insert
Save changes before leaving form (Y/N)? Y<RET>		

## Hepatitis C Antibody Positive (Reference #2)

Hepatitis C is much more prevalent than originally thought at least in certain key patient sub-populations. As new and more sensitive assays come into use, we seem to find more evidence of this pathogen. We are looking for evidence of exposure to Hepatitis C in patients as demonstrated by Hepatitis C antibody positivity. The need for confirmatory testing or demonstration of active disease is not currently necessary in gathering data for this program. Different facilities may use different assays for this test. What we are looking for is evidence of presence of antibody to Hepatitis C, whether it be recorded as “weakly positive”, “strongly positive”, “positive”, or “present”. If other phrases are used to describe a test result, one should be able to differentiate the results upon entry into the program. As an example, the words, “present” “and “not present” would not allow retrieval of only positive cases as both phrases contain the word, “present”.

### Example:

```

                                Lab Search/Extract Primary Menu

ENH    Lab Search/Extract Manual Run (Enhanced)
VR      Print Detailed Verification Report
LO      Local Pathogen Menu ...
PI      Pathogen Inquiry
UP      Lab EPI Parameter Setup
        Lab EPI Protocol Edit
LK      Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS
Select LAB SEARCH/EXTRACT NAME: HEPATITIS C ANTIBODY POS <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 1 of 5
NAME: HEPATITIS C ANTIBODY POS		ACTIVE: YES	
Laboratory Test(s) <b>HEPATITIS C ANTIBODY&lt;RET&gt;</b>	Indicator <b>Contains&lt;RET&gt;</b>	Value <b>Pos&lt;RET&gt;</b>	
<b>Note:</b> Enter the appropriate test for your site, and how the results are reported.			
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>			
ICD Code <RET>	Cd Set	ICD Description	
Exit	Save	Next Page	Refresh
COMMAND: <b>N&lt;RET&gt;</b>		Press <PF1>H for help    Insert	

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: HEPATITIS C ANTIBODY POS		ACTIVE: YES	
Selected Etiology <RET>		Selected Snomed Codes	
Antimicrobial Susceptibility <RET>	NLT Code	NLT Description	
Exit	Save	Next Page	Refresh
COMMAND: <b>N&lt;RET&gt;</b>		Press <PF1>H for help    Insert	

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: HEPATITIS C ANTIBODY POS		ACTIVE: YES
<hr/> <p style="text-align: center;">Topography Selection</p>		
Include <RET>	Exclude <RET>	
<hr/>		
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: HEPATITIS C ANTIBODY POS		ACTIVE: YES
<hr/>		
First Encounter:<RET>	Follow PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: HEPATITIS C ANTIBODY POS		ACTIVE      YES
<hr/>		
Run Date:<RET>	Protocol: LREPI<RET>	
Run Cycle: MONTHLY<RET>	Lag Days: 15<RET>	
General Description:<TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help
Save changes before leaving form (Y/N)? Y<RET>		

### Legionella (Reference #7)

Since the American Legion Convention in Philadelphia in the 1970's, Legionnaires' Disease has been an illness of keen interest to the DVA. Because diagnosis is complex, we have chosen to review for presence of *Legionella* in culture and in ICD DIAGNOSIS file (#80). We will not look at *Legionella* direct fluorescent antibody positivity because of the potential high false positivity of this test. Likewise, serology is not easy to interpret or easily extracted from **VISTA** for our purposes and will **not** be included as a marker in this first iteration of the EPI program. Because it is not yet approved, the newer test of *Legionella* urinary antigen will not be used either. The Selected Etiology screen display has been partially pre-populated.

### Example:

```

                                Lab Search/Extract Primary Menu

ENH      Lab Search/Extract Manual Run (Enhanced)
VR        Print Detailed Verification Report
LO        Local Pathogen Menu ...
PI        Pathogen Inquiry
UP        Lab EPI Parameter Setup
          Lab EPI Protocol Edit
LK        Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: LEGIONELLA<RET>

```

```

                                LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN    Page 1 of 5

NAME: LEGIONELLA                                ACTIVE: YES

Laboratory Test(s)                                Indicator                                Value
<RET>

ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>

```

ICD Code 482.80 <RET>	Cd Set ICD-9	ICD Description LEGIONNARIE'S DISEASE
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: LEGIONELLA		ACTIVE: YES	
Selected Etiology			
<b>Examples:</b> LEGIONELLA BOZEMANII			
LEGIONELLA DUMOFFII			
LEGIONELLA GORMANII			
LEGIONELLA JORDANIS			
LEGIONELLA LONGBEACHAE			
LEGIONELLA MICDAEI			
LEGIONELLA OAKRIDGENSIS			
LEGIONELLA PNEUMOPHILIA			
LEGIONELLA SP			
LEGIONELLA WADSWORTHII			
<RET>			
<p><b>Note:</b> During the post Init, the ETIOLOGY FIELD file (#61.2) was searched to pre-populate the Etiology field (#3) in the EMERGING PATHOGENS file (#69.5). Listed above are examples of etiology entries which may have been populated from your site's file. Additional etiologies may be added or deleted at the <u>Selected Etiology</u> prompt to meet your site-specific needs.</p> <p><b>Note:</b> If spelling differences occur within your ETIOLOGY FIELD file (#61.2) be consistent with your local file and spell the results here, as it is spelled in your file (even if it is spelled differently in the example). We are concerned more importantly with data <u>recovery</u>.</p>			
Antimicrobial Susceptibility	NLT Code	NLT Description	
<RET>			
Exit	Save	Next Page	Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert	

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: LEGIONELLA	ACTIVE: YES	
<hr/> <p style="text-align: center;">Topography Selection</p>		
Include <RET>	Exclude <RET>	
<hr/>		
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: LEGIONELLA	ACTIVE: YES	
<hr/>		
First Encounter:<RET>	Follow PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: E. LEGIONELLA	ACTIVE YES	
<hr/>		
Run Date:<RET>	Protocol:LREPI<RET>	
Run Cycle:MONTHLY<RET>	Lag Days:15<RET>	
General Description:<TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help      Insert
Save changes before leaving form (Y/N)?Y<RET>		



## Leishmaniasis (Reference #14)

Leishmaniasis is a significant tropical disease that can cause serious complications. It is of interest to the Department of Veterans Affairs as Leishmania has caused illness among military personnel for many years. In addition, the Persian Gulf War occurred in an area of the world where the parasite is endemic. Because no simple, straightforward serology exists and no standard culture techniques exist, we have chosen to follow this entity through ICD diagnosis codes.

### **Example:**

```

                                Lab Search/Extract Primary Menu

ENH      Lab Search/Extract Manual Run (Enhanced)
VR        Print Detailed Verification Report
LO        Local Pathogen Menu ...
PI        Pathogen Inquiry
UP        Lab EPI Parameter Setup
          Lab EPI Protocol Edit
LK        Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: LEISHMANIASIS <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 1 of 5
NAME: LEISHMANIASIS		ACTIVE: YES	

---

Laboratory Test(s)	Indicator	Value
<b>&lt;RET&gt;</b>		

ICD Coding System [ICD-9 or ICD-10]? (9/10):9**<RET>**

ICD Code	Cd Set	ICD Description
085.0	<b>ICD-9</b>	VISCERAL LEISHMANIASIS
085.1	<b>ICD-9</b>	CUTAN LEISHMANIAS URBAN
085.2	<b>ICD-9</b>	CUTAN LEISHMANIAS ASIAN
085.3	<b>ICD-9</b>	CUTAN LEISHMANIAS ETHIOP
085.4	<b>ICD-9</b>	CUTAN LEISHMANIAS AMER
085.5	<b>ICD-9</b>	MUCOCUTAN LEISHMANIASIS
085.9	<b>ICD-9</b>	LEISHMANIASIS NOS
<b>&lt;RET&gt;</b>		

---

Exit	Save	Next Page	Refresh
------	------	-----------	---------

COMMAND: **N<RET>** Press <PF1>H for help    Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: LEISHMANIASIS		ACTIVE: YES	

---

Selected Etiology	Selected Snomed Codes
<b>&lt;RET&gt;</b>	

Antimicrobial Susceptibility	NLT Code	NLT Description
<b>&lt;RET&gt;</b>		

---

Exit	Save	Next Page	Refresh
------	------	-----------	---------

COMMAND: **N<RET>** Press <PF1>H for help    Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: LEISHMANIASIS		ACTIVE: YES
Topography Selection		
Include <RET>	Exclude <RET>	
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: LEISHMANIASIS		ACTIVE: YES
First Encounter:<RET>	FOLLOW PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: E. LEISHMANIASIS		ACTIVE YES
Run Date:<RET>	Protocol: LREPI<RET>	
Run Cycle: MONTHLY<RET>	Lag Days: 15<RET>	
General Description:<TAB>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help      Insert
Save changes before leaving form (Y/N)? Y<RET>		

## Malaria (Reference #11)

The plasmodial parasite is responsible for the blood-borne disease of malaria. Malaria can cause acute as well as chronic, relapsing disease. Occasionally, U.S. troops are deployed in malaria endemic areas. This placement could potentially put troops at risk for acquiring this disease. For the Emerging Pathogens Initiative program, we are interested in tracking patients with malaria, either acute or chronic, relapsing, and in either inpatient or outpatient status. No standardized serologic test allows for easy identification. Since not all sites consistently code and record malarial parasites seen histologically or on blood smears (not all of these interpretations are done through the Pathology and Laboratory Service), we have currently decided to track malaria based on ICD coding.

### **Example:**

Lab Search/Extract Primary Menu	
ENH	Lab Search/Extract Manual Run (Enhanced)
VR	Print Detailed Verification Report
LO	Local Pathogen Menu ...
PI	Pathogen Inquiry
UP	Lab EPI Parameter Setup
	Lab EPI Protocol Edit
LK	Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: **UP<RET>** Lab Search/Extract Parameter Setup

Select LAB SEARCH/EXTRACT NAME: **?<RET>**  
 Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER  
 Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? **Y** (Yes)**<RET>**  
 Choose from:

- CANDIDA
- CLOSTRIDIUM DIFFICILE
- CREUTZFELDT-JAKOB DISEASE
- CRYPTOSPORIDIUM
- DENGUE
- E. COLI 0157:H7
- HEPATITIS C ANTIBODY POS
- LEGIONELLA
- LEISHMANIASIS
- MALARIA
- NCH CHOLESTEROL
- NCH PAP SMEAR
- PEN-RES PNEUMOCOCCUS
- STREPTOCOCCUS GROUP A
- TUBERCULOSIS
- VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: **MALARIA <RET>**

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: MALARIA		ACTIVE: YES
<hr/>		
Laboratory Test(s)	Indicator	Value
<RET>		
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>		
ICD Code	Cd Set	ICD Description
084.0	ICD-9	FALCIPARUM MALARIA
084.1	ICD-9	VIVAX MALARIA
084.2	ICD-9	QUARTAN MALARIA
084.3	ICD-9	OVALE MALARIA
084.4	ICD-9	MALARIA NEC
084.5	ICD-9	MIXED MALARIA
084.6	ICD-9	MALARIA NOS
084.7	ICD-9	INDUCED MALARIA
084.8	ICD-9	BLACKWATER FEVER
084.9	ICD-9	MALARIA COMPLICATED NEC
<RET>		
<hr/>		
Exit	Save	Next Page Refresh
COMMAND: N<RET>		Press <PF1>H for help Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: MALARIA		ACTIVE: YES
<hr/>		
Selected Etiology	Selected Snomed Codes	
<RET>		
Antimicrobial Susceptibility	NLT Code	NLT Description
<RET>		
<hr/>		
Exit	Save	Next Page Refresh
COMMAND: N<RET>		Press <PF1>H for help Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: MALARIA		ACTIVE: YES
<hr/> <p style="text-align: center;">Topography Selection</p>		
Include <RET>	Exclude <RET>	
<hr/>		
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: MALARIA		ACTIVE: YES
<hr/>		
First Encounter:<RET>	FOLLOW PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Selected SEX:<RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: MALARIA		ACTIVE YES
<hr/>		
Run Date:<RET>	Protocol: LREPI<RET>	
Run Cycle: MONTHLY<RET>	Lag Days: 15<RET>	
General Description:<TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help
Save changes before leaving form (Y/N)? Y<RET>		

## Penicillin- Resistant Pneumococcus (Reference #3)

The emergence of antibiotic resistance in microbial agents is of great interest and concern for health care. Penicillin (PCN) was once the mainstay of therapy for *Streptococcus pneumoniae* infections but resistance to this agent is becoming more prominent. Different therapeutic strategies need to be developed once the prevalence of PCN-resistant *S. pneumoniae* reaches a critical threshold in a community. In order to monitor this, we are looking for the presence of any resistance in the pneumococci (either “moderate/intermediate” or “frank/high” level resistance). As such, any *S. pneumoniae* which is not fully susceptible to PCN on PCN susceptibility testing should be recorded.

### Example:

```

                                Lab Search/Extract Primary Menu

    ENH   Lab Search/Extract Manual Run (Enhanced)
    VR     Print Detailed Verification Report
    LO     Local Pathogen Menu ...
    PI     Pathogen Inquiry
    UP     Lab EPI Parameter Setup
          Lab EPI Protocol Edit
    LK     Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
    CANDIDA
    CLOSTRIDIUM DIFFICILE
    CREUTZFELDT-JAKOB DISEASE
    CRYPTOSPORIDIUM
    DENGUE
    E. COLI 0157:H7
    HEPATITIS C ANTIBODY POS
    LEGIONELLA
    LEISHMANIASIS
    MALARIA
    NCH CHOLESTEROL
    NCH PAP SMEAR
    PEN-RES PNEUMOCOCCUS
    STREPTOCOCCUS GROUP A
    TUBERCULOSIS
    VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: PEN-RES PNEUMOCOCCUS <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: PEN-RES PNEUMOCOCCUS		ACTIVE: YES
Laboratory Test(s)	Indicator	Value
<RET>		
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>		
ICD Code	Cd Set	ICD Description
<RET>		
Exit	Save	Next Page Refresh
COMMAND: <RET>		Press <PF1>H for help Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: PEN-RES PNEUMOCOCCUS		ACTIVE: YES
Selected Etiology	Selected Snomed Codes	
<b>NOTE:</b> You may enter a new ETIOLOGY, if you wish. <b>STREPTOCOCCUS PNEUMONIAE 12</b> Are you adding 'STREPTOCOCCUS PNEUMONIAE' as a new ETIOLOGY (the 1ST for this EMERGING PATHOGENS)?Y<RET>		
Antimicrobial Susceptibility	NLT Code	NLT Description
<b>Penicillin&lt;RET&gt;</b> Are you adding ' Penicillin ' as a new Antimicrobial Susceptibility (the 1ST for this EMERGING PATHOGENS)?Y <RET>		
Exit	Save	Next Page Refresh
COMMAND: <RET>		Press <PF1>H for help Insert



LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: PEN-RES PNEUMOCOCCUS		ACTIVE: YES
<hr/> <p style="text-align: center;">Topography Selection</p>		
Include <RET>	Exclude <RET>	
<hr/>		
Exit	Save	Next Page      Refresh
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: PEN-RES PNEUMOCOCCUS		ACTIVE: YES
<hr/>		
First Encounter:<RET>	FOLLOW PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Selected SEX:<RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: PEN-RES PNEUMOCOCCUS		ACTIVE YES
<hr/>		
Run Date:<RET>	Protocol: LREPI<RET>	
Run Cycle: MONTHLY<RET>	Lag Days: 15<RET>	
General Description:<TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>		Press <PF1>H for help      Insert
Save changes before leaving form (Y/N)? Y<RET>		

## Streptococcus-Group A (Reference #6)

*Streptococcus*-Group A can be associated with or cause significant disease such as severe fasciitis and streptococcal toxic shock syndrome. We are especially interested to find out how much severe/deep seated disease the VA is experiencing, but other disease entities are of interest also. To this end, we are looking for all episodes of culture positivity for *Streptococcus*-Group A, regardless of site and regardless of inpatient or outpatient status of the person from whom the specimen is obtained. We are aware that some sites may use rapid screenings for *Streptococcus*-Group A, especially from pharyngeal sources. These rapid screens may be difficult to capture, so we are not asking for them on this first iteration of the EPI program.

### **Example:**

```

                                Lab Search/Extract Primary Menu

ENH   Lab Search/Extract Manual Run (Enhanced)
VR     Print Detailed Verification Report
LO     Local Pathogen Menu ...
PI     Pathogen Inquiry
UP     Lab EPI Parameter Setup
       Lab EPI Protocol Edit
LK     Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: STREPTOCOCCUS-GROUP A <RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: STREPTOCOCCUS-GROUP A		ACTIVE: YES
Laboratory Test(s) <RET>	Indicator	Value
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>		
ICD Code <RET>	Cd Set	ICD Description
Exit      Save      Next Page      Refresh		
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: STREPTOCOCCUS-GROUP A		ACTIVE: YES
Selected Etiology <b>STREPTOCOCCUS-GROUP A</b> <RET>	Selected Snomed Codes	
Antimicrobial Susceptibility <RET>	NLT Code	NLT Description
Exit      Save      Next Page      Refresh		
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: STREPTOCOCCUS-GROUP A		ACTIVE: YES
Topography Selection		
Include <RET>	Exclude <RET>	
Exit      Save      Next Page      Refresh		
COMMAND: N<RET>		Press <PF1>H for help      Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: STREPTOCOCCUS-GROUP A		ACTIVE: YES
First Encounter:<RET>	FOLLOW PTF: <b>YES</b> <RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	

Select SEX:<RET>

Exit      Save      Refresh

COMMAND: **E<RET>**

Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN      Page 5 of 5

NAME: STREPTOCOCCUS-GROUP A

ACTIVE YES

Run Date:<RET>

Protocol:**LREPI<RET>**

Run Cycle:**MONTHLY<RET>**

Lag Days:**15<RET>**

General Description:<TAB>

Exit      Save      Refresh

COMMAND: **E<RET>**

Press <PF1>H for help

Insert

Save changes before leaving form (Y/N)?**Y<RET>**

Tuberculosis (Reference #5)

*Mycobacterium tuberculosis* infection is an important public health concern. Recent increases in incidence of disease, and occurrence of multiply-drug resistant strains in outbreak situations along with the increased susceptibility of HIV-infected persons for this disease has generated renewed interest in this entity. Since the national data show that 80-85% of all reported active tuberculosis cases are culture positive (with acid fast bacilli smear-only positive cases increasing the reporting by 2-5% more) we have decided to use culture positivity for *Mycobacterium tuberculosis* to track tuberculosis infections in the current iteration of the EPI software package. Information regarding susceptibility will be tracked as well. For the national EPI program, there will be no need to enter specific antimycobacterial agents to be tracked; it will be done automatically. ICD coding is complex and confusing for many cases of tuberculosis and therefore will **not** be used.

**Example:**

```

                                Lab Search/Extract Primary Menu

    ENH    Lab Search/Extract Manual Run (Enhanced)
    VR      Print Detailed Verification Report
    LO      Local Pathogen Menu ...
    PI      Pathogen Inquiry
    UP      Lab EPI Parameter Setup
           Lab EPI Protocol Edit
    LK      Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
  Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
  Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
  CANDIDA
  CLOSTRIDIUM DIFFICILE
  CREUTZFELDT-JAKOB DISEASE
  CRYPTOSPORIDIUM
  DENGUE
  E. COLI 0157:H7
  HEPATITIS C ANTIBODY POS
  LEGIONELLA
  LEISHMANIASIS
  MALARIA
  NCH CHOLESTEROL
  NCH PAP SMEAR
  PEN-RES PNEUMOCOCCUS
  STREPTOCOCCUS GROUP A
  TUBERCULOSIS
  VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: TUBERCULOSIS<RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: TUBERCULOSIS		ACTIVE: YES
Laboratory Test(s)	Indicator	Value
<RET>		
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>		
ICD Code	Cd Set	ICD Description
<RET>		
Exit	Save	Next Page Refresh
COMMAND: N<RET>		Press <PF1>H for help Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME:TUBERCULOSIS		ACTIVE:YES
Selected Etiology	Selected Snomed Codes	
Mycobacterium tuberculosis<RET>		
Antimicrobial Susceptibility	NLT Code	NLT Description
<RET>		
Exit	Save	Next Page Refresh
COMMAND: N<RET>		Press <PF1>H for help Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: TUBERCULOSIS		ACTIVE:YES
Topography Selection		
Include	Exclude	
<RET>		
Exit	Save	Next Page Refresh
COMMAND: N<RET>		Press <PF1>H for help Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: TUBERCULOSIS		ACTIVE: YES

First Encounter:<RET>	FOLLOW PTF:YES<RET>
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>
Select SEX:<RET>	
<hr/>	
Exit	Save      Refresh
COMMAND: E<RET>	Press <PF1>H for help

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: TUBERCULOSIS	ACTIVE	YES
<hr/>		
Run Date:<RET>	Protocol:LREPI<RET>	
Run Cycle:MONTHLY<RET>	Lag Days:15<RET>	
General Description:<TAB>		
<hr/>		
Exit	Save	Refresh
COMMAND: E<RET>	Press <PF1>H for help	Insert
Save changes before leaving form (Y/N)?Y<RET>		

## Vancomycin-Resistant Enterococcus (VRE) (Reference #1)

Vancomycin-Resistant Enterococcus (VRE) is a pathogen of increasing importance. Not only can it cause significant disease, but also it can be spread within facilities. It is important to capture all positive cultures for VRE (not just disease). As such, all positive cultures for VRE will be reported.

**Note:** This includes cultures positive for prevalence and surveillance review, including specimens of stool and rectal swabs.

Vancomycin-resistant *Enterococcus faecalis* and *E. faecium* are most common, but we wish to look at all vancomycin resistant enterococci whether speciated or not. Therefore, it is important to be sure to list all the places in the Micro Lab package where *Enterococcus* are found, either as *Enterococcus*, *E. (sp.)*, Group D-*Streptococcus*, *E. faecalis*, *E. faecium*, *E. durans*, *E. gallinarum*, *E. casseliflavus*, etc.

**NOTE:** Only a partial pre-populated Etiology list is shown in the screen display example at the Selected Etiology prompt. Please be sure to review the entire Etiology list. If you have other etiology results at your site, they can be added to this Etiology list. Again, if alternate spellings are present in your site's ETIOLOGY FIELD file (#61.2), be certain those spellings assure capture of all data points possible.



**Example:**

```

                                Lab Search/Extract Primary Menu

    ENH      Lab Search/Extract Manual Run (Enhanced)
    VR       Print Detailed Verification Report
    LO       Local Pathogen Menu ...
    PI       Pathogen Inquiry
    UP       Lab EPI Parameter Setup
            Lab EPI Protocol Edit
    LK       Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
    CANDIDA
    CLOSTRIDIUM DIFFICILE
    CREUTZFELDT-JAKOB DISEASE
    CRYPTOSPORIDIUM
    DENGUE
    E. COLI 0157:H7
    HEPATITIS C ANTIBODY POS
    LEGIONELLA
    LEISHMANIASIS
    MALARIA
    NCH CHOLESTEROL
    NCH PAP SMEAR
    PEN-RES PNEUMOCOCCUS
    STREPTOCOCCUS GROUP A
    TUBERCULOSIS
    VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: VANC-RES ENTEROCOCCUS <RET>

```

```

                                LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN    Page 1 of 5

NAME: VANC-RES ENTEROCOCCUS                                ACTIVE: YES

-----
Laboratory Test(s)                Indicator                Value
<RET>

ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>

ICD Code                Cd Set                ICD Description
<RET>

-----
Exit      Save      Next Page      Refresh

COMMAND: N<RET>                                Press <PF1>H for help      Insert

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5						
NAME: VANC-RES ENTEROCOCCUS		ACTIVE: YES						
<div style="display: flex; justify-content: space-between;"> <div style="width: 45%;"> <p>Selected Etiology</p> <p><b>Examples:</b> Enterococcus</p> <p style="margin-left: 20px;">Enterococcus (Strept. faecalis-Group D)</p> <p style="margin-left: 20px;">Streptococcus faecalis</p> <p style="margin-left: 20px;">Streptococcus faecium</p> <p style="margin-left: 20px;">Enterococcus avium</p> <p style="margin-left: 20px;">Enterococcus avium - (Group D)</p> <p style="margin-left: 20px;">Enterococcus casseliflavus</p> <p style="margin-left: 20px;">Enterococcus faecalis</p> <p style="margin-left: 20px;">Enterococcus gallinarum</p> <p style="margin-left: 20px;">Enterococcus malodoratus</p> <p style="margin-left: 20px;">Enterococcus hirae</p> <p style="margin-left: 20px;">Enterococcus mundtii</p> <p style="margin-left: 20px;">Enterococcus raffinosus</p> <p style="margin-left: 20px;">Enterococcus sp.</p> <p style="margin-left: 20px;">Enterococcus species</p> </div> <div style="width: 45%;"> <p>Selected Snomed Codes</p> <p style="margin-left: 20px;">Enterococcus durans</p> <p style="margin-left: 20px;">Streptococcus sp. Group D</p> <p style="margin-left: 20px;">Enterococcus solitarius</p> <p style="margin-left: 20px;">Enterococcus pseudoavium</p> <p style="margin-left: 20px;">Enterococcus faecium</p> <p style="margin-left: 20px;">Enterococcus durans</p> </div> </div> <p>&lt;RET&gt;</p> <p><b>Note:</b> During the post Init, the ETIOLOGY FIELD file (#61.2) was searched to pre-populate the Etiology field (#3) in the EMERGING PATHOGENS file (#69.5). Listed above are examples of etiology entries which may have been populated from your site's file. Additional etiologies may be added or deleted at the <u>Selected Etiology</u> prompt to meet your site specific needs.</p> <p><b>Note:</b> If spelling differences occur within your ETIOLOGY FIELD file (#61.2) be consistent with your local file and spell the results here, as it is spelled in your file (even if it is spelled differently in the example). We are concerned more importantly with data <u>recovery</u>.</p>								
<table style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="text-align: left; width: 40%;">Antimicrobial Susceptibility</th> <th style="text-align: left; width: 20%;">NLT Code</th> <th style="text-align: left; width: 40%;">NLT Description</th> </tr> </thead> <tbody> <tr> <td colspan="3">VANCOMYCIN&lt;RET&gt;</td> </tr> </tbody> </table>			Antimicrobial Susceptibility	NLT Code	NLT Description	VANCOMYCIN<RET>		
Antimicrobial Susceptibility	NLT Code	NLT Description						
VANCOMYCIN<RET>								
<div style="display: flex; justify-content: space-between;"> <span>Exit</span> <span>Save</span> <span>Next Page</span> <span>Refresh</span> </div>								
<div style="display: flex; justify-content: space-between;"> <span>COMMAND: N&lt;RET&gt;</span> <span>Press &lt;PF1&gt;H for help</span> <span>Insert</span> </div>								

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN				Page 3 of 5
NAME: VANC-RES ENTEROCOCCUS		ACTIVE: YES		
<hr/>				
Topography Selection				
Include <RET>			Exclude <RET>	
<hr/>				
Exit	Save	Next Page	Refresh	
COMMAND: N<RET>		Press <PF1>H for help      Insert		

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN				Page 4 of 5
NAME: VANC-RES ENTEROCOCCUS		ACTIVE: YES		
<hr/>				
First Encounter:<RET>			FOLLOW PTF: YES<RET>	
BEFORE DATE OF BIRTH:<RET>			AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>				
<hr/>				
Exit	Save	Refresh		
COMMAND: E<RET>		Press <PF1>H for help      Insert		
Save changes before leaving form (Y/N)?Y<RET>				

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN				Page 5 of 5
NAME: VANC-RES ENTEROCOCCUS		ACTIVE YES		
<hr/>				
Run Date:<RET>			Protocol: LREPI<RET>	
Run Cycle: MONTHLY<RET>			Lag Days: 15<RET>	
General Description:<TAB>				
<hr/>				
Exit	Save	Refresh		
COMMAND: E<RET>		Press <PF1>H for help      Insert		
Save changes before leaving form (Y/N)?Y<RET>				

## Conclusion

Once you have finished entering the information as directed by the national Infectious Diseases Program Office, these fields should **not** again be changed except for the following conditions:

1. As requested nationally via the Veterans Affairs Headquarters (VAHQ) Infectious Disease Program Office to update, modify, add, or delete data from the existing files used by the Laboratory Search/Extract software or an addition of a new entity to be tracked.
2. The yearly review must ensure that the entry is acceptable and to update the EPI files with any changes in etiology, lab tests or results parameters that may have occurred locally at the site during the previous year.

Annually the EPI national program materials should be reviewed by the VAMCs and updated. It is suggested that this review occur in February of each year. If no changes have occurred in lab practices, etiologies, sites, or results parameters leave the information as is until the next review period. If changes did occur, then enter them as appropriate in order to capture the data requested for each EPI national entity (disease/organism) to be tracked.

As entities (diseases/organisms) are no longer to be tracked nationally (“dropped from the list”), or a new entity is to be tracked (“added to the list”), revision will be forwarded to the sites to assist in updating your site files.

## EPI Helpful Hints:

### Clostridium difficile

There are two preferred methods that will make it easy to capture data for *Clostridium difficile* criteria (i.e., as well as several other methods which sites may already employ). **Note:** As long as the designated parameter results being tracked are in a retrievable field (i.e., **not** a “Free Text” or “Comment” field) the method the site chooses is an individual decision.

**Preferred Method:** The first preferred method is to have the site define an etiology of “**Clostridium difficile toxin positive**”. This allows a topography specimen of accession area “**feces/stool**” to be accessioned through the Microbiology accession area. Then, if the stool specimen were indeed positive for *Clostridium difficile* toxin, by any of the known methods of testing, the etiology would be “**Clostridium difficile toxin positive**.” To accomplish this method would require sites to enter three new local etiologies:

- **Clostridium difficile toxin positive**
- **Clostridium difficile toxin negative**
- **Clostridium difficile toxin indeterminant**

These would be different from a culture isolate being positive for *Clostridium difficile*, in that they actually are etiologies/results based on toxin testing. This leaves the etiology of *Clostridium difficile* for actual culture positive specimens for the organism *Clostridium difficile*. The Lab Search/Extract Parameter Setup, the site parameter by which the software will capture a patient diagnosed with proven *Clostridium difficile*-associated colitis, will be by placing “**Clostridium difficile toxin positive**” etiology into the selected etiology entry screen. This has the advantage of being more consistent with other data entry practices in the Microbiology sections of most laboratories.

**Preferred Method:** The second preferred method is having the data in retrievable form would be to enter/accession the specimen for *Clostridium difficile* toxin assay under the chemistry/serology format (regardless of where the test is physically done) with the results being a choice of “positive”, “negative”, or “indeterminate”. This would allow one to enter “*Clostridium difficile* toxin” assay as the test for the EPI software to search in the chemistry/serology format. The result would be retrievable for EPI under a chemistry/serology lab test of “*Clostridium difficile* toxin” with the indicator “contains” and the value of “pos”, as noted in the sample page. If your site does not routinely do *Clostridium difficile* toxin assay testing this way, a different method of accessioning the specimen to get it in chemistry/serology format would be needed.

However, the Chemistry/Serology format would give additional flexibility in placing interpretational guidelines for the test results in the “Comments” field. For the EPI, “positive” or “negative” results can **not** be located in a “Free Text” or “Comments” field as these are **not** retrievable.

Some VAMCs accession the stool specimen for the *Clostridium difficile* toxin assay under the Microbiology format. An etiology is not given under the final culture result, but written into free text or comments section stating the *Clostridium difficile* toxin assay test result. This is not in a retrievable format and therefore not acceptable for the EPI criteria.

Some VAMCs still use cytotoxin assays of cell culture, which are again entered in a “Free Text” or “Comment” field. This again is not acceptable unless it is accessioned and recorded under the chemistry/serology format as a straightforward lab test result of “positive” or “negative” or “indeterminate”.

Some VAMCs choose to report *Clostridium difficile* toxin assay positivity under the Microbiology application. As an etiology/culture result of *Clostridium difficile* (even though culture, was not actually done) this is not a true measure of what is actually being tested (as most sites do not culture the organism but just run the toxin assay test). However, if your site uses this means to represent *Clostridium difficile* toxin assay positivity and there are no exceptions (such as the site reporting an actual positive culture of (*Clostridium difficile* which is toxin assay negative), then this would be acceptable though less desirable for EPI purposes.

## Validating EPI Data Capture

Once the Lab Search/Extract Parameter is setup, the Laboratory Search/Extract software automatically generates an Verification Report message around the 15<sup>th</sup> of each month. The Emerging Pathogen Verification Report should be used to compare the EPI data capture to site-specific data capture. This comparison also helps determine if the Lab Search/Extract Parameter for the EPI criteria is accurately defined. It is recommended that the Lab Search/Extract Manual Run [LREPI (EPI) MANUAL RUN] option be run to evaluate 1-3 months of data (i.e., as determined by the sites).

The Microbiology Laboratory staff, Laboratory Manager, TQI/QI/QA, or other person (i.e., as determined by the sites) may already have data of isolated “organisms of interest”. Several of the nationally defined emerging pathogens may well corresponds with the existing data. Therefore, a quick comparison can be done using the Emerging Pathogens Verification Report message. This comparison also ensures that the Laboratory Search/Extract software is appropriately capturing the EPI cases and numbers

For tests such as Hepatitis C, most Laboratory Managers should be able to generate reports (with patient names) that includes “positive” tests results to use for comparison. Additionally, the Health Information Management Section at each site should be able to generate a report of ICD-9 and ICD-10 Diagnoses by date. This ICD Diagnoses by date Report helps determine if the 14 VAHQ defined Emerging Pathogens data captures concurs with the EPI criterion (i.e., Cryptosporidium-007.8, Legionnaire’s disease--482.80, malaria--084, 084.0, 084.1, 084.2, 084.3, 084.4, 084.5, 084.6, 085.7, 084.8, 084.9, dengue-061, 065.4, Creutzfeldt-Jakob--046.1, and Leishmaniasis--085, 085.0, 085.1, 085.2, 085.3, 085.4, 085.5, 085.9).

Be aware that a number of these emerging pathogens **do not** occur at a high frequency. Sites with previously known cases of emerging pathogens, such as TB, should run the Lab Search/Extract Manual Run [LREPI (EPI) MANUAL RUN] option for the entire month to verify that the TB culture was isolated and to see if it is captured. Additionally, “test patients” known to have these lab results can also be run.

The purpose of this validation is **not** to require extra paperwork for QI monitors and long term document files. The validation should be done at the initial implementation of the Laboratory Search/Extract software to ensure accurate data capture. Thereafter, a review should be done once every 4-6 months to ensure that Lab Search/Extract Parameter Setup [LREPI (EPI) PARAMETER UPDATE] option entries for the EPI criteria remains accurate. Parameter updates may be required if a new lab test/result is to be implemented for one of the Emerging Pathogens Initiative.

## Emerging Pathogens Verification Report Message

The following is an example of the Emerging Pathogens Verification Report mail message. The EPI software automatically transmits the EPI data to the AITC on the 15<sup>th</sup> of each month. The AITC processes the on the 25<sup>th</sup> of each month

### **Example:**

Emerging Pathogens Verification Report [#60004] Page 1

---

REPORTING DATE FROM 12-01-1996 TO 12-31-1996 Message Seq # 1 Auto

JOHN DOEDOE XXX-XX-XXXX 07-07-1913 M WORLD WAR II 45205  
Outpatient Accession Date 12-11-1996@1025

\*\*\*\*\* STREPTOCOCCUS GROUP A \*\*\*\*\*

12-11-1996@1025 BACT 96 10383 MICRO CULTURE LEG  
1 12-13-1996 STREPTOCOCCUS BETA HEMOLYTIC, GROUP A  
2 12-13-1996 STAPHYLOCOCCUS (COAGULASE NEGATIVE)

ORG # 1 12-11-1996@1025 ANTIBIOTIC MIC LEG

ORG # 2 12-11-1996@1025 ANTIBIOTIC MIC LEG

JOHN DOEDOE XXX-XX-XXXX 01-08-1923 M WORLD WAR II 45239  
Inpatient Admission Date 12-19-1996@1125

\*\*\*\*\*4 CLOSTRIDIUM DIFFICILE \*\*\*\*\*  
12-25-1996@1415 MSER 96 418 CHEMISTRY TEST  
FECES

Clostridium Difficile Toxin 12-27-1996@1403  
POSITIVE

Can be verified using standard result reviews for "CH" subscripted tests (e.g., LRRSP, LRRP3, LRSORD, LRSORA, LRGEN)
---

JOHN DOE ###-##-#### 11-05-1910 M WORLD WAR II 45255  
Inpatient Admission Date 12-03-1996@1908  
Discharge Date 12-09-1996@1151 Discharge Disposition REGULAR

250.01 DIABETES MELLI W/O COMP TYP I

276.8 HYPOPOTASSEMIA

427.31 ATRIAL FIBRILLATION

428.0 CONGESTIVE HEART FAILURE

482.30 PNEUM. UNSPEC. STREPTOCOCCUS

PTF data can be verified using several different PTF options: DG PTF ICD DIAGNOSIS SEARCH DG PTF SUMMARY DIAG/OP OUTPUT DG PTF COMPREHENSIVE INQUIRY (most require DGPTFSUP key)
---

\*\*\*\*\*6 STREPTOCOCCUS GROUP A \*\*\*\*\*

12-04-1996 BACT 96 10187 MICRO CULTURE SPUTUM  
1 12-06-1996 STREPTOCOCCUS BETA HEMOLYTIC, GROUP A  
2 12-06-1996 STAPHYLOCOCCUS AUREUS

ORG # 1 12-04-1996 ANTIBIOTIC MIC SPUTUM

ORG # 2 12-04-1996 ANTIBIOTIC MIC SPUTUM



Penicillin	R
Clindamycin	S
Vancomycin	S
TETRCLN	S
TRMSULF	S
Erythromycin	S
Oxacillin	S
Cephalothin	S
Ciprofloxacin	S
AMPICILLIN-SULBACTAM	S

Microbiology subscribed organisms and susceptibilities can be reviewed using LRMIPSZ, LRMIPC, LRMIPOG, LRGEN, LRRSP, and LRRP3.
---

JOHN DOE ###-##-#### 02-23-1920 M PRE-KOREAN 45150  
 Inpatient Admission Date 11-18-1996@2213

\*\*\*\*\*8 CANDIDA \*\*\*\*\*

12-11-1996@0100 BLD 96 3914 MICRO CULTURE BLOOD  
 1 12-15-1996 CANDIDA ALBICANS

ORG # 1 12-11-1996@0100 ANTIBIOTIC MIC BLOOD

JOHN DOE ###-##-#### 12-23-1949 M VIETNAM ERA 45206  
 Outpatient Accession Date 12-20-1996@1309

\*\*\*\*\*2 HEPATITIS C ANTIBODY POS \*\*\*\*\*

12-20-1996@1309 RIA 1220 68 CHEMISTRY TEST SERUM  
 Hepatitis C Ab 01-03-1997@1347 STRONG POSITIVE -

JOHN DOE ###-##-#### 07-06-1919 M WORLD WAR II 41074  
 Outpatient Accession Date 12-19-1996@1007

\*\*\*\*\*1 VANC-RES ENTEROCOCCUS \*\*\*\*\*

12-19-1996@1007 BACT 96 10618 MICRO CULTURE URINE  
 1 12-23-1996 PSEUDOMONAS AERUGINOSA  
 2 12-23-1996 ENTEROCOCCUS FAECIUM

ORG # 1 12-19-1996@1007 ANTIBIOTIC MIC URINE

Gentamicin	S
Cefazolin	R
Ampicillin	R
Tobramycin	S
TRMSULF	R
Amikacin	S
Cefoxitin	R
Cefotaxime	I
Nitrofurantoin	R
Cefoperazone	S
Mezlocillin	S

## Protocols

**LREPI:** This event driver protocol defines the associated parameters needed to build the HL7 Message used to send data to the AITC.

**LREPI CLIENT:** This subscriber protocol defines the parameter required by the HL7 application that determines where to send the HL7 formatted message containing the EPI data.

## Domain

The EPI software Domain name is Q-EPI-MED.GOV. Sites **must** install Patch XM\*DBA\*103. See VA MailMan V. 7.1 manual for instruction on how to set up the domain after the patch has been installed.

## EPI Mail Groups

The EPI software application has two mail groups. These mail groups are use too electronically transmit data (i.e., HL7 format messages) to the AITC and to receive messages from the AITC. Listed below are recommendations for membership and suggestions for managing the EPI mail groups.

### Office of the Director (00)

It is highly recommended that the Office of the Director (00) at each VAMC initially designate the member(s) responsible for overseeing the two EPI mail groups.

The EPI mail group member(s) is responsible for transmitting the EPI data to the AITC.

The **EPI-Report** mail group member(s) will receive the Emerging Pathogens Verification Report Messages, Confirmation Messages and the Processing Report transmitted by the AITC. This individual(s) will also assist in the EPI data validation, data corrections, and re-transmitting the EPI data to the AITC.

**NOTE:** It is highly recommended that a TQI/QI/QA staff, Laboratory Information Manager (LIM), Microbiology director or supervisor, Infection Control Practitioners, or Hospital Epidemiologist, or individual(s) with similar functions be a member(s) of the mail groups. This member(s) is responsible for making EPI data corrections due to the numerous files from which the data is obtained (e.g., PTF, PIMS, Health Information Management, Laboratory, etc.). Once the corrections are made, it is the responsibility of the EPI mail group member(s) to re-transmit the EPI data to the AITC. These members may also be of assistance with the verification and periodic validation processes.

## EPI Mail Groups

The following information explains the functionality of each mail group.

### **EPI:**

This mail group is used for the transmission of EPI HL7 messages derived from the parameters defined in the LAB SEARCH/EXTRACT file (#69.5) to the AITC.

- Transmit EPI data from the VAMCs (i.e., in HL7 format as defined by the Laboratory Search/Extract software) to the AITC.

### **EPI-Report:**

This mail group is used to deliver the Verification Report message (i.e., created in human readable formatted) which is used to review the data sent to the AITC.

Receives the EPI Verification Report that is generated at each VAMC around the 15<sup>th</sup> of each month. The EPI Verification Report is automatically generated monthly from the HL7 formatted mail message (i.e., in a human readable format). The EPI Verification Report can be generated manually using the Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option and re-run whenever deemed necessary.

Receive the Confirmation Messages

Receive Processing Report from AITC. This report displays any errors that occurred and if a message is unacceptable for receiving by the AITC.

## Adding Mail Groups

Add the EPI mail groups to the HL7 APPLICATION PARAMETER file (#771) using VA FileMan V. 21.0:

### **Example:**

```
Select OPTION: ENTER OR EDIT FILE ENTRIES <RET>

INPUT TO WHAT FILE: HL7 APPLICATION PARAMETER file (#771) <RET>
                      (7 entries)
EDIT WHICH FIELD: ALL// [Enter Facility Name field]<RET>

THEN EDIT FIELD:<RET>

Select HL7 APPLICATION PARAMETER NAME: EPI <RET>                ACTIVE
FACILITY NAME: [Enter your facility name or facility number] <RET>

Select HL7 APPLICATION PARAMETER NAME: EPI-Report<RET>          ACTIVE
FACILITY NAME: [Enter your facility name or facility number] <RET>
```

## Starting the Lower Level Protocol of the HL7 V. 1.6 Background Job

### **Example:**

```
Select Systems Manager Menu Option:    HL7 Main<RET> Menu
```

- 1 V1.5 OPTIONS ...
- 2 V1.6 OPTIONS ...
- 3 Activate/Inactivate Application
- 4 Print/Display Menu ...
- 5 Purge Message Text File Entries

```
Select HL7 Main Menu Option: 2<RET>  V1.6 OPTIONS
```

- 1 Communications Server ...
- 2 Interface Workbench
- 3 Message Requeuer

```
Select V1.6 OPTIONS Option: 1<RET>  Communications Server
```

- 1 Edit Communication Server parameters
- 2 Manage incoming & outgoing filers ...
- 3 Monitor incoming & outgoing filers
- 4 Start LLP
- 5 Stop LLP
- 6 Systems Link Monitor
- 7 Logical Link Queue Management
- 8 Report

```
Select Communications Server Option: 4<RET> Start LLP
```

This option is used to launch the lower level protocol for the  
Appropriate device. Please select the node with which you want  
to communicate

```
Select HL LOGICAL LINK NODE: EPI<RET>
```

The LLP was last shutdown on JAN 30, 1997 12:06:19.

Select one of the following:

F	FOREGROUND
B	BACKGROUND
Q	QUIT

```
Method for running the receiver: B//<RET> ACKGROUND
Job was queued as 131225.
```

# EPI HL7 Format Mail Message

The EPI data is transmitted via an HL7 format mail message to the AITC. The HL7 format mail message should to be transmitted to the AITC by the 15<sup>th</sup> of each month. Forward the HL7 format mail message to the [XXX@Q-EPI](#) domain.

**Example:**

Subj: HL7 Message FEB 10,1997@15:56:29 from Station XXX STATION XXX [#63430]  
10 Feb 97 15:56 262 Lines  
From: POSTMASTER (Sender: ANYBODY) in 'IN' basket. Page 1

MSH|^&|EPI-LAB|EPI-LAB|19970210155618|ORU~R01|483|P|2.2||NE|NE|USA

NTE||R~REPORTING DATE FROM 19960101 TO 19970131~1

```
PID|1|000-54-  
3435~4~M10|17~8~M10||DUCK~DONALD||19840426|M|7|~|~|~|~|~|~|~|~|~|000543435
```

[illegible]

NTE | 1 | 1~VANC-RES ENTEROCOCCUS

```
OBR|1|||81121.0000~CHEMISTRY
TEST~VANLT|||19960312123902| |||||SER~~SERUM|||CH
0312 14
```

```
OBX|1|ST|~~~183~CHOLESTEROL~VA60||190|mg/dL|$S(AGE<50:135,1:120)-
288|||1996
0328145221
```

Select MESSAGE Action: IGNORE (IN basket)// **FORWARD<RET>**

Send mail to: **XXX@Q-EPI** <RET>.MED.VA.GOV via FOC-AUSTIN.VA.GOV

## EPI Confirmation Message

An EPI Confirmation mail message is sent to the sending site EPI-Report mail group after the EPI HL7 formatted mail message has been received by the AITC database.

## EPI Processing Report

A Processing Report mail message is sent to the sending facility EPI-Report mail group at the end of the AITC processing cycle (i.e., the 25<sup>th</sup> of each month). The EPI Processing Report mail message confirms that the EPI data has been processed and lists any errors and/or warning codes requiring corrections.

### **Example: EPI Processing Report**

Subj: EPI/LRK #970451447950300 [#1425971] 14 Feb 97 14:55 CST 50 Lines  
From: <POSTMASTER@FOC-AUSTIN.VA.GOV> in 'IN' basket. Page 1 \*\*NEW\*\*

---

2EPI0001 LRK.  
STATION XXX EPI PROCESSING REPORT REPORT DATE 1997/02/11  
PAGE 01

PROCESS DATE	SSN	ENCOUNTER DATE	MESSAGE	ERROR CODES
19970131	000000190	19970131132151	001	NO ERRORS
19970131	003466370	19961002160512	001	500
19970131	012225556	19970121121609	001	NO ERRORS
19970131	026368799	19961229043131	001	NO ERRORS
19970131	061187860	19960917165748	001	NO ERRORS
19970131	062267678	19970131125821	001	NO ERRORS
19970131	062385860	19961230185116	001	NO ERRORS
19970131	062385860	19970107162010	001	NO ERRORS
19970131	075684002	19970109103417	001	NO ERRORS
19970131	084501279	19970128090146	001	NO ERRORS
19970131	097549144	19970108132918	001	NO ERRORS
19970131	099302298	19970108082142	001	NO ERRORS
19970131	099345601	19970113114752	001	NO ERRORS
19970131	103166277	19970122095702	001	NO ERRORS
19970131	110541525	19970106225241	001	NO ERRORS
19970131	113328446	19970114143128	001	NO ERRORS
19970131	133425965	19970131105820	001	NO ERRORS
19970131	135205512	19960607174229	001	NO ERRORS
19970131	138384641	19970114082310	001	NO ERRORS
19970131	186185119	19960220155121	001	NO ERRORS
19970131	193609042	19970129172826	001	NO ERRORS
19970131	202364130	199701141902	001	NO ERRORS
19970131	218565414	19961011133221	001	500
19970131	229749494	19970123203635	001	NO ERRORS
19970131	250461267	19970121115616	001	NO ERRORS
19970131	253769911	19970115085841	001	NO ERRORS
19970131	254783854	19970118162650	001	NO ERRORS
19970131	259387083	19961028154210	001	NO ERRORS

## EPI Table of Reject and Warning Codes

### Example:

ERROR NUMBER	FIELD NAME	EDIT DESCRIPTION	ERROR MESSAGE
<i>000 Series Miscellaneous</i>			
001	Message Control ID	Must not be blank	Message control ID was blank
002	Batch Sending Facility	Sending Station not valid. (Refer to table AA001)	Invalid Batch Sending Facility.
003	Segment Name	PID Segment missing. Do not edit for the existence of PID when NTE segments are present.	PID Segment missing.
004	Segment Name	PV1 Segment missing. Do not edit for the existence of PV1 when NTE segments are present.	PV1 Segment missing.
005	Segment Name	Invalid Segment name.	Invalid HL7 Segment name.
006	Message Creation Date	Must a valid date.	Message Creation Date is invalid.
007	Message Creation Time	Must a valid time.	Message Creation Time is invalid.
<i>100 Series NTE Totals Segment</i>			
100	Action Ind	Currently not being used.	Currently not being used.
105	Totals Total Count	Must be numeric, if Action Ind is 'T'.	NTE Totals Total Count was not numeric.
110	Negative Input Ind	Must be 'N', if Action Ind is not 'T'.	Negative Input Ind was not 'N'.
<i>200 Series PID Segment</i>			
200	Patient Name	Required. Must be alpha numeric. Must not be all numeric. Must not be all blanks.	Patient Name is missing, or not alphanumeric, or all numeric, or all blanks.
205	Patient Date of Birth	Not required. Must be less than the processing year.	Date of Birth is after the Date of transmission. (Also see W03, W04, and W05)
210	Patient Sex	Not required. Must be blank or match table. (Refer to table 0001)	Sex code is not blank or a valid code. (Refer to table 0001)
215	Patient Race	Not required. Must be blank or a valid code. (Refer to table VA07)	Race code is not blank or a valid code. (Refer to table VA07)
220	Patient Address	Must be blank or 'H'.	Patient Address is not blank or 'H'.

ERROR NUMBER	FIELD NAME	EDIT DESCRIPTION	ERROR MESSAGE
224	Patient Zip Code	Not required. Must be Blank or numeric. If numeric, first five digits must not be all zeros. If last four digits exist, then must be numeric.	Address Zip Code is missing or not numeric.
235	Social Security Number	Required. Last byte must be 'P' or blank.	Pseudo SSN is not 'P' or blank.
236	Social Security Number	Required. Must be numeric. Must be greater than zeros.	Social Security Number is missing, or not numeric, or is equal to zeros.
240	Period of Military Service	Must be a valid code. (Refer to table VA011)	Period of Service was invalid. (Refer to table VA011).
<i>300 Series OBR Segment</i>			
300	Universal Service ID	Must be a valid code.	Invalid Universal Service ID
305	Observation Date	Must be numeric date. Must be a valid date. Must be less than processing date.	Observation Date is invalid date or after the date of transmission
307	Observation Time	Not required Must be blank or numeric. If numeric, must be a valid time.	Observation Time is invalid.
310	Specimen Source Code	Not required. If not blank, must be a valid code.	Invalid Specimen Source Code (also see W07)
315	Parent Observation ID	Not required. Must be blank or a valid code.	Invalid Parent Observation ID.
<i>400 Series PV1 Segment</i>			
400	Patient Class	Required. Must be 'I', 'O', or 'U'.	Patient Class is not 'I', 'O', or 'U'.
410	Discharge Date	Not required. Must be blank or a valid date. Must be less than or equal to processing date.	Discharge Date is invalid or after date of transmission.
411	Discharge Time	Not required. Time must be blank or a valid time.	Discharge Time is invalid
420	Admit Date/Time	Required. Must be numeric. date. Must be less than or equal to processing date.	Admit Date is invalid or after date of transmission.
421	Admit Date/Time	Required. Time must be numeric. Must be a valid time.	Admit Time is invalid.



ERROR NUMBER	FIELD NAME	EDIT DESCRIPTION	ERROR MESSAGE
<i>500 Series DG1 Segment</i>			
500	Diagnosis Code	Required. Must be a valid code.	Invalid Diagnosis Code.
<i>600 Series OBX Segment</i>			
600	Observation Nat Lab Num.	If not blank, must be a valid code.	Invalid Observation Nat Lab Num. (Also see W09)
605	Final Result Date	Must be blank or a valid date. Must be numeric date. Must be a less than or equal to the processing date.	Final Result. Date is invalid or after the date of transmission.
<i>W00 Series Warnings</i>			
W03	Patient Date of Birth	Must not be all spaces.	Patient Date of Birth is all spaces. (Also see 205)
W04	Patient Date of Birth	Year must not be all zeros.	Patient Date of Birth Year is all zeros. (see also 205)
W05	Patient Date of Birth	Must be a valid date.	Patient Date of Birth is not in a valid date format. (Also see 205)
W07	Specimen Source Code	Blanks in code.	Specimen Source code is blank. (see also 310)
W09	Observation Nat Lab Num.	Blanks in code.	Observation Nat Lab Num is blank. (Also see 600)



# APPENDIX B

## NCH USER GUIDE



# NCH User Guide

The NCH User Guide provides all the necessary information, instructions, illustrations, and examples required to implement and maintain the **VISTA** Laboratory Search/Extract software application.

## Overview

**National Center for Health Promotion Cholesterol and Pap Screening (NCH):** Pursuant to the Congressional mandates stipulated in its enabling legislation (PL. 102-585, U.S.C. 17, 1704), the Department of Veterans Affairs (DVA) National Center for Health Promotion (NCHP) is tasked with monitoring and improving the prevalence of health promotion screening activities provided to veterans VA-wide. The purpose of the NCH Cholesterol and Pap Screening database is to monitor cholesterol and Pap screening activities at a national level, with the ultimate goal for improving detection and treatment of hyperlipidemia and cervical cancer. The NCH Cholesterol and Pap Screening database will also provide a valuable resource for clinical and health services researchers on screening activities and health outcomes, particularly for high-risk and special emphasis group such as veterans with hyperlipidemia, older veterans, and female veterans.

The NCH database is an ongoing initiative. When completed, the NCH database shall consist of a serial (i.e., annual) patient-level database of Cholesterol and Pap Screening procedures conducted in DVA facilities nationally.

At minimum, the NCH database records the following characteristics of the screening procedures:

1. Records the DVA facility where the procedure was conducted.
2. Records the Social Security Number (SSN) of the patient.
3. Records the date the screening procedure was conducted or the date the specimen taken.
4. Records the type of test provided.
5. Records the result of the screening procedure.

## Mandate

The installation of the **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software should be completed in accordance with PL. 102-585, U.S.C. 17, 1704. It is recommended that all VAMCs installed and implemented the software as soon as possible to help facilitate these important initiatives.

## NCH Database Access

Access to the NCH database is currently restricted to the NCHP staff. Upon completion of the NCH Cholesterol and Pap Screening initiative the NCHP office will initiate a data-sharing protocol that ensures adherence to the Privacy Act (5 U.S.C.552a).

## Impact

Clinicians and clerks in the Laboratory Service are the primary users of the software. Front-line managers, VISN directors, clinical managers, and health researchers may also have the need to use the software. Telecommunications issues are paramount. The **VISTA** Laboratory Search/Extract software automatically generates and transmits the NCH data to the AITC database daily.

## National Roll-Up

The NCH database is stored at the AITC. The database is compiled in SAS format (i.e., serialized by year). Some NCH database routine reports are generated by the AITC. The NCHP may also generate routine and specialized reports as deemed necessary and appropriate. Data elements included in the national roll-up are proportion of unique patients screened by age, type of test and by date range.

The Denominator file criteria is developed by the NCHP for prevalence ratios for screening activities by the DVA facilities and VISNs. The volume of data expected is very large. It is estimated that the number of records in the NCH database may approach 500,000 annually. Funding is provided by the NCHP to the AITC for NCH database storage and report generation.

## NCH Search and Extract Criteria

The **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 software extracts the NCH data directly from the **VISTA** database. The patient social security number, accession date, and time are used as the unique identifier for the search/extract criteria. The **VISTA** Laboratory Search/Extract software does **not** add any additional data to the patient medical record. However, because the unique identifier is appended to the record. The NCH data can be linked to the patient record if deemed appropriate. All NCH data extracted by the **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175 are already routinely collected. No new equipment or staffing is necessary to use the **VISTA** Laboratory Search/Extract software.

## Recommended Users

It is **highly recommended** that the following person(s) jointly participate in reviewing the parameter descriptions:

Laboratory Information Manager (LIM)

Total Quality Improvement/Quality Improvement/Quality Assurance (TQI/QI/QA) staff (or person at the site with similar function)

Representative from the Microbiology section for the Emerging Pathogens Initiative (i.e., director, supervisor, or technologist)

## Periodic Reviews

The NCH parameter descriptions require an ongoing review process (as specified by the NHCP). The person (s) participating in the ongoing review process is responsible for ensuring the following requirements are current:

Periodic reviews of the ICD codes

Periodic reviews of the Lab Search/Extract Parameter Setup

Annual review of the NCHP parameter descriptions

## NCH Data Transmission

### **NCH HL7 formatted mail messages**

The **VISTA** Laboratory Search/Extract software automatically processes the NCH data daily producing HL7 formatted mail messages. A Verification Report mail message is then produced for each HL7 mail message. The Verification Report is a copy of the HL7 formatted mail messages in a human readable format.

The NCH HL7 formatted mail messages are daily transmitted automatically by the **VISTA** Laboratory Search/Extract software to the AITC database.

### **NCH Verification Report mail message**

The NCH Verification Report is sent to the LR NCH-Report mail group. The members of the LR NCH-Report mail group may use this Verification Report to review the NCH data and make correction as deemed necessary (e.g., social security number, valid date of births, period of service, etc.).

The new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option can be used to generate the NCH Verification Report and to transmit NCH data corrections to the AITC database.

**NOTE:** See Appendix B for an example of the NCH Verification Report mail message.

### **NCH Acknowledgment mail message**

The NCH Acknowledgment mail message displays the status of the NCH HL7 formatted mail message transmission to the AITC. The NCH Acknowledgment mail message is sent to the sending facility LR NCH-Report mail group.

### **NCH VA View Alert mail message**

A NCH VA Alert mail message is sent to the LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4), Report Mail Group field (#1), after the NCH data has been processed by the AITC.



## Cholesterol Screening

Cholesterol Screening is a critical mechanism in the detection and treatment of hyperlipidemia and coronary artery disease, particularly in high-risk patients. Hyperlipidemia is defined as an undesirable level of plasma lipids. The major lipids are cholesterol and triglycerides. Evidence suggests that elevated total cholesterol, elevated low-density lipoprotein cholesterol, and reduced high-density lipoprotein cholesterol contribute to the development of coronary artery disease. Treatment of hyperlipidemia involves primary and secondary prevention along with pharmacological intervention.

The NCH database is used to assist VA clinicians, Planners, and Researchers in the prevention and treatment of hyperlipidemia by providing estimates of screening rate, age- and sex-adjusted hyperlipidemia prevalence, and actual values for total low-density and high-density lipoprotein cholesterol. The VHA Health Promotion and Disease Prevention Handbook (VHA Handbook 1101.8) recommends that males ages 35 to 65 and females ages 45 to 65 be screened for hyperlipidemia every five years.

Laboratory Services have tended to identify cholesterol screening using the following test names:

**Cholesterol**  
**Triglyceride**

**HDL Cholesterol**  
**LDL Cholesterol**

### **Cholesterol Screening for Hyperlipidemia**

Hyperlipidemia is a major risk factor for cardiovascular disease. Cholesterol Screening detects abnormal blood lipid levels; and subsequent interventions to normalize these levels are effective strategies in the prevention of coronary artery disease. Abnormal blood lipids that are associated with cardiovascular disease are elevated total cholesterol, triglycerides and low density lipoprotein (LDL) levels, and reduced high density lipoprotein (HDL) levels.

**Example: Lab Search/Extract Parameter Setup for NCH Cholesterol**

```

                                Lab Search/Extract Primary Menu

    ENH    Lab Search/Extract Manual Run (Enhanced)
    VR     Print Detailed Verification Report
    LO     Local Pathogen Menu ...
    PI     Pathogen Inquiry
    UP     Lab EPI Parameter Setup
           Lab EPI Protocol Edit
    LK     Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
  Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
  Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
  CANDIDA
  CLOSTRIDIUM DIFFICILE
  CREUTZFELDT-JAKOB DISEASE
  CRYPTOSPORIDIUM
  DENGUE
  E. COLI 0157:H7
  HEPATITIS C ANTIBODY POS
  LEGIONELLA
  LEISHMANIASIS
  MALARIA
  NCH CHOLESTEROL
  NCH PAP SMEAR
  PEN-RES PNEUMOCOCCUS
  STREPTOCOCCUS GROUP A
  TUBERCULOSIS
  VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: NCH CHOLESTEROL<RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 1 of 5
NAME: NCH Cholesterol		ACTIVE: YES
Laboratory Test(s)	Indicator	Value
CHOLESTEROL		
TRIGLYCERIDE		
HDL CHOLESTEROL		
LDL CHOLESTEROL		
<RET>		
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>		
ICD Code	Cd Set	ICD Description
<RET>		
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		Press <PF1>H for help
		Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: NCH Cholesterol		ACTIVE: YES
Selected Etiology	Selected Snomed Codes	
Antimicrobial Susceptibility	NLT Code	NLT Description
<RET>		
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		Press <PF1>H for help
		Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 3 of 5
NAME: NCH Cholesterol		ACTIVE: YES
Topography Selection		
Include	Exclude	
Exit	Save	Next Page
Refresh		
COMMAND: N<RET>		Press <PF1>H for help
		Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: NCH Cholesterol	ACTIVE: YES	
<hr/>		
FIRST ENCOUNTER:<RET>	Follow PTF: <b>YES</b> <RET>	
BEFORE DATE OF BIRTH:<RET>	AFTER DATE OF BIRTH:<RET>	
Select SEX:<RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: <b>E</b> <RET> Press <PF1>H for help Insert		
Save changes before leaving form (Y/N)? <b>Y</b> <RET>		

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: NCH Cholesterol	ACTIVE: YES	
<hr/>		
Run Date:<RET>	Protocol: <b>LRNCH</b> <RET>	
Run Cycle: <b>DAILY</b> <RET>	Lag Days: <b>10</b> <RET> ( <b>Site Specific</b> )	
General Description:<Tab>		
<hr/>		
Exit	Save	Refresh
COMMAND: <b>E</b> <RET> Press <PF1>H for help Insert		
Save changes before leaving form (Y/N)? <b>Y</b> <RET>		

## Papanicolaou (Pap) Screening

Pap Screening is intended to detect cervical dysplasia and cancer at the earliest time possible. Cervical cancer is a disease associated with considerable morbidity and mortality. Approximately 16,000 new cases of cervical cancer are detected in the United States each year with nearly 5,000 deaths.

### **Pap Smear for Cervical Cancer**

Currently, the U. S. Preventive Service Task Force and the VHA Health Promotion and Disease Prevention Handbook (VHA Handbook 1101.8) recommend that Pap testing ensure with the onset of sexual activity. Pap testing should be repeated every three years with a possible top age limit for testing at 65, although evidence for this age limit is equivocal.

#### **Example: Lab Search/Extract Parameter Setup for NCH Pap smear**

```

                                Lab Search/Extract Primary Menu

ENH      Lab Search/Extract Manual Run (Enhanced)
VR        Print Detailed Verification Report
LO        Local Pathogen Menu ...
PI        Pathogen Inquiry
UP        Lab EPI Parameter Setup
          Lab EPI Protocol Edit
LK        Antimicrobial Link Update

Select Lab Search/Extract Primary menu Option: UP<RET> Lab Search/Extract
Parameter Setup

Select LAB SEARCH/EXTRACT NAME: ?<RET>
Answer with LAB SEARCH/EXTRACT NAME, or REFERENCE NUMBER
Do you want the entire 16-Entry LAB SEARCH/EXTRACT List? Y (Yes)<RET>
Choose from:
CANDIDA
CLOSTRIDIUM DIFFICILE
CREUTZFELDT-JAKOB DISEASE
CRYPTOSPORIDIUM
DENGUE
E. COLI 0157:H7
HEPATITIS C ANTIBODY POS
LEGIONELLA
LEISHMANIASIS
MALARIA
NCH CHOLESTEROL
NCH PAP SMEAR
PEN-RES PNEUMOCOCCUS
STREPTOCOCCUS GROUP A
TUBERCULOSIS
VANC-RES ENTEROCOCCUS

Select LAB SEARCH/EXTRACT NAME: NCH PAP SMEAR<RET>

```

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 1 of 5
NAME: NCH PAP SMEAR		ACTIVE: YES	
Laboratory Test(s)	Indicator	Value	
PAP SMEAR			
<RET>			
ICD Coding System [ICD-9 or ICD-10]? (9/10):<RET>			
ICD Code	Cd Set	ICD Description	
<RET>	<RET>		
Exit	Save	Next Page	Refresh
COMMAND: N<RET>		Press <PF1>H for help	Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: NCH PAP SMEAR		ACTIVE: YES	
Selected Etiology		Selected Snomed Codes	
Antimicrobial Susceptibility	NLT Code	NLT Description	
<RET>			
Exit	Save	Next Page	Refresh
COMMAND: N<RET>		Press <PF1>H for help	Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 3 of 5
NAME: NCH PAP SMEAR		ACTIVE: YES	
Topography Selection			
Include	Exclude		
Exit	Save	Next Page	Refresh
COMMAND: N<RET>		Press <PF1>H for help	Insert

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 4 of 5
NAME: NCH PAP SMEAR	ACTIVE: YES	
<hr/>		
FIRST ENCOUNTER: <b>NO</b> <RET>	Follow PTF: <b>YES</b> <RET>	
BEFORE DATE OF BIRTH: <RET>	AFTER DATE OF BIRTH: <RET>	
Select SEX: <RET>		
<hr/>		
Exit	Save	Refresh
COMMAND: <b>E</b> <RET> Press <PF1>H for help Insert		
Save changes before leaving form (Y/N)? <b>Y</b> <RET>		

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 5 of 5
NAME: NCH PAP SMEAR	ACTIVE: YES	
<hr/>		
Run Date: <RET>	Protocol: <b>LRNCH</b> <RET>	
Run Cycle: <b>DAILY</b> <RET>	Lag Days: <b>10</b> <RET> (Site Specific)	
General Description: <Tab>		
<hr/>		
Exit	Save	Refresh
COMMAND: <b>E</b> <RET> Press <PF1>H for help Insert		
Save changes before leaving form (Y/N)? <b>Y</b> <RET>		

## NCH Mail Messages

The following are examples of mail messages and reports that are produced by the **VISTA** Laboratory Search/Extract Patch LR\*5.2\*175.

### NCH HL7 Formatted Mail Message

The **VISTA** Laboratory Search/Extract software automatically transmits the NCH HL7 formatted mail messages daily to the AITC, formerly AITC, database.

**Example:** NCH HL7 formatted mail messages transmission to the AITC database

```
MSH|^~^&|NCH-LAB|525|NCH-
AITC|200|19980505160420||ORU~R01|1810579|6|2.2||NE|AL|USA
NTE|R~REPORTING DATE FROM 19980422 TO 19980422~1
PID|1|000-10-5466~9~M10|69604~5~M10|EPI~PATFOUR|00041027|M|6|~01603|||||
|034105466|||||2
PV1|1|O|||||||19980422080566
NTE|1|50~NCH CHOLESTEROL
OBR|1||81121.0000~CHEMISTRY
TEST~VANLT||19980422080566|||||SER~~SERUM||R/
CH 0422 340
OBX|1|ST|82466.0000~Cholesterol Total~VANLT~183~CHOLESTEROL~VA60||165|mg/dl|"
SEE TEST DESC"-|||||19980423140447OBX|2|ST|83013.0000~Cholesterol
HDL~VANLT~244~HDL CHOLESTEROL~VA60||30|mg/dl|30-65|||||19980423140447
OBX|3|ST|82350.0000~Calculation~VANLT~901~LDL CHOLESTEROL~VA60||109c|mg/dl|"
```



## NCH Verification Report mail message

The NCH Verification Report mail message may be sent to the LR NCH Report mail group. Members of the LR NCH Report mail may choose to review the NCH data and make corrections (e.g., social security number, valid date of births, period of service, etc.) as deemed appropriate. The new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option automatically transmits the NCH data corrections to the AITC database.

### **Example:** NCH Verification Report mail message.

MailMan message for DOE,DAVID R LAB PROGRAMMER

Printed at BROCKTON.VA.GOV 05 May 98 16:06

Subj: National Center for Health Promotion [#13326619] 05 May 98 16:04  
21 Lines

From: POSTMASTER in 'IN' basket. Page 1

-----  
REPORTING DATE FROM 04-22-1998 TO 04-22-1998 Message Seq # 1 manual

EPI PATFIVE 000-07-5914 10-27-1925 M WORLD WAR II 01603  
Outpatient Accession Date 04-22-1998@0805

\*\*\*\*\* 50 NCH CHOLESTEROL \*\*\*\*\*

04-22-1998@0805 R/CH 0422 340 CHEMISTRY TEST SERUM  
CHOLESTEROL 04-23-1998@1404 165 mg/dl " SEE TEST DESC"-  
HDL CHOLESTEROL 04-23-1998@1404 30 mg/dl 30-65  
LDL CHOLESTEROL 04-23-1998@1404 109c mg/dl " SEE TEST DESC"-  
TRIGLYCERIDE 04-23-1998@1404 160 mg/dl 35-200  
TOTAL 82466.0000 CHOLESTEROL 1  
TOTAL 82466.0000 PATIENTS WITH CHOLESTEROL 1  
TOTAL TRIGLYCERIDE 1  
TOTAL PATIENTS WITH TRIGLYCERIDE 1  
TOTAL 83013.0000 HDL CHOLESTEROL 1  
TOTAL 83013.0000 PATIENTS WITH HDL CHOLESTEROL 1  
TOTAL 82350.0000 LDL CHOLESTEROL 1  
TOTAL 82350.0000 PATIENTS WITH LDL CHOLESTEROL 1

## NCH HL7 Mail Message Status List

The HL7 Message Status List contains the status of each transmission sent by the site to the AITC database.

### **Example: HL7 Message Status List**

```

HL7 MESSAGE STATUS LIST                JUL  8,1998  14:02  PAGE 1
-----
NAME: APPLICATION LEVEL ERROR          CODE: ALE
DESCRIPTION:  An error occurred while the application was processing the
message. The Error Message field of the Message Text file provides a
description of the error.

NAME: AWAITING ACKNOWLEDGEMENT         CODE: AA
DESCRIPTION:  The message has been sent, but has not yet been acknowledged
as having been processed by the receiving application.

NAME: AWAITING PROCESSING               CODE: AP
DESCRIPTION:  The message is awaiting processing by the receiving
application.

NAME: BEING GENERATED                 CODE: BG
DESCRIPTION:  The message is in the process of being generated and stored in
the Message Text file by the Messaging System.

NAME: ERROR DURING GENERATION          CODE: EDG
DESCRIPTION:  An error occurred during the generation of the message.

NAME: ERROR DURING PROCESSING          CODE: EDP
DESCRIPTION:  An error occurred during the processing of the message.

NAME: ERROR DURING TRANSMISSION        CODE: EDT
DESCRIPTION:  An error occurred during the transmission of the message.

NAME: PENDING TRANSMISSION             CODE: PT
DESCRIPTION:  The message is waiting to be transmitted.

NAME: SUCCESSFULLY COMPLETED          CODE: SC
DESCRIPTION:  One of the following actions was accomplished for the message:

The message was successfully generated. The message was successfully sent and
no acknowledgment is required. The message was successfully sent and
acknowledged.

```

## NCH HL7 Formatted (Acknowledgment) mail message

The following is an example of the NCH HL7 formatted Acknowledgment mail message received from the AITC.

**Example:** NCH HL7 formatted Acknowledgment mail message.

```

DATE/TIME ENTERED: MAY 18, 1998@22:13:27
SERVER APPLICATION: NCH-AITC          TRANSMISSION TYPE: INCOMING

MESSAGE ID: 1855276                PARENT MESSAGE: MAY 17,
1998@02:08:05

PRIORITY: IMMEDIATE                  RELATED EVENT PROTOCOL: LRNCH
MESSAGE TEXT:
MSH^~|\&^NCH-AITC^200^NCH-
LAB^525^19980518210822^^ACK~^1855276^P^2.2^^NE^AL|
MSA^AA^1855276
DSP^0001^9999999999::000
DSP^0001^006309421:19980512125558:000 3
DSP^0001^009381088:19980514103659:000

```

## NCH Acknowledgment mail message

The NCH Acknowledgment mail message displays the status of the NCH HL7 formatted Acknowledgment mail message in a human readable format. The NCH Acknowledgment mail message is sent to the LR NCH-Report mail group. If the NCH Acknowledgment mail message displays an error, make corrections (e.g., social security number, valid date of births, period of service, etc.) as deemed necessary. Use the new Lab Search/Extract Manual Run (Enhanced) [LREPI ENHANCED MANUAL RUN] option to automatically transmit the NCH data corrections to the AITC database.

**Example:** NCH Acknowledgment mail messages in a human readable format.

```
Subj: Acknowledgment message from NCH. [#13202750] 22 Apr 98 16:36 3 Lines
From: DOE,DAVID K in 'IN' basket. Page 1
-----
```

```
Flash... This is an acknowledgment message from NCH
Message #: 1 From the status file 771.6 message is categorized as:
Message: 1746675 AWAITING PROCESSING @ Apr 22, 1998@16:32:48
```

```
Subj: Acknowledgement message from Austin. [#13224503] 24 Apr 98 11:32 3 Lines
From: POSTMASTER - POSTMASTER (Sender: DOE,DAVID R - LAB PROGRAMMER) in 'LAB
SEARCH' basket. Page 1
-----
```

```
Flash... This is an acknowledgement message from NHC
Message #: 1 From the status file 771.6 message is categorized as:
Message: 1733492 ERROR DURING TRANSMISSION @ Apr 14, 1998@09:15:44
```

```
Subj: Acknowledgement message from NHC. [#13128237] 14 Apr 98 19:31 3 Lines
From: DOE,EPI R - LAB PROGRAMMER in 'LAB SEARCH' basket. Page 1
-----
```

```
Flash... This is an acknowledgement message from NHC
Message #: 1 From the status file 771.6 message is categorized as:
Message: 1711986 SUCCESSFULLY COMPLETED @ Apr 08,
```

## NCH VA Alert

A NCH VA Alert mail message is sent to the LABORATORY SEARCH/EXTRACT PROTOCOL file (#69.4), Report Mail Group field (#1) after NCH transmissions has been processed by the AITC.

### Example: NCH VA Alert mail message from the AITC

```
1.I National Center for Health Promotion Was processed at JUN 30, 1998@03:
2.I National Center for Health Promotion Was processed at JUN 29, 1998@03:
3.I National Center for Health Promotion Was processed at JUN 28, 1998@02:
4.I National Center for Health Promotion Was processed at JUN 27, 1998@02:
5.I National Center for Health Promotion Was processed at JUN 26, 1998@03:
```

Select from 1 to 5 or enter ?, A I, F, P, M, R, or ^ to exit: A

Processed Alert Number 1

National Center for Health Promotion Was processed at JUN 30,  
1998@03:08:20

Processed Alert Number 2

National Center for Health Promotion Was processed at JUN 29,  
1998@03:11:13

Processed Alert Number 3

National Center for Health Promotion Was processed at JUN 28,  
1998@02:33:00

Processed Alert Number 4

National Center for Health Promotion Was processed at JUN 27,  
1998@02:25:11

Processed Alert Number 5

National Center for Health Promotion Was processed at JUN 26,  
1998@03:19:22



# APPENDIX-C

## EDITING FILES, LINKING DATA, EDITING SCREENS, REQUEST FORM





## Editing Files/Screens, Linking Data, Request Form

This section contains instructions for, editing files, linking data, and a Workload and Suffix Code Request Form (that can be reproduced).

### Editing TOPOGRAPHY file (#61)

Specific HL7 codes **must** be added to the TOPOGRAPHY file (#61). The HL7 Code field (#08) in this file is used to add the entries. Specific HL7 codes that **must** be added to File #61 is located in the HL7 section of this guide, Table 0070 (Specimen Source Codes). The following is an example of how to add the specific HL7 codes to the TOPOGRAPHY file (#61) using VA FileMan - Enter Or Edit File Entries [ ] option.

#### **Example:** How to add specific HL7 codes to TOPOGRAPHY file (#61)

Select OPTION:       **ENTER** OR EDIT FILE ENTRIES

INPUT TO WHAT FILE: TOPOGRAPHY FIELD// **<RET>**

EDIT WHICH FIELD: ALL// **.08** HL7 CODE

THEN EDIT FIELD: **<RET>**

Select TOPOGRAPHY FIELD NAME: ? **<RET>**

Answer with TOPOGRAPHY FIELD NAME, or SNOMED CODE, or ABBREVIATION, or SYNONYM

Do you want the entire 8575-Entry TOPOGRAPHY FIELD List? **NO<RET>**

You may enter a new TOPOGRAPHY FIELD, if you wish

ANSWER MUST BE 2-80 CHARACTERS IN LENGTH

Select TOPOGRAPHY FIELD NAME:       AMNIOTIC FLUID                       8Y300

HL7 CODE: ? **<RET>**

Answer must be 2-4 characters in length.

Enter the two to four character code from the left column:

ABS	ABCs
AMN	Amniotic fluid
ASP	Aspirate
BPH	Basophils
ABLD	Blood arterial
BBL	Blood bag
BON	Bone
BRTH	Breath
BRO	Bronchial
BRN	Burn

Enter RETURN to continue or '^' to exit: ^

HL7 CODE: **AMN<RET>**

## How to Link Antimicrobial Entries to Workload Codes Entries

The post INIT links as many of the ANTIMICROBIAL SUSCEPTIBILITY file (#62.06) entries to the WKLD CODE file (#64) entries that are identified in your site files. However, the ANTIMICROBIAL SUSCEPTIBILITY file (#62.06) entries that were **not** linked (i.e. no match found) to the WKLD CODE file (#64) will require linking. The Antimicrobial Link Update [LREPILK] option contains the following three options that are used to identify and link the entries that were **not** linked by the post INIT.

### Using the Antimicrobial Link Update [LREPILK] options

#### Example:

Select Lab Search/Extract Primary Menu<RET>

```

ENH    Lab Search/Extract Manual Run (Enhanced)
VR      Print Detailed Verification Report
LO      Local Pathogen Menu ...
PI      Pathogen Inquiry
UP      Lab EPI Parameter Setup
        Lab EPI Protocol Edit
LK      Antimicrobial Link Update

```

Select Lab Search/Extract Primary Menu Option: **LK <RET>** Antimicrobial Link Update

This option will allow you to link file '62.06 ANTIMICROBIAL SUSCEPTIBILITY' file with file '64 WKLD CODE'.

Select one of the following:

```

A      AUTO
M      MANUAL
S      SEMI-AUTO

```

#### How to link entries using the AUTO option

The AUTO option identifies and attempts to link entries that are not currently linked. This option also displays linked and non linked entries.

#### Example:

Enter response: **A<RET>**UTO

```

AMIKACN      <----Linked---->    Amikacin
AMPICLN      <----Linked---->    Ampicillin
CLINDAM      <----Linked---->    Clindamycin
POLYMYXIN B  <----Not Linked----> No Match Found
RIFAMPIN     <----Linked---->    Rifampin

```

#### How to add and delete entries to a file using the MANUAL option

The MANUAL option will create or edit the links. Selections are from entries in the ANTIMICROBIAL SUSCEPTIBILITY file (#62.06).

### Example:

```
Enter response: M<RET>ANUAL
Select ANTIMICROBIAL SUSCEPTIBILITY NAME: PENICLIN<RET>          PENICILLIN
NATIONAL VA LAB CODE: Substance P// PEN<RET>
    1  PENFIELD AND CONE STAIN                88010.0000
    2  PENICILLIN Penicillin                   81852.0000
    3  PENTAZOCINE Pentazocine                 81854.0000
    4  PENTOBARBITAL Pentobarbital             81856.0000
CHOOSE 1-4: 2 Penicillin<RET>
Select ANTIMICROBIAL SUSCEPTIBILITY NAME: VANCMCN<RET>          VANCOMYCIN
NATIONAL VA LAB CODE: Shell Vial Technique// VANCOMYCIN<RET> Vancomycin
81485.0000<RET>
Select ANTIMICROBIAL SUSCEPTIBILITY NAME: Ampicillin/sulbactam<RET>
Ampicillin/subalctam
NATIONAL VA LAB CODE: Ampicillin// @<RET>
    SURE YOU WANT TO DELETE? Y (Yes)<RET>
Select ANTIMICROBIAL SUSCEPTIBILITY NAME:
```

### How to add entries using the SEMI-AUTO option.

The SEMI-AUTO option looks for entries that are not currently linked and prompts the user to select the corresponding entry in the WKLD CODE file (#64).

### Example:

```
Enter response: S<RET>EMI-AUTO

AMIKACN          AMIKACIN
NATIONAL VA LAB CODE: AMIK<RET>ACIN  Amikacin          81098.0000
Continue YES/<RET>

AMPICLN          AMPICILLIN
NATIONAL VA LAB CODE: AMP<RET>
    1  AMP CYCLIC                81029.0000
    2  AMPHETAMINE Amphetamine   81528.0000
    3  AMPHOTERICIN B Amphotericin B 81530.0000
    4  AMPICILLIN Ampicillin       81532.0000
CHOOSE 1-4: 4 Ampicillin
Continue YES// <RET>

CLINDAM          CLINDAMYCIN
NATIONAL VA LAB CODE: CLINDAMYCIN  Clindamycin      81676.0000
Continue YES// <RET>

CARBCLN          CARBENICILLIN
NATIONAL VA LAB CODE:
Continue YES// NO <RET>
```

## How to Delete an Entry from the Laboratory Search/Extract Parameters Input Screen

Use the tab key to move the cursor. Highlight the entry that is to be deleted, select the “@” symbol, then press enter/return. You will then receive a deletion warning asking if you are sure.

### **Example:** Deleting an Entry

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN			Page 2 of 5
NAME: CANDIDA		ACTIVE: YES	
Selected Etiology			
CANDIDA PARAPSILOSIS <Tab>			
CANDIDA PSEUDOTROPICALIS <Tab>			
<b>CANDIDA SKIN TEST ANTIGEN @ &lt;Ret&gt;</b>			
CANDIDA STELLATOIDEA			
Antimicrobial Susceptibility	NLT Code	NLT Description	
<Tab>			
Exit	Save	Next Page	Refresh
COMMAND: Press <PF1>H for help			
WARNING: DELETIONS ARE DONE IMMEDIATELY!			
(EXITING WITHOUT SAVING WILL NOT RESTORE DELETED RECORDS.)			
Are you sure you want to delete this entire Subrecord (Y/N)? <b>y &lt;Ret&gt;</b>			

## How to add an entry using the Laboratory Search/Extract Parameters Input Screen

Use the tab key to move the cursor. Highlight a blank line where the entry is to be added.

### Example:

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: CANDIDA	ACTIVE: YES	
Selected Etiology		
CANDIDA		
CANDIDA GUILLIERMONDII		
<b>CAN &lt;Ret&gt;</b>		
Antimicrobial Susceptibility	NLT Code	NLT Description
<b>&lt;Tab&gt;</b>		
1	CAN CANDIDA ALBICANS	4081
2	CANARYPOX VIRUS	3604
3	CANDICIDIN	7328
4	CANDIDA, NOS	4080
5	CANDIDA GUILLIERMONDII	4082
Choose 1-5 or '^' to quit: <b>1 &lt;Ret&gt;</b>		

LABORATORY SEARCH/EXTRACT PARAMETERS INPUT SCREEN		Page 2 of 5
NAME: CANDIDA	ACTIVE: YES	
Selected Etiology		
CANDIDA GUILLIERMONDII		
CANDIDA KRUSEI		
CANDIDA ALBICANS      <- The entry will appear after answering yes to the adding a new ETIOLOGY prompt.		
Antimicrobial Susceptibility	NLT Code	NLT Description
<b>&lt;Tab&gt;</b>		
CAN CANDIDA ALBICANS		
Are you adding 'CANDIDA ALBICANS' as a new ETIOLOGY? <b>Y &lt;Ret&gt;</b>		

## Additional Workload and Suffixes Codes Request Form

This Request Form may be reproduced and used for requesting additional Workload and Suffixes codes. Please submit the Request Form to the address located at bottom of form.

### **Additional Workload and Suffixes Codes Request Form**

Site Name:\_\_\_\_\_ Site Number:\_\_\_\_\_ Date:\_\_\_\_\_

Contact Person:\_\_\_\_\_ Commercial PH#\_\_\_\_\_ EXT\_\_\_\_\_

FTS PH#:\_\_\_\_\_ EXT\_\_\_\_\_

Procedure Name \_\_\_\_\_ Lab Section \_\_\_\_\_

Abbreviations: \_\_\_\_\_

Procedure Name \_\_\_\_\_ Lab Section \_\_\_\_\_

Abbreviations: \_\_\_\_\_

Procedure Name \_\_\_\_\_ Lab Section \_\_\_\_\_

Abbreviation: \_\_\_\_\_

Method: \_\_\_\_\_ Lab Section \_\_\_\_\_

Abbreviation: \_\_\_\_\_

Method: \_\_\_\_\_ Lab Section \_\_\_\_\_

Abbreviation: \_\_\_\_\_

Method: \_\_\_\_\_ Lab Section \_\_\_\_\_

Abbreviation: \_\_\_\_\_

Instrument Name: \_\_\_\_\_ Manufacturer's Name: \_\_\_\_\_

Instrument Name: \_\_\_\_\_ Manufacturer's Name: \_\_\_\_\_

Instrument Name: \_\_\_\_\_ Manufacturer's Name: \_\_\_\_\_

Instrument Name: \_\_\_\_\_ Manufacturer's Name: \_\_\_\_\_

Instrument Name: \_\_\_\_\_ Manufacturer's Name: \_\_\_\_\_

**Submit Request Forms to:** *Frank Stalling, P&LMS Informatics Manager*  
*1901 North Highway 360, Suite 351*  
*Grand Prairie, Texas 75050*