

RADIOLOGY/NUCLEAR MEDICINE ADPAC GUIDE

(Package Implementation and Maintenance)

Version 5.0 December 1999

(Revised June 2019)

Department of Veterans Affairs
VistA Health System Design and Development

Revision History

Date	Page	Change
July 2002	V-8, V-13, V-16	Added a Division Parameter for HL7 Applications
July 2002	IX-41	Added HL7 Radiology Menu
September 2002	VII-45	New entries and modifications to existing entries are prohibited
		without approval from Radiology Service VACO.
September 2002	VII-89	Added new option VistaRad Category Entry/Edit
September 2002	VIII-39	Added new option VistaRad Category Print
May 2003	V-2	Patch RA*5*34: Inserted paragraphs regarding if parameter is set
		to NO or not answered at all.
May 2003	VII-6	Patch RA*5*34: Moved 'PRINT ON ABNORMAL REPORT'
		entry up to align with top of row.
May 2003	VII-64	Patch RA*5*34: Inserted text regarding semicolon restriction on
		procedure name.
May 2003	XIV-5	Patch RA*5*34: Inserted text regarding several cases sharing one
		printset.
March 2004	VII-2, 3, 68, 69, VIII-21,	Patch RA*5*45: Inserted text regarding no longer using AMIS
	32, XI-1, XII-5, XVII-1,	codes to identify contrast medial use, and associating
	XVII-6	parent/descendant procedures that use contrast media.
December 2004	VIII-22-28, 32, 35, 37,	Patch RA*5*45. Inserted documentation for new Procedure File
	40, 43	Listings: HIST Display Rad/NM Procedure Contrast Media
		History and List of Procedures with Contrast.
		Inserted a Note for several Procedure File Listings that contrast
		media associations will display, when appropriate.
December 2016	V-3,4,10	Patch RA*5*133. Document changes to Reason Enter/Edit and
	VII-82	Division parameter Set-up.
June 2017	IX 7-8	Patch RA*5*137. Documentation added for new item on the RA
		HL7 Menu.
	V-35	Add Registered Request Printer to 'Location Parameter List' [RA
		SYSLOCLIST] output.
	VII-28	Add pregnancy screen to Exam Status Tracking example.
August 2018	III-4, IV-1	Patch RA*5*124. Documentation added for new RA
August 2016	111-4, 1 V - 1	UNVERIFY security key.
June 2019		Patch RA*5.0*158: Update the 'Procedure Modifier Entry' [RA
Julic 2017	I-3, II-9, VII-1,81-83,	MODIFIER] description, remove the 'Reason Edit' option [RA
	XV-2	REASON EDIT] and update of the FAQ (section XV) pertaining
	2X V - 2	to procedure edits.
		to procedure cures.

Table of Contents

I.	INTRODUCTION	
	Implementation Checklist	
	Supervisor Menu	
	Division, Location, Imaging Type Definitions	
	■ Parent/Descendent Exam and Printsets	I-7
	■ Radiopharmaceuticals	I-7
II.	IMPLEMENTATION CHECK LIST	II-1
	Update from V. 4.5 to V. 5.0	
	Virgin Installation (no previous version installed or implemented)	II-10
III.	SCREENING METHODS	III-1
IV.	SECURITY KEYS	IV-1
	RA ALLOC	IV-1
	RA MGR	IV-1
	RA UNVERIFY	IV-1
	RA VERIFY	IV-1
V.	System Definition Menu	V-1
	Division Parameter Set-up	V-2
	Division-wide Order Entry Parameters	V-2
	Exam Entry/Edit Parameters	V-3
	Films Reporting Parameters	V-4
	Miscellaneous Division Parameters	V-7
	HL7 Applications Associated with this Division	V-8
	Imaging Locations Associated with this Division	
	Print Division Parameter List	V-15
	Camera/Equip/Rm Entry/Edit	V-19
	List of Cameras/Equip/Rms	
	Location Parameter Set-up	V-22
	Flash Card, Jacket Labels, and Exam Label Param	
	Order Entry Parameters	
	Report Parameters	
	Cameras/Equip/Rooms Used by this Location	
	Default CPT Modifiers used by this Location	V-25
	Recipients of the 'Stat' Alert for this Location	
	Imaging Type for this Location	
	Location Parameter List	V-34
VI.	RAD/NUC MED PERSONNEL MENU	VI-1
	Classification Enter/Edit	VI-2
	Clerical List	VI-6

i

	Interpreting Resident List	VI-6
	Interpreting Staff List	VI-6
	Technologist List	VI-6
VII.	UTILITY FILES MAINTENANCE MENU	VII-1
	Complication Type Entry/Edit	VII-2
	Diagnostic Code Enter/Edit	VII-4
	Examination Status Entry/Edit	
	Exported Setup	
	Example Status Tracking Sessions	
	Film Type Entry/Edit	
	Label/Header/Footer Formatter	
	Major AMIS Code Entry/Edit	VII-45
	☞Nuclear Medicine Setup Menu	VII-47
	■ Lot (Radiopharmaceutical) Number Enter/Edit	VII-48
	Route of Administration Enter/Edit	VII-49
	Site of Administration Enter/Edit	VII-50
	■ Vendor/Source (Radiopharmaceutical) Enter/Edit	VII-51
	Order Entry Procedure Display Menu	VII-53
	Common Procedure Enter/Edit	VII-54
	Create OE/RR Protocol from Common Procedure	VII-60
	Display Common Procedure List	
	Procedure Edit Menu	VII-62
	■ Cost of Procedure Enter/Edit	VII-63
	Procedure Enter/Edit	VII-64
	Procedure Message Entry/Edit	VII-79
	Procedure Modifier Entry	VII-80
	Reports Distribution Edit	
	Sharing Agreement/Contract Entry/Edit	
	Standard Reports Entry/Edit	
	Valid Imaging Stop Codes Edit	
	VistaRad Category Entry/Edit	VII-89
VIII.	MAINTENANCE FILES PRINT MENU	VIII-1
	Complication Type List	VIII-2
	Diagnostic Code List	VIII-3
	Examination Status List	
	Film Sizes List	VIII-
	Label/Header/Footer Format List	VIII-
	11 M: AMEC 1 1:4	* 7***
	Major AMIS Code List	VIII-

Modifier List	VIII-
™Nuclear Medicine List Menu	VIII-
14	
■ Lot (Radiopharmaceutical) Number List	VIII-
15	
■Route of Administration List	VIII-
16	
Site of Administration List	VIII-
17	
■ Vendor/Source (Radiopharmaceutical) List	VIII-
18	
Procedure File Listings	VIII-
■ HIST Display Rad/NM Procedure Contrast Media History	VIII-
20	
Sample Contrast Media Edit History by Procedure Report	VIII-
Active Procedure List (Long)	VIII-
24 Active Procedure List (Short)	VIII-
25	V 111
Alpha Listing of Active Procedures	VIII-
■Barcoded Procedure List	VIII-
27	
Inactive Procedure List (Long)	VIII-
30 Invalid CPT/Stop Code List	VIII-
31	
List of Inactive Procedures (Short)	VIII-
EList of Procedures with Contrast	VIII-
34	
■ Parent Procedure List	VIII-
38	
Procedure Message List	VIII-
Series of Procedures List	VIII-
Report Distribution Lists	VIII-

	Sharing Agreement/Contract List	VIII-
	Standard Reports Print	VIII-
	44	, , , , , , , , , , , , , , , , , , , ,
	VistaRad Category Print	VIII-
IX.	RADIOLOGY HL7 MENU	IX-1
	Rad/Nuc Med HL7 Voice Reporting Errors	IX-2
	Resend Radiology HL7 Message	
	Send ADT HL7 message to PACS	IX-7
X.	Worksheets	X-1
	Division Parameters	X-1
	Camera/Equip/Room	X-3
	Location Parameters	X-4
	Classification	X-7
	Blank Examination Status Worksheet	X-8
	Procedure Setup	X-12
	■ Radiopharmaceuticals	X-13
	Example Radiopharmaceutical Drug Classifications	X-15
	™ Nuclear Medicine Setup	X-17
XI.	₽PARENT/DESCENDENT EXAM AND PRINTSETS	XI-1
	Definition	
	Common Data for Printsets	
	Non-Common Data for Printsets	
	Single Report Setup for Printsets	
	Special Character Designation for Printsets	
	HL7 Interface	
	Registration	
	Editing Data	
	Labels, Headers, Footers	
	VA Fileman Prints/Inquiries	
	Add Exams to Last visit	XI-5
	CPRS Interface for Printsets	XI-5
	Imaging/Multi-Media Interface Effects on Printsets	XI-5
XII.	CASE EDITS AND STATUS TRACKING	XII-1
	Overview	
	Taking Control	
	Stuffing Data and Default Responses	
	Radiopharmaceutical Logic	XII-2
	■ Medication Logic	XII-3

V V / / / / / / / / / / / / / / / / / / /	Typery	1
XVII. 1	GLOSSARY	XVII-
XVI.	₽ STOP TASK (TASKMAN)	XVI-1
	General Editing Helpful Hints	XV-14
	Research	
	Contracting Out	
	Complications	
	Effect of Portable and Operating Room Modifiers	
	Doubling Workload Credit	
	Hold/Cancelled/Pending/Scheduled	
	Verifying Reports	
	Purging Records	
	Printing/Batching/Etc	
	Requesting/Scheduling Requests	
	Case Numbers	
	Procedures	
XV.	FREQUENTLY ASKED QUESTIONS Imaging Types	
	•	
	Reports	
AIV.	TROUBLESHOOTING	
XIV.	Though echooting	VIV 1
	Errors	
	Crediting	
	Options that can Potentially Trigger Crediting	
	Location Parameter Edit	
	Invalid CPT/Stop Code List	
	Imaging Stop Codes File 71.5	
XIII.	OUTPATIENT WORKLOAD CREDITING	
	Population of Imaging Exam Data	
	Radiopharm Prescription & Prescribing Physician/Dose Logic	XII-4

I. Introduction

As application coordinator, you are required to build, customize, and maintain the Radiology/Nuclear Medicine module files to meet the specific needs of your site. This document is provided to assist you in this function.

This icon designates something as being new or changed in some way from the previous version.

Implementation Checklist

We suggest using the checklist to make sure you do not miss any necessary steps in implementation. If your site implemented version 4.5, use the Update V. 4.5 to 5.0 Installation Checklist. If this is a first time installation or setup of this program, use the Virgin Installation Checklist.

Supervisor Menu

Options are shown using their external name (the name users see) with the internal name in brackets. This is an example of what you might see:

Division Parameter Set-up [RA SYSDIV]

There are four main menus used for implementation and maintenance of this software. The IRM Menu which Information Resource Management handles also plays a roll in implementation, but it is addressed in the technical manual. Each option in the four menus is explained in the following chapters. They are exercised during implementation and throughout the use of the software to maintain or update the files.

- 1. The first menu is the System Definition Menu. It contains options to set up and view Division and Imaging Location parameters and to add new cameras, equipment, or rooms to the system.
- 2. The second is the Rad/Nuc Med Personnel Menu in which you can enter or edit the classification and access of clerks, technologists, residents, and staff. It also contains print options to obtain a listing of each of these categories.

3. The third is the Utility Files Maintenance Menu which contains options to add and edit:

Complication types

Diagnostic codes

Exam statuses activation and parameters

Film types

Label/header/footer formats

Major AMIS codes (use only if authorized by VA HEADQUARTERS)

Nuclear Medicine parameters

Order entry procedures

Procedures and their characteristics

Procedure messages

Procedure modifiers

Reasons for cancellation


Report distribution queues

Sharing agreements and contracts

Standard reports

Valid imaging stop codes

4. The fourth, the Maintenance Files Print Menu is used in conjunction with these add/edit options to display the contents of the files.

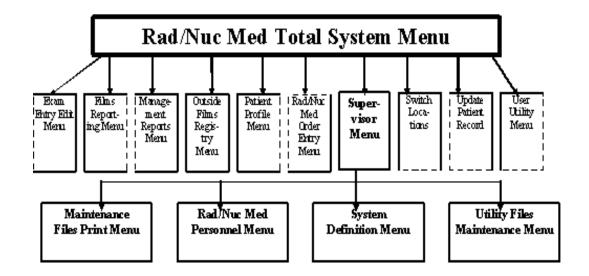

Complication Type List Classification Enter/Edit Camera/Equip/Rm Entry/Edit **Complication Type** Entry/Edit Diagnostic Code List **Clerical List Division Parameter** Diagnostic Code Enter/Edit Set-up **Examination Status List Interpreting Resident List** List of Cameras/ **Examination Status** Entry/Edit Equip/Rms Film Sizes List **Interpreting Staff List Location Parameter List** Film Type Entry/Edit Label/Header/ **Technologist List Location Parameter** Label/Header/ **Footer Format List** Set-up **Footer Formatter** Major AMIS Code List **Print Division Parameter List** Major AMIS Code Entry/Edit Nuclear Medicine Setup Menu **Modifier List** Nuclear Medicine List Menu **Order Entry Procedure Display Menu Procedure File Listings Procedure Edit Menu Report Distribution Lists Sharing Agreement Reports Distribution Edit** /Contract List Standard Reports Print **Sharing Agreement** /Contract Entry/Edit Standard reports Entry/Edit **Valid Imaging Stop Codes**

June 2019

Division, Location, Imaging Type Definitions

There are three terms you must be familiar with in order to understand how this package functions. They are Division, Imaging Type, and Imaging Location.

A **Division** is defined as a major unit of an institution which requires separate management. In most instances, there is only one Division at a given medical center which is also referred to as the site. Integrated sites usually require at least one Division defined per site. The Division Parameter Set-up and Location Parameter Set-up options allow you to tailor the package to fit the particular needs of your department. All parameters are specific to the Division, which allows you, as the coordinator, to select the parameters you want to utilize with either a YES or NO answer.

An **Imaging Location** is defined as a physical location where a procedure can be performed. Each location has its own set of parameters defining location-specific criteria, e.g., label/header/footer formats, printer devices for flash cards, jacket labels requests and reports and how many labels to print. An Imaging Location may be associated with **only one** Division and **only one** Imaging Type.

An **Imaging Type** must be assigned to each procedure and each Imaging Location. Therefore, many Imaging Types can be active within a single Division. Imaging Types other than General Radiology and Nuclear Medicine should only be "activated" if they indeed are Imaging Types used at your facility and report procedures and/or workload separately from the other Imaging Types. In other words, if you want CT scans to be on separate workload reports rather than on General Radiology workload reports, then you should activate the CT Scan Imaging Type (by associating locations that do CT Scans with the Imaging Type of CT scan and assigning CT scan procedures an Imaging Type of CT Scan). Imaging Types are limited in this package to the following:

ABBREVIATION

ANI
CARD
CT
RAD
MRI
NM
US
VAS
MAM


IMAGING TYPE

You should **not** attempt to add any additional Imaging Types to the above list.

It is important to know how to activate an Imaging Type. When at least one Imaging Location and one or more procedures have been assigned a new Imaging Type, that Imaging Type becomes "active".

Below is a schematic showing relationships:

Division, Imaging Location, Imaging Type, and Camera/Equipment/Room Relationships

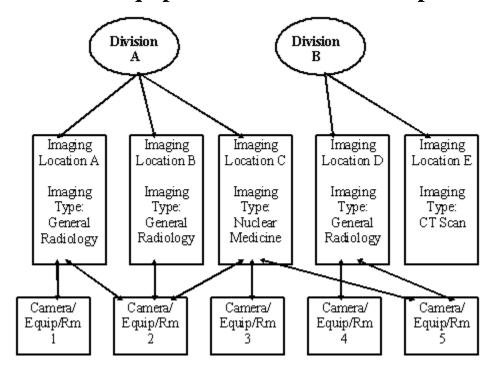

Note: The Imaging Type of any given procedure will determine in which Imaging Locations it may be registered and performed. Example: If the Imaging Type for procedure X is General Radiology, it can be registered in Imaging Locations A, B, and D.

Parent/Descendent Exam and Printsets

A major enhancement to this version is a single combined report for a set of related procedures. This is a "printset" mechanism for entering a single report for all descendent cases registered from a parent order (See Procedure Enter/Edit for more information on parent procedures). For more detailed information see Parent/Descendent Exams and Printsets. The ability to create a separate report for each procedure ordered under a single parent procedure still exists.

Radiopharmaceuticals

Another important addition with this version is the ability to enter and edit information specific to radiopharmaceuticals for Nuclear Medicine. A new menu, Nuclear Medicine Setup Menu under the Utility Files Maintenance Menu, allows the site to define parameters for radiopharmaceuticals concerning lot number, route and site of administration, and source/vendor. The addition of radiopharmaceutical fields has a major effect on case and status edits for Nuclear Medicine and Cardiology Studies Imaging Types. For more complete information, refer to the chapter Case Edits and Status Tracking.

II. Implementation Check List

Update from V. 4.5 to V. 5.0

This implementation checklist assumes that you have fully implemented version 4.5 and now want to review for implementation the new functionality available in this version.

Before IRM loads this version into production:

- a. Review the Release Notes to acquaint yourself with any changes that were made to this version of the software and update the users to those changes.
- b. Review all new items in this manual (those designated by) and read the chapters on Single Reports for Parent/Descendent Exams and Printsets and Case Edits and Status Tracking. Also, Procedure Enter/Edit in this manual and Report Entry/Edit and Diagnostic Code and Interpreter Edit by Case No. in the user manual are affected by printsets.
- c. Radiology/Nuclear Medicine accesses the Pharmacy Drug file #50 for radiopharmaceuticals and their associated VA Classifications; therefore, the Drug file must contain all the radiopharmaceuticals used at your site. Provide the Pharmacy ADPAC with a list of the radiopharmaceuticals and work with the ADPAC to enter them into the Drug file. When entering Radiopharmaceuticals into the Drug file, it is recommended to enter the Nuclide then the Pharmaceutical as the Drug Name and also a synonym to make it easier for Radiopharmaceuticals to be retrieved. (See Example Radiopharmaceutical Drug Classifications under the chapter Worksheets.)

The Radiopharmaceuticals should be assigned one of the following VA Classifications under Diagnostic Agents:

DX200 Radiopharmaceuticals, Diagnostic DX201 Imaging Agents (in vivo), Radiopharmaceuticals DX202 Non-Imaging Agents Radiopharmaceuticals

- d. If a status's Order is removed through the Examination Status Entry/Edit option but the Appear on Status Tracking field is set to YES, exams in that status will, after the installation of version 5.0, be available for edit in Status Tracking of Exams. The inactivated status will be processed after the active statuses. The purpose of this feature is to allow processing of exams to complete if they are in an inactive status and the case edit options do not prompt for the required fields to allow data entry to complete the cases. If your site has a large volume of historic cases that fall into this category, you may want to make sure the Appear on Status Tracking field is set to NO prior to the installation of this version.
- e. Print the Active Procedure List (Short) to get a list of Nuclear Medicine procedures that can be used to fill out your worksheets for procedure set-up.
- f. It will be necessary to train the technologists and transcriptionists in handling case and report entry/edits for printsets. This should be done in your training account prior to installing live.

The technologists must be made aware of the following:

- 1) If an exam is added through the Add Exam to Last Visit option, and the last visit was a printset/parent, the added exam becomes part of the printset. The report will be a single one even if the added exam is quite different from the substance of the parent.
- 2) If there are two parent requests of the same Imaging Type, do not use the Add Exam to Last Visit option to add the procedures of the second parent to the first. Doing this will create a single report for all procedures from both parents, and the second request will never be completed because all descendants will apply to completing the first request.

While you do the setup required after installation, Imaging Service users must be off the system. You should prepare as much as possible to minimize the downtime for Imaging Service users.

The following is a listing and discussion of all the new functionality that you may want to edit in setting up this version. It is strongly suggested that you review the following checklist and complete any worksheets **prior to the actual setup**. Then use the completed checklist and worksheets during the setup.

than the	RM DOSE WARNIN default message	G MESSAGE - If	vou want to o	
than the		G MESSAGE - If	you want to o	
in the I	ser when the radioph Pose" and the "High A	armaceutical dose Adult Dose", enter ceuticals multiple	ed directiv falls outside the message in the Rad/Nu	e by a physician. the range of the "Low here. These fields re c Med Procedures file.
	e:			

a. ALLOW 'RELEASED/NOT VERIFIED' - If you want to be able to display unverified reports to users outside the Imaging service, this Imaging Location parameter should be set to YES. This field was moved from the Division parameters. If as a Division parameter it contained a YES, then all the Imaging Locations for that Division will contain a YES in this field as a default after this version is installed. Otherwise, the default will be no.

Location	Send Alert	Allow Released/ Unverified

b. STAT REQUEST ALERT RECIPIENTS - List any Radiology/Nuclear Medicine personnel that should receive an alert whenever a STAT request is submitted to the Imaging Location. Bypass this step if you don't want alerts sent.

Note: To receive a STAT request alert, recipients must have a Radiology/Nuclear Medicine personnel classification. To use this feature, the Division parameter, "Ask Imaging Location" must be set to YES, and CPRS (OE/RR V. 3.0 or higher) must be installed.

Locations:	Recipients:	
_		

- c. URGENT REQUEST ALERTS? If you want to send alerts to selected personnel whenever there is a STAT request, enter YES in this field. Personnel designated to receive STAT request alerts will also receive URGENT alerts. Note: To use this feature, the Division parameter "Ask Imaging Location" must be set to YES, and CPRS (OE/RR V. 3) or higher must be installed.
- 3) Diagnostic Code Enter/Edit.

INACTIVE - Print out a list of your diagnostic codes using the option Diagnostic Code List. If you want to inactivate any of the codes, so they do not appear as selections to the users, highlight them and refer to that list when doing this step of the setup.

4) Examination Status Entry/Edit.

Editing this option takes a great deal of up front preparation. Make sure you read the information in the chapters Examination Status Entry/Edit and Case Edits and Status Tracking. Use the worksheet (See chapter Worksheets) when defining your changes.

A new Imaging Type, Mammography, was exported with this version. If you intend to activate Mammography as an Imaging Type, then be sure to include it when defining your changes to status tracking.

Note: This installation will **turn on** the "Ask Radiopharmaceuticals and Dosages" parameter for all nuclear medicine and cardiology studies procedures. All other new parameters will be off. In order for nuclear medicine technologists to accurately enter radiopharmaceuticals, all radiopharmaceuticals used at your site must be entered into the Pharmacy Drug file.

The new fields to edit for all Imaging Types are as follows:

- a. GENERATE EXAMINED HL7 MESSAGE
- b. ASK MEDICATIONS & DOSAGES
- c. ASK MEDICATION ADMIN DT/TIME & PERSON

There are a number of new fields for Nuclear Medicine. If your facility wants to record data in most or all of the new radiopharmaceutical data fields, you will have to use Status Tracking to edit cases. Only through Status Tacking will you be able to control which prompts are asked at which statuses. In order to provide simplicity and avoid extra undesirable prompts at facilities who want to enter a minimum set of radiopharmaceutical data, the case editing options (Case No. Exam Edit and Edit Exam by Patient) do not prompt for all the new Nuclear Medicine fields. Instead, only the radiopharmaceutical, dose administered, date/time of dose administration, and person who administered dose will be consistently asked for nuclear medicine or cardiology studies cases. However, if data has already been entered via Status Tracking in other nuclear medicine related fields, the case edits will prompt for those fields to allow correction of erroneous radiopharmaceutical data already entered. In addition, if the procedure parameter so specifies, the case edits will prompt for radiopharmaceutical dose prescribed, prescribing physician, and witness to dose administration. The new fields to edit for Nuclear Medicine and Cardiology Studies Imaging Types are as follows:

- a. RADIOPHARMS & DOSAGES REQUIRED
- b. ACTIVITY DRAWN REQUIRED (RADIOPHARM)
- c. DRAWN DT/TIME AND PERSON REQUIRED (RADIOPHARM)
- d. ADMIN DT/TIME/PERSON REQUIRED (RADIOPHARM)
- e. ROUTE/SITE REQUIRED (RADIOPHARM ADMINISTERED)
- f. LOT NO. REQUIRED (RADIOPHARM)
- g. VOLUME/FORM REQUIRED (RADIOPHARM)
- h. ASK RADIOPHARMACEUTICALS AND DOSAGES
- i. ASK ACTIVITY DRAWN (RADIOPHARMACEUTICAL)
- j. ASK DRAWN DT/TIME & PERSON (RADIOPHARMACEUTICAL)
- k. ASK ADMIN DT/TIME & PERSON (RADIOPHARMACEUTICAL)
- 1. ASK ROUTE/SITE OF ADMIN (RADIOPHARM)
- m. ASK LOT NO. (RADIOPHARM)
- n. ASK VOLUME/FORM (RADIOPHARM)
- o. PRINT DOSAGE TICKET WHEN EXAM REACHES THIS STATUS

 Note: Printing of the dosage ticket requires that a valid print device assigned to
 the location (see IRM Menu, Device Specifications for Imaging Locations option,
 Dosage Ticket Printer prompt) and YES answered for the Radiopharms/Dosages
 Required field at the status you are editing (or lower).

Even if you do not make changes to the exam status parameters, select the option for each Imaging Type. This will trigger the Data Inconsistency Report and identify all problems. You can up-arrow at the Select an Examination Status prompt as shown here:

Select an Imaging Type: NUCLEAR MEDICINE

Select an Examination Status: ^	
Checking order numbers used for status progression	
within : NUCLEAR MEDICINE	
Required order numbers are in place.	
Data Inconsistency Report For Exam Statuses	
DEVICE: HOME// (Enter a printer name)	

6) Label/Header/Footer Formatter. There are a number of new fields available that could be added to labels, headers or footers:

Long Case Number Barcode Resident Signature Name Staff Signature Name Verifying Signature Name Verifying Signature Title Patient Address Line 1 Patient Address Line 2 Patient Address Line 3

- 7) Nuclear Medicine Setup Menu contains four options. Routes and sites of administration of radiopharmaceuticals were exported with this version. Review the data in those files, using the print options Route of Administration List and Site of Administration List, before adding more entries. Worksheets for all four options are available in the chapter Worksheets. Add sites and vendors/sources before adding routes and lot numbers. Sites should be entered before routes so they will be selectable as valid sites for a route and lot numbers can be associated with vendors/sources. Also, routes and sites must be available for the next step, Procedure Enter/Edit. You may also want to request that IRM add the Lot Number Enter/Edit as a secondary menu option for the Nuclear Medicine technologists.
- 8) Procedure Enter/Edit. This is another complex setup that requires considerable up front preparation. If you plan to use the nuclear medicine features, you will need to enter radiopharmaceutical parameters on the appropriate procedures. Use the worksheet, making as many copies as necessary. The new fields to edit for all Imaging Types are as follows:
 - a. SINGLE REPORT
 - b. PROMPT FOR MEDS
 - c. DEFAULT MEDICATION
 - d. EDUCATIONAL DESCRIPTION
 - e. DISPLAY EDUCATIONAL DESCRIPTION

Note: If you turn parents into printsets by answering YES to the "Single Report" prompt, the printset features will only take effect on cases registered **after** V. 5.0 is installed.

The new fields to edit for Nuclear Medicine and Cardiology Studies Procedures are as follows:

- a. SUPPRESS RADIOPHARM PROMPT
- b. PROMPT FOR RADIOPHARM RX
- c. DEFAULT RADIOPHARMACEUTICAL

LOW DOSE HIGH DOSE USUAL DOSE ROUTE FORM

If your site uses Parent procedures and you want a single report for all descendents, you can use the Print File Entries option to obtain a list of Parent procedures. You need not print the descendents as shown in this example. Or, you can use one of the existing procedure listing options to get complete data for all active procedures.

```
Select Supervisor Menu Option: Print File Entries

OUTPUT FROM WHAT FILE: (file name)// RAD/NUC MED PROCEDURES

SORT BY: NAME// TYPE OF PROCEDURE
```

- 9) If you activate the new Imaging Type Mammography, be sure to enter any modifiers you want for that Imaging type using the option Procedure Modifier Entry.
- 10) If you want Request Cancellations to print, give your IRM support staff a list showing which printer for each imaging location should be used to print them. They will enter these via their IRM menu using the Device Specifications for Imaging Location option.
- 11) This release is year 2000 compliant. The Date of Birth and Current Patient Location (which is date/time stamped) fields for labels, headers and footers are affected because the year is now printed as four digits instead of two. Retest any labels or report headers/footers that contain these fields to make sure they don't line wrap.

Virgin Installation (no previous version installed or implemented)

1)	☐ Before IRM loads this system into production, you should read this manual, read over the
	installation instructions listed here, and complete any related worksheets. The software should be installed in a Training environment first so you can familiarize yourself with the setup procedure and train users as much as possible before going live.
2)	☐ Before you begin the process of building your files, you should consult with your IRM
	Service to be sure you have access to the security keys (RA MGR, RA ALLOC and RA VERIFY) and have been assigned the Rad/Nuc Med Total System Menu.
3)	☐ Provide the Pharmacy ADPAC with a list of the radiopharmaceuticals and work with the
	ADPAC to enter them into the Drug file. See an example of Radiopharmaceutical entries in the Drug file under the chapter Worksheets. They should be assigned one of the following VA Classifications under Diagnostic Agents:
	(a) DX200 Radiopharmaceuticals, Diagnostic(b) DX201 Imaging Agents (in vivo), Radiopharmaceuticals(c) DX202 Non-Imaging Agents Radiopharmaceuticals
4)	☐ Edit parameters for each Division at your facility using the option, Division Parameter
	Set-up. (See chapter Worksheets)
5)	☐ Edit parameters for each Imaging Location by using the option Location Parameter Set-
	up. (See chapter Worksheets)
6)	☐ Use the option Examination Status Entry/Edit for each Imaging Type that you activate at
	your facility. Read about examination status tracking first, including the information in chapter Case Edits and Status Tracking so that you understand the implications of activating an Imaging Type. You should have the exam status worksheets filled out ahead of time so that this part of the setup goes smoothly.

7)	Review the rest of the options in the System Definition Menu and the Utility Files
	Maintenance Menu and make any changes that are appropriate for your site. Look for worksheets at the back of this manual.
8)	☐ Use the print options in the System Definition Menu, Rad/Nuc Med Personnel Menu, and
	the Maintenance Files Print Menu to obtain dated hard copies of the data for your records.
9)	☐ Determine the menus and/or options you need to assign to each of your users. There are
	a number of work related menus (such as clerk, file room, radiologist, ward, technologist, etc.) exported with this package. See the technical manual for more information. Give the list to IRM. You may want to review the Mailman bulletins provided by the Radiology/Nuclear Medicine software with IRM so mail groups with the appropriate members can be set up.
10)	This item is extremely important and should be done as soon as the above steps
	are complete. No user will be able to sign-on to the Radiology/Nuclear Medicine package unless s/he is assigned at least one Imaging Location. Use the option Classification Enter/Edit in the Rad/Nuc Med Personnel Menu now and give each user access to whatever Imaging Locations they work with. It may be desirable to assign the RA ALLOC key to transcriptionists if they need to enter reports for all Imaging Locations.
11)	Continue training your users.

III. Screening Methods

Most options are screened by a combination of Imaging Type, Division, and Imaging Location. Others are screened by "ownership". Division, Imaging Type, and Location screening methods are based on the following:

- 1. The set of locations the user is privileged to access (entered by the ADPAC or IRM through the Classification Enter/Edit option). Anyone who owns the RA ALLOC key can access all Imaging Locations so it isn't necessary to assign Imaging Locations to owners of this key. (Appropriate owners are ADPACs and transcriptionists.)
- 2. The sign-on location, which must be a location to which the user has access.
- 3. The user's Imaging Type access "set" and Division access "set" are calculated by the system, based on the Imaging Locations to which the user has access. They include all the Imaging Types and Divisions of all the locations in the user's access set. For example, if a user only has access to the General Radiology "X-RAY 1" Imaging Location which is in Division A, access will be calculated by the system as follows:
 - a. If the user selects an option that is screened by location, s/he will only be able to access exam data within X-RAY 1.
 - b. If the user selects an option that is screened by Imaging Type, s/he will be able to access data that is within General Radiology regardless of location or Division. This means the user could access more than one location if there are multiple locations of that Imaging Type.
 - c. If the user selects an option that is screened by Imaging Type and Division, s/he will be able to access data within Division A and General Radiology Imaging Type regardless of location.
 - d. If the user selects a non-screened option or an option screened by "ownership" of a report, batch, etc., the location/Imaging Type/Division screening will have no effect.

So far we have the following variations of screening:

Method No. and Description

- 1. Ownership (usually of reports or report batches).
- 2. Some variation of single location screening Used for the option Register Patient for Exams because the Imaging Type of the procedure selected must match the Imaging Type of the sign-on location. Used for Add Exams to last Visit because the sign-on location must match the location of the exams previously entered for the visit.
- 3. Selected location(s) Locations selected must be in the user's location access set, and must have an Imaging Type that matches the Imaging Type of the sign-on location. This method is used for Exam Status Display/Tracking options, and will be used for some of the reports that have to do with Orders/Requests. (Exam status options must work with one Imaging Type at a time since the exam statuses are now imaging-type specific.) Default response at location selection prompt is "ALL", which means all locations to which the user has access.
- 4. Sign-on Imaging Type User is only allowed to select patients and/or cases with an Imaging Type that matches the Imaging Type of the sign-on location. This is used for the remainder of the case edit/view/cancel options.
- 5. Selected Imaging Types If user can only access one Imaging Type, the system will default to that Imaging Type without prompting. If the user can access multiple Imaging Types, the user is prompted to select one or more Imaging Types. Default response is "ALL", which means all Imaging Types to which the user has access.
- 6. Division and Imaging Type selection User first selects one or more Divisions from his/her Division access list (if only one Division is accessible, the system will default to it rather than prompt the user). Next, user selects one or more Imaging Types from his/her Imaging Type access list (which may now be further restricted, since they may have chosen a Division in which not all Imaging Types to which they have access are active). This method was used for Daily Report, Complication Report, Abnormal Report so far, and will be used for most workload and management reports on exams that have at least progressed to the "registration" stage. It is also used for transcriptionist options. Default responses to the Division and Imaging Type selection prompts are "ALL", which means all to which the user has been given access privileges.
- 7. Menu option allocation This is always at the discretion of the site.
- 8. Division selection.

Here's a list of options and their screening methods:

~	•
Cre	ening
\mathcal{L}	CIIIII

Menu Option Method Key Needed

Access Erroneous Reports None (On Supervisor Menu only)

Cancel a Request None

Create a Batch None

Delete a Report None RA MGR
Delete Printed Batches by Date None RA MGR

Display Rad/Nuc Med Procedure Information None

Display a Rad/Nuc Med Report None
Display Patient Demographics None
Distribution Queue Menu (all options) None

Enter Last Past Visit Before VistA None

Exam Deletion None RA MGR

Exam Profile (selected sort) None

Hold a Request None

Inquire to File Entries (FileMan) N/A (On Supervisor Menu only)

Jacket Labels None

List Exams with Inactive/Invalid Statuses None (On Supervisor Menu only)

List Reports in a Batch None

Maintenance Files Print Menu (all options)

None (On Supervisor Menu only.

If appropriate, printout will

show Imaging Type)

Nuclear Medicine List Menu (all options)

None
Nuclear Medicine Setup Menu (all options)

Outside Films Profile

Outside Films Registry Menu (all options)

None
Pending/Hold Rad/Nuc Med Request Log

None

Print a Batch of Reports

None

Print File Entries (FileMan) N/A (On Supervisor Menu only)

Print Rad/Nuc Med Requests by Date
Print Selected Requests by Patient
Print Worksheets
None
Profile of Rad/Nuc Med Exams
None

Rad/Nuc Med Personnel Menu (all options) None (On Supervisor Menu only)

Request an Exam None

Schedule a Request None

Search File Entries (FileMan) N/A (On Supervisor Menu only)

Select Report to Print by Patient None

Menu Option	Scree:	•	Key Needed
System Definition Menu (all options)	None	_	isor Menu only. ate, edits will be g Type)
Test Label Printer	None		
Transcription Report	None		
Unverify a Report for Amendment	None		RA UNVERIFY
Update Patient Record	None		
Verify Batch	None		RA VERIFY
View Exam by Case No	None		
Ward/Clinic Scheduled Request Log	None		
Utility Files Maintenance Menu (all options)	None	(On Superv	isor Menu only)
Add/Remove Report from Batch	1		
Delete Printed Batches	1		
On-line Verifying of Reports	1		RA VERIFY
Resident On-line Pre-Verification	1		RA VERIFY
Add Exams to Last Visit	2		
Register Patient for Exams	2		
Detailed Request Display	3		
Exam Status Display	3		
Log of Scheduled Requests by Procedure	3		
Status Tracking of Exams	3		
Switch Locations	None		
Cancel an Exam	4		
Case No. Exam Edit	4		
Diagnostic Code and Interpreter Edit			
by Case No.	4		
Duplicate Flash Card	4		
Draft Report (Reprint)	4		
Edit Exam by Patient	4		
Indicate No Purging of an Exam/Report	4		
Mass Override Exam Status	4		RA MGR
Override a Single Exam Status to 'complete'	4		RA MGR
Update Exam Status	4 (O	n Supervisor	Menu only)

Menu Option	Screening <u>Method</u>	Key Needed
Abnormal Exam Report	6	
Complication Report	6	
Daily Log Report	6	
Delinquent Status Report	6	
Examination Statistics	6	
Functional Area Workload Reports (all rots)	6	
Incomplete Exam Report	6	
Personnel Workload Reports (all options	6	
except the Transcription Report)		
Report Entry/Edit	6	
Special Reports (all except AMIS reports)	6	
Unverified Reports	6	
Verify Report Only	6	RA VERIFY
AMIS Code Dump by Petient	0	
AMIS Code Dump by Patient	8	
AMIS Report	8	

Note: The RA ALLOC key gives a user access to the broadest possible set of locations, Imaging Types, and Divisions possible with any given option. It is the equivalent of having access to all Imaging Locations.

IV. Security Keys

RA ALLOC

Ownership of the RA ALLOC key overrides the location access security given to personnel through Classification Enter/Edit. Owners of the RA ALLOC key have expanded access to Imaging Locations, Imaging Types, and Divisions. In the case of most workload reports, this means they can select from a list of all Divisions and Imaging Types to include on the report. In the case of various edit and ordering functions, it means they can select from all locations within the Imaging Type to which they are currently signed on through the "Select sign-on location" prompt.

RA MGR

The RA MGR key gives the user access to supervisor-type functions. Those functions are the following:

- 1. editing completed exams
- 2. adding an exam to a visit that is older than yesterday
- 3. showing the user all non-completed exams, not just those associated with the user's current Division, during execution of the "status tracking" function
- 4. updating the exam status of an exam to complete
- 5. deleting exams
- 6. deleting reports
- 7. entering a report on a cancelled exam if site parameters don't allow it
- 8. deleting printed batches by date

RA UNVERIFY

The RA UNVERIFY key allow users to unverify a report for amendment.

RA VERIFY

The RA VERIFY key allows users to verify reports.

V.System Definition Menu

Supervisor Menu ... System Definition Menu ...

> Camera/Equip/Rm Entry/Edit Division Parameter Set-up List of Cameras/Equip/Rms Location Parameter List Location Parameter Set-up Print Division Parameter List

Defining and Maintaining Your System Parameters

In this chapter, we will show you how to define your system parameters using the above options, Division Parameter Set-up, Location Parameter Set-up, and Camera/Equip/Rm Entry/Edit. Once you have defined and entered or edited your system parameters, you should use the print options in the System Definition Menu to provide a dated hard copy of the parameters.

New Divisions, locations or cameras/equipment/rooms may be added or their parameters edited at any time throughout the use of the software. These options are used for file maintenance as well as implementation. If this is not a new implementation of the software, you should take this time to review your system parameters for any necessary updating of the files. In particular,

check those items that appear with the ricon.

Division Parameter Set-up [RA SYSDIV]

The Division Parameter Set-up option is used to define Division unique information. The Radiology/Nuclear Medicine module was designed to handle the functions of multiple Divisions within a medical center even though most medical centers have only one Division. As a result, the program is structured around the Division concept. Therefore, **even if you are a single Division medical center you must enter at least one Division** using this option.

All Imaging Types and locations associated with this Division are affected by the parameters in this option. For instance, if a clinical history message is entered (number 3 below), the message will appear for any type of procedure ordered whether for general radiology, nuclear medicine, ultrasound, etc. If you require the impression on the report, it will be required on reports from all Imaging Types.

An existing Division may be edited or a new Division may be entered. All Divisions must be in the Institution file (# 4). The following is a listing of Division parameters and a brief explanation for each.

Division-wide Order Entry Parameters

1. ASK 'IMAGING LOCATION' - This field determines whether or not the "Submit Request To" prompt appears when a procedure is requested. If set to YES, the user will be **required** to enter the name of the Imaging Location where the request is to be sent and printed at the time they request an exam. During the requesting process, the division of the requesting location (i.e., patient location) is used to determine whether to prompt for a "Submit Request To" location. If the requesting location is not assigned a division in the Hospital Location file #44, the primary division will be used. The user will only be allowed to choose from Imaging Locations whose Imaging Type matches the Imaging Type of the procedure. It is strongly recommended that you answer "YES" to this question.

If you do not answer this prompt, the system default is NO and the prompt, "Submit Request To:" will not appear. Setting it to NO will make order tracking more difficult, and is not recommended. If it is set to YES, the Division will be able to implement alerts on STAT and URGENT requests once CPRS (OE/RR V. 3) is installed.

¹If it is set to NO or not answered at all, and there is only one imaging location with the same imaging type as the ordered procedure's imaging type, then this lone imaging location will be assigned to the order.

If it is set to NO or not answered at all, and there are more than one imaging location with the same imaging type as the ordered procedure's imaging type, then the order will not

¹ Patch RA*5*34 Inserted paragraphs regarding parameter being set to NO or not answered at all.

have an imaging location and will show up in the UNKNOWN section of the 'Detailed Request Display' [RA ORDER DISPLAY] option.

2. TRACK REQUEST STATUS CHANGES - If set to YES, this parameter establishes the collection of the date and time of request status changes. Request status changes pertain solely to the physician's order and whether it is in one of the following statuses: discontinued, complete, hold, pending, active, scheduled, unreleased. Request status changes can be used to find the status of an order at any given time in the past.

This prompt is now bypassed and automatically set to YES.

3. CLINICAL HISTORY MESSAGE - This is an **optional**, free text message, 5-70 characters in length that appears just prior to the prompt for the clinical history field. This is often used to remind the person entering the request that specific information is needed in the patient's clinical history.

If you bypass this field, no message will appear prior to the user seeing the prompt for the clinical history.

Exam Entry/Edit Parameters

4. DETAILED PROCEDURE REQUIRED - If set to YES, entry of a Parent, Series or Detailed procedure type, will be required. If set to NO, procedures with the procedure type Broad will be accepted during order entry as a valid entry. This parameter in the division of the requesting location is used to determine if a broad procedure can be ordered. At the time the procedure is registered, broad types won't be allowed if either the division of the imaging location or the division of the requesting location requires a detailed procedure. If no division is entered on the requesting (or imaging) location in the Hospital Location file #44, the primary division will be used. Depending on parameters in the Examination Status file, you can require that the procedure be changed to either a Series or Detailed procedure by the radiology or nuclear medicine staff in order for the examination to acquire a given status. Again think about those entering orders; will they be able to discriminate between the more specific Detailed and Series type procedure, or will they be more comfortable selecting a general procedure name from the Broad types?

If you choose to bypass this prompt, the system default is NO and users will be able to request Broad type procedures.

5. ASK 'CAMERA/EQUIP/RM' - Setting this parameter to YES indicates that the question "Primary Camera/Equip/Rm" will be asked when the user is entering data about the examination in status tracking. If you want to track camera, equipment or room use for each procedure, answer YES to this prompt, keeping in mind that the user may still bypass the question if s/he chooses to do so. (The Examination status

parameters can be set to **require** this data before an exam can progress to a Complete (or other) status.)

If you accept the default of NO, the user will not see the Primary Camera/Equip/Rm prompt.

- 6. AUTO USER CODE FILING The access code entered at sign-on is used automatically for the activity logs if this parameter is set to YES. If NO, the user(s) will be asked to enter access codes during various processes. This feature was used in older versions to minimize security risks of unattended logged-in terminals and multiple users sharing one terminal. It should not be necessary at most sites now, and should therefore be set to YES.
- 7. TRACK EXAM STATUS CHANGES Setting this parameter to YES results in the recording of the dates/times of examination status changes. These statuses pertain to the exam, not the request. Statuses are defined by the site using the option Examination Status Entry/Edit. Examples of statuses that might be defined are as follows: Called for Exam, Cancelled, Complete, Examined, Transcribed, Waiting for Exam, etc.

This prompt is now bypassed and automatically set to YES.

8. ASK EXAM STATUS TIME - If the Track Exam Status Changes parameter is set to YES, then this parameter will appear. Setting this parameter to YES indicates that the user should be asked the time of the status change. The Date/Time is added to the Status Change Dt/Time subfile and to the Activity Log.

If set to NO or bypassed, the system automatically uses the current date and time and the prompt for time will not appear. This data is added only to the Activity Log.

9. TIME LIMIT FOR FUTURE EXAMS - This parameter sets a specific number of hours in the future during which a patient can be registered for an exam. If this prompt is bypassed, the user will not be able to register a patient for a future date or time. The number of hours cannot exceed 168 (one week). Generally, registering a patient is done when the patient arrives at the Imaging Location. Registering a patient is not the same as scheduling a patient for a procedure, which takes any future date you want.

Note: Registering patients in advance opens the possibility for exam updating before the exam really takes place. This can cause failure of workload crediting. Registration more than 24 hours in the future is not recommended.

Films Reporting Parameters

- 10. ALLOW STANDARD REPORTS Transcriptionists will be able to select a standard report when using the Report Entry/Edit option if this parameter is set to YES. There is an option, Standard Reports Entry/Edit, which can be used to build a standard report with a report text that is often used. An example of a standard/routine exam might be, Normal Chest with the text reading "Chest exam is within normal limits". If there are a number of such exams reported, allowing the transcriptionists this ability saves time. These standard reports are discussed in a later chapter.
- 11. ALLOW BATCHING OF REPORTS If set to YES, the transcriptionist can group reports, generally by physician name, for easy, organized access or printing. If changes must be made, the transcriptionist can look for the report within a batch rather than by keeping track of record numbers and patient names. It provides a time-saving way to organize work.
- 12. ALLOW COPYING OF REPORTS Contents of one report can be copied to another if this parameter is set to YES. If a patient has had more than one test and the report text and impression of one report are needed for another, this handy ability helps save time. It places the entire report text and impression text of the selected report into the one you are editing. The text can then be further edited. The Clinical History will be moved from the patient's current exam record to the report record, not from the copied report.

- 13. IMPRESSION REQUIRED ON REPORTS Set this parameter to YES if the impression is required on a report in order for it to become verified and for the exam to progress to a status of Complete. (Also, see Examination Status parameter instructions for setting up at least the Complete status to require Impression if this parameter is set to YES.)

 Note: It is highly recommended that this parameter be set to YES. When reports are purged from the system, the clinical history and report text are deleted. Only the impression is retained on line.
- 14. ALLOW VERIFYING BY RESIDENTS This parameter can be set to YES to allow interpreting residents to verify their own reports while doing On-line Verifying of Reports. Further, if this parameter is set to YES and the resident's **individual** parameter, Allow Verifying of Others (in the New Person file), is **also** set to YES, the resident will be allowed to verify other interpreting physicians' reports in addition to their own.

It is recommended that residents NOT be allowed to verify reports because of the legal implications of report verification.

15. ALLOW RPTS ON CANCELLED CASES? - If this field is set to YES, any user can enter a report on a cancelled exam.

If this parameter is set to NO, only users who own the RA MGR key can edit or enter a report on a Cancelled case.

16. WARNING ON RPTS NOT YET VERIF? - If this field is set to YES, results reports in all statuses except "Verified" will contain the status surrounded by asterisks under the body of the printed report. This can be useful in making sure that unverified report copies are not placed in patient charts or mistakenly interpreted as final results. (This is not a new feature. Only the name was changed.)

17. AUTO E-MAIL TO REQ. PHYS? - If this field is set to YES, the system will automatically send the Radiology/Nuclear Medicine findings report to the Requesting Physician via e-mail at the time it is verified. If this field is left blank, verified reports will not be e-mailed to the requesting physician.

When a report is verified, the Auto E-Mail to Req. Phys. Field of the division for the requesting location is used to determine whether to send results to the requesting physician. If the requesting location is not assigned to a division in the Hospital Location file #44, the primary division will be used. **Note: If the requesting location is assigned to a division in the Hospital Location file #44, make sure that division is also defined in the Rad/Nuc Med Division file #79.**

18. ¹ALLOW E-SIG ON COTS HL7 RPTS – If this field is set to YES, the system will automatically add the Electronic Signature Block printed name of the Verifying Physician that signed the report being transmitted to Radiology/Nuclear Medicine via an HL7 interface to a COTS voice reporting system. This parameter is only checked for COTS HL7 interfaces that use the Generate/Process Routine 'RAHLTCPB' in the subscriber protocol. Note: The HL7 COTS interfaces supported by this option are the IBM MedSpeak interface released with patch RA*5*4 and the TalkStation and PowerScribe interfaces released with patch RA*5*12. See the appropriate appendix in the Rad/Nuc Med Technical Manual for additional information on these interfaces.

Miscellaneous Division Parameters

- 19. PRINT FLASH CARD FOR EACH EXAM To have **one** flash card automatically printed with each registration of a **case** for all Imaging Locations, set this parameter to YES. To allow the printing of more than one flash card, set this parameter to NO and use the Location Parameter Set-up option to enter the number you want printed for each location.
- 20. PRINT JACKET LBLS W/EACH VISIT Indicate to the system that a jacket label be printed with each visit (may include many cases/exams) by setting this parameter to YES. The number of labels printed will be the same as the entry in the location parameter How Many Jacket Lbls Per Visit.

-

¹ Patch RA*5*12

21. CONTRAST REACTION MESSAGE - This is an **optional**, free text field to be used for entering a warning message that will be displayed each time a patient who has had a reaction to contrast medium is registered for a procedure that uses a contrast medium.

If you bypass this field, no message will be displayed.

22. RPHARM DOSE WARNING MESSAGE - This field contains the message

displayed to the user when the dose falls outside the range of the "Low Adult Dose" and the "High Adult Dose" for the radiopharmaceutical. These fields reside in the Default Radiopharmaceuticals multiple in the Rad/Nuc Med Procedures file. Use the Procedure Enter/Edit option to adjust these values. If no values exist for a procedure, this message will have no effect.

¹HL7 Applications Associated with this Division

23. HL7 RECEIVING APPLICATION – this field should only be set if using ROUTING LOGIC with the Radiology HL7 subscriber PROTOCOLs (File #101). (See the Radiology/Nuclear Medicine V. 5.0 HL7 Manual for further details.)

If using ROUTING LOGIC then only the HL7 RECEIVING APPLICATIONs defined here will receive Radiology HL7 messages for this Division. If ROUTING logic IS BEING USED, BUT NO hl7 receiving applications ARE DEFINED, then Radiology HL7 messages for the Division will not be created.

If no applications are entered and ROUTING LOGIC is not set up (this is the normal setup), then the system assumes that all HL7 applications should receive HL7 messages created in this Division.

Note: If in doubt, do not set this field.

Imaging Locations Associated with this Division

24. IMAGING LOCATION - This field should contain all Imaging Locations associated with this Division. To determine the locations available for your selection before editing this field, you can use the option Print File Entries within the Supervisor Menu to obtain a list.

Only those locations in the Hospital Location file #44 whose TYPE field is **CLINIC** and whose PATTERN DATE field is NULL will appear in this list. Also, print the Inactivate Date so you can avoid locations that have been inactivated.

OUTPUT FROM WHAT FILE: RAD/NUC MED PROCEDURES// HOSPITAL LOCATION

_

¹ Patch RA*5*25 Added a Division Parameter for HL7 Applications.

```
(148 entries)
SORT BY: NAME// 2 TYPE
START WITH TYPE: FIRST// CLINIC USES INTERNAL CODE: C
GO TO TYPE: LAST// CLINIC USES INTERNAL CODE: C
 WITHIN TYPE, SORT BY: PATTERN
                                     (multiple)
 PATTERN SUB-FIELD: .01 PATTERN DATE
 START WITH PATTERN DATE: FIRST// @
 GO TO PATTERN DATE: LAST// @
   WITHIN PATTERN DATE, SORT BY: <RET>
STORE IN 'SORT' TEMPLATE: <RET>
FIRST PRINT FIELD: NAME
THEN PRINT FIELD: INACTIVATE DATE
THEN PRINT FIELD: <RET>
Heading (S/C): HOSPITAL LOCATION LIST Replace <RET>
START AT PAGE: 1// <RET>
DEVICE: (Printer name)
```

If you are in the process of editing this field in the Division Parameter Set-up option, entering one or two question marks will list your location choices.

New locations may also be entered in the Imaging Locations file #79.1 at this step. Entry of **new** Imaging Locations should be coordinated with Medical Administration Service and/or Information Resource Management. An Imaging Location must be defined in the Division parameters through entry at this step in order to process examinations from that location. Once an Imaging Location is associated with a Division, all the Division parameters become those of the location, unless stated otherwise.

At the completion of your edits, the program checks each Imaging Location associated

with the Division for invalid locations. A valid Imaging Location has a "Type" of CLINIC and has no appointment patterns. An error message will appear for any invalid Imaging Location. Any problem not fixed will result in failure to credit of procedures in that location. Correct the problem, and then run through the Division Parameter Set-up option again to make sure the location returns OK.

Two fields from the Division parameters directly impact the location parameters: Print Flash Card for Each Exam and Print Jacket Lbls W/Each Visit. If the flag for the first parameter is set to YES, then one flash card will be printed automatically for each examination and you will not be asked to enter an amount to print through the Location Parameters Set-Up option. If the Division parameter Print Jacket Lbls W/Each Visit is set to YES, jacket labels will be printed for each examination. The number of labels printed will be determined by the entry in the How Many Jacket Lbls Per Visit location parameter.

Note: Whenever changes are made to the Division parameters, any users currently signed on to the system who want to use the new changes must

exit and reenter the package or use the Switch Locations option to make the new parameters take effect.

Here is an example of how the Division parameters might be set up. Enter one or two question marks at the prompt if you need more information. There are a good number of default answers for your convenience. They are suggestions as to what is a usual answer but they can be changed by you. Striking the return key accepts the default as it is.

Prompt/User Response

Discussion

request tracking purposes.

Division Parameter Set-up

Select Division: NEW DIVISION
VAMC 234
Are you adding 'NEW DIVISION' as a new RAD/NUC MED
DIVISION (the 2ND)? Y (YES)

Division-wide Order Entry Parameters:

Ask 'IMAGING LOCATION': yes// <RET>

Answering "yes" is recommended so that all requests will be associated with an Imaging Location for report sorting and

CLINICAL HISTORY MESSAGE: Enter clinical history relevant to the procedure.

Exam Entry/Edit Parameters:
-----DETAILED PROCEDURE REQUIRED: no// <RET>

ASK 'CAMERA/EQUIP/RM': no// yes

AUTO USER CODE FILING: yes// <RET>

ASK EXAM STATUS TIME: no// Y yes

TIME LIMIT FOR FUTURE EXAMS: 12

We will allow cases to be registered 12 hours in advance of the current date/time.

Films Reporting Parameters:

ALLOW STANDARD REPORTS: no// Y yes

ALLOW BATCHING OF REPORTS: yes// <RET>

ALLOW COPYING OF REPORTS: yes// <RET>

IMPRESSION REQUIRED ON REPORTS: yes// <RET>

The impression is saved even though report texts and clinical histories may be purged. Respond YES to this prompt to prevent verifying and completing reports before an impression is entered.

ALLOW VERIFYING BY RESIDENTS: no// <RET>

ALLOW RPTS ON CANCELLED CASES?: no// <RET>

If you answer NO to this prompt, only users with the RA MGR key may enter reports on cancelled cases.

WARNING ON RPTS NOT YET VERIF?: yes// <RET>

If this field is set to YES, the status of the results report will be printed on the output.

AUTO E-MAIL TO REQ PHYS? yes// <RET>

If this field is set to YES, the report will be automatically sent to the Requesting Physician via email at the time it is verified.

¹ALLOW E-SIG ON COTS HL7 RPTS: yes// <RET>

If this field is set to YES, the electronic signature of the verifying physician is printed on COTS HL7 reports.

Miscellaneous Division Parameters:

PRINT FLASH CARD FOR EACH EXAM: no// \langle RET \rangle

PRINT JACKET LBLS W/EACH VISIT: no// <RET>

CONTRAST REACTION MESSAGE: ***WARNING*** THIS PATIENT HAS HAD CONTRAST MEDIA REACTIONS

This message cannot exceed 70 characters.

RPHARM DOSE WARNING MESSAGE: THIS DOSE DEVIATES FROM THE PRESCRIBED RANGE AS OUTLINED IN THE DEPARTMENTAL PROCEDURES MANUAL, AND REQUIRES A SIGNED, DATED, AND WRITTEN DIRECTIVE BY AN AUTHORIZED PHYSICIAN.

This is a word processing field and may be as long as you like.

²HL7 Applications Associated with this Division:

Select HL7 RECEIVING APPLICATION: RA-TALKLINK-TCP ACTIVE

Example shows how to restrict HL7 messages for this Division to TalkStation Voice Recognition tool only. If no HL7 applications are being used, or only one HL7 application, then leave this field blank.

Select HL7 RECEIVING APPLICATION: <RET>

¹ Patch RA*5*12

² Patch RA*5*25 Added a Division Parameter for HL7 Applications.

Imaging Locations Associated with this Division: -----

Select IMAGING LOCATION: ?

Answer with IMAGING LOCATIONS

You may enter a new IMAGING LOCATIONS, if you wish Enter a location for this Division.

Answer with IMAGING LOCATIONS, or TYPE OF IMAGING

To obtain a list of locations available to you, enter a ? at this prompt. The entries in the Imaging Location file display with their Imaging Types in parentheses to the right. Entering ?? will display a list of all possible selections. If you do not see the location that you want, it is either not in the Hospital Location file or the Type field of the Hospital Location file does not contain CLINIC.

Choose from:

X-RAY(GENERAL RADIOLOGY-499)

NUC MED LOC (NUCLEAR MEDICINE-639)

ULTRASOUND (ULTRASOUND-499)

MAGNETIC RESONANCE IMAGING (MAGNETIC RESONANCE IMAGING-499)

NUC(NUCLEAR MEDICINE-499) FLUORO (GENERAL RADIOLOGY-499) RAD(GENERAL RADIOLOGY-639)

CAT SCAN (CT SCAN-499) US (ULTRASOUND-578)

You may enter a new IMAGING LOCATIONS, if you wish

Select the appropriate location.

The entry selected must be set up as an

'Imaging' type in the Hospital Location file.

Select IMAGING LOCATION: X-RAY (GENERAL

RADIOLOGY-499)

...OK? yes// **<RET>** (YES)

GENERAL RADIOLOGY

Location is already assigned to the HINES IRMFO

Division! ??

Remember a location can only belong to one Division. This is telling us X-RAY already belongs to Hines IRMFO.

```
(ULTRASOUND-578)
Select IMAGING LOCATION: US
     ...OK? yes// <RET> (YES)
     ULTRASOUND
Are you adding 'US' as a new IMAGING LOCATIONS (the
1ST for this RAD/NUC MED DIVISION)? Y (YES)
Select IMAGING LOCATION: RAD (GENERAL
RADIOLOGY-639)
     ...OK? yes// <RET> (YES)
     GENERAL RADIOLOGY
Are you adding 'RAD' as a new IMAGING LOCATIONS (the
2ND for this RAD/NUC MED DIVISION)? Y (YES)
Select IMAGING LOCATION: <RET>
Division Parameters have been set!
Imaging Location X-RAY is OK.
Imaging Location ULTRASOUND is OK.
Imaging Location CAT SCAN is OK.
```

Print Division Parameter List [RA SYSDIVLIST]

Once all of the parameters are set for each of your Divisions, use this option to print a dated, hard copy of the settings for reference. Any time you make a change to the parameters, it is useful to keep a copy for historical reference. Here we show the settings for two Divisions, New Division and Hines CIO Field Office.

Prompt/User Response

Discussion

```
<sup>1</sup>Print Division Parameter List

DEVICE: (Printer Name) RIGHT MARGIN: 80// <RET>
```

```
Rad/Nuc Med Division List
                                           AUG 17,1994 14:44 PAGE 1
Division: NEW DIVISION
 Track Request Status Changes? : NO
Track Exam Status Changes? : yes
Ask for Exam Status Time? : yes
  Print Flash Card For Each Exam? : no
 Ask for Imaging Location? : NO
 Ask for Requested Date?
 Require Detail Procedure...
       ...upon Initial Exam Entry?: no
  Print Jacket Labels...
                ...with each Visit?: no
 Ask for Camera/Equipment/Room? : yes
  Allow Rpts on Cancelled Exams? : no
  Auto e-mail to Req. Physicians? : yes
  Time Limit For Future Exams...
                   ...in Hours?: 0
 Allow Standard reports? : yes Allow Batching of reports? : yes
```

-

¹ Patch RA*5*25 Added new fields and reformatted parameters for better display.

```
AUG 17,1994 14:44 PAGE 2
Rad/Nuc Med Division List
       Allow Copying of reports?
 Impression Required on reports? : yes
 Allow Verifying by Residents? : no Interpreting Staff Required? : NO
 Warning on reports not yet...
                   ... Verified? : yes
 Allow e-sig on COTS HL7 reports? :yes
 Contrast Message:
   ***WARNING*** THIS PATIENT HAS HAD CONTRAST MEDIA REACTIONS
 Clinical History Message:
 _____
   ***Enter clinical history relevant to procedure***
 Resource Device:
 _____
```

```
RAD/Nuc Med Division List

AUG 17,1994 14:44 PAGE 3

RAD RESOURCE DEVICE

Radiopharmaceutical Dose Warning Message:

THIS DOSE DEVIATES FROM THE PRESCRIBED RANGE AS OUTLINED IN THE DEPARTMENTAL PROCEDURES MANUAL, AND REQUIRES A SIGNED, DATED, AND WRITTEN DIRECTIVE BY AN AUTHORIZED PHYSICIAN.

HL7 Receiving Applications for Division:

RA-PSCRIBE-TCP

Imaging Locations:

US
RAD
NUC
```

```
AUG 17,1994 14:44 PAGE 4
Rad/Nuc Med Division List
_____
Division: HINES CIO FIELD OFFICE
 Track Request Status Changes? : YES
 Track Exam Status Changes? : yes
Ask for Exam Status Time? : yes
 Print Flash Card For Each Exam? : no
 Ask for Imaging Location? : YES
 Ask for Requested Date?
                                : yes
 Auto User Code Filing?
                                : yes
 Require Detail Procedure...
       ...upon Initial Exam Entry?: no
 Print Jacket Labels...
             ... with each Visit?: no
 Ask for Camera/Equipment/Room? : yes
 Allow Rpts on Cancelled Exams? : no
 Auto e-mail to Req. Physicians? : yes
 Time Limit For Future Exams...
                 ...in Hours?: 96
 Allow Standard reports? : yes Allow Batching of reports? : yes
```

```
AUG 17,1994 14:44 PAGE 5
Rad/Nuc Med Division List
 Allow Copying or reports? : yes
 Impression Required on reports? : yes
 Allow Verifying by Residents? : yes
 Interpreting Staff Required?
                               : NO
 Warning on reports not yet...
                    ...Verified?: yes
 Allow e-sig on COTS HL7 reports? : yes
 Contrast Message:
 _____
   ***WARNING*** THIS PATIENT HAS HAD CONTRAST MEDIA REACTIONS
 Clinical History Message:
   **CLINICAL HISTORY IS MANDATORY**
 Resource Device:
 _____
```

Rad/Nuc Med Division List

AUG 17,1994 14:44 PAGE 6

Radiopharmaceutical Dose Warning Message:

HL7 Receiving Applications for Division:

RA-TALKLINK-TCP

Imaging Locations:

X-RAY
ULTRASOUND
MAGNETIC RESONANCE IMAGING
FLUORO
CAT SCAN

Camera/Equip/Rm Entry/Edit [RA SYSEXROOM]

This option is used to enter cameras, equipment, and/or rooms that are used by the department. It will provide the list of choices for a user when entering a patient's procedure if, within the Division Parameters Set-up option, the field "Ask Camera/Equip/Room" contains a YES. If you have already decided that tracking this information is not important to your service, then skip this option.

When defining the parameters in each Imaging Location (Location Parameter Set-up), you will be asked to enter any cameras, equipment and/or rooms associated with the location. Each can be associated with one or more locations, since some sites use the same room across locations, Imaging Types, and Divisions. Cameras/Equipment/Rooms cannot be deleted through this option. However, you can "inactivate" them by removing them from the Imaging Location through the Location Parameter Set-up option.

- 1. CAMERA/EQUIP/ROOM Enter new cameras, equipment, or rooms at this prompt. This field is directly associated with one or more Imaging Locations.
- 2. CAMERA/EQUIP/RM DESCRIPTION Enter a description 3 40 characters in length to describe the type of camera, equipment, room or procedures performed.

Prompt/User Response

Camera/Equip/Rm Entry/Edit Select Camera/Equip/Room: ??

Choose from:

CAMERA 1 Triple Head Spect Scanner
CAT SCAN PICKER PQ200 SCANNER
INP1 X-RAY GASTRO X-RAYS
INP2 X-RAY BROKEN LEG X-RAYS

MAMMOGRAPHIC Used in mammography exams. MRI ROOM MAGNETIC RESONANCE IMAGING ROOM

NUC1 HEART IMAGING NUC2 BRAIN IMAGING

OUTP1 X-RAY X-RAYS ONLY IN THIS ROOM OUTP2 X-RAY ACCIDENT X-RAYS

PORTABLE PORTABLE

This field contains the name of this camera/equipment/room. Names must be between 1 and 30 characters in length.

Discussion

Enter one or two? to obtain a list of selections already in the file.

System Definition Menu

Select Camera/Equip/Room: **PORTABLE II**Are you adding 'PORTABLE II' as a new CAMERA/EQUIP/RM (the 12TH)? **Y** (YES)

CAMERA/EQUIP/RM DESCRIPTION: ?
Enter a brief description, between 3 and 40 characters in length, to define the types of procedures performed.

CAMERA/EQUIP/RM DESCRIPTION: Portable used in

Surgery

CAMERA/EQUIP/RM: PORTABLE II// <RET>

DESCRIPTION: Portable used in Surgery// <RET>

Select Camera/Equip/Rm: <RET>

This field can be used to describe this camera/equipment/room and define what types of procedures may be performed there. "Triple Head Spect Scanner" is an example of a description.

List of Cameras/Equip/Rms [RA SYSEXLIST]

Once you complete your entries for the Camera/Equip/Rm file #78.6, use this option to print a dated, hard copy of the contents of the file.

Prompt/User Response

Discussion

List of Cameras/Equip/Rms

DEVICE: (Printer Name) RIGHT MARGIN: 80// <RET>

Imaging Camera/Equ	ipment/Room List	AUG 18,1994 08:08 PAGE 1
CAMERA/EQUIP/RM	RAD/NUC MED DIVISION	LOCATION
_		
CAMERA 1	HINES INFORMATION SY	Y-RAY
CAT SCAN	HINES INFORMATION SY	
INP1 X-RAY	HINES INFORMATION SY	X-RAY
INP1 X-RAY	HINES INFORMATION SY	FLUORO
INP2 X-RAY	HINES INFORMATION SY	X-RAY
MRI ROOM	HINES INFORMATION SY	MAGNETIC RESONANCE I
NUC1	NEW DIVISION	NUC MED LOC
NUC2	NEW DIVISION	NUC MED LOC
NUC2	HINES INFORMATION SY	NUC
OUTP1 X-RAY	HINES INFORMATION SY	X-RAY
OUTP1 X-RAY	HINES INFORMATION SY	FLUORO
OUTP2 X-RAY	HINES INFORMATION SY	X-RAY
PORTABLE	NEW DIVISION	NUC MED LOC
PORTABLE	HINES INFORMATION SY	X-RAY
PORTABLE	HINES INFORMATION SY	ULTRASOUND
PORTABLE II	**UNASSIGNED**	**UNASSIGNED**

Note that "PORTABLE II", created on the previous page, appears on this list as "**UNASSIGNED**". This is because "PORTABLE II" has not been assigned an Imaging Location. All camera/equipment/rooms that have been entered in the Camera/Equip/Rm file, but not assigned to any Imaging Location will appear at the end of this list.

A camera/equipment/room can belong to more than one division/location.

Location Parameter Set-up [RA SYSLOC]

The Radiology/Nuclear Medicine program is designed to handle multiple locations within a Division. Each location must have its own set of parameters. The Location Parameter Set-up option is used by the coordinator to define Imaging Location parameters for new locations and edit location parameters for existing locations. **At least one location must be entered for each Division**. The following is a listing of location parameters and brief explanations.

1. LOCATION - An existing location may be edited or a new location may be entered at this prompt. Before entering a new location, use the option Print File Entries within the Supervisor Menu to see what locations are already in the file. Print the Imaging Type entries in the Hospital Location file #44 (sort by Type (from CLINIC to CLINIC) and print Name and Inactivate Date). Any location whose type is CLINIC is available for you to select at this prompt. Any locations you need but are not listed, you will need to create using this Location Parameter Set-up option. For the purpose of workload crediting, when a new location is selected, it will automatically be flagged as an "occasion of service" clinic in File 44. If there are appointment patterns, a duplicate clinic of the same name will automatically be added in File 44 and marked "occasion of service." Note: Always contact MAS and coordinate with them the addition of new locations.

When you and MAS add a new location that is not in the Hospital Location file, you must also enter a Hospital Location Type. The "HOSPITAL LOCATION TYPE:" prompt should always be answered with "CLINIC".

Flash Card, Jacket Labels, and Exam Label Param

Note: Labels can be printed from the option Duplicate Flash Card even if the site parameters are not set up to print labels.

Note: Printers for Flash Cards, Jacket Labels, Requests, and Reports for each Imaging Location are set up via the IRM Menu option Device Specifications for Imaging Locations. If printers have not been assigned, you will get a message "No default flash card printer has been assigned. Contact IRM."

If under the Division parameters, you answered YES to "PRINT FLASH CARD FOR EACH EXAM:" and/or "PRINT JACKET LBLS W/EACH VISIT:", you **will not** see the "HOW MANY..." prompt.

2. HOW MANY (FLASH CARDS, JACKET LABELS) PER VISIT - This is a whole number between 0 and 10 to establish a parameter for printing the number of cards and/or labels needed for each patient visit at this location.

This question will **not** appear if you entered YES at the Division parameter, "PRINT FLASH CARD FOR EACH EXAM:" and this location is associated with the Division. Bypassing each of these prompts or entering 0 will suppress the printing of the item.

HOW MANY EXAM LABELS PER EXAM - Exam labels are printed for each registered case. This is a whole number between 0 and 9.

3. DEFAULT (FLASH CARD, JACKET LABEL, EXAM LABEL) FORMAT - A default format can be entered for each card or label type. Each time cards or labels are printed for an exam registered to this location, they will be printed using the specified format. Formats for each may be designed using the option Label/Header/Footer Formatter.

Order Entry Parameters

This parameter identifies the printer used to print request forms (form 519) as procedures are ordered. It is set up via the IRM Menu option Device Specifications for Imaging Locations, at the prompt, "REQUEST PRINTER NAME:". If a printer has not been assigned to this prompt, a message will tell you to contact IRM.

Report Parameters

- 4. DEFAULT REPORT HEADER FORMAT A report header format can be selected that will be printed at the top of all patient result reports printed from this location. Header formats can be designed using the option Label/Header/Footer Formatter.
- 5. DEFAULT REPORT FOOTER FORMAT This is a report footer format that will be printed at the bottom of all patient result reports printed from this location. Footer formats can be designed using the option Label/Header/Footer Formatter.

- 6. REPORT LEFT MARGIN This is a column number which indicates the left margin of patient result reports. This only affects the left margin for the impression, report text, and clinical history.
- **7.** REPORT RIGHT MARGIN The number entered here is used as the Report Right Margin. This affects the right margin for the impression, report text, and clinical history.
- 8. PRINT DX CODES IN REPORT? If this field is set to YES, all primary and secondary diagnostic codes associated with the exam will be printed on the exam report. Since the primary diagnostic code may be printed in exam report headers and footers, you may also need to change report header and footer setups to get the diagnostic codes printed exactly as you wish. If this field is set to NO, you can include the primary diagnostic code in headers or footers by using the Label/Header/Footer Formatter.

Note: The device for this parameter is set up via the IRM Menu option Device Specifications for Imaging Locations. If a printer has not been assigned to this prompt, a message will tell you to contact IRM.

9. VOICE DICTATION AUTO-PRINT - If this field is set to YES, and voice dictation is used to enter a report, a copy of the report will be automatically printed on the device defined in the field Report Print Name.

Cameras/Equip/Rooms Used by this Location

10. CAMERA/EQUIP/RM - This multiple field will contain all cameras, equipment, or rooms associated with this Imaging Location. Cameras, equipment, and rooms already entered through the Camera/Equip/Rm Entry/Edit option may be used or new ones can be entered at this prompt.

¹Default CPT Modifiers used by this Location

Default CPT modifier(s) entered for this location will be automatically entered into each exam record during exam registration when the following 2 conditions are met:

- There is no default CPT Modifier for this exam's procedure.
- This location is the current sign-on (or switched-to) location.

Values entered for this field are not screened, because there is no CPT CODE for this location with which to compare it. The default values entered here may be rejected later during exam registration, when the default CPT Modifier(s) are screened, using the exam procedure's CPT CODE. If so, the following message will appear:

** Unable to enter default CPT Modifier __ (Modifier name) **

Recipients of the 'Stat' Alert for this Location

- 11. STAT REQUEST ALERT RECIPIENTS To receive a stat request alert, a recipient must have a Radiology/Nuclear Medicine personnel classification. They will receive an alert whenever a STAT request is submitted to this Imaging Location. To use this feature, the Division parameter, "Ask Imaging Location" must be set to YES, and OE/RR Version 3.0 or higher must be installed.
- 12. **URGENT REQUEST ALERTS? This field will determine if alerts will be sent to the Imaging Personnel defined in the "Stat Request Alert Recipients" field. If this field is set to YES, then all the users in the "Stat Request Alert Recipients" field will be sent an alert for urgent requests. To use this feature, the Division parameter "Ask Imaging Location" must be set to YES, and OE/RR Version 3.0 or higher (CPRS) must be installed.
- 13. **ALLOW 'RELEASED/NOT VERIFIED' This Imaging Location parameter can be set to YES to allow reports to have the status "RELEASED/NOT VERIFIED" which permits the displaying and printing of the report to users outside of the Imaging service even if it is not yet verified. This field was moved from the Division parameters.
- 14. INACTIVE This field should only be used if an Imaging Location is no longer used. If a date is entered in this field, the Imaging Location will not be available for selection at exam registration or request submission (i.e. the "SUBMIT REQUEST TO" prompt will not allow ordering clinicians to choose this location for request printing.) However, users can still select an inactive location to process exams previously registered to it.

_

¹ Patch RA*5*19 May 2000 New field added to Imaging Locations file #79.1.

Imaging Type for this Location

15. TYPE OF IMAGING - Type of imaging done at this location. Assign **one** of the following Imaging Types to your location:

Angio/Neuro/Interventional CT Scan General Radiology Magnetic Resonance Imaging Cardiology Studies (Nuc Med) Nuclear Medicine Ultrasound Vascular Lab Mammography

This Location parameter should contain the name of the type of imaging from the Imaging Type file (#79.2). If this field is filled in, the system automatically adds the proper Imaging Type on the registration record. Only procedures assigned the same Imaging Type can be registered to the location.

Assigning a new Imaging Type (i.e. one not previously selected for any other location) implies that procedures of the new Imaging Type are performed under a separate service or section at your site, and separate workload report is desired. For instance, at some sites Ultrasound may be done separately from Radiology, and it would be appropriate to assign the Ultrasound Imaging Type to one or more locations. Once the Ultrasound Imaging Type is activated in this way, and procedures designated as Ultrasound Imaging Type have been registered, the system will report Ultrasound workload separately from Radiology workload.

Note: The Cardiology Studies Imaging Type should NOT be assigned to a location unless nuclear cardiology studies are done under a service other than Nuclear Medicine.

The software was designed to handle nuclear cardiology alone, and therefore other types of cardiology examinations such as echocardiography should not be entered into the system. Further, local modifications to add Imaging Types are strongly discouraged due to the complexity of the underlying code and associated workload implications.

Be cautious about assigning an Imaging Type not previously used to a location. This "activates" the new Imaging Type. Imaging workload reports, enter/edit options, and screen displays are all separated by their Imaging Types. Be sure this is how your service wants to handle the package. Once this package is loaded, users sign on to a specific location, which is associated with an Imaging Type. Procedures are also associated with an Imaging Type. The best bet is to look at workload for each Imaging Type and each location. Are each now treated almost separately from all others?

- 16. CREDIT METHOD This field indicates the type of credit this location receives for examinations. The default is REGULAR (code 0), which means that for this Imaging Location both the exam and the interpretation are done at your facility. If this Imaging Location has been set up to register exams done in another facility, with the imaging interpretation only done at this facility, enter code 1 (Interpretation Only). If exams are done in this location, but interpretation is done elsewhere, enter code 3 (Technical Component Only). If both the exam and the interpretation are done elsewhere, and exams entered for this Imaging Location are entered only for record-keeping purposes, enter code 2 (No Credit). (This field was added to the software with the installation of patch RA*4.5*4.) Until changes are made in the other VistA software packages and in the Austin system, codes 1 and 3 will remain inactive and not selectable.
- 17. DSS ID This field is used to identify the general outpatient area where an imaging exam was done. It is used only by DSS, and only when processing outpatients. It is roughly equivalent to the "Stop Codes" previously used for crediting outpatient exams. Only entries also in the Imaging Stop Codes file #71.5 are valid. (This field was added to the software with the installation of patch RA*4.5*4.)

Note: There are three things you must remember to do now:

- 1) For every location you set up, be sure each is assigned to a Rad/Nuc Med Division in file #79, using the option Division Parameter Set-up.
- 2) Printers may not be set up for every location. If you see the following message when editing this option, you will need to provide IRM with the name of the location and determine with IRM which printer(s) should be used to print labels when exams are registered in this location and which printer should print requests submitted to this location. Also, a default results report printer may be assigned.

```
No default flash card printer has been assigned. Contact IRM.
```

3) Whenever changes are made to the location parameters, any users currently signed on to the system who want to use the new location must use the Switch Locations option or exit and reenter the package to make the new parameters take effect.

Here is an example of editing a new location:

Prompt/User Response

Discussion

Location Parameter Set-up
Select Location: ??

Enter two question marks to obtain a list of Imaging Locations already set up. The Imaging Type of each location appears to the right of the location name.

```
Choose from: X-RAY
```

FLUORO

X-RAY
NUC MED LOC
ULTRASOUND
MAGNETIC RESONANCE IMAGING
NUC

(GENERAL RADIOLOGY-449)
(NUCLEAR MEDICINE-639)
(ULTRASOUND-449)
(MAGNETIC RESONANCE IMAGING-449)
(NUCLEAR MEDICINE-449)

(GENERAL RADIOLOGY-449)

This field points to the Hospital Location File, and contains the name of this location.

You also get a list of all Choose from: available locations that are CAT SCAN in the Hospital Location file FLUORO #44 whose "Type" field is MAGNET set to "CLINIC". MAGNETIC RESONANCE IMAGING MAMMOGRAPHY NUC MED LOC NUCLEAR MEDICINE ULTRASOUND WESTSIDE XRAY X-RAY Name given by user to any ward, clinic, fileroom, operating room or other location within a VA facility. Select Location: MAMMOGRAPHY Are you adding 'MAMMOGRAPHY' as a new IMAGING LOCATIONS (the 10TH)? Y (Yes) IMAGING LOCATIONS TYPE OF IMAGING: ? Enter an Imaging Type for this location. Answer with IMAGING TYPE, TYPE OF IMAGING, or ABBREVIATION Choose from: ANGIO/NEURO/INTERVENTIONAL CARDIOLOGY STUDIES (NUC MED) CT SCAN GENERAL RADIOLOGY MAGNETIC RESONANCE IMAGING NUCLEAR MEDICINE ULTRASOUND VASCULAR LAB MAMMOGRAPHY IMAGING LOCATIONS TYPE OF IMAGING: MAMMOGRAPHY Since you have added a new Imaging Location, remember to assign it to a Rad/Nuc Med Division through Division Parameter Set-up. LOCATION: MAMMOGRAPHY// <RET>

Flash Card Parameters:

HOW MANY FLASH CARDS PER VISIT: 1

System Definition Menu

DEFAULT FLASH CARD FORMAT: ??

This 'locations' parameter can contain the default format (from the 'Flash Cards Formats' file) to be used when queuing a flash card to print. This format is the default 'flash card' for the current location.

Choose from:

EXAM LABEL1
FLASH1 FORMAT
FOOTER1
HEADER1
JACKET1 FORMAT
TEST FORMAT

DEFAULT FLASH CARD FORMAT: FLASH1 FORMAT

No default flash card printer has been assigned. Contact IRM.

Whenever you have a question about what to enter, use one or two question marks to give you information about the prompt. Here we get a list of available formats. If there are no selections, build some formats using the option Label/Header/Footer Formatter.

This is the notice telling you that no printer has been assigned by IRM to this location for use as a flash card printer.

Jacket Label Parameters:

HOW MANY JACKET LBLS PER VISIT: 1

DEFAULT JACKET LABEL FORMAT: JACKET1 FORMAT

No default jacket label printer has been assigned. Contact IRM.

Exam Label Parameters:

HOW MANY EXAM LABELS PER EXAM: 1// <RET>

DEFAULT EXAM LABEL FORMAT: EXAM LABEL1

Order Entry Parameters:

No default request printer has been assigned. Contact TRM .

Report Parameters:

DEFAULT REPORT HEADER FORMAT: **HEADER1**DEFAULT REPORT FOOTER FORMAT: **FOOTER1**

All the labels and formats were built using the option Label/Header/Footer Formatter and reside in the LBL/HDR/FTR Formats file #78.2. (Device settings may override these margins.) For typical usage or $8 \Omega \times 11$ paper, margins of 10 and 70 are common.

REPORT LEFT MARGIN: 10// <RET>
REPORT RIGHT MARGIN: 70// <RET>

PRINT DX CODES IN REPORT?: ??

If this field is set to Yes, all primary and secondary diagnostic codes associated with the exam will be printed on the exam report.

(Note: Since the primary diagnostic code may be printed in exam report headers and footers, you may also need to change report header and footer set-ups to get the diagnostic codes printed exact as you wish.)

Choose from:

Y yes N no

PRINT DX CODES IN REPORT?: Y yes

No default report printer has been assigned. Contact $\ensuremath{\mathsf{IRM}}$.

```
Cameras/Equip/Rooms Used by this Location:
______
Select CAMERA/EQUIP/RM: ??
This field contains a list of cameras/equip/rms (from
the 'Camera/Equip/Rm' file) that can be chosen when
entering data for exams done in this 'location'.
Cameras/equip/rms can be used by more than one
'location'.
Choose from:
  CAMERA 1 Triple Head Spect Scanner
CAT SCAN PICKER PQ200 SCANNER
INP1 X-RAY GASTRO X-RAYS
INP2 X-RAY BROKEN LEG X-RAYS
MAMMOGRAPHIC Used in mammography exams.
MRI ROOM MAGNETIC RESONANCE IMAGING ROOM
  NUC1 HEART IMAGING NUC2 BRAIN IMAGING
  OUTP1 X-RAY X-RAYS ONLY IN THIS ROOM OUTP2 X-RAY ACCIDENT X-RAYS
   PORTABLE PORTABLE
  PORTABLE II Portable used in Surgery PORTABLE III PORTABLE IN ECC
This field contains the name of this
camera/equipment/room. Names must be between 1 and
30 characters in length.
Select CAMERA/EQUIP/RM: MAMMOGRAPHIC Used in
mammography exams.
 Are you adding 'MAMMOGRAPHIC' as a new
CAMERAS/EQUIP/RMS (the 1ST for this IMAGING
LOCATIONS)? Y (YES)
Select CAMERA/EQUIP/RM: <RET>
<sup>1</sup>Default CPT Modifiers used by this Location:
                                                      More than one modifier
_____
                                                      may be entered.
+----+
| Your entry cannot be compared with a CPT CODE, so be very sure |
| that this is the Default CPT Modifier that you want to stuff |
| into every registered exam from this imaging location.
+-----+
Select DEFAULT CPT MODIFIERS (LOC): 26 PROFESSIONAL COMPONENT
 Are you adding '26' as a new DEFAULT CPT MODIFIERS (LOC) (the 1ST for this
```

IMAGING LOCATIONS)? No// Y (Yes)

Select DEFAULT CPT MODIFIERS (LOC): <RET>

_

¹ Patch RA*5*19 May 2000 New field for Imaging Locations file #79.1

Recipients of the 'Stat' Alert for this Location:

Select STAT REQUEST ALERT RECIPIENTS: PROVIDER1, ONE

Select STAT REQUEST ALERT RECIPIENTS: <RET>

URGENT REQUEST ALERTS?: y yes

ALLOW 'RELEASED/UNVERIFIED': no// <RET>

INACTIVE: <RET>

TYPE OF IMAGING: GENERAL RADIOLOGY// <RET>

CREDIT METHOD: 0

DSS ID: 703

Select Location: <RET>

Only persons with Radiology/Nuclear Medicine personnel classification can be entered at this prompt.

If this field is set to YES, then all the users in the STAT REQUEST ALERT RECIPIENTS field be sent an alert for urgent requests.

Entering YES in this field permits the displaying of any report from this location to users outside of the Imaging service even if it is not verified.

If this location is now inactive, enter the date when it became inactive

Enter:

- 0 Regular Credit
- 1 Interpretation Only
- 2 No Credit
- 3 Technical

Component Only

Only entries also in the Imaging Stop Codes file #71.5 are valid as selections.

Location Parameter List [RA SYSLOCLIST]

After setting up all your locations, you should print out a hard copy of the data. **Make sure printers have been designated for each location where appropriate and each location has been given a Division. Any location not assigned to a Division will not appear on this report.** Note the sort order of this listing: Imaging Type, Division, Imaging Location.

Prompt/User Response

Discussion

```
Location Parameter List
Do you wish to include inactive Imaging Locations? No// <RET> NO
Select Imaging Location: All// <RET>
                                                           In this example, we are
                                                           eliminating locations that
Another one (Select/De-Select): -nuc
                                                           belong to Nuclear
     1 NUC
                                                           Medicine. We accepted the
     2 NUC MED LOC
                                                           default ALL because it was
Choose 1-2: 2
                            (NUCLEAR MEDICINE-449)
                                                           easier to eliminate the 2
Another one (Select/De-Select): -nuc
                                                           Nuclear Medicine locations
(NUCLEAR MEDICINE)
                                                           than it was to enter the 5
Another one (Select/De-Select): <RET>
                                                           Radiology locations.
                                                           The following is just one
DEVICE: (Enter printer) RIGHT MARGIN: 80// <RET>
                                                           page from a printout.
```

```
OCT 10,1996 09:14
Imaging Location List
                                                              PAGE 7
               Location: X-RAY
               Division: HINES INFORMATION SYSTEMS CTR
               Inactivation Date:
           Imaging Type: GENERAL RADIOLOGY
          Credit Method: Regular Credit
                 DSS ID: X-RAY
               Allow 'Released/Unverified': no
   Flash Card Parameters:
                                            Jacket Label Parameters:
   Default Printer: P-DOT MATRIX BACK
                                            Default Printer: P-DOT MATRI
   Cards Per Visit:
                                            Labels Per Visit: 0
   Default Format : TESTALL
                                            Default Format : JACKET1 FOR
   Exam Label Parameters:
                                            Order Parameters:
   Default Printer: [Flash Card Printer]
Labels Per Exam: 0
                                           Default Printer: LINE
                                           Cancellation Ptr: P-DOT MATRI
   Default Format : EXAM LABEL1
                                            Registered Ptr: RAD DIAG
                                            Radiopharmaceutical
   Reporting Parameters:
                                            Dosage Ticket Parameters
   -----
   Default Printer: P-DOT MATRIX BACK
                                            Default Printer:
   Default Header: HEADER1
   Default Footer: FOOTER1
   DX Codes in Rpt: yes
   Voice Dictation auto-print: yes
   Allowable Cam/Equip/Rms For This Location:
       OUTP1 X-RAY
       OUTP2 X-RAY
       INP1 X-RAY
       INP2 X-RAY
       PORTABLE
       CAMERA 1
       ROOM 1
| 1 Default CPT Modifiers for this location:
   ______
       26
             PROFESSIONAL COMPONENT
   'Stat' Alert Recipients For This Location: Urgent Request Alerts?: yes
   _____
       PROVIDER, ONE
       PROVIDER, TWO
       PROVIDER, THREE
```

_

¹ Patch RA*5*19 May 2000 New field for Imaging Locations file #79.1

VI. Rad/Nuc Med Personnel Menu ...

Supervisor Menu ... Rad/Nuc Med Personnel Menu ...

> Classification Enter/Edit Clerical List Interpreting Resident List Interpreting Staff List Technologist List

Overview:

This menu is used to assign and track classifications of personnel who work in Diagnostic Radiology and Nuclear Medicine Services. There is only one edit option in this menu, Classification Enter/Edit. The other four options provide dated, hard copies of the data.

Classification Enter/Edit [RA PNLCLASS]

This function allows you to enter or edit personnel and their corresponding classifications. An employee will usually have only one user classification assigned to him or her at a given time; however, more than one classification per employee is allowed, as long as "Interpreting Resident" and "Interpreting Staff" are not assigned to the same employee.

The Imaging Location Access prompt was changed to Rad/Nuc Med Location Access.

There are currently four classifications the personnel can be placed under. They are as follows:

- 1) Technologist
- 2) Interpreting Resident
- 3) Interpreting Staff
- 4) Clerk

All users requiring an electronic signature code (generally physician staff) should be assigned the XUSESIG option of KERNEL by the site manager.

- 1. All users of this package must be assigned Radiology/Nuclear Medicine location access. If they do not have at least one Radiology/Nuclear Medicine location, they will not be able to access the system.
- 2. Whenever a user is on the system and their Radiology/Nuclear Medicine location access is edited via this option, that user must exit the package and reenter or use the Switch Locations option to have the new Radiology/Nuclear Medicine location access take effect.

The prompts for technologist and clerk are identical, so only one example is shown here:

Prompt/User Response

Discussion

```
Classification Enter/Edit
Select Personnel: ??

Choose from:
   PROVIDER, FOUR
   PROVIDER, FIVE
   ...
   '^' TO STOP: ^
```

Select Personnel: PROVIDER, SIX

Select RAD/NUC MED CLASSIFICATION: ??

Defines the role of this individual in the Radiology/Nuclear Medicine world.

Choose from:

T technologist
R resident
S staff
C clerk

Multiple classifications for one person are not usually necessary. Certain exam field edits are screened by classification. For example, only people who have a "technologist" classification are valid choices at the technologist prompt during case edits. Occasionally, residents are also classified as technologists because they act as technologists. ADPACs sometimes give themselves multiple classifications for testing purposes.

Select RAD/NUC MED CLASSIFICATION: ${f T}$ (technologist) Are you adding 'technologist' as a new RAD/NUC MED CLASSIFICATION (the 1ST for this NEW PERSON)? ${f Y}$ (YES)

Select RAD/NUC MED CLASSIFICATION: <RET>

Select RAD/NUC MED LOCATION ACCESS: ??

This field will allow the user to access all data specific to a particular Radiology/Nuclear Medicine Location.

Choose from:

X-RAY

NUC MED LOC

ULTRASOUND

MAGNETIC RESONANCE IMAGING

NUC

FLUORO

WESTSIDE XRAY

(GENERAL RADIOLOGY-449)

(NUCLEAR MEDICINE-639)

(MAGNETIC RESONANCE IMAGING-449)

(NUCLEAR MEDICINE-449)

(GENERAL RADIOLOGY-449)

(GENERAL RADIOLOGY-639)

CAT SCAN (CT SCAN-449)
US (ULTRASOUND-578)

RAD (GENERAL RADIOLOGY-578)

Rad/Nuc Med Personnel Menu

Select IMAGING LOCATION ACCESS: MAMMOGRAPHY GENERAL RADIOLOGY

Are you adding 'MAMMOGRAPHY' as a new IMAGING LOCATION ACCESS (the 1ST for this NEW PERSON)? Y (YES)

Select IMAGING LOCATION ACCESS: <RET>

RAD/NUC MED INACTIVE DATE: <RET>

Since this employee can access Mammography and Mammography is within General Radiology type of imaging, this employee will have access to the exam data belonging to the Mammography location and, in options where access is determined by Imaging Type, this employee will also be able to view data for all of General Radiology.

Select Personnel: PROVIDER, SEVEN

Select RAD/NUC MED CLASSIFICATION: **R** (resident)
Are you adding 'resident' as a new RAD/NUC MED
CLASSIFICATION (the 1ST for this NEW PERSON)? **Y**(YES)

Select RAD/NUC MED CLASSIFICATION: <RET>

Select IMAGING LOCATION ACCESS: MAMMOGRAPHY GENERAL RADIOLOGY

Are you adding 'MAMMOGRAPHY' as a new IMAGING LOCATION ACCESS (the 1ST for this NEW PERSON)? Y (YES)

Select IMAGING LOCATION ACCESS: <RET>

RAD/NUC MED INACTIVE DATE: <RET>

STAFF REVIEW REQUIRED: ??

This field applies to 'Interpreting Resident' personnel. If it contains a 'yes', an interpreting staff is required to review this resident's report results. Choose from:

Y yes N no

STAFF REVIEW REQUIRED: Y yes

The same prompts appear for both resident and staff. Here is an example of data entry for resident. ALLOW VERIFYING OF OTHERS: ??

If this field is set to 'YES' and the 'ALLOW VERIFYING BY RESIDENTS' Division parameter is also set to 'YES' then this resident is allowed to verify reports associated with another interpreting physician. (If both parameters are set to 'YES' the 'On-line Verifying of Reports' option will prompt the user to 'Select Interpreting Physician:' allowing the user to select an interpreting physician other than him/herself.) If this field is set to 'NO' then this resident is only allowed to verify his/her own reports. If the Division parameter 'ALLOW VERIFYING BY RESIDENT' is set to 'NO' then regardless of how this field is set, the resident will not be allowed to verify other interpreting physicians' reports.

If the user is classified as Interpreting Staff, s/he will always be allowed to select another interpreting physician's name and reports if this field is set to 'YES'.

Choose from:

y YES

ALLOW VERIFYING OF OTHERS: NO

Select Personnel: <RET>

Clerical List
[RA PNLCLERK]
Interpreting Resident List
[RA PNLRES]
Interpreting Staff List
[RA PNLSTAFF]
Technologist List
[RA PNLTECH]

Each personnel type is printed in the same manner. Here is an example of the Interpreting Resident List. Keep in mind that ownership of the RA ALLOC key means the user has access to every Imaging Location, even though none may be assigned to him/her.

Prompt/User Response

Discussion

INTERPRETING R	ESIDENT LI	ST	AUG 18,	1994 15:18 PAGE	1
NAME	-		,	IMAGING LOC ACCESS	
_					
PROVIDERA, ONE Rad/Nuc Med	l Keys:		resident	X-RAY	
-					
PROVIDER, SEVEN Rad/Nuc Med		yes	resident	MAMMOGRAPHY	
-					
PROVIDERC, ONE		no	technolog	rist	
			resident	ULTRASOUND X-RAY	
Rad/Nuc Med	l Keys: RA	ALLOC		A-KAI	
_					
PROVIDERD, THRE Rad/Nuc Med		no	resident	X-RAY	

VII. Utility Files Maintenance Menu...

Supervisor Menu ... Utility Files Maintenance Menu ...

Complication Type Entry/Edit
Diagnostic Code Enter/Edit
Examination Status Entry/Edit
Film Type Entry/Edit
Label/Header/Footer Formatter
Major AMIS Code Entry/Edit
Nuclear Medicine Setup Menu ...
Order Entry Procedure Display Menu ...
Procedure Edit Menu ...
Reports Distribution Edit
Sharing Agreement/Contract Entry/Edit
Standard Reports Entry/Edit
Valid Imaging Stop Codes Edit

Overview:

New data may be added or edited via these options at any time throughout the use of the software. These options are therefore used for file maintenance as well as implementation. If this is not a new implementation of the software, you should take this time to review their contents for any necessary updating of the files. In particular, check those items that appear with

the icon.

This menu has corresponding options in the Maintenance Files Print Menu. As with the System Definition Menu edit options, you may want to use the Maintenance Files Print Menu to print a dated hard copy of the data entered in each file to keep your own record of changes.

Complication Type Entry/Edit [RA COMPEDIT]

This function allows the coordinator to enter or edit complication types that may occur during an exam. This option will produce the selection list you see when prompted for a complication type in the Case No. Exam Edit and Edit Exam by Patient options. These two options are used to associate complication types with specific examinations.

¹When there is a contrast medium reaction, the patient's record in the Adverse Reaction Tracking package will be automatically updated to indicate the patient has had a contrast media reaction. As a result, every time the patient is registered for a procedure which may involve contrast media, the user will be warned and s/he should then take the appropriate action. (The system uses CPT codes to determine whether contrast media may be used for a given procedure.)

Two complications are exported with the system, "Contrast Reaction" and "No Complication". Other site specific entries can be added using this function. You should not delete or change complication type names.

A list of complication types can be obtained through the Complication Type List option.

Prompt/User Response

Complication Type Entry/Edit

Select COMPLICATION TYPES: ??

Choose from:
 CONTRAST REACTION
 NAUSEA AND VOMITING
 NO COMPLICATION
 SEVERE MORBIDITY

This field contains the name of the type of complication that may be associated with an exam.

Initially, the system contains two complication types: 'Contrast Reaction', and 'No Complication'. Other complication types may be entered by the rad/nuc med coordinator at each site and can be 3 to 30 characters in length.

Discussion

Before adding a new complication, check to see what your site already has in the file. In this example, we see four complications. Pregnancy is not there so let's add that to the file.

¹Patch RA*5*45 March 2004 AMIS Code

```
Select COMPLICATION TYPES: PREGNANCY
Are you adding 'PREGNANCY' as a new COMPLICATION
TYPES (the 5TH)? Y (Yes)

COMPLICATION: PREGNANCY// <RET>

CONTRAST MEDIUM REACTION?: ?
Enter 'Y' for YES if the complication is a contrast medium reaction.
Choose from:
Y YES
N NO

CONTRAST MEDIUM REACTION?: N NO

Select COMPLICATION TYPES: <RET>
```

If the Contrast Medium Reaction field is set to YES, when a patient's exam complication field is edited to include this complication, the Adverse Reaction Tracking package interface will be triggered so that contrast media can be recorded as an allergy/reaction for this patient.

The Complication Report also uses this field to determine how many patient reactions were due to contrast media.

Note: All contrast media allergies and reactions are stored in the Adverse Reaction Tracking (ART) package rather than in the Radiology/Nuclear Medicine package. The automatic interface to the ART package will add contrast media to the patient's allergy records when a Rad/Nuc Med user enters this reaction through the Update Patient Record option or when entering complications marked as contrast media type reactions. When a procedure that uses a contrast medium is ordered or registered for a patient who has had a previous adverse reaction to the contrast medium, the contrast reaction message is displayed to the user.

There are other data items that most facilities require to complete the record or to submit an FDA report. It is suggested that one person from the Imaging department be responsible for following up on contrast media reactions. This person should be a member of the GMRA VERIFY DRUG ALLERGY and GMRA P&T COMMITTEE FDA mail groups so that s/he receives bulletins on patient reactions to contrast media. This person should be given the necessary ART key(s) and menu option(s) to enter additional information in a consistent manner and verify the data. Talk to IRM and your Clinical Coordinator to set this up.

Diagnostic Code Enter/Edit [RA DIAGEDIT]

This function allows the coordinator to enter/edit diagnostic codes used in reports. This file builds the selection list the user sees when prompted for a Diagnostic Code.

The following is the list of diagnostic codes that are loaded as part of the initialization process:

- 1. Normal
- 2. Minor Abnormality
- 3. Major Abnormality, No Attn. Needed
- 4. Abnormality, Attn. Needed
- 5. Major Abnormality, Physician Aware
- 6. Undictated Films Not Returned, 3 Days
- 7. Unsatisfactory/Incomplete Exam
- 8. Possible Malignancy, Follow-up Needed

An Inactive field was added so that a code could be inactivated and no longer appear as a selection for the user. Entries should not be deleted. The names of entries should also not be changed, once in production, since this could affect their meaning in historical data.

Entering a diagnostic code for each case serves a number of purposes:

- 1. Diagnostic codes can be used to generate alerts to the requesting clinician whenever an examination has an abnormal result. In this option, you may designate which diagnostic codes trigger these alerts. More than one code may be designated to trigger an alert message.
- 2. If the site decides to enter a diagnostic code for exams, database searches using VA FileMan can be performed using the Diagnostic Code field of the exam record in the search criteria. For example, the database can be searched for all CHEST PA&LAT procedures which have a diagnostic code of "Normal", or all that are "Unsatisfactory/Incomplete Exam".

- 3. Each case can be assigned **multiple** diagnostic codes. This allows you to not only indicate the **urgency** of the problem, but the **nature** of the problem as well, assuming you create the appropriate additional diagnostic codes. Diagnostic codes can also be used to flag teaching cases or other interesting cases if you create the appropriate diagnostic codes (example: "Teaching Case").
- 4. Diagnostic codes that have a YES in the "Print on Abnormal Report" field will result in exams with those codes appearing on an abnormal report that can be used for follow up and tracking.

Prompt/User Response

Diagnostic Code Enter/Edit

Select DIAGNOSTIC CODES: ??

Choose from:

1	NORMAL
2	MINOR ABNORMALITY
3	MAJOR ABNORMALITY, NO ATTN. NEEDED
4	ABNORMALITY, ATTN. NEEDED
5	MAJOR ABNORMALITY, PHYSICIAN AWARE
6	UNDICTATED FILMS NOT RETURNED, 3 DAYS
7	UNSATISFACTORY/INCOMPLETE EXAM
8	POSSIBLE MALIGNANCY, FOLLOW-UP NEEDED
9	POSSIBLE MALIGNANCY, FOLLOW-UP CRITICAL

This field contains a short diagnostic code to indicate the results of the interpreting physician's analysis of the images. Diagnostic codes can be between 3 and 40 characters in length and can be used as a basis for database searches. (i.e.. How many 'Normal' chest exams were performed during a specific time period?)

Eight diagnostic codes are included with the original package. Other site-specific diagnostic codes may be entered by the rad/nuc med coordinator at each site.

Discussion

We want to add a diagnostic code for overexposure. Let's first see what the file contains by entering "??" at the Diagnostic Codes prompt. Codes 1 - 8 were exported with the software. Code #9 was added by the site.

We want to add a code so we enter the name of the code in 3 - 40 characters.

Select DIAGNOSTIC CODES: **OVEREXPOSURE, MECH ERROR**Are you adding 'OVEREXPOSURE, MECH ERROR' as
a new DIAGNOSTIC CODES (the 10TH)? **Y** (Yes)

It will be the 10th code.

DIAGNOSTIC CODES NUMBER: 10// <RET>

DIAGNOSTIC CODE: OVEREXPOSURE, MECH ERROR

Replace <RET>

We have an opportunity to change the code number if we want and edit the name. Some facilities use special code # ranges for each Imaging Type (e.g., 50-59 for CT, 90-99 for MRI).

DESCRIPTION: OVEREXPOSURE DUE TO A MECHANICAL ERROR

We can give the code a longer description, 1 - 80 characters.

¹PRINT ON ABNORMAL REPORT: Y YES

Should exams given this code print on the Abnormal Exam Report? If yes, then enter a YES at this prompt.

GENERATE ABNORMAL ALERT?: N NO

Do you want this code to generate an Alert to the requesting physician? For this example, an alert is not needed.

INACTIVE: <RET>

For more information about alerts, refer to OE/RR and Kernel documentation.

Select DIAGNOSTIC CODES: <RET>

¹ Patch RA*5*34 Moved 'PRINT ON ABNORMAL REPORT' entry up to align with top of row.

Examination Status Entry/Edit [RA EXAMSTATUS]

The Examination Status file #72 contains various statuses for each Imaging Type. An examination status is a designated point from the time a case is registered in an Imaging Location (Waiting for Exam) through its completion (Complete). These statuses **do not** include request statuses such as Pending and Hold or results reporting statuses such as Draft and Verified. If your site decides to activate any Imaging Type other than General Radiology, then you must use this option to define how the statuses will work for each Imaging Type activated. If this is a virgin installation of the software, General Radiology should be treated no differently than the other Imaging Types.

When the system is initialized, the statuses will be set up, **but not activated**, for each exported Imaging Type (excluding General Radiology which is activated). **Note:** 6 statuses will exist in the file for each of the 9 Imaging Types (i.e., 54 total statuses, 6 per Imaging Type). This makes it possible for the status parameter setup to be tailored differently for each Imaging Type.

The fields you edit in this option generally fall into three basic categories:

- 1. Fields that indicate what data is mandatory for an exam to progress to the given status.
- 2. Which prompts should be asked of the user when trying to move an exam into the given status (applies to the Status Tracking option only).
- 3. Fields that indicate to the various workload report routines whether an exam in a particular exam status should be used in the compilation of that report.

Except for the Imaging Type, Examination Status, User Key Needed, and Default Next Status fields, all responses are answered with a YES/NO or bypassed.

This version of the software changes the appearance of a number of the prompts by using the Title for the prompt instead of the actual field Name.

Note: Only the Status Tracking of Exams option will follow the status parameters for the "Ask" prompts. Case No. Exam Edit or Edit Exam by Patient will not do this. All three options will follow the "Require" parameters to progress to the next status.

IMAGING TYPE - Statuses are defined for each Imaging Type. You begin in this option by telling the system which activated Imaging Type you want to edit. Remember, an Imaging Type is activated when one or more Imaging Locations and one or more procedures have been assigned to it.

EXAMINATION STATUS - This field is the current exam status being defined. When the software is installed, the required entries of the Examination Status file are created for each Imaging Type:

- a. ¹There is a status with an order of "1" (WAITING FOR EXAM) for each Imaging Type. This is the status that the module will place an exam in upon registration regardless of which data has been entered. **The name and order of the status must not be changed.**
- b. There is a status with an order of "9" (COMPLETE) for each Imaging Type. When an exam reaches this status, then the module considers the exam complete and the exam can no longer be called up by case number. The case number can then be reused. The name and order of the status must not be changed.
- c. There is a status with an order of "0" (CANCELLED) for each Imaging Type. Cancelled case numbers are also reused. The name and order of the status must not be changed.

If you intend to assign statuses other than the above three, a numbering system of 2 through 8 is used to assign each status. The higher the number assigned to a status, the more complete the exam should be.

The statuses Called for Exam, Examined, and Transcribed are also installed with the software and appear as part of the selection list. No order numbers are installed for Called for Exam, Examined, or Transcribed. They are available for your use but are not required for status tracking to function. If you apply an order number to any one of these three statuses, you will activate that status. We suggest you use the worksheet provided at the back of this manual to plan your strategy when implementing new Imaging Types and/or statuses.

-

¹ Patch RA*5*15

ORDER IN SEQUENCE OF STATUS PROGRESSION - This field contains the logical order that this exam status falls into. The order of the progression is set up by your site and must be a number from 2 through 8. If this field is left blank, then the Exam Status Tracking function will not use the status you selected.

SHOULD THIS STATUS APPEAR IN STATUS TRACKING - For any statuses that you activate, consider answering YES at this prompt. The Appears on Status Tracking field, if set to YES, will cause exams of that status to come up for edit in the Status Tracking option. If left blank or NO, exams of that status will be bypassed. If, for example, you answer NO or bypass this field for every status, you would get nothing on Status Tracking.

Unless there is some site-specific reason to do otherwise, all sites should have Waiting for Exam appear on status tracking, because a newly registered exam is automatically set to it. Examined is a status that sites often want to appear in Status Tracking. The system does not allow Complete and Cancelled statuses to appear on status tracking.

An example of when this field would be set to NO is for a status called Transcribed. Many times after a report is transcribed it may not be verified for a day or two and the number of Transcribed exams may get very large. If Transcribed was displayed in the status tracking function, then the system performance would suffer because all these exams would have to be processed and displayed.

USER KEY NEEDED TO MOVE AN EXAM TO THIS STATUS - When extra protection is desired on an exam status, use this field to define the security key. If this field is set to a valid security key, then the user will need that key in order to move an exam into this status.

If a key is desired, IRM will have to create it and allocate it to the appropriate users. It is the coordinator's responsibility to enter this key into this field after IRM creates it. The key name must start with "RA".

DEFAULT NEXT STATUS FOR EXAM - This field contains the name of the next status that the current status can advance to, after the parameters of that next status are met. The next status must be active, i.e., it must have data for the parameter Order in Sequence of Status Progression. The Exam Status Tracking function uses this field in determining which prompts should be displayed for edit by the user. For example, if the "Default Next Status" is Examined for a case in the Waiting for Exam status, the Status Tracking option will prompt for all the "Ask..." fields you have set to YES for the Examined status.

Note: If you do not enter a "Default next status for exam", the user will be asked to choose a next status during status tracking. Most sites choose a default next status that requires the most possible data entry at one session. For example, the default next status for Waiting for Exam might be Examined, and the default next status for Examined might Complete.

The more default next statuses you use, the more times a user will have to use the status tracking option to force the exam to go to the next status. However, if you go to the other extreme and make Waiting for Exam's default next status Complete, you may not get a fine enough granularity of status differentiation for status tracking.

A simple, but still viable setup would be the following:

Appear on Status Tracking:

Waiting for Exam

Examined

Complete

CAN AN EXAM BE CANCELLED WHILE IN THIS STATUS - While an exam is in this status, should the user be allowed to cancel an exam? Answer with YES to allow canceling and NO to disallow canceling. If no entry is made in this field for the current status, then the user will **not** be allowed to cancel the exam.

GENERATE EXAM ALERT FOR REQUESTING PHYSICIAN - This field contains a YES if the exam status should generate an alert to the requesting physician that the patient has been examined.

GENERATE EXAMINED HL7 MESSAGE - If an HL7 message should be generated

when a patient is examined, enter a YES at this prompt. The usual purpose of this field is to identify the status where images have been collected. The HL7 message recipient is usually a vendor-supplied PACS or other system that needs this HL7 message to trigger some action. There are other events that automatically trigger different HL7 messages to be sent by Rad/Nuc Med. These are fully discussed in the technical manual.

¹VistARAD CATEGORY – This field will only appear if your site is running the Imaging package and is only needed by sites using VistaRad for soft-copy reading of images. **Vistarad Category Values:** If this Examination Status is to be used for exams that will be read with VistaRad, then enter a value that corresponds to it from the following selections:

- W Waiting for Exam
- E Examined (assigned when the tech has completed image acquisition)
- D Dictated/Interpreted (assigned when the radiologist has completed the diagnostic reading)
- T Transcribed (assigned when the report has been transcribed)

Note: Not all status codes should be assigned a Vistarad Category value. Exam statuses other than the above selections should not be entered into this field.

Vistarad Types of Lists: The Vistarad software produces the following types of lists of Radiology exams. Each type includes exams based on the value of the Examination Status at the time the list is compiled according to the Vistarad Category that is assigned.

Vistarad Type of List	Rad/Nuc Med Exam Vistarad Category
Unread	Includes all exams with Category "E"
Dogant	Includes Cotogories "D" and "T"

Recent Includes Categories "D" and "T"
All Active Includes Categories "D", "T", and "E"
Pending Includes all exams with Category "W"

Exam Locking: When a radiologist reads exams from the Vistarad workstation, the software uses an exam "lock" to guarantee that the same exam is not read by two or more radiologists at the same time. In order for an exam to be locked, it must pass these tests:

- 1. It is not already locked by another radiologist, and
- 2. Its current Status must be "Examined"

The "Examined" Status is assigned at the time the technologist completes the image acquisition and performs the corresponding Case Edit operation in the Radiology package. The Vistarad Category field is used to let Vistarad know which Status values mean "Examined" and therefore which ones are eligible to be locked for reading.

•

¹ Patch RA*5*23 new field in the Examination Status file #72. Field added to worksheet (see Worksheets section).

Status Change Requirements

TECHNOLOGIST REQUIRED - This field may be set to YES to require that the name of the Technologist must be entered in order for an exam to be moved into the current status being defined.

RESIDENT OR STAFF REQUIRED - This field may be set to YES to require that an Interpreting Resident or Staff Physician be entered in order for an exam to be moved into the current status being defined.

DETAILED PROCEDURE REQUIRED - This field may be set to YES to require that a Detailed Procedure type be entered in order for an exam to be moved into the current status being defined. If the Division parameter requiring that a Detailed Procedure type be entered when ordering an exam is set to YES, then the Division requirement takes precedence over this requirement.

FILM ENTRY REQUIRED - This field may be set to YES to denote that Film Entry is required in order for an exam to be moved into the current status being defined.

DIAGNOSTIC CODE REQUIRED - This field may be set to YES to require that the Primary Diagnostic Code be entered in order for an exam to be moved into the current status being defined.

CAMERA/EQUIP/RM REQUIRED - This field may be set to YES to require that the Camera/Equip/Rm be entered in order for an exam to be moved into the current status being defined.

¹PROCEDURE MODIFIERS REQUIRED – This field should be set to YES if an exam requires at least 1 procedure modifier be entered in order for an exam to be moved into the current status being defined.

²CPT MODIFIERS REQUIRED - This field should be set to YES if an exam requires at least 1 CPT modifier be entered in order for an exam to be moved into the current status being defined.

¹ Patch RA*5*10 April 2000

² Patch RA*5*10 April 2000

(Radiopharmaceuticals requirements)

The following Radiopharmaceutical prompts only appear if you are editing the Imaging Types of Nuclear Medicine or Cardiology Studies.

RADIOPHARMS & DOSAGES REQUIRED - This field may be set to YES to require that the Radiopharmaceuticals and Dosages be entered in order for an exam of this Imaging Type to be moved into the current status being defined.

To override the requirement for radiopharmaceutical on an individual procedure, set the "Suppress Radiopharm Prompt" field of the Rad/Nuc Med Procedure file to "1" for that procedure using Procedure Enter/Edit option.

- ACTIVITY DRAWN REQUIRED (RADIOPHARM) This field may be set to YES to require that the Activity Drawn be required in order for an exam of this Imaging Type to be moved into the current status being defined.
- DRAWN DT/TIME AND PERSON REQUIRED (RADIOPHARM) This field may be set to YES to require that the Date/Time of the measure and the Person who measured the dose be entered in order for an exam of this Imaging Type to be moved into the current status being defined.
- ADMIN DT/TIME/PERSON REQUIRED (RADIOPHARM) This field may be set to YES to require that the Date/Time of Dose Administration and the Person who administered the dose be entered in order for an exam to be moved into the current status being defined.
- ROUTE/SITE REQUIRED (RADIOPHARM ADMINISTERED) This field may be set to YES to require that the Route of Administration and Site of Administration be entered in order for an exam of this Imaging Type to be moved into the current status being defined.

Note: If a selected route has no associated site, this requirement is ignored.

LOT NO. REQUIRED (RADIOPHARM) - This field may be set to YES to require that the Lot No. be entered in order for an exam of this Imaging Type to be moved into the current status being defined.



VOLUME/FORM REQUIRED (RADIOPHARM) - This field may be set to YES to

require that the Volume and Form be entered in order for an exam of this Imaging Type to reach the current status being defined.

WARNING: You must use the reporting feature of the system as a prerequisite for the following requirements:

REPORT ENTERED REQUIRED - This field may be set to YES to require that an exam have an entry in the Rad/Nuc Med Reports file #74. This does **not** mean that the report text field of the report record must be populated.

IMPRESSION REQUIRED - This field may be set to YES to require that a report Impression be entered in order for an exam to be moved into the current status being defined. For each Imaging Type, the Impression should be required at some point (probably for the Transcribed or Complete status). Since the report text can be purged, it is strongly advised that an Impression be mandatory with each exam.

VERIFIED REPORT REQUIRED - This field may be set to YES to require that the report be Verified in order for an exam to be moved into the current status being defined.

Status Tracking Functions

Depending on how each site sets up the following fields, the user may or may not be prompted for the data while using the exam status tracking function.

ASK FOR TECHNOLOGIST - If this field is set to YES, then the user will be prompted for the Technologist in order to update an exam to this status.

ASK FOR INTERPRETING PHYSICIAN - If this field is set to YES, then the user will be prompted for both the Resident and Staff Interpreting Physicians in order to update an exam to this status.

ASK FOR PROCEDURE - If this field is set to YES, then the user will be prompted for the Procedure in order to update an exam to this status.

ASK FOR FILM DATA - If this field is set to YES, then the user will be prompted for the Film Data in order to update an exam to this status.

ASK FOR DIAGNOSTIC CODE - If this field is set to YES, then the user will be prompted for the Primary Diagnostic Code in order to update an exam to this status. If a Primary Diagnostic code is entered, another prompt will appear asking for a Secondary Diagnostic Code.

ASK FOR CAMERA/EQUIP/RM - If this field is set to YES, then the user will be prompted for the Camera/Equip/Rm in order to update an exam to this status.

¹ASK FOR PROCEDURE MODIFIERS - If this field is set to YES, then the user will be prompted for Procedure Modifiers in order to update an exam to this status.

ASK FOR CPT MODIFIERS - If this field is set to YES, then the user will be prompted for CPT Modifiers in order to update an exam to this status.

ASK FOR USER CODE - If this field is set to YES, then the user will be prompted for the User Code in order to update an exam to this status.

ASK MEDICATIONS & DOSAGES - This parameter works in conjunction with a parameter on the Rad/Nuc Med Procedure file #71. If this field is set to YES, **and** the "Prompt for Meds" field of the procedure (in the Rad/Nuc Med Procedure file #71) is set to YES, then the user will be prompted for the Medication and Dose in order to update an exam of this Imaging Type to this status.

ASK MEDICATION ADMIN DT/TIME & PERSON If this field is set to YES, then the user will be prompted for the Date and Time the medication was administered and the Person who administered the dose in order to update an exam of this Imaging Type to this status. However, this will only happen if medications and dosages were entered.

(Radiopharmaceuticals functions)

The following Radiopharmaceutical prompts only appear if you are editing the Imaging Types of Nuclear Medicine or Cardiology Studies.

ASK RADIOPHARMACEUTICALS AND DOSAGES - If this field is set to YES, then the user will be prompted for the Radiopharmaceutical and Dose in order to update an exam of this Imaging Type to this status.

To override this parameter, and turn OFF the prompt for radiopharmaceuticals on an individual procedure, set the "Suppress Radiopharm Prompt" field of the Rad/Nuc Med Procedure file to "1" for that procedure using the Procedure Enter/Edit option.

_

¹ Patch RA*5*10 April 2000

- ASK ACTIVITY DRAWN (RADIOPHARMACEUTICAL) If this field is set to YES and a Radiopharmaceutical has been entered, then the user will be prompted for the Activity Drawn in order to update an exam of this Imaging Type to this status.
- ASK DRAWN DT/TIME & PERSON (RADIOPHARMACEUTICAL) If this field is set to YES and a Radiopharmaceutical has been entered, then the user will be prompted for the Date/Time of measure and Person who measured the dose when updating an exam of this Imaging Type to this status.
- ASK ADMIN DT/TIME & PERSON (RADIOPHARMACEUTICAL) If this field is set to YES and a Radiopharmaceutical has been entered, then the user will be prompted for the Date/Time of dose administration and the Person who administered the dose in order to update an exam of this Imaging Type to this status.
- ASK ROUTE/SITE OF ADMIN (RADIOPHARM) If this field is set to YES and a Radiopharmaceutical has been entered, then the user will be prompted for the Route of administration and Site of administration in order to update an exam of this Imaging Type to this status.
- ASK LOT NO. (RADIOPHARM) If this field is set to YES and a Radiopharmaceutical has been entered, then the user will be prompted for the Lot No. in order to update an exam of this Imaging Type to this status.
- ASK VOLUME/FORM (RADIOPHARM) If this field is set to YES and a

 Radiopharmaceutical has been entered, then the user will be prompted for Volume and
 Form in order to update an exam of this Imaging Type to this status.

Management Report Criteria

When the workload report routine is processing an exam that is in the exam status currently being defined, it will check the following associated parameters to determine whether to include the exam in the report compilation.

Cancelled and/or Complete exams will be included on any report addressed in this option if the field for Cancelled and/or Complete contains a YES. The only exception to this is the Delinquent Status Report. Cancelled and Complete exams will not be included on the

Delinquent Status Report even if you do enter YES in the "Delinquent Status Report?" field for the Cancelled and/or Complete statuses.

CLINIC REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Clinic Workload Report.

PTF BEDSECTION REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the PTF Bedsection Workload Report.

SERVICE REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Service Workload Report.

SHARING/CONTRACT REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Sharing/Contract Workload Report.

WARD REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Ward Workload Report.

FILM USAGE REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Film Usage Workload Report.

TECHNOLOGIST REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES to indicate that exams with the exam status currently being defined should be included in the Technologist Workload Report.

AMIS REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the AMIS Workload Report.

DETAILED PROCEDURE REPORT (WORKLOAD) SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Detailed Procedure Workload Report.

CAMERA/EQUIP/RM REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Camera/Equip/Rm Workload Report.

PHYSICIAN REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Physician Workload Report.

RESIDENT REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Resident Workload Report.

STAFF REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Staff Workload Report.

DELINQUENT STATUS REPORT SHOULD INCLUDE EXAMS OF THIS STATUS - Set this parameter to YES, to indicate that exams with the exam status currently being defined should be included in the Delinquent Status Report.

R	ac	lie	ol	oł	12	ır	n	12	ıc	e	u	ti	ca	al	I)(DS	a	g	e	T	i	ck	æ	t	P	r	ir	ıt	(20	r	ıt	r	0	l					
		_	-	_	-	_	_	-	_	_	-	_	_	-	_	_	-	_	-	_	_	-	_	_	_	_	_	_	_	_	_	_	_	_	_	_	_	-	_	_	_

Note: The following radiopharmaceutical prompt only appears if you are editing the Imaging Type Nuclear Medicine or Cardiology Studies.

PRINT DOSAGE TICKET WHEN EXAM REACHES THIS STATUS - If an exam of this Imaging Type that uses radiopharmaceuticals reaches this status and this field is set to YES, a dosage ticket will be printed if the following are also true:

- 1. A valid print device has been assigned to the Imaging Location.
- 2. YES has been answered for the Radiopharms/Dosages Required field at this (or lower) status for this Imaging Type.

Exported Setup

The following setup is exported for all Imaging Types. Any setup already done by your site will be left untouched. Recommended answers are shown for all fields, but you will have to fill in the unanswered fields. The fields shaded, only apply to Imaging Type file entries whose "Radiopharmaceuticals Used?" field is set to YES (Nuclear Medicine and/or Cardiology Studies). This field is not editable.

Note: The Statuses are populated, but not activated.

If Imaging Types other than General Radiology are activated at your medical center, you will have to assign order numbers to any additional statuses you will be using, indicate whether the statuses should appear on Status Tracking, assign a default next status for each, and if you feel it is necessary, review all the other parameters associated with the statuses.

Field: Statuses:

	2000000					
		Waiting	Called			
	Can-	for	for Exam	Exam-	Tran-	Com-
	celled	Exam		ined	scribed	plete
Order in sequence of status progression	0	1				9
Should this Status appear in Status Tracking		Y				
User Key needed to move an exam to this status						
Default next status for exam						
Can an exam be cancelled while in this status	Y	Y	Y			
Generate exam alert for requesting physician						
Generate Examined HL7 Message						
VISTARAD Category						

		Waiting	Called			
	Can-	for	for Exam	Exam-	Tran-	Com-
	celled	Exam		ined	scribed	plete
Data Requirements for this Status						
Technologist Required			Y	Y	Y	Y
Resident or Staff Required					Y	Y
Detailed Procedure Required				Y	Y	Y
Film Entry Required				Y	Y	Y
Diagnostic Code Required					Y	Y
Camera/Equip/RM Required					Y	Y
Radiopharms & Dosages required						
Activity Drawn required (Radiopharm)						
Drawn Dt/Time and Person required (Radiopharm)						
Admin Dt/Time/Person required (Radiopharm)						
Route/Site required (Radiopharm administered)						
Lot No. required (Radiopharm)						
Volume/Form required (Radiopharm)						
Report Entered required					Y	Y
Impression required						Y
Verified Report required						Y
To Move to this Status						
Ask for Technologist			Y	Y		
Ask for Interpreting Physician					Y	
Ask for Procedure				Y		
Ask for Film Data				Y		
Ask for Diagnostic Code					Y	
Ask for Camera/Equip/RM				Y		
Ask for User Code						
Ask Medications & Dosages						
Ask Medication Admin Dt/Time & Person						
Ask Radiopharmaceuticals and Dosages				Y		
Ask Activity Drawn (Radiopharmaceutical)						
Ask Drawn Dt/Time & Person (Radiopharm)						
Ask Admin Dt/Time & Person (Radiopharm)						
Ask Route/Site of Admin (Radiopharm)						
Ask Lot No. (Radiopharm)						
Ask Volume/Form (Radiopharm)						

		Waiting	Called			
	Can-	for	for Exam	Exam-	Tran-	Com-
	celled	Exam		ined	scribed	plete
Exams in this Status should appear on						
these Reports						
Clinic Report		Y	Y	Y	Y	Y
PTF Bedsection Report		Y	Y	Y	Y	Y
Service Report (Workload)		Y	Y	Y	Y	Y
Sharing/Contract Report (Workload)		Y	Y	Y	Y	Y
Ward Report (Workload)		Y	Y	Y	Y	Y
Film Usage Report (Workload)		Y	Y	Y	Y	Y
Technologist Report (Workload)		Y	Y	Y	Y	Y
AMIS Report (Workload)		Y	Y	Y	Y	Y
Detailed Procedure Report (Workload)		Y	Y	Y	Y	Y
Camera/Equip/Rm Report		Y	Y	Y	Y	Y
Physician Report		Y	Y	Y	Y	Y
Resident Report		Y	Y	Y	Y	Y
Staff Report		Y	Y	Y	Y	Y
Delinquent Status Report		Y	Y	Y	Y	
Radiopharmaceutical Dosage Ticket Print						
Control						
Print Dosage Ticket when exam reaches this status						

The following example takes you through entries for Complete when setting up the Imaging Type, Nuclear Medicine. Remember, the prompts for radiopharmaceuticals will only appear when the Imaging Type is Nuclear Medicine or Cardiology Studies.

```
Examination Status Entry/Edit
Select an Imaging Type: NUCLEAR MEDICINE
Select an Examination Status: COMPLETE
                                          Imaging Type: NUCLEAR MEDICINE
                                  Order: 9
        ...OK? Yes// <RET> (Yes)
Name of Current Exam Status: COMPLETE// <RET>
Order in sequence of status progression: 9// <RET>
Should this Status appear in Status Tracking ?: NO// <RET>
User Key needed to move an exam to this status: <RET>
Default next status for exam: <RET>
Can an exam be cancelled while in this status ?: NO// <RET>
Generate exam alert for requesting physician ?: NO// <RET>
Generate Examined HL7 Message:
    Status Change Requirements
    Please indicate which of the following is required
     in order to place an exam into the 'COMPLETE' status:
```

Utility Files Maintenance Menu

```
Technologist Required ?: YES// <RET>
Resident Or Staff Required ?: YES// <RET>
Detailed Procedure required ?: YES// <RET>
Film Entry Required ?: NO// <RET>
Diagnostic Code Required ?: NO// <RET>
Camera/Equip/Rm Required ?: NO// <RET>
     (Radiopharmaceutical requirements)
Radiopharms & Dosages required ?: YES// <RET>
Activity Drawn required (Radiopharm) ?: YES// <RET>
Drawn Dt/Time and Person required (Radiopharm) ?: YES// <RET>
Admin Dt/Time/Person required (Radiopharm) ?: YES// <RET>
Route/Site required (Radiopharm administered) ?: YES// <RET>
Lot No. required (Radiopharm) ?: YES// <RET>
Volume/Form required (Radiopharm) ?: YES// <RET>
   WARNING: You must use the reporting feature of the system as
           a prerequisite for the following requirements:
Report Entered required ?: YES// <RET>
Impression required ?: YES// <RET>
Verified Report required ?: YES// <RET>
     Status Tracking Functions
     ______
     Please indicate which of the following should be asked when
     changing an exam's status to 'COMPLETE' while using
     the 'status tracking' feature:
Ask for Technologist ?: YES// <RET>
Ask For Interpreting Physician ?: YES// <RET>
Ask for Procedure ?: YES// <RET>
Ask For Film Data ?: NO// <RET>
Ask For Diagnostic Code ?: YES// <RET>
Ask For Camera/Equip/Rm ?: NO// <RET>
Ask for User Code ?: NO// <RET>
Ask Medications & Dosages ?: YES// <RET>
Ask Medication Admin Dt/Time & Person ?: YES// <RET>
     (Radiopharmaceutical functions)
Ask Radiopharmaceuticals and Dosages ?: YES// <RET>
Ask Activity Drawn (Radiopharmaceutical) ?: YES// <RET>
Ask Drawn Dt/Time & Person (Radiopharmaceutical) ?: YES// <RET>
Ask Admin Dt/Time & Person (Radiopharmaceutical) ?: YES// <RET>
Ask Route/Site of Admin (Radiopharm) ?: YES// <RET>
Ask Lot No. (Radiopharm) ?: YES// <RET>
Ask Volume/Form (Radiopharm) ?: YES// <RET>
     Management Report Criteria
     Please indicate which of the following workload reports should
     include exams with the 'COMPLETE' status:
```

Clinic Report should include exams of this status ?: YES// <ret> PTF Bedsection Report should include exams of this status ?: YES// <ret> Service Report (Workload) should include exams of this status ?: YES// <ret> Sharing/Contract Report (Workload) should include exams of this status ?: YES// <ret> Ward Report (Workload) should include exams of this status ?: YES// <ret> Film Usage Report (Workload) should include exams of this status ?: YES// <ret> Technologist Report (Workload) should include exams of this status ?: YES// <ret> AMIS Report (Workload) should include exams of this status ?: YES// <ret> Detailed Procedure Report (Workload) should include exams of this status ?: YES// <ret> Camera/Equip/Rm Report should include exams of this status ?: YES// <ret> Physician Report should include exams of this status ?: YES// <ret> Staff Report should include exams of this status ?: YES// <ret> Delinquent Status Report should include exams of this status ?: YES// <ret></ret></ret></ret></ret></ret></ret></ret></ret></ret></ret></ret></ret></ret>
Radiopharmaceutical Dosage Ticket Print Control
Print Dosage Ticket when exam reaches this status ?: NO// <ret></ret>
Select an Examination Status: <ret></ret>
Checking order numbers and Default Next Status used for status progression within : NUCLEAR MEDICINE
Required order numbers are in place.
Data Inconsistency Report For Exam Statuses
DEVICE: HOME// <ret> (enter a printer or <ret> to display on screen)</ret></ret>
Page: 1 Date: Sep 17, 1997
Data Inconsistency Report For Exam Statuses
Checking order numbers and Default Next Status used for status progression within : NUCLEAR MEDICINE
Required order numbers are in place.
Exam Status Data Inconsistencies Not Found.

The following example shows a setup with inconsistencies as demonstrated in the Data Inconsistency Report that follows the setup:

```
Select an Imaging Type: NUCLEAR MEDICINE
Select an Examination Status: COMPLETE Imaging Type: NUCLEAR MEDICINE
                                  Order: 9
Name of Current Exam Status: COMPLETE// <RET>
Order in sequence of status progression: 9// <RET>
Should this Status appear in Status Tracking ?: NO// <RET>
User Key needed to move an exam to this status: <RET>
Default next status for exam: <RET>
Can an exam be cancelled while in this status ?: NO// <RET>
Generate exam alert for requesting physician ?: NO// <RET>
Generate Examined HL7 Message: <RET>
    Status Change Requirements
    _____
    Please indicate which of the following is required
    in order to place an exam into the 'COMPLETE' status:
Technologist Required ?: YES// <RET>
Resident Or Staff Required ?: YES// <RET>
Detailed Procedure required ?: YES// <RET>
Film Entry Required ?: YES// <RET>
Diagnostic Code Required ?: YES// <RET>
Camera/Equip/Rm Required ?: YES// <RET>
     (Radiopharmaceutical requirements)
Radiopharms & Dosages required ?: NO// <RET>
Activity Drawn required (Radiopharm) ?: NO// <RET>
Drawn Dt/Time and Person required (Radiopharm) ?: NO// <RET>
Admin Dt/Time/Person required (Radiopharm) ?: NO// <RET>
Route/Site required (Radiopharm administered) ?: NO// <RET>
Lot No. required (Radiopharm) ?: NO// <RET>
Volume/Form required (Radiopharm) ?: NO// <RET>
  WARNING: You must use the reporting feature of the system as
           a prerequisite for the following requirements:
Report Entered required ?: YES// <RET>
Impression required ?: NO// <RET>
Verified Report required ?: YES// <RET>
    Status Tracking Functions
    Please indicate which of the following should be asked when
     changing an exam's status to 'COMPLETE' while using
     the 'status tracking' feature:
```

```
Ask for Technologist ?: YES// <RET>
Ask For Interpreting Physician ?: YES// <RET>
Ask for Procedure ?: NO// <RET>
Ask For Film Data ?: YES// <RET>
Ask For Diagnostic Code ?: YES// <RET>
Ask For Camera/Equip/Rm ?: NO// <RET>
Ask for User Code ?: NO// <RET>
Ask Medications & Dosages ?: NO// <RET>
Ask Medication Admin Dt/Time & Person ?: NO// <RET>
     (Radiopharmaceutical functions)
Ask Radiopharmaceuticals and Dosages ?: NO// <RET>
Ask Activity Drawn (Radiopharmaceutical) ?: NO// <RET>
Ask Drawn Dt/Time & Person (Radiopharmaceutical) ?: NO// <RET>
Ask Admin Dt/Time & Person (Radiopharmaceutical) ?: NO// <RET>
Ask Route/Site of Admin (Radiopharm) ?: NO// <RET>
Ask Lot No. (Radiopharm) ?: NO// <RET>
Ask Volume/Form (Radiopharm) ?: NO// <RET>
    Management Report Criteria
     _____
     Please indicate which of the following workload reports should
     include exams with the 'COMPLETE' status:
Clinic Report should include exams of this status ?: YES// <RET>
PTF Bedsection Report should include exams of this status ?: YES// <RET>
Service Report (Workload) should include exams of this status ?: YES// <RET>
Sharing/Contract Report (Workload) should include exams of this status ?:
YES// <RET>
Ward Report (Workload) should include exams of this status ?: YES// <RET>
Film Usage Report (Workload) should include exams of this status ?: YES//
<RET>
Technologist Report (Workload) should include exams of this status ?: YES//
AMIS Report (Workload) should include exams of this status ?: YES// <RET>
Detailed Procedure Report (Workload) should include exams of this status ?:
YES// <RET>
Camera/Equip/Rm Report should include exams of this status ?: YES// <RET>
Physician Report should include exams of this status ?: YES// <RET>
Resident Report should include exams of this status ?: YES// <RET>
Staff Report should include exams of this status ?: YES// <RET>
Delinquent Status Report should include exams of this status ?: NO// <RET>
     Radiopharmaceutical Dosage Ticket Print Control
```

Print Dosage Ticket when exam reaches this status ?: NO// <RET>

Select an Examination Status: <RET>

When exiting this option, you receive a report on the consistency of the data you entered. In the example shown here, several inconsistencies are noted and may or may not require correction for the status tracking function to work properly.

Checking order numbers
and Default Next Status used for status progression _____
within : NUCLEAR MEDICINE

Required order numbers are in place.

Data Inconsistency Report For Exam Statuses

DEVICE: HOME// (Enter a Device)

		Page: 1 Date: May 20, 1997
	Data Inconsistency Report For Exam S	
	on reaching Complete : NUCLEAR MEDICINE	
	g are permissible, but could lead to es when prompts are not answered in l	
STATUS	PROMPT	DATA
EXAMINED	ASK FOR INTERPRETING PHYS? RESIDENT OR STAFF REQUIRED?	Y N
EXAMINED	ASK FOR FILM DATA? FILM ENTRY REQUIRED?	Y N
EXAMINED	ASK FOR DIAGNOSTIC CODE? DIAGNOSTIC CODE REQUIRED?	Y N
EXAMINED	ASK FOR CAMERA/EQUIP/RM? CAMERA/EQUIP/RM REQUIRED?	Y N
TRANSCRIBED	ASK FOR CAMERA/EQUIP/RM? CAMERA/EQUIP/RM REQUIRED?	Y N
COMPLETE	ASK FOR CAMERA/EQUIP/RM? CAMERA/EQUIP/RM REQUIRED?	N Y
	ng order numbers	
	Fault Next Status used for status pro : NUCLEAR MEDICINE	gression
Required orde	er numbers are in place.	
	Exam Status Data Inconsistencies Not	Found.

In the above sample report, if status tracking were used to enter data, the parameters would potentially allow users to leave the Camera/Equip/Rm field blank in every Status Tracking edit session and never get a prompt at the final edit session. This would result in the exam staying in the Transcribed status. A better solution would be to change the "Ask For Camera/Equip/Rm" field to YES on the Complete status or change the "Camera/Equip/Rm Required" field on the Transcribed and Examined statuses to YES.

Example Status Tracking Sessions

The following sessions show how four statuses are affected by changes in data. These are the assigned statuses for the examples:

Status	Order	Default Next Status
WAITING FOR EXAM	1	TRANSCRIBED
CALLED FOR EXAM	3	TRANSCRIBED
TRANSCRIBED	7	COMPLETE
COMPLETE	9	

Backing down to lower status: Deleting data previously entered may cause an exam status to back down to a lower status.

```
Enter Case #, Status, (N)ext status, '^' to Stop: NEXT// 521
Case # being tracked: 521
Name: TEST, FEMALE Case # : 20
Division: HINES CIO FIELD OFFICE Location: ULTRASOUND A
Procedure: ECHOGRAM ABDOMEN COMPLETE (Detailed)
   ***** Old Status: CALLED FOR EXAM
   ***** New Status: TRANSCRIBED
Do you wish to continue? YES// <ret>
    PRIMARY INTERPRETING RESIDENT: <ret>
    PRIMARY INTERPRETING STAFF: <ret>
    Select FILM SIZE: <ret>
    PRIMARY DIAGNOSTIC CODE: <ret>
    PRIMARY CAMERA/EOUIP/RM: <ret>
    Select TECHNOLOGIST: PROVIDER, NINE// @
   SURE YOU WANT TO DELETE? Y (Yes)
    Select TECHNOLOGIST: <ret>
    PREGNANT AT TIME OF ORDER ENTRY: NO
    PREGNANCY SCREEN: u Patient is unable to answer or is unsure
    PREGNANCY SCREEN COMMENT: Patient sent to lab for test
   No 'technologist' entered for this exam.
No 'film data' entered for this exam.
No 'camera/equip/room' entered for this exam.
No 'report' entered for this exam.
   No 'report'
                                       entered for this exam.
       ...exam status backed down to 'WAITING FOR EXAM'
       STATUS CHANGE DATE/TIME: JUN 10,1997@10:06// <ret>
... Status backed down for case #: 521
```

Upgrading status but not to Default Next Status: A case may be upgraded to a status that is not its DEFAULT NEXT STATUS, if that case's data satisfy all required data for a status that's higher than the old status, but lower than the default next status.

```
Case # being tracked: 521
Name:PATIENT,ONE Case # : 521
Division: HINES CIO FIELD OFFICE Location: ULTRASOUND A
Procedure: ECHOGRAM ABDOMEN COMPLETE (Detailed)
                    **** Old Status: WAITING FOR EXAM
   ***** New Status: TRANSCRIBED default next is not CALLED FOR EXAM
Do you wish to continue? YES// <ret>
     PRIMARY INTERPRETING RESIDENT: <ret>
    PRIMARY INTERPRETING STAFF: <ret>
     Select FILM SIZE: <ret>
    PRIMARY DIAGNOSTIC CODE: <ret>
    PRIMARY CAMERA/EQUIP/RM: <ret>
    Select TECHNOLOGIST: PROVIDER, NINE
    Select TECHNOLOGIST: <ret>
   No 'film data' entered for this exam.
No 'camera/equip/room' entered for this exam.
No 'report' entered for this exam.
No 'camera/equip/room' entered for this exam.
No 'report' entered for this exam.
       ...though upgraded, new status level (CALLED FOR EXAM)
      is not as high as the desired level (TRANSCRIBED)
       STATUS CHANGE DATE/TIME: JUN 10,1997@10:06// <ret>
...Status successfully changed for case #: 521
```

Status unchanged: A case's status will remain the same, if no data had been added to, nor deleted from it.

```
Name: PATIENT, ONE
                                         Case # : 521
Name: PATIENT,ONE Case # : 521
Division: HINES CIO FIELD OFFICE Location: ULTRASOUND A
Procedure: ECHOGRAM ABDOMEN COMPLETE (Detailed)
   ***** Old Status: CALLED FOR EXAM
   ***** New Status: TRANSCRIBED
Do you wish to continue? YES// <ret>
    PRIMARY INTERPRETING RESIDENT: <ret>
    PRIMARY INTERPRETING STAFF: PROVIDER, TEN// <ret>
    Select SECONDARY INTERPRETING STAFF: <ret>
    Select FILM SIZE: <ret>
    PRIMARY DIAGNOSTIC CODE: NORMAL// <ret>
    Select SECONDARY DIAGNOSTIC CODE: <ret>
    PRIMARY CAMERA/EQUIP/RM: <ret>
    Select TECHNOLOGIST: PROVIDER, NINE // <ret>
   No 'film data' entered for this exam.
No 'camera/equip/room' entered for this exam.
      ...exam status not changed??
...Status unchanged for case #: 521
```

Qualifies for higher status than Default Next Status: A case may have enough data to qualify for a higher status than its DEFAULT NEXT STATUS, but Status Tracking will not allow the case to be upgraded to a status higher than the DEFAULT NEXT STATUS. If this should occur, a message will be displayed to remind the user to re-edit this case.

```
Exam Status Tracking Module
Division: HINES CIO FIELD OFFICE
Date : 06/10/97 11:16 AM
Locations: ULTRASOUND A

Case # Date Patient Procedure Equip/Rm
----- 9:55 AM PATIENT, ONE ECHOGRAM ABDOMEN COMPLETE

Enter Case #, Status, (N) ext status, '^' to Stop: NEXT// 521

Case # being tracked: 521

Name: PATIENT, ONE Case # : 521
Division: HINES CIO FIELD OFFICE Location: ULTRASOUND A
Procedure: ECHOGRAM ABDOMEN COMPLETE (Detailed)
***** Old Status: CALLED FOR EXAM
***** New Status: TRANSCRIBED
```

```
Do you wish to continue? YES// <ret>
    PRIMARY INTERPRETING RESIDENT: <ret>
    PRIMARY INTERPRETING STAFF: PROVIDER, TEN// <ret>
    Select SECONDARY INTERPRETING STAFF: <ret>
    Select FILM SIZE: 2
     1 2X12
     2 2X4
     3 2X6
CHOOSE 1-3: 1
     AMOUNT (#films or cine ft): 1
    Select FILM SIZE: <ret>
    PRIMARY DIAGNOSTIC CODE: NORMAL// <ret>
    Select SECONDARY DIAGNOSTIC CODE: <ret>
    PRIMARY CAMERA/EQUIP/RM: PORTABLE PORTABLE
    Select TECHNOLOGIST: PROVIDER, NINE // <ret>
This case also qualifies for higher status(es):
    COMPLETE
Since Status Tracking can only upgrade one status level at a time,
please edit this exam again.
      STATUS CHANGE DATE/TIME: JUN 10,1997@11:18//
...Status successfully changed for case #: 521
 Exam Status Tracking Module Division: HINES CIO FIELD OFFICE Date : 06/10/97 11:18 AM Status : CALLED FOR EXAM
Exam Status Tracking Module
Locations: ULTRASOUND A
                                           Procedure
Case # Date Patient
                                                                    Equip/Rm
Enter Case #, Status, (N) ext status, '^' to Stop: NEXT// N
 Exam Status Tracking Module Division: HINES CIO FIELD OFFICE Date : 06/10/97 11:18 AM Status : TRANSCRIBED
Locations: ULTRASOUND A
Case # Date Patient Procedure
                                                                  Equip/Rm
-----
     06/20/94 PATIENT, THREE ZZLINE ULTRASOUND PORTABLE 01/19/95 PATIENT, FOUR ECHOGRAM ABDOMEN COMPLETE ROOM1 9:55 AM PATIENT, ONE ECHOGRAM ABDOMEN COMPLETE PORTABLE
20
169
521
Enter Case #, Status, (N) ext status, '^' to Stop: NEXT// 521
Case # being tracked: 521
```

```
Name: PATIENT, ONE
                                    Case # : 521
Division: HINES CIO FIELD OFFICE Location: ULTRASOUND A
Procedure: ECHOGRAM ABDOMEN COMPLETE (Detailed)
   ***** Old Status: TRANSCRIBED
   ***** New Status: COMPLETE
Do you wish to continue? YES// <ret>
    PRIMARY INTERPRETING RESIDENT: <ret>
    PRIMARY INTERPRETING STAFF: PROVIDER, TEN // <ret>
    Select SECONDARY INTERPRETING STAFF: <ret>
    Select FILM SIZE: 2X12// <ret>
     FILM SIZE: 2X12// <ret>
     AMOUNT (#films or cine ft): 1// <ret>
    Select FILM SIZE: <ret>
    PRIMARY DIAGNOSTIC CODE: NORMAL// <ret>
    Select SECONDARY DIAGNOSTIC CODE: <ret>
    PRIMARY CAMERA/EQUIP/RM: PORTABLE// <ret>
    Select TECHNOLOGIST: PROVIDER, NINE // <ret>
      STATUS CHANGE DATE/TIME: JUN 10,1997@11:18// <ret>
...Status successfully changed for case #: 521
   ...will now designate request status as 'COMPLETE'...
    Visit credited.
```

... request status successfully updated.

Film Type Entry/Edit [RA FILMEDIT]

This function allows you to edit existing film types and enter new film types. These film types build the list of choices which the technologist can choose from when selecting a Film Size during the Status Tracking of Exams, Edit Exam by Patient, and Case No. Exam Edit options to indicate which films were used during a procedure.

An entry should be an identifying name or size. To speed up processing, assign each film type a name with unique characters at the beginning of the entry. If your medical center wants to get a report showing exams that did not use films, an entry of NO FILMS on this list is an appropriate choice. Once a film size has been entered through this option, it cannot be deleted; however, it can be inactivated. Also, enter a corresponding wasted film for each film size.

Prompt/User Response

Discussion

Film Type Entry/Edit

Select FILM SIZES: 10X12
Are you adding '10X12' as a new FILM SIZES (the 46TH)? Y (Yes)

FILM: 10X12// <RET>

IS THIS A 'CINE' FILM SIZE?: NO

Enter ?? if you want a list of film sizes already entered at your site.
The entry should be 3-30 characters in length.

If this field contains a YES, it indicates that this is cine film. Cine film is treated differently in compiling the AMIS reports and the Film Usage Report.

If this field contains a YES, it indicates that this film is used for fluoro. Exams where this film size is entered are assumed to be fluoro not requiring the use of films. Leaving this field blank is the same as

answering No.

FLUORO ONLY?: <RET>

INACTIVATION DATE: <RET>

If this film size entry should no longer be needed, enter a date here to inactivate it. That will remove it from the user selections. If the film type is active, there should be **no entry** in this field.

WASTED FILM?: <RET>

Bypass this field, as shown here, if this film size is NOT a wasted type. If this field contains a YES, it indicates that the film will be considered a wasted film type and will be excluded from the "Film Usage Report". Instead, it will be included on the "Wasted Film Report". If a film size is already associated with a wasted film size, the former cannot also be set to a wasted film size by answering "YES" at this prompt.

Let's enter a wasted film for the 10X12.

Select FILM SIZES: W-10X12

Are you adding 'W-10X12' as a new FILM SIZES (the 46TH)? \mathbf{Y} (Yes)

FILM: W-10X12// <RET>

IS THIS A 'CINE' FILM SIZE?: NO

FLUORO ONLY?: <RET>
INACTIVATION DATE: <RET>

We suggest you distinguish all wasted film in the same way, such as using "W-" before the name of its counterpart. WASTED FILM?: Y YES

Select FILM SIZES: <RET>

ANALOGOUS UNWASTED FILM SIZE: 10x12

Entering "10X12" here means that when a technician enters "W-10X12" on a case, it is the "10X12" size that has been wasted.

Enter the corresponding

film.

When you exit this option, something similar to the following is displayed:

Relationship between 'Unwasted Film Size' and 'Corresponding Wasted FilmSize':

'Unwasted Film Size' 'Corresponding Wasted Film Size'

10X12 W-10X12
11X14 W-11X14
14X17 W-14X17
2X12 W-2X12
2X4 W-2X4
2X6 unassociated with a 'Wasted Film' type
....

Enter RETURN to continue or '^' to exit: <RET>

Label/Header/Footer Formatter [RA FLASHFORM]

The Label/Header/Footer Formatter option is used to enter and edit flash card, jacket and examination label formats, and headers and footers. A format specifies the data you want printed and how you want it printed on each type of card or label. Headers are the line or lines of information at the top of a result report, usually including the patient name, social security number and similar information. Footers are the information printed at the bottom of a report, for example, patient name, examination date, page number, etc.

The header, footer, and label format specifications that you enter by using the Label/Header/Footer Formatter option are stored in the Label Print Fields file #78.7. Various programs use these specifications when printing exam results reports and the labels that are printed at the time an exam is registered.

On the following page is a list of the choices you have for formatting flash cards, exam labels, jacket labels, headers and footers. When the data depends on the availability of a specific requirement, we've also included what that requirement is.

Label Print Fields: Requirement:

Label Print Fields:	Requirement:		
Age of Patient	Date of Birth entered in MAS		
Attending Physician at Order	Must be entered in MAS		
Attending Physician Current	Must be entered in MAS		
Case Number	Case must be registered		
Case Number Barcode	Case must be registered, printer hardware		
	must be capable of printing barcodes, printer		
	device definition must be set up by IRM to		
	support barcode printing		
Current Date/Time			
Current Patient Location			
Date of Birth	Date of Birth entered in MAS		
Date of Exam	Case must be registered		
Date Report Entered	After report is entered		
Date Reported	After report is entered		
Diagnostic Code	After diagnostic codes are entered		
Exam Modifiers	Modifiers must be entered		
Free Text			
Hospital Division (Institution)			
Imaging Location	Exam must be registered		
Last Time this Exam Performed			
Last Visit for any Exam			
Long Case Number	Case must be registered		
Long Case Number Barcode			
Name of Patient			
Page Number			
	First line of street address must be entered in		
Patient Address Line 1	MAS		
	Second and third lines of address must be		
Patient Address Line 2	entered in MAS		
1 aucht Address Lille 2			
	City, State, and Zip must be entered in MAS		
Patient Address Line 3			
Primary Physician at Order	Must be entered in MAS		
Primary Physician Current	Must be entered in MAS		
Procedure			
	·		

Procedure Barcode, 30 chars	Printer hardware must be capable of printing		
,	barcodes, printer device definition must be set		
	up by IRM to support barcode printing		
Procedure Barcode, 60 chars	Printer hardware must be capable of printing		
	barcodes, printer device definition must be set		
	up by IRM to support barcode printing		
Report Status	Report must be entered		
Request Entered Date/Time			
Requesting Area			
Requesting Physician			
Resident	Must be entered		
Resident Signature	Report must be entered and verified by resident		
	Must be verified by resident, then looks for		
Resident Signature Name	Signature Block Printed Name or defaults to		
	name from the .01 field of File 200.		
Service			
Sex of Patient	Must be entered in MAS		
SSN of Patient	Must be entered in MAS		
SSN of Patient Barcode	Printer hardware must be capable of printing		
	barcodes, printer device definition must be set		
	up by IRM to support barcode printing		
	Printer hardware must be capable of printing		
¹ SSN of Patient Barcode-No Dash	barcodes, printer device definition must be set		
	up by IRM to support barcode printing		
Staff Radiologist	Primary Staff must be entered		
Staff Signature	Report must be entered and verified by staff		
Staff Signature Name	Must be verified by a staff member, then looks		
	for Signature Block Printed Name or defaults		
	to name from the .01 field of File 200.		
Technologist	Must be entered in Exam Edit option		
Usual Patient Category	Must be entered via Update Patient Record		
Verified Date	Report must be verified		
Verifying Electronic Signature	Report must be verified		

_

¹ Patch RA*5*15 May 2000 NOIS: CIH-0999-43119 Added label print field.

Verifying Radiologist	Report must be verified		
Verifying Signature	Report must be verified		
Verifying Signature Block	Report must be verified		
Verifying Signature Name	Must be verified, then uses Signature Block Name.		
Verifying Signature Title	Must be verified, then uses Signature Block Title from the New Person file #200.		

When you specify formats you can also enter a "TITLE" to be printed for each field. When the field is printed, the title will precede the actual data. The positioning of each field on the page is site-specific. If no title is entered at this step, a default title is filled in by the system. If no title is desired, then the word "NONE" should be entered when you are prompted for "TITLE". This might be useful when printing the patient's address since the default would be PATIENT ADDRESS LINE 1, etc.

You will be prompted for row and column positions. These positions determine how your flash card or label is formatted. There can be a maximum of 66 rows from top to bottom of a page in a format but you must limit the number of rows used to fit on your labels and/or results report headers and footers. There can be a maximum of 80 columns across a page in a format, but labels usually allow less than 80 columns.

After you have selected all of the fields you need, you will be asked the number of rows in the format. Be sure to answer with the highest row number you have assigned. This is necessary so that the print program will know how much space to allow on a page for headers and footers.

After a format has been entered or edited, the system will display a test version of the format, so that you can see what it will actually look like and make any necessary changes.

When you enter Location Set-up Parameters, the Imaging Locations file #79.1 stores information telling which specifications in the Label Print Fields file #78.7 should be used to print the following:

- 1. Default Flash Card Format
- 2. Default Exam Label Format
- 3. Default Jacket Label Format
- 4. Default Report Header Format
- 5. Default Report Footer Format

These parameters are set using the Location Parameter Set-up option.

Note: When designing a card or label for a specific purpose, keep in mind which data will be available at the time the card or label is printed. For instance, if the card is going to be printed at the time of registration, the system will not know the name of the technologist. Therefore, the technologist will not print out on the form even if you include the technologist field when defining the format. See the chart on the preceding page for data requirements for the fields.

Fields that are to be printed twice in any one format should be enclosed in double quotes (e.g., a second Free Text field would be entered "FREE TEXT"). The Free Text field can be used to enter a statement which will serve as a heading or to indicate an area on the flash card where a comment or description is to be entered. When entered as free text, this statement will be interpreted by the system in its literal sense; i.e., exactly as you have entered it. You must then specify what will appear by entering the free text through the sub-field, "Literal Text". For example, if you make Free Text a part of the jacket label format and enter "***RETURN TO FILE ROOM***" at the "Literal Text:" prompt, then "***RETURN TO FILE ROOM***" will print on all jacket labels using the format.

Before starting, design your format on paper, adding spacing, counting characters for alignment, etc. Using graph paper makes it easier to count characters, lines, columns and spaces. You can use the <Space> key to adjust spacing. We'll show you how that works in our example.

Let's design the following flash card format. Note our plan to line up data by using spaces between the field and the colon.

NAME :		SSN:
EXAM DATE:		CASE:
PROCEDURE:		
COMMENTS :		
DECLIECHTNO	DIIVC.	
REQUESTING	PHIS:	

Prompt/User Response

Label/Header/Footer Formatter

>>> Label/Header/Footer Formatter <<<

Select flash card formats format name: Flash & ${\bf COMMENT}$

Are you adding 'FLASH & COMMENT' as a new FLASH CARD FORMATS (the 7TH)? \mathbf{Y} (Yes)

FORMAT NAME: FLASH & COMMENT// <RET>

Select FIELD: FREE TEXT

Are you adding 'FREE TEXT' as a new PRINT FIELDS (the 1ST for this FLASH CARD FORMATS)? Y (Yes)

ROW: 5

COLUMN: 2

Discussion

Note: The system refers to anything you design as a Flash Card Format even if it is a label or report header.

Enter "??" to obtain a list of formats already designed by your site. We'll call our format, FLASH & COMMENT.

This is the 7th one we have on file. They can be selected and edited at any time.

We'll make our first field a free text field.

We'll have it print on row 5 and begin printing in column 2.

LITERAL TEXT: COMMENTS<sp>:<sp>

We want it printed as "COMMENTS:". We added one space between the S and the colon to line up colons and data that will be printed to the right of the colons.

Select FIELD: REQUESTING PHYSICIAN

Are you adding 'REQUESTING PHYSICIAN' as a new PRINT FIELDS (the 2ND for this FLASH CARD FORMATS)? Y (Yes)

ROW: 10 COLUMN: 2

TITLE (OPTIONAL): REQUESTING PHYS:<sp>

Select FIELD: NAME OF PATIENT

Are you adding 'NAME OF PATIENT' as a new PRINT FIELDS (the 3RD for this FLASH CARD FORMATS)? ${\bf Y}$ (Yes)

ROW: 1 COLUMN: 2

TITLE (OPTIONAL): NAME<sp><sp><sp><sp><sp><sp>

Select FIELD: SSN OF PATIENT

Are you adding 'SSN OF PATIENT' as a new PRINT FIELDS (the 4TH for this FLASH CARD FORMATS)? ${\bf Y}$

(Yes)
ROW: 1
COLUMN: 40

TITLE (OPTIONAL): SSN<sp>:<sp>

Requesting Physician will be another field.

It will print in row 10, column 2 as "REQUESTING PHYS: ". We added a space after the colon before hitting the <RET> key.

We are using five spaces between NAME and the colon and one space after the colon to line up colons and data.

There's one space between SSN and the colon and one after the colon.

```
Select FIELD: DATE OF EXAM
 Are you adding 'DATE OF EXAM' as a new PRINT FIELDS
  (the 5TH for this FLASH CARD FORMATS)? Y (Yes)
 ROW: 2
 COLUMN: 2
 TITLE (OPTIONAL): EXAM DATE: < sp>
Select FIELD: PROCEDURE
 Are you adding 'PROCEDURE' as a new PRINT FIELDS
  (the 6TH for this FLASH CARD FORMATS)? Y (Yes)
 ROW: 3
 COLUMN: 2
 TITLE (OPTIONAL): PROCEDURE: < sp>
Select FIELD: CASE NUMBER
  Are you adding 'CASE NUMBER' as a new PRINT FIELDS
  (the 7TH for this FLASH CARD FORMATS)? Y (Yes)
 ROW: 2
 COLUMN: 40
 TITLE (OPTIONAL): CASE:<sp>
Select FIELD: <RET>
                                                       key.
NUMBER OF ROWS IN FORMAT: 10
```

When we finish adding fields, we hit the <RET> key.

Find the field with the highest row number and enter that number here. In this case it was Requesting Physician.

Then the format of the flash card is displayed.

Comments:

The new patient address label print fields are useful for printing mailing labels for reminder postcards, mammography follow-up, etc.

The order of field entry does not matter. For example, the items on the 7th row may be entered before rows 1-6 because the compiler reorders them per your numbering instructions.

Major AMIS Code Entry/Edit [RA MAJORAMIS]

This function only allows you to edit existing AMIS codes. AMIS (Automated Medical Information System) codes are used in the compilation of various workload reports, including the AMIS Report.

When the system is first installed, all necessary entries and data associated with each entry will be part of the module. This option should only be used if VA Headquarters changes the AMIS code description or weight.

Prompt/User Response	Discussion
Major AMIS Code Entry/Edit	Enter the code you want to edit
New entries and modifications to existing entries are prohibited without approval from Radiology Service VACO.	
Do you have approval from Radiology Service VACO to add/change any AMIS codes and weight? No// YES	
Select MAJOR RAD/NUC MED AMIS CODES DESCRIPTION: GASTROINTESTINAL	
DESCRIPTION: GASTROINTESTINAL// <ret></ret>	
WEIGHT: 6//	Remember; only make changes that have been made by VA Headquarters.

Select MAJOR RAD/NUC MED AMIS CODES DESCRIPTION: <RET>

-

¹ Patch RA*5*31 September 2002 Need VACO approval for new entries and modifications.

Supervisor Menu ... Utility Files Maintenance Menu ...

Nuclear Medicine Setup Menu

Lot (Radiopharmaceutical) Number Enter/Edit Route of Administration Enter/Edit Site of Administration Enter/Edit Vendor/Source (Radiopharmaceutical) Enter/Edit

Overview:

This menu contains the options to build and manage the Nuclear Medicine radiopharmaceutical lot numbers, routes and sites of administration, and vendor source.

A number of Routes of Administration and Sites of Administration were exported with the software. Print the contents of their respective files (Route of Administration #71.6 and Site of Administration #71.7) using the options Route of Administration List and Site of Administration List before you attempt to add data. You may edit the names of the routes and sites, but you may not delete them.

Add sites (Site of Administration Enter/Edit) before Routes of Administration since the sites must be in place so that you can select them during route edits and associate them with routes. Likewise, vendors/sources must be in place before editing lots.

Lot (Radiopharmaceutical) Number Enter/Edit

[RA NM EDIT LOT]

This option allows you to enter lot numbers for the radiopharmaceuticals used at your facility. Before entering new lots, you may need to first use the Vendor/Source (Radiopharmaceutical) Enter/Edit option to enter sources for use with this option.

Prompt/User Response

Discussion

Lot (Radiopharmaceutical) Number Enter/Edit

Select RADIOPHARMACEUTICAL LOT NUMBER/ID: 048

Are you adding '048' as a new RADIOPHARMACEUTICAL LOT (the 26TH)? ${\bf y}$ (Yes)

RADIOPHARMACEUTICAL LOT SOURCE: UNIVERSITY LABS

NUMBER/ID: 048// <RET>

SOURCE: UNIVERSITY LABS// <RET>

EXPIRATION DATE:

KIT #: <RET>

RADIOPHARM: SULFUR COLLOID TC-99M DX201

Select RADIOPHARMACEUTICAL LOT NUMBER/ID:

You can enter a new Lot number or ID here.

Select a source from the list that you built using the option Vendor/Source (Radiopharmaceutical) Enter/Edit.

When a lot is no longer being used, enter an expiration date.

Enter a kit number if it applies to this lot.

Entering one or two "?" here will give you a list of active radiopharmaceuticals that you can associate with this lot.

Route of Administration Enter/Edit

[RA NM EDIT ROUTE]

Use this option to build and manage a list of routes of administration and their corresponding sites of administration. You may want to use the Site of Administration Enter/Edit option first to define sites that will be associated with the routes you enter here.

Note: If no sites are associated with a particular route and status tracking parameters require/ask a route/site, the status parameter is ignored.

```
Route of Administration Enter/Edit
Select ROUTE OF ADMINISTRATION NAME: INTRAVENOUS
Are you adding 'INTRAVENOUS' as a new ROUTE OF ADMINISTRATION (the 19TH)? Y
NAME: INTRAVENOUS// <RET>
Select SYNONYMS: <RET>
   Current parameters for entry of sites for this route :
     PROMPT FOR FREE TEXT SITE?
    VALID SITES OF ADMINISTRATION =
                     -- NOTE --
         If 'PROMPT FOR FREE TEXT SITE?' is 'Y',
          then users will not be given a selection
          of predefined 'VALID SITES'
     Select one of the following:
          Р
                  PROMPT FOR FREE TEXT SITE?
                   VALID SITES OF ADMINISTRATION
Edit which field: VALID SITES OF ADMINISTRATION
Select VALID SITES OF ADMINISTRATION: ??
    Valid sites of drug administration for this route of administration
     should be entered here. During case edits, users will be restricted to
     selecting only sites entered here when this route of administration is
    used.
Choose from:
   CERVICAL SUBARACHNOID SPACE
  ENDOTRACHEAL TUBE
  GALL BLADDER
  LEFT ANTECUBITAL FOSSA
Select VALID SITES OF ADMINISTRATION: LEFT JUGULAR
Select VALID SITES OF ADMINISTRATION: RIGHT JUGULAR
Select VALID SITES OF ADMINISTRATION: (etc.) <RET>
```

Site of Administration Enter/Edit

[RA NM EDIT SITE]

Use this option to build a list of sites for the administration of radiopharmaceuticals. You can also enter a synonym for the site, so users have a choice of terms for selection. This list is used by the option Route of Administration Enter/Edit in the Valid Sites of Administration field.

Prompt/User Response

Discussion

Site of Administration Enter/Edit

Select SITE OF ADMINISTRATION NAME: LEFT JUGULAR

Are you adding 'LEFT JUGULAR' as a new SITE OF

ADMINISTRATION (the 36TH)? Y (Yes)

NAME: LEFT JUGULAR//

Select SYNONYM: <RET>

Select SITE OF ADMINISTRATION NAME:

Vendor/Source (Radiopharmaceutical) Enter/Edit

[RA NM EDIT SOURCE]

Use this option to build a list of vendors, pharmacies or other sources for the radiopharmaceuticals used in your facility. This list is used by the Lot (Radiopharmaceutical) Number Enter/Edit option in the Radiopharmaceutical Lot Source field.

Prompt/User Response

Discussion

```
Vendor/Source (Radiopharmaceutical) Enter/Edit
Select RADIOPHARMACEUTICAL SOURCE VENDOR/PHARMACY:
CHEMICALS INC.
  Are you adding 'CHEMICALS INC.' as
     a new RADIOPHARMACEUTICAL SOURCE (the 11TH)? Y
(Yes)
VENDOR/PHARMACY: CHEMICALS INC.// <RET>
Select RADIOPHARMACEUTICAL SOURCE VENDOR/PHARMACY:
```

Supervisor Menu ... Utility Files Maintenance Menu ... Order Entry Procedure Display Menu ...

> Common Procedure Enter/Edit Create OE/RR Protocol from Common Procedure Display Common Procedure List

Overview:

This menu contains the options to manage the common procedure display used in either OE/RR or in the Request an Exam option.

Common Procedure Enter/Edit [RA COMMON PROCEDURE]

The most frequently ordered procedures can be set up in a Common Procedure Display to make ordering a Rad/Nuc Med procedure easier for the ordering clinician. You may enter as many procedures into the Common Procedure file as you wish, however, you may only use sequence numbers from 1 to 40 per Imaging Type. Only procedures with sequence numbers are displayed, therefore the maximum number of common procedures that can be displayed is 40 per Imaging Type.

This menu option may be used to add and/or remove (inactivate) procedures from the display (Common Procedure List) users see when entering requests. Procedures should not be deleted from this file, just inactivated. Active common procedure entries require a sequence number to tell the system the order in which they should be displayed on the screen. When you enter "YES" to the Inactivate from List prompt field it will automatically remove the procedure's sequence number. The system will automatically resequence active common procedures upon exit from the option.

If you want to assign a new entry a particular sequence number on the Common Procedure List, you should reassign the sequence numbers of all procedures on the list that follow the insertion point and assign the selected sequence number to the new entry. For example, if you have 10 entries on your Common Procedure List and you wish to enter a procedure as sequence number 8, you could reassign sequence numbers 8, 9 and 10 to numbers 11, 12, and 13 and then assign the new entry to number 8. This manual realignment of sequence numbers should all be done before exiting the option, because this tells the system the order in which the procedures should be displayed.

The system will automatically resequence the numbers to remove gaps in numbering upon exit. For example, say you exit the option leaving 10 active common procedures for Vascular Lab Imaging Type and they happen to have sequence number 3, 5, 6, 7, 10, 20, 35, 37, 38, 40. The automatic re-sequencing will replace their order numbers with 1, 2, 3, 4, 5, 6, 7, 8, 9, 10 respectively.

The Radiology/Nuclear Medicine Request an Exam option will display these procedures in the order in which you sequence them. Version 2.5 of OE/RR will also display an identically numbered common procedure list. Version 3.0 of OE/RR will display the same common procedures; however, they will always appear in alphabetical order. To accommodate this, you may wish to number the procedures so that their names are alphabetized if the same clinicians will be using both OE/RR version 3.0 and the Request and Exam option.

If you are running version 2.5 of OE/RR, you may also establish a procedure as an orderable item in the OE/RR package through this option. You will be prompted for the name the procedure can be ordered by through OE/RR, Imaging Location where the procedure will be performed, category of exam, requested date of exam, mode of transport of the patient, urgency of the exam, whether isolation procedures are required and any modifiers associated with the procedure. You **must** enter a Name and Imaging Location or the procedure will not be accepted as an orderable item. After establishing a procedure as an orderable item in OE/RR, you should utilize the Create OE/RR Protocol from Common Procedure option to actually set up the orderable item (protocol) in the OE/RR package. You may not edit a protocol once it has been created through the Create OE/RR Protocol from Common Procedure option. If you are running a version of OE/RR higher than 2.5, none of the above applies, and you will not see any prompts for creating protocols. Creating quick order protocols for OE/RR version 3.0 and beyond will be done through OE/RR and will no longer be possible through Radiology/Nuclear Medicine.

There is no limit to the number of entries that may be added to the Rad/Nuc Med Common Procedure file #71.3. Any procedure entered into this file (with or without a sequence number) can be converted into an OE/RR protocol, thus making the procedure an orderable item selectable by any user requesting an exam through OE/RR. When you establish an orderable item through this option and then create a protocol through the Create OE/RR Protocol from Common Procedure option, you are building a quick order. This means you may select the protocol through OE/RR and have certain information already associated with the procedure.

It is suggested that before you add a new entry to the Common Procedure List, you use the option Display Common Procedure List to obtain a list of procedures and their sequence numbers. That will help you determine where you want the new entry placed and what sequence numbers can be used.

Prompt/User Response

Discussion

Common Procedure Enter/Edit

INACTIVATE FROM LIST: <RET>

```
Select one of the following Imaging Types:
   GENERAL RADIOLOGY
  NUCLEAR MEDICINE
  ULTRASOUND
  CT SCAN
  CARDIOLOGY STUDIES (NUC MED)
Select IMAGING TYPE: NM NUCLEAR MEDICINE
Select RAD/NUC MED COMMON PROCEDURE: ??
  PARATHYROID SCAN (3)
  GASTROESOPHAGEAL REFLUX STUDY
                                  (2)
  GI BLOOD LOSS-ACUTE GI BLEED
                                 (5)
  WHOLE BLOOD VOLUME (4)
  LYMPHATIC/LYMPH NODE IMAGING (no sequence number)
  HEPATOBILIARY SCAN (HIDA SCAN) (6)
  ESOPHAGEAL MOTILITY (7)
  CARDIAC SHUNT DETECTION
                            (8)
  CEREBRAL BLOOD FLOW (1)
   You may enter a new RAD/NUC MED COMMON PROCEDURE
   if you wish. Enter a procedure name. Only
   active procedures may be selected.
Answer with RAD/NUC MED PROCEDURES NAME
Do you want the entire RAD/NUC MED PROCEDURES
             List? NO
Select RAD/NUC MED COMMON PROCEDURE: BONE SCAN
                                   (NM Broad
 Are you adding 'BONE SCAN' as
   a new RAD/NUC MED COMMON PROCEDURE (the 19TH)? Y
                                               (Yes)
PROCEDURE: BONE SCAN// <RET>
```

Note that the procedures already entered as Common Procedures appear with their sequence number to the right. If there is no sequence number, that procedure has been "inactivated" from the common procedure list.

When adding a new procedure, simply bypass this prompt. To remove a procedure from the Common Procedure list, enter YES at this prompt. To reactivate a procedure, enter an @ at this prompt to remove the YES response.

SEQUENCE NUMBER: 12// <RET>

Do you wish to establish this procedure as a quick order item in OE/RR? No// ${\bf YES}$

Let's make this an orderable protocol in OE/RR. all prompts below this are for the purpose of defining enough data to establish a quick order protocol. If you answer NO to this prompt, none of the prompts below will appear. (This prompt and the ones below will not appear if your site has installed a version of OE/RR higher than 2.5.)

NAME OF ORDERABLE ITEM: BONE SCAN

This is the text that will appear to the user ordering the exam.

Select MODIFIERS: <RET>

There are no specific modifiers that go with this exam so we can bypass this prompt.

IMAGING LOCATION: ?

Enter the Imaging Location for this procedure.
The Imaging Type of the procedure and the Imaging
Type of the location must be the same.
Answer with IMAGING LOCATIONS TYPE OF IMAGING
Choose from:

NUC MED LOC (NUCLEAR MEDICINE-449)
NUC (NUCLEAR MEDICINE-639)

IMAGING LOCATION: NUC

The exam request will be printed at the location entered here.

```
CATEGORY OF EXAM: ?
                                                         In this example, the default
     If an entry is not made for this field, the
                                                         of Inpatient is what we
     category of 'Inpatient' will
                                                         want, so we strike the
     automatically be entered for each request
                                                         <RET> key to bypass this
     generated from the protocol that
     is built for this procedure.
                                                         prompt.
     Choose from:
       Т
               INPATIENT
       Ω
               OUTPATIENT
       С
               CONTRACT
                SHARING
       F.
               EMPLOYEE
               RESEARCH
CATEGORY OF EXAM: <RET>
REQUESTED DATE: ?
                                                         The default of Today is fine
     If an entry is not made in this field, the date
                                                         so we strike the <RET> key
     of today (TODAY) will automatically be
                                                         to bypass this prompt.
     entered for each request generated from the
     protocol that is built for this procedure.
     Choose from:
       TODAY
              TODAY
       TODAY+1 TOMORROW
REQUESTED DATE: <RET>
MODE OF TRANSPORT: ?
                                                         Ambulatory is acceptable,
     If an entry is not made in this field, the mode
                                                         so we can bypass this
     of transport of 'ambulatory' will
                                                         prompt.
     automatically be entered for each request
     generated from the protocol that
     is built for this procedure.
     Choose from:
       а
               AMBULATORY
                PORTABLE
               STRETCHER
               WHEEL CHAIR
MODE OF TRANSPORT: <RET>
ISOLATION PROCEDURES: ?
                                                         Isolation is not required.
     If an entry is not made in this field, a value
     of 'NO' (isolation procedures not required)
     will automatically be entered for each request
     generated from the protocol that is built for
     this procedure.
     Choose from:
                YES
                NO
       n
ISOLATION PROCEDURES: <RET>
```

REQUEST URGENCY: ?

Routine is acceptable.

If an entry is not made in this field, a request urgency of 'ROUTINE' will automatically be entered for each request generated from the protocol that is built for this procedure. Choose from:

1 STAT
2 URGENT
9 ROUTINE
REQUEST URGENCY: <RET>

Select RAD/NUC MED COMMON PROCEDURE: <RET>

Create OE/RR Protocol from Common Procedure [RA CREATE OE/RR PROTOCOL]

Note: This option is functional only if OE/RR Version 2.5 is installed. If a higher version of OE/RR is installed, this option is not needed and will display a message to that effect.

This option is used to create orderable items (protocols) in the Order Entry/Results Reporting (OE/RR) package Version 2.5 based upon entries in the Rad/Nuc Med Common Procedures file. When you choose the option, the entire list of procedures in the Rad/Nuc Med Common Procedure file is presented for your selection as shown below. Press the <RET> key until your choice comes up. Then select the number of that choice. We added Bone Scan to the file in the previous chapter. Now a protocol will be created for it.

Prompt/User Response

Discussion

```
Create OE/RR Protocol from Common Procedure
         Create an OE/RR Protocol from a Radiology/Nuclear Medicine
         Common Procedure
     ACUTE GI BLOOD LOSS IMAGING
     ANGIO CAROTID CEREBRAL UNILAT S&I
     ANGIO CERVICOCEREBRAL CATH CP
4
    ARTHROGRAM KNEE CP
5
    BONE AGE
6
    BRAIN IMAGING COMPLETE STUDY WITH VASCULAR FLOW, PLANAR
7
    CHEST STEREO PA
    ECHOCARDIOGRAM M-MODE COMPLETE
8
    MAMMOGRAM BILAT
9
10 RADIONUCLIDE THERAPY, HYPERTHYROIDISM
11
   ABDOMEN 1 VIEW
12
   ARTHROGRAM ANKLE S&I
13 ANGIO CORONARY BILAT INJ S&I
14
    ECHOGRAM ABDOMEN COMPLETE
15
    RADIONUCLIDE THERAPY, THYROID SUPPRESSION
Type '^' to STOP, or
Choose from 1-15: <RET>
16 BRAIN IMAGING, PLANAR ONLY
17
    BONE SCAN
     ACROMIOCLAVICULAR J BILAT
Type '^' to STOP, or
Choose from 1-18: 17
Processing the BONE SCAN procedure.
 Located in the RA (RADIOLOGY/NUCLEAR MEDICINE) namespace.
```

RA584 BONE SCAN PROTOCOL Created

Display Common Procedure List [RA DISPLAY COMMON PROCEDURES]

This option allows you to look at what the common procedure list will look like to the user, when using the request a procedure either through OE/RR version 2.5 or the Radiology/Nuclear Medicine option. OE/RR version 3.0 will also display the same procedures, but they will be in alphabetical order.

Prompt/User Response

Discussion

Display Common Procedure List Select one of the following Imaging Types: GENERAL RADIOLOGY NUCLEAR MEDICINE ULTRASOUND CT SCAN CARDIOLOGY STUDIES (NUC MED)

Select IMAGING TYPE: GENERAL RADIOLOGY

COMMON RADIOLOGY/NUCLEAR MEDICINE PROCEDURES (GENERAL RADIOLOGY) ______

- 1) ACUTE GI BLOOD LOSS IMAGING
 2) ANGIO CAROTID CEREBRAL UNILAT S&
 3) ANGIO CORONARY BILAT INJ S&I
 4) ARTHROGRAM ANKLE S&I
 5) BRAIN IMAGING, PLANAR ONLY
 6) CHEST STEREO PA
 7) ECHOGRAM ABDOMEN COMPLETE
 8) RADIONUCLIDE THERAPY, THYROID SU
 9) ABDOMEN 1 VIEW
 10) MAMMOGRAM BILAT
 11) BONE AGE
 12) BONE SCAN

Press RETURN to continue... <RET>

Utility Files Maintenance Menu

Supervisor Menu ... Utility Files Maintenance Menu ... Procedure Edit Menu ...

> Cost of Procedure Enter/Edit Procedure Enter/Edit Procedure Message Entry/Edit Procedure Modifier Entry

Overview:

This menu contains the options to define and manage data for procedures.

Cost of Procedure Enter/Edit

[RA PROCOSTEDIT]

This function allows you to apply a cost to procedures. Unit cost and Total cost appear on the Procedure/CPT Statistics Report [RA CPTSTATS].

You are asked for an imaging type, the procedure you want to edit, and the cost. Only procedures matching the selected imaging type(s) can have their cost edited. Select procedures by entering either the CPT code for the procedure or the procedure name.

Cost of Procedure Enter/Edit

Select Imaging Type: All// <ret>

Another one (Select/De-Select): <ret>

Select RAD/NUC MED PROCEDURES NAME: 71020 CHEST X-RAY CHEST 2 VIEWS

PA&LAT (RAD Detailed) CPT:71020

COST OF PROCEDURE: 99.98

Procedure Enter/Edit [RA PROCEDURE]

This function allows you to enter new procedures into the system and to edit existing procedures. Entries in this file are used as the allowable selections for users at any Procedure prompt. This function must be used to associate the CPT (Current Procedural Terminology) codes for applicable descendent procedures with contrast media. This prompt appears as part of the Request an Exam, Register an Exam, Add Exam to Last Visit, Case No. Exam Edit, and the Exam Status Tracking options.

1. RAD/NUC MED PROCEDURES NAME - You can only edit the name of the procedure if the procedure was added by your site; if the procedure was distributed by the package, editing the name field is not allowed. Entries are related to Imaging Types. If the same procedure is done by two different locations of differing Imaging Types then you need to enter them with slightly different names but with the same CPT and AMIS codes:

US - Carotid Doppler VAS - Carotid Doppler

¹**Note**: In order to comply with the CPRS procedure naming convention, the name must not contain a semicolon (;).

2. TYPE OF IMAGING - This field is used to associate a particular Imaging Type to a procedure. Each procedure can only be associated with one Imaging Type. You may enter a new procedure with a slightly different name to make it show up for a second Imaging Type.

Example: You want Carotid Doppler for Ultrasound and for Vascular Lab. One procedure can be called Carotid Doppler and assigned to the Ultrasound Imaging Type and the other VAS - Carotid Doppler with the Vascular Lab Imaging Type.

If you have activated a new Imaging Type and are changing the Imaging Type of procedures to match the new Imaging Type, you will also have to deactivate these procedures from the Common Procedure list (if they are common procedures). They can later be put on a separate new common procedure list under the new Imaging Type. Refer to Common Procedure Enter/Edit for more information.

_

¹ Patch RA*5*34 Inserted text regarding semicolon restriction on procedure name.

3. TYPE OF PROCEDURE - There are four types of procedures:

Detailed - Procedure is associated with an AMIS code and must have a CPT code assigned to it.

Series - Procedure is associated with more than one AMIS code and must have a CPT code assigned to it. It is recommended that the Type of Procedure field be set to Series for all procedures with more than one AMIS code. In addition to being reported under each AMIS code, the procedures with multiple AMIS codes are also tallied in the "Series of AMIS Codes" line item of the AMIS Report. Setting these procedures as a Series type serves as a reminder to users. The AMIS Report counts each Series procedure the same as the Detailed procedure in the body of the report, but an additional line tallies the Series procedures by themselves. The Bilateral, Operating Room, and Portable modifiers are not allowed for Series procedures as that would cause the AMIS Report to be incorrect.

Broad - Procedure is mainly used by the receptionist, clerk, or ordering clinician when scheduling a patient for an exam. It is used when s/he is not exactly sure which Detailed or Series procedure will be performed. However, before the exam can be considered Complete, the procedure should be changed to a Detailed or Series in order to ensure adequate workload credit to your facility. Broad procedures are not associated with a CPT or AMIS codes. These prompts will not be seen if the user selects Broad as the procedure type. If no CPT code is entered for a Detailed or Series procedure, the procedure type is automatically changed to Broad.

If the Division parameter requiring a Detailed or Series procedure upon initial exam registration is set to YES, then the ordering party will not be able to select a Broad procedure.

Parent - Procedure is used for ordering purposes only. It must have descendents (other procedures of Detailed or Series type) that are actually registered. Parent procedures are used to simplify the ordering process when a group of related procedures must be done. (See Parent/Descendent Exam and Printsets) Similar to Broad procedures, they act as temporary placeholders. When a parent order is requested, the system automatically finds all the descendants and checks for the use of contrast media. The parent procedure is never actually registered. If one or more of the registered descendent procedures associated with the parent procedure uses contrast media, the parent procedure will be flagged as being associated with contrast media.

When creating a Parent procedure, make the parent name slightly different than any of its descendents so the system will accept it as a new procedure. **Remember, when creating a Parent, no CPT code is entered therefore, it will not be registered**. All descendents must be either Detailed or Series Type with a CPT code.

Example:

Parent: Arthrogram of the Shoulder (or Arthrogram -

Shoulder)

Descendents: 73040

23350

4. MODALITY - This is a new field that only appears when you are defining either Detailed or Series type procedures. Each procedure can have more than one modality. The modality data is used by the Imaging package to produce a modality work list that will help to correctly identify a patient. If you are not running the Imaging package at your site, you may bypass this field.

- 5. SINGLE REPORT This is a new field that only appears when you define a procedure as a Parent procedure. If this field is left unanswered, then the reporting of the procedure's descendents will be put into separate reports, as in the previous version. If this field is answered YES, then the reporting of the procedure's descendents will be combined into one report, where the case records for these descendents will have the same data for report text, impression, pre-verifier, verifier, staff (primary and secondary), residents (primary and secondary), diagnoses, and all other report-related fields.
- 6. HEALTH SUMMARY WITH REQUEST This field points to the Health Summary Type file. When a procedure is requested, this field is checked for an associated Health Summary Type and, if one exists, it is printed along with the request. You can select an appropriate Health Summary Type here. If necessary, you may need to check with your clinical coordinator for help in setting up useful health summary formats.
- 7. SYNONYM One or more synonyms may be associated with a procedure. At any Procedure prompt, a synonym can be entered to select a procedure if you have predefined the synonym here.

¹ Patch RA*5*3

- 8. PROMPT FOR MEDS This is a new field. It only appears for procedures that are Detailed or Series. If this field is set to YES, case edits will prompt for medications administered for this procedure. To make Status Tracking prompt for medications administered, this field must be set to YES and the Examination Status parameter "Ask Medications & Dosages" must be set to YES for one or more statuses of the procedure's Imaging Type.
- 9. DEFAULT MEDICATION This is a new field. It only appears for procedures that are Detailed or Series. Use this field to enter any medications when they are usually used for the procedure. Radiopharmaceuticals may not be entered in this field. Medications, if defined here, are entered automatically into the patient exam record by the system during registration.
 - DEFAULT MED DOSE will only appear if you enter a Default Medication. If a default dose is entered, it will be automatically entered on the patient's exam record at the time the procedure is registered. Acceptable format: 99999.999 unit_of_measure. Unit_of_measure may be between 1-30 characters in length. Example: 2.5 mg
- 10. SUPPRESS RADIOPHARM PROMPT This is a new prompt that only applies to Detailed or Series procedures with an Imaging Type that uses radiopharmaceuticals (i.e., nuclear medicine or cardiology studies). If the exam status setup for the affected Imaging Types specifies that radiopharmaceutical data should be asked and/or required to progress to various statuses, this field can be used to over-ride that requirement for a particular procedure.

Available categories are defined as follows:

- 0 or blank Ask for Radiopharmaceutical if and when exam status parameters specify (this is the default).
- Suppress the Radiopharmaceutical prompt. Never ask or require Radiopharmaceutical even if exam status parameters specify that a Radiopharmaceutical should be entered. Even if default Radiopharmaceuticals are entered for the procedure they will not be automatically entered on the exam record by the system when the case is registered.

11. PROMPT FOR RADIOPHARM RX - This is a new prompt that only applies to

Detailed or Series procedures with an Imaging Type that uses radiopharmaceuticals (nuclear medicine or cardiology studies). It controls whether case edit and status tracking options prompt for prescribing physician, prescribed radiopharmaceutical dose, and witness to the dose administration. If you answer YES to this field, prescriber and witness prompts will appear for any radiopharmaceuticals entered for this procedure. (For example, a witness is necessary for a Total Body Iodine procedure so you would answer YES at this field for that procedure.) If this field is left blank, prescriber and witness prompts will not appear. If this field is set to YES, Status Tracking will prompt during each edit session until the data is entered; case edits will always prompt, even if the data has already been entered.

12. DEFAULT RADIOPHARMACEUTICAL - This is a new prompt that only applies to

Detailed or Series procedures with an Imaging Type that uses radiopharmaceuticals (nuclear medicine or cardiology studies). It will be automatically entered as a radiopharmaceutical on patients' nuclear medicine exam records during registration. This information can be changed, if necessary, during case edit. If nothing is entered in this field, no radiopharmaceuticals will be automatically entered during registration for this procedure.

Note: The radiopharmaceutical is the only data item that is actually automatically entered on the patient exam record during registration. No dosage amounts or routes of administration are entered automatically; however, they will appear as default responses.

LOW ADULT DOSE (in mCi) - This is the lowest usual adult dosage commonly associated with this radiopharmaceutical when administered for this procedure. The system will display this value and the high adult dose when asking for the actual dosage during case editing. If the dosage entered for an exam falls outside of the usual high/low range, the person editing the case will see a warning message. If you do not want a warning to appear no matter how low the dosage entered, do not enter a low range dosage here.

- HIGH ADULT DOSE (in mCi) This is the highest usual adult dosage commonly associated with this radiopharmaceutical when administered for this procedure. The system will display this value and the low adult dose when asking for the actual dosage during case editing. If the dosage entered for an exam falls outside of the usual high/low range, the person editing the case will see a warning message. This is a mandatory data item.
- USUAL DOSE (in mCi) This is the usual dosage (in mCi) associated with this radiopharmaceutical. It will appear as an informational item on the screen just before the user is prompted for the actual Dose. If nothing is entered in this field, no usual dose information will be displayed. The usual dose you enter here must fall within the Low Adult Dose and High Adult Dose described above.
- DFLT ROUTE OF ADMINISTRATION This is the route of administration commonly associated with this radiopharmaceutical for this procedure. It will appear as the default response when the route of administration is asked during case editing. If nothing is entered in this field, no default route of administration will be displayed.
- DEFAULT SITE OF ADMINISTRATION This is the bodily site for administering the radiopharmaceutical for the procedure being defined. The system will use this for the default response during case editing. If nothing is entered here, there will be no default response to the Site of Administration question during case editing.
- DEFAULT FORM This is the form of this Radiopharmaceutical when used in this procedure. If entered, it will be used as the default answer to the "Form" prompt during case edits.

June 2019

13. AMIS CODE -This field contains the valid AMIS code(s) from the Major AMIS Codes file associated with this exam procedure. There may be more than one AMIS code associated with each procedure. If the procedure type is Detailed, there will be one AMIS code. If the procedure type is Series, there will be more than one AMIS code. Broad and Parent types of procedures have no AMIS code.

AMIS WEIGHT MULTIPLIER - You are not likely to change an AMIS Weight Multiplier unless VA HEADQUARTERS rules the change. However, if you add a new procedure, you must determine from VA HEADQUARTERS guidelines what should be entered at this prompt.

This field contains a number (0-99) to indicate to the various workload report routines how many times to multiply the weighted work units associated with the AMIS code. The weight for each AMIS code is stored in the Major Rad/Nuc Med AMIS Code file.

Most multipliers will be 1. However, there are some that are greater than 1. For example, the procedure called Upper GI + Small Bowel has the AMIS code 9-Gastrointestinal which has a weight of 6 work units. But, the AMIS multiplier is 2. Therefore on the workload reports, the site will get credit for 12 weighted work units each time it is performed.

Radiology Central Office has directed that the weighted work unit for Ultrasound procedures involving multiple organs in the same visit can be counted only once if the organs are located close to each other (i.e. Liver and Gallbladder).

BILATERAL - This field contains a YES to indicate that this is a bilateral procedure or a NO if it is not a bilateral procedure. If it contains a YES, then the system multiplies the AMIS weight by 2. No harm is done by indicating the Multiplier as 2 and Bilateral as YES. The system will only use 2 as the multiplier, not (2*2).

- 14. CPT CODE Editing the CPT Code for existing procedures is no longer allowed because of possible contrast media associations. If the CPT Code of the procedure in the RAD/NUC MED PROCEDURE (#71) file must be changed, the procedure should be INACTIVATED and a new procedure, with the correct CPT Code, should be created.
- 15. ¹DEFAULT CPT MODIFIERS(PROC) Default CPT Modifier(s) entered for this procedure will be automatically entered into the exam record during exam registration.
- 16. STAFF REVIEW REQUIRED If the procedure requires an interpreting staff physician's review of the resident's report, even though the division parameter is set to allow residents to verify reports, enter a YES at this prompt. (Refer to Personnel Classification and Division Parameter Set-up for other factors that control whether residents can verify reports.) This field does not appear with Broad or Parent procedures.
- 17. REQUIRED FLASH CARD PRINTER This field contains the name (1-20 characters) of the location where the flash cards for this procedure are to be printed. Normally, the flash cards are printed at the reception desk, however, for particular procedures the flash cards may be printed at a different location such as where the exam is to be performed. This field does not appear with Parent procedures since Parent procedures are not registered.
- 18. **Note:** The exam labels will print out on the flash card printer associated with the imaging location parameter.
- 19. REQUIRED FLASH CARD FORMAT If a flash card format has been defined, you can associate a specific format with this procedure. Otherwise the format defined for the location will be used. This field does not appear with Parent procedures, since flash cards are only printed during registration and parent procedures are never registered.
- 20. FILM TYPE This field contains a default value of the type of film from the Film Sizes file that is usually needed for this procedure. During the exam registration process, the system checks this field for a film type. If film type exists, then it is automatically entered into the Film Type field of the exam record. The entry in the exam record can be edited. Only procedures that are performed a relatively high number of times, such as Chest PA&LAT, should have data entered for this field. This field does not appear with Parent procedures since parent procedures are never registered.

-

¹ Patch RA*5*19 May 2000 New field added to the Rad/Nuc Med Procedures file #71.

- 21. NORMAL AMOUNT NEEDED FOR EXAM This field contains the default value (a number from 1 through 10) of the normal amount of film usually required for this procedure. During the exam registration process, the system checks this field for data. If data is present, then the information in this field is stuffed into the "Number Used" field of the exam record. The technologist is allowed to modify this data after performing the exam. Only procedures that are performed a relatively high number of times should have data entered for this field. The procedure CHEST PA&LAT is a good example.
- 22. RAD/NM PHYS APPROVAL REQUIRED If YES is entered at this prompt, when a request for this procedure is entered, the request will require the name of the approving physician. This name will print on the request form.
- 23. MESSAGE You can select a message from the Rad/Nuc Med Procedure Message file for this procedure that will automatically appear at the time of the request and will be printed on the request form.
- 24. DESCENDENTS This field only appears with Parent procedures. Enter any Detailed or Series procedures that make up the set. Remember that the Parent procedure is a placeholder for ordering and will not be registered, so all procedures to be registered must be entered as descendents.
- 25. **Note**: To add a second procedure that is the same as one already entered, enter either the CPT code or the procedure name in quotes as shown here:

```
Select DESCENDENTS: FOOT 2 VIEWS// ?
Answer with DESCENDENTS
Choose from:
ANKLE 2 VIEWS
FOOT 2 VIEWS
TOE(S) 2 OR MORE VIEWS

You may enter a new DESCENDENTS, if you wish
Enter a procedure name which is related to the parent procedure.
Select only active detailed or series procedures of the same imaging type as the parent.
Answer with RAD/NUC MED PROCEDURES NAME, or CPT CODE, or SYNONYM
Do you want the entire RAD/NUC MED PROCEDURES List? n (No)
Select DESCENDENTS: FOOT 2 VIEWS//"FOOT 2 VIEWS"
```

26. EDUCATIONAL DESCRIPTION - This is a new field used to describe the

procedure and give information about its purpose in confirming or eliminating various diagnoses. It may be used for educational purposes by clinicians who are trying to familiarize themselves with imaging procedures that are available. If you enter an educational description, you are also asked if you want to "Display Educational Description When Ordered". Entering YES at that field displays the description along with any Procedure Message during the ordering process.

Note: If you are using the Screen Editor to add an Educational Description and you want a list of items, to keep the list from wrapping add a space (<sp>) before each item as shown in the following example:

```
When ordering this procedure, we recommend the following steps be taken:
<sp>1) Do this first.
<sp>2) Next do this.
<sp>3) Etc.
```

27. DISPLAY EDUCATIONAL DESCRIPTION WHEN ORDERED - This field

determines whether or not the Educational Description of a procedure should be displayed along with the procedure message (if any) at the time the procedure is ordered. Enter YES to display the description with the procedure message, NO to not display the description.

28. INACTIVATION DATE - If the procedure has an inactivation date prior to the current date, then the procedure is not a valid choice and will not appear on the list of active procedures when displayed to the user. Entering an inactivation date makes the procedure invisible to users so it cannot be selected. This should only be done if the procedure becomes obsolete or should not be ordered or performed at your site. You can enter a future date. This is useful if for example a CPT becomes inactive on January 1 and you want to enter it earlier than that date. Inactivation takes effect at **midnight** of the date you enter.

Here is an example of a Nuclear Medicine, Detailed procedure. The prompts for radiopharmaceuticals only appear when the procedure has an Imaging Type of Nuclear Medicine or Cardiology Studies:

Prompt/User Response

```
Select RAD/NUC MED PROCEDURES NAME: ACUTE GI BLOOD LOSS IMAGING
                                                                        (NM
   Detailed) CPT:78278
   NAME: ACUTE GI BLOOD LOSS IMAGING Replace <RET>
   TYPE OF IMAGING: NUCLEAR MEDICINE// <RET>
   TYPE OF PROCEDURE: DETAILED// <RET>
  <sup>1</sup>Select MODALITY: <RET>
   HEALTH SUMMARY WITH REQUEST: INPATIENT// <RET>
   Select SYNONYM: INPT// <RET>
PROMPT FOR MEDS: YES//<RET>
🖙 Select DEFAULT MEDICATION: LIDOCAINE 0.5% W/EPI INJ MDV// <RET>
   DEFAULT MEDICATION: LIDOCAINE 0.5% W/EPI INJ MDV// <RET>
DEFAULT MED DOSE: <RET>
     Reminder --
      Selection of 'Suppress' for 'SUPPRESS RADIOPHARM PROMPT' will result in
      the deletion of all default radiopharmaceutical information for this
     procedure.
SUPPRESS RADIOPHARM PROMPT: Always// <RET>
PROMPT FOR RADIOPHARM RX: YES// <RET>
Select DEFAULT RADIOPHARMACEUTICAL: SULFUR COLLOID TC-99M// PERCHLORACAP 250MG
   CAPS
                   DX201
            ...OK? Yes// <RET> (Yes)
     DEFAULT RADIOPHARMACEUTICAL: PERCHLORACAP 250MG CAPS
           // <RET>
     LOW ADULT DOSE (in mCi): 2// <RET>
     HIGH ADULT DOSE (in mCi): 5// <RET>
        Input a 'Usual Dose' value within the range of: 2-5
     USUAL DOSE (in mCi): 3// <RET>
     DFLT ROUTE OF ADMINISTRATION: ORAL// <RET>
     DEFAULT SITE OF ADMINISTRATION: MOUTH// <RET>
     DEFAULT FORM: Solid (Other)// <RET>
```

¹ Patch RA*5*3

```
Select DEFAULT RADIOPHARMACEUTICAL: < RET>
   Select AMIS CODE: NUCLEAR MEDICINE// <RET>
    AMIS CODE: NUCLEAR MEDICINE// <RET>
    AMIS WEIGHT MULTIPLIER: 1// <RET>
   Select AMIS CODE: <RET>
   CPT CODE: 78278// (no editing) <>>> NO EDITING ALLOWED
   STAFF REVIEW REQUIRED: NO// <RET>
   RAD/NM PHYS APPROVAL REQUIRED: NO// <RET>
   REQUIRED FLASH CARD PRINTER: <RET>
   REQUIRED FLASH CARD FORMAT: ER FORMAT// <RET>
   Select FILM TYPE: 11X14// <RET>
    FILM TYPE: 11X14// <RET>
    NORMAL AMOUNT NEEDED FOR EXAM: 1// <RET>
   Select FILM TYPE: <RET>
   Select MESSAGE: CALL DR. SMITH BEFORE ORDERING THIS PROCEDURE// <RET>
EDUCATIONAL DESCRIPTION:
   This test is designed to LOCATE the site of KNOWN GI bleeding, NOT to
   determine whether there IS GI bleeding.
     Edit? NO// <RET>
📂 DISPLAY EDUCATIONAL DESCRIPTION WHEN ORDERED: YES// <RET>
   INACTIVATION DATE: <RET>
   Here is an example of a Parent procedure entry:
   Prompt/User Response
   Select RAD/NUC MED PROCEDURES NAME: AORTAL FEMORAL RUNOFF
    Are you adding 'AORTAL FEMORAL RUNOFF' as
       a new RAD/NUC MED PROCEDURES (the 737TH)? Y (Yes)
     RAD/NUC MED PROCEDURES TYPE OF PROCEDURE: P PARENT
   NAME: AORTAL FEMORAL RUNOFF Replace <RET>
   TYPE OF IMAGING: GENERAL RADIOLOGY
   TYPE OF PROCEDURE: PARENT// <RET>
SINGLE REPORT: YES
   HEALTH SUMMARY WITH REQUEST: <RET>
   Select SYNONYM: <RET>
   RAD/NUC MED PHYS APPROVAL REQUIRED: NO// <RET>
   Select MESSAGE: <RET>
   Select DESCENDENTS: 75625 CONTRAST X-RAY EXAM OF AORTA AORTO ABD TRANS L
  W/SERIAL FILMS S&I (Detailed) CPT:75625
    Are you adding 'AORTO ABD TRANS L W/SERIAL FILMS S&I' as
       a new DESCENDENTS (the 1ST for this RAD/NUC MED PROCEDURES)? Y (Yes)
```

Select DESCENDENTS: **36200** PLACE CATHETER IN AORTA, CATHETER (FEMORAL OR AXILLARY APPROACH) (Detailed) CPT:36200

Are you adding 'AORTA, CATHETER (FEMORAL OR AXILLARY APPROACH)' as

a new DESCENDENTS (the 2ND for this RAD/NUC MED PROCEDURES)? \mathbf{Y} (Yes) Select DESCENDENTS: **75716** ARTERY X-RAYS, ARMS/LEGS ANGIO EXTREMITY BILAT S&I (Detailed) CPT:75716

Are you adding 'ANGIO EXTREMITY BILAT S&I' as

a new DESCENDENTS (the 3RD for this RAD/NUC MED PROCEDURES)? ${\bf Y}$ (Yes) Select DESCENDENTS: <a href="https://doi.org/10.1007/j.nep-10.1007

EDUCATIONAL DESCRIPTION:

No existing text Edit? NO// <RET>

INACTIVATION DATE: <RET>

Here is an example of editing an existing CT SCAN, Detailed, procedure. In this edit session, we are adding contrast media to the procedure.

Select RAD/NUC MED PROCEDURES NAME: CT ABDOMEN W/CONT (CT Detailed) CPT:74160

TYPE OF IMAGING: CT SCAN// <RET>
TYPE OF PROCEDURE: DETAILED// <RET>

CONTRAST MEDIA USED: No// Yes

Select CONTRAST MEDIA: ?

You may enter a new CONTRAST MEDIA, if you wish Enter the contrast agent(s) associated with this procedure.

Choose from:

- I Ionic Iodinated
- N Non-ionic Iodinated
- L Gadolinium
- C Cholecystographic
- G Gastrografin
- B Barium
- M unspecified contrast media

Select CONTRAST MEDIA: N (N Non-ionic Iodinated)

Are you adding 'Non-ionic Iodinated' as

a new CONTRAST MEDIA (the 1ST for this RAD/NUC MED PROCEDURES)? No// Y

(Yes)

Select CONTRAST MEDIA: <RET>
Select MODALITY: CT// <RET>

HEALTH SUMMARY WITH REQUEST: <RET>

Select SYNONYM: **<RET>**PROMPT FOR MEDS: **<RET>**

Select DEFAULT MEDICATION: <RET>

Select AMIS CODE: COMPUTED TOMOGRAPHY// <RET>
AMIS CODE: COMPUTED TOMOGRAPHY// <RET>

AMIS WEIGHT MULTIPLIER: 1// <RET>

BILATERAL?: <RET>

CT HEAD OR BODY?: BODY// <RET>

Select AMIS CODE: <RET>

CPT CODE// 74160 (no editing) < < < < NO EDITING ALLOWED

Select DEFAULT CPT MODIFIERS(PROC): <RET>

STAFF REVIEW REQUIRED: NO// <RET>

RAD/NM PHYS APPROVAL REQUIRED: NO// <RET>

REQUIRED FLASH CARD PRINTER: <RET>
REQUIRED FLASH CARD FORMAT: <RET>
Select FILM TYPE: DRY LASER FILM// <RET>
FILM TYPE: DRY LASER FILM// <RET>

NORMAL AMOUNT NEEDED FOR EXAM: 1// <RET>

Select FILM TYPE: <RET>

Select MESSAGE: CALL DR. SMITH BEFORE ORDERING THIS PROCEDURE

// **<RET>**

EDUCATIONAL DESCRIPTION: <RET>

No existing text Edit? NO// **<RET>**

INACTIVATION DATE: <RET>

Utility Files Maintenance Menu

CT ABDOMEN W/CONT is a descendent procedure, updating parent(s)... Updating parent: CT ABDOMEN

Updating parent: CT CHEST & ABDOMEN & PELVIS Updating parent: CT ABDOMEN & PELVIS Updating parent: CT CHEST & ABDOMEN

Select RAD/NUC MED PROCEDURES NAME:

Procedure Message Entry/Edit [RA PROCMSGEDIT]

This option is used to enter, edit and delete messages related to procedures. These messages are used at the time of request to alert the requester of special requirements related to the procedure. The requirements may be other tests that need to be done before the procedure or special approval needed to request the procedure.

Messages entered through this option are associated with procedures in the Rad/Nuc Med Procedures file using the option Procedure Enter/Edit. Multiple messages can be associated with one procedure.

Note: Procedure messages cannot be deleted from this file. Instead, use the Procedure Enter/Edit option to remove the message from the procedure.

Prompt/User Response

Procedure Message Entry/Edit

Select RAD/NUC MED PROCEDURE MESSAGE TEXT: ??

Choose from:

Patient must be NPO after midnight before test Patient must have RUQ pain, fever, elevated WBC Patient must not have eaten within 3 hours This test is not useful for assessing CHRONIC cholecystitis!

This field contains messages concerning special requirements when ordering a procedure. One or more of these messages can be tied to a procedure in the RAD/NUC MED PROCEDURES file (#71).

Select RAD/NUC MED PROCEDURE MESSAGE TEXT: Make sure the patient has not eaten for 3 hours prior to the test to avoid a false positive result!

Are you adding 'Make sure the patient has not eaten for 3 hours prior to the test to avoid a false positive result!' as a new RAD/NUC MED PROCEDURE MESSAGE (the 17TH)? Y (Yes)

TEXT: Make sure the patient has not eaten for 3 hours prior to the test to avoid a false positive result! Replace <RET>

Select RAD/NUC MED PROCEDURE MESSAGE TEXT: <RET>

Discussion

In this example, we show you how to add a new message. It must be 3 to 240 characters in length. Do not begin the message with any special characters such as "*" or "#". Use letters or numbers only.

Note: It is strongly recommended that the 1st 30 characters of the message text be unique to facilitate selection of entries.

Procedure Modifier Entry [RA MODIFIER]

This option allows the entry of new modifiers to the Procedure Modifier file (#71.2). Procedure modifiers are used to further describe the procedure. Frequently used modifiers include Right, Left, Bilateral, Portable, etc. This file provides the list of selections for the prompt "Select MODIFIERS" found in the exam edit options. An examination report will show any modifiers that have been entered for that examination.

You cannot delete a modifier. To make it inactive for an Imaging Type, type an at sign (@) at the prompt "Select TYPE OF IMAGING" for the Imaging Type you want.

If all Imaging Types are deleted from the "TYPE OF IMAGING:" field for a modifier, users will not be able to select it for any exams.

Prompt/User Response

Select TYPE OF IMAGING: <RET>

Discussion

Procedure Modifier Entry Select Procedure Modifier: CORONAL Are you adding 'CORONAL' as a new PROCEDURE MODIFIERS (the 11TH)? Y (Yes) Select TYPE OF IMAGING: ? Answer with TYPE OF IMAGING You may enter a new TYPE OF IMAGING, if you wish Enter the Imaging Type(s) for which this modifier is valid. Answer with IMAGING TYPE TYPE OF IMAGING, or ABBREVIATION Choose from: ANGIO/NEURO/INTERVENTIONAL CARDIOLOGY STUDIES (NUC MED) CT SCAN GENERAL RADIOLOGY MAGNETIC RESONANCE IMAGING NUCLEAR MEDICINE ULTRASOUND VASCULAR LAB MAMMOGRAPHY Select TYPE OF IMAGING: GENERAL RADIOLOGY

Enter a modifier 2 - 30 characters in length.

Enter one or more Imaging Types that will use this modifier.

AMIS CREDIT INDICATOR1: ?

If this modifier receives extra AMIS credit, enter the type of credit indicator it requires. Choose from:

b BILATERAL o OPERATING ROOM

p PORTABLE
AMIS CREDIT INDICATOR: <RET>

Select Procedure Modifier: <RET>

If this modifier would mean that the procedure is done bilaterally or in the operating room or with portable equipment, you should select the appropriate response here. For example, the modifier whose name is "Bilateral" should have "b" entered in this field. AMIS credit is doubled when a procedure or exam has a bilateral modifier. Operating room and portable exams are counted under both the exams AMIS category and the operating room or portable AMIS category.

June 2019

¹ AMIS CREDIT INDICATOR is no longer used in the calculation of workload. RVUs are the new metric to capture workload. Any checks for AMIS CREDIT INDICATOR at the time the order is signed in CPRS are ignored. (RA*5.0*158)

Reports Distribution Edit [RA DISTEDIT]

The distribution queue functionality of this package makes it possible to print results reports by ward, clinic, file room, etc. If a distribution queue is activated, at the time a report is verified, it will automatically be placed in the appropriate queue(s). This option allows you to edit specific report distribution queues. The types of distribution queues which can be edited are the following:

Clinic Reports
File Room
Medical Records
Other than Ward or Clinic
Requesting Physician
Ward Reports

Information which may be edited for the distribution queues include: a top of page message, the category of reports to be included in the queue, and the date after which the distribution queue should become inactivated.

Prompt/User Response

Discussion

```
Reports Distribution Edit

Select REPORT DISTRIBUTION QUEUE NAME: ??

Select the queue you want to edit.

CLINIC REPORTS
FILE ROOM
MEDICAL RECORDS
OTHER THAN WARD OR CLINIC
REQUESTING PHYSICIAN
WARD REPORTS

Select REPORT DISTRIBUTION QUEUE NAME: FILE ROOM

TOP OF PAGE MESSAGE: ***** Radiology File Room Distribution *****
Replace <RET>
```

CATEGORY OF REPORTS: ALL REPORTS// ??

This field contains the categories available for rad/nuc med reports.

Choices are: 'I' for inpatient, 'O' for outpatient, 'N' for non ward or clinic, 'A' for all reports.

Choose from:

I INPATIENT
O OUTPATIENT
N NON WARD OR CLINIC
A ALL REPORTS

CATEGORY OF REPORTS: ALL REPORTS// I INPATIENT

INACTIVATION DATE: T (AUG 29, 1994)

As long as a queue is active, when results reports are verified, the Report Distribution file will be populated with the reports that belong in the active queues. For example, if the patient is an inpatient and his report is verified, it would be placed in the Ward Reports queue, the File Room queue, and the Medical Records queue if all three are active. As long as no inactivation date is entered here, the distribution queue is active.

Note: If your facility is not using Distribution Queues, you should enter an Inactivation Date. This will avoid an accumulation of unwanted reports in the distribution queues.

Select REPORT DISTRIBUTION QUEUE NAME: <RET>

Sharing Agreement/Contract Entry/Edit [RA SHARING]

This function allows you to enter new sharing agreements/contracts into the system and edit existing ones.

A Sharing Agreement between two sites might be a situation where both facilities share the cost and maintenance of a particular piece of equipment. An example would be a mobile MRI which goes from site to site.

In Contract Agreement, one site is solely responsible for the equipment's cost and maintenance and is reimbursed for each examination performed on a patient from the other facility.

Note: A patient whose category of exam is "Contract" or "Sharing", must have an entry from this file associated with the exam.

When the category of an exam is "Contract" or "Sharing" agreement and an exam's Contract/Sharing Source field is entered, that exam will be included in the Sharing Agreement/Contract Report. Users must choose one of the agreements you enter here. To make the exam category "Contract" or "Sharing", the category of exam must be edited when the exam is ordered or during registration or case edit.

Prompt/User Response

Sharing Agreement/Contract Entry/Edit Select Agreement/Contract: ?? Choose from: LFL CONTRACTORS MEMORIAL HOSPITAL Select Agreement/Contract: UNIVERSITY HOSPITAL Are you adding 'UNIVERSITY HOSPITAL' as a new CONTRACT/SHARING AGREEMENTS (the 3RD)? Y (Yes) AGREEMENT NAME: UNIVERSITY HOSPITAL// <RET> TYPE OF AGREEMENT: ?? Choose from: S SHARING C CONTRACT TYPE OF AGREEMENT: S SHARING INACTIVATION DATE: <RET>

Discussion

Check to make sure your entry is not already in the file by entering "??" at the first prompt.

Standard Reports Entry/Edit [RA STANDRPTS]

This function allows you to enter new standard reports or to edit existing ones. These are generic reports (e.g., normal chest) that may be inserted into a procedure report. They are site specific, meaning the Imaging Services decide what reports are to be standard and what the text of the report should be. Any Division or Imaging Location may use any standard report. A standard report may be used as the entire report or may serve as a starting point for a report.

These standard reports are used by the transcriptionist in the Report Entry/Edit option of the Films Reporting Menu. If the transcriptionist selects a standard report for a specific exam, the report and impression fields will be filled in for that exam.

The Transcriptionist will only be given the "Select 'Standard' Report to Copy" prompt during the Report Enter/Edit function, if the Division parameter Allow Standard Reports is set to "YES".

Standard Report Strategies:

- 1. May be an **entire** report (e.g., "Normal Chest")
- 2. May include standard information which is in **every** report for a given procedure (e.g., how an UGI, B.E., or Nuc Med procedure was performed, what was injected, what images were obtained). Then specifics of the case will be added in the dictation. (If you are using this strategy, you do not necessarily have to put in a standard impression.)
- 3. Could have a standard report for each test for each attending, if you want, since each may have his/her own wording preference.

Prompt/User Response

Discussion

Standard Reports Entry/Edit

Select STANDARD REPORTS: NORMAL CHEST Are you adding 'NORMAL CHEST' as a new STANDARD REPORTS (the 3RD)? Y (Yes)

Enter the name of a report 3 - 30 characters in length.

STANDARD REPORTS REPORT NUMBER: 3// <RET> STANDARD REPORT: NORMAL CHEST// <RET>

If you are using the VA FileMan screen editor, this is what will appear next. A line editor will be different. Please contact IRM if you have questions concerning screen or line editing.

REPORT TEXT:

No existing text Edit? NO// YES

______ =[WRAP][INSERT]=< REPORT TEXT >[<PF1>H=Help TChest exam is within normal limits. Saving text ...

IMPRESSION:

No existing text Edit? NO// YES

______ =[WRAP][INSERT]=< IMPRESSION >[<PF1>H=Help TNormal Chest X-Ray. Saving text ...

Valid Imaging Stop Codes Edit [RA VALID STOP CODES]

The Imaging Stop Codes file #71.5 must be updated every year around October 1 when VA Headquarters MAS issues revised stop codes (also known as DSS ID) and after your IRM installs the yearly MAS maintenance patch that updates the stop codes in the Clinic Stop file #40.7. You may add and delete entries from the Imaging Stop Codes file using this option.

Only stop codes entered here will be valid selections when editing the DSS ID field in the option Location Parameter Set-Up (Imaging Locations file).

Prompt/User Response

Discussion

```
Valid Imaging Stop Codes Edit
Select IMAGING STOP CODES VALID CODE: ??
Choose from:
 NUCLEAR MEDICINE 109
      146
 RADIONUCLIDE THERAPY 144
 ULTRASOUND 115
 VASCULAR LABORATORY 421
 X-RAY 105
    Only stop codes valid for imaging should be
     entered in this file.
    When a stop code becomes invalid, it should be
    deleted from this file.
    When a new stop code becomes valid for imaging
    procedures, it should be added to this file.
    Entries are used as a selection list when
    editing procedure stop codes, and to screen the
    principal clinic edit of Imaging Locations.
```

Utility Files Maintenance Menu

Choose from:

ADMITTING/SCREENING 102

ADULT DAY HEALTH CARE 190

ALCOHOL SCREENING 706

ALCOHOL TREATMENT-GROUP 556

Your list of choices contains the entries in the Clinic Stop file #40.7.

Select IMAGING STOP CODES VALID CODE: MAMMOGRAM 703

Are you adding 'MAMMOGRAM' as a new IMAGING STOP CODES (the 5TH)? \mathbf{Y} (Yes) VALID CODE: MAMMOGRAM//

Select IMAGING STOP CODES VALID CODE: <RET>

¹VistaRad Category Entry/Edit [RA VISTARAD CATEGORY E]

This option may be displayed under the 'Utility Files Maintenance Menu' ONLY IF the VistaRad component of the Imaging package is being used; otherwise, this option will not be displayed.

This option allows the user to update only the VISTARAD CATEGORY field that is in the EXAMINATION STATUS file (#72), without having to go through all the other fields that are seen in the 'Examination Status Entry/Edit' option.

In this option, you will be asked to select one Imaging type and an examination status, then you will be prompted to enter/edit the value of the VISTARAD CATEGORY. The option will continue to ask for another examination status, until you are done with that Imaging type.

Prompt/User Response

```
VistaRad Category Entry/Edit
Select an Imaging Type: GENERAL RADIOLOGY

Select EXAMINATION STATUS: EXAMINED Imaging Type: GENERAL RADIOLOGY
Order: 3
...OK? Yes// <RET> (Yes)

VISTARAD CATEGORY: ?
```

This field is only needed for sites that will be using VistaRad for soft-copy reading of images. This information is used by VistaRad software to prepare the various types of exam lists that are displayed on the VistaRad workstation, and to properly manage exam locking for the radiologists.

If this Examination Status is to be used for exams that will be read with VistaRad, then enter a value that corresponds to it from the following list. Note that not all status codes should be assigned a VistaRad Category value, but only those that apply.

All other Exam Status codes that may be defined in the Radiology Exam Status file should NOT be entered into this field.

```
Choose from:
W Waiting for Exam

E Examined

D Dictated/Interpreted
T Transcribed

VISTARAD CATEGORY: Exami Examined

Select EXAMINATION STATUS: ^
```

¹ Patch RA 5*31 September 2002 New option added.

VIII. Maintenance Files Print Menu ...

Supervisor Menu ... Maintenance Files Print Menu ...

Complication Type List
Diagnostic Code List
Examination Status List
Film Sizes List
Label/Header/Footer Format List
Major AMIS Code List
Modifier List
Nuclear Medicine List Menu...
Procedure File Listings
Report Distribution Lists
Sharing Agreement/Contract List
Standard Reports Print
Rad/Nuc Med HL7 Voice Reporting Errors

Overview:

This menu corresponds to the Utility Files Maintenance Menu. Each option allows you to print out a dated, hard copy of the data from the related file.

Complication Type List [RA COMPRINT]

This option provides a list of the complications from the Complication Types file #78.1 and whether or not each is considered a Contrast Medium Reaction.

For more information on complication types, refer to the Index and the chapter on Complication Type Entry/Edit.

Prompt/User Response

Discussion

Complication Type List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

Complication Types	CONTRAST MEDIUM	SEP	30,1996	08:52	PAGE 1
COMPLICATION	REACTION?				
CONTRAST REACTION	YES				
NAUSEA AND VOMITING	NO				
NO COMPLICATION	NO				
PREGNANCY	NO				
SEVERE MORBIDITY	NO				

Diagnostic Code List [RA DIAGP]

This option lists the contents of the Diagnostic Codes file #78.3.

For more information diagnostic codes, refer to the Index and the chapter on Diagnostic Code Enter/Edit.

Prompt/User Response

Discussion

Diagnostic Code List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen. See the next page for an example of a list.

Diagnostic Codes SEP 30,1996 08:53 PAGE 1 PRINT					
NUMBER DIAGNOSTIC CODE	ON ABNORMAL REPORT	GENERATE ABNORMAL ALERT?			
1 NORMAL There were no problems for this exam!	NO				
2 MINOR ABNORMALITY There was a slight abnormality in this exam!	NO				
3 MAJOR ABNORMALITY, NO ATTN. NEEDED There was a major problem, but no attention w	NO ras needed!				
4 ABNORMALITY, ATTN. NEEDED There was a major problem, and attention is n	YES ecessary!	YES			
5 MAJOR ABNORMALITY, PHYSICIAN AWARE There is a major problem, and the doctor is a	NO ware of it!				
6 UNDICTATED FILMS NOT RETURNED, 3 DAY	S NO				
7 UNSATISFACTORY/INCOMPLETE EXAM This exam was either unsatisfactory or incomp					
8 POSSIBLE MALIGNANCY, FOLLOW-UP NEEDE A possible malignancy was noted. Follow-up a	D YES ction should be taken.				
9 POSSIBLE MALIGNANCY, FOLLOW-UP CRITI Follow-up must be done, and documented					
10 OVEREXPOSURE, MECH ERROR OVEREXPOSURE DUE TO A MECHANICAL ERROR	YES	NO			

Note: A blank in the "Print on Abnormal Report" column or the "Generate Abnormal Alert?" column is the same a "NO".

Examination Status List [RA EXAMSTATUSP]

This option lists the contents of the Examination Status file #72. It lists all contents of the Examination Status file whether or not the Imaging Types and statuses are activated. Our example shows only one of the Imaging Types that has been activated.

For more information on exam statuses, refer to the index and the chapter on Examination Status Entry/Edit.

Prompt/User Response

Examination Status List

Select Imaging Type: NUCLEAR MEDICINE

Another one (Select/De-Select): <RET>

Select Device: (Printer Name)

Discussion

This report now allows printing by selected Imaging Types. Enter the name of a printer or press the <RET> key to display the output on your screen.

See the next three pages for an example of the output.

Examination Statuses Page: 1

Run Date: Jan 27, 1997 8:45:12 am

Type Of Imaging: NUCLEAR MEDICINE

Status: ***CANCELLED*** Order: 0

Default Next Status: User Key Needed:

Generate Examined HL7 Message:

Generate Exam Alert: Allow Canceling?: YES
Appear On Status Tracking?: NO Print Dosage Ticket: NO

ASK ON STATUS TRACKING: REQUIRED FOR CHANGE TO THIS STATUS:

WORKLOAD REPORTS THAT USE THIS STATUS IN ITS COMPLETION:

AMIS

Status: ***WAITING FOR EXAM*** Order: 1

Default Next Status: EXAMINED User Key Needed:

Generate Examined HL7 Message:

Generate Exam Alert: NO Allow Canceling?: YES
Appear On Status Tracking?: YES Print Dosage Ticket: NO

ASK ON STATUS TRACKING: REQUIRED FOR CHANGE TO THIS STATUS:

WORKLOAD REPORTS THAT USE THIS STATUS IN ITS COMPLETION:

CLINIC

PTF BEDSECTION

SERVICE

SHARING/CONTRACT

WARD

FILM USAGE

TECHNOLOGIST

AMIS

DETAILED PROCEDURE

CAMERA/EQUIP/RM

PHYSICIAN

RESIDENT

STAFF

DELINQUENT STATUS

Examination Statuses Page: 2

Run Date: Jan 27, 1997 8:45:12 am

Status: ***EXAMINED*** Order: 2

Default Next Status: TRANSCRIBED User Key Needed:

Generate Examined HL7 Message: YES

Generate Exam Alert: NO Allow Canceling?: NO Appear On Status Tracking?: YES Print Dosage Ticket: YES

ASK ON STATUS TRACKING: REQUIRED FOR CHANGE TO THIS STATUS:

-----______

TECHNOLOGIST TECHNOLOGIST

PROCEDURE DETAILED PROCEDURE FILM DATA FILM ENTRY

CAMERA/EQUIP/RM CAMERA/EQUIP/RM MEDICATIONS & DOSAGES RADIOPHARMS/DOSAGES

MED ADM DT/TIME/PERSON RADIOPHARMS AND DOSAGES

ADM DT/TIME/PERSON LOT NO.

WORKLOAD REPORTS THAT USE THIS STATUS IN ITS COMPLETION:

CLINIC

PTF BEDSECTION

SERVICE

SHARING/CONTRACT

WARD

FILM USAGE

TECHNOLOGIST

DETAILED PROCEDURE

CAMERA/EQUIP/RM

PHYSICIAN

RESIDENT

STAFF

DELINQUENT STATUS

Status: ***TRANSCRIBED*** Order: 3

Default Next Status: COMPLETE User Key Needed:

Generate Examined HL7 Message:

Generate Exam Alert: YES Allow Canceling?: NO Print Dosage Ticket: NO Appear On Status Tracking?: NO

ASK ON STATUS TRACKING: REQUIRED FOR CHANGE TO THIS STATUS:

_____ ______

TECHNOLOGIST TECHNOLOGIST INTERPRETING PHYS REPORT ENTERED

PROCEDURE IMPRESSION FILM DATA RESIDENT OR STAFF

DIAGNOSTIC CODE DETAILED PROCEDURE

CAMERA/EQUIP/RM FILM ENTRY MEDICATIONS & DOSAGES CAMERA/EQUIP/RM VOLUME/FORM

Examination Statuses Page: 3

Run Date: Jan 27, 1997 8:45:12 am

MED ADM DT/TIME/PERSON
RADIOPHARMS AND DOSAGES
ACTIVITY DRAWN
DRAWN DT/TIME/PERSON
ADM DT/TIME/PERSON
ROUTE/SITE OF ADM
LOT NO.

RADIOPHARMS/DOSAGES

WORKLOAD REPORTS THAT USE THIS STATUS IN ITS COMPLETION:

CLINIC

PTF BEDSECTION

SERVICE

SHARING/CONTRACT

WARD

FILM USAGE

TECHNOLOGIST

AMIS

DETAILED PROCEDURE

CAMERA/EQUIP/RM

PHYSICIAN

RESIDENT

STAFF

DELINOUENT STATUS

Status: ***COMPLETE***

Default Next Status:

Generate Examined HL7 Message:

Generate Exam Alert: NO

Appear On Status Tracking?: NO

ASK ON STATUS TRACKING:

TECHNOLOGIST
INTERPRETING PHYS

PROCEDURE

FILM DATA

DIAGNOSTIC CODE

CAMERA/EQUIP/RM

MEDICATIONS & DOSAGES

MED ADM DT/TIME/PERSON

RADIOPHARMS AND DOSAGES

ACTIVITY DRAWN

DRAWN DT/TIME/PERSON

ADM DT/TIME/PERSON

ROUTE/SITE OF ADM

LOT NO.

VOLUME/FORM

Order: 9

User Key Needed:

Allow Canceling?: NO

Print Dosage Ticket: NO

REQUIRED FOR CHANGE TO THIS STATUS:

REQUIRED FOR CHANGE TO THIS STATUS:

TECHNOLOGIST

REPORT ENTERED VERIFIED REPORT

IMPRESSION

RESIDENT OR STAFF

DETAILED PROCEDURE

FILM ENTRY

CAMERA/EQUIP/RM

RADIOPHARMS/DOSAGES

Page: 4

Examination Statuses Run Date: Jan 27, 1997 8:45:12 am

WORKLOAD REPORTS THAT USE THIS STATUS IN ITS COMPLETION:

CLINIC

PTF BEDSECTION

SERVICE

SHARING/CONTRACT

WARD

FILM USAGE

TECHNOLOGIST

AMIS

DETAILED PROCEDURE

CAMERA/EQUIP/RM

PHYSICIAN

RESIDENT

STAFF

DELINQUENT STATUS

Film Sizes List [RA FILMP]

This option lists the contents of the Film Sizes file #78.4.

For more information on films, refer to the Index and the chapter on Film Type Entry/Edit.

Prompt/User Response

Discussion

Film Sizes List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

Film Size Specif	ficatio	ns			SEP 30,1996	08:55	PAGE 1
-		F		W	·		
		1		а			
	С	u		S			
	i	0		t	Analogous		
	n	r		е	Unwasted Film		
Film	е	o Ina	ctivation	Dt d	Size		
10X12	N						
11X14	N						
14X17	N						
2X12	N						
2X4	N	JUL	5,1996				
2X6	N						
3X FILMS	N						
CHROMATIC	N						
CINE	Y						
DENTAL	N						
FLUORO ONLY	N						
NO FILMS	N						
W-10X12	N			Y	10X12		
W-11X14	N			Y	11X14		
W-14X17	N			Y	14X17		
W-2X12	N			Y	2X12		
W-2X4	N	JUL	5,1996	Y	2X4		
W-2X6	N			Y	2X6		
W-3X FILMS	N			Y	3X FILMS		
W-CHROMATIC	N			Y	CHROMATIC		
W-CINE	Y			Y	CINE		
W-DENTAL	N			Y	DENTAL		
W-FLUORO ONLY	N			Y	FLUORO ONLY		

Label/Header/Footer Format List [RA FLASHFORMP]

This option prints the contents of the Lbl/Hdr/Ftr Formats file #78.2.

For more information on formats, refer to the Index and the chapter on Label/Header/Footer Formatter.

Prompt/User Response

Discussion

Label/Header/Footer Format List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

```
Exam Label/Report Header/Report Footer/Flash Card Formats

SEP 30,1996 08:55 PAGE 1

FIELD ROW COL TITLE

TEXT

FORMAT NAME: ER FORMAT

NAME OF PATIENT 1 2 NAME :
SSN OF PATIENT 1 40 SSN :
DATE OF EXAM 2 1 EXAM DATE:
TECHNOLOGIST 2 40 TECH:
PROCEDURE 3 1 PROCEDURE:
FREE TEXT 5 2 COMMENTS:
VERIFYING RADIOLOGIST 10 2 INTERPRETER:
```

```
Exam Label/Report Header/Report Footer/Flash Card Formats
SEP 30,1996 08:55 PAGE 2
FIELD ROW COL TITLE

FORMAT NAME: EXAM HEADER
FREE TEXT
1 1 EXAM REPORT
NAME OF PATIENT
CURRENT DATE/TIME
SSN OF PATIENT
2 1 SSN:
ATTENDING PHYSICIAN AT ORDER 2 40 Att at Order:

REQUESTING PHYSICIAN 0 3 1 REQ PHYS:
ATTENDING PHYSICIAN CURRENT 3 40 Att Current:

REQUEST ENTERED DATE/TIME 4 1
CURRENT PATIENT LOCATION 5 1
```

Major AMIS Code List [RA MAJORAMISP]

This option lists the contents of the Rad/Nuc Med AMIS Codes file #71.1.

For more information on AMIS codes, refer to the Index and the chapter on Major AMIS Code Entry/Edit.

Prompt/User Response

Discussion

Major AMIS Code List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

Major MAJOR	AMIS Codes	SEP 30,1996 08:56 PAGE 1
AMIS		
	DESCRIPTION	WEIGHT
CODE	DESCRIPTION	WEIGHI
1	SKULL, INC. SINUS, MASTOID, JAW, ETC	3
2	CHEST-SINGLE VIEW	1
3	CHEST MULTIPLE VIEW	2
4	CARDIAC SERIES	3
5	ABDOMEN-KUB	2
6	OBSTRUCTIVE SERIES	3
7	SKELETAL-SPINE & SACROILIAC	3
8	SKELETAL-BONE & JOINTS	2
9	GASTROINTESTINAL	6
10	GENITOURINARY	6
11	CHOLECYSTOGRAM, ORAL	5
12	CHOLANGIOGRAM	10
13	LAMINOGRAM	5
14	BRONCHOGRAM	5
15	DIGITAL SUBTRACTION ANGIOGRAPHY	15
16	ANGIOGRAM, CATH- CEREBRAL	15
17	ANGIOGRAM, CATH- VISCERAL	20
18	ANGIOGRAM, CATH- PERIPHERAL	10
19	VENOGRAM	15
20	MYELOGRAM	10
21	COMPUTED TOMOGRAPHY	8
22	INTERVENTIONAL RADIOGRAPHY	20
23	ULTRASOUND, ECHOENCEPHALOGRAM	7
24	OTHER	5
25	EXAMS IN OPER.SUITE AT SURGERY	1
26	PORTABLE (BEDSIDE) EXAMINATIONS	1
27	NUCLEAR MEDICINE	1

Modifier List [RA MODIFIERP]

This option lists the contents of the Procedure Modifiers file #71.2.

For more information on modifiers, refer to the Index and the chapter on Procedure Modifier Entry.

Prompt/User Response

Discussion

Modifier List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

Procedure Modifiers NOV 18,1996 13:52 PAGE 1 AMIS CREDIT NAME INDICATOR TYPE OF IMAGING: CARDIOLOGY STUDIES (NUC MED) RIGHT TYPE OF IMAGING: CT SCAN RIGHT TYPE OF IMAGING: GENERAL RADIOLOGY BILATERAL BILATERAL EXAM PORTABLE EXAM PORTABLE OPERATING ROOM EXAM OPERATING ROOM LEFT RIGHT

PORTABLE

BILATERAL

OPERATING ROOM

OPERATING ROOM

TEST

TEST OR

NEW ONE

OBLIQUE

TEST BILATERAL

Maintenance Files Print Menu

Supervisor Menu ... Maintenance Files Print Menu ...

Nuclear Medicine List Menu ...

Lot (Radiopharmaceutical) Number List Route of Administration List Site of Administration List Vendor/Source (Radiopharmaceutical) List

Use this menu in conjunction with the Nuclear Medicine Setup Menu.

™ Lot (Radiopharmaceutical) Number List

[RA NM PRINT LOT]

This option prints a list of entries in the Radiopharmaceutical Lot file #71.9. The file is exported without data. It prints in 132 column format.

RADIOPHARMACEUTICAL LOT NUMBER/ID	RADIOPHARMACEUTICAL KIT #		MAY 30,1997 EXPIRATION DATE	12:59 PAGE 1 INACTIVE	
SOURCE: CHEMICALS INC.		123598Y34512548			
1234	TECHNETIUM Tc99m SESTAMIBI		JAN 20,1997		
SOURCE: DUPONT					
PER 1234	PERCHLORACAP 250MG CAPS	235			
456	SULFUR COLLOID TC-99M	445	DEC 31,1997	YES	
MAA 7136 7635	TECHNETIUM Tc99m ALBUMIN AGGREGATED TECHNETIUM Tc99m MEDRONATE	1-5			
MDP 6136P	TECHNETIUM Tc99m MEDRONATE	100		YES	
MDP 6137P	TECHNETIUM Tc99m MEDRONATE	100	JAN 20,1997		
MDP 6137P	TECHNETIUM Tc99m MEDRONATE	101	JAN 13,1997		
DTPA 1234	TECHNETIUM Tc99m PENTETATE				
SOURCE: SQUIBB					
12345	TECHNETIUM Tc99m ALBUMIN AGGREGATED	333	JAN 22,1997		
MDP 6137P	TECHNETIUM Tc99m MEDRONATE	1			

Route of Administration List

[RA NM PRINT ROUTE]

This option prints a list of entries in the Route of Administration file #71.6. A number of routes were exported with this version of the software.

ROUTE OF ADMINISTRATION		JAN 21,1997 08:36 PAGE 1
NAME		
SITES OF ADMIN		PROMPT
CATHETER		
URINARY BLADDER		
G-TUBE	(synonym)	GASTROSTOMY TUBE
		FEEDING TUBE
INHALATION	(synonym)	VENTILATOR
MOUTH		
TRACHEAL TUBE		
NOSE		
ENDOTRACHEAL TUBE		
TRACHEOSTOMY		
INTRA-ARTICULAR	(synonym)	JOINT
		ARTICULAR
LEFT HIP		
RIGHT HIP		
LEFT KNEE		
RIGHT KNEE		
LEFT WRIST		
RIGHT WRIST		
LEFT SHOULDER		
RIGHT SHOULDER		Yes
INTRADERMAL	(synonym)	
INTRALESIONAL	(-11-17	Yes
INTRAMUSCULAR	(synonym)	
111111111111111111111111111111111111111	(synonym)	
RIGHT FOREARM	(0]11011]111)	
LEFT FOREARM		
LEFT SHOULDER		
RIGHT SHOULDER		
LEFT HIP		
RIGHT HIP		
INTRAOCULAR	(synonym)	EYE
RIGHT EYE	(Symonym)	
LEFT EYE		
INTRAPERITONEAL	(synonym)	PERITONEAL
INTRAPLEURAL	(synonym)	
LEFT PLEURA	(Symonym)	LUDOIM
RIGHT PLEURA	(aimonim)	CDIMAI
INTRATHECAL LUMBAR SUBARACHNOID SPACE	(synonym)	SLIMMT
CERVICAL SUBARACHNOID SPACE		

Site of Administration List

[RA NM PRINT SITE]

This option prints a list of entries in the Site of Administration file #71.7. A number of sites were exported with this version of the software.

SITE OF ADMINISTRATION	JAN 21,1997 08:36 PAGE 1
NAME	(SYNONYM)
CERVICAL SUBARACHNOID SPACE	SUBARACHNOID SPACE
ENDOTRACHEAL TUBE	ETT
	E.T.T.
GALLBLADDER	GB
	G.B.
INTRAOCULAR	EYE
LEFT ANTECUBITAL FOSSA	L. ANTECUBITAL FOSSA
	L ANTECUBITAL FOSSA
	ANTECUBITAL FOSSA
LEFT ARM	L ARM
	L. ARM
	ARM
LEFT BUTTOCK	L BUTTOCK
	L. BUTTOCK
	BUTTOCK
LEFT DELTOID	L DELTOID
	L. DELTOID
	DELTOID
LEFT EYE	L. EYE
	L EYE
	EYE
LEFT FEMORAL	L. FEMORAL
	L FEMORAL
	FEMORAL
	L. GROIN
	L GROIN
	GROIN
LEFT FOOT	L. FOOT
	L FOOT
	FOOT
LEFT FOREARM	L. FOREARM
	L FOREARM
	FOREARM
LEFT GLUTEUS	L GLUTEUS
	L. GLUTEUS
	GLUTEUS
LEFT HAND	L. HAND
	L HAND
	HAND

Vendor/Source (Radiopharmaceutical) List

[RA NM PRINT SOURCE]

This option prints a list of entries in the Radiopharmaceutical Source file #71.8. The file is exported without data.

RADIOPHARMACEUTICAL SOURCE	JAN 21,1997	08:37	PAGE 1
VENDOR/PHARMACY			
CHEMICALS INC.			
DUPONT			
NATURAL CHEMICALS LABORATORY			
PHARM-A-LOGIC (OFF SITE)			
SQUIBB			
UNIVERSITY LABS			
UNIVERSITY NUCLEAR PHARMACY			

Supervisor Menu ... Maintenance Files Print Menu ... Procedure File Listings ...

¹HIST Display Rad/NM Procedure Contrast Media History

Active Procedure List (Long)
Active Procedure List (Short)
Alpha Listing of Active Procedures
Barcoded Procedure List
Display Rad/NM Procedure Contrast Media History
Inactive Procedure List (Long)
Invalid CPT/Stop Code List
List of Inactive Procedures (Short)

List of Procedures with Contrast

Parent Procedure List Procedure Message List Series of Procedures List

All of these reports access data in the Rad/Nuc Med Procedures file #71. Many of the following printouts are lengthy. "Q"ueuing these reports to be printed after hours should be considered.

-

¹ Patch RA*5*45 March 2004 Associating parent/descendant that use contrast media.

¹HIST Display Rad/NM Procedure Contrast Media History

[RA CMAUDIT HISTORY]

This option allows the user to look back at the prior contrast medium/media data associated with a Rad/Nuc Med Procedure.

1. From the Menu below, select the Supervisor Menu by typing 'Supervisor Menu' next to Select Rad/Nuc Med Total System Menu Option:

```
Exam Entry/Edit Menu ...
Films Reporting Menu ...
Management Reports Menu ...
Outside Films Registry Menu ...
Patient Profile Menu ...
Radiology/Nuclear Med Order Entry Menu ...
Supervisor Menu ...
Switch Locations
Update Patient Record
User Utility Menu ...
Select Rad/Nuc Med Total System Menu Option: Supervisor Menu
```

2. From the Supervisor Menu, select the Maintenance Files Print Menu by typing 'Maintenance Files Print Menu' next to Select Supervisor Menu Option:

```
Select Rad/Nuc Med Total System Menu Option: Supervisor Menu
   HL7
          Radiology HL7 Menu
          Access Uncorrected Reports
          Delete a Report
          Delete Printed Batches By Date
          Exam Deletion
          Inquire to File Entries
          List Exams with Inactive/Invalid Statuses
          Maintenance Files Print Menu
          Mass Override Exam Status
          Override a Single Exam's Status to 'complete'
          Print File Entries
          Rad/Nuc Med Personnel Menu ...
          Search File Entries
          Switch Locations
          System Definition Menu ...
          Unverify a Report for Amendment
          Update Éxam Status
          Utility Files Maintenance Menu ...
Select Supervisor Menu Option: Maintenance Files Print Menu 🤦
```

¹ Patch RA*5*45 December 2004 New Procedure File Listings: HIST Display Rad/NM Procedure Contrast Media History

3. From the Maintenance Files Print Menu, select Procedure File Listings by typing '**Procedure File Listings**' next to Select Maintenance Files Print Menu Option:

```
Select Supervisor Menu Option: Maintenance Files Print Menu

Complication Type List
Diagnostic Code List
Examination Status List
Film Sizes List
Label/Header/Footer Format List
Major AMIS Code List
Modifier List
Nuclear Medicine List Menu ...
Procedure File Listings ...
Report Distribution Lists
Sharing Agreement/Contract List
Standard Reports Print
VistaRad Category Print

Select Maintenance Files Print Menu Option: Procedure File Listings
```

4. From the Procedure file Listings Menu, select HIST Display Rad/NM Procedure Contrast Media History by typing 'HIST Display Rad/NM Procedure Contrast Media History' next to Select Procedure File Listings Option:

5. The Select Procedure Prompt displays. Enter a specific procedure, (for example, CT PELVIS W/CONT) or type 'ALL' to select all available procedures.

```
HIST
         Display Rad/NM Procedure Contrast Media History
         Active Procedure List (Long)
         Active Procedure List (Short)
         Alpha Listing of Active Procedures
         Barcoded Procedure List
         Inactive Procedure List (Long)
         Invalid CPT/Stop Code List
         List of Inactive Procedures (Short)
         List of Procedures with Contrast
         Parent Procedure List
         Procedure Message List
         Series of Procedures List
Select Procedure File Listings Option: hist Display Rad/NM Procedure Contrast M
edia History
Select Procedure: CT PELVIS W/CONT
```

- 6. At the prompt: 'Another one (Select/De-Select),' enter another procedure or <RET> to continue.
- 7. At the appropriate prompt, enter the start and end dates for the search.
- 8. At the **DEVICE HOME**// prompt, enter the device on which to print the data or **<RET>** for the default (screen).

```
Select Procedure: CT PELVIS W/CONT (CT Detailed) CPT:72193

Another one (Select/De-Select):

Enter the start date for the search: 1/11/11 (JAN 11, 1911)

Enter the ending date for the search: 11/03/04 (NOV 03, 2004)

DEVICE: HOME// TELNET
```

Sample Contrast Media Edit History by Procedure Report

Contrast Media Edit History By Procedure Run Date: Dec 16, 2004 From: JAN 11,1919 To: NOV 3,2004 Page 1 Date/Time Changed User Contrast Media ______ CT PELVIS W/CONT Oct 27, 2004 1:28 pm VISTPATIENT, ONE Gastrografin Oct 27, 2004 1:29 pm VISTPATIENT, TWO CT PELVIS W/CONT **User deleted all contrast media data** Oct 27, 2004 2:39 pm VISTPATIENT, THREE CT PELVIS W/CONT Gastrografin, Barium CT PELVIS W/CONT Oct 27, 2004 2:40 pm VISTPATIENT, FOUR **User deleted all contrast media data** CT PELVIS W/CONT Oct 27, 2004 2:41 pm VISTPATIENT, FIVE Gastrografin, Barium CT PELVIS W/CONT Oct 27, 2004 2:41 pm VISTPATIENT, SIX Enter RETURN to continue or '^' to exit:

Active Procedure List (Long) [RA PROCLONG]

This option prints a listing of all active procedures, their CPT and AMIS codes, AMIS weight multiplier, whether it is a bilateral procedure or a CT of the head or body, the type of procedure (broad, series, detailed or parent), flash card information, and whether review of the report and approval of the request is required. The report is sorted by AMIS category, then procedure name. It is a lengthy printout, so queuing the output should be considered.

¹If you assign one or more modalities to a procedure, the abbreviated and long name of each modality will print as the last item for the procedure.

²Note: Contrast media associations will display, when appropriate.

Prompt/User Response

Discussion

```
Active Procedure List (Long)

Select Imaging Type: All// NUCLEAR MEDICINE

Another one (Select/De-Select): <RET>

This report requires a 132 column output device.

DEVICE: (Printer Name or "Q")
```

This report now allows printing by selected Imaging Types.

If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

```
Active Radiology/Nuclear Medicine Procedures (Long)
                                                                                          DEC 13,1996 08:43
                                                                                                                  PAGE 1
         PROCEDURE / 3CPT MODIFIER
                                                                                                                 HEAD/BODY
CODE
                                                      AMIS CODE
                                                                                MULTIPLIER
                                                                                                BILATERAL
        ACUTE GI BLOOD LOSS IMAGING
                                              27 - NUCLEAR MEDICINE
     Type of Procedure
                                        DETAILED
     Required Flash Card Printer:
                                         ER FORMAT
     Required Flash Card Format :
     Cost of Procedure
    Health Summary With Request:
Staff Review of Reports Req:
                                         TNPATTENT
     Rad Approval of Request Req:
                                         NUCLEAR MEDICINE
     Type of Imaging
     Suppress Radiopharm Prompt:
                                         Always
     Prompt for Radiopharm RX
     Default Radiopharmaceutical:
                                         PERCHLORACAP 250MG CAPS
     Usual Dose :
Dflt Administration Route :
Dflt Administration Site :
                                          ORAL
     High Adult Dose
      Low Adult Dose
      Default Form
                                          Solid (Other)
     Prompt for Medication?
     Descendents
     Procedure Message
     Educational Description
                                         This test is designed to LOCATE the site of KNOWN GI
                                            bleeding, NOT to determine whether there IS GI bleeding.
```

¹ Patch RA*5*3

² Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

³ Patch RA*5*19 May 2000 "/ CPT MODIFIER" added to report header.

Active Procedure List (Short) [RA PROCSHORT]

This is a list of active procedures, sorted by AMIS category and procedure name. Since the report can be lengthy, consider queuing the output. Only procedures with an AMIS code are included.

¹Note: Contrast media associations will display, when appropriate.

Prompt/User Response

Discussion

Active Procedure List (Short)

Select Imaging Type: All// GENERAL RADIOLOGY

Another one (Select/De-Select): <RET>

This report requires a 132 column output device.

DEVICE: (Printer Name or "Q")

This report now allows printing by selected Imaging Types.

Enter the name of a printer that prints 132 columns. If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

	Radiology/Nuclear Medicine Procedures	(Short)			30,1996 10:5	4 PAGE 1
CPT CODE	PROCEDURE	AMIS CODE	MULT	BI- LAT	IMG CT TYPE	MODALITY
76062	BONE SURV COMP (INCL APPENDIC 1	-SKULL, INC.SINUS, MASTOID, JAW, ET	1		RAD	
76061	BONE SURV LMTD (E.G. METASTATI 1	-SKULL, INC.SINUS, MASTOID, JAW, ET	1		RAD	

-

¹ Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

Alpha Listing of Active Procedures [RA ALPHALIST]

As you can see by the example below, this is a shorter, alphabetical listing of the active procedures. Since the report can be lengthy, consider queuing the output.

¹Note: Contrast media associations will display, when appropriate.

Prompt/User Response

Discussion

Alpha Listing of Active Procedures

Select Imaging Type: All// GENERAL RADIOLOGY

Another one (Select/De-Select): <RET>

This report requires a 132 column output device.

DEVICE: (Printer Name or "Q")

This report now allows printing by selected Imaging Types.

Enter the name of a printer that prints 132 columns. If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

Alpha Listing of Active Procedures SEP 30,1996 10:54 PAGE 1							
NAME			AMIS CODE	MULT		TYPE MODALITY	
ABDOMEN 1 VIEW	74000	5	ABDOMEN-KUB	1		RAD	
ABDOMEN 2 VIEWS	74010	5	ABDOMEN-KUB	1		RAD	
ABDOMEN FOR FETAL AGE 1 VIEW	74720	5	ABDOMEN-KUB	1		RAD	
ABDOMEN FOR FETAL AGE MULT VIEWS	74725	5	ABDOMEN-KUB	1		RAD	
ABDOMEN MIN 3 VIEWS+CHEST	74022	2 6	CHEST-SINGLE VIEW OBSTRUCTIVE SERIES	1 1		RAD	
ACROMIOCLAVICULAR J BILAT	73050	8	SKELETAL-BONE & JOINT	S 1	YES	RAD	

-

¹ Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

Barcoded Procedure List

[RA BARPROCPRINT]

If your site uses barcodes, use this option to print a list of barcoded CPTs and/or Procedures for a selected Imaging Type.

Prompt/User Response

Discussion

Barcoded Procedure List

If you don't know the height of the barcode, you need to print a test version and count the lines as shown here:

To print barcoded procedure list, you will need to know the height (in vertical lines) of the barcode output on the printer to be used.

Do you wish to print a sample barcode for the purpose of determining the height (in vertical lines) of the barcode? No// ${\bf Y}{\rm ES}$

If you don't know the height of the barcodes printed at your site, enter YES to obtain a sample barcode. You may need to contact IRM for the name of a printer for barcodes.

Select a printer with barcode setup: (Enter the name of the device) DO YOU WANT YOUR OUTPUT QUEUED? NO// \mathbf{Y} (YES)

Requested Start Time: NOW// <RET> (DEC 11, 1996@12:46:25)
Request Queued, Task #: 35507

Enter the height of the barcode: ^

Exit the option and get the printout of the test barcode. Here's an example:

	Barcode test print					
			Printed on: Dec 11, 1996	Page 1		
LINE 1	 1					
LINE 2						
LINE 3						
LINE 4						
LINE 5	5					
LINE 6	6					
LINE 7	7					
LINE 8	8					
LINE 9	9					
LINE 1	10					
TEST E	BARCODE	PRINT	BarcodeBarcodeBarcodeBarcode			
			BarcodeBarcodeBarcodeBarcode			
			BarcodeBarcodeBarcodeBarcode			
			TEST BARCODE PRINT			
LINE 1						
LINE 2						
LINE 3	_					
LINE 4						
LINE 5						
LINE 6						
LINE 7	'					
LINE 8						
LINE 9						
LINE 1	TU					

Using the line counts printed above and below the example barcode to determine how many vertical lines the barcode occupies (include the procedure name line printed below the barcode), it appears to be 4 lines. Now you can print a list of barcodes used by your Imaging Type. Return to the option.

To print barcoded procedure list, you will need to know the height (in vertical lines) of the barcode output on the printer to be used.

Do you wish to print a sample barcode for the purpose of determining the height (in vertical lines) of the barcode? No// $\langle RET \rangle$

Enter the height of the barcode: 4

Select one of the following:

- C CPT barcode
- P Procedure barcode
- B Both barcodes

Which of the above would you like to print?: Both

Select the data element you wish to sort by: ??

Enter 'C' to sort by the CPT Code, or 'P' to sort by the procedure.

Select one of the following:

- C CPT Code
- P Procedure

Select the data element you wish to sort by: Procedure

Select Imaging Type: GENERAL RADIOLOGY

Another one (Select/De-Select): <RET>

Select Procedure: ALL

Another one (Select/De-Select): <RET>

Select a printer with barcode setup: (Enter the barcode printer name) DO YOU WANT YOUR OUTPUT QUEUED? NO// \mathbf{Y} (YES)

Requested Start Time: NOW// <RET> (DEC 11, 1996@12:47:58)

Request Queued, Task #: 35508

Inactive Procedure List (Long) [RA INACPRCLONG]

This output is the same as that for Active Procedure List (Long) only these procedures have all been inactivated. Since the report can be lengthy, consider queuing the output.

¹If you assigned one or more modalities to a procedure, the abbreviated and long name of each modality will print as the last item for the procedure.

²Note: Contrast media associations will display, when appropriate.

Prompt/User Response

Discussion

Inactive Procedure List (long)
Select Imaging Type: All// GENERAL RADIOLOGY
Another one (Select/De-Select): <RET>
This report requires a 132 column output device.
DEVICE: (Printer Name or "Q")

This report now allows printing by selected Imaging Types.

Enter the name of a printer that prints 132 columns. If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

```
Inactive Radiology/Nuclear Medicine Procedures (Long)
                                                                                    SEP 30,1996 10:54 PAGE 1
                                                                             MULTIPLIER BILATERAL HEAD/BODY
CODE
        PROCEDURE / <sup>3</sup>CPT MODIFIER
                                              AMIS CODE
70350 CEPHALOGRAM ORTHODONTIC
                                         1 -SKULL, INC. SINUS, MASTOID, JAW, ET
      Type of Procedure
                                     DETAILED
      Required Flash Card Printer:
      Required Flash Card Format :
      Cost of Procedure
      Health Summary With Request:
      Staff Review of Reports Req:
      Rad Approval of Request Req:
      Type of Imaging
                                       GENERAL RADIOLOGY
      Prompt for Medications
      Descendents
      Procedure Message
      Educational Description
```

¹ Patch RA*5*3

² Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

³ Patch RA*5*19 May 2000 "/ CPT MODIFIER" added to report header.

Invalid CPT/Stop Code List [RA INVALID CPT/STOP]

This option prints a list of Invalid CPT or Stop Codes. When generating the report, you can choose to include procedures that already have been marked Inactive (by entering an inactive date through the option Procedure Enter/Edit). It should be noted that the system sees procedure and CPTs as two separate entities. In other words, an active procedure can have a valid or invalid CPT associated with it. This report will print a warning if a Parent procedure has no Descendents, and if a detailed or series procedure has no AMIS code(s). It will also tell you if an Imaging Location setup doesn't meet criteria for workload crediting.

Prompt/User Response

Discussion

Invalid CPT/Stop Code List

This option prints a list of Radiology/Nuclear Medicine Procedures with missing or invalid CPT's and DSS ID's, and Imaging Locations pointing to Hospital Locations with inappropriate set-up parameters. Broad, parent and inactive procedures are not required to have codes. To be valid, DSS ID's must be in the Imaging Stop Codes file 71.5; CPT's must be nationally active.

Include Inactive procedures? NO// <RET>

Select Imaging Type: All// GENERAL RADIOLOGY

Another one (Select/De-Select): <RET>

DEVICE: (Printer Name or "O")

This report now allows printing by selected Imaging Types.

Enter the name of a printer or press the <RET> key to bring the output to the screen.

If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

```
RADIOLOGY/NUCLEAR MEDICINE INVALID CPT AND STOP CODES
Run Date: Sep 30, 1996@13:41:04
                                                                      Page: 1
ABDOMEN FOR FETAL AGE 1 VIEW (Detailed, RAD)
    Invalid CPT 74720
ABDOMEN FOR FETAL AGE MULT VIEWS (Detailed, RAD)
    Invalid CPT 74725
ANGIO CERVICOCEREBRAL SELECT 1 VESSEL S&I (Detailed, RAD)
    No CPT entered.
    No stop code(s) entered.
ANGIO CERVICOCEREBRAL SELECT 2 VESSEL S&I (Type missing, Imaging Type
missing)
    Invalid CPT 75654
     . . . Etc.
Imaging Type(s): GENERAL RADIOLOGY
TOTAL INVALID CPT CODES: 53
TOTAL MISSING CPT CODES: 1
TOTAL BROKEN POINTERS TO CPT CODES: 0
TOTAL PARENT PROCEDURES W/O DESCENDENTS: 0
TOTAL SERIES AND/OR DETAILED PROCEDURES W/O AMIS CODES: 3
Note: Remember to review the entries in the Imaging Stop Codes file #71.5
every year in October to make sure they agree with VA HQ rules. If the
entries in file 71.5 are not accurate, this report will not be accurate.
Imaging Location: MAMMOGRAPHY
  700-series noncredit Stop Code being used
```

List of Inactive Procedures (Short) [RA INACPRCSHORT]

This output is the same as Inactive Procedure List (Long) excluding the Type of Procedure, Required Flash Card Printer, Required Flash Card Format, Staff Review of Reports Req, and Rad Approval of Request Req fields. Since the report can be lengthy, consider queuing the output.

¹Note: Contrast media associations will display, when appropriate.

Prompt/User Response

List of Inactive Procedures (Short)

Select Imaging Type: All// GENERAL RADIOLOGY

Another one (Select/De-Select): <RET>

This report requires a 132 column output device. DEVICE: (Printer Name or "Q")

Discussion

This report now allows printing by selected Imaging Types.

Enter the name of a printer that prints 132 columns. If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

Inactive CPT	Radiology/Nuclear Medicine Proced	dures (Short)		SEP BI-	30,1	1996 10):54 PAGE	1
CODE	PROCEDURE	AMIS CODE	MULT	LAT	CT	TYPE	MODALITY	
70350	CEPHALOGRAM ORTHODONTIC	1 -SKULL, INC. SINUS, MASTOID, JAW, ET	1			RAD		
70355	PANOREX	1 -SKULL, INC. SINUS, MASTOID, JAW, ET	1			RAD		

-

¹ Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

1 List of Procedures with Contrast

[RA PROCMEDIA]

This option allows the supervisor to generate a list of procedures that are associated with contrast media or medium. The option allows the user to select procedures by imaging type and whether the procedures are active or inactive.

The format of the report includes: Procedure Name, CPT Code, CPT Modifier(s), Procedure Type, Imaging Type, and Inactivation Date (if applicable).

1. From the Menu below, select the Supervisor Menu by typing 'Supervisor Menu' next to Select Rad/Nuc Med Total System Menu Option:

```
Exam Entry/Edit Menu ...
Films Reporting Menu ...
Management Reports Menu ...
Outside Films Registry Menu ...
Patient Profile Menu ...
Radiology/Nuclear Med Order Entry Menu ...
Supervisor Menu ...
Switch Locations
Update Patient Record
User Utility Menu ...
Select Rad/Nuc Med Total System Menu Option: Supervisor Menu
```

2. From the Supervisor Menu, select the Maintenance Files Print Menu by typing 'Maintenance Files Print Menu' next to Select Supervisor Menu Option:

_

¹ Patch RA*5*45 December 2004 New Procedure File Listing: List of Procedures with Contrast.

```
Select Rad/Nuc Med Total System Menu Option: Supervisor Menu
   HL7
          Radiology HL7 Menu ...
          Access Uncorrected Reports
          Delete a Report
          Delete Printed Batches By Date
          Exam Deletion
          Inquire to File Entries
          List Exams with Inactive/Invalid Statuses
          Maintenance Files Print Menu .
Mass Override Exam Status
          Override a Single Exam's Status to 'complete'
          Print File Entries
          Rad/Nuc Med Personnel Menu ...
          Search File Entries
          Switch Locations
          System Definition Menu ..
          Unverify a Report for Amendment
          Update Exam Status
          Utility Files Maintenance Menu ...
Select Supervisor Menu Option: Maintenance Files Print Menu
```

3. From the Maintenance Files Print Menu, select Procedure File Listings by typing 'Procedure File Listings' next to Select Maintenance Files Print Menu Option:

```
Select Supervisor Menu Option: Maintenance Files Print Menu

Complication Type List
Diagnostic Code List
Examination Status List
Film Sizes List
Label/Header/Footer Format List
Major AMIS Code List
Modifier List
Nuclear Medicine List Menu ...
Procedure File Listings ...
Report Distribution Lists
Sharing Agreement/Contract List
Standard Reports Print
VistaRad Category Print

Select Maintenance Files Print Menu Option: Procedure File Listings
```

4. From the Procedure file Listings Menu, select List of Procedures with Contrast by typing 'List of Procedures with Contrast' next to Select Procedure File Listings Option:

```
HIST Display Rad/NM Procedure Contrast Media History
Active Procedure List (Long)
Active Procedure List (Short)
Alpha Listing of Active Procedures
Barcoded Procedure List
Inactive Procedure List (Long)
Invalid CPT/Stop Code List
List of Inactive Procedures (Short)
List of Procedures with Contrast
Parent Procedure List
Procedure Message List
Series of Procedures List
Series of Procedures List
```

5. The Select Imaging Type Prompt displays. Enter a specific, Type of Imaging (see Imaging Types below) (for example, CT SCAN) or type 'ALL' to select all available imaging types.

Imaging Types

ANGIO/NEURO/INTERVENTIONAL CARDIOLOGY STUDIES (NUC MED) CT SCAN GENERAL RADIOLOGY MAGNETIC RESONANCE IMAGING MAMMOGRAPHY NUCLEAR MEDICINE ULTRASOUND VASCULAR LAB

```
Select Imaging Type: All// ALL

Another one (Select/De-Select):

Select one of the following:

A Active
I Inactive
B Both

Select Procedure Status: A// Both
DEVICE: HOME// TELNET
```

- 6. At the prompt: 'Another one (Select/De-Select),' enter another procedure or <RET> to continue.
- 7. At the Select Procedure Status prompt, enter the desired status.

8. At the **DEVICE: HOME**// prompt, enter the device on which to print the data or **<RET>** for the default (screen).

Rad/Nuc Med Procedures with Contrast Media/Medium

Date: Dec 17, 2004 Page 1

CPT Procedure

CPT Modifiers Contrast Media
Procedure Type Imaging Type Inactivation Date

74010 ABDOMEN 2 VIEWS

RESIDENT/TEACHING PHYS SERV Barium

'OPT OUT' PHYS/PRACT EMERG OR URGENT SERVICE Cholecystographic Ionic Iodinated Non-ionic Iodinated

Gadolinium Gastrografin

unspecified contrast media

DETAILED GENERAL RADIOLOGY

Press RETURN to continue...

Parent Procedure List

[RA PARENT PROCEDURE LIST]

The Parent Procedure List includes only parent type procedures and all fields applicable to parent procedures that are included on the Active Procedure long listing. It shows whether they are printsets (i.e., Single Report = Yes), whether a health summary is associated with it for printing with the request, whether Radiology/Nuclear Medicine approval is required to order it, its imaging type, all descendants and their CPT codes, whether the descendents use contrast media, and the educational description.

¹If you assigned one or more modalities to a procedure and the procedure is a descendant of the parent, then the abbreviated modality displays to the right of the descendant.

²Note: Contrast media associations will display, when appropriate.

```
Parent Procedure List

Select Imaging Type: All// <RET>

Another one (Select/De-Select): <RET>

Select one of the following:

A Active
I Inactive
B Both

Select Procedure Status: A Active

DEVICE: HOME// (Printer Name)
```

1

¹ Patch RA*5*3

² Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

MAY 29,1997 08:02 Active Parent Procedure List PAGE 1

Parent Name

C&P SPINAL SURVEY Single Report : YES Hlth Summary w Req. : Rad Approval of Req.: NO

Type of Imaging : GENERAL RADIOLOGY

Descendants : (72040) SPINE CERVICAL MIN 2 VIEWS (72070) SPINE THORACIC 2 VIEWS

(72100) SPINE LUMBOSACRAL MIN 2 VIEWS

Procedure Message :

Educational Desc : THIS STUDY IS USED FOR C&P WHEN A COMPLETE SPINE

IS NEEDED FOR RATINGS

BONE SCAN PARENT Single Report : YES Hlth Summary w Req. : Rad Approval of Req.: NO

Type of Imaging : NUCLEAR MEDICINE
Descendants : (78306) BONE IMAGING, WHOLE BODY Descendants

(78990) PROVISION OF DIAGNOSTIC RADIONUCLIDES

Procedure Message :

Educational Desc : Patient should drink 32 ounces of fluid after

Radionuclide injection and before imaging. Patient

is imaged 2 to 3 hours after the radionuclide

injection.

Procedure Message List [RA PROCMSGPRINT]

This option prints the contents of the Rad/Nuc Med Procedure Message file #71.4. If there is a large number of entries in this file, consider queuing the output.

Prompt/User Response

Discussion

Procedure Message List

DEVICE: (Printer Name or "Q")

Enter the name of a printer or press the <RET> key to bring the output to the screen.

If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

RAD/NUC MED PROCEDURE MESSAGE LIST TEXT

SEP 30,1994 13:42

PAGE 1

TEXT

_

If the patient has not eaten for more than 36 hours or is on TPN, please inform Nuclear Medicine Service so CCK can be administered.

If the patient is on morphine, you cannot assess for CBD obstruction since it will cause spasm of the Sphincter of Oddi.

Make sure the patient has not eaten for 3 hours prior to the test, to avoid a false positive result!

Patient must be NPO after midnight the night before the test!

Patient must have RUQ pain, fever, elevated WBC

Patient must not have eaten within 3 hours of this test!

This test is useful in assessing for ACUTE cholecystitis only, not chronic cholecystitis! Try US for chronic cholecystitis.

Series of Procedures List [RA PROCSERIES]

This list is the same as the Active Procedure List (Long) except it contains only Series procedures.

¹If you assigned one or more modalities to a procedure, then the abbreviated and long name of each modality will print as the last item for the procedure.

²Note: Contrast media associations will display, when appropriate.

Prompt/User Response

Discussion

Series of Procedures List

Select Imaging Type: All// GENERAL RADIOLOGY

Another one (Select/De-Select): <RET>

This report requires a 132 column output device.

DEVICE: (Printer Name or "Q")

This report now allows printing by selected Imaging Types.

Enter the name of a printer that prints 132 columns. If you enter "Q" instead of a printer name, you will also see prompts for a device and a time to print.

```
SEP 30,1994 10:54 PAGE 1
Radiology/Nuclear Medicine Procedures (Series Only)
      PROCEDURE / <sup>3</sup>CPT MODIFIERS
                                                                          MULTIPLIER BILATERAL HEAD/BODY
                                                 AMIS CODE
76062 BONE SURV COMP (INCL APPENDIC 1 -SKULL, INC. SINUS, MASTOID, JAW, ET 1
      Type of Procedure
      Required Flash Card Printer:
      Required Flash Card Format :
      Staff Review of Reports Req:
      Rad Approval of Request Req:
      Type of Imaging : GENERAL RADIOLOGY
Descendents :

Proceeding Manager
     Procedure Message
76061 BONE SURV LMTD (E.G. METASTATI 1 -SKULL, INC. SINUS, MASTOID, JAW, ET 1
      Type of Procedure
                                      SERIES
      Required Flash Card Printer:
      Required Flash Card Format :
      Staff Review of Reports Req:
      Rad Approval of Request Req:
      Type of Imaging :
                                      GENERAL RADIOLOGY
      Descendents
      Procedure Message
```

¹ Patch RA*5*3

² Patch RA*5*45 December 2004 Inserted Note about contrast media associations.

³ Patch RA*5*19 May 2000 "/ CPT MODIFIERS" added to report header.

Report Distribution Lists [RA DISTP]

This option displays data in the Report Distribution Queue file #74.3.

For more information on report distribution, refer to the Index and the chapter on Reports Distribution Edit.

Prompt/User Response

Discussion

Report Distribution Lists

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

Report Distribution Queue List	FEB 10,1997 12:51 PAGE 1 CATEGORY OF INACTIVATION
NAME	REPORTS DATE
TOP OF PAGE MESSAGE	
CLINIC REPORTS	OUTPATIENT
***** Clinic Distribution *****	
FILE ROOM	INPATIENT
***** Radiology File Room Distr	ibution ****
VEDTOIL DECORDS	ALL DEDODES
MEDICAL RECORDS	ALL REPORTS
***** MAS Medical Records Room	Distribution ****
MAS Medical Records Room	=======================================
OTHER THAN WARD OR CLINIC	NON WARD OR CLINIC
OTHER THAN WARD OR CHINIC	NON WAILD OIL CHINIC
**** Non Ward or Clinic Report:	s ****
_	
REQUESTING PHYSICIAN	ALL REPORTS
~	
WARD REPORTS	INPATIENT
**** Ward Distribution ****	

Sharing Agreement/Contract List [RA SHARINGP]

This option prints data from the Contract/Sharing Agreement file #34.

For more information on sharing agreements and contracts, refer to the Index and the chapter on Sharing Agreement/Contract Entry/Edit.

Prompt/User Response

Discussion

Sharing Agreement/Contract List

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

	OCT 3,1996	09:40	PAGE 1
TYPE OF	TNACTIVATION		
_			
AGREEMENT	DATE		
$C \cap V \cap D \setminus C \cap D$			
SHARING			
SHARING			
	TYPE OF AGREEMENT CONTRACT SHARING SHARING	TYPE OF INACTIVATION AGREEMENT DATE CONTRACT SHARING	AGREEMENT DATE CONTRACT SHARING

Standard Reports Print [RA STANDPRINT]

This option prints data from the Standard Reports file #74.1.

For more information on standard reports, refer to the Index and the chapter on Standard Reports Entry/Edit.

Prompt/User Response

Discussion

Standard Reports Print

DEVICE: (Printer Name)

Enter the name of a printer or press the <RET> key to display the output on your screen.

Here is an example of what the report contains:

Standard Reports List OCT 3,1996 09:40 PAGE 1
REPORT
NUMBER STANDARD REPORT

1 NORMAL CHEST - PA/LAT
Report Text:
PA and Lateral views of the chest were obtained. The bones and soft tissues are unremarkable. The cardiac silhouette is normal in size. The lung fields are clear bilaterally.
Impression:
Normal Chest X-ray.

¹VistaRad Category Print [RA VISTARAD CATEGORY P]

This option may be displayed under the 'Maintenance Files Print Menu' ONLY IF the VistaRad component of the Imaging package is being used; otherwise, this option will not be displayed.

This option allows the user to list only the VISTARAD CATEGORY field that is in the EXAMINATION STATUS file (#72), without having to list all the other fields that are seen in the 'Examination Status List' option.

In this option, you will be asked to select one or more Imaging types and then a device. The current VISTARAD CATEGORY value for all exam statuses of the selected Imaging types will be displayed/printed.

Prompt/User Response

```
VistaRad Category Print
Select Imaging Type: GENERAL RADIOLOGY
Another one (Select/De-Select): CT SCAN
Another one (Select/De-Select): <RET>
Select Device: HOME// <RET> TELNET Right Margin: 80// <RET>
                         VistaRad Categories
                                                            Page: 1
                  Run Date: Jul 22, 2002 1:56:26 pm
       Type Of Imaging: CT SCAN
                     Order VistaRad Category
Status
CANCELLED
WAITING FOR EXAM
                         1 Waiting for Exam
                         2 Examined
EXAMINED
                         3 Transcribed
TRANSCRIBED
COMPLETE
        Type Of Imaging: GENERAL RADIOLOGY
          Order VistaRad Category
Status
CANCELLED
                          Ω
WAITING FOR EXAM
                          1
                   CALLED FOR EXAM
EXAMINED
                          3
TRANSCRIBED
Enter RETURN to continue or '^' to exit: <RET>
                       VistaRad Categories
                                                            Page: 2
                   Run Date: Jul 22, 2002 1:56:26 pm
```

_

¹ Patch RA*5*31 September 2002 Added new option.

Maintenance Files Print Menu

COMPLETE 9
Enter RETURN to continue or '^' to exit:

IX. Radiology HL7 Menu

Supervisor Menu ...
¹Radiology HL7 Menu ...

Rad/Nuc Med HL7 voice Reporting Errors Resend Radiology HL7 Message

Overview:

This menu is used to manage HL7 messages created and received by the Radiology/Nuclear Medicine V. 5.0 package.

¹ Patch RA*5*25 Added Radiology HL7 Menu.

Rad/Nuc Med HL7 Voice Reporting Errors¹ [RA HL7 VOICE REPORTING ERRORS]

This option displays HL7 messages received from Voice Recognition Reporting Tools (MedSpeak, PowerScribe, TalkStation, etc.) which have been rejected by VistA Radiology. These messages are stored in the HL7 Message Exceptions file #79.3.

Further options are available to allow the list to be printed and individual messages can be resent and/or deleted.

Prompt/User Response	Discussion
Rad/Nuc Med HL7 Voice Reporting Exception List	This report allows printing by selected Voice Reporting Application, Date and Radiology User.
Select Voice Reporting Application: ALL// RA-CLIENT-TCP	Enter the name of a Voice
Another one (Select/De-Select): RA-TALKLINK-TCP	Reporting Application (File #771) or press the <ret></ret>
Another one (Select/De-Select): <ret></ret>	key to select all applications.
Exception starting date/time: JAN1 (JAN 01, 1999)	Enter a date to begin reporting from.
Exception ending date/time: NOW// <ret> (JUN 11, 1999@</ret>	12:18:26)
	Enter the most recent
	date/time to report on or press the <ret> key to</ret>
	accept the default of now.

_

¹ Patch RA*5*12

Select Radiology User: ALL// ${\tt TAY}{\tt LOR}$, CAMERON

Another one (Select/De-Select): <RET>

Individual users can be selected for reporting. (Entries associated with unknown users will always be displayed). Alternatively, press the <RET> key to display all users.

HL7	Voice Rep	oorting Errors	Jun 11, 1999 12	2:18:38	Page:	1 of 1
	Date	Time	Patient	Case#	User	
	HL7 Voic	ce Reporting A	pplication: MEDSPEA	AK		
1.	•		NAMEB, PATIENT Pointer (field #7		USER, ONE	
2.	•	.999 14:53:16 Missing Exam D	Not known ate and/or Case Nur		Not Known	
3.	•	.999 15:08:38 Invalid Provid	NAMEA, PATIENT er Name	1	Not Known	
	HL7 Void	ce Reporting A	pplication: TALKSTA	ATION		
4.	•	.999 13:45:28 Case Number In	NAMEC, PATIENT valid	2022	USER, ONE	
En	ter ?? for	more actions				
RS DE	Resend HI Delete HI	17 Message	FS First Screen NS Next Screen LS Last Screen //	DN Next	HL7 Exception	

June 2019

Prompt/User Response

Select Action: Next Screen// <RET>

Select Action: Next Screen// LS

Select Action: Quit// FS

Select Action: Next Screen// <↓>

Select Action: Next Screen// <1>

Select Action: Next Screen// PL

Select Device: (enter a printer)

Do you want your output QUEUED? NO// Y (YES)

Requested Start Time: NOW// <RET>
 Request Queued, Task #: 81303

Discussion

List Manager will display the next page of HL7 Exceptions in the center screen.

List Manager will jump to the last screen of entries in the file.

List Manager will re-display the first screen of entries shown above.

Scroll down through each entry individually.

Scroll back up through each entry individually.

The entire list will be sent to a printer. (Each Sending Application will start on a separate page).

Print immediately or request the job to start later.

______ HL7 Voice Reporting Errors (MEDSPEAK - RADIOLOGY/NUCLEAR MEDICINE) Printed: JUN 11, 1999 14:03

Exception Date: MAR 12,1999 at 11:17:20 User: USER,ONE Patient Name: NAMEB, PATIENT Case: 1111 Reason Rejected: Event Protocol Pointer (field #773,8) missing

Exception Date: MAR 23,1999 at 14:53:16 User: Not Known Patient Name: Not known Case: ?????

Reason Rejected: Missing Exam Date and/or Case Number

Exception Date: MAR 23,1999 at 15:08:38
Patient Name: NAMEA, PATIENT User: Not Known

Case: 1

Reason Rejected: Invalid Provider Name

HL7 Voice Reporting Errors Page: 2 (TALKSTATION - RADIOLOGY/NUCLEAR MEDICINE) Printed: JUN 11, 1999 14:03 ______

Exception Date: MAR 16,1999 at 13:45:28 User: USER,ONE Patient Name: NAME, PATIENT1
Reason Rejected: Case Number Invalid Case: 2022

Exception Date: MAR 19,1999 at 15:12:25 User: USER,ONE Patient Name: NAME,PATIENT2 Case: 2098
Reason Rejected: Missing Exam Date

Exception Date: MAR 23,1999 at 14:38:57 User: USER,ONE Patient Name: NAME,PATIENT3 Case: 1111 Patient Name: NAME, PATIENT3
Reason Rejected: Invalid Procedure

Exception Date: MAR 23,1999 at 15:07:53 User: USER,ONE

Patient Name: NAME, PATIENT4 Case: 22

Reason Rejected: Invalid Diagnostic Code

Exception Date: APR 5,1999 at 15:37:09 User: USER, ONE Case: 278

NAME, PATIENT3 Patient Name:

Reason Rejected: Report already on file

User: USER, ONE

Exception Date: APR 6,1999 at 14:39:14
Patient Name: NAME, PATIENT3 Case: 278

Reason Rejected: Report already on file

Prompt/User Response

Select Action: Next Screen// RS
Select HL7 Exception (1-34) :9
 Re-sending Message #3749...

Select Action: Next Screen// DE
Select HL7 Exception (1-34) :10
 Deleting Exception...

Discussion

To re-send a message from the list, Select an entry from the list. Once resent it cannot be re-selected. A message similar to "Message Re-submitted on JUN 11, 1999@14:15:51" is displayed instead of the rejection reason.

To permanently delete an exception report entry, select an entry from the list.

Once deleted it cannot be re-selected. A message similar to "Reported Exception Deleted on JUN 11, 1999@14:22:40" is displayed instead of the rejection reason.

Resend Radiology HL7 Message

[RA HL7 MESSAGE RESEND]

There are many reasons why an HL7 message may not reach its attempted destination. This option allows you to re-send HL7 messages for such cases. The option will prompt for a Radiology case and will then build a new HL7 message in the hope that it will reach the subscribing applications this time.

Prompt/User Response

Discussion

Enter Case Number: 120400-111

Enter an active case number, or a completed case as "MMDDYY-999", or a patient's name, or a patient's 9 digit SSN, or the first character of the patient's last name and the last four digits of the patient's SSN.

Choice	Case No.	Procedure		Name	Pt ID
1	120400-111	ABDOMEN FOR	R FETAL AGE 1 V	PATIENT, TEST	9999
===					=======
Re-Broad	cast the 'EXA	AM REGISTERED'	HL7 message for	this case???	YES// y YES

A confirmation prompt is displayed.
Enter "Y" to Continue or "N" to cancel.
The prompt will either read:
EXAM REGISTERED,
EXAM CANCELLED, OT
REPORT VERIFIED depending on the
current exam and report status.

Send ADT HL7 message to PACS

[RA ADT]

This option forces an ADT-A08 patient demographic update message to be sent from VistA to PACS. If the patient selected had a patient identifier (SSN) update, an ADT-A47 message will also be sent. Note that this option requires patch MAG*3*183 be installed and configured. Please see the MAG*3*183 patch description for details on setting up ADT messages to PACS.

```
Select Patient(s): TEST, PATIENT TRES III 5-2-53 100001113
NO

NSC VETERAN

Another one (Select/De-Select):
Would you like to task this Job? Yes// n (No)
Demographic updates sent for 1 patients
```

Prompt/User Response

Discussion

Select Patient(s):

Enter a patient or 'ALL' to send demographic updates for all radiology

patients.

Another one (Select/De-Select):

Enter additional patients

Would you like to task this Job? Yes// n

(No)

Enter Yes or No

Select Patient(s): TEST, PATIENT TRES III 5-2-53 100001113

NO

NSC VETERAN

Another one (Select/De-Select):

Would you like to task this Job? Yes// (Yes)

Requested Start Time: NOW// (OCT 05, 2017@09:13:26)

Task# 148673

Prompt/User Response

Discussion

Would you like to task this Job? Yes// (Yes)

Requested Start Time: NOW//

Task# 148673

Enter Yes or No

If you choose to task it, you will be

prompted for a start time

Confirmation TASK number displayed

X. Worksheets

Division Parameters

Division:			
Parameter	Parameter Set Equal to:	YES	NO
Division-wide Order Entr	y Parameters:		===
ASK IMAGING LOCATION	ON:		
TRACK REQUEST STAT	US CHANGES:		
CLINICAL HISTORY ME	ESSAGE:		
Exam Entry/Edit Parame	ters:		
DETAILED PROCEDURE	E REQUIRED:		
ASK 'CAMERA/EQUIP/R	M':		
AUTO USER CODE FILI	NG:		
TRACK EXAM STATUS	CHANGES:		
ASK EXAM STATUS TIN	ΛE:		
TIME LIMIT FOR FUTUR	RE EXAMS:		
Films Reporting Paramet	ers:		
ALLOW STANDARD RE	PORTS:		
ALLOW BATCHING OF	REPORTS:		
ALLOW COPYING OF R	EPORTS:		
IMPRESSION REQUIRE	O ON REPORTS:		
ALLOW VERIFYING BY	RESIDENTS:		

DIVISION PARAMETERS (CONTINUED)

Division:			
Parameter	Parameter Set Equal to:	YES	NO
ALLOW RPTS ON CAN	CELLED CASES?:		
WARNING ON RPTS NO	OT YET VERIF?:		
AUTO E-MAIL TO REQ	. PHYS?:		
¹ ALLOW E-SIG ON CO	ΓS HL7 RPTS		
Miscellaneous Division I	Parameters:		
PRINT FLASH CARD F	OR EACH EXAM:		
PRINT JACKET LBLS V	V/EACH VISIT:		
CONTRAST REACTION	NMESSAGE:		
RPHARM DOSE WARN	ING MESSAGE:		
Imaging Locations Asso	ciated with this Division:		

¹ Patch RA*5*12

Camera/Equip/Room

NAME:	DESCRIPTION:

Location Parameters

Imaging Location:		_								
Imaging Type:			Pr	inci	pal	Clir	nic:			
Parameter	Pa	aran	nete	r Se	et Ec	qual	to			
Flash Card Parameters										
HOW MANY FLASH CARDS PER VISIT: 0	1	2	3	4	5	6	7	8	9	10
DEFAULT FLASH CARD FORMAT:										
Jacket Label Parameters										
HOW MANY JACKET LBLS PER VISIT: 0	1	2	3	4	5	6	7	8	9	10
DEFAULT JACKET LABEL FORMAT:										
Exam Label Parameters										
HOW MANY EXAM LABELS PER VISIT: 0	1	2	3	4	5	6	7	8	9	
DEFAULT EXAM LABEL FORMAT:										
Report Parameters										
DEFAULT REPORT HEADER FORMAT:										
DEFAULT REPORT FOOTER FORMAT:										
REPORT LEFT MARGIN:										
REPORT RIGHT MARGIN:									_	
PRINT DX CODES IN REPORT?:			Y	ES:		_	1	NO:		
VOICE DICTATION AUTO-PRINT	YF	ES:			N	IO:				

LOCATION PARAMETERS (CONTINUED)

Imaging Location:	
Parameter	Parameter Set Equal to
Cameras/ Equipment/Rooms Used	by this Imaging Location:
¹ Default CPT Modifiers used by th	
Recipients of the 'Stat' Alert for the	is Location
STAT REQUEST ALERT RECIPIE	NTS:

¹ Patch RA*5*19

LOCATION PARAMETERS (CONTINUED)

Imaging Location:		
Parameter ===================================	Parameter Set Equal to	==
URGENT REQUEST ALERTS?:	YES: NO:	
ALLOW 'RELEASED/NOT VERIFIED':	YES: NO:	
Credit and DSS Parameters		
CREDIT METHOD:		
DSS ID:		

Classification

NAME	CLASSIFICATION (Resident/Staff/Clerk/Technologist)	LOCATION ACCESS
		======

Blank Examination Status Worksheet

Field: Statuses:

rieia:	Statuses:					
		Waiting	Called			
	Can-	for	for	Exam-	Tran-	Com-
	celled	Exam	Exam	ined	scribed	plete
Order in sequence of status progression						
Should this Status appear in Status Tracking						
User Key needed to move an exam to this status						
Default next status for exam						
Can an exam be cancelled while in this status						
Generate exam alert for requesting physician						
Generate Examined HL7 Message						
VISTARAD Category						
Data Requirements for this Status						
Technologist Required						
Resident Or Staff Required						
Detailed Procedure Required						
Film Entry Required						
Diagnostic Code Required						
Camera/Equip/RM Required						
Procedure Modifiers Required						
CPT Modifiers Required						
Radiopharms & Dosages required						
Activity Drawn required (Radiopharm)						
Drawn Dt/Time and Person required (Radiopharm)						
Admin Dt/Time/Person required (Radiopharm)						
Route/Site required (Radiopharm Administered)						
Lot No. required (Radiopharm)						
Volume/Form required (Radiopharm)						
Report Entered Required						
Impression Required						
Verified Report Required						

To Move to this Status	 	 	
Ask for Technologist			
Ask for Interpreting Physician			
Ask for Procedure			
Ask for Film Data			
Ask for Diagnostic Code			
Ask for Camera/Equip/RM			
Ask for Procedure Modifiers			
Ask for CPT Modifiers			
Ask for User Code			
Ask Medications & Dosages			
Ask Medication Admin Dt/Time & Person			
Ask Radiopharmaceuticals and Dosages			
Ask Activity Drawn (Radiopharmaceutical)			
Ask Drawn Dt/Time & Person (Radiopharm)			
Ask Admin Dt/Time & Person (Radiopharm)			
Ask Route/Site of Admin (Radiopharm)			
Ask Lot No. (Radiopharm)			
Ask Volume/Form (Radiopharm)			

Exams in this Status should appear on	 	 	
these Reports			
Clinic Report			
PTF Bedsection Report			
Service Report (Workload)			
Sharing/Contract Report (Workload)			
Ward Report (Workload)			
Film Usage Report (Workload)			
Technologist Report (Workload)			
AMIS Report (Workload)			
Detailed Procedure Report (Workload)			
Camera/Equip/RM Report			
Physician Report			
Resident Report			
Staff Report			
Delinquent Status Report			
Radiopharmaceutical Dosage Ticket Print	 	 	
Control			
Print Dosage Ticket when exam reaches this status			

Examination Statuses

Below is the set of exam statuses as released for each Imaging Type. **Note: They are populated, but not activated.** To activate, assign an Order number, a Default Next Status, and indicate that it should Appear on Status Tracking.

Field: Statuses:

rieia:	Statuses:					
		Waiting	Called			
	Can-	for	for Exam	Exam-	Tran-	Com-
	celled	Exam		ined	scribed	plete
Order in sequence of status progression	0	1				9
Should this Status appear in Status Tracking		Y				
User Key needed to move an exam to this status						
Default next status for exam						
Can an exam be cancelled while in this status	Y	Y	Y			
Generate exam alert for requesting physician						
Generate Examined HL7 Message						
Data Requirements for this Status						
Technologist Required			Y	Y	Y	Y
Resident Or Staff Required					Y	Y
Detailed Procedure Required				Y	Y	Y
Film Entry Required				Y	Y	Y
Diagnostic Code Required					Y	Y
Camera/Equip/RM Required					Y	Y
Radiopharms & Dosages required						
Activity Drawn required (Radiopharm)						
Drawn Dt/Time and Person required (Radiopharm)						
Admin Dt/Time/Person required (Radiopharm)						
Route/Site required (Radiopharm administered)						
Lot No. required (Radiopharm)						
Volume/Form required (Radiopharm)						
Report Entered required					Y	Y
Impression required						Y
Verified Report required						Y
To Move to this Status						
Ask for Technologist			Y	Y		
Ask for Interpreting Physician					Y	
Ask for Procedure				Y		
Ask for Film Data				Y		
Ask for Diagnostic Code					Y	
Ask for Camera/Equip/RM				Y		
Ask for User Code						
Ask Medications & Dosages						
Ask Medication Admin Dt/Time & Person						

Ask Radiopharmaceuticals and Dosages			Y		
Ask Activity Drawn (Radiopharmaceutical)					
Ask Drawn Dt/Time & Person (Radiopharm)					
Ask Admin Dt/Time & Person (Radiopharm)					
Ask Route/Site of Admin (Radiopharm)					
Ask Lot No. (Radiopharm)					
Ask Volume/Form (Radiopharm)					
Exams in this Status should appear on	 				
these Reports					
Clinic Report	Y	Y	Y	Y	Y
PTF Bedsection Report	Y	Y	Y	Y	Y
Service Report (Workload)	Y	Y	Y	Y	Y
Sharing/Contract Report (Workload)	Y	Y	Y	Y	Y
Ward Report (Workload)	Y	Y	Y	Y	Y
Film Usage Report (Workload)	Y	Y	Y	Y	Y
Technologist Report (Workload)	Y	Y	Y	Y	Y
AMIS Report (Workload)	Y	Y	Y	Y	Y
Detailed Procedure Report (Workload)	Y	Y	Y	Y	Y
Camera/Equip/Rm Report	Y	Y	Y	Y	Y
Physician Report	Y	Y	Y	Y	Y
Resident Report	Y	Y	Y	Y	Y
Staff Report	Y	Y	Y	Y	Y
Delinquent Status Report	Y	Y	Y	Y	
Radiopharmaceutical Dosage Ticket Print	 				
Control					
Print Dosage Ticket when exam reaches this status					

Procedure Setup

(This worksheet addresses **only the new fields**, except Single Report which is only seen when editing Parent Procedures. Remember, the radiopharmaceutical prompts are only used with Nuclear Medicine/Cardiology Studies Imaging Type procedures.)

Procedure:	
PROMPT FOR MEDS: YES NO	
DEFAULT MEDICATION:	DOSE:
SUPPRESS RADIOPHARM PROMPT: ASK (O)	NEVER ASK (1)
PROMPT FOR RADIOPHARM RX: YES NO _	
*DEFAULT RADIOPHARMACEUTICAL:	
LOW DOSE: HIGH DOSE: ROUTE: SITE:	· · · · · · · · · · · · · · · · · · ·
EDUCATIONAL DESCRIPTION:	
DISPLAY EDUCATIONAL DESCRIPTION: YES	NO

^{*} Default Radiopharmaceutical is the only field automatically "stuffed" (added to the record) in the case data. Low, High, and Usual doses are displayed for information during exam edits, and route, site, and form will appear as defaults during exam edits.

Radiopharmaceuticals

X201, or DX202. e following this v	See an example of Rad vorksheet.	iopharmaceutical entrie	es with synonyms in the	e Dr
		· 		-
		·		-
				-
		· 		-
		· 		-
				-
		· 		-
				-
				
				-
	 - :	· <u></u>	· 	_
				-
		· 		-

Provide a list of radiopharmaceuticals and their synonyms to your Pharmacy ADPAC so that the

Worksheets		

Example Radiopharmaceutical Drug Classifications

The following is an example of the radiopharmaceutical entries in the Drug file from one of our test sites.

	ADIOPHARMS IN DRUG FILE GENERIC NAME
DX200	Tc99m SODIUM PERTECHNETATE Synonym: SODIUM PERTECHNETATE Synonym: TECHNETIUM Synonym: TC99M
DX200	Tc99m MEDRONATE Synonym: MDP Synonym: TC99M
DX200	Tc99m PENTETATE (DTPA) Synonym: DTPA Synonym: TC99M
DX200	TL-201 THALLOUS CHLORIDE Synonym: TL201 Synonym: THALLIUM Synonym: T1
DX200	Tc99m MEBROFENIN Synonym: TC99M Synonym: HIDA Synonym: MEBROFENIN
DX200	Tc99m ALBUMIN COLLOID Synonym: TC99M Synonym: ALBUMIN COLLOID Synonym: AC
DX200	Tc99m SULFUR COLLOID Synonym: TC99M
DX200	Tc99m SESTAMIBI Synonym: TC99M Synonym: SESTAMIBI
DX200	Ga-67 GALLIUM CITRATE Synonym: GA67
DX200	I-123 SODIUM IODIDE Synonym: SODIUM IODIDE Synonym: 1123
DX200	I-131 SODIUM IODIDE Synonym: SODIUM IODIDE Synonym: 1131
DX200	Xe-133 XENON GAS Synonym: XE133 Synonym: XENON
DX200	In-111 PENTETATE (DTPA)

Worksheets

	Synonym: IN111
	Synonym: INDIUM
DX200	In-111 OXINE
	Synonym: IN111
	Synonym: INDIUM
DX200	Co-57 CYANCOBALAMIN
	Synonym: CO57
	Synonym: CO-57
	Synonym: COBALT 57
	Synonym: CYANCOBALAMIN
DX200	Tc99m LABELLED WBC's
	Synonym: TC99M
	Synonym: WBC
DX200	In-111 LABELLED WBC's
	Synonym: WBC
	Synonym: INDIUM
DX200	Tc99m ALBUMIN AGGREGATED
	Synonym: TC99M
	Synonym: MAA
DX200	Tc99m EXAMETAZIME (CERETEC)
	Synonym: TC99M
	Synonym: CERETEC
DX200	Tc99m BICISATE (NEUROLITE)
	Synonym: TC99M
	Synonym: NEUROLITE
DX200	Tc99m LABELLED RBC's
	Synonym: TC99M
	Synonym: RBC
	Synonym: SODIUM PERTECHNETATE
DX200	Tc99m SODIUM PERTECHNETATE/Cold Pyrolite
	Synonym: TC99M
	Synonym: PYP
DX200	In-111 OctreoScan
	Synonym: INDIUM
	Synonym: OCTREOSCAN

Nuclear Medicine Setup

Sites of Administration Syn	nonym(s)
Vendor/Source	
	
	
	<u></u>

NUCLEAR MEDICINE SETUP (CONTINUED)

Routes of Administration

ROUTE:	(P)rompt for Free Text or (V)alid Site:
Synonyms:	Valid Sites:
ROUTE:	(P)rompt for Free Text or (V)alid Site:
Synonyms:	Valid Sites:
ROUTE:	
Synonyms:	Valid Sites:
ROUTE:	(P)rompt for Free Text or (V)alid Site:
Synonyms:	Valid Sites:

XI. Parent/Descendent Exam and Printsets

¹Definition

An exam set or printset contains a Parent procedure and its Detailed or Series descendent procedures. Requesting a Parent will automatically cause each descendent to be presented for registration as separate cases under a single visit date and time. If the parent is defined to be a printset, the collection and printing of all common report related data between the descendents is seen as one entity.

Common Data for Printsets

For a printset, a single report is entered regardless of how many descendent cases are registered. Common data for all cases in a printset includes all fields recorded in the Rad/Nuc Med Report file #74 and a few fields in the Rad/Nuc Med Patient file #70: report text, impression, diagnostic codes, primary and secondary residents and staff, verifier, pre-verifier, clinical history, verification and pre-verification date/time and physician, transcriptionist, no-purge indicator, reported date, report status, and exam date/time. For printsets, case-related fields that are common data (residents, staff, diagnostic codes, etc.) are now editable only through report editing functions, and are no longer editable via case editing options. This is necessary to ensure that there is no opportunity for users to enter different data in these fields for each case in a printset.

Non-Common Data for Printsets

Non-common data still exists for each procedure in a printset. Non-common data includes modifiers, case number, exam status, complications, film, and most other fields that can be edited separately for each case.

Single Report Setup for Printsets

When setting up Parent type procedures through the Procedure Enter/Edit option, you will see a new prompt, "SINGLE REPORT:", that should be set to YES to indicate that reporting for the registered descendents of the parent procedure will be combined into a single report, making it a printset. If the field is left blank, a separate report will be required for each registered descendent and the cases will NOT be treated as a printset. When descendents are registered in the

¹ Patch RA*5*45 March 2004 Parent/descendent procedures that use contrast media

Rad/Nuc Med Patient file #70, the "MEMBER OF SET" field on the case record is set automatically by the system as follows:

- 1 = descendent procedure, separate reports required
- 2 = descendent procedure, single report for all cases under this exam date/time

blank = not a descendent, separate report required

Selection of any member of a printset in report related functions will bring up certain identifying information (case numbers, procedures) for all cases in the printset.

Report printouts and displays are modified for printsets to show all non-common data (such as procedure names and modifiers) on one report. The report print format for printsets has been changed to include each procedure and its modifiers. Residents, staff, and diagnoses are printed once for each printset because they are the same for cases in the printset. The same changes have been made to screen displays of reports, to email messages containing reports, and to the retained erroneous report that is stored on the report record before amending it.

If a report for a printset is deleted, each exam in the printset will show that its report was deleted.

If a user tries to create a report for a cancelled case, and the cancelled case belongs to a printset where the other cases of the printset already have a report, then the cancelled case will refer to the same report as the rest of the printset if:

the Division parameter allows reports on cancelled cases, or

the user owns the RA MGR security key.

If neither of these conditions is true, the cancelled case will NOT have a report associated with it. If a cancelled case is selected via a voice recognition system, and the case is part of a printset, the user will not be able to enter a report. Instead, an error message will be transmitted notifying him to use VistA for reporting on this case.

Special Character Designation for Printsets

All options that do case lookups and displays will now show special characters to the left of the case number column for cases belonging to a printset. **Note: These will only appear on cases that are entered after the installation of V. 5.0.** This also includes the List Reports in a Batch option and the Distribution Queue Menu options: Clinic Distribution List, Ward Distribution List, and Unprinted Reports List. The Register an Exam and Display Patient Demographics options also display special characters for cases in a printset when the list of the last five procedures is displayed. Options such as Status Tracking do NOT display the

special characters because cases of the same printset may be in different statuses and therefore would not appear on the same screen.

In the example below, "+" denotes the first case registered in a set of descendents, and "." appears to the left of the rest of the descendents in the set. Only printsets are displayed this way. So, if the parent procedure's "SINGLE REPORT" field is not set (i.e., if it is left blank), there will be no special characters to the left of the case number. The sequence of case numbers/procedures within a printset are shown in the same order as they were registered regardless of whether the lookup is from an exam edit or report edit/print options.

Sample look-up display:

	Case No.	Procedure	Exam Date
1	122	HAND 3 OR MORE VIEWS	03/23/96
2	+113	ANKLE 2 VIEWS	04/30/96
3	.114	FOOT 2 VIEWS	04/30/96
4	.115	TOE(S) 2 OR MORE VIEWS	04/30/96
5	81	NECK SOFT TISSUE	05/03/96

HL7 Interface

The HL7 interface to the voice recognition systems is changed so that the voice recognition system screen showing active cases for a selected patient will display cases belonging to a printset. For more complete information on the interface, refer to the Rad/Nuc Med Technical Manual.

The PACS/Imaging interface is changed to include multiple OBR and OBX segments for each individual procedure when a single HL7 message broadcasting an imaging report must account for all the cases in a printset.

Registration

During registration, a single Detailed, Series, or Broad procedure request can be replaced by a parent-printset from the same Imaging Type and the predefined descendant(s) can either be registered or discarded for other procedures by doing the following:

- 1. Enter "Pn" at the prompt, where "P" indicates that you want to trigger the parent-printset registration feature and "n" is the request number. The request must NOT be a parent. Only one "n" value is allowed. You will then see a prompt for a parent procedure.
- 2. Enter a parent procedure of the same imaging type as the requested procedure. The parent must be predefined as a printset.
- 3. Then proceed to register the predefined descendant(s) OR, discard (^) its descendant(s) and register any descendants you choose when it asks you for more procedures to add on.

The list of requests displayed will have a "+" in front of procedures that are printset parents. (A "+" will not display in front of non-printset parents.) After the replacement printset parent is registered, all outstanding orders that have matching procedures to any just registered will be displayed as a reminder that these may be duplicates and might have to be cancelled.

See the session capture that shows the use of Pn in the chapter Register Patient for Exams, Radiology/Nuclear Medicine User Manual.

Editing Data

During exam edit, each individual case can be retrieved and edited separately. When editing a case or Diagnostic Code on a case that is a member of a printset, no entering or editing will be allowed on common data. Common data may only be entered and edited through Report Entry/Edit since it must be the same for all cases in the printset. The Report Entry/Edit option now displays all case numbers, exam statuses and procedures that belong to the same printset as the case selected.

When common data is entered through Report Entry/Edit all printset cases' statuses will be updated accordingly. A screen display will warn transcriptionists that all cases are being updated with the same common data. Case editing options only update the status of the individual case being edited; case editing no longer allows editing of report-related fields that are common to all cases in the printset, and, therefore, case editing cannot affect the status of other cases in the printset.

Labels, Headers, Footers

Label Print field output on printed labels and report headers and footers is modified to print a "+" next to data that is non-common to indicate that the data is part of a related set. This was done because multiple values cannot be printed in the restricted space available in report headers and footers, or on labels. Label print fields affected (case number, procedure, last time this exam performed, long case number, technologist, resident, staff radiologist, and exam modifiers) will have a YES in the new "UNIQUE" field #6 on the Label Print Fields file #78.7. This uneditable field is released with data and should NOT need further tailoring by the site.

VA Fileman Prints/Inquiries

The same "+" is used for computed fields on the Rad/Nuc Med Report file #74 when output is generated using VA FileMan Inquire or Print options. An example of this would be using VA FileMan to print File #74's Primary Camera/Equip/Rm computed field. Since the data for this field could be different for each case in a printset, only the data from the first case registered is printed, and a "+" is printed next to it to indicate there could be different Primary Camera/Equip/Rms in other members of the printset.

Any case number in a printset can be used in VA FileMan lookups of the Rad/Nuc Med Report file #74 to retrieve a report for the printset. A new "OTHER CASE#" field #4.5 has been added to this file, and contains all case numbers for a printset report other than the case number already in the "DAY/CASE#" field #.01.

Add Exams to Last visit

The Add Exams to Last Visit option now allows cases to be added to exam sets **as long as no report has been entered yet for any of the cases**. No new order will be created; the added case will point to the same order as the rest of the cases in the exam set. This is true for printsets, as well as exam sets that require separate reports for each case.

CPRS Interface for Printsets

The Rad/Nuc Med OE/RR V. 3.0 (CPRS) interface includes functionality to make it plain to the user that a single report applies to multiple procedures.

Imaging/Multi-Media Interface Effects on Printsets

If your facility is running the Imaging interface, a placeholder ("stub") report will be created to hold the image IDs as soon as the images are collected. From that point on, a lowercase "i" will appear next to the procedure on various lookup screens, but the stub report is otherwise transparent to the user. If you need to delete a printset exam member, you will not be allowed to since a stub report exists. However, you can cancel the case instead.

XII. Case Edits and Status Tracking

Overview

The case edit and status tracking options can vary considerably due to the way the system is set up at your site. By defining Division parameters using Division Parameter Set-up [RA SYSDIV], procedure parameters using the Procedure Enter/Edit [RA PROCEDURE] option, and status parameters using the Examination Status Entry/Edit [RA EXAMSTATUS] option you can control a great deal of what the user sees when doing case edits or updating the status of a record. The flow of prompts you see in the case edit and status tracking sections of the user manual may not reflect what your user sees. This chapter is dedicated to providing you some insight into the system response to variations in setup.

Taking Control

The case edit options **always** follow the same logic; you have little control over what is asked and when it is asked. Case No. Exam Edit and Edit Exam by Patient will prompt for data items that are automatically entered by the system during registration (such as inpatient, ward, bed section, etc.) and data already entered previously. This gives the opportunity to change incorrect data. Complications can also be entered in the case edit options.

To gain control of what and when specific fields appear for editing, you should set up statuses via Examination Status Enter/Edit and use Status Tracking of Exams to enter case information. In other words, always use Status Tracking of Exams for entering case data if you want to control the data flow. Use Case Edit options when you want to edit complications and fields already populated. If your facility does not use Status Tracking, you must make sure Exam Status Parameters don't require any fields that are not asked in the Case Edit options or else you will have problems getting cases to Complete.

Stuffing Data and Default Responses

Stuffed data is data that automatically becomes an entry in the record because certain conditions are met. Depending on the data that is stuffed, it may or may not be editable. For instance, the sign-on user name may be stuffed but not editable, whereas the requested detailed procedure may be stuffed but it is editable at registration. Default responses may be stuffed data, or they may be provided as

the most logical response to a prompt (e.g., Default Response//) and only become part of the record when you strike the <RET> key to accept it. Some data items such as Default Radiopharmaceutical, can be set up through Procedure Enter/Edit to be automatically entered (stuffed) on a patient exam record. Other data items, such as radiopharmaceutical dose, form etc., are not stuffed into the record, but will appear as default responses if defined for a procedure.

Only detailed and series procedures are capable of having stuffed and default responses. Therefore, a parent does not have default responses, but its descendents can have them. A broad procedure cannot have default responses. However, if the requested broad procedure is changed to a detailed procedure at registration, then the user can benefit from the detailed procedure's default answers. If the procedure is changed later on, during case edits or status tracking, the user will need to check the data in the default fields to verify accuracy.

The fields that can have responses stuffed when entering data are Film Size,

Radiopharmaceuticals, and Med Administered. The default data for these fields are entered via the option Procedure Enter/Edit. If the procedure is not a Nuclear Medicine or Cardiology Studies procedure, then radiopharmaceuticals cannot be associated with it.

If a radiopharmaceutical and/or medication was stuffed during registration and then the procedure was later changed, the program will delete the stuffed data and display information messages on the screen.

Radiopharmaceutical Logic

The following radiopharmaceutical logic pertains only to those procedures that use radiopharmaceuticals. There are three parameters that work in conjunction to determine what prompts the user sees:

Whether the user is **asked** for radiopharmaceuticals, When the user is **asked** for radiopharmaceuticals, and If radiopharmaceuticals are **required**.

Suppress Radiopharm Prompt, edited via the option Procedure Enter Edit, determines **whether** or not the user is **asked** for a radiopharmaceutical. It can contain Always (0 or Null (no answer entered)) or Never (1) using the following definitions:

- O Always (same as Null) ask for Radiopharmaceutical if and when exam status parameters specify (this is the default).
- Never ask (suppress the radiopharmaceutical prompts) or require
 Radiopharmaceutical even if exam status parameters specify that a
 Radiopharmaceutical should be entered. If default Radiopharmaceuticals are
 entered for the procedure they will not be automatically entered on the exam
 record by the system when the case is registered.

Ask Radiopharms and Dosages? (while in a specific status), edited via the option Examination Status Entry/Edit, can contain a YES (**when** trying to progress to this status, **ask** it) or NULL (don't ask it).

Radiopharms/Dosages Required?, edited via the option Examination Status Entry/Edit, determines **if** entry of a radiopharmaceutical is **required**. It can contain either a YES (it is required) or NULL (it is not required). This determines whether it is required to progress to the status.

To suppress asking/requiring radiopharmaceuticals on an individual procedure, you must enter NEVER on the procedure field, Suppress Radiopharm Prompt. Status Tracking and Case Edit both skip radiopharmaceutical prompts if NEVER is entered. Status update logic will also skip over the radiopharmaceutical requirement if NEVER is in the field, even if the Exam Status, Radiopharms/Dosages Required?, parameters are set to YES. Also, Registration will **not** stuff default radiopharmaceuticals if NEVER is in that field, even if default radiopharmaceuticals have been entered on the procedure file. The Procedure Enter/Edit option will delete any default radiopharmaceuticals previously entered for this procedure if this field is NEVER.

Medication Logic

Medication logic pertains to all types of procedures. There are two parameters that work in conjunction to determine what prompts the user sees:

Whether the user is **asked** for medications, and When the user is **asked** for medications.

Prompt for Meds (Rad/Nuc Med Procedure file #71, field #5), edited via the option Procedure Enter Edit, determines **whether** or not the user is **asked** for a medication.

Ask Medications & Dosages? (Examination Status file #72, field #.27), edited via the option Examination Status Entry/Edit, can contain a YES (**when** in this status, **ask** it) or NULL (don't ask it).

Medication prompts will only appear if both Ask Medications & Dosages and Prompt for Meds are set to YES. Case Edits and Status Tracking both use Prompt for Meds, but Case Edits will always ask (and re-ask on additional edit sessions), whereas Status Tracking will only present the medication prompt for the appropriate status as set up in File #72. Medications do not affect status updates because they are **never** required.

Radiopharm Prescription & Prescribing Physician/Dose Logic

This logic affects only those procedures that use radiopharmaceuticals. There is only one determining factor:

Whether the user is asked for the prescribing physician and prescribed dose.

Examination Status parameter fields concerning prescriptions for radiopharmaceuticals do not exist. A single field on the Rad/Nuc Med Procedure file #71, Prompt for Radiopharm RX, field #19, is used for determining prompting logic. During the Status Tracking session where the dose administered has been entered, Status Tracking will ask for prescription data only if Field #19 is set to YES. Once prescription data is entered, Status Tracking will no longer prompt for prescription data in subsequent edits; however, if only one prescription field (physician or dose) is entered, Status Tracking will continue to ask both fields at every subsequent edit session. Case edits will always ask for prescription data if Field #19 is set to YES. In both Case Edits and Status Tracking, a warning will continue to display if the dose administered is outside of the high/low range. Prescription data do not affect status updates because it cannot be required.

Population of Imaging Exam Data

Patient Exam Data Field	Automatically Entered	Dx Code Edit	Asked in Case Edits	Asked in Case Edits only if already populated	Asked in Status Tracking if "Ask" param is set to YES	Cannot be asked in Status Trking
Visit:						
Exam Date						X
Type of Imaging	When registered					X
Hospital Division	When registered					X
Imaging Location	When registered					X
Examinations:						
Case Number	When registered		**		**	X
Procedure	A () 1 11 ()		X		X	37
Exam Status	After each edit of case and report.					X
Category of Exam	When registered		X			X
Ward	When registered if			X		X
	inpatient			**		**
Service	When registered if inpatient			X		X
Bedsection	When registered if inpatient			X		X
Principal Clinic	When registered if outpatient		If category = opt	X		X
Contract/Sharing Source	When registered if entered on request			Also asked if category = contract/sharing		X
Research Source	When registered if entered on request			Also asked if category = research		X
¹ Contrast Media Used	When registered procedure is associated with contrast media, this field will default to "YES."		X	Always Asked	Always Asked	
² Contrast Media	When Contrast Media Used is YES, the user can edit the contrast media used for the exam.		Only Asked if 'Contrast Media Used' ="Yes"	Only Asked if 'Contrast Media Used' ="Yes"	Only Asked if 'Contrast Media Used' ="Yes"	
Imaging Order	When registered					X
Primary Interpreting		If NOT a			If NOT a]
Staff		printset			printset	
Secondary Interpreting		If NOT a			If NOT a	
Staff	1	printset			printset	
Primary Interpreting		If NOT a			If NOT a	
Resident Secondary Interpreting		printset If NOT a			printset If NOT a	-
Residents		printset			printset	
Primary Diagnostic		If NOT a			If NOT a	
Code		printset			printset	
Secondary Diagnostic		If NOT a			If NOT a	1
Codes		printset			printset	
Requesting Physician	When registered		X			X
Complication	,		X			X

Patch RA*5*45 March 2004 Parent/descendent procedures that use contrast media
 Patch RA*5*45 March 2004 Parent/descendent procedures that use contrast media

Case Edits and Status Tracking

Complication Text		When a		X
		complication is		
		entered		

Patient Exam Data Field	Automatically Entered	Dx Code Edit	Asked in Case Edits	Asked in Case Edits only if already populated	Asked in Status Tracking if "Ask" param is set to YES	Cannot be asked in Status Trking
Primary			If division		X	
Camera/Equip/Rm			param = YES			
Requested Date	When registered					X
Requesting Location	When registered					X
Visit	When complete					X
Dosage Ticket Printed?	When printed					X
HL7 Examined Msg Sent?	When sent					X
Film Size Data:						
Film Size	When registered if set up on procedure		X		X	
Amount (#films)	When registered if set up on procedure		X		X	
Exam Status Times						
Status Change	After each status				If division	
Date/Time	change if division param = NO				param = YES	
New Status	X					X
Computer User	X					X
Activity Log						
Log Date	If division param = YES					X
Type of Action	X					X
Computer User	X					X
1 Technologist			X			
Comment						
Modifiers	When registered, from request		X			X
Technologist		X	X		X	
Medications						
Med Administered	If procedure defaults set up		If procedure param for meds = YES		If procedure param for meds = YES	
Med Dose			If procedure param for meds = YES		If procedure param for meds = YES	
Date/Time Med Administered				X	If procedure param for meds = YES	
Person Who Administered Med				X	If procedure param for meds = YES	
Clinical History for Exam	When registered, from request					

⁻

¹ Patch RA*5*18 August 2000 Added field for comments by the technologist.

Patient Exam Data Field Radiopharmaceuticals	Automatically Entered	Dx Code Edit	Asked in Case Edits	Asked in Case Edits only if already populated	Asked in Status Tracking if "Ask" param is set to YES	Cannot be asked in Status Trking
Radiopharms	When registered, if		If procedure		If procedure	
	procedure defaults set up		param Suppress Radiopharm = NO		param Suppress Radiopharm = NO	
Fields below are asked only if a radiopharm is entered. Radiopharms and related fields only apply to procedures of Nuclear Medicine or Cardiology Studies.						
Prescribed Dose by MD Override			If procedure Rx param = YES		If not already entered, and if procedure Rx param = YES	
Prescribing Physician			If procedure Rx param = YES		If not already entered, and if procedure Rx param = YES	
Activity Drawn				X	X	
Date/Time Drawn				X	X	
Person Who Measured Dose				X	X	
Dose Administered			X		X	
Date/Time Dose Administered			X		X	
Person Who Administered Dose			X		X	
Witness to Dose Administration			If procedure Rx param = YES		If not already entered, and if procedure Rx param = YES	
Route of Administration				X	X	
Site of Administration				X	X (unless no predefined sites for route, and Free Text site param is OFF)	
Site of Admin Text				X	X	
Lot No.				X	X	
Volume Form				X	X	

XIII. Outpatient Workload Crediting

Imaging Stop Codes File 71.5

This file contains the stop codes that can be used by imaging services. They are determined by MAS, VA Headquarters. You can edit these stop codes via the option Valid Imaging Stop Codes Edit and delete any not used at your site. It provides the list of selections for stop codes assignment to Imaging Locations. Upkeep of this file is entirely up to you, the ADPAC. It is not a pointed-to file (meaning preservation of its data does not affect reports, etc.) therefore entries should be deleted if they don't need to be on the selection list.

Invalid CPT/Stop Code List

The report option Invalid CPT/Stop Code List under the Procedure File Listings menu should be used as a tool to review the validity of procedure CPTs, Imaging Location setup, and several other "sanity checks". Warnings and errors on this report must be addressed in order to properly credit imaging workload at your facility.

This report shows any CPTs that are nationally inactive, missing, or have broken pointers. It prints warnings for Parent type procedures with no Descendents, and for detailed or series procedures with no AMIS code(s). The report will also tell you if there are any Imaging Locations with incorrect DSS IDs (i.e., Stop codes not in the Imaging Stop Codes file), and if an Imaging Location is not pointing to a Hospital Location file #44 entry marked as Occasion of Service.

Valid Imaging Stop Codes Edit

All stop codes are determined by the Medical Administration Office in VA Headquarters which updates the contents of the Clinic Stop file #40.7 every October 1 via a software release. After IRM installs the software that updates that file, you can use the option Inquire to File Entries to get a listing of the entries in File #40.7 and check them against entries in your Imaging Stop Codes file #71.5. If there are any changes, you need to update File #71.5 by deleting the ones that have become obsolete and adding the codes that correspond to the imaging services done at your site via the option Valid Imaging Stop Codes Edit. If you add new Stop Codes, it is your responsibility to first verify that the codes are recognized as valid for Imaging Services by MAS, VA Headquarters.

Location Parameter Edit

Whenever you use the Location Parameter Edit option to edit or set up an Imaging Location, the software automatically checks the "pointed to" entry in the Hospital Location file #44. You should see a short message such as "Location is OK" after editing. If the software detects a problem, it will attempt to correct the problem, possibly by creating a new File #44 entry for the Imaging Location to point to. This will happen if someone has erroneously entered appointment patterns for the File #44 location. If the Stop Codes don't agree, the Stop Code entry in File #44 will be overwritten to match the DSS ID entered by the Rad/Nuc Med ADPAC. If an error cannot be resolved by the software, you will see a message telling you to take some action. If the software can fix the problem, you will see a message telling you what was done.

Options that can Potentially Trigger Crediting

Case No. Exam Edit
Edit Exam by Patient
Report Entry/Edit
On-line Verifying of Reports
Verify Report Only
Status Tracking of Exams
Update Exam Status
Diagnostic Code Entry by Case No.
Resident On-line Pre-Verification
Verify Batch
Report Entry through vendor supplied Voice Recognition Systems

Crediting

- 1. Rad/Nuc Med sends data to the PCE (Patient Care Encounter) package whenever an exam's status progresses to Complete.
- 2. Using the option Mass Override Exam Status or Override a Single Exam Status to Complete does **not** result in any stop code or CPT crediting.
- 3. Crediting is only filed for outpatients.

4. The following data is mandatory for each case:

Detailed Procedure with a valid CPT code
Primary Interpreting Staff or Primary Interpreting Resident
Patient
Exam Date/Time
DSS ID
Requesting Location

- 5. The PCE package sends the Scheduling Package the crediting data that needs to be transmitted to Austin.
- 6. If PCE accepts the Rad/Nuc Med crediting data, a value will be written to the Visit field #27 of the Rad/Nuc Med Patient file #70 on that case. If not, the failure bulletin is sent to recipients.
- 7. If descendent cases are registered for a given exam date/time, each case may be either Complete or Cancelled before data is sent to PCE.

Errors

There is a bulletin (RAD/NUC MED CREDIT FAILURE) that can be used to notify a mail group of credit errors generated specifically through Radiology/Nuclear Medicine. Your IRM should set up a mail group that includes the ADPAC, or whoever will be responsible for researching and resolving crediting problems.

Note: For more information on outpatient crediting, please see the PCE and Scheduling documentation.

XIV. Troubleshooting

On-line Verifying of Reports

The following is a chart using various parameter set-ups in the options Division Parameter Set-up and Classification Enter/Edit to show how On-line Verifying of Reports is affected. If selecting another physician's reports in on-line verifying is not working the way you feel it should, check this chart to see what you need to change.

Division Parameter Set-up	Classifi- cation Enter/Edit	Classification Enter/Edit	Depending on parameters, the physician using on-line verification may or may not be allowed to select another physician's reports.
Allow Verifying by Residents	Allow Verifying of Others	Classification of User logged on	First prompts displayed when On-line Verifying of Reports is selected
Y	Y	Resident	Select Interpreting Physician*
N	N	Resident	Interpreting Residents are not allowed to verify reports.
Y	N	Resident	Select one of the following: 1 Pre-verified 2 Draft 3
N	Y	Resident	Interpreting Residents are not allowed to verify reports.
Y	Y	Staff	Select Interpreting Physician
N	N	Staff	Select one of the following: 1 Pre-verified 2 Draft 3
Y	N	Staff	Select one of the following: 1 Pre-verified 2 Draft 3
N	Y	Staff	Select Interpreting Physician

^{*} The Select Interpreting Physician prompt allows the user to choose another physician and review that physician's reports. If this isn't the first prompt displayed, the user cannot edit another physician's reports.

Physician Title on Reports: The title that prints to the right of the physicians' names on the results report is taken from the New Person File #200 Signature Block Title field #20.3. So, if you find that incorrect titles are printing, have the user or IRM edit this field.

Reports

Workload Reports: Please note that the separation of workload figures by imaging type supported in version 4.5 will take effect only on exams entered after the ADPAC activates the new imaging types. For example, if an exam was entered before the new imaging type was activated, the exam will be counted under its old imaging type. To include only exams entered after the new imaging types were activated, the beginning date entered at the date range prompts when generating the report must be after the date when the new imaging types were activated.

For most workload reports that are sortable/selectable by one-many-all, division and imaging type, the division totals page will only print if there is more than one imaging type in the division. If there is only one imaging type in the division, the imaging type total page is used for the division total.

The following reports calculate workload counts (i.e., exam counts, patient visit counts, and weighted work units) in a similar way:

Functional Area Workload Reports...

Clinic Report PTF Bedsection Report Service Report Sharing Agreement/Contract Report

Ward Report

Special Reports... Camera/Equip/Rm Report Personnel Workload Reports...

Physician Report Resident Report Staff Report Technologist Report Transcription Report

- 1. Data must fall within the date range selected.
- 2. Exam Status file (#72) parameters are checked to see if the case's exam status is one that should cause the exam to appear on the report.
- 3. Exam status must have an Imaging Type the user has selected or defaulted to, and the exam's division must be one the user selected or defaulted to.
- 4. Checks the Category of Exam to identify the patient category (inpatient, outpatient, research, or other). The personnel reports identify the patient as an inpatient or outpatient. The functional reports identify the patient as an inpatient, outpatient, research, or other. The clinic report will ignore a case if the category is inpatient. The ward, bedsection, and service reports will ignore a case if the category is outpatient or other.
- 5. Ignores the case if the procedure does not exist in the Rad/Nuc Med Procedure file (#71). (This would mean there is a broken pointer type of data corruption problem in the record.) Theoretically, this shouldn't happen. If it does, it means that someone completely deleted a procedure.
- 6. Ignores the case if there are no AMIS codes associated with the procedure. (AMIS codes can be entered using the Procedure Enter/Edit option.)
- 7. Ignores the case on the Sharing Agreement/Contract report if the pointed to entry in the Contract/Sharing Agreements file (#34) does not exist. ("Does not exist" means there is a broken pointer type of data corruption, which happens when someone inadvertently deletes an entry from the Contract/Sharing Agreements file.) If no contract/sharing source was entered, the case is reported under other.
- 8. Calculates the weighted work units (WWUs). The AMIS Weight Multiplier field of the Rad/Nuc Med Procedures file contains a number (0-99) to indicate to the various workload report routines how many times to multiply the weighted work units associated with the AMIS code. The Weight for each AMIS code is stored in the Weight field of the Major Rad/Nuc Med AMIS Code file.

Most multipliers will be 1. However, there are some that are greater than 1. For example, a procedure called Upper GI and Small Bowel might have AMIS code 9-Gastrointestinal which has a Weight of 6 work units, and an AMIS Weight Multiplier of 2. Therefore, on the workload reports, the site will get credit for 12 weighted work units (WWUs) each time it is performed. If there are multiple AMIS codes for the procedure, each AMIS Weight Multiplier is multiplied by the AMIS weight, then the results are summed.

Below is a sample of the exam/procedure/AMIS data relationship. Using this sample, the WWUs for the exam would be 12, and the exam would count as 2 exams.

Rad/Nuc Med exam data stored in Rad/Nuc Med Patient file #70:	Procedures stored in Rad/Nuc Med Procedures file #71:	AMIS categories stored in Rad/Nuc Med Major AMIS Codes file #71.1:
Patient: Joe Veteran Procedure:> Exam Modifiers: (none)	Procedure: Upper GI AMIS codes data: AMIS Code(s):> AMIS Weight Multiplier: 2 Bilateral: CT Head or Body:	Desc: Gastrointestinal Weight: 6

In the above example, this would be the calculation for Weighted Work Units: $6 \times 2 = 12$.

The modifiers (operating room, portable, and bilateral) have special meaning and cause increased counts. If the AMIS Weight Multiplier is less than 2, "bilateral" affects the WWUs by multiplying the AMIS Weight Multiplier by 2. The result is multiplied by the AMIS code's "Weight". In the sample above, if the procedure had been bilateral, WWUs would never have been affected. Here's step-by-step how it works:

It finds all the AMIS codes associated with each procedure. If no AMIS weight has been entered, the default is 1.

If bilateral (i.e., the procedure's "Bilateral" field is set to "YES" **or** a bilateral modifier was entered during exam ordering or case editing), the multiplier value (AMIS Weight Multiplier) is doubled.

For each AMIS code on the procedure, it multiplies the weight times the weight multiplier of each code, then sums the products (if more than one AMIS code on the procedure).

- 9. Increments the counts for division, Imaging Type report, procedure etc.
- 10. If more than one AMIS code exists for the procedure, the appropriate exam category count is incremented by one. If only one AMIS code for the procedure, increment the patient count by the weight multiplier.
- 11. Cases with more than one technologist, resident, and staff will be incremented accordingly with a message warning that the total number of exams and weighted work units will be higher in these personnel reports than the other workload reports.
- 12. WWUs are printed on all functional area reports, the Technologist Report, and the Camera/Equip/Rm Report. WWUs do not apply to the Physician, Resident, Staff and Transcription Reports.

Abnormal Exam Report: This report includes only those exams that have a Diagnostic Code whose "PRINT ON ABNORMAL REPORT" field in the Diagnostic Codes file #78.3 is set to "YES".

¹If several cases share one report (printset), then they will be printed together under the same patient and exam date/time.

Delinquent Status Report: There are times when an exam seems like it should be in the Complete status when it is actually in a Transcribed or lower status. The following steps describe how this might happen.

- . A case is registered; the status is Waiting for Exam.
- . The exam is completed quickly and images are sent to the radiologist or nuclear medicine physician before status tracking or case edit is used.
- . The transcriptionist enters the report, but the status either cannot be upgraded because the case has not been edited by the technician or cannot be upgraded to the status immediately before Complete.
- . The report is verified, but the technician still has not done case edit or status tracking, so the exam remains in some status lower than Complete.
- . The technician now uses status tracking (not Case No. Edit or Edit Exam by Patient). Status tracking pushes the exam status to Examined (or whatever the next ordered status is) because it is designed to push the status to the next ordered status rather than the highest possible status.
- . Status Tracking will display all higher status(es) for which this exam now qualifies. However, if the technician does not follow up with additional status tracking edits of the same case, the case stays in Examined or Transcribed (or whatever) and goes into the "black hole" until it shows up on the Delinquent Status report.

Another scenario circumstance that may cause the same finding:

If status tracking cannot account for incomplete cases, the other possibility is that an error occurred on a tasked job that was supposed to update the status. This is likely if there are very few occurrences. Since status update tasks call so many

¹ Patch RA*5*34 Inserted text regarding several cases sharing one printset.

Troubleshooting

other programs, the error might not even show evidence to IRM that Radiology/Nuclear Medicine is responsible. This could be reconstructed based on log subfiles that are populated (if the logs are turned on) together with careful scrutiny of the error log that IRM can access.

Use the Delinquent Status report regularly to catch this infrequent happening.

XV. Frequently Asked Questions

Many of the following questions were originally asked over E-Mail. Whenever you see quotes around a statement (some editing was necessary), it was lifted from actual questions and answers from the sites. We wanted to provide the user's viewpoint along with items of particular importance that we feel need to be discussed as often as possible. The following should be viewed as possible solutions to handle various conditions, but not necessarily the only way or the best solution at every site. When in doubt, consult with your hospital IRM Service.

Imaging Types

Question: How do I deactivate an Imaging Type that we have been using?

Answer: Assigning an Imaging type not previously used to locations and procedures "activates" that Imaging Type. That causes workload reports for that Imaging Type to be separately reported. Imaging Types also affect the enter/edit options and screen displays. When a user signs on to the computer, s/he signs on under a specific location and Imaging Type which affects the data the user can see or edit.

Once an Imaging Type has been "activated" and used by registering cases to the associated location, it takes a great deal of effort and time to deactivate the Imaging Type.

If an Imaging Type is to be deactivated, these are the steps that must be taken. We'll use the example of wanting to deactivate Ultrasound (there may be one or more Ultrasound Imaging locations) and send everything to General Radiology.

- a. First and foremost, **DO NOT DELETE THE IMAGING TYPE**, **LOCATION(S)**, **OR ASSOCIATED PROCEDURES**. Whenever you have the urge to delete, remember that deletion causes loss of historical data and possible database corruption.
- b. All procedures are associated with an Imaging Type. In our example scenario, any procedure with an Imaging Type of Ultrasound will need to be changed to an Imaging Type of General Radiology by using the option Procedure Enter/Edit. This must be done to prevent additional procedures from being registered in the Ultrasound location(s).

- c. Sign on under each Ultrasound location (or use the Switch Locations option) and check for all pending cases for that location using the option Pending/Hold Rad/Nuc Med Request Log. Print a list of these for tracking purposes.
- d. As patients arrive, sign on under the location(s) you want to send the pending cases to and register each of the cases at the selected location (in this case General Radiology locations). You can either use existing General Radiology locations or add a new one using the option Location Parameter Set-up. Over time, all the requests that were submitted to the Ultrasound location(s) should be registered in a General Radiology location.
- e. You can use the Exam Status Display or Status Tracking of Exams options to find all cases in progress for the Ultrasound location(s). Close out, gradually, all the cases in the locations for the Imaging Type that will be deactivated.

Procedures

1. **Question**: Can I edit the procedures?

Answer: Yes, you can edit the procedure. If the procedure has a CPT Code you cannot edit that procedure's CPT Code, you can only edit the name field if the procedure was added by the site.

If no CPT Code is entered, the procedure type is automatically changed to "Broad". "Series" procedures may have multiple AMIS codes. A "Parent" procedure must have descendants.

If you want to use the same procedure in two different Imaging Locations with different Imaging Types, enter procedure names slightly differently but use the same CPT and AMIS codes.

2. **Question**: Explain the concept behind broad procedures.

"I know that they do not have associated CPT or AMIS codes. I gather that they are to be used PRELIMINARILY, and then changed to detailed procedures. Is this correct? This would seem to make sense only if a site is using remote on-line order requests (either OE/RR or from within the Radiology package)."

Answer: "...when ward/clinic clerks are transcribing written chart orders to VISTA, they are not always able to correlate the provider's written order to a

Radiology Procedure found in the computer. Rather than halting the process and trying to find the provider for a further clarification, etc., a 'Broad' procedure may be ordered with an explanation in the Clinical History field."

3. **Question**: Is there any way to restrict users to only being able to request broad procedures? "The main problem is that there are many similar procedures and requesting physicians usually will not know the difference between them. In fact, they CANNOT know what we will ultimately decide to do for the patient." "...I think the best option for Nuclear Medicine is to list only the broad procedures "bone scan", "renal scan", etc. Then WE will choose which detailed procedure to convert it to. This is one of the ways in which Nuclear Medicine needs differ from those of Radiology, since, as was mentioned, many of their procedures are straight forward enough so that the requester CAN be given a choice of the detailed procedure."

Answer: You can use the Rad/Nuc Med Common Procedure file to implement this. "Enter broad procedures as the only "active" ones in the Rad/Nuc Med Common Procedure file. That way, only the broad ones will be displayed when the ordering party sees the order screen."

"In addition, if your site is using OE/RR, then you should make sure the Quick Order protocols (if there are any set up) only use broad procedures. Although that wouldn't completely prohibit ordering of series and detailed procedures, people are much more likely to choose from what they can see on the screen."

4. **Question**: "When and how does the conversion from a request for a broad procedure to detailed or series take place?"

Answer: The conversion to a detailed or series procedure takes place at registration or during case editing and/or Status Tracking of Exams.

"When the patient is registered, (the Imaging Service) can clear up the murkiness by switching the order to Series or Detailed procedure."

"I would also assign a msg. (procedure message) to each of the Broad procedures that requests a Clinical History that spells out the exact exam they are looking for. This Clinical History could then be used, in conjunction with the expertise of the imaging department, to decide which Detailed procedures to switch to."

5. **Question**: When a physician orders a procedure, I want the physician to also see what the preps are for that procedure so the preps can be ordered at the same time. For instance, an angiogram should "show you how many IV units are needed for that test." How do I do that?

Answer: A message can be displayed with any radiology/nuclear medicine procedure. Build your message(s) using the option Procedure Message Entry/Edit. Then use the option Procedure Enter/Edit to assign the message(s) to the procedure. When ordering a procedure that has one or more messages, the message is displayed as soon as the procedure is selected.

Case Numbers

1. **Question**: I have two active cases with the same case number. Will this be a problem?

Answer: A case number can be assigned to two different cases, although it normally should not happen. Case numbers for completed cases are reused so this is what has taken place:

- Case X is completed and its case number is reused by Case Y.
- Case Y is not completed yet when Case X is "Unverified", which reopens it with the same case number. Both cases are active and identified by the same case number.

This will not cause problems but may be confusing when attempting to pull up a record by its case number. In this situation, be sure to verify the name of the patient whose report you want to work with, to make sure you have the correct case.

2. **Question**: "Case numbers are reaching the ten thousands. Does this mean anything". In the past case numbers were usually two or three digits.

Answer: Case numbers for Completed cases are deleted each week and are then available for reuse. It appears that your cases are not reaching the status of Complete.

Using the option Examination Status Entry/Edit in the Utility Files Maintenance menu, check your exam statuses for Complete within each Imaging Type. Any item under Status Change Requirements that contains a YES is required data for a case to reach the status of Complete. Also check the Status Tracking Functions items; if data must be present for a case to reach Complete, then the user must be asked for that data. Here is an example of what you will find in that option. Note the data in **bold** print.

```
Examination Status Entry/Edit
Select an Imaging Type: nucLEAR MEDICINE

Select an Examination Status: COMPLETE Imaging Type: NUCLEAR MEDICINE Order: 9
...OK? Yes// (Yes)

Name of Current Exam Status: COMPLETE//
Order in sequence of status progression: 9//
```

```
Should this Status appear in Status Tracking ?: NO
        //
User Key needed to move an exam to this status:
Default next status for exam:
Can an exam be cancelled while in this status ?: NO
Generate exam alert for requesting physician ?: NO
         //
Generate Examined HL7 Message:
     Status Change Requirements
     Please indicate which of the following is required
     in order to place an exam into the 'COMPLETE' status:
Technologist Required ?: YES//
Resident Or Staff Required ?: YES//
Detailed Procedure required ?: YES//
Film Entry Required ?: YES//
Diagnostic Code Required ?: NO//
Camera/Equip/Rm Required ?: YES//
     (Radiopharmaceutical requirements)
Radiopharms & Dosages required ?: YES//
Activity Drawn required (Radiopharm) ?: NO//
Drawn Dt/Time and Person required (Radiopharm) ?: NO//
Admin Dt/Time/Person required (Radiopharm) ?: NO//
Route/Site required (Radiopharm administered) ?: NO//
Lot No. required (Radiopharm) ?: NO//
Volume/Form required (Radiopharm) ?: NO//
   WARNING: You must use the reporting feature of the system as
            a prerequisite for the following requirements:
Report Entered required ?: YES//
Impression required ?: YES//
Verified Report required ?: YES//
     Status Tracking Functions
     _____
     Please indicate which of the following should be asked when
     changing an exam's status to 'COMPLETE' while using
     the 'status tracking' feature:
Ask for Technologist ?: YES//
Ask For Interpreting Physician ?: YES//
Ask for Procedure ?: YES//
Ask For Film Data ?: YES//
Ask For Diagnostic Code ?: YES//
Ask For Camera/Equip/Rm ?: YES//
Ask for User Code ?: NO//
Ask Medications & Dosages ?: YES//
Ask Medication Admin Dt/Time & Person ?: YES//
     (Radiopharmaceutical functions)
Ask Radiopharmaceuticals and Dosages ?: YES//
```

Ask Activity Drawn (Radiopharmaceutical) ?: YES//
Ask Drawn Dt/Time & Person (Radiopharmaceutical) ?: YES//

```
Ask Admin Dt/Time & Person (Radiopharmaceutical) ?: YES//
Ask Route/Site of Admin (Radiopharm) ?: YES//
Ask Lot No. (Radiopharm) ?: YES//
Ask Volume/Form (Radiopharm) ?: YES//

Management Report Criteria

Please indicate which of the following workload reports should include exams with the 'COMPLETE' status:

Clinic Report should include exams of this status ?: YES// ^
```

One site keeps on top of their Incomplete cases by running the Incomplete Exam Report within the Daily Management Reports menu several times a day. It helps them "to track down lost or misplaced films in order to get them read and also to see if a Technician accidentally forgot to case edit." They correct any errors daily.

Requesting/Scheduling Requests

1. **Question**: Could someone provide a step-by-step procedure used for requesting/scheduling requests, complete with who does what?

Answer: The following is an example of the work flow at one VAMC:

- a. Request comes to Radiology via VISTA.
- b. Request goes into Waiting to be Scheduled bucket. (i.e., Request status = Pending)
- c. Techs are responsible daily for scheduling their area (e.g. CT tech gets all CT requests, US gets US requests, etc.) All other requests go to the regular staff techs for scheduling.
- d. Inpatients get called down the same day, unless they are for Fluoro or some other exam needing prep. Then they are given a date in the computer using the option Schedule a Request. **Note**: we do NOT give a time for inpatient appointments. They simply are ON CALL. As a rule, they are called first before our outpatients show up.
- e. Outpatient requests are given a date and a time in the Schedule a Request option and also we go in to the Scheduling package under Make an Appointment and give it the same date and time. (Use of the Scheduling package is optional.)
- f. All outpatients are printed appointment letters by the techs and our letters have the preps on them also. (Done through the Scheduling package)
- g. If a patient calls and cannot come in, the tech puts the request on hold using the Hold a Request option and you can choose whatever reasons you want from your file. They also write on the request that the patient was unable to come. It then goes back to the Waiting to be scheduled bucket, and gets rescheduled as above with letters and all. If a patient fails to report, we do the same thing and give him another appointment.

h. If a patient fails to report two times, we use the option Cancel a Request and return the request to the attending physician.

"This makes it great because seldom are there any requests in a Pending status for very long, and you can easily see if you have missed someone. Also, once it is in a Scheduled status you can easily tell your workload for any given day, and also find out by date or by patient when a test is scheduled."

2. **Question**: How are routine requests handled for No Shows of patients that have called in for a re-appointment? Is there a way for the user to edit the desired/requested date of the request if the x-ray is desired on the day of the appointment?

Answer: "Under the Rad/Nuc Med Order/Entry menu you can schedule a request. It will list pending requests and you choose which ones you want to schedule, then it will prompt you for date and time desired. One thing, if these requests already have a schedule date, you must place them on hold first, then move to schedule."

3. **Question**: For those who are using the Schedule a Request option, do your "MAS people have it and are you keeping a suspense file for your hard copies of rescheduled requests?"

Answer: One site stated they use a suspense file: "The only reason we do is if the computer system should happen to go down, we can enter from the printed requests in our back-up system and continue business as usual without having to try to get patient charts for routine x-rays.

Another site had this to say: "We put the hard copies in an accordion type file for the day of the month it was scheduled. If a patient is a no show, we put the request on hold and reschedule it once if they fail to report a second time, we d/c the request with reason: Patient failed to report."

4. **Question**: "Once the X-ray request has been scheduled, the scheduled date on the request cannot be changed. Correct?"

Answer: "Incorrect. You **can** change the scheduled date again, but it takes a few extra steps. If you have a request for which the scheduled date must be changed, use the Rad/Nuc Med Order/Entry menu and place the request on Hold (Hold a Request option). Then use the Schedule a Request option and enter the new date.

Printing/Batching/Etc.

Question: "Is it possible to have Angiograms and Venograms queued to a printer in the Vascular Lab after the report has been transcribed?"

Answer: One view: "Are angios and venos done in a separate location? If so, they will print to the report printer assigned for that location. Or, you can Batch all the angios and venos in the same Batch, and then print this Angio/Veno Batch anytime to any printer.

Distribution Queues are for distribution of **verified reports** and are broken up by inpatient/outpatient/file room/etc. You wouldn't be able to get angios/venos broken out by using Distribution Queues."

Another view: "Your transcriptionists could use Batch names such as 'transcriptionist name/radiology section/date' or whatever scheme your folks feel comfortable with to distinguish the batches. Read up on Batches in the ADPAC manual, especially the management of same so you can set up a proper print/purge scenario at your VAMC that will get the job done but not allow the Report Batches file to grow too large."

In addition to the above, we would like to stress the fact that verifying a report triggers its placement into the Distribution Queue. If a site leaves reports Released/Unverified, it will not be able to use the Distribution Queue.

Purging Records

1. **Question**: "Can anyone doing transcription within Radiology Service give me some input reference time frame for purging as follows:

Activity log cut off.

Report log cut off.

Clinical history cut off.

Tracking time cut off. (Exam status changes)

After the report goes the impression remains forever -- IRMS needs to get some purging done and I would sure like some input if possible from other stations."

Answer: There are many ways to handle purging. Here are a few of the examples from the sites:

- a. We "have the report purge set to coincide with when we purge the film files. That is 5 years."
- b. Activity log cut off 999 days

Report cut off 999 days Clinical history 999 days Tracking cut off 999 days Order data cut off 30 days

c. Activity log cut off 90days

Report cut-off 1825 days
Clinical History 90 days
Tracking time cut-off 90 days
Order data cut-off 30 days

- d. "We started full purges of everything over a year old but have now moved to doing it every six months. Do be careful that you have ensured that impressions are supplied. We were getting a lot of 'Refer to report such-and-such' and, of course, it would have been purged too. You will only be able to purge those exams which are Complete, so there may be some clean up to do before you can proceed. If you're really timid, you could archive the global on tape, until you're confident of your procedure."
- 2. **Question**: What do you do with report text such as "See above" or "Refer to..." when purging?

Answer: One view: "A very important issue was raised. If you plan to purge report text after a period of time, leaving only the impression, make sure that the radiologists know this (and keep reminding them). I have seen many dictation where the impression read something like "As above", which obviously is not that helpful if "the above" has been deleted! In fact, you must specifically warn those who dictate not to do this."

Another view: "We use the Indicate No Purge option when the physician dictates the impression as "see above report". " **Note**: Transcriptionists need to be trained to do this.

Verifying Reports

1 **Question**: "I would like some input about who can Verify a report. Can it be only the radiologist who dictated and signed the report? Or can a radiologist sign other radiologists' reports?"

Answer: "Radiologists can verify their own reports and, if the site parameters are turned on, a radiologist may also be given the ability to verify other radiologists' reports."

2. **Question**: "How can I get the titles of contract radiologists to appear on the report without local modifications?"

Answer: "The Rad/Nuc Med software uses the Signature Block Title field in the New Person file #200 as the title printed to the right of each primary and secondary staff and resident's name. You can edit this field so that it shows that a person classified as Radiology/Nuclear Medicine resident or staff is a contract radiologist."

(See the Troubleshooting chapter for more detail)

Hold/Cancelled/Pending/Scheduled

1. **Question**: Does anyone use the Pending/Hold Rad/Nuc Med Request Log (used to be Log of Pending Radiology Requests)?

Answer: "We use this list extensively. We also print it to each ward every evening at 8 pm. The nurses love it. They can run down the list and make sure that patients that are to be NPO, are. We haven't had a miss in the last year. We also print to the lab a 5 am each morning. They look for the CTs and Specials and draw these patients early. By the time the patient is ready to be done, the lab work is ready. It has really helped."

2. **Question**: "What do you do in the event the patient is Scheduled. These don't show up on a Pending report. We have so many patients scheduled and the floors "forget" to give the patients preps.

Answer: One site: "There is a option Log of Scheduled Requests by Procedure. We print one for us and each ward prints one daily (Ward/Clinic Scheduled Request Log) to prep the next day's patients. This will group GIs, CTs, BEs, Chests, etc. together and gives the patient name, SSN, procedure, patient

location and time of appointment. For inpatients when we schedule a request we give date only with no time since they are On Call. We do give outpatients a time of course. Since we have instituted this, we have had little or no repeat preps, since all shifts have this option and it saved us work all around."

Another site: "We are on full ward and clinic order entry which includes scheduling. The wards just call for special things like BEs, IVPs, etc. and we give them a date, and they enter the request in the pending file for that date.

3. **Question**: "After you schedule a request, there should be a way of looking at it and seeing the time and date of the appointment. As it stands now, you can see that it is scheduled but unless you print out days forward, you cannot find the appointment."

Answer: "Go under Detailed Request Display either by patient or date. If you have ... given it a time and date, it will appear on the Detailed Request Display. Just select the # of the request (it will have an "s" in the column for status)". Enter Yes at the prompt "Do you wish to display request status tracking log?".

4. **Question**: How is workload credited for the Transcription Report?

Answer: The transcriptionist that initiates/creates the report is credited with the total number of lines regardless of how many transcriptionists may have contributed to the report. The length of each line is summed and divided by 75 to determine the total number of lines. The number of lines counted is always the current number of lines in the report. So, if lines have been added, changed, or deleted, only the final number of lines are counted at the time the Transcription Report is run.

Doubling Workload Credit

Question: What is the effect of the modifier on special reports? When exactly is the workload credit doubled in the system?

Answer: The following reports do a calculation that will double workload counts if the Bilateral modifier is used and the AMIS weight multiplier associated with the procedure AMIS pair is 1 or "blank":

Physician Report [RA WKLPHY] Resident Report [RA WKLRES] Exam Room Report [RA WKLROOM] Staff Report [RA WKLSTAFF]

Technologist Report [RA WKLTECH]

PTF Bedsection Report [RA LWKLBEDSEC]

Clinic Report [RA LWKLCLINIC]

Service Report [RA LWKLSERVICE]

Sharing Agreement/Contract Report [RA LWKLSHARING]

Ward Report [RA LWKLWARD]

Cost Distribution Report [RA CDR REPORT]

Film Usage Report [RA FILMUSE]

AMIS Code Dump by Patient [RA AMISDUMP]

AMIS Report [RA AMIS]

Detailed Procedure Report [RA WKLPROCEDURE]

Effect of Portable and Operating Room Modifiers

Question: Which workload reports are affected and how?

Answer: The AMIS Report and the AMIS Code Dump by Patient are affected. In the AMIS Code Dump by Patient when an AMIS Code is selected all procedures not in that AMIS category are usually bypassed. However, if the user selects either AMIS category 25 (Operating Room) or 26 (Portable) then all exams that meet the other criteria, regardless of the AMIS code of the procedure performed, will be included if the exam contains an Operating Room or Portable type of modifier. Therefore, if an ankle x-ray (AMIS code 8) is done, and a modifier of the "Operating Room" type was entered for the exam, then that exam will show up under and AMIS Code Dump by Patient when AMIS code 25 is selected and when AMIS Code 8 is selected.

In the AMIS Report and Detailed Procedure Report (Detailed Procedure Workload Report) the same procedure described above would be counted under both AMIS code 8 and AMIS code 25.

Complications

Question: "How would you delete a complication done months ago on an exam? This patient who has never had a reaction to contrast dye has a contrast reaction message come up when a CT is done."

Answer: Use the Update Patient Record option. Enter NO at the question Contrast Medium Allergy. The Adverse Reaction Tracking package will automatically create an 'Entered in Error record.

Contracting Out

Question: "Is anyone contracting out mammograms and how are they handled through VISTA?"

Answer:

One site: "Once the appointment is made, our scheduling clerk enters a request through VistA from the handwritten consult. When the results are mailed back, the study is registered in to get a case no. Under case edit, No Room used is entered along with No Films with a 0. After this is completed, the clerk then enters the report into the system from the copy stating at the beginning of the report that THIS STUDY WAS PERFORMED AT ... BY DRS., etc. The report is then verified. MRIs are handled the same way in order to make the reports for these exams readily available to the medical staff."

Another site: "We also indicate under category (IP, OP, etc.) that they are a contract source and the contractor's name. That way they are separated from the actual workload performed on site."

Research

Question: "Is there any facility out there involved in X-raying body parts for research purposes?"

Answer: "We register any patient involved in research as a Research patient in the Exam Category field. We get credit for this in our CDR report."

General Editing Helpful Hints

Question: "How can I keep text from wrapping when using the Screen Editor?"

Answer: On each new line, enter a space <sp> before entering the text as in the following example:

```
When ordering this procedure, we recommend the following steps be taken:
<sp>1) Do this first.
<sp>2) Next do this.
<sp>3) Etc.
```

Question: "How do I add a second identical entry in a field that accepts multiple entries?"

Answer: Enter the second, third, etc. entry in quotes as shown here:

```
Select DESCENDENTS: FOOT 2 VIEWS// ?

Answer with DESCENDENTS

Choose from:

ANKLE 2 VIEWS

FOOT 2 VIEWS

TOE(S) 2 OR MORE VIEWS

You may enter a new DESCENDENTS, if you wish
Enter a procedure name which is related to the parent procedure.

Select only active detailed or series procedures of the same imaging type as the parent.

Answer with RAD/NUC MED PROCEDURES NAME, or CPT CODE, or SYNONYM
Do you want the entire RAD/NUC MED PROCEDURES List? n (No)

Select DESCENDENTS: FOOT 2 VIEWS//"FOOT 2 VIEWS"
```

June 2019

XVI. Stop Task (TaskMan)

One of the new features in this version is the ability of the user to stop tasks. All reports that can consume a significant amount of CPU/print time are now stoppable through the Stop Task action of the Taskman User [XUTM USER] option under the User's Toolbox [XUSERTOOLS] menu. The enhanced report logic checks for a stop flag during processing that is done before printing actually begins as well as during printing. You can get to this menu from any menu prompt in the Radiology/Nuclear Medicine software simply by typing TBOX. Here is an example of how it works. Note that TBOX is a Common Option under the menu.

```
Exam Entry/Edit Menu ...
         Films Reporting Menu ...
         Management Reports Menu ...
         Outside Films Registry Menu ...
         Patient Profile Menu ...
         Radiology/Nuclear Med Order Entry Menu ...
          Supervisor Menu ...
          Switch Locations
          Update Patient Record
          User Utility Menu ...
Select Rad/Nuc Med Total System Menu Option: ??
          Exam Entry/Edit Menu ... [RA EXAMEDIT]
          Films Reporting Menu ... [RA RPT]
          Management Reports Menu ... [RA MGTRPTS]
          Outside Films Registry Menu ... [RA OUTSIDE]
          Patient Profile Menu ... [RA PROFILES]
          Radiology/Nuclear Med Order Entry Menu ... [RA ORDER]
          Supervisor Menu ... [RA SUPERVISOR]
          Switch Locations [RA LOC SWITCH]
          Update Patient Record [RA PTEDIT]
          User Utility Menu ... [RA USERUTL]
You can also select a secondary option:
          Rad/Nuc Med Total System Menu ... [RA OVERALL]
          Radiology/Nuclear Med Order Entry Menu ... [RA ORDER]
          Supervisor Menu ... [RA SUPERVISOR]
          Utility Files Maintenance Menu ... [RA MAINTENANCE]
Or a Common Option:
  TBOX User's Toolbox ... [XUSERTOOLS]
   Δ7.7
         View Alerts [XQALERT]
              Press 'RETURN' to continue, '^' to stop: <RET>
          Continue [XUCONTINUE]
            **> Reverse lock ZZLUKE
          Halt [XUHALT]
```

```
MailMan Menu ... [XMUSER]
         Restart Session [XURELOG]
         Time [XUTIME]
         Where am I? [XUSERWHERE]
Select Rad/Nuc Med Total System Menu Option: tbox User's Toolbox
         Display User Characteristics
         Edit User Characteristics
         Electronic Signature code Edit
         Menu Templates ...
         Spooler Menu ...
         Switch UCI
         TaskMan User
         User Help
Select User's Toolbox Option: taskMan User
Select TASK: ??
Please wait while I find your tasks...searching...finished!
1: (Task #35437) DQ^RAFLH, DQ^RAFLH. Device LINE. POC, POC.
  From 12/09/96 at 13:58, By you.
  Waiting for device LTA1707:
2: (Task #35591) ENTASK^RAMAINP, Rad/Nuc Med Active Procedures (Long).
  Device LINE. POC, POC. From 12/13/96 at 8:42, By you. Unscheduled by you.
______
3: (Task #35680) DQ^RAFLH, DQ^RAFLH. Device LINE. POC, POC.
  From 12/16/96 at 13:08, By you. Unscheduled by you.
4: (Task #35823) DQ^RAFLH, DQ^RAFLH. Device LINE. POC, POC.
  From 12/20/96 at 8:26, By you. Unscheduled by you.
5: (Task #35825) PRTORD^RAORDQ, PRTORD^RAORDQ. Device FIELD SUPPORT. POC, POC.
  From 12/20/96 at 8:30, By you. Completed 12/20/96 at 8:31.
______
6: (Task #35826) DQ^RAFLH, DQ^RAFLH. Device LINE. POC, POC.
  From 12/20/96 at 8:31, By you. Unscheduled by you.
Press RETURN to continue or '^{\prime} to exit: ^{\bullet}
Select TASK: 1 DQ^RAFLH
             Taskman User Option
                  Display status.
                  Stop task.
                  Edit task.
                  Print task.
                  List own tasks.
                  Select another task.
             Select Action (Task # 35437): stop Stop task.
Task unscheduled and stopped.
```

XVI-2

```
Taskman User Option
                    Display status.
                    Stop task.
                   Edit task.
                    Print task.
                    List own tasks.
                    Select another task.
               Select Action (Task # 35437): <RET>
          Display User Characteristics
         Edit User Characteristics
         Electronic Signature code Edit
         Menu Templates ...
         Spooler Menu ...
         Switch UCI
         TaskMan User
         User Help
Select User's Toolbox Option: <RET>
         Exam Entry/Edit Menu ...
         Films Reporting Menu ...
         Management Reports Menu ...
         Outside Films Registry Menu ...
         Patient Profile Menu ...
         Radiology/Nuclear Med Order Entry Menu ...
         Supervisor Menu ...
          Switch Locations
          Update Patient Record
          User Utility Menu ...
```

Select Rad/Nuc Med Total System Menu Option:

XVII. Glossary

Active An order status that occurs when a request to perform a

procedure on a patient has been registered as an exam, but

before it has reached a status of Complete.

Activity log A log of dates and times data was entered and/or changed.

The Radiology/Nuclear Medicine system is capable of maintaining activity logs for reports, exam status changes, imaging type parameter changes, purge dates, outside

film registry activity, and order status changes.

Alert Alerts consist of information displayed to specific users

triggered by an event. For example, alerts pertaining to Rad/Nuc Med include the Stat Imaging Request alert, an

Imaging Results Amended alert, and an Abnormal

Imaging Results alert. The purpose of an alert is to make a user aware that something has happened that may need attention. Refer to Kernel and OE/RR documentation for

more information about alerts.

¹AMIS code For imaging, one of 27 codes used to categorize

procedures, calculate workload crediting and weighted work units. AMIS codes are determined by VA Central Office and should not be changed at the medical centers.

AMIS weight multiplier A number associated with a procedure-AMIS code pair

that is multiplied by the AMIS code weighted work units. If the multiplier is greater than 1, a single exam receives

multiple exam credits.

Attending physician The Radiology/Nuclear Medicine software obtains this

data from the MAS package, which is responsible for its entry and validity. Refer to the VistA MAS package documentation for more information and a description of

the meaning of this term as it applies to VISTA.

Batch In the Radiology/Nuclear Medicine system, a batch is a

set of results reports. Transcriptionists may create

batches to keep similar reports together and cause them to print together. One possible purpose might be to print all

reports dictated by the same physician together.

¹ Patch RA*5*45 March 2004 AMIS Code

Glossary

Bedsection

See PTF Bedsection.

Bilateral

A special type of modifier that can be associated with an exam, a procedure, or an AMIS code. When an exam is bilateral due to one of the aforementioned associations, workload credit and exam counts are doubled for that exam on most workload and AMIS reports.

Broad procedure

A non-specific procedure that is useful for ordering when the ordering party is not familiar enough with imaging procedures to be able to specify the exact procedure that is to be performed. Before an exam status can progress to Complete, the imaging service must determine a more specific procedure and change the exam procedure to reflect the actual Detailed or Series procedure done. Depending on site parameters, broad procedures may or may not be used at a given facility. Also see Detailed and Series procedure.

Bulletin

A special type of mail message that is computergenerated and sent to a designated user or members of a mail group. Bulletins are usually created to inform someone of an event triggered by another user's data entry, or exam and request updates that require some action on the part of the bulletin recipient.

CPT code

See Current Procedural Terminology

Camera/Equipment/Room

The specific room or piece of equipment used for a patient's imaging exam. Each is associated with one or more imaging locations. The Radiology/Nuclear Medicine system supports, but does not require users to record the camera/equipment/room used for each exam.

Cancelled

A status that can be associated with an exam. Also see Discontinued.

Case number

A computer-generated number assigned to the record generated when one patient is registered for one procedure at a given date/time. Note that when multiple procedures are registered for a patient at the same date/time, each procedure will be given a different case number. Case numbers will be recycled and reused by a new patient/procedure/date instance when the exam attains a Complete status.

Category of exam For the purposes of this system, category of exam must be

Outpatient, Inpatient, Contract, Sharing, Employee, or Research. Several workload and statistical reports print exam counts by category. Others use the category to determine whether exams should be included on the

report.

Clinic Hospital locations where outpatients are cared for. In

VISTA, clinics are represented by entries in the Hospital

Location file (#44). Radiology/Nuclear Medicine

Imaging Locations, represented by entries in the Imaging

Location file (#79.1), are a subset of the Hospital

Location file.

Clinical history Data entered in the Radiology/Nuclear Medicine system

during exam ordering and exam edit. Usually entered by the requesting party to inform the imaging service why the exam needs to be done and what they hope to find out

by doing the exam.

Clinical history message Text that, if entered in system parameter setup, will

always display when the user is prompted for clinical history. Generally used to instruct the user on what they

should enter for clinical history.

Common procedure A frequently ordered procedure that will appear on the

order screen. Up to 40 per imaging type are allowed by the system. Other active Rad/Nuc Med procedures are selectable for ordering, but only the ones designated as common procedures and given a display sequence

number will be displayed prior to selection.

Complete A status that can be attained by an order or an exam.

Complication A problem that occurs during or resulting from an exam,

commonly a contrast medium reaction.

Contract A possible category of exam when imaging services are

contracted out.

Contrast medium A radio-opaque injectable or ingestible substance that

appears on radiographic images and is helpful in image interpretation. It is used in many imaging procedures.

Contrast reaction message A warning message that will display when a patient who

has had a previous contrast medium reaction is registered for a procedure that uses contrast media. The message text is entered during Rad/Nuc Med division parameter

setup.

Credit stop code See Stop Code. Also see MAS package documentation.

Current Procedural Terminology

(CPT)

A set of codes published annually by the American Medical Association which include Radiology/Nuclear Medicine procedures. Each active detailed or series procedure must be assigned a valid, active CPT code to facilitate proper workload crediting. In VISTA, CPT's

are represented by entries in the CPT file #81.

Descendent A type of Rad/Nuc Med procedure. One of several

associated with a 'Parent' type of procedure. The descendent procedures are actually registered and

performed. Also see Parent.

Desired date (of an order)

The date the ordering party would like for the exam to be

performed. Not an appointment date. The imaging service is at liberty to change the date depending on their

availability.

Detailed procedure A procedure that represents the exact exam performed,

and is associated with a CPT code and an AMIS code.

Diagnostic code Represented, for the purposes of this system, by entries in

the Diagnostic Codes file (#78.3). Diagnostic codes describe the outcome of an exam, such as normal, abnormal. A case may be given one or more (or no)

diagnostic code(s).

Discontinued An order status that occurs when a user cancels an order.

Distribution queue A mechanism within the Radiology/Nuclear Medicine

system that facilitates printing results reports at various hospital locations, such as the patient's current ward or

clinic, the file room, and medical records.

Division, Rad/Nuc Med A subset of the VistA Institution file (#4). Multi-

divisional sites are usually sites responsible for imaging at

more than one facility.

Draft A report status that is assigned to all Rad/Nuc Med results

reports as soon as they are initially entered into the

system, but before they are changed to a status of Verified

or (if allowed) Released/Not Verified.

DSS ID Formerly Stop Code associated with each procedure, now

DSS ID associated with each Imaging Location.

Electronic signature code A security code that the user must enter to identify

him/herself to the system. This is required before the user is allowed to electronically verify Rad/Nuc Med reports.

This is not the same as the Access/Verify codes.

Exam label One of 3 types of labels that can be printed at the time

exam registration is done for a patient. Also see jacket

label, flash card.

Exam set An exam set contains a Parent procedure and its Detailed

or Series descendent procedures. Requesting a Parent will automatically cause each descendent to be presented for registration as separate cases under a single visit date and

time.

Exam status The state of an exam that describes its level of progress.

Valid exam statuses are represented in this system by entries in the Examination Status file (#72). Examples are ordered, cancelled, complete, waiting for exam, called for exam, and transcribed. The valid set of exam statuses can, to a degree, be tailored by the site. There are many parameters controlling required data fields, status tracking

and report contents that are determined when the

parameters of this file are set up. There are separate and

different set of statuses for requests and reports.

Exam status parameter setup See Exam status.

Exam status time The date/time when an exam's status changes, triggered

by exam data entry that can be done through over a dozen

different options.

Film size Represented by entries in the Film Sizes file (#78.4) in

this system. Used to facilitate film use/waste tracking.

Glossary

Flash card One of 3 labels that can be generated at the time an exam

is registered for a patient. The flash card was named because it can be photographed along with an x-ray, and its image will appear on the finished x-ray. Helpful in marking x-ray images with the patient's name, SSN, etc.

to ensure that x-rays are not mixed up.

Footer The last lines of the results report, the format of which

can be specified using the Lbl/Hdr/Ftr formatter.

Format The specification for print locations of fields on a printed

page. In this system, print formats can be specified using

the Lbl/Hdr/Ftr formatter.

Header The top lines of the results report, the format of which can

be specified using the Lbl/Hdr/Ftr formatter.

Health Summary Refers to a report or VistA software package that

produces a report showing historical patient data. Can be configured to meet various requirements. Refer to the Health Summary documentation for more information.

Hold An order status occurring when a users puts an order on

hold, indicating that a study should not yet be done or scheduled, but that it will likely be needed in the future.

Hospital location Represented in VistA by entries in the Hospital Location

file (#44). Rad/Nuc Med locations are a subset of the

Hospital Location file.

Imaging location One of a subset of Hospital Locations (See Hospital

location) that is represented in the Imaging Location file

#79.1, and is a location where imaging exams are

performed.

Imaging type

For the Rad/Nuc Med software, the set of valid imaging types is:

ANGIO/NEURO/INTERVENTIONAL

GENERAL RADIOLOGY MAMMOGRAPHY

NUCLEAR MEDICINE

ULTRASOUND VASCULAR LAB

CARDIOLOGY STUDIES (NUC MED)

CT SCAN

MAGNETIC RESONANCE IMAGING

These are the imaging types that are supported by this version of the software. Each imaging location and procedure may be associated with only one imaging type.

Impression

A short description or summary of a patient's exam results report. Usually mandatory data to complete an exam. The impression is not purged from older reports even though the lengthier report text is.

Inactivate

The process of making a record in a file inactive, usually by entering an inactive date on that record or deleting a field that is necessary to keep it active. When a record is inactive, it becomes essentially "invisible" to users. Procedures, common procedures, modifiers, and exam statuses can be inactivated.

Inactivation date

A date entered on a record to make it inactive. See Inactivate.

Information Resources Management The service within VA hospitals that supports the installation, maintenance, troubleshooting, and sometimes local modification of VistA software packages and the hardware that they run on.

Interpreting physician (also Interpreting Resident, Interpreting Staff)

The resident or staff physician who interprets exam images.

IRM

See Information Resources Management.

Jacket label

One of the 3 types of labels that can be generated at the time an exam is registered for a patient. Usually affixed to the x-ray film jacket. (See also exam label, flash card.) Glossary

Key See security key.

Label print fields Fields that are selectable for printing on a label, report

header, or report footer on formats that are designed using

the Lbl/Hdr/Ftr formatter.

Lbl/Hdr/Ftr formatter

The name given to the option/mechanism that allows

users to define formats for labels and for headers and footers on results reports. Users can specify which fields to print at various columns and lines on the label or report

header/footer.

Mode of transport The patient's method of moving within the hospital,

(ambulatory, wheelchair, portable, stretcher) designated at

the an exam is ordered.

Modifier Additional information about the characteristics of an

exam or procedure (such as bilateral, operating room, portable, left, right). Also see bilateral, operating room,

portable.

No purge indicator A flag that can be set on the exam record to force the

purge process to bypass the record. Guarantees that the record will not be purged when a historic data purge is

scheduled by IRM. Also see Purge.

Non-credit stop code Certain stop codes, usually for health screening, that do

not count toward workload credit. If a non-credit stop code is assigned to a procedure, another credit stop code

must also be assigned. Also see Stop code.

OE/RR See Order Entry/Results Reporting.

On-line verification The option within the Radiology/Nuclear Medicine

package that allows physicians to review, modify, and

electronically sign patient result reports.

Operating room A special type of procedure modifier, that, when assigned

to an exam will cause the exam to be included in

workload/AMIS reports under both the AMIS code of the

procedure and under the AMIS code designated for

Operating Room.

Order The paper or electronic request for an imaging exam to be

performed.

Order Entry The process of requesting that one or more exams be

performed for a patient. Order entry for

Radiology/Nuclear Medicine procedures can be

accomplished through a Rad/Nuc Med software option or through a separate VistA package, Order Entry/Results

Reporting (OE/RR).

Order Entry/Results Reporting A VistA package that performs that functions of

ordering for many clinical packages, including

Radiology/Nuclear Medicine.

Outside films registry A mechanism in this system that allows users to track

films done outside of the medical center. This can also be accomplished through the VistA Records Tracking package. Refer to Records Tracking documentation for

more information.

Parent procedure A type of Rad/Nuc Med procedure that is used for

ordering purposes. It must be associated with Descendent procedures that are of procedure type Detailed and/or Series. At the time of registration the descendent procedures are actually registered. Setting up a parent procedure provides a convenient way to order multiple related procedures on one order. Parent/descendent procedure relationships must be set up ahead of time during system definition and file tailoring by the ADPAC.

Pending An order status that every Rad/Nuc Med order is placed

in as soon as it is ordered through this system's ordering option. This system also receives orders from the VistA

OE/RR system and places them in a pending status.

Portable A special type of procedure modifier, that, when assigned

to an exam will cause the exam to be included in

workload/AMIS reports under both the AMIS code of the

procedure and under the AMIS code designated for

Portable.

Pre-verification The process whereby a resident reviews a report and

affixes his/her electronic signature to indicate that the report is ready for staff (attending) review, facilitated through an option in this system for Resident Pre-

verification.

Primary Interpreting Staff/ Resident The attending or resident primarily responsible for the

interpretation of the case. (See also Secondary

Interpreting Staff/Resident.)

Primary physician The Radiology/Nuclear Medicine software obtains this

data from the MAS package, which is responsible for its entry and validity. Refer to the VistA MAS package documentation for more information and a description of

the meaning of this term as it applies to VISTA.

Principal clinic For the purposes of the Radiology/Nuclear Medicine

system, this term is usually synonymous with 'referring clinic'. However, for the purposes of crediting, it is defined as the DSS (clinic/stop) code that is associated with the imaging location where the exam was performed.

Printset A printset contains a Parent procedure and its Detailed or

Series descendent procedures. If the parent is defined to be a printset, the collection and printing of all common report related data between the descendents is seen as one

entity.

Problem draft A report status that occurs when a physician identifies a

results report as having unresolved problems, and designates the status to be Problem Draft. Depending on site parameters, a report may be designated as a Problem Draft due to lack of an impression. Also see Problem

statement.

Problem statement When a results report is in the Problem Draft status, the

physician or transcriptionist is required to enter a brief statement of the problem. This problem statement

appears on report displays and printouts.

Procedure For the purposes of this system, a medical procedure done

with imaging technology for diagnostic purposes.

Procedure message Represented in this system by entries in the Rad/Nuc Med

Procedure Message file (#71.4). If one or more procedure messages are associated with a procedure, the text of the messages will be displayed when the procedures is ordered, registered, and printed on the request form. Useful in alerting ordering clinicians and imaging personnel of special precautions, procedures, or

requirements of a given procedure.

Procedure type A characteristic of a Rad/Nuc Med procedure that affects

exam processing and workload crediting. See Detailed,

Series, Broad, and Parent.

PTF Bedsection See MAS documentation.

Purge The process that is scheduled at some interval by IRM to

purge historic computer data. In this system, purges are done on results report text, orders, activity logs, and

clinical history.

Registration (of an exam)

The process of creating a computer record for one or

more patient/procedure/visit date-time instances. Usually done when the patient arrives at the imaging service for

an exam.

Released/not verified A results report status that may or may not be

implemented at a given medical center. Reports in this status may be viewed or printed by hospital staff outside

of the imaging service.

Report batch See Batch.

Report status The state of a report that describes its progress level.

Valid report statuses in this system are Draft, Problem draft, Released/not verified (if the site allows this status), and Verified. The status of a report may affect the status of an exam. Also see Exam status. Exams and requests

each have a separate and different set of statuses.

Request Synonymous with order. See Order.

Request status The state of a request (order) that describes its progress

level. Valid request statuses in this system are

Unreleased(only if created through OE/RR), Pending, Hold, Scheduled, Active, Discontinued, and Complete. Reports and exams each have a separate and different set

of statuses.

Request urgency Data entered at the time an exam is ordered to describe

the priority/criticality of completing the exam quickly

(i.e. Stat, Urgent, Routine).

Staff/Resident

Requesting location Usually the location where the patient was last seen or

treated (an inpatient's ward, or an outpatient's clinic). All requesting locations are represented by an entry in the VistA Hospital Location file (#44). The requesting location may be, but is usually not an imaging location.

Requesting physician The physician who requested the exam.

Research source A research project or institution that refers a patient for a

Radiology/Nuclear Medicine exam.

Scheduled An order status that occurs when imaging personnel enter

a date when the exam is expected to be performed.

Secondary Interpreting This generally refers to an attending/resident that assisted

or sat in on review of the case, but is not primarily responsible for it. It may also be used to indicate a second reviewer of the case, for quality control or peer

review purposes.

Security key Represented by an entry in the VistA Security Key file.

Radiology/Nuclear Medicine keys include RA MGR, RA

ALLOC, and RA VERIFY. Various options and

functionalities within options require that the user "own"

a security key.

Staff Imaging Attending.

Staff review (of reports) The requirement where an attending imaging physician is

required to review the reports written by a resident

imaging physician.

Standard report Represented in this system by entries in the Standard

Reports file (#74.1). Standard reports can be created by the ADPAC during system definition and set-up. If the

division setup specifies that they are allowed,

transcriptionists will be offered a selection of standard report text and impressions to minimize data entry effort.

Status tracking

The mechanism within this system that facilitates exam tracking from initial states to the complete state. ADPACs must specify during exam status parameter setup which statuses will be used, which data fields will be required to progress to each status, which data fields will be prompted, and exams of which statuses will be included on various management reports.

Stop code

Member of a coding system designed by VA Central Office to aid in determining workload and reimbursement of the medical centers. Stop codes are controlled by VA Central Office MAS. Imaging stop codes are represented by entries in the Valid Imaging Stop Code file #71.5. Imaging stop codes are a subset of the VistA Clinic Stop file #40.7. See MAS documentation for more information.

Synonym

In the Radiology/Nuclear Medicine package, synonyms are alternate terms that can be associated with procedures for the purposes of convenient look-up/retrieval. A given procedure may have more than one synonym, and a given synonym may be used for more than one procedure.

Technologist

Radiology/Nuclear Medicine personnel who usually are responsible for performing imaging exams and entering exam data into the system.

Time-out

The amount of time allowed before a user is automatically logged out of the system if no keystrokes are entered. This is a security feature, to help prevent unauthorized users from accessing your $V \mbox{IST} A$ account in case you forget to log off the system or leave your terminal unattended.

Transcribed

An exam status that may occur when a results report is initially entered into the system for an exam. If this status is not activated at the site, it will not occur.

Unreleased

An order status that occurs when an exam order is created, but no authorization to carry out the order has been given. This is possible only if the order is created through the OE/RR software.

Glossary

Verification For the purposes of this system, the process of causing a

results report to progress to the status of Verified. This happens when a physician affixes his/her electronic signature to the report, or when a transcriptionist, with the proper authorization, enters the name of a physician who has reviewed and approved a report. This is analogous to

a physician signing a paper report.

Verified A results report status that occurs at the time of

verification. A report is verified when the interpreting physician electronically signs the report or gives his/her authorization that the report is complete and correct. Also

see Verification.

VISTA Veterans Health Information Systems and Technology

Architecture. Formerly known as DHCP.

Waiting for exam An exam status that occurs as soon as the exam is first

registered. The system automatically places all exams in

this status upon registration.

Ward The hospital location where an inpatient resides. In

VISTA, wards are a subset of the Hospital Location file

(#44).

Weighted work unit The number that results from multiplying the weight of a

procedure's AMIS code with the procedure's AMIS weight multiplier for that AMIS code. If a procedure has more than one AMIS code, the multiplication is done for

each and the results are summed.

Workload credit A general term that can refer to either the CPT type of

workload credit that is used in the VA to calculate reimbursement to medical centers for work done, or the AMIS crediting used by the AMIS workload system.

XVIII. Index

Ask 'Camera/Equip/Rm'V-3 Ask Drawn Dt/Time & Person (Radiopharmaceutical). II-6, VII-16 Ask Exam Status TimeV-4 Abnormal Alert......VII-6 Ask for Camera/Equip/Rm......VII-15 Abnormal Exam Report.....VII-6, XIV-5 Abnormal Report III-2, III-5 Ask for Film Data......VII-14 Access codeV-4 Ask for Interpreting Physician......VII-14 Access Erroneous Reports III-3 Ask for Procedure VII-14 Activation of exam statusVII-19 Ask for Technologist......VII-14 Active Procedure List (Long)VIII-19, VIII-24 Ask for User Code......VII-15 exampleVIII-24 Ask Imaging Location II-4 Active Procedure List (Short)VIII-19, VIII-25 Ask 'Imaging Location'V-2 exampleVIII-25 Ask Lot No. (Radiopharm).....II-6, VII-16 Activity Drawn Ask Medication Admin Dt/Time & PersonII-5, VII-15 exam statusVII-16 Ask Medications & DosagesII-5, VII-15 Activity Drawn Required (Radiopharm)......II-6, VII-13 case edits/status tracking......XII-4 Activity logsV-4 Ask Radiopharmaceutical Add Exams to Last Visit...... III-4, VII-64, XI-5 case edits/status tracking......XII-2 Add/Remove Report from Batch III-4 field......VII-13, VII-15 Admin Dt/Time/Person Required (Radiopharm)II-6, VII-13 Ask Radiopharmaceuticals and Dosages II-5, II-6, VII-15 Adverse reaction trackingVII-2 Ask Radiopharms and Dosages exam statusVII-10 case edits/status tracking......XII-3 Ask Route/Site of Admin (Radiopharm)II-6, VII-16 stat request......II-4 Ask Volume/Form (Radiopharm)......II-6, VII-16 Allergy Auto E-Mail to Req. Phys?.....II-3, V-7 contrast media.....VII-2 Auto User Code FilingV-4 Allow Batching of Reports V-5 Allow Copying of Reports......V-5 Allow E-Sig on COTS HL7 Rpts......V-7 B Allow Released/Not Verified.....II-3, V-5, V-25 Barcoded Procedure ListVIII-19, VIII-27 Allow Rpts on Cancelled Cases V-6 Allow Standard Reports......V-4, VII-85 example......VIII-27 Allow Verifying by Residents V-6, XIV-1 Batch namesXV-9 Allow Verifying of Others V-6, XIV-1 Batching reportsV-5 Alpha Listing of Active Procedures....... VIII-19, VIII-26 Bilateral VII-81 exampleVIII-26 effect on workload reportsXIV-3 Alphabetical listing of active procedures......VIII-26 AMISVII-2 Broad procedures VII-65, XV-2 categoryXV-13 suggested use ofXV-3 code(s).....XIII-1 contrast media use determined by......VII-2 crediting failure.....XIII-3 credit indicator.....VII-81 Bulletins II-11 VACOVII-45 weight.....XIV-3 \mathbf{C} weight.....VII-45 weight multiplierXIV-3 Called for ExamVII-8 workload reportXIV-3 Camera/Equip/Rm.....V-24 AMIS Code......VII-70, XV-13 required.....V-3 AMIS Code Dump by Patient III-5, XV-13 worksheetX-3 AMIS ReportIII-5, VII-45, VII-65, XV-13 Camera/Equip/Rm DescriptionV-19 AMIS reportsVII-33 Camera/Equip/Rm Entry/EditV-1, V-19, V-24 AMIS Weight MultiplierVII-70 example.....V-19 AMIS Workload Report Camera/Equip/Rm ReportXIV-2 exam statusVII-17 Camera/Equip/Rm Required VII-12 Appears on Status Tracking Camera/Equipment/Room.....V-19, V-21 fieldVII-9 Cameras/equip/rooms used by this location V-24 Ask Activity Drawn (Radiopharmaceutical)..... II-6, VII-16 Can an Exam be Cancelled while in this Status...... VII-10 Ask Admin Dt/Time & Person (Radiopharmaceutical) .II-6, Cancel a Request III-3, XV-8

VII-16

Cancel an Exam III-4	CPT
CancelledVII-8, XV-11	nationally activeVIII-3
Cancelled casesV-6	CPT codes
Cancelled examsVII-16	invalid/missingXIII-
Case editXII-1, XV-3	Create a BatchIII-
Case Edits and Status TrackingII-5	Create OE/RR ProtocolVII-5
Case No. Exam Edit III-4, VII-2, VII-7, VII-33, VII-64,	Create OE/RR Protocol from Common Procedure VII-53
XIII-2	VII-55, VII-60
Case number	exampleVII-6
exceptionally highXV-4	Credit clinic stop codesVIII-3
processing to completionXV-4	Credit failureV-
reuse/recycling ofXV-4	Credit MethodV-2
same on two active casesXV-4	Interpretation OnlyV-2
Category of exam	No CreditV-2
effect on workload reportsXIV-3	RegularV-2
CDR reportXV-14	Technical Component OnlyV-2
Check list	CreditingV-2
Cine filmVII-33	required dataXIII-
Classification	workloadVII-45, VII-81, VII-8
worksheet X-7	workloadXIII-
Classification Enter/Edit II-11, III-1, VI-1, VI-2, XIV-1	Crediting workloadVII-8
exampleVI-2	110
Classification of userXIV-1	D
Clerical ListVI-1, VI-6	D
Clerk	Daily Log ReportIII-
Clinic ReportXIV-2, XV-13	
Clinic Stop fileXIII-1	Daily Management ReportsXV- Daily ReportIII-
Clinic Workload Report	
exam statusVII-17	Data Inconsistency Report for Exam Statuses VII-2
Clinical history	Date/Time of dose administration
for broad procedure requestXV-3	exam status
Clinical History Message	Date/Time of measure
Common Data for PrintsetsXI-1	exam status
Common procedureXV-3	Default (Flash Card, Jacket Label, Exam Label) Format . V
imaging type, byVII-54	23
Common Procedure Enter/EditVII-53, VII-54	Default Exam Label Format
exampleVII-56	Default Flash Card Format
Common Procedure ListVII-54, VII-64	Default FormVII-6
Complete examsVII-16	Default Jacket Label Format
Complete status	Default Med Dose
ComplicationXV-13	Default MedicationII-8, VII-6
Complication Report III-2, III-5, VII-3	associated with a procedure
Complication Type Entry/EditVII-1, VII-2, VIII-2	Default Next Status
exampleVII-2	Default Next Status for Exam
Complication Type ListVII-2, VIII-1, VIII-2	Default Radiopharmaceutical II-
	associated with a procedure
example	formII-
Complication types	high doseII-
Contract radiologistsXV-11	low dose II-
Contract/Sharing	routeII-
Contract/Sharing Agreements fileXIV-3	usual doseII-
Contract/Sharing Source	Default Report Footer Format
Contracting out XV-14	Default Report Header Format
Contrast Medium AllergyXV-13	Default ResponsesXII-
Contrast medium reaction	Delete a ReportIII-
Contrast Medium Reaction	Delete Printed BatchesIII-
Contrast reaction	Delete Printed Batches by DateIII-
correction of erroneousXV-13	Deletion of stuffed dataXII-
Contrast ReactionVII-2	Delinquent Status Report III-5, VII-1
Contrast Reaction Message	exam statusVII-1
Cost Distribution ReportXV-13	Descendent proceduresVII-66, VII-72, XIII-
Cost of Procedure Enter/Edit	DetailedVII-6
CPRS Interface for PrintsetsXI-5	Detailed procedure entry

Detailed Procedure Required V3, VII-12	exampleVII-74	DX202 Non-Imaging Agents Radiopharmaceuticals II-1, II-
Detailed Procedure Workload Report		10
Exam status		
Exam status		E .
Stuffied /default responses		£
Detailed Request Display.	•	Edit Exam by Patient III-4, VII-2, VII-7, XIII-2
Device Specifications for Imaging Locations V-14 Diff Route of Administration V-11-69 DIff Site of Administration V-11-61 Diagnostic Code Enter/Edit V-11-1 Diagnostic Code Entry by Case No III-1 V-11-2 Diagnostic Code Required V-11-1 V-11-2 Diagnostic Code Required V-11-2 Diagnostic Code Required V-11-3 Diagnostic Code Required V-11-4 Diagnostic Code Required V-11-4 Diagnostic Code Required V-11-4 Diagnostic Code Required V-11-4 Diagnostic Code Required V-11-5 Diagnostic C		
Dit Rite of Administration		
Dit Site of Administration		
Diagnostic Code Enter/Edit VII-1, VII-3		
Diagnostic Code Enter/Edit VII-1, VII-3, vII-1, vIII-3 inactive II-1, vIII-3 vIII-1, vIII-3 vIII-1, vIII-3 vIII-1, vIII-3		
default		
Diagnostic Code Entry by Case No. III-4, XIII-2 Exam labels V-22, V-23, Diagnostic Code Entry by Case No. III-4, XIII-2 Exam profile (selected sort) III-3 Exam Room Report XV-12 Exam Room Report SV-12 Exam Retistant Status Flate VII-14 Cancelled VII-7, VII-8, VII-9 Complete VII-14 Cancelled VII-7, VII-8, VII-9 Complete VII-16 Cancelled VII-7, VII-15, VII-15, VII-15, VII-16 VII-17 VI		
Diagnostic Code Entry by Case No. III-4, XIII-2		
Diagnostic Code List VIII-1, VIII-3		,
Exam Som Report XV-12		
Diagnostic Code Required		
Diagnostic codes		-
abnormal results VII-4 triggering alerts VII-5 VII-6 triggering alerts VII-7 VII-8 VII-9 triggering alerts VII-10 triggering alerts VII-10 triggering alerts VII-10 triggering alerts VII-10 triggering alerts VII-11 triggering alerts VII-12 triggering alerts VII-13 triggering alerts VII-13 triggering alerts VII-13 triggering alerts VII-14 triggering alerts VII-15 trigge		
System, include in. VII-4 VII-2 VII-2, VII-8, VII-9 VII-7 VII-8, VII-9 VII-7, VII-8, VII-9 VII-19 VII-7, VII-8, VII-9 VII-10, VIII-10, VII-10, VII-10, VII-10, VII-10, VII-10, VII-10, VII-10, VII		
Complete		
Diagnostic Codes file	•	
Inactive field. VII-4 Display a Rad/Nuc Med Report III-3 Display Common Procedure List VII-51, VII-55, VII-51, VII-55, VII-51 example VII-61 Display Educational Description III-8 Display Educational Description When Ordered VII-73 Display Patient Demographics III-3 Display Patient Demographics III-3 Display Rad/Nuc Med Procedure Information III-3 Display Rad/Nuc Med Procedure Information III-3 Distribution queue XV-9 multiple active XVII-83 Reports Distribution Edit option VV-9 multiple active XVII-83 Exam Status Display Tracking options III-2 Distribution Queue Menu III-3 Division Queue Menu III-3 Division Parameter Set-up.I-4, III-10, V-1, V-2, V-28, XIV-1 changes to V-9 changes, when they take effect V-9 exam entry/edit, affecting XV-9 changes, when they take effect V-9 exam entry/edit, affecting V-9 worksheet X-1 Dosage ticket XVII-18 Dosage ticket XVII-19 Downtime III-2 Draft Report III-4 Draft Report III-5 Draft Report III-6 Draft Report III-7 Draft Report III-8 Draft Report III-9 Draft Report III-9 Draft Report III-1		
Display a Rad/Nuc Med Report	8	
Display Common Procedure List VII-53, VII-51 Display Educational Description		
example VII-15 Display Educational Description III-8 Display Educational Description When Ordered VII-73 Display Patient Demographics III-3 Display Rad/Nuc Med Procedure Information III-3 Disribution queue XV-9 multiple active VII-83 Reports Distribution Edit option VII-82 Distribution Queue XVII-82 Exam Status Display III-4, XV-2 Exam Status Display III-4, XV-2 Exam Status Display III-4, XV-2 Exam Status Tracking VII-16 Exam status Tracking VII-17 Division III-1, III-1, V-2, V-22 Division Parameter Set-up .I-4, III-10, V-1, V-2, V-28, XIV-1 1 affecting case edits/status tracking XII-1 Division parameters I-1 changes to V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam entry/edit affecting V-19 worksheet X-1 Dosage ticket VIII-18 Doubling workload credit XV-12 Downtime III-2 Downtime III-3 Division V-27, XIII-1, XIII-2 workload crediting VII-18 Doubling workload credit XV-12 Downtime III-2 Downtime III-3 Division V-27, XIII-1, XIII-2 workload crediting V-27, XIII-1, III-10 Division V-10 Division III-3 Division parameters I-1 III-3 Division parameters I-1 changes to V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam status VIII-18 Examination Status Entry/Edit II-5, II-6, II-10, V-4, VII-1 VII-7, VIII-5, XV-4 affecting case edits/status tracking XII-1 Sexamination Status bear vivial to the control of the	Display Common Procedure List VII-53, VII-61	• •
Display Educational Description		
Display Educational Description When Ordered VII-73 Display Patient Demographics. III-3 Display RadivNuc Med Procedure Information III-3 Display RadivNuc Med Procedure Information III-3 Distribution queue XV-9 multiple active VII-83 Reports Distribution Edit option VII-82 Distribution Queue Report Distribution Lists VIII-42 report verification triggering entry into XV-9 Distribution Queue Menu III-3 Division — I-4, III-1, V-2, V-22 Division Parameter Set-up.I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking XII-1 changes to V-9 exam entry/edit, affecting V-9 exam status Subjealy III-4, XV-2 Exam Status Display IIII-4, XV-2 Exam Status Display I		
Display Patient Demographics. III-3 Display Rad/Nuc Med Procedure Information. III-3 Distribution queue		•
Display Rad/Nuc Med Procedure Information III-3 Distribution queue		
Distribution queue XV9- multiple active VIII-82 Reports Distribution Edit option VIII-82 Distribution Queue Report Distribution Lists VIII-42 report verification triggering entry into XV-9 Distribution Queue Menu III-3 Division — I-4, III-1, V-2, V-22 Division Parameter Set-up.I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking XII-1 Division parameters — I-1 changes to V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam status VIII-8 Dosage ticket VIII-8 exam status VIII-8 Dosage ticket VIII-8 Downtime III-2 Downtime III-2 Downtime III-2 Downtime III-2 Dox200 Radiopharmaceuticals, Diagnostic III-1, III-10 DX200 Radiopharmaceuticals, Diagnostic III-1, III-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals III- Division Queue Report Distribution Edit option VIII-8 Exam Status Display IIII-4, V-2, V-4, V-5, Exam Status Display/Tracking options IIII-4, V-4, V-5 Exam Status Display IIII-4, V-2, V-2 Exam Status Display IIII-4, V-2, V-3 Exam Status Display IIII-4, VIII-6 Exam Status Display IIII-10 Exam Status III-1, III-6 Exam Status III-4, III-6 Exam Status III-1, III-6 Exam Status III-4, I		
multiple active	Distribution queueXV-9	
Reports Distribution Edit option. VII-82 Distribution Queue Report Distribution Lists VIII-42 report verification triggering entry into XV-9 Distribution Queue Menu III-3 Division — .I-4, III-1, V-2, V-22 Division Parameter Set-up. I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking XII-1 Division parameters — .I-1 changes to V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam status Division Status Entry/Edit .II-5, II-6, II-10, V-4, VII-1, VIII-5, XV-4 Examination Status Entry/Edit .II-5, II-6, II-10, V-4, VII-1, VIII-5, XV-4 affecting case edits/status tracking XII-1 Division parameters V-9 exam entry/edit, affecting V-9 exam entry/edit, affecting V-9 exam status VII-18 exam status Status Display/Tracking options III-2 Exam Status Tracking VIII-1 Exam Status Sactivation and parameters 1-1-2 Exam status Tracking VIII-5 Examination Status Sta	multiple activeVII-83	
Distribution Queue Report Distribution Lists VIII-42 report verification triggering entry into XV-9 Distribution Queue Menu III-3 Division — I-4, III-1, V-2, V-22 Division Parameter Set-up I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking XII-1 Changes to V-9 changes when they take effect V-9 exam entry/edit, affecting V-3 example V-9 worksheet X-1 Dosage ticket VIII-8 Dosage ticket VIII-8 Doubling workload crediti XV-1 Downtime III-2 Draft Report VIII-9 Day On Radiopharmaceuticals, Diagnostic III-1, II-10 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals III-1 DX201 Imaging Agents (in vivo), Radiopharmaceuticals III-1 DIIII-3 Exam Status Tracking VIII-6 Exam statuss Tracking III-2 Examination Status Tracking III-2 Examination Status Is Examination Status III-3 Examination Status Entry/Edit .II-5, II-6, II-10, V-4, VII-1, VIII-5 Examination Status Entry/Edit .II-5, II-6, II-10, V-4, VII-1, VIII-5 Examination Status Intervent III-4 Examination Status Intervent III-4 Examination Status Intervent III-4 Examination Status Intervent III-5 Examinat	Reports Distribution Edit optionVII-82	
report verification triggering entry into XV-9 Distribution Queue Menu III-3 Division — I-4, III-1, V-2, V-22 Division Parameter Set-up I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking XII-1 Changes to V-9 changes, when they take effect V-9 exam entry/edit, affecting V-3 example V-9 worksheet X-1 Dosage ticket VII-18 Dosage ticket VII-18 Doubling workload credit XV-12 Downtime II-2 Downtime II-2 Downtime Arabic Argonic Medical Computing Will-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-19 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- Examination Status istaus Sustaus V-3 worksheet X-2- worksheet X-3 worksheet X-3 worksheet X-4 Examination Status Entry/Edit .II-5, II-6, II-10, V-4, VII-1, VII-5 worksheet X-1- VII-7, VIII-5, XV-4 affecting sate edits/status tracking XIII-1 Examination Status Entry/Edit III-5, II-6, II-10, V-4, VII-1, VII-5 worksheet X-1- VII-7, VIII-5, XV-4 affecting sate edits/status tracking XIII-1 Examination Status Entry/Edit III-5, II-6, II-10, V-4, VII-1, VII-5 worksheet X-3 affecting case edits/status tracking XII-1 Examination Status Entry/Edit III-5, II-6, II-10, V-4, VII-1 Examination Status Entry/Edit III-5, III-6 Examination Status Entry/Edit III-5, III-6 Examination Status Entry/Edit III-5, III-6 Examination Status Entry/E		
Distribution Queue Menu III-3 Division — I-4, III-1, V-2, V-22 Division Parameter Set-up. I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking XII-1 Changes to V-9 changes, when they take effect V-9 exam entry/edit, affecting V-3 example V-9 worksheet X-1 Dosage ticket VII-18 Dosage ticket VII-18 Doubling workload credit XV-12 Downtime III-2 Downtime III-2 Downtime III-3 Division Parameter Set-up. I-4, II-10, V-1, V-2, V-28, XIV-1 Examination Status Sentry/EditII-5, II-6, II-10, V-4, VII-1, VII-7, VIII-5, XV-4 affecting case edits/status tracking XII-1 suggested entries VIII-19 Examination Status Entry/EditII-5, II-6, II-10, V-4, VII-1, VII-7, VIII-5, XV-4 affecting case edits/status tracking XIII-1 suggested entries VIII-19 Examination Status List VIII-1, VIII-5 Examination Status List VIII-1, VIII-5 Examination Status for Nuclear Medicine example VIII-12 Film Entry Required VIII-10 Example VIII-10 Examination Status Status tracking XIII-1 Examination Status Sentry/EditII-5, II-6, II-10, V-4, VII-1, VIII-10 Examination Status Sentry/EditII-5, II-6, II-10, V-4, VIII-1, VIII-10 Examination Status Entry/EditIII-5, II-6, II-10, V-4, VIII-1, VIII-10 Examination Status Entry/EditIII-5, II-6, III-10, V-4, VIII-1, VIII-10 Examination Status Entry/EditIII-5, II-6, II-10, V-4, VIII-1, VIII-10 Examination Status Entry/EditIII-5, II-6, II-10, V-4, VIII-1, VIII-10 Examination Status Entry/EditIII-5, II-6, II-10, V-4, VIII-1 Examination Status Entry/EditIII-5, II-6, II-10, V-4, VIII-1 Examination Status Entry/EditIII-5, II-6, II-10, V-4, VIII-1 Examination Status Entry/EditIII-5, II-6, III-10, V-4, VIII-1 Examination Status Entry/EditIII-5, II-6, III-10, V-4, VIII-1 Examination Status Entry/EditIII-5, II-6, III-10, V-4, VIII-1 Examination Status Entry/EditIII-5, III-6, III-10, VIII-10, VIII-10 Examination Status Entry/EditIII-5, II-6, III-10, VIII-10 Examination Status Entry/EditIII-5, III-6, III-10, VIII-10 Examination Status Entry/EditIII-5, III-6,		Exam statuses activation and parametersI-2
Division ————————————————————————————————————		Examination StatisticsIII-5
Division Parameter Set-up .I-4, II-10, V-1, V-2, V-28, XIV-1 affecting case edits/status tracking. XII-1 Division parameters		
affecting case edits/status tracking		parametersV-3
affecting case edits/status trackingXII-1 Division parameters	Division Parameter Set-up .I-4, II-10, V-1, V-2, V-28, XIV-	
Division parameters I-I changes to V-9 changes, when they take effect V-9 exam entry/edit, affecting V-3 example V-9 worksheet X-1 Dosage ticket VIII-18 exam status VIII-18 exam status VIII-18 Doubling workload credit XV-12 Downtime III-2 Draft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) III-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VIII-8 Dx200 Radiopharmaceuticals, Diagnostic III-1, III-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II-1, III-10 Dx201 Imaging Agents (in vivo), Radiopharmaceuticals II-1, III-10 IIII-10 III-10 III-10 IIII-10 IIII-	1	
changes to V-9 changes, when they take effect V-9 exam entry/edit, affecting V-3 example V-9 worksheet X-1 Dosage ticket VIII-18 Doubling workload credit XV-12 Downtime III-4 Drawn Dt/Time and Person Required (Radiopharm) III-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VIII-8 Dx200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- d affecting case edits/status tracking XIII-1 suggested entries VIII-19 Examination Status List VIII-19 Examination Status List VIII-1, VIII-5 example VIII-5 Examination status parameters V-3, V-6 Examination status parameters V-3, V-6 Examination status parameters V-3, V-6 Examination status for Nuclear Medicine example VIII-12 Film Entry Required VIII-12 Film Sizes List VIII-1, VIII-10 example VIII-10 Film Type Entry/Edit VIII-1, VIII-33, VIII-10		
changes, when they take effect V-9 exam entry/edit, affecting V-3 example V-9 worksheet X-1 Dosage ticket VIII-18 exam status VIII-18 Doubling workload credit XV-12 Downtime III-2 Doraft Report IIII-4 Drawn Dt/Time and Person Required (Radiopharm)II-6, VIII-13 DSS ID. V-27, XIII-1, XIII-2 workload crediting. VIII-8 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic III-1, III-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals III-1, III-10 DX201 Imaging Agents (in vivo), Radio		
exam entry/edit, affecting V-3 example V-9 example V-9 worksheet X-1 Dosage ticket VII-18 exam status VII-18 Doubling workload credit XV-12 Downtime III-2 Draft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) II-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-87 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic II-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- L H 10 Examination Status List VIII-1, VIII-5 example V-3, VII-7 Examination Status List VIII-1, VIII-5 example VIII-5 Examination Status List VIII-1, VIII-5 example VIII-7 Examination Status List VIII-1, VIII-5 example VIII-8 Examination Status List VIII-1, VIII-5 example VIII-8 Examination Status List VIII-1, VIII-5 example VIII-8 Examination Status List VIII-1 Examination Status List VIII-1, VIII-5 example VIII-8 Examination Status List VIII-1 Examination Status List VIII-1, VIII-5 example VIII-5 Examination Status List VIII-1, VIII-5 Examination Status List VIII-1 Examination Status List VIII-1, VIII-5 example VII-3 Film Entry Required VIII-1, VIII-10 example VIII-1, VIII-10 Film Type VII-1, VIII-33, VIII-10 Film Type Entry/Edit VII-1, VII-33, VIII-10		
example V-9 worksheet X-1 Dosage ticket VIII-18 exam status VIII-18 Doubling workload credit XV-12 Downtime III-2 Draft Report IIII-4 Drawn Dt/Time and Person Required (Radiopharm) III-6, VIII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VIII-8 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals		
worksheet X-1 Dosage ticket VII-18 exam status VII-18 Doubling workload credit XV-12 Downtime III-2 Draft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) II-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-8 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II-1, II-10 DX201 I		
Dosage ticket VII-18 exam status VII-18 Doubling workload credit XV-12 Downtime III-2 Doraft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) II-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-87 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- I M 10 Examination status parameters V-3, V-6 Examined Status for Nuclear Medicine example example VIII-23 Film Entry Required VIII-12 Film Sizes List VIII-1, VIII-10 example VIII-10 Film Type Intry/Edit VII-1, VII-33, VIII-10	•	
exam status VII-18 Doubling workload credit XV-12 Downtime III-2 Doraft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) III-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-87 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- Examined VII-8 Examined Status for Nuclear Medicine example example VII-8 Film Entry Required VII-12 Film Sizes List VII-11 Examined Status for Nuclear Medicine example example VII-12 Film Sizes List VIII-11 Film Type VII-11 Film Type Entry/Edit VII-12 VII-13 VII-13 VII-13 VII-14 Film Type Entry/Edit VII-15 VII-16 VII-17 Film Type Entry/Edit VII-17 VII-18 Examined Status for Nuclear Medicine example example VII-12 Film Sizes List VIII-13 VIII-10 Film Type Entry/Edit VII-13 VII-14 VII-15 VII-15 VII-15 VII-16 VII-16 VII-17 Film Type Entry/Edit VII-18 VII-19 VII-19		
Doubling workload credit XV-12 Downtime III-2 Draft Report IIII-4 Drawn Dt/Time and Person Required (Radiopharm) III-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-87 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II-1, II-10 DX201 Imaging Agents (in vivo)		
Downtime III-2 example VII-23 Draft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) III-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-87 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II-1, II-10 The state of the sample VII-12 Film Entry Required VIII-10 Example VIII-10 Film Sizes List VIII-1, VIII-10 Film Type VIII-10 Film Type Intry/Edit VIII-1, VII-33, VIII-10		
Draft Report III-4 Drawn Dt/Time and Person Required (Radiopharm) II-6, VII-13 DSS ID V-27, XIII-1, XIII-2 workload crediting VII-87 Duplicate Flash Card III-4 DX200 Radiopharmaceuticals, Diagnostic III-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- Telm Type Intry/Edit VII-1, VII-33, VIII-10		
Drawn Dt/Time and Person Required (Radiopharm)II-6, VII-13 DSS ID		Crampic
DSS ID	Drawn Dt/Time and Person Required (Radiopharm)II-6,	F
workload crediting		Film Entry Required VII 12
Duplicate Flash Card		
DX200 Radiopharmaceuticals, DiagnosticII-1, II-10 DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- Tilm TypeVII-1, VII-33, VIII-10		
DX201 Imaging Agents (in vivo), Radiopharmaceuticals II- Film Type Entry/EditVII-1, VII-33, VIII-10		
	1, II-10	

Film typesI-2	Imaging locations	
Film Usage ReportVII-33, VII-34, XV-13	data corruption due to deletion of	XV-1
Film Usage Workload Report	new	
exam statusVII-17	Imaging locations associated with division	
Films Reporting MenuVII-85	Imaging Locations file	
Films reporting parameters	Imaging stop codes	
Flash card format	Imaging Stop Codes file	
defaultVII-40	Imaging type I-4, III-1,	
Flash card formatter	activating	
example of set-upVII-41	activation of	
Flash card, jacket labels, exam label parametersV-22	associating a procedure with	
Flash cardsI-4, V-7, V-9, V-22, V-23	common procedure, by	
Footers	data corruption due to deletion of	
results reportsVII-36	deactivating of	
Frequently asked questionsXV-1	displays, effect on	
Functional Area Workload Reports III-5	edits, effect on	
Future RegistrationV-4	exam status	
1 uture Registration V-4	examination status	
~	modifier screening by	
${f G}$	Nuclear Medicine	
Generate Exam Alert for Requesting PhysicianVII-10	procedures associated with	
Generate Examined HL7 MessageII-5, VII-10	reports, effect on	
	separation of workload	
H	sign-on definition of	
**	Imaging Type file	
Header/Footer formatter	Imaging type of this location	
example of set-upVII-41	Implementation	
Headers	Implementation check list	
results reportsVII-36	Impression	
Health screening stop codesXIII-1	common problem in data entry of	
Health Summary Type fileVII-66	Impression Required	
Health Summary with RequestVII-66	Impression Required on Reports	
High Adult DoseVII-69	Impressions	
HIST Display Rad/NM Procedure Contrast Media History	Inactivate	VII-54
VIII-20	Inactivation Date	
HL7 InterfaceXI-3	associated with a procedure	
HL7 messageVII-10	Inactivation of Imaging Location	V-25
HoldXV-11	Inactive	
Hold a Request III-3, XV-7, XV-8	Diagnostic Code Enter/Edit	II-5
Hospital Location file V-8, V-13, V-22, V-29, XVII-3,	Inactive Procedure List (long)	VIII-30
XVII-6, XVII-12, XVII-14	example	VIII-30, VIII-41
Occasion of Service entryXIII-1	Inactive Procedure List (Long)	
Pattern Date fieldV-8	Inactive proceduresVI	
Type field	Incomplete cases	
	Incomplete Exam Report	
Hospital Location Type, prompt	Indicate No Purge	
How Many (Flash Cards, Jacket Labels, Exam Labels) per	Indicate No Purging of an Exam/Report	
VisitV-23	Inquire to File Entries	
How Many Jacket Lbls Per Visit	Inquire to File Options	
	Installation check	
I		
	virgin install	
Imaging location II-11, V-2, V-8, V-19, V-22	virgin install	
accessII-11	Institution file	
assignment to usersII-11	Interpreting physician	
Imaging LocationI-4, VII-55	Interpreting Resident	
assignment to usersVI-2	Interpreting Resident List	
Imaging location access	Interpreting residents	
when changes take effectVI-2	Interpreting Staff	
Imaging Location AccessVI-2	Interpreting Staff List	
Imaging location parametersI-1	Invalid CPT/Stop Code ListVIII-19	
Imaging location, inactivation ofV-25	example	VIII-31
6	Invalid locations	V-9

IRM MenuI-1	Low Adult DoseVII-68	
J	\mathbf{M}	
Jacket label format	Maintenance Files Print MenuI-2, II-11, III-3, VII-1, VIII-	
defaultVII-40	1, VIII-18	
Jacket labelsI-4, V-7, V-9, V-22, V-23	Major AMIS Code Entry/Edit VII-1, VII-45, VIII-12	
Jacket Labels	exampleVII-45	
THE STATE OF THE S	Major AMIS Code ListVIII-1, VIII-12	
-	exampleVIII-12	
${f L}$	Major AMIS codes (use only if authorized by VACO)I-2	
Lab workXV-11	Major AMIS Codes (use only if authorized by VACO)1-2 Major AMIS Codes fileVII-70	
Label and Report header/footer formatsI-2	Major Rad/Nuc Med AMIS Code file	
Label formatter	Management report criteria exam statusVII-16	
example of set-up		
Label Print Fields file VII-40	status tracking	
Label/Header/Footer Format ListVIII-11	Margin, left	
exampleVIII-11	Mass Override Exam StatusIII-4	
Label/header/footer formatsI-4	Medication Logic	
Label/Header/Footer FormatterII-7, V-23, V-30, V-31,	case edits/status tracking	
VII-1, VII-36, VIII-11	Medications	
exampleVII-41	associated with a procedure	
new fieldsII-7	exam statusVII-15	
Label/Header/Footer ListVIII-1	MessageVII-72	
Labels, Headers, FootersXI-4	Miscellaneous Division ParametersV-7	
Lbl/Hdr/Ftr formatsXVII-8	Missing or invalid CPTsVIII-31	
flash cardXVII-6	Modifier ListVIII-1, VIII-13	
footerXVII-6	exampleVIII-13	
formatXVII-6	ModifiersVII-55	
headerXVII-6	effect on workload reportsXIV-3	
List Exams with Inactive/Invalid Statuses III-3	Multiple AMIS codes	
List of active procedures	effect on workload reportsXIV-5	
shortVIII-25	Multipliers	
List of Cameras/Equip/RmsV-1, V-21	effect on workload reportsXIV-3	
List of inactive proceduresVIII-30		
List of Inactive Procedures (Short)VIII-19, VIII-33	N	
exampleVIII-33	14	
List of invalid CPT or stop codesVIII-31	No ComplicationVII-2	
	No showsXV-7	
List of series of procedures	Non-Common Data for PrintsetsXI-1	
List Reports in a Batch	Non-count stop codesXIII-1	
Listing of all active proceduresVIII-24	NPOXV-11	
Location	Nuclear Medicine List MenuVIII-1, VIII-14	
Location Parameter ListV-1, V-34	Nuclear Medicine List Menu (all options)III-3	
exampleV-34	Nuclear Medicine parameters	
Location Parameter Set-up I-4, II-3, II-10, V-1, V-22, VII-	Nuclear Medicine setup	
40, XV-2	worksheetX-17	
camera/equip/rm, associatedV-19	Nuclear Medicine Setup Menu I-3, I-7, II-8, III-3, VII-1,	
exampleV-28	VII-47, VIII-14	
Location Parameter Set-UpV-9	routesII-8	
Location parametersV-9	sites II-8	
when changes take effectV-28	Sites11-0	
worksheetX-4		
Location Set-up ParametersVII-40	0	
Log of Scheduled Requests by Procedure III-4, XV-11		
Long Case Number BarcodeII-7	OE/RRVII-55, XV-3	
Lot (Radiopharmaceutical) Number Enter/Edit VII-47, VII-	using common procedures to create protocols VII-60	
48, VII-51	On-line Verifying of Reports III-4, V-6, XIII-2, XIV-1	
exampleVII-48	Operating room	
Lot (Radiopharmaceutical) Number ListVIII-14, VIII-15	effect on workload reportsXIV-3	
Lot No. Required (Radiopharm)II-6, VII-13	Order entryXV-12	
Lot Number Enter/EditII-8	Order Entry ParametersV-2, V-23	

Order Entry Procedure Display MenuVII-1, VII-53	exampleV-15
Order entry proceduresI-2	Print Dosage Ticket when Exam Reaches this Status II-6
Order Entry/Results ReportingV-2	Print DX Codes in Report?V-24
Order in Sequence of Status ProgressionVII-9	Print File EntriesIII-3, V-8
Order numberVII-8	optionII-8
Order of the progressionVII-9	Print Flash Card for Each ExamV-7, V-9
Ordering	promptV-22, V-23
broad proceduresXV-3	Print Jacket Labels with Each VisitV-7, V-9, V-22
common proceduresXV-3	Print Rad/Nuc Med Requests by DateIII-3
OE/RRXV-3	Print Selected Requests by PatientIII-3
reschedulingXV-7	Print WorksheetsIII-3
sample workflowXV-7	Printer set-up V-23, V-28, V-34
user problem in selecting proceduresXV-3	Printing/batchingXV-9
Orders	Printsets
accounting for all pendingXV-2	Add Exams to Last VisitXI-5
Outside Films ProfileIII-3	cancelled casesXI-2
Outside Films Registry Menu III-3	common dataXI-1
Override a Single Exam Status to 'complete' III-4	CPRS interfaceXI-5
	editing dataXI-4
P	HL7 interfaceXI-3
•	labels/headers/footersXI-4
Parameters	non-common dataXI-1
changes, when they take effect V-9	RA MGR keyXI-2
when changes take effectV-28	single reportXI-1
Parent procedureVII-65	special character designationXI-2
separate reportsVII-66	VA Fileman prints/inquiriesXI-5
single reportVII-66	ProcedureV-2, V-3
Parent procedure entryVII-75	AMIS weight multiplierXIV-3
Parent Procedure ListVIII-19, VIII-38	broadV-3
Parent proceduresXIII-1	cost, entering ofVII-63
single report all descendentsII-8	detailedV-3
Parent/Descendent Exam and PrintsetsI-7, XI-1	done by locations of different imaging types VII-64
Patient Care EncounterXIII-2	seriesV-3
Patient chartsXV-8	Procedure message
PendingXV-11	for broad typeXV-3
Pending orders	Procedure Enter/Edit . II-1, II-3, II-8, VII-1, VII-62, VII-64,
accounting forXV-2	VII-79, VIII-31, XV-1, XV-4
Pending requestsXV-8	affecting case edits/status trackingXII-1
Pending/Hold Rad/Nuc Med Request Log III-3, XV-2, XV-	medication parameters II-8
11	parent procedure entry example VII-75
Person who administered the dose	radiopharmaceutical parametersII-8
exam statusVII-16	Procedure File ListingsVIII-1, VIII-18, XIII-1
Person who measured	Procedure Message Entry/EditVII-62, VII-79, XV-4
exam statusVII-16	exampleVII-79
Personnel classification	Procedure Message ListVIII-19, VIII-40
Personnel classification lists	exampleVIII-40
exampleVI-6	Procedure messagesI-2
Personnel Workload Reports III-5	Procedure modifier
Pharmacy Drug file	bilateralVII-81
radiopharmaceuticalsII-1, II-10	operating roomVII-81
VA ClassificationsII-1	portableVII-81
Physician ReportXIV-2, XV-12	screened by imaging typeVII-80
Physician Title on ReportsXIV-2	Procedure Modifier EntryVII-62, VII-80, VIII-13
Physician Workload Report	exampleVII-80
exam statusVII-18	Procedure modifiersI-2, VII-80
Population of Imaging Exam DataXII-5	Procedure setup
PortableVII-81, XV-13	worksheetX-12
effect on workload reportsXIV-3	Procedures
PrepsXV-11	broadXV-2
Primary Camera/Equip/Rm promptV-3	changing broad to detailed/seriesXV-3
Print a Batch of ReportsIII-3	changing imaging type ofXV-1
Print Division Parameter List	CPTXV-2

data corruption due to deletion of	XV-1	field	II-c
detailed	XV-2	Radiopharms/Dosages Required?	
edit	XV-2	case edits/status tracking	XII-3
imaging type of	XV-1	Reason for hold/DC of order	XV-8
message	XV-3	Reasons for cancellation	I-2
name edit	XV-2	Register an Exam	VII-64
originally distributed	XV-2	Register Patient for Exams	III ²
Procedures and their characteristics		Registration	
Profile of Rad/Nuc Med Exams		changing a broad to detailed or seri-	
Prompt for Meds		future dates	V-4
case edits/status tracking	XII-3	use of Pn	
Prompt for Radiopharm Rx		Released/not verified	
PTF Bedsection Report	XIV-2, XV-13	Released/Unverified	
PTF Bedsection Workload Report		Report	
exam status		Report Distribution Lists	
Purging		example	
Purging records	XV-9	Report distribution queues	
		Report Entered Required	VII-14
Q		Report Entry/Edit II-1, III-5	
· ·		Report footer	V-23
Quick order	VII-55	Report footer format	T.T.T. 47
		default	
R		Report header	V-23
		Report header format	VIII 40
RA ALLOC II-10, II-11, III-1, III-5, IV	-1, VI-6, XVII-12	default	
RA HL7 Voice Reporting Errors	IX-2	Report Left Margin	
RA MGR II-10, III-3, III-4, IV-1, V-		Report Parameters	
RA VERIFY II-10, III-4, III-		Report purge Report Right Margin	
Rad/NM Phys Approval Required		Report text	
Rad/Nuc Med HL7 voice Reporting Errors		Reports Distribution EditV	v-0 /II_1 VII_82 VIII_//
Rad/Nuc Med HL7 Voice Reporting Error		example	
Rad/Nuc Med Location Access		Request	
Rad/Nuc Med Personnel Menu I-1, II-11		Request an Exam	
Rad/Nuc Med Procedures file		Request Printer Name	
Rad/Nuc Med Procedures Name		Request status	
Rad/Nuc Med Total System Menu		Requesting physician	
Rad/Nuc Personnel Menu		Requests	,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,
Radiopharm Prescription & Prescribing Ph	•	Display of	XV-11
Logic		lab work	
Radiopharmaceutical Dosage Ticket Print	Collubi VII-18	NPO	XV-11
Radiopharmaceutical dose	II 2	preps	XV-11
High Adult DoseLow Adult Dose		scheduled	XV-11
Radiopharmaceutical Drug Classifications		status tracking log	XV-11
example		tracking	XV-11
Radiopharmaceutical functions	21 13	Required Flash Card Format	
exam status	VII-15	Required Flash Card Printer	
Radiopharmaceutical Logic		Rescheduling requests	
in case edits/status tracking	XII-2	Research	
Radiopharmaceutical Lot Source		Research patient	
field	VII-51	Resend Radiology HL7 Message	
Radiopharmaceuticals	I-7	Resident	
associated with a procedure	VII-67	Resident On-line Pre-Verification	
requirements		Resident Penart	
setup/definition of		Resident Report	
suppress asking/requiring		Resident Signature Name	11-
VA Drug Classifications	II-1, II-10	Resident Workload Report	1711 10
worksheet	X-13	exam status	V11-18
Radiopharmaceuticals Used?		Results reports batch printing	VVI
field		distribution queues	
Radiopharms & Dosages Required	II-6, VII-13	headers/footers	
Radiopharms/Dosages Required		11Cauci 5/100tc15	v11-30

print location ofXV-9	imaging location	III-2
queuing to printersXV-9	imaging type	III-2
title of physicianXIV-2	Incomplete Exam Report	III-5
Review RequiredVII-71	Indicate No Purging of an Exam/Report	III4
Route of Administration Enter/Edit . VII-47, VII-49, VII-50	Inquire to File Options	
exampleVII-49	Jacket Labels	
Route of Administration fileVII-47	List Exams with Inactive/Invalid Statuse	sIII-3
Route of Administration List . II-8, VII-47, VIII-14, VIII-16	List Reports in a Batch	
Route/Site of administration	Log of Scheduled Requests by Procedure	
exam statusVII-16	Maintenance Files Print Menu	
Route/Site Required (Radiopharm Administered) . II-6, VII-	Mass Override Exam Status	
13	Nuclear Medicine List Menu (all options	s)III-3
Routes	Nuclear Medicine Setup Menu (all optio	
radiopharmaceuticalII-8	On-line Verifying of Reports	
Rpharm Dose Warning MessageII-3, V-8	option allocation	
	Outside Films Profile	
C	Outside Films Registry	
\mathbf{S}	Override a Single Exam Status to 'compl	
Sample Contrast Media Edit History by Procedure Report	Pending/Hold Rad/Nuc Med Request Lo	
VIII-23	Personnel Workload Reports	
Schedule a Request III-3, XV-7, XV-8	Print a Batch of Reports	
ScheduledXV-11	Print File Entries	
Scheduling XV-11 Scheduling XIII-3	Print Rad/Nuc Med Requests by Date	
Scheduling packageXV-7	Print Selected Requests by Patient	
Scheduling requests	Print Worksheets	
reschedulingXV-7	Profile of Rad/Nuc Med Exams	
step-by-step exampleXV-7	Rad/Nuc Med Personnel Menu	
Screening Screening	Register Patient for Exams	III-4
Abnormal ReportIII-5	Report Entry/Edit	
Access Erroneous Reports III-3	report ownership	
Add Exams to Last Visit	reports	
Add/Remove Report from Batch III-4	Request an Exam	
AMIS Code Dump by Patient III-5	Resident On-line Pre-Verification	
AMIS Report III-5	Schedule a Request	
Cancel a Request III-3	Search File Entries	
Cancel an ExamIII-4	security keys	
Case No. Exam Edit III-4	security keys	
Complication Report III-5	Select Report to Print by Patient	
Create a Batch III-3	selected imaging type	
Daily Log ReportIII-5	Special Reports	
Delete a Report III-3	Status Tracking of Exams	
Delete Printed Batches	Switch Locations	
Delete Printed Batches by Date III-3	System Definition Menu	III-4
Delinquent Status Report	Test Label Printer	
Detailed Request Display III-4	Transcription Report	
Diagnostic Code Entry by Case No III-4	Unverified Reports	
Display a Rad/Nuc Med Report III-3	Unverify a Report	
Display Patient Demographics III-3	Update Exam Status	
Display Rad/Nuc Med Procedure Information III-3	Update Patient Record	
distribution queue III-3	Utility Files Maintenance Menu	
division III-2	Verify Batch	
Draft ReportIII-4	Verify Report Only	
Duplicate Flash Card III-4	View Exam by Case No	
Edit Exam by Patient III-4	Ward/Clinic Scheduled Request Log	
effects of RA ALLOC key III-5	Screening stop codes	
Enter Last Past Visit Before DHCP III-3	Search File Entries	
Exam Deletion III-3	Security key	
Exam Profile (selected sort)	status tracking	VII-9
Exam Status Display III-4	Security keys	
Examination Statistics III-5	RA ALLOC	
Functional Area Workload Reports III-5	RA MGR	
Hold a Request III-3	RA VERIFY	

used in options III-5	effect of outpatient statusXIII-2
Select Interpreting PhysicianXIV-1	error tracking/recovery using MAS mail group bulletins
Select Report to Print by Patient III-3	XIII-3
Select sign-on locationIV-1	options that triggerXIII-2
SeriesVII-65	Stop code crediting, error tracking/recovery using Rad/Nuc
Series of Procedures ListVIII-19, VIII-41	Med bulletinXIII-3
Service ReportXIV-2, XV-13	Stop codesXIII-
Service Workload Report	background errorsXIII-3
exam statusVII-17	crediting ofXIII-2
Seven Hundred series stop codesXIII-1	how they are passed to MASXIII-2
Sharing Agreement/ContractXIV-3	invalid/missingXIII-
Sharing Agreement/Contract Entry/EditVII-1, VII-84, VIII-	yearly updateVII-87
43	Stop Codes
exampleVII-84	assignment to proceduresXIII-1
Sharing Agreement/Contract ListVIII-1, VIII-43	Stopping tasked jobsXVI-
exampleVIII-43	Stuffing DataXII-
Sharing Agreement/Contract Report.VII-84, XIV-2, XV-13	Submit Request To promptV-2, V-25
Sharing agreements and contractsI-2	Supervisor Menu. V-1, VI-1, VII-1, VII-53, VIII-1, VIII-18
Sharing/Contract Workload Report	Suppress Radiopharm PromptII-8, VII-67
exam statusVII-17	Switch LocationsIII-4, V-10, V-28, XV-2
Should this Status Appear in Status TrackingVII-9	Synonym
Signature block title, use in results reportXIV-2	associated with a procedure
Sign-on imaging type III-2	System Definition MenuI-1, II-11, III-4, V-1, VII-
Single Report II-8, VII-66	
parent procedureVII-66	T
Single Report Setup for PrintsetsXI-1	1
Site of Administration Enter/EditVII-47, VII-50	Tasked jobs
exampleVII-50	stoppingXVI-
Site of Administration fileVII-47	TechnologistVI-2
Site of Administration List II-8, VII-47, VIII-14, VIII-17	Technologist ListVI-1, VI-0
Sites	Technologist ReportXIV-2, XV-13
radiopharmaceuticalII-8	Technologist RequiredVII-12
Special Character Designation for PrintsetsXI-2	Technologist Workload Report
Special Reports III-5, XIV-2	exam statusVII-17
StaffXIV-1	Test Label Printer III-4
Staff ReportXIV-2, XV-13	Time Limit for Future ExamsV-2
Staff Signature NameII-7	Title of interpreting physician on results reports XV-1
Staff Workload Report	Title of physician on results reportsXIV-2
exam statusVII-18	Track Exam Status ChangesV
Standard reports I-2, V-4, XVII-12	Track Request Status ChangesV-3
Standard Reports Entry/Edit V-5, VII-1, VII-85, VIII-44	Tracking requestsXV-1
exampleVII-86	TranscribedVII-8
Standard Reports PrintVIII-1, VIII-44	Transcription
exampleVIII-44	standard reportsVIII-44
Stat Request Alert RecipientsII-4, V-25	Transcription ReportIII-4, XIV-2
Status	Transcriptionist
examVII-7	Trouble shootingXIV-
orderVII-7	Type of ImagingV-26, VII-64
reportVII-7	Type of ProcedureVII-65
Status Change RequirementsXV-4	Type of Procedure fieldVII-65
Status parameter	••
route/siteVII-49	${f U}$
Status trackingVII-8, XII-1	U
Status TrackingXV-4	Unverified ReportsII-3, III-5, V-6
security keyVII-9	Unverify a ReportIII-2
Status tracking functionsVII-14	Update Exam Status III-4, XIII-2
Status Tracking of ExamsIII-4, VII-7, VII-33, XIII-2, XV-2	Update Patient RecordIII-4, XV-13
Stop codes	Urgent Request Alerts
VA HeadquartersVII-87	Urgent Request Alerts?V-25
Stop code crediting	User Key Needed to Move an Exam to this Status VII-9
effect of division onXIII-2	Usual DoseVII-69
effect of exam status overrideXIII-2	Utility Files Maintenance XV-
	Carry 1 nos manicolanec

Utility Files Maintenance Menu .I-2, I-7, II-11, III-4, VII-1,	Ward Workload Report	
VII-47, VII-53, VIII-1	exam status	
	Ward/Clinic Scheduled Request Log	gIII-4, XV-11
${f V}$	Warning on Rpts Not Yet Verified	V-6
Y	Wasted Film Report	VII-34
VA Fileman Prints/InquiriesXI-5	Wasted film sizes	VIII-10
VA HeadquartersXIII-1	Wasted Film, prompt	VII-34
Valid imaging stop codesI-2	Weight	VII-45, XIV-3
Valid Imaging Stop Codes Edit VII-1, VII-87, XIII-1	Weighted work unit	
exampleVII-87	Weighted work units	XIV-3, XIV-4
Valid Sites of Administration	Workload	
fieldVII-50	multiple AMIS codes	XIV-5
Vendor/Source (Radiopharmaceutical) Enter/EditVII-47,	Workload credit	
VII-51	Workload crediting	
exampleVII-51	doubling of	XV-12
Vendor/Source (Radiopharmaceutical) List VIII-14, VIII-18	effect of modifiers	XV-12
Verified Report RequiredVII-14	options that trigger	XIII-2
Verified reportsXV-9	Workload report	
Verify a reportXV-11	Workload reports	
Verify Batch III-4, XIII-2	effect of contract/sharing on	
Verify Report OnlyIII-5, XIII-2	effect of exam status on	XIV-3
Verifying others' reportsXV-11	effect of multiple residents	XIV-5
Verifying reportsV-6, XV-11	effect of multiple staff	XIV-5
Verifying Signature NameII-7	effect of multiple technologists	
View Exam by Case No III-4	effects of AMIS code on	
Voice Dictation Auto-PrintV-24	multipliers	XIV-3
Voice Recognition SystemsXIII-2	troubleshooting	
Volume/Form Required (Radiopharm)II-6, VII-14	WWU	
	WWUs	XIV-4
\mathbf{W}		
V V	X	
Waiting for examVII-8		
Ward ReportXIV-2, XV-13	XUSESIG	VI-2