



Update_2_0_92

Clinical Reminders

VA-BFA TEMPLATE UPDATES

Install Guide

September 2019

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION	1
INSTALL DETAILS	3
INSTALL EXAMPLE.....	3
POST INSTALL	6

Introduction

Changes to this version of the BFA treatment note include: adding the option to document Battlefield Auricular Acupressure (BAA), creating the ability to pull what type of acupuncture device was used to assess effectiveness of the different needle types including the effectiveness of acupressure devices. We have also added the category of device “ear magnet” as an option for device used. Several small formatting options were corrected/changed.

UPDATE_2_0_92 contains 1 Reminder Exchange entry:
UPDATE_2_0_92 VA-BFA TEMPLATE UPDATE

The exchange file contains the following components:

TIU TEMPLATE FIELD

VA-BFA R-L-BL
WORD PROCESSING SHORT
BFA SCALE HYPERLINK
BFA SUPPLEMENTAL QUESTIONS HYPERLINK

HEALTH FACTORS

VA-REMINDER UPDATES
VA-UPDATE_2_0_92
VA-BFA CATEGORY
VA-BFA PT AGREED BAA
VA-BFA WITH OTHER TREATMENT
VA-BFA ONLY
VA-BFA NEEDLE TYPE EAR MAGNET
VA-BFA NEEDLE TYPE OTHER
VA-BFA NEEDLE TYPE STANDARD NEEDLES
VA-BFA NEEDLE TYPE EAR SEED
VA-BFA NEEDLE TYPE EAR PRESS TACKS
VA-BFA NEEDLE TYPE SEMI PERM ASP
VA-BFA PAIN POST SCALE 10
VA-BFA PAIN POST SCALE 9
VA-BFA PAIN POST SCALE 8
VA-BFA PAIN POST SCALE 7
VA-BFA PAIN POST SCALE 6
VA-BFA PAIN POST SCALE 5
VA-BFA PAIN POST SCALE 4
VA-BFA PAIN POST SCALE 3
VA-BFA PAIN POST SCALE 2
VA-BFA PAIN POST SCALE 1
VA-BFA PAIN POST SCALE 0
VA-BFA PAIN SCALE 0

VA-BFA PAIN SCALE 10
VA-BFA PAIN SCALE 9
VA-BFA PAIN SCALE 8
VA-BFA PAIN SCALE 7
VA-BFA PAIN SCALE 6
VA-BFA PAIN SCALE 5
VA-BFA PAIN SCALE 4
VA-BFA PAIN SCALE 3
VA-BFA PAIN SCALE 2
VA-BFA PAIN SCALE 1
VA-BFA PAIN ASSESS STRESS 24H CATEGORY
VA-BFA PAIN ASSESS STRESS 10
VA-BFA PAIN ASSESS STRESS 9
VA-BFA PAIN ASSESS STRESS 8
VA-BFA PAIN ASSESS STRESS 7
VA-BFA PAIN ASSESS STRESS 6
VA-BFA PAIN ASSESS STRESS 5
VA-BFA PAIN ASSESS STRESS 4
VA-BFA PAIN ASSESS STRESS 3
VA-BFA PAIN ASSESS STRESS 2
VA-BFA PAIN ASSESS STRESS 1
VA-BFA PAIN ASSESS STRESS 0
VA-BFA PAIN AFFECTS MOOD 24H CATEGORY
VA-BFA PAIN ASSESS AFFECT MOOD 2
VA-BFA PAIN ASSESS AFFECT MOOD 1
VA-BFA PAIN ASSESS AFFECT MOOD 0
VA-BFA PAIN ASSESS AFFECT MOOD 10
VA-BFA PAIN ASSESS AFFECT MOOD 9
VA-BFA PAIN ASSESS AFFECT MOOD 8
VA-BFA PAIN ASSESS AFFECT MOOD 7
VA-BFA PAIN ASSESS AFFECT MOOD 6
VA-BFA PAIN ASSESS AFFECT MOOD 5
VA-BFA PAIN ASSESS AFFECT MOOD 4
VA-BFA PAIN ASSESS AFFECT MOOD 3
VA-BFA PAIN INTERFERE SLEEP 24H CATEGORY
VA-BFA PAIN ASSESS INTERFERE SLEEP 10
VA-BFA PAIN ASSESS INTERFERE SLEEP 9
VA-BFA PAIN ASSESS INTERFERE SLEEP 8
VA-BFA PAIN ASSESS INTERFERE SLEEP 7
VA-BFA PAIN ASSESS INTERFERE SLEEP 6
VA-BFA PAIN ASSESS INTERFERE SLEEP 5
VA-BFA PAIN ASSESS INTERFERE SLEEP 4
VA-BFA PAIN ASSESS INTERFERE SLEEP 3
VA-BFA PAIN ASSESS INTERFERE SLEEP 2
VA-BFA PAIN ASSESS INTERFERE SLEEP 1
VA-BFA PAIN ASSESS INTERFERE SLEEP 0

VA-BFA PAIN INTER/ACTIV. 24H CATEGORY
VA-BFA INTERFERE ACTIVITY 10
VA-BFA PAIN INTERFERE ACTIVITY 9
VA-BFA PAIN INTERFERE ACTIVITY 8
VA-BFA PAIN INTERFERE ACTIVITY 7
VA-BFA PAIN INTERFERE ACTIVITY 6

REMINDER SPONSOR

Office of Patient Care Services

REMINDER TERM

VA-REMINDER UPDATE_2_0_92

REMINDER DIALOG

VA-BFA BATTLEFIELD ACUPUNCTURE NOTE

Install Details

This update is being distributed as a web host file. The address for the host file is:

http://vista.med.va.gov/reminders/UPDATE_2_0_92.PRD

The file will be installed using Reminder Exchange, programmer access is not required.

Installation:


=====

This update can be loaded with users on the system. Installation will take less than 20 minutes.

Install Example

If your site is unable to use LWH functionality, you will use vista mail to install the exchange file.

To Load the Web Host File. Navigate to Reminder exchange in Vista

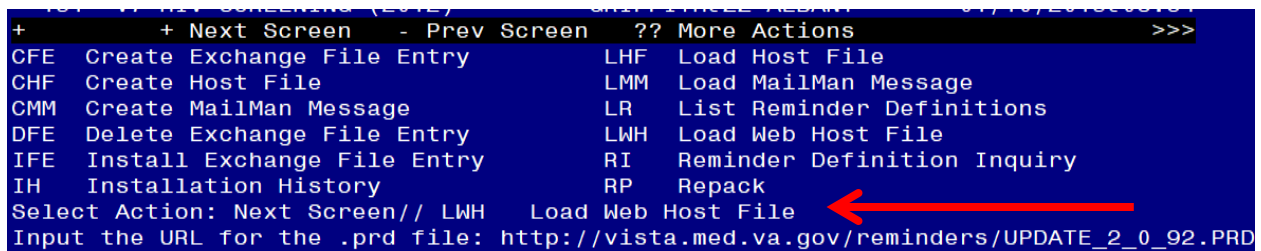

```
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                LMM  Load MailMan Message
CMM  Create MailMan Message          LR   List Reminder Definitions
DFE  Delete Exchange File Entry      LWH  Load Web Host File
IFE  Install Exchange File Entry     RI   Reminder Definition Inquiry
IH   Installation History             RP   Repack
Select Action: Next Screen// LWH    Load Web Host File
Input the URL for the .prd file: http://vista.med.va.gov/reminders/UPDATE_2_0_92.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

http://vista.med.va.gov/reminders/UPDATE_2_0_92.PRD

```
http://vista.med.va.gov/reminders/UPDATE_2_0_92.PRD successfully loaded.
```

You should see a message at the top of your screen that the file successfully loaded. Search and locate an entry titled **UPDATE_2_0_92 VA-BFA TEMPLATE UPDATE** in reminder exchange

+Item	Entry	Source	Date Packed
194	UPDATE_2_0_92 VA-BFA TEMPLATE UPDATE	GRIFFITH@SALT LAKE CI	08/06/2019@07:14
195	V7-HIV SCREENING (2012)	MURPHY@ATLANTA VAMC	01/09/2019@14:51
196	V7-HIV SCREENING (2012)	GRIFFITH@ZZ ALBANY	01/10/2019@06:25
197	V7-HIV SCREENING (2012)	GRIFFITH@ZZ ALBANY	01/10/2019@08:54
198	VA BRANCHING LOGIC REMINDER UPDATES OEF/OIF	VOLPP@NORTHERN CAL	09/13/2010@11:42
199	VA MH SCREENING REMINDERS UPDATE	VOLPP@NORTHERN CAL	09/13/2010@11:36
200	VA VACO TRANSPLANT REFERRAL	ROZELLE@N. FLORIDA/S	11/14/2011@13:37
+ + Next Screen - Prev Screen ?? More Actions >>>			
CFE	Create Exchange File Entry	LHF	Load Host File
CHF	Create Host File	LMM	Load MailMan Message
CMM	Create MailMan Message	LR	List Reminder Definitions
DFE	Delete Exchange File Entry	LWH	Load Web Host File
IFE	Install Exchange File Entry	RI	Reminder Definition Inquiry
IH	Installation History	RP	Repack
Select Action: Next Screen// IFE Install Exchange File Entry			
Enter a list or range of numbers (1-292): 194			

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry
Enter the number that corresponds with your entry titled **UPDATE_2_0_92 VA-BFA TEMPLATE UPDATE** (in this example it is entry 194, it will vary by site)

Component	Category	Exists
Source: GRIFFITH,ELIZABETH A at SALT LAKE CITY		
Date Packed: 08/06/2019@07:14:21		
Package Version: 2.0P35		
Description:		
The following Clinical Reminder items were selected for packing:		
REMINDER DIALOG		
VA-BFA BATTLEFIELD ACUPUNCTURE NOTE		
REMINDER TERM		
VA-REMINDER UPDATE_2_0_92		
Keywords:		
Components:		
+ Enter ?? for more actions >>>		
IA	Install all Components	IS Install Selected Component
Select Action: Next Screen// IA		

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// **IA Install all Components**

During the installation, you will see several prompts: *Note that you will **Install** all new components*

You will be prompted to install the dialog component – you will **INSTALL ALL** even though all components may have an X to indicate that they already exist, the components have been updated and should be installed:

Packed reminder dialog: VA-BFA BATTLEFIELD ACUPUNCTURE NOTE [NATIONAL DIALOG]

Item	Seq.	Dialog Findings	Type	Exists
1		VA-BFA BATTLEFIELD ACUPUNCTURE NOTE	dialog	X
2	1	VA-BFA BATTLEFIELD ACUPUNCTURE NOTE VERSION Finding: *NONE*	element	X
3	5	VA-BFA GP INTIAL/FOLLOW UP Finding: *NONE*	group	X
4	5.5	VA-BFA GP INITIAL ONLY Finding: *NONE*	group	X
5	5.10	VA-BFA GP FOLLOW UP ONLY Finding: *NONE*	group	X
6		PXRM COMMENT	prompt	X
7	8	VA-BFA GP TYPE OF TREATMENT Finding: *NONE*	group	X

+ + Next Screen - Prev Screen ?? More Actions

DD	Dialog Details	DT	Dialog Text	IS	Install Selected
DF	Dialog Findings	DU	Dialog Usage	QU	Quit
DS	Dialog Summary	IA	Install All		

Select Action: Next Screen// IA

At the **Select Action** prompt, type **IA** to install the dialog – **VA-BFA BATTLEFIELD ACUPUNCTURE NOTE**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: **Y// ES**

Packed reminder dialog: VA-BFA BATTLEFIELD ACUPUNCTURE NOTE [NATIONAL DIALOG]

VA-BFA BATTLEFIELD ACUPUNCTURE NOTE (reminder dialog) installed from exchange fi

Item	Seq.	Dialog Findings	Type	Exists
1		VA-BFA BATTLEFIELD ACUPUNCTURE NOTE	dialog	X
2	1	VA-BFA BATTLEFIELD ACUPUNCTURE NOTE VERSION Finding: *NONE*	element	X
3	5	VA-BFA GP INTIAL/FOLLOW UP Finding: *NONE*	group	X
4	5.5	VA-BFA GP INITIAL ONLY Finding: *NONE*	group	X
5	5.10	VA-BFA GP FOLLOW UP ONLY Finding: *NONE*	group	X
6		PXRM COMMENT	prompt	X
7	8	VA-BFA GP TYPE OF TREATMENT Finding: *NONE*	group	X

+ + Next Screen - Prev Screen ?? More Actions

DD	Dialog Details	DT	Dialog Text	IS	Install Selected
DF	Dialog Findings	DU	Dialog Usage	QU	Quit
DS	Dialog Summary	IA	Install All		

Select Action: Next Screen// Q

When the dialog has completed installation, you will then be returned to this screen. At the **Select Action** prompt, type **Q**.

Component	Category	Exists
Source: GRIFFITH,ELIZABETH A at SALT LAKE CITY		
Date Packed: 08/06/2019@07:14:21		
Package Version: 2.0P35		
Description:		
The following Clinical Reminder items were selected for packing:		
REMINDER DIALOG		
VA-BFA BATTLEFIELD ACUPUNCTURE NOTE		
REMINDER TERM		
VA-REMINDER UPDATE_2_0_92		
Keywords:		
Components:		
+ Enter ?? for more actions >>		
IA	Install all Components	IS Install Selected Component
Select Action: Next Screen// Q		

You will be returned to this screen. Type Q to quit.

Install Completed

Post Install

When you install the updated template, it should replace the existing template in your system, no post install setup is required. To confirm that the install was successful, if you open the note or template, you should see a version number of 3.1 post install.

NOTE: If you require further technical assistance, if there are any questions or problems during the installation, please notify your local IT support to log a national help desk ticket Desk and have them submit a national ticket to:

Category: Enterprise Applications

Enterprise Application: VistA - Clinical Reminders

Assignment Group: NTL SUP Clin 2