



Update _2_0_45

Clinical Reminders

VA-Maternity Care Coordinator Note

Install Guide

March 2019

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION	1
PRE-INSTALLATION	4
INSTALL DETAILS	4
INSTALL EXAMPLE	5
POST INSTALLATION	8

Introduction

Description:

This national dialog is needed to standardized documentation of Maternity Care Coordination, as well as to capture health factor data of maternity care for veteran patients.

The dialog template is for maternity care coordinators (MCCs), or those acting in that capacity, to document and track pregnant Veterans whose maternity care is being paid for by the VA, according to the VHA Directive 1330.03.

The dialog template matches the current call templates being used which exist as Word templates. Calls are made at specified gestational weeks and touch on topics MCCs are required to review with the Veterans. The MCCs are currently documenting using the MCC Telephone Care Program generated Word templates, copying and pasting into CPRS, or not using the templates at all thus risking missing topics needed to be covered. In addition - to capture required data on maternal and fetal outcomes, additional work is required to gather data on those outcomes-facilities have various ways they are doing this including using SharePoint forms and Excel spreadsheets – none of that data is being collected at a national level.

Our goal is to provide an electronic template for this required documentation, to capture data and to ease the overall workload of the MCCs. The format will be familiar to the MCCs.

One reminder dialog is being distributed with this update:

UPDATE_2_0_45 contains 1 Reminder Exchange entry:

UPDATE_2_0_45 VA-MATERNITY CARE COORDINATOR NOTE

The exchange file contains the following components:

TIU TEMPLATE FIELD

VA-MCC DIRECTIVE
TEXT (1-60 CHARACTERS) REQ
TEXT 3 SPACES
DATE
TEXT (6 CHAR)
TEXT (1-60 CHAR)1
MCC GESTATE AGE
DATE (*)
BLANK SPACE1
TEXT (1-60 CHAR)

HEALTH FACTORS

VA-REMINDER UPDATES
VA-UPDATE_2_0_45
VA-MCC CALLS (C)

VA-MCC CALL7 MADE
VA-MCC CALL6B MADE
VA-MCC CALL6A MADE
VA-MCC CALL5 MADE
VA-MCC CALL4 MADE
VA-MCC CALL3 MADE
VA-MCC CALL2 MADE
VA-MCC CALL1 MADE
VA-MCC INTIMATE PARTNER VIOLENCE (C)
VA-MCC VIOLENCE YES
VA-MCC VIOLENCE NO
VA-MCC DEPRESSION (C)
VA-MCC DEPRESSION EDINBURG SCREEN GIVEN
VA-MCC DEPRESSION SCREEN POSITIVE
VA-MCC DEPRESSION SCREEN INTERMEDIATE
VA-MCC DEPRESSION SCREEN NEGATIVE
VA-MCC DELIVERY INFORMATION (C)
VA-MCC BABY HOSP STAY DOCUMENTED
VA-MCC PATIENT HOSP STAY DOCUMENTED
VA-MCC COMPLICATION CLAMPSIA
VA-MCC COMPLICATION INFECTIONS
VA-MCC COMPLICATION HEMORRAGE
VA-MCC LABOR HOURS DOCUMENTED
VA-MCC BIRTH WEIGHT DOCUMENTED
VA-MCC DELIVERED CSECTION
VA-MCC DELIVERED ASSIST VAGINAL
VA-MCC DELIVERED VAGINAL
VA-MCC BIRTH MULTIPLES
VA-MCC BIRTH TWINS
VA-MCC BIRTH SINGLETON
VA-MCC DELIVERY DATE
VA-MCC FETAL DEMISE
VA-MCC LIVE BIRTH
VA-MCC TOBACCO USE (C)
VA-MCC TOBACCO SMOKE 2ND HAND EXPOSE NO
VA-MCC TOBACCO SMOKE 2ND HAND EXPOSE YES
VA-MCC TOBACCO SMOKES DOES NOT
VA-MCC TOBACCO SMOKES RECENTLY QUIT
VA-MCC TOBACCO SMOKES YES
VA-MCC ALCOHOL USE (C)
VA-MCC ALCOHOL DOES CONSUME
VA-MCC ALCOHOL DOES NOT CONSUME
VA-MCC HISTORY OF OBSTETRIC PROBLEMS (C)
VA-MCC HISTORY OF PRETERM BIRTH
VA-MCC HISTORY OF PRECLAMPIA
VA-MCC HISTORY OF POSTPARTUM DEPRESSION

VA-MCC HISTORY OF GESTATIONAL DIABETES
VA-MCC HISTORY OF CESAREAN
VA-MCC HISTORY OF HEALTH PROBLEMS (C)
VA-MCC HISTORY OF OBESITY
VA-MCC HISTORY OF THYROID
VA-MCC HISTORY OF STI
VA-MCC HISTORY OF HYPERTENSION
VA-MCC HISTORY OF RENAL DISEASE
VA-MCC HISTORY OF SEIZURE DISORDER
VA-MCC HISTORY OF NACOTIC OBESITY
VA-MCC HISTORY OF LUPUS
VA-MCC HISTORY OF HYPERCOAGULABLE
VA-MCC HISTORY OF HIV
VA-MCC HISTORY OF DIABETES
VA-MCC HISTORY OF DEPRESSION/MH DISORDER
VA-MCC HISTORY OF CHRONIC PAIN
VA-MCC HISTORY OF ASTHMA

REMINDER SPONSOR

WOMEN VETERANS HEALTH PROGRAM

REMINDER TERM

VA-REMINDER UPDATE_2_0_45

REMINDER DIALOG

VA-MATERNITY CARE COORDINATION

HEALTH SUMMARY TYPE

VA-MCC CALLS MADE
VA-MCC DELIVERY INFORMATION
PALLI PHQ2 PH9
VA-MCC TOBACCO
VA-MCC ALCOHOL
VA-MCC IPV
VA-MCC DEPRESSION SCREEN
VA-MCC HEALTH PROBLEMS

HEALTH SUMMARY OBJECTS

VA-MCC CALLS MADE (TIU)
VA-MCC DELIVERY INFORMATION (TIU)
PALLI PHQ2 PH9 (TIU)
VA-MCC TOBACCO (TIU)
VA-MCC ALCOHOL (TIU)
VA-MCC IPV (TIU)
VA-MCC DEPRESSION SCREEN (TIU)
VA-MCC HEALTH PROBLEMS (TIU)

TIU DOCUMENT DEFINITION

VA-MCC CALLS MADE
VA-MCC DELIVERY INFORMATION
PALLI PHQ2 PHQ9 (TIU)
VA-MCC TOBACCO
VA-MCC ALCOHOL
VA-MCC INTIMATE PARTNER VIOLENCE
VA-MCC DEPRESSION SCREEN
VA-MCC HEALTH PROBLEMS

Pre-Installation

Setup of Note title

This note title may already exist at your site, If not, create the following and add it to the same document class as your other women's health notes.

Name: **VA-MATERNITY CARE COORDINATION**
VHA Enterprise Standard Title: **WOMENS HEALTH PROGRESS NOTE** OR
WOMENS HEALTH NOTE

Comprehensive information on Creating Document Definitions can be found beginning on page 47 of the [TIU/ASU Implementation Guide](#)

Install Details

This update is being distributed as a web host file. The address for the host file is:
http://vista.med.va.gov/reminders/UPDATE_2_0_45.PRD

The file will be installed using Reminder Exchange, programmer access is not required.
The total time for installation should be less than 30 minutes.

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

```
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                LMM  Load MailMan Message
CMM  Create MailMan Message          LR   List Reminder Definitions
DFE  Delete Exchange File Entry      LWH  Load Web Host File
IFE  Install Exchange File Entry      RI   Reminder Definition Inquiry
IH   Installation History
Select Action: Next Screen// lwh  Load Web Host File
Input the URL for the .prd file: http://vista.med.va.gov/reminders/UPDATE_2_0_45.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

http://vista.med.va.gov/reminders/UPDATE_2_0_45.PRD

http://vista.med.va.gov/reminders/UPDATE_2_0_45.PRD successfully loaded.

You should see a message at the top of your screen that the file successfully loaded.

Search and locate an entry titled **UPDATE_2_0_45 VA-MATERNITY CARE COORDINATOR NOTE** in reminder exchange

```
+Item  Entry                                     Source                                Date Packed
118  UPDATE_2_0_45 VA-MATERNITY CARE COORDINATOR NOTE  PLOTT@SALT LAKE CI  02/08/2019@10:51
119  UPDATE_2_0_46 VA-MH TOBACCO USE SCREEN          GRIFFITH@SALT LAKE CI  08/14/2018@09:34
120  UPDATE_2_0_47 VA-HEPATITIS B IMMUNIZATIONS      MONTGOMERY@SALT LAKE CI  08/15/2018@13:48
121  UPDATE_2_0_48 MISC UPDATE 6                     GRIFFITH@SALT LAKE CI  05/09/2018@13:39
122  UPDATE_2_0_49 VA-SUICIDE BEHAVIOR AND OVERDOSE REPORT  GRIFFITH@SALT LAKE CI  01/17/2019@06:50
123  UPDATE_2_0_5 MISC UPDATE 1                     GRIFFITH@SALT LAKE CI  05/03/2016@09:42
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE  Create Exchange File Entry      LHF  Load Host File
CHF  Create Host File                LMM  Load MailMan Message
CMM  Create MailMan Message          LR   List Reminder Definitions
DFE  Delete Exchange File Entry      LWH  Load Web Host File
IFE  Install Exchange File Entry      RI   Reminder Definition Inquiry
IH   Installation History
Select Action: Next Screen// IFE  Install Exchange File Entry
Enter a list or range of numbers (1-237): 118
```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry
Enter the number that corresponds with your entry titled **UPDATE_2_0_45 VA-MATERNITY CARE COORDINATOR NOTE** (*in this example it is entry 118 it will vary by site*)

```

Component                                     Category   Exists
Source:      PLOTT,JOSEPH G at SALT LAKE CITY
Date Packed: 02/08/2019@10:51:13
Package Version: 2.0P35

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
    VA-MATERNITY CARE COORDINATION

REMINDER TERM
    VA-REMINDER UPDATE_2_0_45
Non-exchangeable TIU object(s):
-----
TIU Object: PATIENT AGE
Object Method: S X=$$AGE^TIULO(DFN)

+      Enter ?? for more actions      >>>
IA  Install all Components             IS  Install Selected Component
Select Action: Next Screen// IA

```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// **IA Install all Components**

You will see several prompts, for all new entries you will choose **I to Install**

For components that already exists but the packed component is different, you will choose **O to Overwrite**. **One Exception You May Be Prompted to address:** Choose **S to Skip** if prompted for Health Summary Type “PALLI PHQ2 PHQ9”

HEALTH SUMMARY TYPE entry named PALLI PHQ2 PH9 already exists but the packed component is different, what do you want to do?

Select one of the following:

- C Create a new entry by copying to a new name
- O Overwrite the current entry
- U Update
- Q Quit the install
- S Skip, do not install this entry

Enter response: O//**S**

You will be prompted to install the reminder dialog component:

```
Packed reminder dialog: VA-MATERNITY CARE COORDINATION [NATIONAL DIALOG]
```

+Item	Seq.	Dialog Findings	Type	Exists
15	45	BLANK LINE NOTE ONLY Finding: *NONE*	element	
16	50	VA-GP MCC PREGNANCY STATUS Finding: *NONE*	group	
17	50.10	VA-MCC PREGNANCY STATUS PREGNANT Finding: *NONE*	element	
18	50.20	VA-MCC PREGNANCY STATUS DELIVERED Finding: *NONE*	element	
19	50.30	VA-GP MCC PREGNANCY LOSS Finding: *NONE*	group	
20	50.30.10	VA-MCC PREGNANCY LOSS DISCUSS Finding: *NONE*	element	
21	50.30.20	VA-MCC PREGNANCY LOSS REFER Finding: *NONE*	element	
+ + Next Screen - Prev Screen ?? More Actions				
DD	Dialog Details	DT	Dialog Text	IS Install Selected
DF	Dialog Findings	DU	Dialog Usage	QU Quit
DS	Dialog Summary	IA	Install All	

Select Action: Next Screen// IA

At the **Select Action** prompt, type **IA** to install the dialog – **VA-MATERNITY CARE COORDINATOR**. Select Action: Next Screen// **IA Install All**
Install reminder dialog and all components with no further changes: Y// **Yes**

```
Packed reminder dialog: VA-MATERNITY CARE COORDINATION [NATIONAL DIALOG]
```

```
VA-MATERNITY CARE COORDINATION (reminder dialog) installed from exchange file.
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-MATERNITY CARE COORDINATION	dialog	X
2	10	VA-MCC PREVIOUS CALLS INTRO TEXT Finding: *NONE*	element	X
3	20	VA-GP MCC PREVIOUS CALLS Finding: *NONE*	group	X
4	20.10	VA-GP MCC REVIEW PREVIOUS CALLS MADE Finding: *NONE*	group	X
5	20.10.10	VA-MCC DATA IMPORT PREVIOUS CALLS MADE Finding: *NONE*	element	X
6	20.20	VA-GP MCC REVIEW PREVIOUS CALL NOTES Finding: *NONE*	group	X
7	20.20.10	VA-MCC DATA IMPORT PREVIOUS CALL NOTES Finding: *NONE*	element	X
+ + Next Screen - Prev Screen ?? More Actions				
DD	Dialog Details	DT	Dialog Text	IS Install Selected
DF	Dialog Findings	DU	Dialog Usage	QU Quit
DS	Dialog Summary	IA	Install All	

Select Action: Next Screen// 0

When the dialog has completed installation, you will then be returned to this screen. At the **Select Action** prompt, type **Q**.

You will then be prompted to install the health summary components. **Choose I to install all components. One Exception You May Be Prompted to address:** Choose **S to Skip** if prompted for Health Summary Type “PALLI PHQ2 PHQ9”

HEALTH SUMMARY TYPE entry named PALLI PHQ2 PH9 already exists but the packed component is different, what do you want to do?

Select one of the following:

- C Create a new entry by copying to a new name
- O Overwrite the current entry
- U Update
- Q Quit the install
- S Skip, do not install this entry

Enter response: O// **S**

Component	Category	Exists
Source: PLOTT,JOSEPH G at SALT LAKE CITY		
Date Packed: 02/08/2019@10:51:13		
Package Version: 2.0P35		
Description:		
The following Clinical Reminder items were selected for packing:		
REMINDER DIALOG		
VA-MATERNITY CARE COORDINATION		
REMINDER TERM		
VA-REMINDER UPDATE_2_0_45		
Non-exchangeable TIU object(s):		

TIU Object: PATIENT AGE		
Object Method: S X=\$\$AGE*TIULO(DFN)		
+ Next Screen - Prev Screen ?? More Actions >>>		
IA	Install all Components	IS Install Selected Component
Select Action: Next Screen// Q		

You will then be returned to this screen. At the **Select Action** prompt, type **Q**.
Install complete

Post Installation

1. Make the Dialog able to be attached to a Personal/Shared template or to a progress note title.

From the Reminders Manager Menu:

Select Reminder Managers Menu Option: **CP CPRS Reminder Configuration**

CA Add/Edit Reminder Categories

CL CPRS Lookup Categories

CS CPRS Cover Sheet Reminder List

MH Mental Health Dialogs Active

PN Progress Note Headers

RA Reminder GUI Resolution Active

TIU TIU Template Reminder Dialog Parameter

DL Default Outside Location

PT Position Reminder Text at Cursor

Select CPRS Reminder Configuration Option: **TIU TIU Template Reminder Dialog Parameter**

Reminder Dialogs allowed as Templates may be set for the following:

1 User USR [choose from NEW PERSON]

3 Service SRV [choose from SERVICE/SECTION]

4 Division DIV [choose from INSTITUTION]

5 System SYS [NATREM.FO-SLC.MED.VA.GOV]

Enter selection: **5** Choose the appropriate number for your site. Your site may do this by System or other levels. The example below uses SYSTEM level

Setting Reminder Dialogs allowed as Templates for System:

Select Display Sequence: **?**

148 VA-VETERANS CHOICE NOTE

149 VA-TBI SCREENING

When you type a question mark above, you will see the list of #'s (with dialogs) that are already taken. Choose a number **NOT** on this list. For this example, looking above I see the number 150 is not present, so I will use 150.

Select Display Sequence: **150**

Are you adding 150 as a new Display Sequence? Yes// **Y YES**

Display Sequence: 150// **<Enter>** 150

Clinical Reminder Dialog: **type in the name of the dialog**

VA-MATERNITY CARE COORDINATION then **<enter>** reminder dialog NATIONAL

OK? Yes// **<Enter>** (Yes)

2. Associate the reminder dialog with the note title or shared template in CPRS


To do this, go to Template Editor, then Document Titles, select the appropriate folder (if applicable), and choose New Template (top right-hand corner).

1. Type in the following information:


Name: **VA-MATERNITY CARE COORDINATION**

Template Type: **Reminder Dialog**

Reminder Dialog: **Va-Maternity Care Coordination**

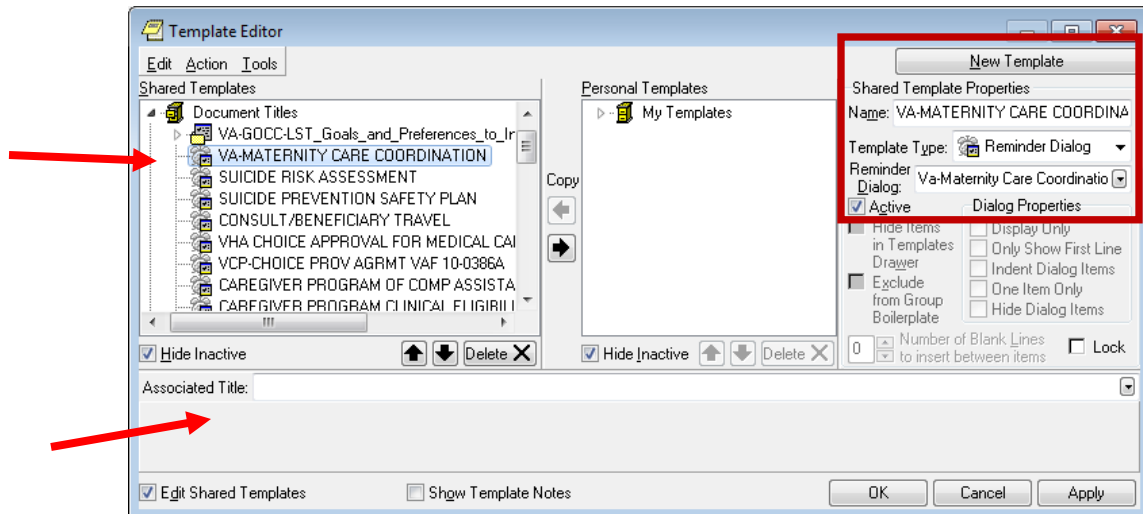

Next, type in the Associated Title, which in this case is:

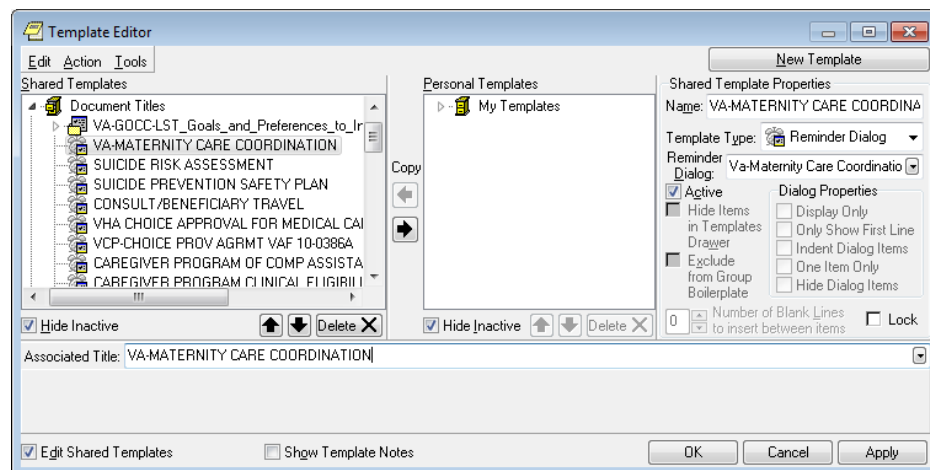

Hit Apply and OK

- Complete section **3a** (edit object) if the VA-MCC (TIU) object exists in your system.
Complete section **3b** (create object) if the VA-MCC (TIU) object does NOT exist in your system.

3a - Edit the VA-MCC (TIU) Health Summary Type: (step by step is highlighted in yellow)

Log into VISTA and go to the Document Definitions menu:
^document definitions (Manager)

```

--- Manager Document Definition Menu ---
<TEST ACCOUNT> 1  Edit Document Definitions
<TEST ACCOUNT> 2  Sort Document Definitions
<TEST ACCOUNT> 3  Create Document Definitions
<TEST ACCOUNT> 4  Create Objects
<TEST ACCOUNT> 5  Create TIU/Health Summary Objects

```

Select Document Definitions (Manager) <TEST ACCOUNT> Option:
5 Create TIU/Health Summary Objects

TIU Health Summary Object Sep 07, 2018@11:10:54 Page: 1 of 28

TIU Object Name	Health Summary Type
1 ABD IMAGES LAST 5Y	ABD IMAGES LAST 5Y
2 ABDOMINAL ULTRASOUND	ABDOMINAL ULTRASOUND
3 ACTIVE INPT MEDS	ACTIVE INPT MEDS
4 ACTIVE OUTPT MEDS	ACTIVE OUTPT MEDS
5 ACTIVE PROBLEMS W/COMMENT	ACTIVE PROBLEMS W/COMMENT
+ Enter ?? for more actions	
Create New TIU Object	Find
Detailed Display/Edit TIU Object	Detailed Display/Edit HS Object
Quit	

Select Action: Next Screen// sl SL
Search for: **VA-MCC (TIU)**

TIU Health Summary Object Sep 07, 2018@10:48:54 Page: 19 of 28

TIU Object Name	Health Summary Type
253 VA-MCC (TIU)	VA-MCC (TIU)
254 VA-MCC ALCOHOL	VA-MCC ALCOHOL
255 VA-MCC CALLS MADE	VA-MCC CALLS MADE
256 VA-MCC DELIVERY INFORMATION	VA-MCC DELIVERY INFORMATION
+ Enter ?? for more actions	
Create New TIU Object	Find
Detailed Display/Edit TIU Object	Detailed Display/Edit HS Object
Quit	

Select Action: Next Screen// **253**

Note: the sequence number for VA-MCC (TIU) may be different in your system

TIUHS Detailed Display/Edit Sep 07, 2018@10:48:56 Page: 1 of 1

TIU Object Name: VA-MCC (TIU)
Owner: xxxxxx
Status: ACTIVE

HS Object: VA-MCC (TIU) (TIU)
HS Type: VA-MCC (TIU)
Technical Field: S X=\$\$TIU^GMTSOBJ(DFN,6580193)

Enter ?? for more actions
Change HS Object Detail Display/Edit HS Object
Change Health Summary Type
Select Item(s): Quit// **det** Detail Display/Edit HS Object

HS OBJECT DISPLAY Sep 07, 2018@10:49:02
Detailed Display for VA-MCC (TIU)

HS Object: VA-MCC (TIU) (TIU)
Health Summary Type: VA-MCC (TIU)
Report Period:
Creator:
HS Object

Print Label: NO Print Report Date and Time: NO
Print Blank Line after Label: NO Print Confidentiality Banner: NO
Customized Header: YES Print Report Date and Time: NO
Suppress Components w/o Data: YES Print Component Header: NO
Print Deceased Information: NO Print Time-Occurrence Limits: NO
National Object: NO Underline Component Header: NO
Blank Line After Header: NO
Overwrite No Data:

Enter ?? for more actions
Select Action: Quit// **ed**
1 Edit HS Object
2 Edit HS Type
CHOOSE 1-2: **2** Edit HS Type

Editing Health Summary Type 'VA-MCC (TIU)'

NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)//
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: yes//
LOCK:
OWNER: <ENTER YOUR NAME HERE>
Select COMPONENT: SPN

Select COMPONENT: spn
1 SPN1 SPN1
2 SPN2 SPN2
3 SPN3
4 SPN4
5 SPN5

Press <Enter> to see more, '^' to exit this list, OR
CHOOSE 1-5:
6 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-6: 6 PROGRESS NOTES SELECTED SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME: Selected Prog Notes//

No selection items chosen.

Select new items one at a time in the sequence you want them displayed.
You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH PROGRESS NOTE
...OK? Yes// (Yes)

Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Type <Enter> to continue or '^' to exit:


Done!

3b. Create the VA-MCC (TIU) Health Summary Type:


You only need to complete this Post-Installation step if the VA-MCC (TIU) object does not exist in your system. If you completed step 3a, you do not need to complete this step.

Step 1: Log into VISTA


Step 2: Go to the DOCUMENT DEFINITION menu and select CREATE DOCUMENT DEFINITION;
Option 5

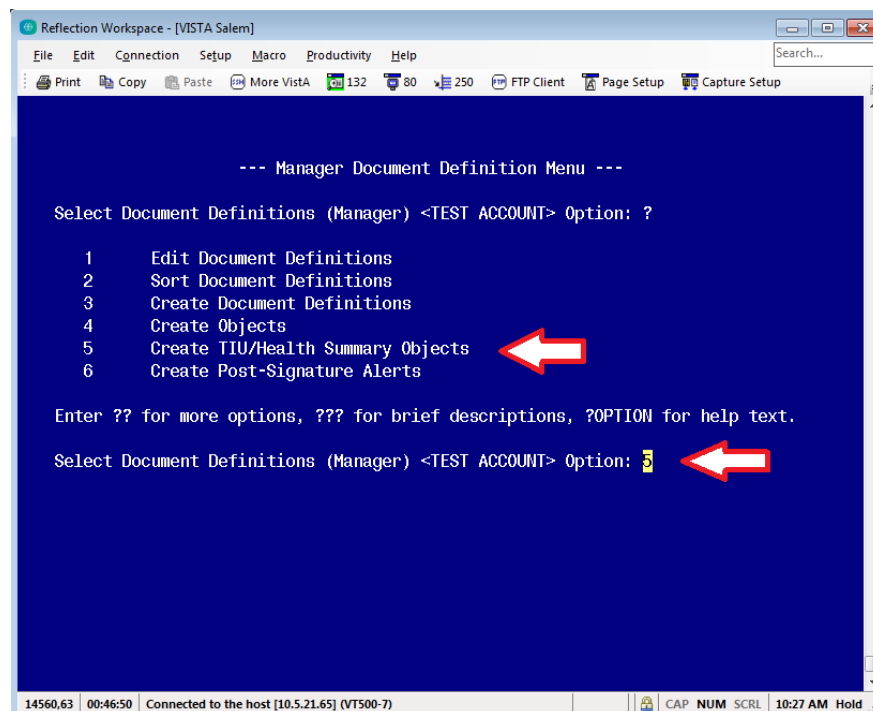


Step 3: Type “Create New TIU Object” and press ENTER

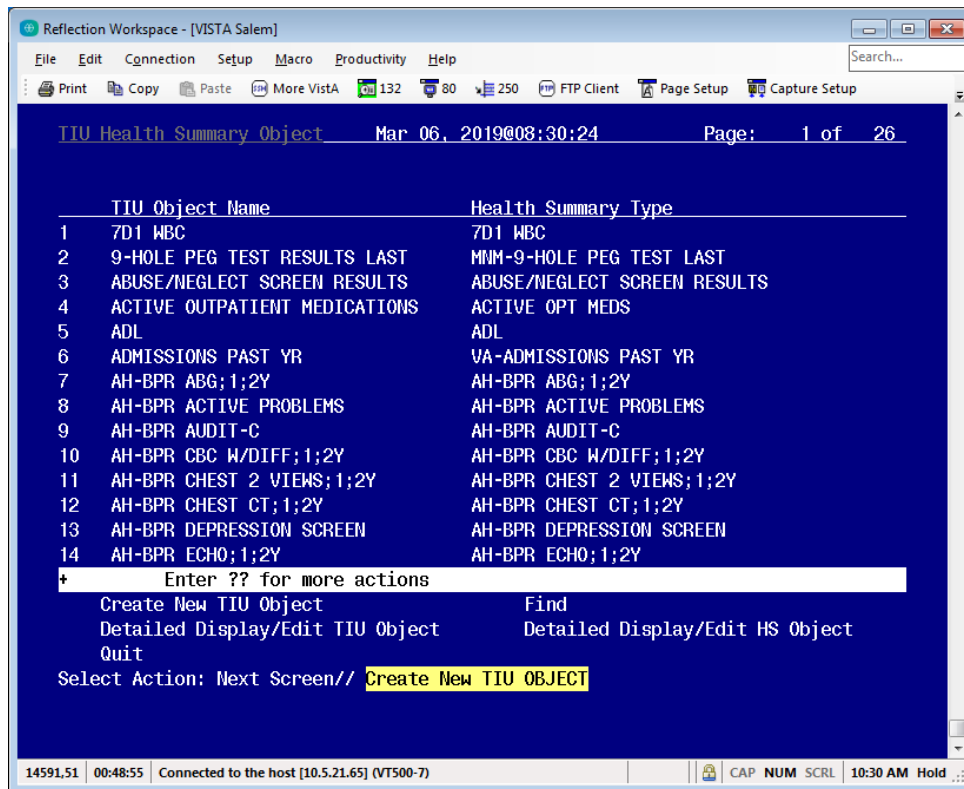


Step 4: At the ‘Enter a new TIU OBJECT NAME’ prompt, type “VA-MCC (TIU)”

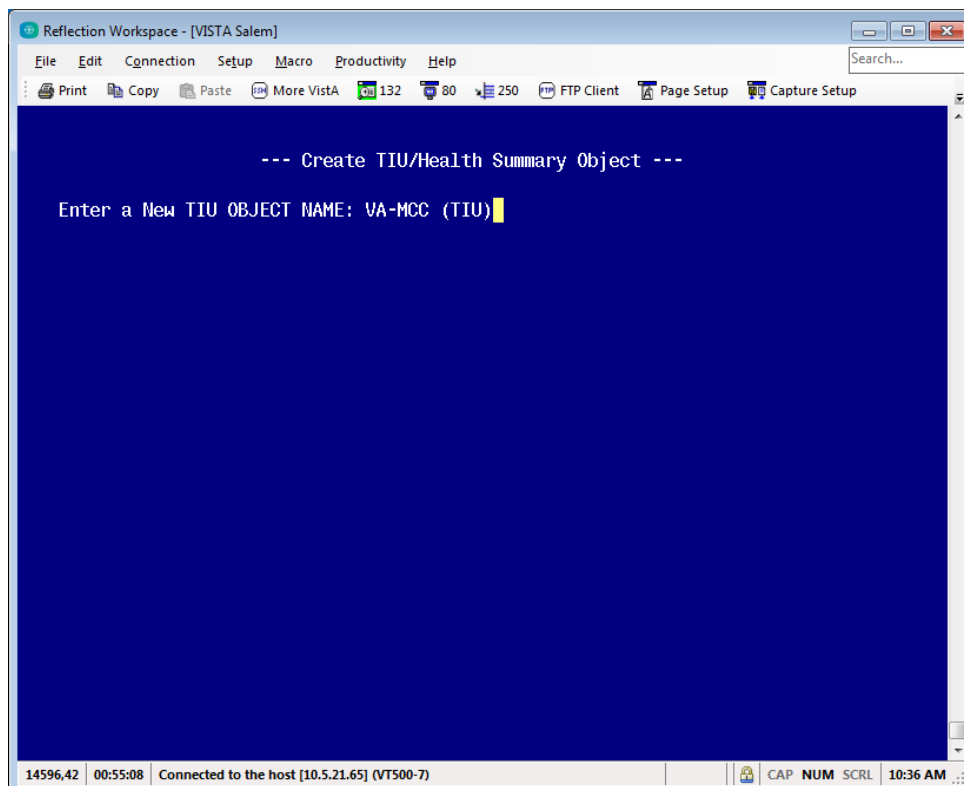


Step 5: At the ‘Is this correct?’ prompt, press ENTER for the default “YES” ... and at the ‘Use a pre-existing Health Summary Object?’ prompt, press ENTER for the default “NO”

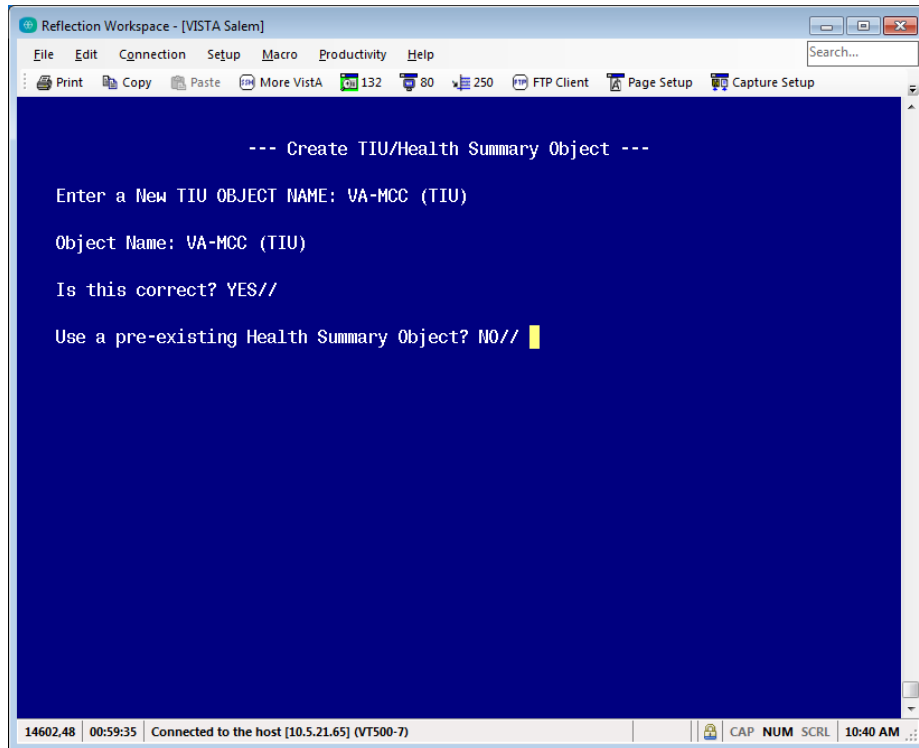

Step 6: At the ‘Select Health Summary Type:’ prompt, type “VA-MCC (TIU)”

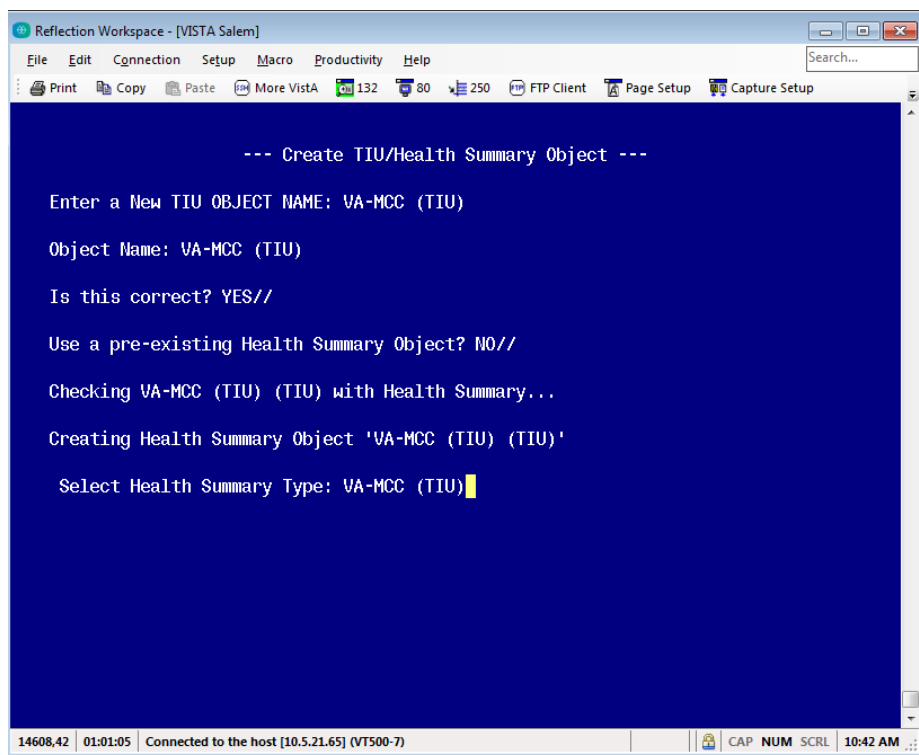

Step 7: At the ‘Are you adding VA-MCC (TIU) ...’ prompt, type “YES”

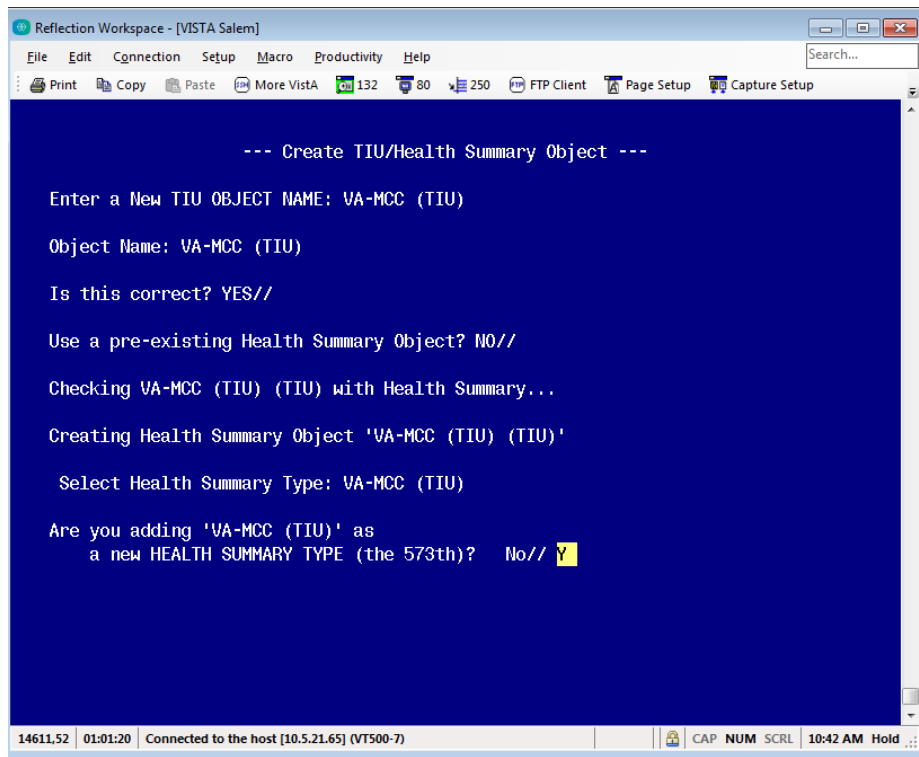


The screenshot shows a terminal window titled 'Reflection Workspace - [VISTA Salem]'. The terminal displays the following text:

```
--- Create TIU/Health Summary Object ---  
Enter a New TIU OBJECT NAME: VA-MCC (TIU)  
Object Name: VA-MCC (TIU)  
Is this correct? YES//  
Use a pre-existing Health Summary Object? NO//  
Checking VA-MCC (TIU) (TIU) with Health Summary...  
Creating Health Summary Object 'VA-MCC (TIU) (TIU)'  
Select Health Summary Type: VA-MCC (TIU)  
Are you adding 'VA-MCC (TIU)' as  
a new HEALTH SUMMARY TYPE (the 573th)? No// Y
```

The status bar at the bottom shows: 14611.52 | 01:01:20 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL | 10:42 AM Hold ...

Step 8: At the ‘NAME’ prompt, press ENTER

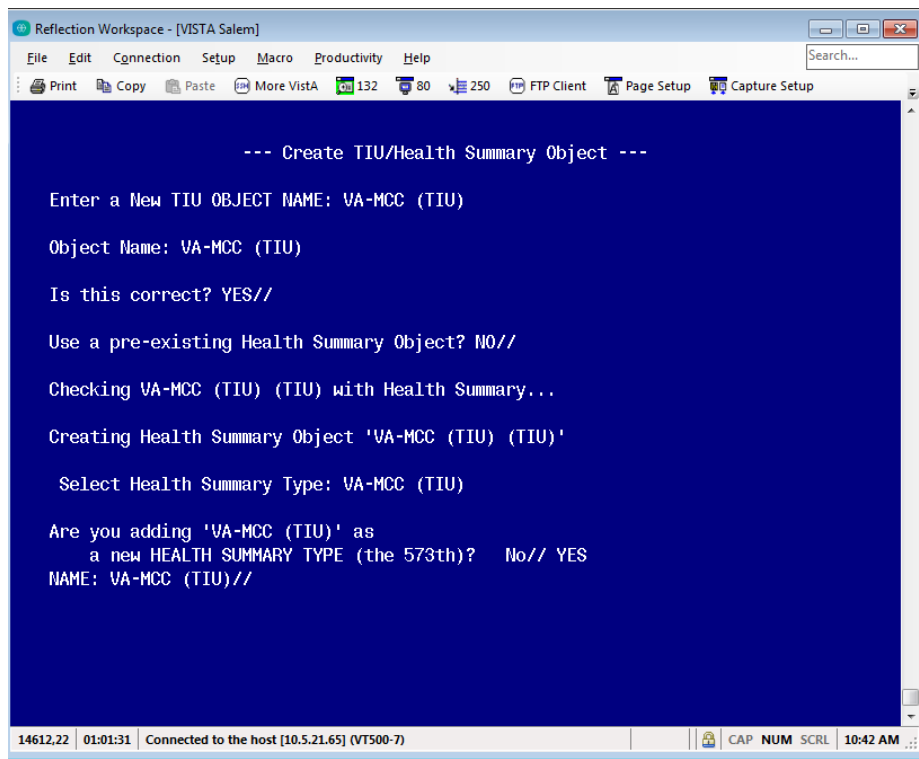


The screenshot shows the same terminal window as in Step 7, but with additional text at the bottom of the prompt sequence:


```
--- Create TIU/Health Summary Object ---  
Enter a New TIU OBJECT NAME: VA-MCC (TIU)  
Object Name: VA-MCC (TIU)  
Is this correct? YES//  
Use a pre-existing Health Summary Object? NO//  
Checking VA-MCC (TIU) (TIU) with Health Summary...  
Creating Health Summary Object 'VA-MCC (TIU) (TIU)'  
Select Health Summary Type: VA-MCC (TIU)  
Are you adding 'VA-MCC (TIU)' as  
a new HEALTH SUMMARY TYPE (the 573th)? No// YES  
NAME: VA-MCC (TIU)//
```

The status bar at the bottom shows: 14612.22 | 01:01:31 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL | 10:42 AM ...


Step 9: At the 'TITLE' prompt, type "VA-MCC (TIU)" and press ENTER

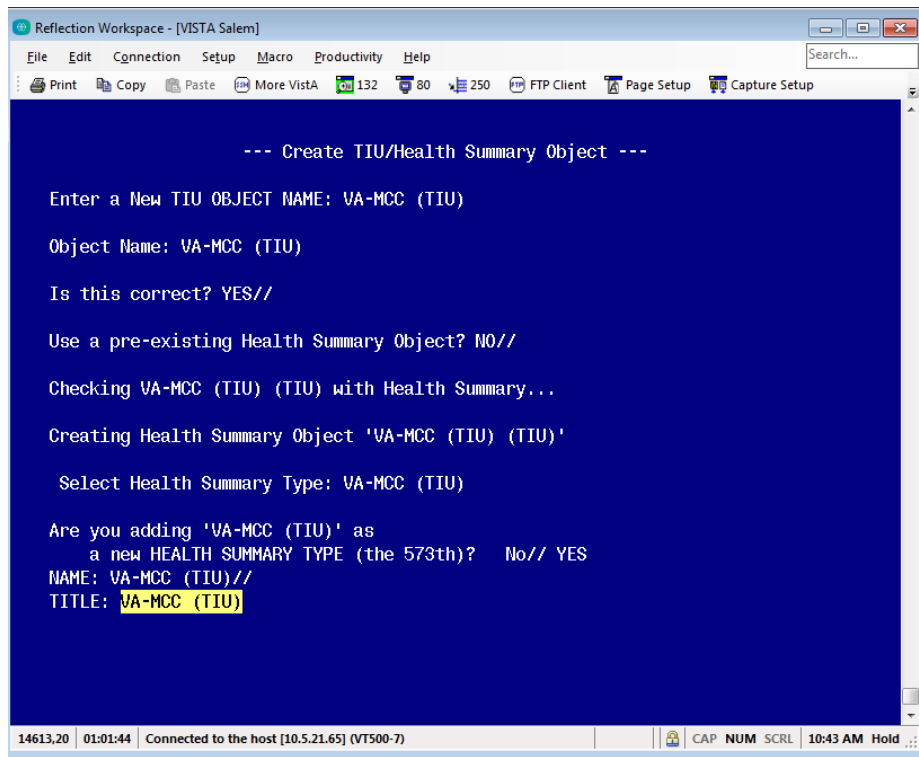

Step 10: At the 'SUPPRESS PRINT OF COMPONENTS ...' prompt, type "YES" and press ENTER

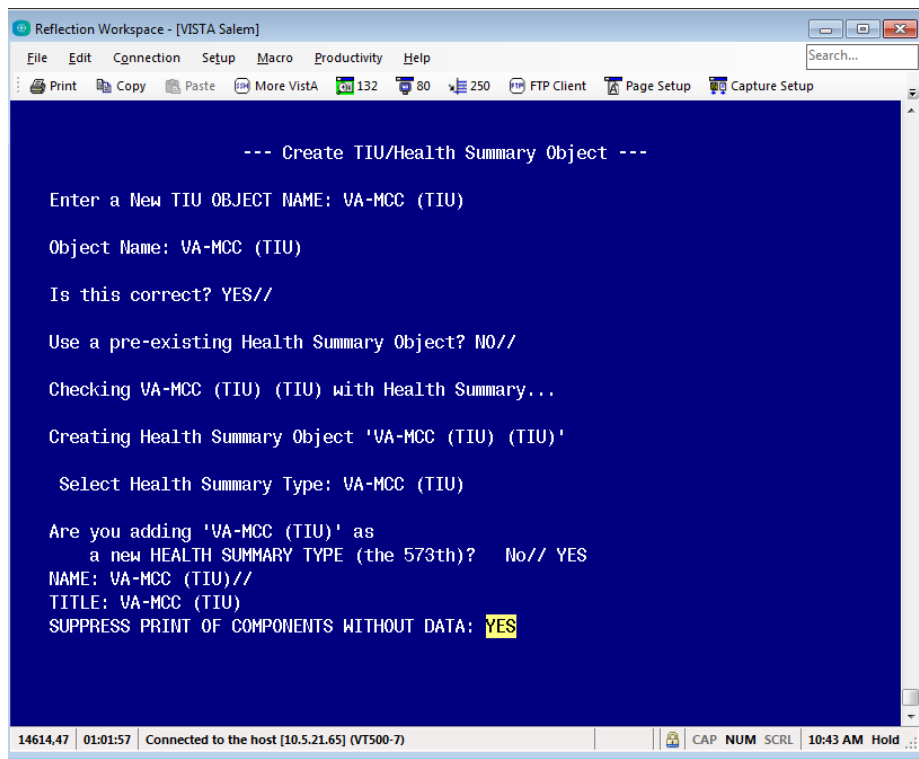

Step 11: At the 'SUPPRESS SENSITIVE PRINT DATA:' prompt, type "NO SSN" and press ENTER

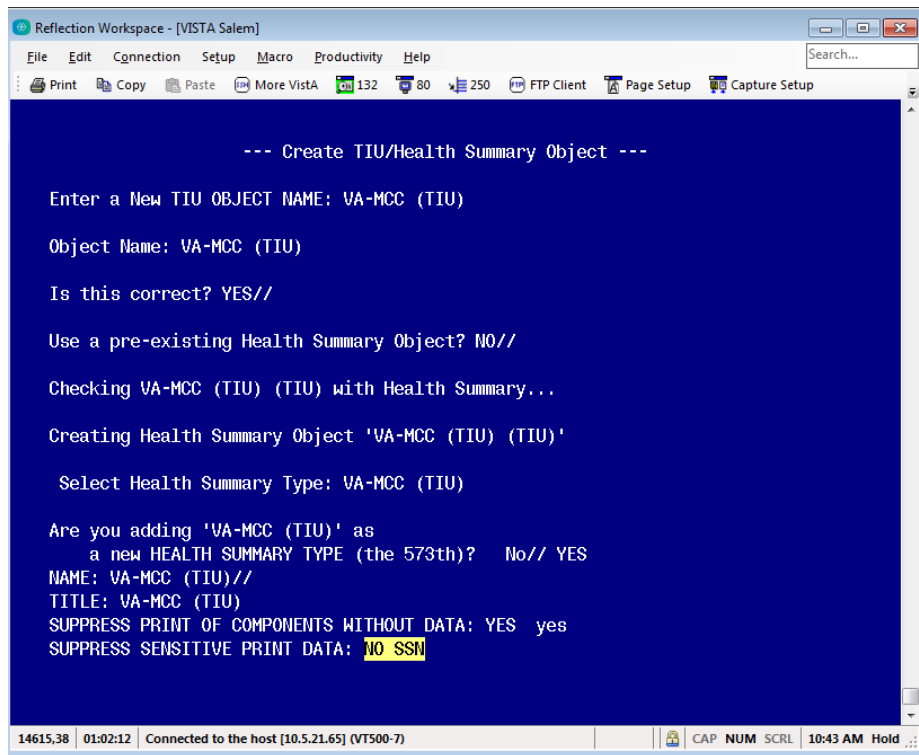

Reflection Workspace - [VISTA Salem]

File Edit Connection Setup Macro Productivity Help

Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

```
--- Create TIU/Health Summary Object ---

Enter a New TIU OBJECT NAME: VA-MCC (TIU)

Object Name: VA-MCC (TIU)

Is this correct? YES//

Use a pre-existing Health Summary Object? NO//

Checking VA-MCC (TIU) (TIU) with Health Summary...


Creating Health Summary Object 'VA-MCC (TIU) (TIU)'

Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN
```

14615,38 01:02:12 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:43 AM Hold

Step 12: At the 'LOCK' prompt, press ENTER

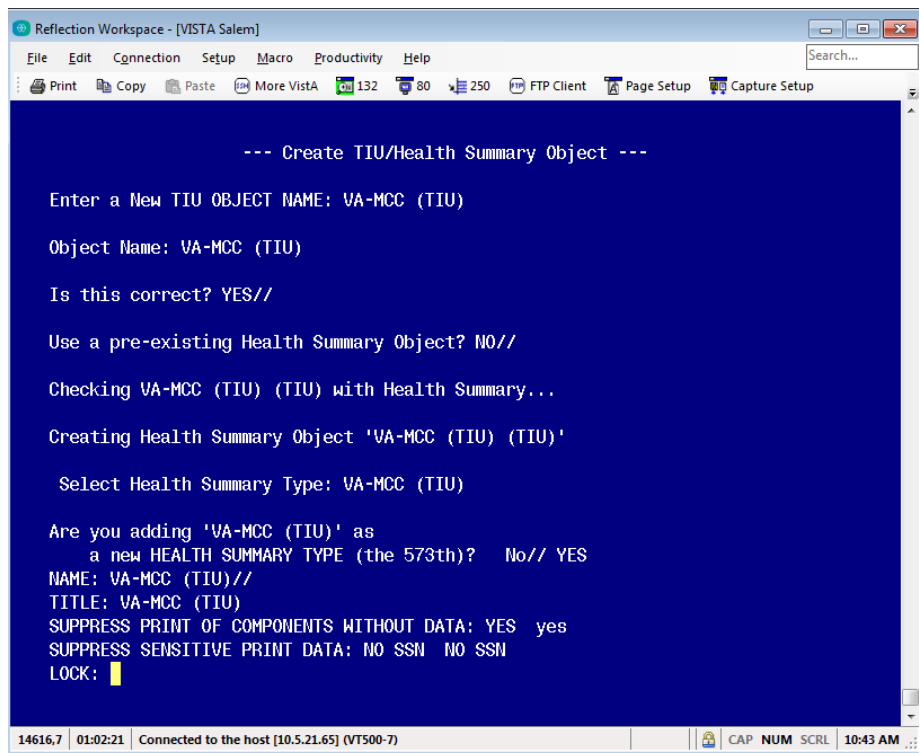

Reflection Workspace - [VISTA Salem]

File Edit Connection Setup Macro Productivity Help

Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

```
--- Create TIU/Health Summary Object ---

Enter a New TIU OBJECT NAME: VA-MCC (TIU)

Object Name: VA-MCC (TIU)

Is this correct? YES//

Use a pre-existing Health Summary Object? NO//

Checking VA-MCC (TIU) (TIU) with Health Summary...


Creating Health Summary Object 'VA-MCC (TIU) (TIU)'

Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
```

14616,7 01:02:21 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:43 AM

Step 13: At the 'OWNER' prompt, press **ENTER**

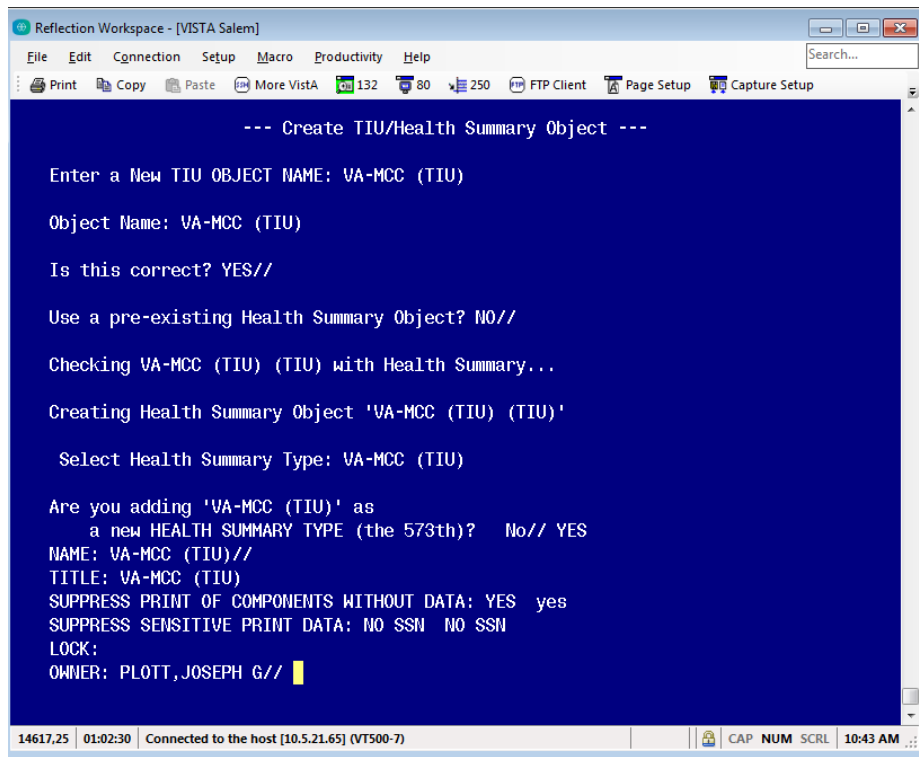

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

--- Create TIU/Health Summary Object ---

Enter a New TIU OBJECT NAME: VA-MCC (TIU)

Object Name: VA-MCC (TIU)

Is this correct? YES//

Use a pre-existing Health Summary Object? NO//

Checking VA-MCC (TIU) (TIU) with Health Summary...


Creating Health Summary Object 'VA-MCC (TIU) (TIU)'

Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//
```

14617,25 01:02:30 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:43 AM

Step 14: At the 'Do you wish to copy ...' prompt, type **NO** and press ENTER

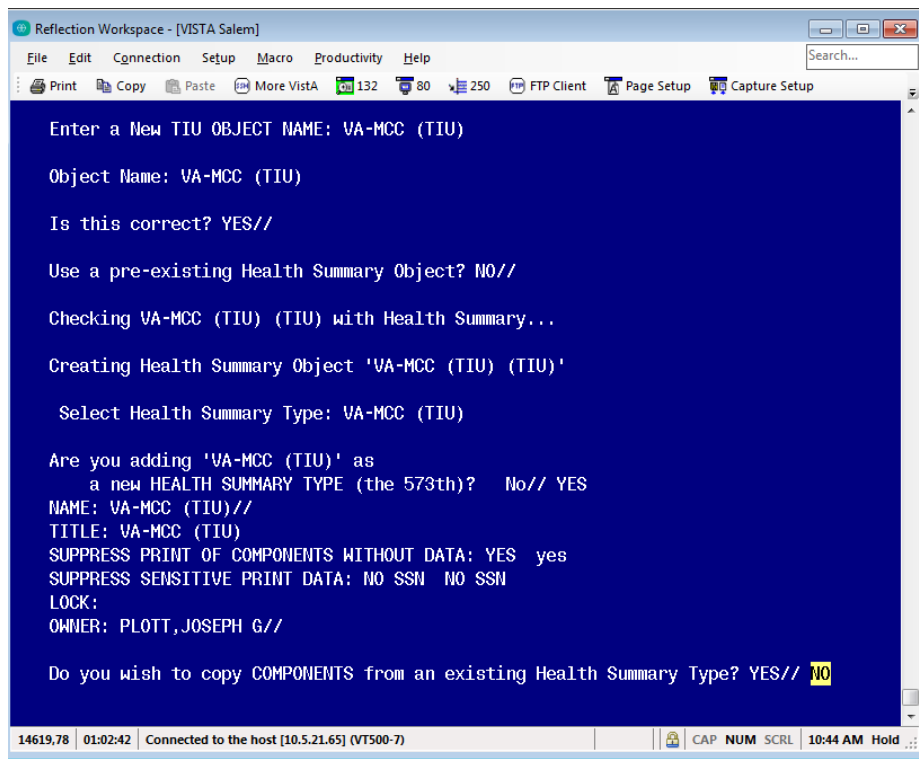

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Enter a New TIU OBJECT NAME: VA-MCC (TIU)

Object Name: VA-MCC (TIU)

Is this correct? YES//

Use a pre-existing Health Summary Object? NO//

Checking VA-MCC (TIU) (TIU) with Health Summary...

Creating Health Summary Object 'VA-MCC (TIU) (TIU)'


Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//

Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
```

14619,78 01:02:42 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:44 AM Hold

Step 15: At the 'Select COMPONENT:' prompt, type "SPN" and press ENTER

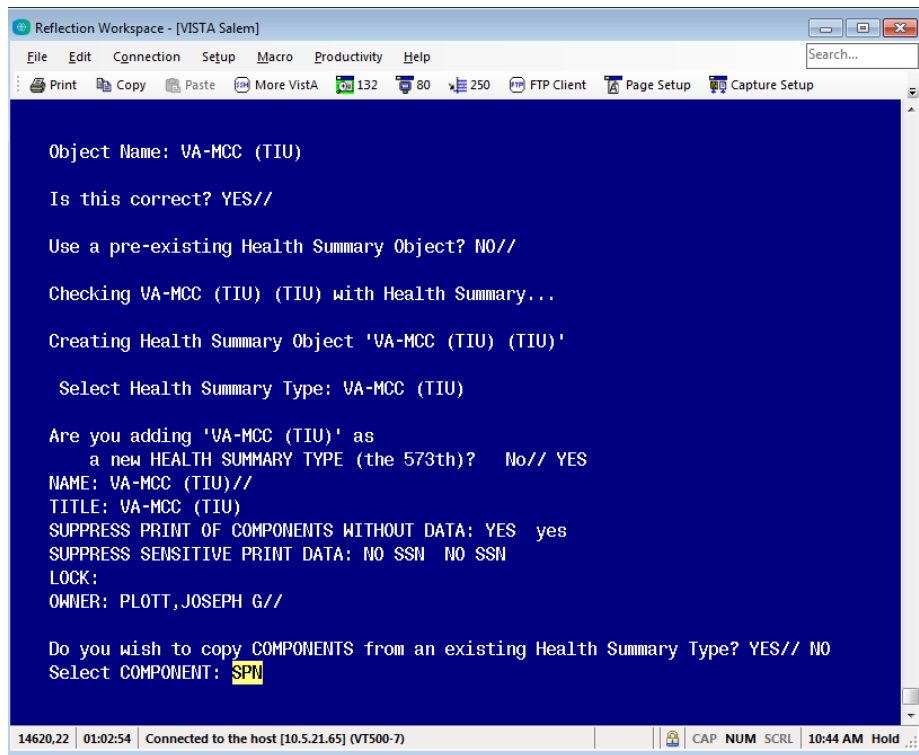


```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Object Name: VA-MCC (TIU)
Is this correct? YES//
Use a pre-existing Health Summary Object? NO//
Checking VA-MCC (TIU) (TIU) with Health Summary...
Creating Health Summary Object 'VA-MCC (TIU) (TIU)'
Select Health Summary Type: VA-MCC (TIU)
Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//
Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN

14620,22 01:02:54 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:44 AM Hold
```

Step 16: Select # 1 (i.e. SPN) and press ENTER

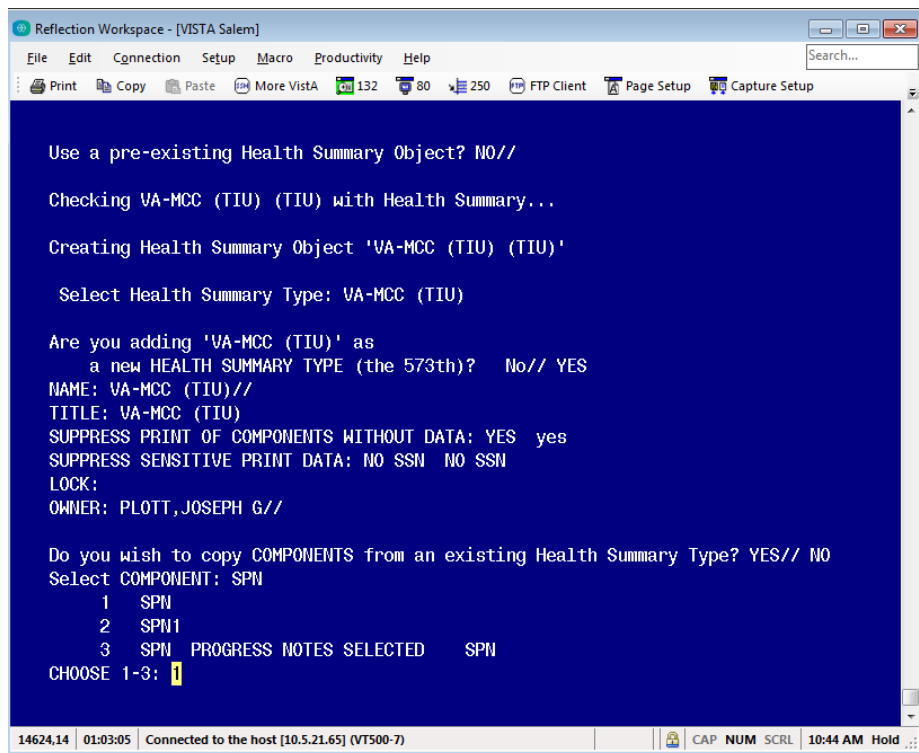


```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Use a pre-existing Health Summary Object? NO//
Checking VA-MCC (TIU) (TIU) with Health Summary...
Creating Health Summary Object 'VA-MCC (TIU) (TIU)'
Select Health Summary Type: VA-MCC (TIU)
Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//
Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN
1 SPN
2 SPN I
3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1

14624,14 01:03:05 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:44 AM Hold
```

Step 17: At the 'SUMMARY ORDER:' prompt, leave the default as-is and press **ENTER**

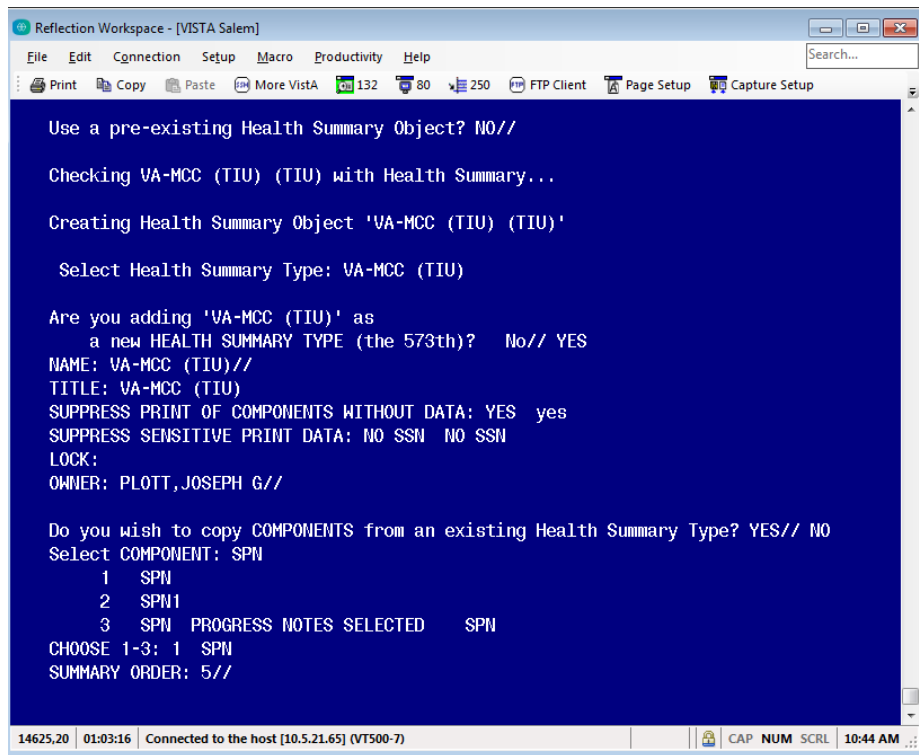

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Use a pre-existing Health Summary Object? NO//

Checking VA-MCC (TIU) (TIU) with Health Summary...

Creating Health Summary Object 'VA-MCC (TIU) (TIU)'


Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//

Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN
1 SPN
2 SPN1
3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5//

14625,20 01:03:16 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:44 AM
```

Step 18: At the 'OCCURRENCE LIMIT:' prompt, type **8** and press ENTER

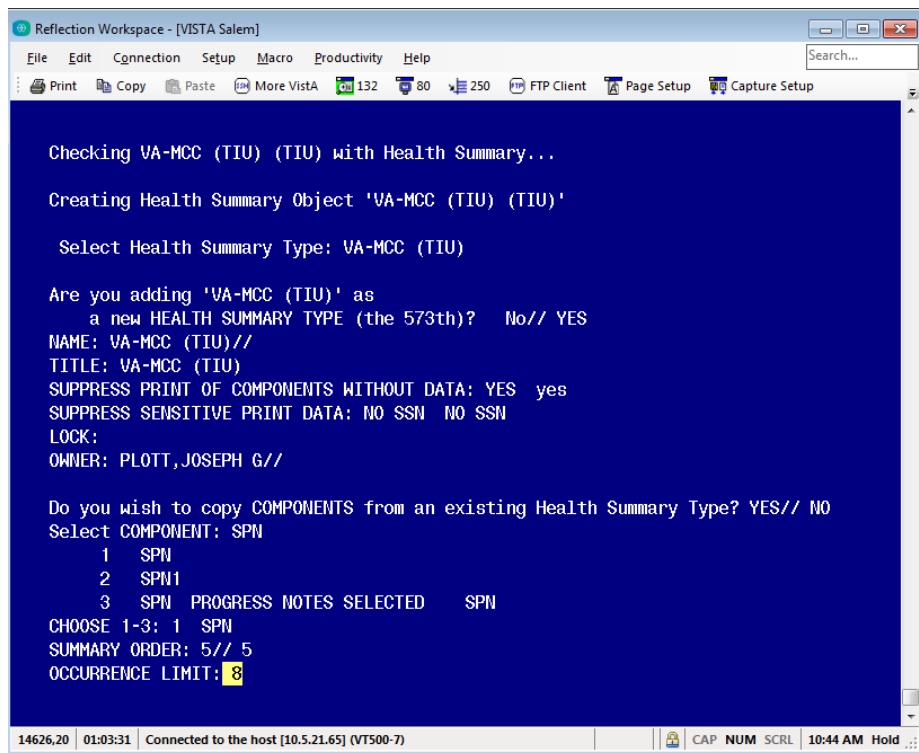

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Checking VA-MCC (TIU) (TIU) with Health Summary...

Creating Health Summary Object 'VA-MCC (TIU) (TIU)'


Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//

Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN
1 SPN
2 SPN1
3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8

14626,20 01:03:31 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:44 AM Hold
```

Step 19: At the 'TIME LIMIT:' prompt, just press **ENTER**

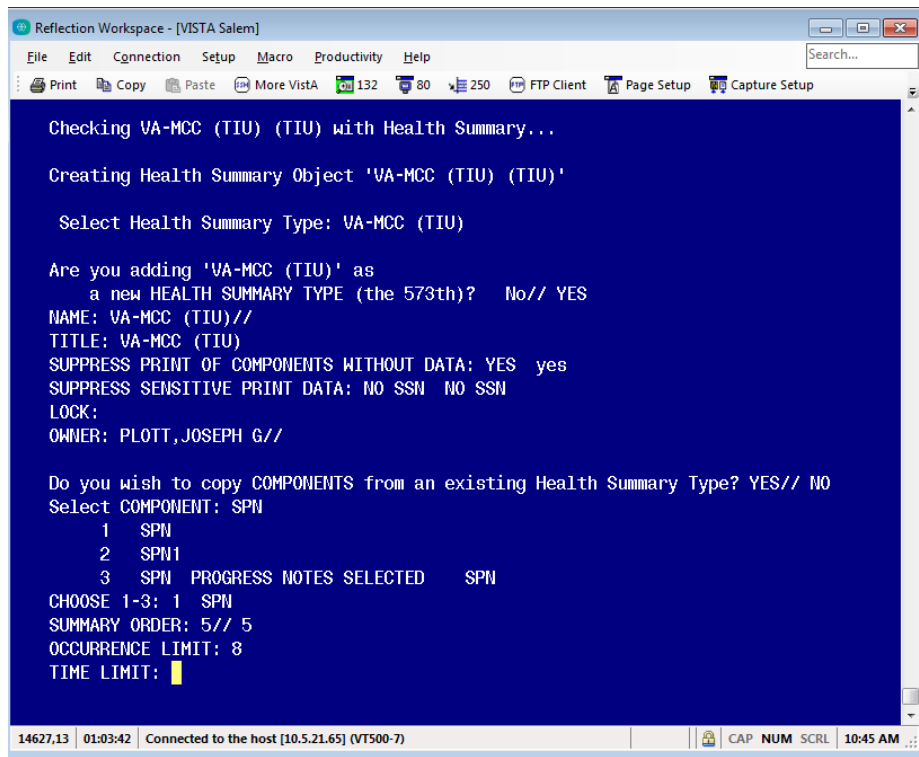

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Checking VA-MCC (TIU) (TIU) with Health Summary...


Creating Health Summary Object 'VA-MCC (TIU) (TIU)'

Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
  a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//

Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN
  1 SPN
  2 SPN1
  3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
14627,13 01:03:42 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:45 AM
```

Step 20: At the 'HEADER NAME:' prompt, just press **ENTER**

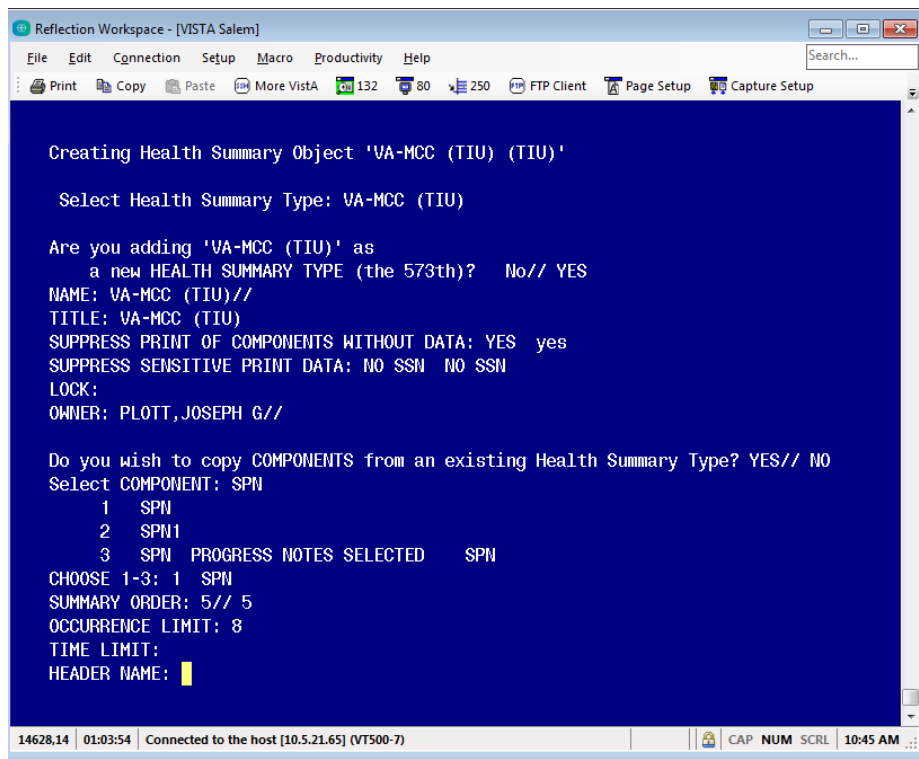

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup


Creating Health Summary Object 'VA-MCC (TIU) (TIU)'

Select Health Summary Type: VA-MCC (TIU)

Are you adding 'VA-MCC (TIU)' as
  a new HEALTH SUMMARY TYPE (the 573th)? No// YES
NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//

Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN
  1 SPN
  2 SPN1
  3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME:
14628,14 01:03:54 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:45 AM
```

Step 21: At the 'SELECTION ITEM:' prompt, type "**VA-MATERNITY CARE COORDINATION**" and press ENTER

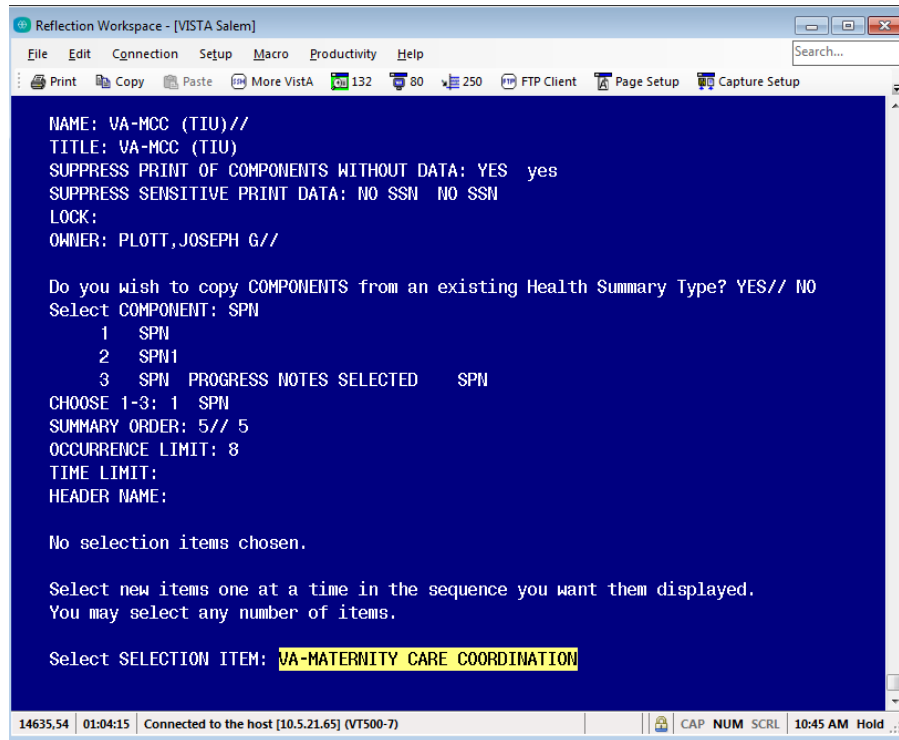

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More Vista 132 80 250 FTP Client Page Setup Capture Setup

NAME: VA-MCC (TIU)//
TITLE: VA-MCC (TIU)
SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: YES yes
SUPPRESS SENSITIVE PRINT DATA: NO SSN NO SSN
LOCK:
OWNER: PLOTT,JOSEPH G//

Do you wish to copy COMPONENTS from an existing Health Summary Type? YES// NO
Select COMPONENT: SPN
  1 SPN
  2 SPN1
  3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME:


No selection items chosen.

Select new items one at a time in the sequence you want them displayed.
You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

14635.54 01:04:15 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:45 AM Hold
```

Step 22: Leave the default at "**YES**" and press ENTER

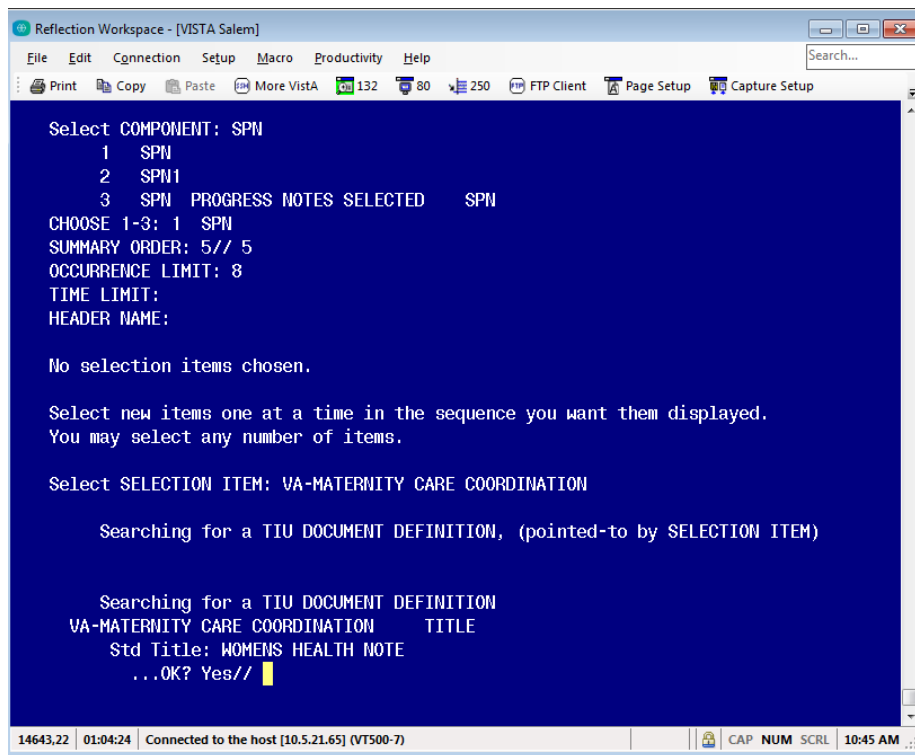

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More Vista 132 80 250 FTP Client Page Setup Capture Setup

Select COMPONENT: SPN
  1 SPN
  2 SPN1
  3 SPN PROGRESS NOTES SELECTED SPN
CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME:

No selection items chosen.

Select new items one at a time in the sequence you want them displayed.
You may select any number of items.


Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes//

14643.22 01:04:24 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:45 AM
```


Step 23: At the 'SELECTION ITEM:' prompt, just press **ENTER**

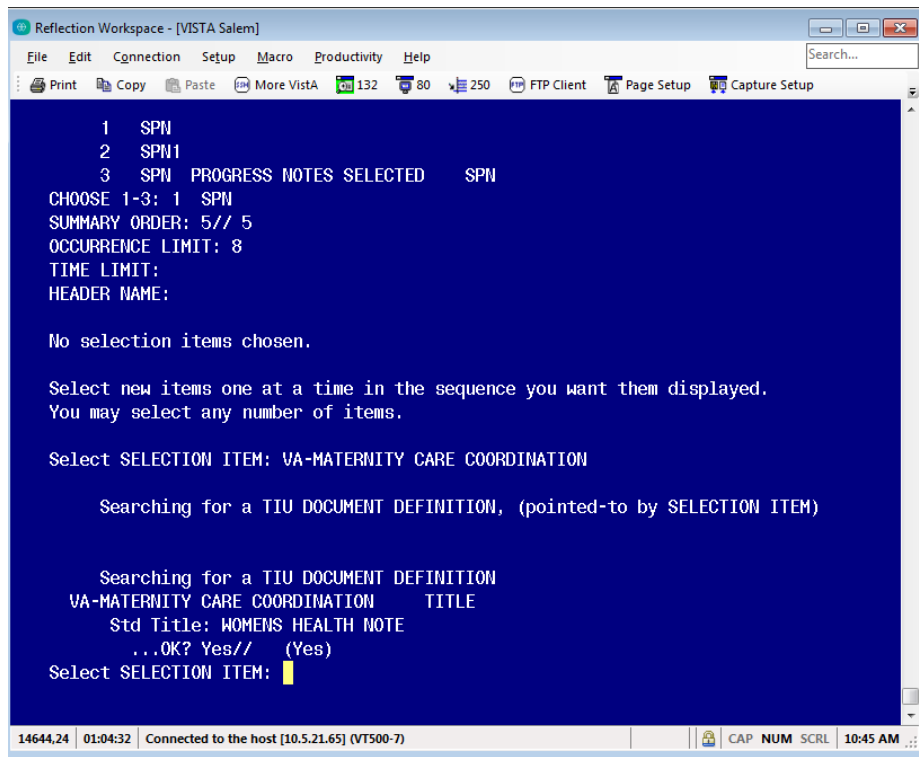

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

1  SPN
2  SPN1
3  SPN  PROGRESS NOTES SELECTED  SPN
CHOOSE 1-3: 1  SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME:

No selection items chosen.


Select new items one at a time in the sequence you want them displayed.
You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION  TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM: 
```

Step 24: At the 'COMPONENT:' prompt, just press **ENTER**

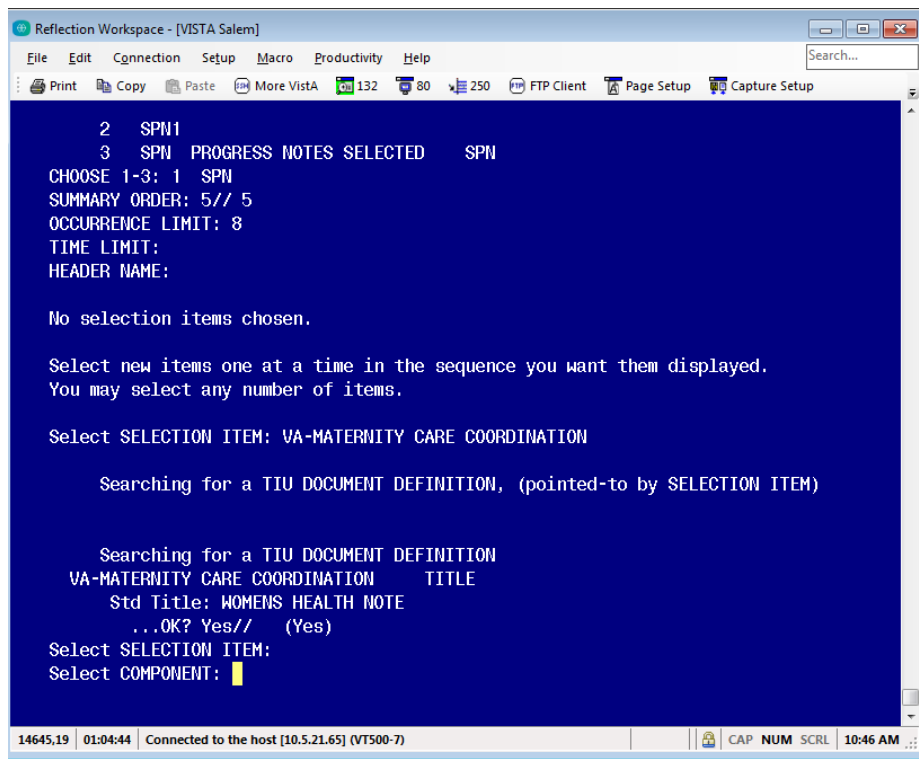

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

2  SPN1
3  SPN  PROGRESS NOTES SELECTED  SPN
CHOOSE 1-3: 1  SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME:

No selection items chosen.


Select new items one at a time in the sequence you want them displayed.
You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION  TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT: 
```

Step 25: At the ‘Do you wish to review ...?’ prompt, leave at “NO” and press ENTER

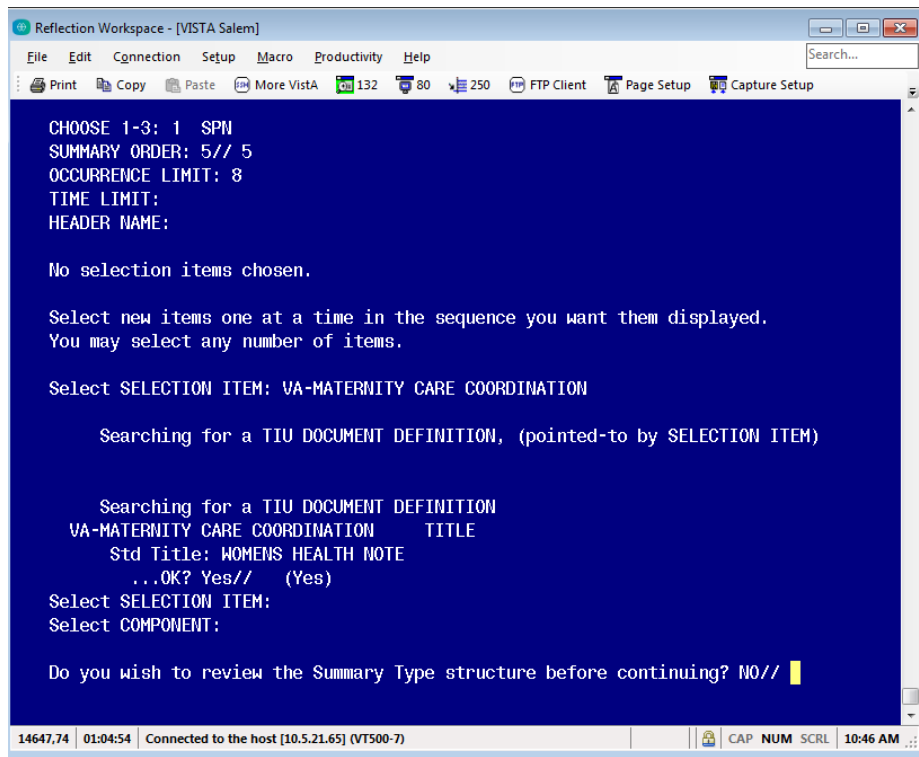

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

CHOOSE 1-3: 1 SPN
SUMMARY ORDER: 5// 5
OCCURRENCE LIMIT: 8
TIME LIMIT:
HEADER NAME:

No selection items chosen.

Select new items one at a time in the sequence you want them displayed.
You may select any number of items.


Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
```

Step 26: At the ‘...TIME LIMITS ...’ prompt, leave at “NO” and press ENTER

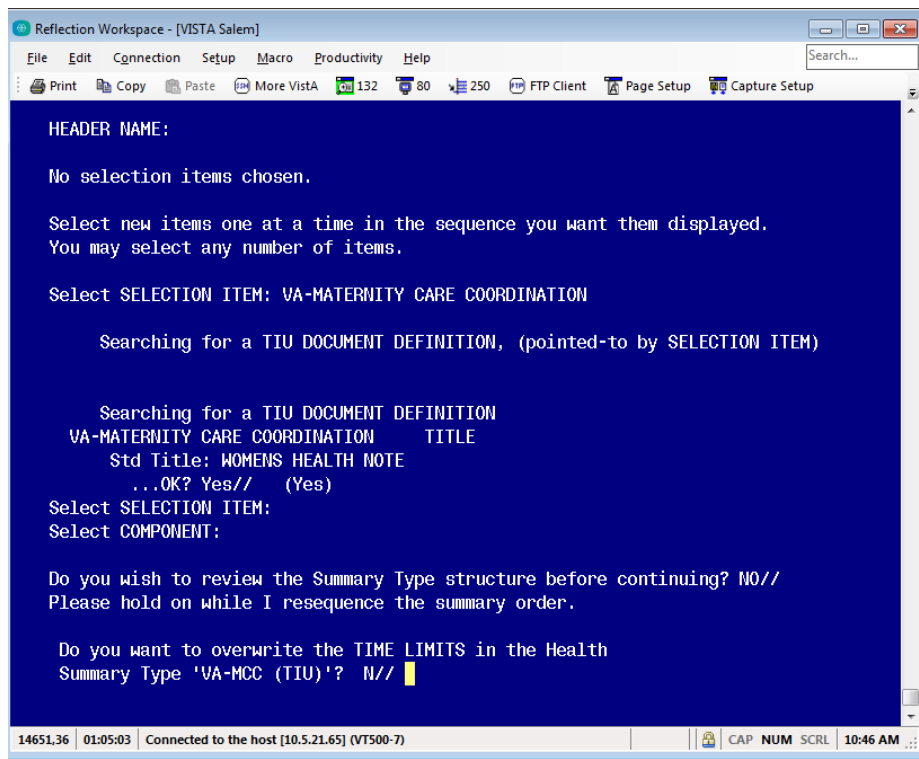

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

HEADER NAME:

No selection items chosen.

Select new items one at a time in the sequence you want them displayed.
You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION


Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N//
```

Step 27: At the ‘... Health Summary Header ...’ prompt, leave at “NO” and press ENTER

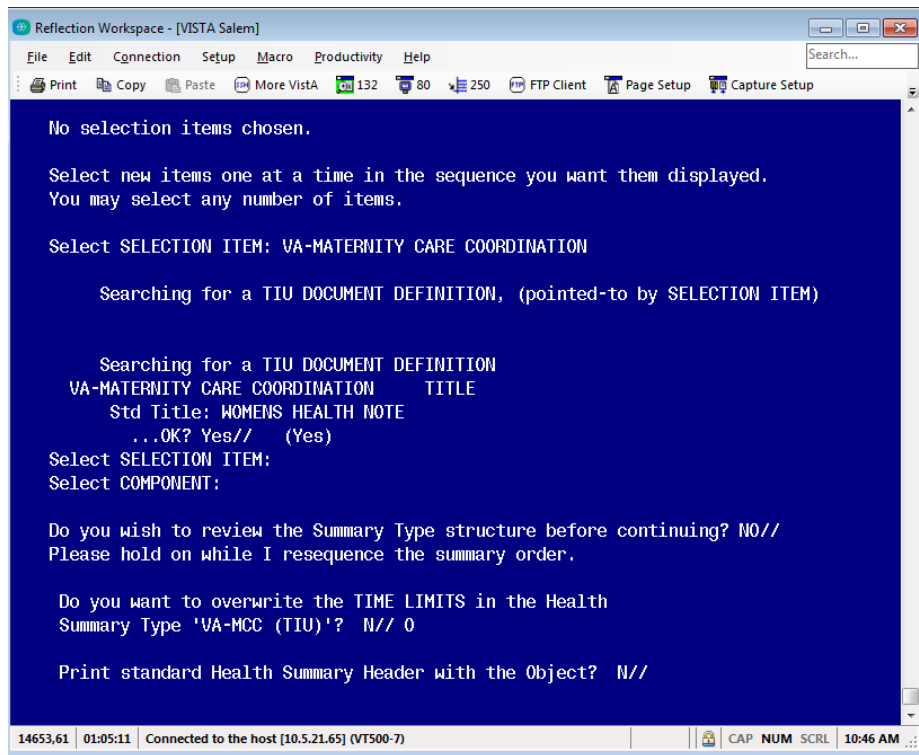

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

No selection items chosen.

Select new items one at a time in the sequence you want them displayed.
You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

    Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

        Searching for a TIU DOCUMENT DEFINITION
        VA-MATERNITY CARE COORDINATION    TITLE
        Std Title: WOMENS HEALTH NOTE
        ...OK? Yes//    (Yes)
Select SELECTION ITEM:
Select COMPONENT:


Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'?  N// 0

Print standard Health Summary Header with the Object?  N//

14653,61 | 01:05:11 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL 10:46 AM
```

Step 28: At the ‘Print Report Date?’ prompt, leave at “NO” and press ENTER

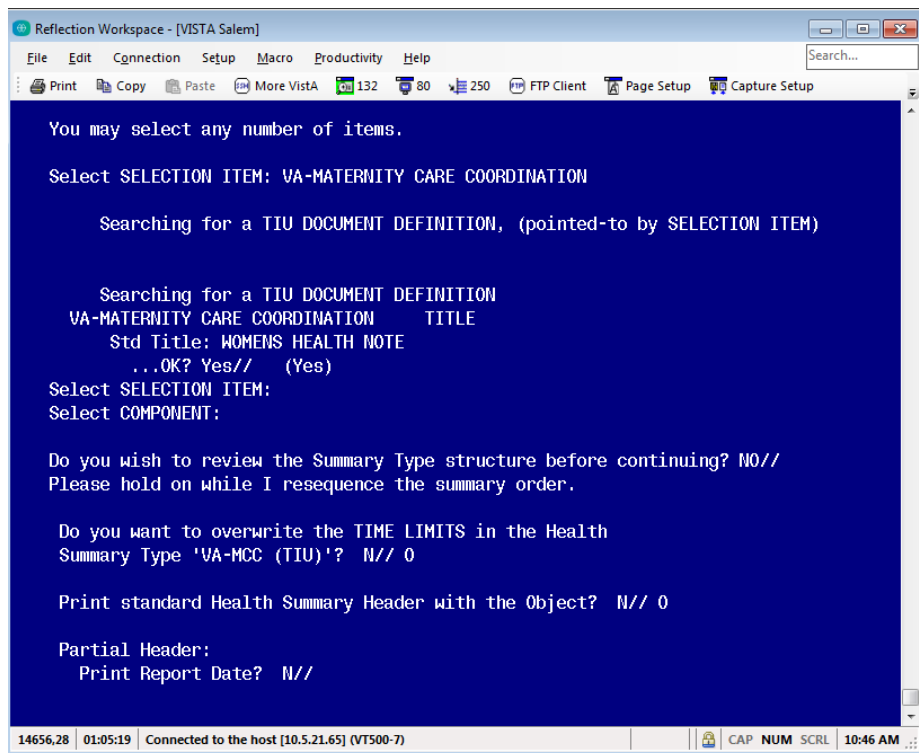

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

You may select any number of items.

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

    Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

        Searching for a TIU DOCUMENT DEFINITION
        VA-MATERNITY CARE COORDINATION    TITLE
        Std Title: WOMENS HEALTH NOTE
        ...OK? Yes//    (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.


Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'?  N// 0

Print standard Health Summary Header with the Object?  N// 0

Partial Header:
Print Report Date?  N//

14656,28 | 01:05:19 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL 10:46 AM
```

Step 29: At the ‘Print Confidentiality Banner?’ prompt, leave at “NO” and press ENTER

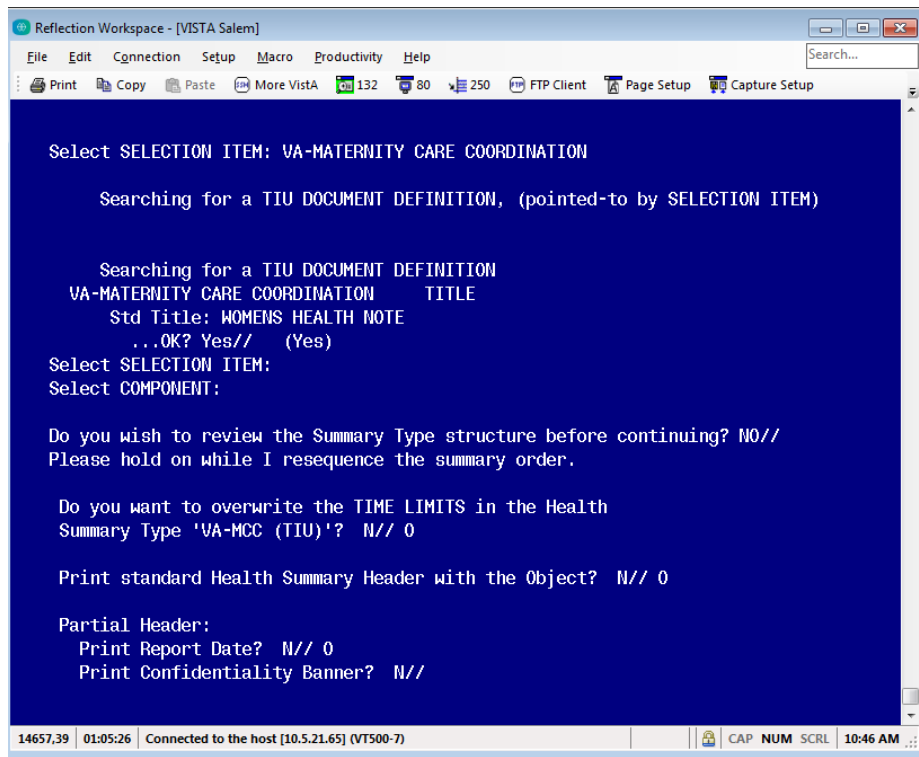

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.


Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N//

14657,39 | 01:05:26 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL 10:46 AM
```

Step 30: At the ‘Print Report Header?’ prompt, leave at “NO” and press ENTER

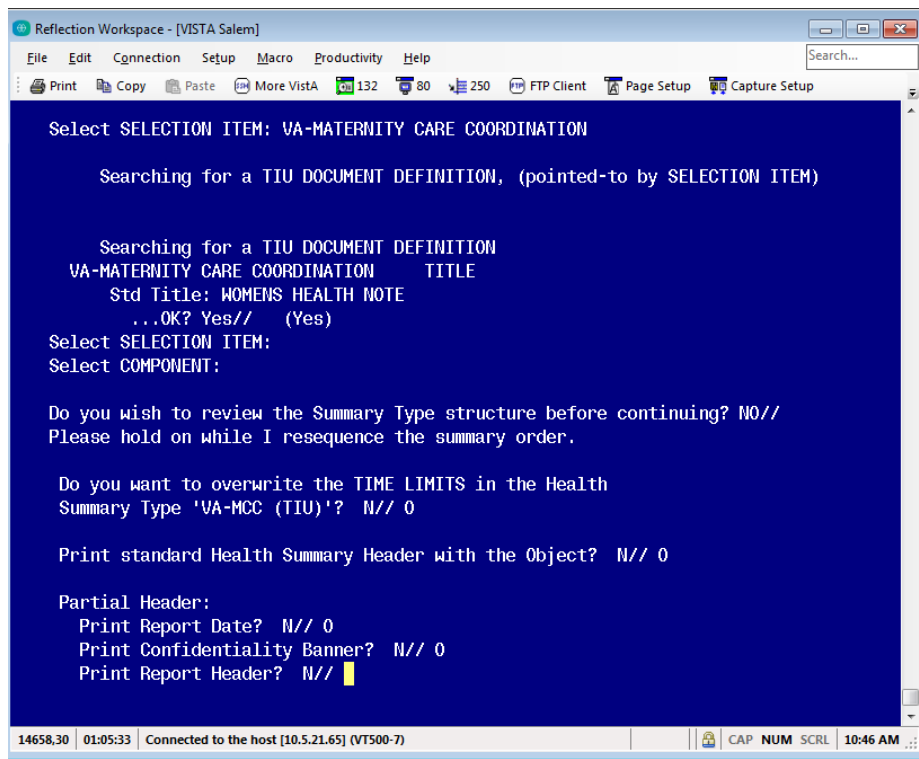

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Select SELECTION ITEM: VA-MATERNITY CARE COORDINATION

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.


Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N//

14658,30 | 01:05:33 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL 10:46 AM
```

Step 31: At the ‘Print the standard Component Header?’ prompt, change to “NO” and press ENTER

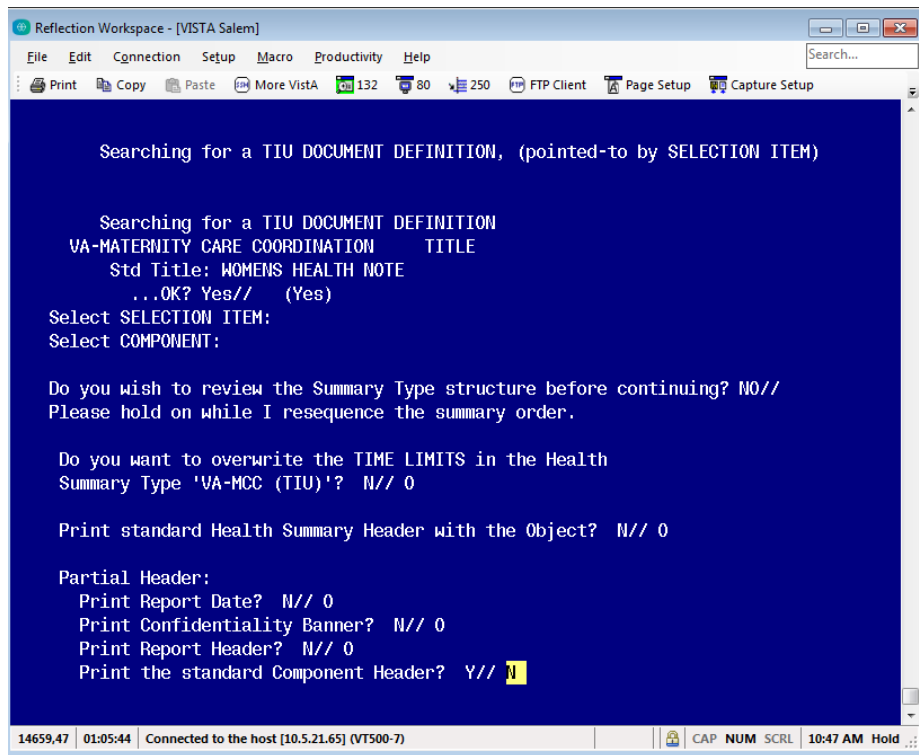

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.


Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N// 0
Print the standard Component Header? Y// 1

14659,47 01:05:44 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:47 AM Hold
```

Step 32: At the ‘Print the date a patient was deceased?’ prompt, leave at “NO” and press ENTER

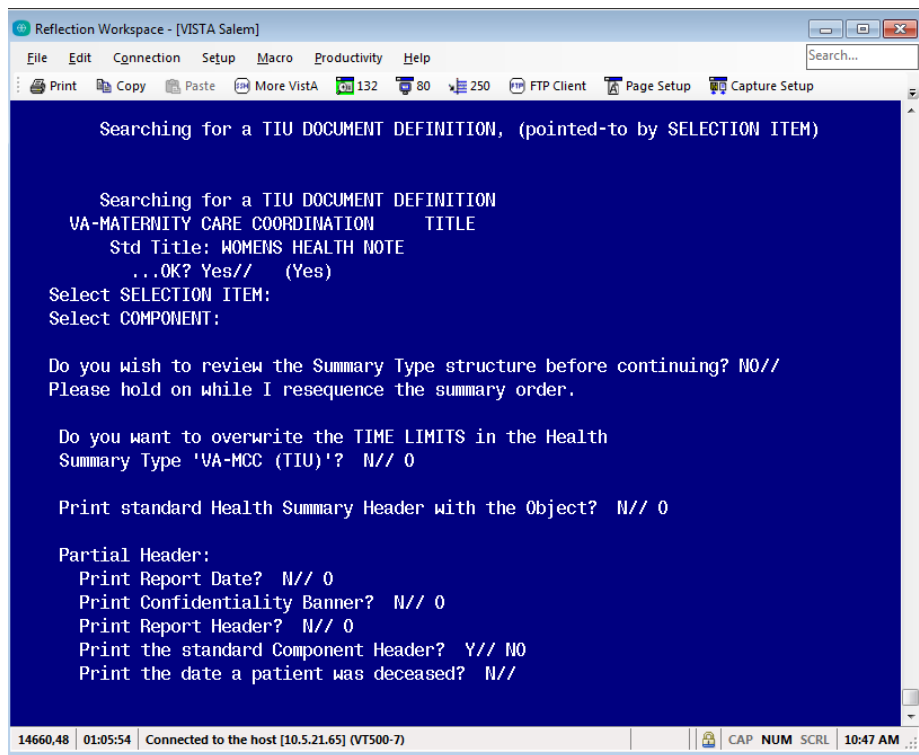

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Searching for a TIU DOCUMENT DEFINITION, (pointed-to by SELECTION ITEM)

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.


Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N// 0
Print the standard Component Header? Y// NO
Print the date a patient was deceased? N//

14660,48 01:05:54 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:47 AM
```

Step 33: At the ‘Print a LABEL ...?’ prompt, leave at “NO” and press ENTER

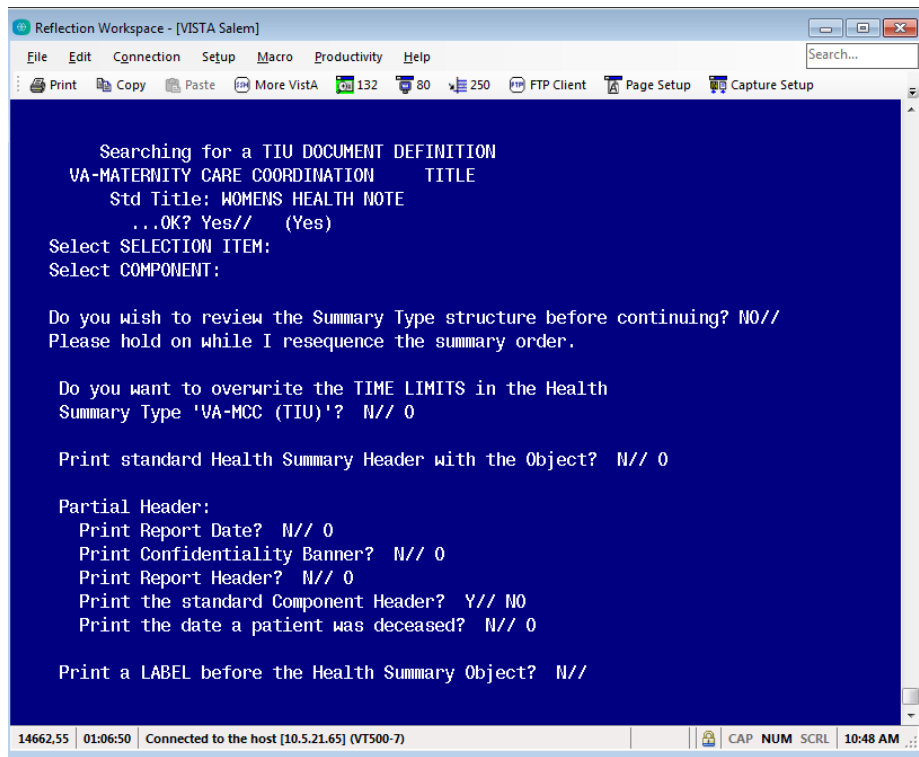

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Searching for a TIU DOCUMENT DEFINITION
VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0


Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N// 0
Print the standard Component Header? Y// NO
Print the date a patient was deceased? N// 0

Print a LABEL before the Health Summary Object? N//

14662,55 01:06:50 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:48 AM
```

Step 34: At the ‘Suppress Components without data?’ prompt, type “YES” and press ENTER

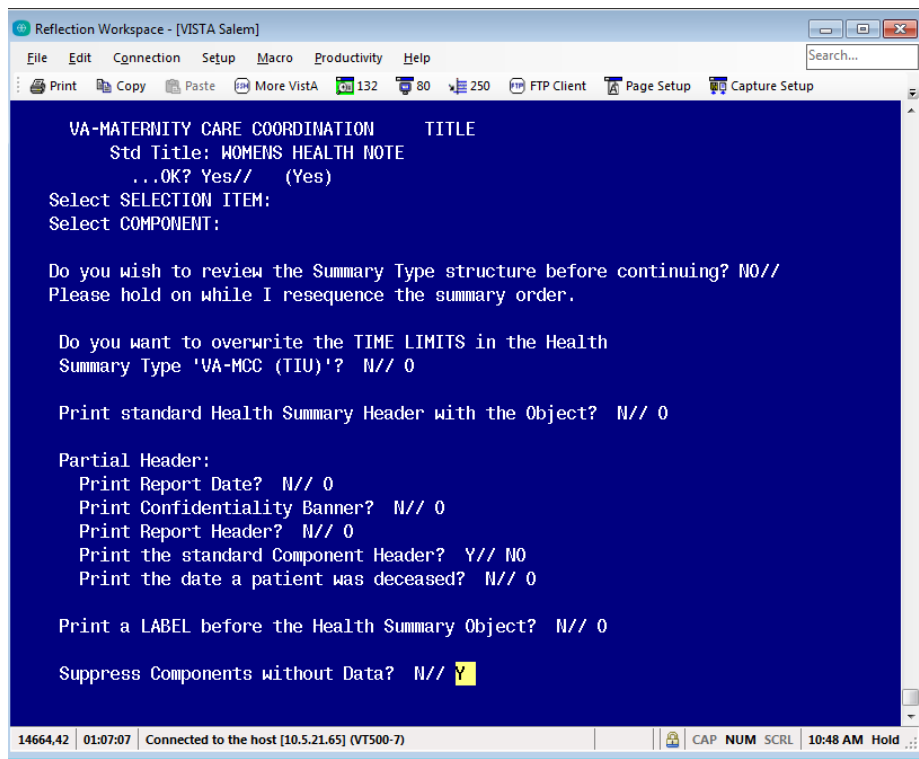

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

VA-MATERNITY CARE COORDINATION TITLE
Std Title: WOMENS HEALTH NOTE
...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0


Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N// 0
Print the standard Component Header? Y// NO
Print the date a patient was deceased? N// 0

Print a LABEL before the Health Summary Object? N// 0

Suppress Components without Data? N// Y

14664,42 01:07:07 Connected to the host [10.5.21.65] (VT500-7) CAP NUM SCRL 10:48 AM Hold
```

Step 35: At the 'Overwrite No Data Available Message' prompt, just press **ENTER**

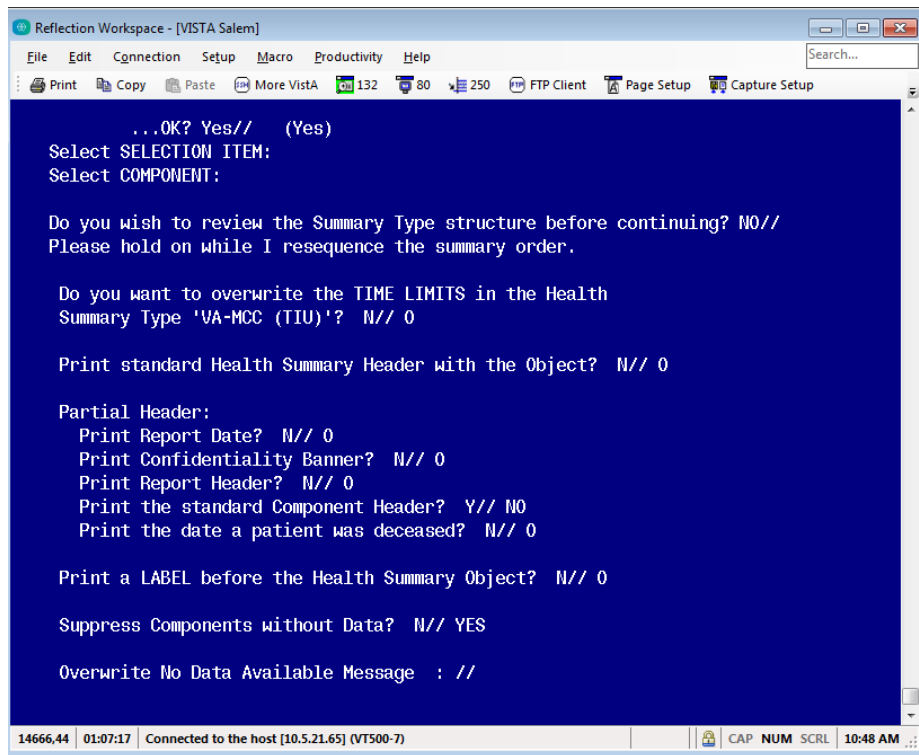

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

...OK? Yes// (Yes)
Select SELECTION ITEM:
Select COMPONENT:

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N// 0
Print the standard Component Header? Y// NO
Print the date a patient was deceased? N// 0


Print a LABEL before the Health Summary Object? N// 0

Suppress Components without Data? N// YES

Overwrite No Data Available Message : //

14666,44 | 01:07:17 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL | 10:48 AM
```

Step 36: At the 'OBJECT DESCRIPTION' prompt, type **YES** and press **ENTER**

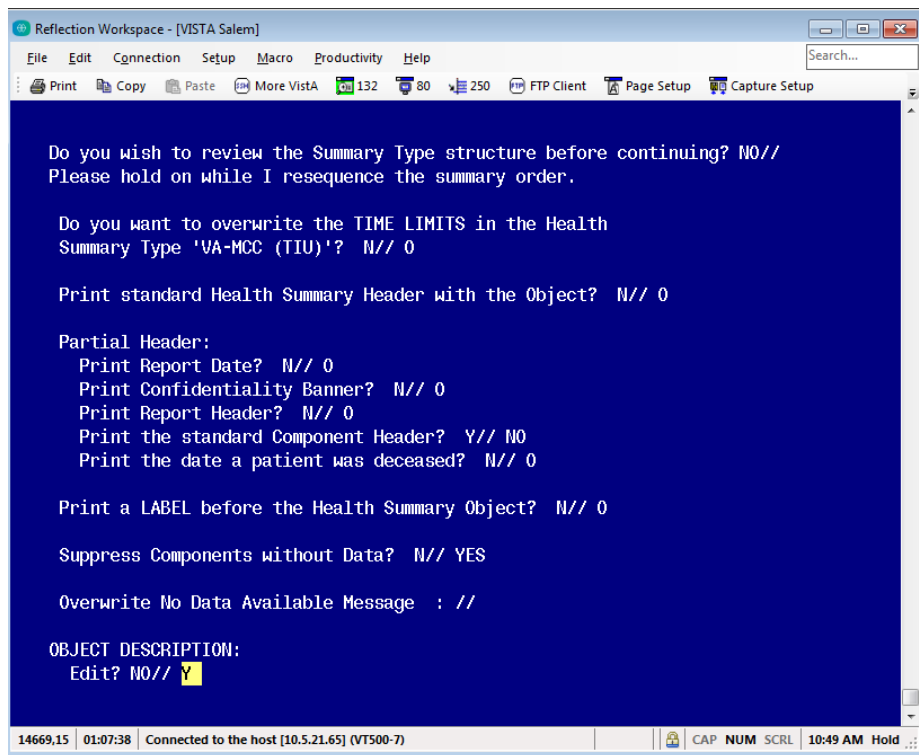

```
Reflection Workspace - [VISTA Salem]
File Edit Connection Setup Macro Productivity Help
Print Copy Paste More VistA 132 80 250 FTP Client Page Setup Capture Setup

Do you wish to review the Summary Type structure before continuing? NO//
Please hold on while I resequence the summary order.

Do you want to overwrite the TIME LIMITS in the Health
Summary Type 'VA-MCC (TIU)'? N// 0

Print standard Health Summary Header with the Object? N// 0

Partial Header:
Print Report Date? N// 0
Print Confidentiality Banner? N// 0
Print Report Header? N// 0
Print the standard Component Header? Y// NO
Print the date a patient was deceased? N// 0

Print a LABEL before the Health Summary Object? N// 0


Suppress Components without Data? N// YES

Overwrite No Data Available Message : //

OBJECT DESCRIPTION:
Edit? NO// Y

14669,15 | 01:07:38 | Connected to the host [10.5.21.65] (VT500-7) | CAP NUM SCRL | 10:49 AM Hold
```

Step 37: Enter the text you see below and press F1+E

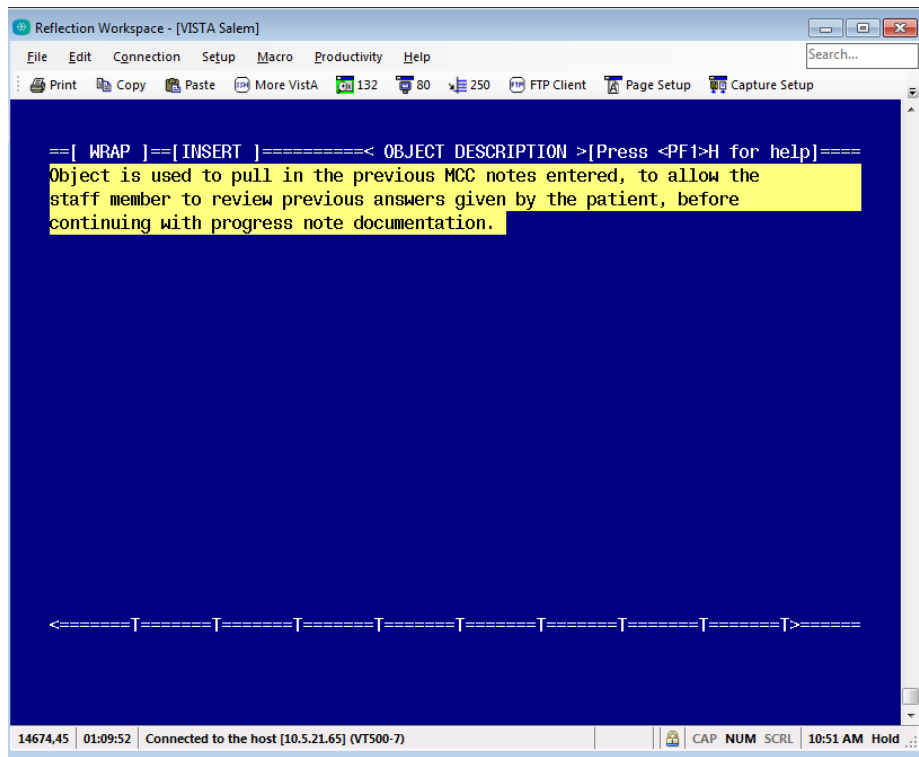

Step 38: At the 'Ok?' prompt, just press **ENTER**

Step 39: Done! Press **ENTER**

NOTE: If you require further technical assistance, if there are any questions or problems during the installation, please notify your local IT support to log a national help desk ticket or contact the VA Service Desk and have them submit a national ticket to:
Category: Enterprise Applications
Subcategory: Applications
Enterprise Application: VistA - Clinical Reminders